provided by St George's Online Research Arch PLOS | NEGLECTED | TROPICAL DISEASES

From Innovation to Application

Detection of Viable *Mycobacterium ulcerans* in Clinical Samples by a Novel Combined 16S rRNA Reverse Transcriptase/IS*2404* Real-Time qPCR Assay

Marcus Beissner¹**, Dominik Symank¹*, Richard Odame Phillips², Yaw Ampem Amoako², Nana-Yaa Awua-Boateng³, Fred Stephen Sarfo², Moritz Jansson¹, Kristina Lydia Huber¹, Karl-Heinz Herbinger¹, Florian Battke⁴, Thomas Löscher¹, Ohene Adjei², Gisela Bretzel¹

1 Department of Infectious Diseases and Tropical Medicine (DITM), University Hospital, Ludwig-Maximilians University, Munich, Germany, 2 Komfo Anokye Teaching Hospital (KATH), Kwame Nkrumah University of Science and Technology (KNUST), Kumasi, Ghana, 3 Kumasi Centre for Collaborative Research (KCCR), Kwame Nkrumah University of Science and Technology (KNUST), Kumasi, Ghana, 4 Dr. Battke SCIENTIA GmbH, Taufkirchen, Germany

Introduction

Buruli ulcer disease (BUD) caused by Mycobacterium ulcerans involves the skin and soft tissue. If left untreated, extensive destruction of tissue followed by scarring and contractures may lead to severe functional limitations. Following the introduction of standardized antimycobacterial chemotherapy with rifampicin and streptomycin, recurrence rates of less than 2% were reported. However, treatment failures occur and a variety of secondary lesions necessitating customized clinical management strategies have been reported. True recurrences by definition occur more than three months after completion of antibiotic treatment, are characterised by the presence of viable bacilli, and require a second course of antibiotics. "Non-healers" may harbour viable, possibly drug-resistant M. ulcerans strains and may benefit from surgical intervention. Early-onset immune-mediated paradoxical reactions emerging during or shortly after treatment do not contain viable bacilli and may heal under conventional wound care and/or minor surgery; lateonset secondary lesions presumably attributable to secondary infection foci may clear spontaneously through enhanced immune responses primed by initial treatment. None of the current diagnostic techniques is applicable to rapidly address the pivotal question of the presence of viable bacilli in non-healers and patients with secondary BUD lesions, and optimal time points for collection of follow-up samples have not yet been investigated. Therefore, to date treatment monitoring is mainly based on clinical observation [1-5]. Reverse transcriptase assays targeting 16S rRNA and mRNA were successfully applied for the rapid detection of viable mycobacteria in clinical samples from patients with tuberculosis and leprosy [6,7]. To employ this technique for classification of BUD lesions and monitoring of treatment success we developed a *M. ulcerans*—specific RNA-based viability assay combining a 16S rRNA reverse transcriptase real-time PCR (RT-qPCR) to determine bacterial viability with an IS2404 quantitative real-time PCR (qPCR) for increased specificity and simultaneous quantification of bacilli.

Development and Validation Ethical Approval

The study was approved by the Committee of Human Research Publication and Ethics, School of Medical Sciences, Kwame Nkrumah University of Science and Technology, Kumasi, Ghana (CHRPE/28/09). Written informed consent was obtained from all study participants, or their legal representatives.

Bacterial Strains, DNA Extracts, and Clinical Samples

Technical validation of the assay was performed with 29 *M. ulcerans* strains originating from Cameroon [8] and Ghana (Table 1), as well as DNA extracts from 18 closely related human pathogenic mycobacterial species and five bacterial species frequently colonizing human skin (Table 2).

Clinical validation was conducted on pre-treatment swab samples in PANTA (BD, Heidelberg, Germany) from 24 suspected BUD cases from Agogo Presbyterian Hospital (n=14) and Tepa Government Hospital (n=10), Ghana (Protocol S1). In addition, post-treatment swab samples from seven IS2404 PCR confirmed BUD patients with incomplete wound healing were collected at week nine (Figures 1 and 2).

All clinical samples were subjected to routine diagnostics (microscopy and IS2404 dry-reagent-based [DRB] PCR) at the Kumasi Centre for Collaborative Research (KCCR) [3].

Primers and Probes

Primers and a hydrolysis probe (TibMol-Biol, Berlin, Germany) for specific amplification of *M. ulcerans* 16S rRNA were designed using DNAsis Max (MiraiBio, San Francisco, USA) by alignment of 16S rRNA gene sequences (GenBank, National Center for Biotechnology Information [NCBI]) from closely related mycobacteria and other bacteria potentially contaminating the human skin (Table 2).

For simultaneous quantification by IS2404 qPCR, the primers described by Fyfe et al. [9] were used in combination

Citation: Beissner M, Symank D, Phillips RO, Amoako YA, Awua-Boateng N-Y, et al. (2012) Detection of Viable *Mycobacterium ulcerans* in Clinical Samples by a Novel Combined 16S rRNA Reverse Transcriptase/IS2404 Real-Time qPCR Assay. PLoS Negl Trop Dis 6(8): e1756. doi:10.1371/journal.pntd.0001756

Editor: Pamela L. C. Small, University of Tennessee, United States of America

Published August 28, 2012

Copyright: © 2012 Beissner et al. This is an open-access article distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

Funding: The research leading to these results has received funding from the European Community's Seventh Framework Programme (FP7/2007–2013) under grant agreement N° 241500. The funders had no role in study design, data collection and analysis, decision to publish, or preparation of the manuscript.

Competing Interests: The authors have declared that no competing interests exist.

- * E-mail: beissner@lrz.uni-muenchen.de
- 9 These authors contributed equally to this work.

1

Table 1. M. ulcerans cultures subjected to the 16S rRNA RT/IS2404 gPCR assay.

M. ulcerans Strain	Source	Origin ^a	16S rRNA RT-qPCR ^b	IS <i>2404</i> qPCR ^c	IS <i>2404</i> qPCR – Wipeout Control ^d
K4s-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K4s-C2	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K4s-C3	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5d-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5d-C2	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5d-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5d-C2	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5d-C3	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5d-C4	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5s-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5s-C2	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5s-C3	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5s-C4	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K5s-C5	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K7b-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K7b-C2	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K7b-C3	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K7b-C4	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K7s-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K7s-C2	DITM	Human isolate – Kamerun	Positive	Positive	Negative
K12S-C1	DITM	Human isolate – Kamerun	Positive	Positive	Negative
941328-C1	DITM	Human isolate – Ghana	Positive	Positive	Negative
07-C1	DITM	Human isolate – Ghana	Positive	Positive	Negative
DS1-C1	DITM	Human isolate – Ghana	Positive	Positive	Negative
97680-C1	DITM	Human isolate – Ghana	Positive	Positive	Negative
G.A.P.001-C1	KCCR	Human isolate – Ghana	Positive	Positive	Negative
G.A.P.033-C1	KCCR	Human isolate – Ghana	Positive	Positive	Negative
G.A.P.071-C1	KCCR	Human isolate – Ghana	Positive	Positive	Negative
G.A.P.078-C1	KCCR	Human isolate – Ghana	Positive	Positive	Negative

Table 1 shows 29 *M. ulcerans* cultures that were available at the Department of Infectious Diseases and Tropical Medicine (DITM) and the Kumasi Centre for Collaborative Research (KCCR) for development and technical validation of the 16S rRNA RT/IS2404 qPCR viability assay and the corresponding test results. Sequence analysis of 16S rRNA genes from the listed strains revealed 100% nucleotide concordance of the corresponding genomic regions amplified by the 16S rRNA RT-qPCR; no SNPs or mutations were detected, suggesting a high selectivity of the assay. Sequencing primers are described in Table 3 [11].

with a hydrolysis probe (Table 3) that was re-designed by DNAsis Max for thermodynamic reasons.

Combined RNA/DNA Extraction, Reverse Transcription, and Real-Time qPCR

Culture suspensions and swab samples were stabilized by RNA protect (Qiagen, Hilden, Germany) and subjected to AllPrep DNA/RNA extraction kit (Qiagen) (Protocols S1 and S2).

M. ulcerans whole transcriptome RNA from cultures and swab samples was transcribed to cDNA by QuantiTect Reverse Transcription Kit (Qiagen) including genomic DNA (gDNA) wipeout (Protocol S2). DNA and cDNA were subjected to IS2404 qPCR and 16S rRNA RT-qPCR, respectively, with corresponding controls (Table 4, Protocols S3 and S4).

Intra- and Inter-Assay Variability

Intra- and inter-assay variability was assessed by testing of each sample in quadruplicate within one 96-well plate, repeated on three different days (Table 5).

Sensitivity

The analytical sensitivity was determined as lower limit of detection (LOD, lowest template concentration rendering amplification of 95% of samples) [10] for

^aM. ulcerans cultures were available from previous studies from Kamerun (n = 21) and Ghana (n = 4) at DITM [8] or were available at KCCR (n = 4) from the present study. All strains were of human origin (BUD patients) and confirmed by conventional IS2404 PCR and sequencing of rpoB- and rpsL-genes that revealed the M. ulcerans Agy99 wild-type sequences (GenBank accession no. CP000325.1) [11,12].

^bResults of the 16S rRNA RT-qPCR of mycobacterial RNA extracts.

^cResults of the IS2404 qPCR of mycobacterial DNA extracts.

^dResults of the IS2404 qPCR of genomic DNA (gDNA) wipeout controls (see Protocols S2 and S3); a positive result indicates gDNA contamination of RNA extracts following DNAse digestions, and a negative result indicates RNA extracts free of gDNA. doi:10.1371/journal.pntd.0001756.t001

Table 2. Specificity of 16S rRNA and IS2404 gPCR assays.

Bacterial Species	Source ^a	Origin ^b	16S rRNA ^d	IS <i>2404</i> e
M. abscessus	NRZ	Human isolate ^p	-	-
M. africanum	NRZ	Human isolate ^p	-	_
M. avium	NRZ	Human isolate ^p	-	-
M. bovis	NRZ	Cattle isolate ^p	-	_
M. chelonae	NRZ	Human isolate ^p	-	-
M. fortuitum	NRZ	Human isolate ^c	-	_
M. gordonae	NRZ	Human isolate ^c	-	-
M. gordonae	DITM	Human isolate ^c	-	_
M. kansasii	NRZ	Human isolate ^p	-	-
M. leprae	DITM	Human isolate ^p	-	_
M. malmoense	NRZ	Human isolate ^c	-	_
M. marinum	NRZ	Human isolate ^p	+	_
M. microti	NRZ	Mouse isolate ^p	-	-
M. scrofulaceum	NRZ	Human isolate ^p	-	-
M. smegmatis	NRZ	Human isolate ^p	-	-
M. szulgai	NRZ	Human isolate ^p	-	-
M. tuberculosis	NRZ	Human isolate ^p	-	-
M. ulcerans	DITM	Human isolate ^p	+	+
M. xenopi	NRZ	Human isolate ^c	-	-
E. coli	MVP	Human isolate ^c	-	_
P. acnes	MVP	Human isolate ^p	-	-
Staph. aureus	MVP	Human isolate ^c	_	_
Staph. epidermidis	MVP	Human isolate ^c	-	-
Str. pyogenes	MVP	Human isolate ^p	_	_

Table 2 shows DNA extracts from closely related mycobacterial species and bacteria potentially contaminating the human skin subjected to the combined 16S rRNA RT/ IS2404 qPCR viability assay and the corresponding test results. Mycobacterial species were selected according to their respective genetic contiguousness to *M. ulcerans* Agy99 (GenBank accession no. CP000325.1) within the 16S rRNA gene sequences as determined by BLASTN analysis (GenBank, NCBI) [13]. *M., Mycobacterium; E., Escherichia; P., Propionibacterium; Staph., Staphylococcus; Str., Streptococcus.* While in-silico analysis revealed that the combined 16S rRNA RT/IS2404 assay will also amplify mycolactone-producing mycobacteria (MPM) other than *M. ulcerans* (e.g., *M. pseudoshottsii, M. liflandii*, and the environmental *M. marinum* [GenBank accession No. NR_042988.1, AY500838.1, and AF456241.1, respectively]), these MPM species were not included in specificity testing.

^aDNA extracts that were not available at the DITM were provided by the National Reference Center (NRZ) for Mycobacteria, Borstel, Germany, and the Max von Pettenkofer-Institute (MVP), Ludwig-Maximilians University, Munich, Germany.

^bThe respective primary patient isolates were considered as ^ppathogenic bacteria or as ^ccommensals/contaminants of clinical samples.

dResults of the 16S rRNA RT-qPCR of DNA extracts; "+" indicates a positive and "-" a negative test result.

 $^{\mathrm{e}}$ Results of the IS2404 qPCR of DNA extracts; "+" indicates a positive and "-" a negative test result.

doi:10.1371/journal.pntd.0001756.t002

Inclusion criteria

- Written informed consent provided
- Participant >5 years of age
- New case (no history of BUD in the past)
- Ulcerative lesion with undermined edge
- Duration of disease < 3 months
- Antibiotic treatment not started

Exclusion criteria

- Denial to participate
- Participant <5 years of age
- Recurrence or secondary lesion
- Ulcer with scarred edge
- Duration of disease > 3 months
- Antibiotic treatment started

Figure 1. Enrolment criteria for the pre-treatment study population. Figure 1 describes enrolment criteria for clinically suspected BUD patients presenting at Agogo Presbyterian Hospital (n = 14) and Tepa Governmental Hospital (n = 10), Ghana, respectively. None of the eligible study participants was excluded. doi:10.1371/journal.pntd.0001756.g001

3

Inclusion criteria

- Written informed consent provided
- Participant >5 years of age
- New case
- Wound has not completely healed
- Antibiotic treatment completed (56 doses)
- No dose missed within 56 days

Exclusion criteria

- Denial to participate
- Participant <5 years of age
- Recurrence or secondary lesion
- Completely healed wound
- Antibiotic treatment not completed
- compliance not ensured

Figure 2. Enrolment criteria for the post-treatment study population. Figure 2 describes enrolment criteria for IS2404 PCR confirmed BUD patients with incomplete wound healing (collection of swab samples feasible) who presented at Agogo Presbyterian Hospital, Ghana (n=7), following completion of 56 doses of rifampicin and streptomycin administered within eight weeks. None of the eligible study participants was excluded.

doi:10.1371/journal.pntd.0001756.g002

both qPCR components using 10-fold serial dilutions of cloned IS2404 templates (GenExpress, Berlin, Germany) with known copy numbers (IS2404 qPCR) and exactly quantified *M. ulcerans* whole genome DNA extracts from cultures (16S rRNA RT-qPCR). The LOD was two (IS2404) and six templates (16S rRNA gene), respectively (Figures 3 and 4).

M. ulcerans DNA and rRNA was detected in all culture extracts. Out of 24 pre-treatment swab samples, 18 (75.0%; 95%-CI: 57.7%-92.3%) had a positive IS2404 qPCR result, 12 out of those were also positive in routine DRB PCR, and rRNA was detected in 15 out of these 18 samples (83.3%; 95%-CI: 66.1%-100%); quantification of the three negative samples revealed a bacillary load below the LOD of the 16S rRNA RT-qPCR (Table 6).

All seven post-treatment swab samples were IS2404 qPCR positive and 16S rRNA negative.

Specificity

Analysis of DNA extracts revealed 100% specificity for the combined assay. *M. marinum* (human isolate) was amplified by 16S rRNA RT-qPCR; however, simultaneous IS2404 qPCR was negative (Table 2).

Bacillary Survival Times

To investigate the effect of sample transport on bacillary survival, mycobacteriological transport media (PANTA and LTM) [3] were spiked with viable M. ulcerans and stored at 4° C and 31° C. RNA was detectable in both media for >4 weeks $(4^{\circ}$ C and 31° C).

After heat-inactivation of *M. ulcerans*—spiked PANTA-samples, RNA positivity

decreased significantly within 12 h, whereas DNA was still detectable after seven days.

Future Application

The assay will support clinicians in classification of secondary lesions and selection of adequate clinical management strategies and provides a powerful tool for clinical research evaluating novel treatment regimens (Box 1).

Through analysis of sequential samples collected during antimycobacterial treatment, the assay will be employed to determine the proportional decrease of bacterial viability over time and to establish laboratory-based evidence for optimal time-points to collect follow-up samples for treatment monitoring.

Whereas the current format of the assay is restricted to reference laboratories,

Table 3. Primers and probes.

Primer/Probe ^a	Sequence (5'-3')	Target Gene ^b	Nucleotide Position ^c	Amplicon Size ^d
MU16S TF MU16S TR MU16S TP	CGA TCT GCC CTG CAC TTC CCA CAC CGC AAA AGC TT 6 FAM-CAC AGG ACA TGA ATC CCG TGG TC-BBQ ^e	16S rRNA	4414800–4414817 4414718–4414734 4414740–4414762	100 bp
IS2404 TF IS2404 TR IS2404 TP2	AAA GCA CCA CGC AGC ATC T AGC GAC CCC AGT GGA TTG 6 FAM-CCG TCC AAC GCG ATC GGC A-BBQ°	IS2404	96685–96667 96627–96644 96664–96646	59 bp
T13 ^f T39 ^f	TGC ACA CAG GCC ACA AGG GA CG AAC GGG TGA GTA ACA CG	16S rRNA	4413906–4413925 4414822–4414840	935 bp

Table 3 indicates primers and probes designed for the 16S rRNA RT-qPCR, the primers described by Fyfe et al., and a re-designed hydrolysis probe used for the amplification, detection, and quantification of IS2404 [9].

^fPrimers T13 (forward) and T39 (reverse) were used for the amplification of a 935-bp region of the *M. ulcerans* 16S rRNA gene, encompassing the region amplified by qPCR primers MU16S TF and MU16S TR, to generate single copy replicates. Furthermore, these primers were used for sequencing of the *M. ulcerans* 16S rRNA gene (Table 1).

doi:10.1371/journal.pntd.0001756.t003

^aTF, forward primer; TR, reverse primer; TP2, hydrolysis probe (TibMolBiol, Berlin, Germany).

^b16S rRNA, gene for the ribosomal 16S RNA detected as 16S cDNA; IS2404, insertion sequence 2404.

^cNucleotide positions are provided for the first (IS2404) or single (16S rRNA) copy of the respective amplicon in *M. ulcerans* Agy99 (GenBank accession no. CP000325.1) as determined by BLASTN analysis within GenBank (NCBI) [13].

dbp, base pairs.

^e6 FAM, 6-Carboxyfluorescein fluorescent dye; BBQ, BlackBerry Quencher.

Table 4. Controls applied in 16S rRNA RT/IS2404 gPCR.

Control	Purpose	Material				
		16S rRNA RT-qPCR ^a	IS <i>2404</i> qPCR ^b			
gDNA wipeout control ^c	To exclude DNA contamination of RNA extracts	Aliquot of each RNA extract following gDNA wipeout before reverse transcription	NA			
Internal positive control	To exclude false negative results due to inhibition	TaqMan exogenous internal positive control (IPC) ^d	TaqMan exogenous internal positive control (IPC) ^d			
Positive run control	To ensure adequate performance of PCR	M. ulcerans cDNA ^e	Cloned IS2404 standard			
Negative no template control	To exclude contamination during PCR set up	H ₂ O	H ₂ O			
Negative extraction control	To exclude contamination during extraction procedure	NA	Extract treated in the same way as samples			

Table 4 indicates controls applied in 16S rRNA RT/IS2404 qPCR. NA, not applicable.

Table 5. Intra- and inter-assay variability of the 16S rRNA RT/IS2404 gPCR assay.

qPCR Target ^a	Standard No.	Run No. 1		Run No. 2		Run No. 3		Intra-Assay Variability		Inter-Assay Variability			
		Ct-range ^b	CV ^c	Ct-range ^b	CV^c	Ct-range ^b	CV^c	Δ Ct max. ^d	CV max.e	Ct-range ^f	CN_a	∆Ct max.h	CV max.
16S rRNA	1	0.23	0.50	0.12	0.48	0.17	0.42	0.23	0.49	0.55	1.33	0.75	1.33
	2	0.09	0.19	0.16	0.30	0.19	0.35			0.24	0.53		
	3	0.12	0.18	0.06	0.20	0.20	0.32			0.31	0.55		
	4	0.15	0.22	0.17	0.25	0.12	0.22			0.75	1.15		
	5	0.07	0.10	0.15	0.20	0.16	0.20			0.71	0.92		
IS2404	1	0.12	0.53	0.13	0.54	0.10	0.42	0.35	0.65	0.61	2.67	0.80	2.66
	2	0.18	0.65	0.15	0.48	0.18	0.57			0.71	2.35		
	3	0.02	0.07	0.23	0.60	0.11	0.28			0.80	2.13		
	4	0.18	0.39	0.14	0.28	0.10	0.22			0.80	1.76		
	5	0.31	0.58	0.25	0.42	0.22	0.38			0.58	1.09		
	6	0.15	0.23	0.31	0.47	0.20	0.32			0.31	0.58		
	7	0.35	0.48	0.15	0.33	0.08	0.29			0.74	1.10		

Table 5 shows intra- and inter-assay variability of the 16S rRNA RT/IS2404 qPCR assay.

16S rRNA RT-qPCR: 16S rRNA gene standards (935 bp) were generated by conventional PCR according to Talaat et al. [12]. Quantification of PCR products was conducted by Picogreen fluorometry (Invitrogen) and copy numbers were calculated based on the known mass of one amplicon. Serial standards were prepared from PCR products in 5 Log dilutions ranging from 3E+6 (standard no. 1) to 3E+2 copies (standard no. 5) of the 16S rRNA amplicon (PCR template: 2 µl) and were subjected to the 16S rRNA RT-qPCR in quadruplicate on one 96-well plate to assess intra-assay variability. The runs were repeated on three days to determine the inter-assay variability between runs 1 through 3. The intra- and inter-assay variability of the 16S rRNA RT-qPCR was low with maximum coefficients of variation (CV) of 0.49 (intra-assay) and 1.33 (inter-assay).

IS2404 qPCR: Cloned IS2404 replicates (1,047 bp, complete sequence; *M. ulcerans* Agy99) were used as standards. Quantification of IS2404 templates was conducted by Picogreen fluorometry (Invitrogen) and copy numbers were calculated based on the known mass of one template. Serial standards were prepared in 7 Log dilutions ranging from 2E+8 to 2E+2 copies of the IS2404 (PCR template: 2 µl) and were subjected to the IS2404 qPCR in quadruplicate on one 96-well plate to assess intra-assay variability. The runs were repeated on three days to determine the inter-assay variability between runs 1 through 3. The intra- and inter-assay variability of the IS2404 qPCR was low with maximum CV of 0.65 (intra-assay) and 2.66 (inter-assay).

^a16S rRNA RT PCR, reverse transcriptase real-time PCR targeting the 16S ribosomal RNA of *M. ulcerans*.

^bIS2404 qPCR, quantitative real-time PCR targeting the insertion sequence (IS) 2404 of M. ulcerans.

^cgDNA, genomic DNA wipeout was conducted using DNAses provided in the QuantiTect Reverse Transcription Kit (Qiagen).

^dTaqMan exogenous internal positive control (Applied Biosystems, Carlsbad, CA).

^ecDNA, complementary DNA obtained through reverse transcription of *M. ulcerans* RNA by QuantiTect Reverse Transcription Kit (Qiagen).

doi:10.1371/journal.pntd.0001756.t004

^a16S rRNA, target of the 16S rRNA RT-qPCR; IS2404, target of the IS2404 qPCR.

^bCt-range, range of Ct-values of samples tested in the same dilution.

^cCV, coefficient of variation of copy numbers from samples tested in quadruplicate of the same dilution.

 $^{^{\}text{d}}\Delta\text{Ct}$ max., maximum Ct-variation of all samples tested within one run.

^eCV max., maximum CV of all samples tested within one run.

^fCt-range, range of Ct-values of samples tested in the same dilution within three runs.

⁹CV of samples in the same dilution tested within three runs.

^h∆Ct max. of all samples tested within three runs.

ⁱCV max. of all samples tested within three runs.

doi:10.1371/journal.pntd.0001756.t005

Standard curve and limit of detection of the 16S rRNA RT-qPCR

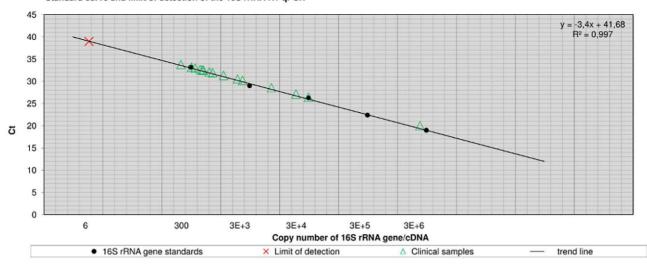


Figure 3. Standard curve and limit of detection of the 16S rRNA RT-qPCR. Figure 3 shows Ct-values of clinical samples plotted versus quantified 16S rRNA copy numbers. Standards for the 16S rRNA RT-qPCR were generated by conventional PCR amplification (Table 5). Log 10 fold serial dilutions (n = 5) were prepared ranging from 3E+6 to 300 copies of the 16S rRNA gene (PCR template: 2 μ I) and were subjected to the assay in quadruplicate to generate a calibration curve. The regression line was y = -3.4x + 41.68 with a coefficient of correlation > 0.99 and the efficiency was E = 0.97. *M. ulcerans* whole genome extracts were quantified by means of IS2404 qPCR and the analytical sensitivity was determined as limit of detection (LOD) by subjecting 10 aliquots of a dilution series containing 30, 15, 10, 8, 6, 3, or 2 copies of the 16S rRNA gene to the assay. The LOD was 6 copies of the target sequence. The copy number (n = 1) of the 16S rRNA gene per *M. ulcerans* genome was determined by copy number variation assay (unpublished data).

doi:10.1371/journal.pntd.0001756.g003

Standard curve and limit of detection of the IS2404 qPCR

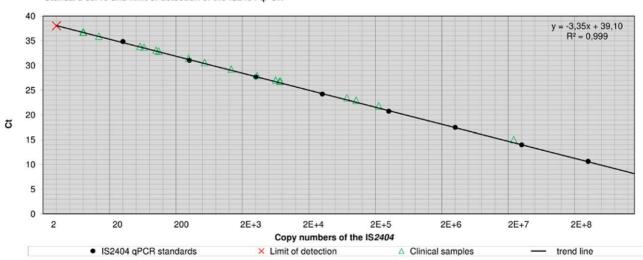


Figure 4. Standard curve and limit of detection of the IS 2404 **qPCR.** Figure 4 shows mean Ct-values of calibration standards and clinical samples plotted versus the quantified copy number of IS2404. Cloned IS2404 templates were used as standards (Table 5). Log 10 fold serial dilutions (n = 8) were prepared ranging from 2E+8 to 20 copies of the IS2404 (PCR template: 2 μ I) and were subjected to the IS2404 qPCR in quadruplicate to generate a calibration curve. The regression line was y = -3.35x+39.10 with a coefficient of correlation >0.99 and the efficiency was E = 0.97. The analytical sensitivity was determined as limit of detection (LOD) by subjecting 10 aliquots of a dilution series containing 10, 5, 4, 3, 2, or 1 copy of the IS2404 to the assay. The LOD was 2 copies of the target sequence. doi:10.1371/journal.pntd.0001756.g004

Table 6. Study participants, clinical information, and diagnostic results.

Clinical Data			Molecular Viabil	Routine Diagnostics ^b				
No. ^c	BUD Patient ^d	Duration (Weeks) ^e	Category of Lesion ^f	IS 2404 [Ct] ⁹	Bacillary Load ^h	16S rRNA ⁱ	MICk	PCRI
1	No	NA	NA	Neg [NA]	NA	Neg	0	Neg
2	Yes	6	III	Pos [15,04]	>1000	Pos	+1	Pos
3	Yes	4	III	Pos [26,80]	584	Pos	+1	Pos
4	Yes	9	III	Pos [32,93]	6–10	Pos	0	Neg
5	Yes	4	1	Pos [35,94]	1–5	Neg	0	Neg
6	Yes	8	II	Pos [36,72]	1–5	Neg	0	Neg
7	Yes	2	1	Pos [36,74]	1–5	Neg	0	Neg
8	Yes	10	1	Pos [27,05]	497	Pos	+1	Pos
9	No	NA	NA	Neg [NA]	NA	Neg	0	Neg
10	Yes	3	I	Pos [30,61]	42	Pos	+1	Pos
11	Yes	8	II	Pos [33,89]	6–10	Pos	0	Neg
12	Yes	9	1	Pos [33,68]	6–10	Pos	0	Neg
13	Yes	3	III	Pos [29,27]	106	Pos	+1	Pos
14	Yes	3	I	Pos [27,98]	261	Pos	+1	Pos
15	Yes	1	1	Pos [26,85]	571	Pos	+1	Pos
16	Yes	2	I	Pos [33,07]	6–10	Pos	0	Pos
17	Yes	2	II	Pos [31,44]	24	Pos	+1	Pos
18	Yes	3	II	Pos [21,85]	>1000	Pos	+2	Pos
19	Yes	4	III	Pos [22,98]	>1000	Pos	+1	Pos
20	Yes	3	I	Pos [23,47]	>1000	Pos	+2	Pos
21	No	NA	NA	Neg [NA]	NA	Neg	0	Neg
22	No	NA	NA	Neg [NA]	NA	Neg	0	Neg
23	No	NA	NA	Neg [NA]	NA	Neg	0	Neg
24	No	NA	NA	Neg [NA]	NA	Neg	0	Neg

Table 6 shows suspected BUD cases with ulcerative lesions enrolled in the pre-treatment cohort (Figure 1), clinical information, and diagnostic results. Swab samples from 24 suspected BUD cases were subjected to 165 rRNA RT/IS2404 qPCR viability assay (swab 1 in PANTA), microscopic examination and enumeration of acid fast bacilli (AFB) following Ziehl-Neelsen staining (swab 2, direct smear), and conventional IS2404 dry-reagent-based (DRB) PCR (swab 3 in Cell Lysis Solution [Qiagen]). 18 patients were laboratory confirmed by IS2404 qPCR and 15 out of those were RNA positive; the quantification by IS2404 qPCR revealed a bacillary load (1–2 bacilli per sample) below the lower limit of detection of the RNA assay for samples from three RNA negative patients. All samples from six IS2404 qPCR negative study participants were also RNA negative. Direct correlation of AFB enumeration with IS2404 qPCR quantification is not feasible due to inhomogeneous distribution of M. ulcerans in different clinical samples. NA, not applicable; Neg, negative test result; Pos, positive test result.

^aResults of the 16S rRNA RT/IS2404 qPCR viability assay. Clinical swab samples in PANTA were directly processed at KCCR, and *M. ulcerans* DNA and cDNA were transported to DITM and subjected to qPCR.

doi:10.1371/journal.pntd.0001756.t006

^bRoutine diagnostics were conducted following standardized procedures at KCCR [3].

^cNo., consecutive number of study participants.

^dYes, IS2404 qPCR confirmed BUD patients; No, IS2404 negative study participants.

^eDuration of disease before presentation of study participants in weeks.

^fCategory of lesion according to the World Health Organization's clinical criteria [1].

⁹Results of the IS2404 qPCR with corresponding cycle threshold (Ct)-values.

^hThe bacillary load in the respective swab samples (No. 2) was estimated on the basis of IS2404 quantification given an IS2404 copy number of 209 copies per *M. ulcerans* genome [9]. For bacterial numbers <10 ranges were estimated.

Results of the 16S rRNA RT-qPCR.

^{*}MIC, microscopic detection and enumeration of AFB was conducted at KCCR including external quality assurance by DITM. The following scale was applied: 0 = negative, +1 = 10 - 99 AFB/100 fields, +2 = 1 - 10 AFB/1 field, +3 = more than 10 AFB/1 field.

¹PCR, conventional, single target gel-based IS2404 DRB PCR.

Box 1. Advantages and Disadvantages of the Molecular Viability Assay

Advantages

- Provides a rapid, sensitive, and specific tool to detect viable bacilli in clinical samples of BUD patients, thus offering an alternative to cultures.
- Supports classification of secondary BUD lesions and monitoring of treatment success.

Disadvantages

- Current test format requires well equipped laboratory with real-time PCR facilities.
- Costs per test (approximately 14 €) may limit the applicability.

sample collection on FTA cards in combination with isothermal dry-reagent-based reverse transcription and amplification formats would facilitate processing of samples also at a peripheral level and at lower costs.

Conclusions

The novel combined 16S rRNA RT/IS2404 qPCR assay proved to be highly sensitive, specific, and efficient in detecting viable *M. ulcerans* in clinical samples under field conditions. The assay is applicable for classification of secondary lesions and

monitoring of treatment success and provides a powerful tool for clinical research.

GenBank Accession Numbers

Genes or DNA sequences of mycobacterial strains used in this study were retrieved from GenBank (NCBI) [13]. The respective sequences and accession numbers are summarized in Table S1.

Supporting Information

Protocol S1 Preparation of PANTA transport medium and stabilisation of M.

vides a powerful tool for clinical research.

- clinical and environmental samples. Appl Environ
- Microbiol 73(15): 4733–4740.

 10. Taylor S, Wakem M, Dijkman G, Alsarraj M, Nguyen M (2010) A practical approach to RT-qPCR-Publishing data that conform to the MIQE
- qPCR-Publishing data that conform to the MIQE guidelines. Methods 50(4): S1–S5.

 11. Beissner M, Awua-Boateng NY, Thompson W, Nienhuis WA, Klutse E, et al. (2010) A genotypic approach for detection, identification, and char

acterization of drug resistance in Mycobacterium

ulcerans in clinical samples and isolates from

- Ghana. Am J Trop Med Hyg 83(5): 1059–1065.

 12. Talaat AM, Reimschuessel R, Trucksis M. (1997)
 Identification of mycobacteria infecting fish to the species level using polymerase chain reaction and restriction enzyme analysis. Vet Microbiol 58(2–4): 200–227.
- Benson DA, Karsch-Mizrachi I, Lipman DJ, Ostell J, Wheeler DL (2008) GenBank. Nucleic Acids Res 36 (Database issue): D25–D30.

References

- World Health Organization (2008) Buruli ulcer: progress report, 2004–2008. Wkly Epidemiol Rec 83: 145–156.
- World Health Organization (2008) Meeting of the WHO Technical Advisory Group on Buruli ulcer, 3 April 2008, Geneva, summary report. Geneva: World Health Organization.
- Beissner M, Herbinger KH, Bretzel G (2010)
 Laboratory diagnosis of Buruli ulcer disease.
 Future Microbiol 5(3): 363–370.
- Ruf MT, Chauty A, Adeye A, Ardant MF, Koussemou H, et al. (2011) Secondary Buruli ulcer skin lesions emerging several months after completion of chemotherapy: paradoxical reaction or evidence for immune protection? PLoS Negl Trop Dis 5(8): e1252. doi:10.1371/journal. pntd.0001252
- Nienhuis WA, Stienstra Y, Abass KM, Tuah W, Thompson WA, et al. (2011) Paradoxical responses after start of antimicrobial treatment in

- Mycobacterium ulcerans infection. Clin Infect Dis 54(4): 519–526.
- Martinez AN, Lahiri R, Pittman TL, Scollard D, Truman R, et al. (2009) Molecular determination of Mycobacterium leprae viability by use of realtime PCR. J Clin Microbiol 47(7): 2124–2130.
- Desjardin LE, Perkins MD, Wolski K, Haun S, Teixeira L, et al. (1999) Measurement of sputum Mycobacterium tuberculosis messenger RNA as a surrogate for response to chemotherapy. Am J Respir Crit Care Med 160: 203–210.
- Bretzel G, Siegmund V, Racz P, van Vloten F, Ngos F, et al. (2005) Post-surgical assessment of excised tissue from patients with Buruli ulcer disease: progression of infection in macroscopically healthy tissue. Trop Med Int Health 10(11): 1199–1206.
- Fyfe JA, Lavender CJ, Johnson PD, Globan M, Sievers A, et al. (2007) Development and application of two multiplex real-time PCR assays for the detection of Mycobacterium ulcerans in

ulcerans RNA/DNA in swab samples and culture suspensions.
(PDF)

Protocol S2 Simultaneous RNA/DNA extraction from swab samples and reverse transcription of whole transcriptome RNA from *M. ulcerans*. (PDF)

Protocol S3 Combined 16S rRNA RT/IS2404 qPCR assay. (PDF)

Protocol S4 16S rRNA RT/IS2404 qPCR run protocol. (XLS)

Table S1 GenBank accession numbers. (DOC)

Acknowledgments

The authors thank Erna Fleischmann, Carolin Mengele and Kerstin Helfrich (DITM), and Mabel Peprah and Michael Frimpong (KCCR) for excellent technical assistance. The authors thank Dr. Sabine Rüsch-Gerdes and Dr. Elvira Richter (National Reference Center for Mycobacteria, Borstel, Germany), as well as Dr. Soeren Schubert (Max von Pettenkofer-Institute, Ludwig-Maximilians University, Munich, Germany) for providing (myco-) bacterial DNA extracts. The manuscript contains parts of the doctoral thesis of Dominik Symank.