

IOWA STATE UNIVERSITY
Digital Press

Animal Industry Report

Animal Industry Report

AS 663

ASL R3130

2017

Desmin and Peroxiredoxin-2 are Potential Biomarkers for Pork Tenderness

Kelsey Carlson
Iowa State University

Kenneth Prusa
Iowa State University, kprusa@iastate.edu

Chris Fedler
Iowa State University, cfedler@iastate.edu

Ed Steadham
Iowa State University, steadham@iastate.edu

Elisabeth Huff-Lonergan
Iowa State University, elonerga@iastate.edu

See next page for additional authors

Follow this and additional works at: https://lib.dr.iastate.edu/ans_air

 Part of the **Meat Science Commons**

Recommended Citation

Carlson, Kelsey; Prusa, Kenneth; Fedler, Chris; Steadham, Ed; Huff-Lonergan, Elisabeth; and Lonergan, Steven (2017) "Desmin and Peroxiredoxin-2 are Potential Biomarkers for Pork Tenderness," *Animal Industry Report*: AS 663, ASL R3130.

DOI: https://doi.org/10.31274/ans_air-180814-313

Available at: https://lib.dr.iastate.edu/ans_air/vol663/iss1/3

This Animal Products is brought to you for free and open access by the Animal Science Research Reports at Iowa State University Digital Repository. It has been accepted for inclusion in Animal Industry Report by an authorized editor of Iowa State University Digital Repository. For more information, please contact digirep@iastate.edu.

Desmin and Peroxiredoxin-2 are Potential Biomarkers for Pork Tenderness

Authors

Kelsey Carlson, Kenneth Prusa, Chris Fedler, Ed Steadham, Elisabeth Huff-Lonergan, and Steven Lonergan

Desmin and Peroxiredoxin-2 are Potential Biomarkers for Pork Tenderness

A.S. Leaflet R3130

Kelsey Carlson, Graduate Student, Department of Animal Science; Kenneth Prusa, Professor, Department of Food Science and Human Nutrition; Chris Fedler, Research Associate, Department of Food Science and Human Nutrition; Ed Steadham, Research Associate, Department of Animal Science; Elisabeth Huff-Lonergan, Professor, Department of Animal Science; Steven Lonergan, Professor, Department of Animal Science

Summary and Implications

Aged pork loins were selected to have similar ultimate pH, color, lipid content, and purge. The star probe values (kg) obtained from these loins were evaluated and loins were sorted into either a low star probe (LSP) group or a high star probe (HSP) group. Two-dimensional difference in gel electrophoresis (2D-DIGE) was used to identify differences in sarcoplasmic proteomes between LSP and HSP groups. Metabolic, stress response, structural, and regulatory protein spots were identified to be significantly different between classification groups. The structural protein desmin was identified to be more abundant in LSP samples and the stress response protein peroxiredoxin-2 was identified to be more abundant in HSP samples. Because these proteins were identified as being significantly different between classification groups in the soluble, sarcoplasmic fraction of samples, these proteins are potential candidates to be biomarkers to differentiate between tough and tender aged pork loins.

Introduction

Consumers, purveyors, and chefs continue to be presented with inconsistent fresh pork quality, specifically pork tenderness. Currently, there is no method to rapidly differentiate between tough and tender fresh pork in processing facilities. Pork loins that exhibit differences in measured tenderness have shown differences in abundance of metabolic, stress response, structural, and regulatory proteins. Rapidly identifying proteins that differ between tough and tender pork is fundamental to identify biomarkers that can be used to differentiate between tough and tender fresh pork products.

Materials and Methods

Commercial pork loins ($n = 159$) were collected from Duroc-sired crossbred pigs and aged for 9 to 11 days. Chops (2.54 cm thick) were collected from loins and data were collected for ultimate pH, visual color and marbling, Hunter L, a, b color, sensory, star probe (kg), and total lipid. Loins

that had extreme star probe measurements were selected if they were also within the following ranges for the following traits: ultimate pH (5.54 – 5.86), visual marbling score (1.0 – 3.0), and total lipid (1.61 – 3.37%). These ranges were selected because previous research showed ultimate pH and lipid content can influence measured tenderness in fresh pork products. Selected samples were grouped into either a LSP ($n = 12$) or HSP ($n = 12$) group.

Sarcoplasmic protein samples were prepared from the selected loins to be used for 2D-DIGE. Proteins are separated by charge and size using this method. In order to have a direct in-gel comparison between LSP and HSP samples, samples were labeled with fluorescent dyes that excite at different wavelengths. DeCyder statistical software (version 6.5) was used to identify significantly different spots between classification groups. Identities of these protein spots was determined using mass spectrometry. In order to confirm the identity of the proteins from 2D-DIGE and mass spectrometry, sarcoplasmic protein samples were prepared for SDS-PAGE and Western blot analysis. Samples from both classification groups were run in duplicate for 2D-DIGE, SDS-PAGE, and Western blots. Data for SDS-PAGE and Western blots were analyzed using the mixed procedure in SAS (version 9.4), with a fixed effect of star probe force and random effect of gel.

Results and Discussion

Samples from both classification groups had similar ultimate pH, color, and loin/chop purge. The average star probe for the LSP samples was 4.95 kg and the average star probe force for the HSP samples was 7.75 kg. These values demonstrate there was an extreme difference in instrumental tenderness between classification groups. There was a total of 102 protein spots significantly different between LSP and HSP groups in the sarcoplasmic protein fraction ($P < 0.01$). Of these spots, 16 of the most prevalent were picked for identification. The protein spots identified included metabolic, stress response, structural, and regulatory proteins.

The structural protein desmin was found to be 90% more abundant in LSP samples (Figure 1). After further analysis of the peptides identified from this protein spot, it was found that only a fragment of desmin (rod portion) was identified. This demonstrates desmin was potentially degraded into the sarcoplasmic protein fraction. This is a novel result because desmin is typically thought to be associated with the insoluble myofibrillar protein fraction due to the role it plays with structure of myofibrils. The stress response protein peroxiredoxin-2 (two spots identified) was found to be 37% and 51% more abundant in

HSP samples (Figure 1). Peroxiredoxin-2 plays a role in protecting cells from oxidative stress.

Protein spot identities were confirmed with SDS-PAGE and Western blots. Western blots confirmed the presence, size, and abundance of these proteins in the sarcoplasmic protein fraction (Figure 2). The sarcoplasmic protein fraction is soluble and is found in meat purge. The

....

ease of solubility of these proteins, coupled with the extreme differences in abundance between LSP and HSP samples, demonstrates the potential for use of desmin and peroxiredoxin-2 as biomarkers to differentiate between tough and tender pork samples.

Acknowledgements

Figure 1. Representative 2D-DIGE gel of the 16 protein spots identified from the sarcoplasmic protein fraction of aged pork loins. Desmin and Peroxiredoxin-2 are labeled.

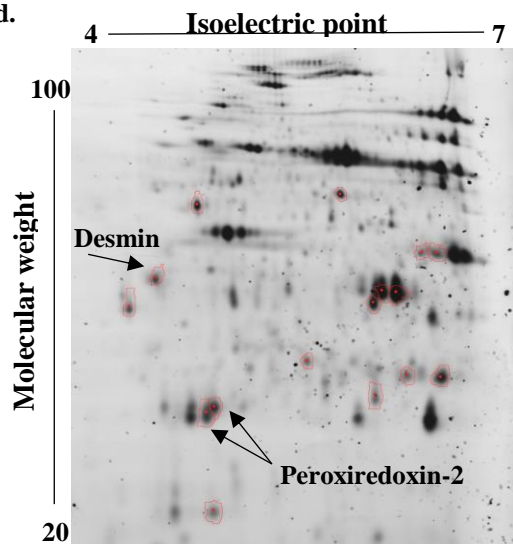


Figure 2. Representative Western blots images of desmin and peroxiredoxin-2 from the sarcoplasmic protein fraction of aged pork loins. Western blots were used to confirm protein identity, size, and abundance from protein spots identified from 2D-DIGE and mass spectrometry. Low star probe (LSP) and high star probe (HSP) samples are labeled.

