



DiSciacca, Jack and Gomez, Carlos and Thompson, Adam and Lawes, Simon and Leach, Richard and Colonna de Lega, Xavier and de Groot, Peter (2017) True-color 3D surface metrology for additive manufacturing using interference microscopy. In: Joint Special Interest Group meeting between euspen and ASPE Dimensional Accuracy and Surface Finish in Additive Manufacturing, 10-12 October 2017, Leuven, Belgium.

Access from the University of Nottingham repository:

<http://eprints.nottingham.ac.uk/48598/1/True-color%203D%20surface%20metrology%20for%20additive%20manufacturing%20using%20interference.pdf>

Copyright and reuse:

The Nottingham ePrints service makes this work by researchers of the University of Nottingham available open access under the following conditions.

This article is made available under the University of Nottingham End User licence and may be reused according to the conditions of the licence. For more details see: http://eprints.nottingham.ac.uk/end_user_agreement.pdf

A note on versions:

The version presented here may differ from the published version or from the version of record. If you wish to cite this item you are advised to consult the publisher's version. Please see the repository url above for details on accessing the published version and note that access may require a subscription.

For more information, please contact eprints@nottingham.ac.uk

True-color 3D surface metrology for additive manufacturing using interference microscopy

Jack DiSciacca¹, Carlos Gomez², Adam Thompson², Simon Lawes², Richard Leach², Xavier Colonna de Lega¹, Peter de Groot¹

¹Zygo Corporation, USA

²Manufacturing Metrology Team, Faculty of Engineering, University of Nottingham, NG7 2RD, UK

jack.disciacca@ametec.com

Abstract

Coherence scanning interferometry (CSI) is widely used for surface topography characterization. With the ability to measure both rough surfaces with high slopes and optical finishes, CSI has made contributions in fields from industrial machining to optical fabrication and polishing [1,2]. While the low coherence sources for CSI are typically broadband and suitable for color imaging, the metrology is usually performed without regard for the color information [3]. We present color surface topography measurements from a CSI instrument designed to provide true color images in addition to areal surface topography of additive manufactured samples. The addition of color measurements enables contamination and defect detection along with blemish and discoloration identification.

Coherence scanning interferometry, metal additive manufacturing, surface topography, metrology, true color

1. Introduction

Additive manufacturing (AM) holds substantial promise for the fields of engineering and production. By providing the ability to selectively add material to a part as it is formed, the technique can produce novel part geometries that were previously beyond the capabilities of conventional, subtractive manufacturing. The new geometries and parts have applications in fields from biomedical implants to aerospace engineering [4,5]. However, while there are substantial advantages with AM, there are open challenges in improving the process. One such challenge is the ability to accurately assess the surface finish and identify defects; both of which are important areas for improving AM [6]. For this work we present novel measurements of AM samples that combine true color with surface topography. While foundational work has been done demonstrating the efficacy of coherence scanning interferometry (CSI) measurements of surface topography of additive manufactured samples [7], we present the first color topographies. These provide examples of discoloration and defects, providing more complete metrology feedback to the manufacturing process.

In the sections below, we describe the methodology used to acquire the measurements, along with the analysis and acquisition. We also describe the samples and the results from the measurements and possible future directions.

2. Methodology

2.1. Measurement System

We used a ZYGO Nexview™ system to generate the true color surface topographies. This microscope maps one million height points in a few seconds using CSI and specialized hardware to synthesize color images [1]. This instrument is capable of measuring topographies of highly complex surfaces

characteristic of AM samples, overcoming many of the perceived limitations of interferometry for these surfaces [7-9]. The measurements relied on a 50× objective and 0.5× tube lens to resolve the sample defects and discoloration, and also to capture data from the high slope and rough surfaces. The single measurement field of view (FOV) with this configuration was 0.34 mm by 0.34 mm, with a spatial sampling of 0.33 μm and an optical resolution of 0.55 μm (given by the Sparrow criteria at an imaging wavelength of 570 nm).

2.2. Acquisition and Analysis

The main principle behind a coherence scanning optical microscope is the finite coherence length of a broad illumination spectrum. The light that illuminates the sample has a mean wavelength of about 570 nm and a bandwidth of 100 nm. This broad source bandwidth translates to a coherence length of a few micrometers. To measure the topography of a sample, the objective is scanned along an axis orthogonal to the nominal plane of the surface, and analysis of the coherence-limited interference fringes at each image pixel provides the height topography. Example measurements of AM samples are illustrated in Figures 4a and 5a.

In addition to surface topography, true color measurements are possible using CSI. As noted above, the 100 nm wavelength bandwidth covers much of the visible spectrum. Sequential illumination of the sample for each color channel (red, green and blue) provides the reflected intensity for these composite colors. The color information combines with the surface topography to give the true color results at best focus [3].

2.3. Samples

The two samples are respectively made of Ti-6Al-4V and Al-Si-10Mg. These materials both have good strength-to-weight ratios and high resistance to fatigue and corrosion [10,11]. This leads to possible application in both automotive and aerospace

manufacturing [10, 11]. Furthermore, Ti-6Al-4V is a biocompatible material which enables biomedical applications [12]. These materials were selected to build two samples using laser powder bed fusion (LPBF). The Al-Si-10Mg and Ti-6Al-4V LPBF samples are shown in Figure 1a and 1b; both are cubes with a 20 mm lateral extent.

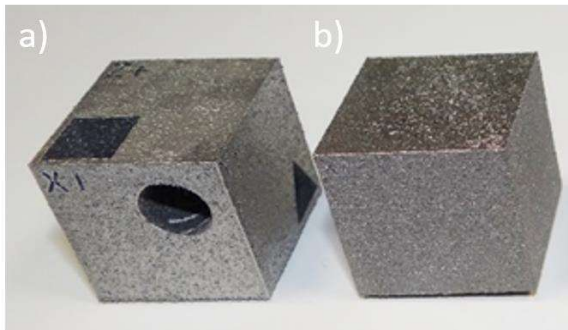


Figure 1. Pictures of the samples: (a) Al-Si-10Mg, (b) Ti-6Al-4V LPBF cubes. Both samples are 20 mm in lateral extent.

3. Results

The results for the true color measurements of the Ti-6Al-4V and Al-Si-10Mg LPBF samples are shown in figures 2 and 3. In figure 2, the data for the Ti-6Al-4V sample shows clear evidence of weld tracks and ripples [7-9]. These represent both micrometer level structures as well as more global features at the hundred micrometer scale. One of the benefits of color measurements is the clear discoloration as the weld ripples are tightly spaced, indicating a possible signal to feedback on the manufacturing process. In addition to the weld ripples, there are also several surface pores and protrusions from spatter or unmelted particles [7-9]. These are well resolved and can also help inform the manufacturing process.

Figure 3 shows the color measurement of the Al-Si-10Mg sample. In contrast to the Ti-6Al-4V sample, this surface is far less ordered and exhibits very little evidence of spatter or pore defects. The overall color is more uniform, with differences largely from surface roughness and topography. Furthermore, there is no evidence of discoloration or localized blemishes.

In addition to the single FOV measurements presented in figures 4 and 5, we also present the combination of multiple measurements along the sample surface to create larger scale surface topographies. The results of both the surface topography and the true color measurements are displayed in figures 4 and 5. These are useful for larger sampling statistics in the context of defect identification and blemish inspection.

4. Summary

The challenges of metrology for AM samples come from the diverse surface topographies and structures. In addition to high slopes and large departures, it is necessary to capture additional information from blemishes and discoloration to help improve and control the manufacturing process. The combination of high-quality interferometric topographic measurements with true color imaging provides a unique metrology opportunity for AM materials.

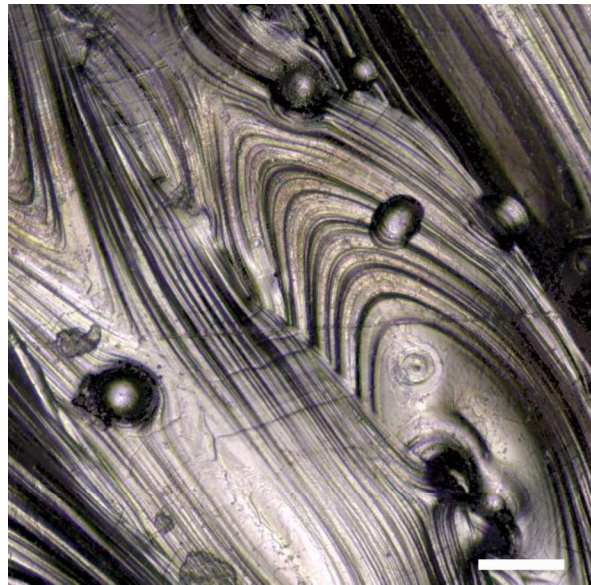


Figure 2. True color measurement for the Ti-6Al-4V LPBF sample. The measurement represents a single field of view for a 50x objective lens (NA 0.55, 0.34 mm × 0.34 mm field of view). The scale bar is 50 μm.

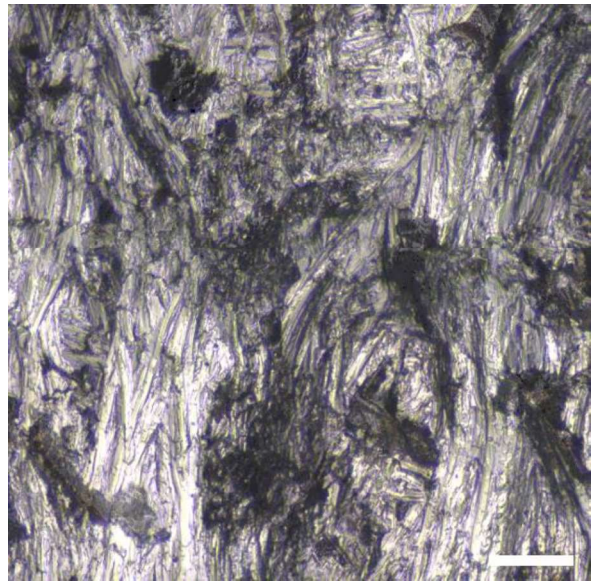


Figure 3. True color measurement for the Al-Si-10Mg LPBF sample. The measurement represents a single field of view for a 50x objective lens (NA 0.55, 0.34 mm × 0.34 mm field of view). The scale bar is 50 μm.

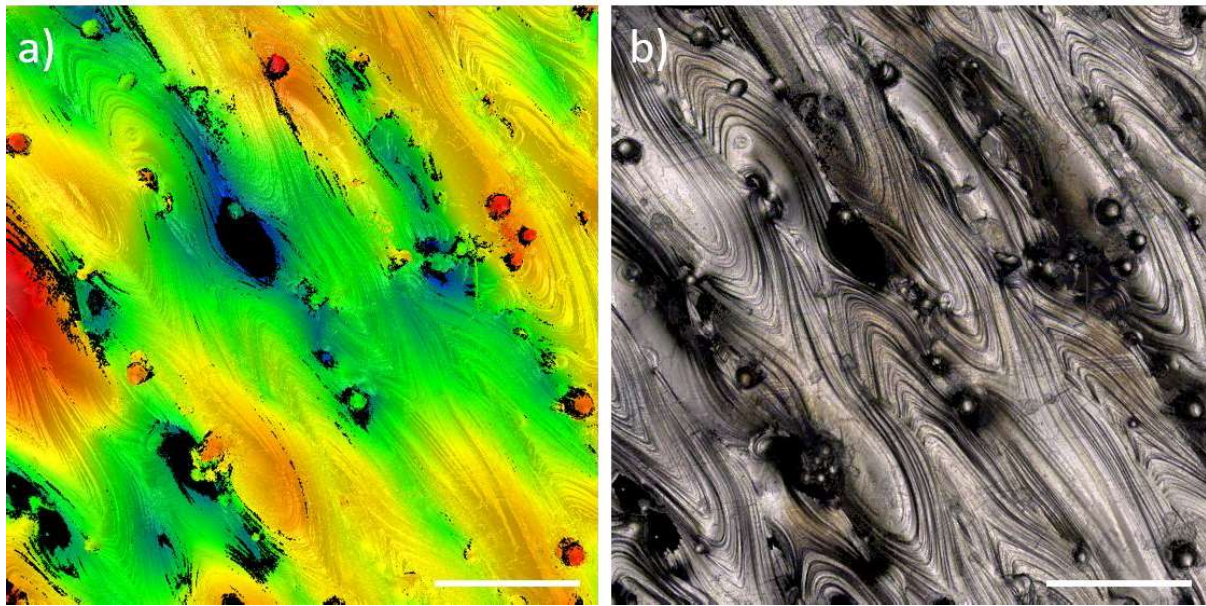


Figure 4. a) Surface topography the Ti-6Al-4V LPBF cube displayed in false color. The color scale gives the surface height, the scale is linear from $-60\ \mu\text{m}$ to $+60\ \mu\text{m}$ b) True color measurement combined with the surface topography. This is a stitched image with a 1 mm by 1 mm field of view with a 50x objective lens (NA 0.55, $0.34\ \text{mm} \times 0.34\ \text{mm}$ single measurement field of view). The lateral scale bar is $250\ \mu\text{m}$.

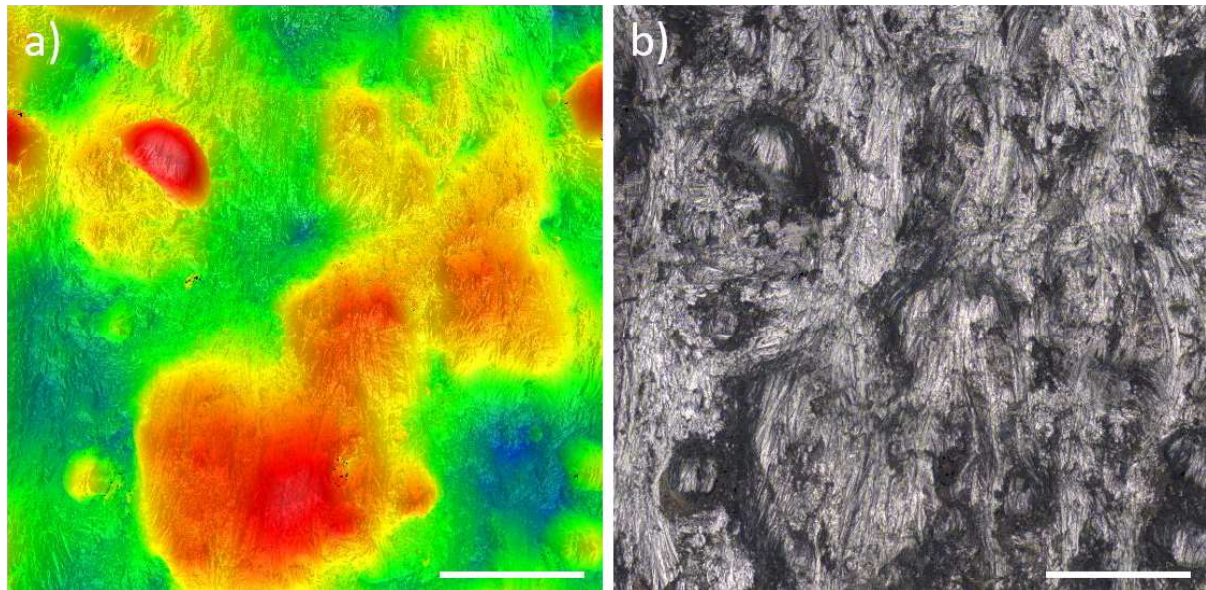


Figure 5. a) Surface topography the Al-Si-10Mg LPBF cube displayed in false color. The color scale gives the surface height, the scale is linear from $-60\ \mu\text{m}$ to $+60\ \mu\text{m}$ b) True color measurement combined with the surface topography. This is a stitched image with a 1 mm by 1 mm field of view with a 50x objective lens (NA 0.55, $0.34\ \text{mm} \times 0.34\ \text{mm}$ single measurement FOV). The lateral scale bar is $250\ \mu\text{m}$.

References

- [1] de Groot P (2015). Principles of interference microscopy for the measurement of surface topography, *Adv. Opt. Photon.*, **7**, 1-65.
- [2] Fay M, Colonna de Lega X, de Groot P (2014). Measuring high-slope and super-smooth Optics with high-dynamic-range coherence scanning interferometry, *Classical Optics*, OW1B.3.
- [3] Beverage, J, Colonna de Lega X, Fay, M (2014). Interferometric microscope with true color imaging. *Proc. SPIE 9203*, *Interferometry XVII: Techniques and Analysis*
- [4] Townsend A, Senin N, Blunt L, Leach R K and Taylor J S 2016 Surface texture metrology for metal additive manufacturing: a review *Precis. Eng.* **46** 34–47
- [5] Gu D 2015 *Laser additive manufacturing of high-performance materials* (Berlin: Springer-Verlag)
- [6] Leach R K 2014 Surface topography measurement instrumentation *Fundamental principles of engineering nanometrology* (Oxford, UK: Elsevier) pp 115–175
- [7] Gomez C, Thompson A, Su R, DiSciaccia J, Lawes S, Leach R K 2017 Optimisation of coherence scanning interferometry for powder bed fusion surfaces *Opt. Engineering* under review
- [8] Senin N, Thompson A, Leach R K 2017 Characterisation of the topography of metal additive surface features with different measurement technologies *Measurement Science and Technology* accepted
- [9] Thompson A, Senin N, Giusca C, Leach R K 2017 Topography of selectively laser melted surfaces: A comparison of different measurement methods *CIRP Annals - Manufacturing Technology*

- [10] ASM Aerospace Specification Metals 2015 *Titanium Ti-6AL-4V-AMS-4911*
- [11] Brandl E, Heckenberger U, Holzinger V, Buchbinder D 2012 Additive manufactured AlSi10Mg samples using selective laser melting (SLM): microstructure, high cycle fatigue, and fracture behavior *Mater. Des.* 34 pp 159–169
- [12] Sidambe A T 2014 Biocompatibility of advanced manufactured titanium implants – a review *Materials* 7 pp 8168–8188