

Technical University of Denmark



Rapid and Simultaneous Determination of Np and Pu in Environmental Samples Using Sequential Injection Anion Exchange Chromatography and ICP-MS

Qiao, Jixin

Publication date:
2009

[Link back to DTU Orbit](#)

Citation (APA):

Qiao, J. (2009). Rapid and Simultaneous Determination of Np and Pu in Environmental Samples Using Sequential Injection Anion Exchange Chromatography and ICP-MS [Sound/Visual production (digital)]. NKS-B RadWorkshop, Risø, Denmark, 16/11/2009

DTU Library
Technical Information Center of Denmark

General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

Rapid and Simultaneous Determination of Np and Pu in Environmental Samples Using Sequential Injection Anion Exchange Chromatography and ICP-MS

PhD student: Jixin Qiao

Supervisor: Xiaolin Hou

Co-supervisor: Per Roos, Manuel Miró

Risø-DTU, Technical University of Denmark

BACKGROUND

Plutonium isotopes ($^{238,239,240,241}\text{Pu}$) and Neptunium (^{237}Np) are highly hazardous radioactive pollutants in the environment due to:

- 1) long radioactive half-lives;
- 2) high radiological toxicities;
- 3) long-term persistence in environment.

Table 1. Nuclear Properties of Important Plutonium Isotopes

Isotope	Half-life	Specific activity (Bq/g)	Principal decay mode	Decay energy (MeV)
^{238}Pu	87.7yr	6.338×10^{11}	α	α 5.499 (70.9%)
^{239}Pu	$2.411 \times 10^4\text{yr}$	2.296×10^9	α	α 5.157 (70.77%)
^{240}Pu	$6.561 \times 10^3\text{yr}$	8.401×10^9	α	α 5.168 (72.8%)
^{241}Pu	14.35yr	3.825×10^{12}	β ->99.99%	α 4.896 (83.2%)
^{237}Np	$2.411 \times 10^6\text{yr}$	2.603×10^7	α	α 4.788 (51%)

BACKGROUND

The determination of plutonium isotopes and Neptunium in the environment is important for:

- 1) environmental risk assessment and monitoring of sites around nuclear facilities;
- 2) emergency preparedness;
- 3) surveys for the contaminated area resulting from nuclear weapon tests, nuclear accidents, and the discharge of nuclear waste.


BACKGROUND

- 1) The **levels of plutonium** isotopes and **neptunium** in the environment are very **low** and depending of the location.
- 2) **Plutonium** and **neptunium** often coexist with **matrix elements** (Ca, Mg, Al, V...) and **other radionuclides** (Th, U, Am, Cm...).

Table 2. Environmental level of ^{238}Pu and $^{239,240}\text{Pu}$

Sample	^{238}Pu	$^{239,240}\text{Pu}$
Soil, Bq/kg	0.07	0.1-7
Herbaceous plants ,Bq/kg	4.5×10^{-4}	0.3-2
Lichen, Bq/kg	-	4-10
Grain, vegetables, Bq/kg	$(0.2-14) \times 10^{-4}$	$(4-89) \times 10^{-4}$
Lake water ,Bq/L	-	$(0.1-29) \times 10^{-6}$
Sea water, Bq/L	-	$(0.7-52) \times 10^{-6}$

BACKGROUND

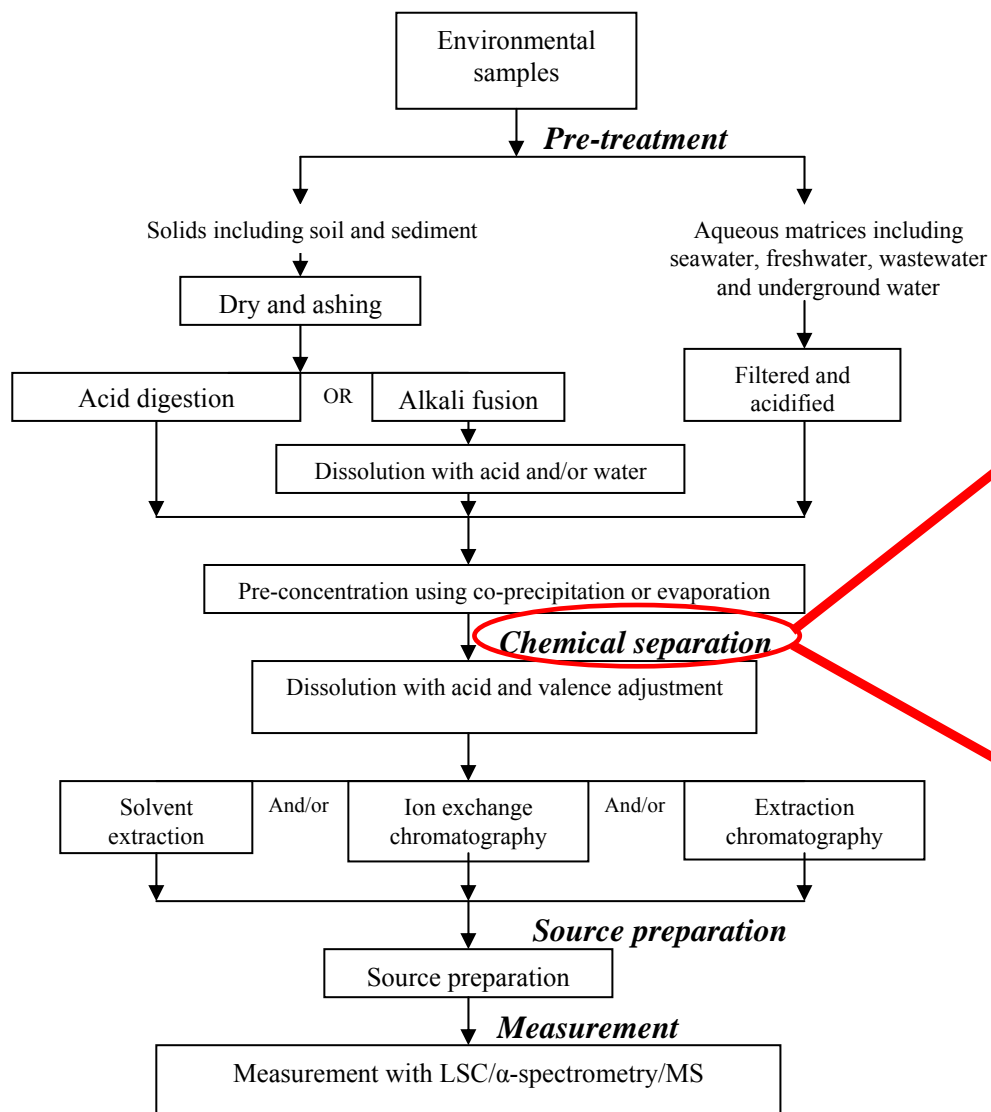

♣ Advantages

- ♣ sensitive;
- ♣ precise;
- ♣ accurate.

♣ Disadvantages

- ♣ time-consuming;
- ♣ labour intensive;
- ♣ generate hazardous liquid and solid waste.

Fig. 1 Analytical procedure for the determination of Pu and Np in environmental samples

OBJECTIVE

Objective:

To develop a **new** analytical method for determination of plutonium isotopes and neptunium in environmental samples.

Main Points: 1)Automatic
2)Rapid
3)Simultaneous

MAIN CHALLENGES

- **Small column size**
- **Same behavior of Pu and Np on the column**
- **Valence adjustment**
- **High chemical yields**
- **Good decontamination factors (U, Th, Pb)**

STRATEGY

Ion-exchange
Chromatography



Sequential
Injection (SI)

ICP-MS

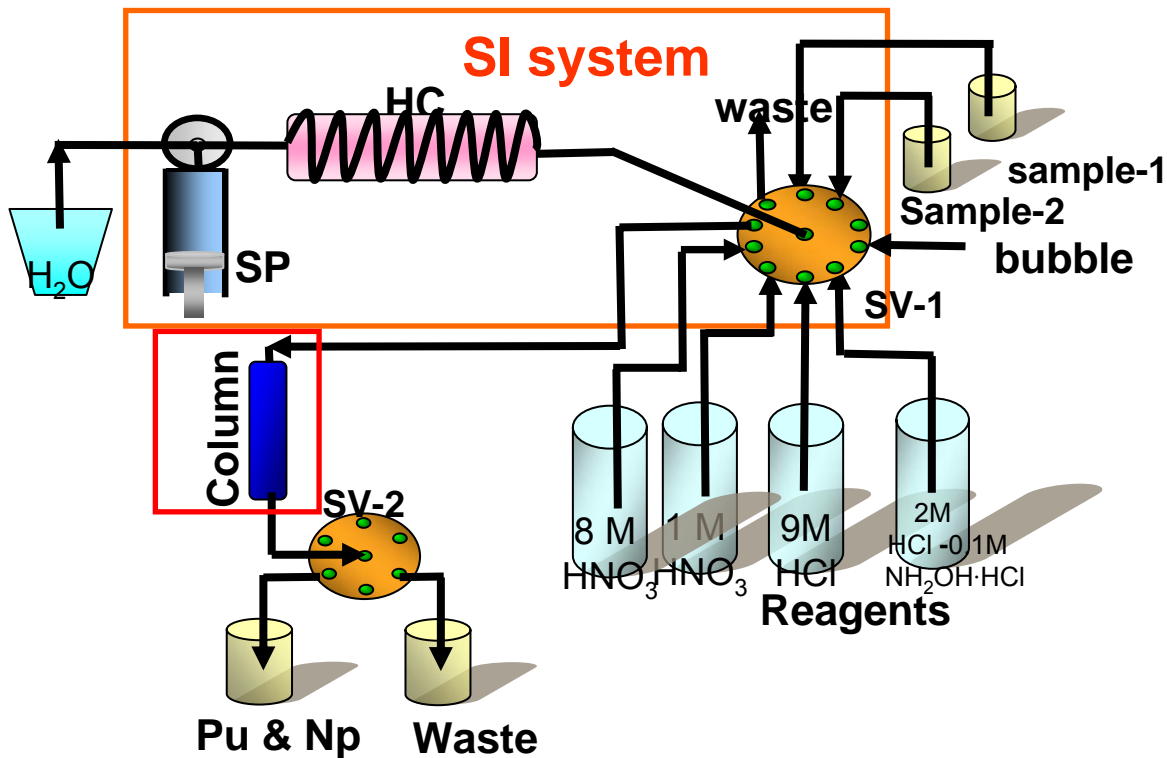


Fig. 2 Scheme of a SI system

Samples

Soil: Danish soil, reference material from a laboratory round-robin intercomparison. The reference values of ^{239}Pu and ^{240}Pu are 0.140 ± 0.008 and 0.098 ± 0.006 Bg/kg.

Sediment, plants, seawater...

Anion exchange chromatographic column

Column size: 16mL (1.0 x 20 cm)

8mL (0.7 x 20 cm)

4mL (0.7 x 10 cm)

2mL (0.5 x 10 cm)

2mL (0.7 x 5.0 cm)

**Resin: AG 1x2 (50-100mesh), AG 1x4 (50-100mesh),
AG 1x4(100-200mesh), AG 1x8(50-100mesh).**

Instrumentation

1) FIA lab system 3500

- ♣ Syringe/peristaltic pump
- ♣ 10-port selection valve



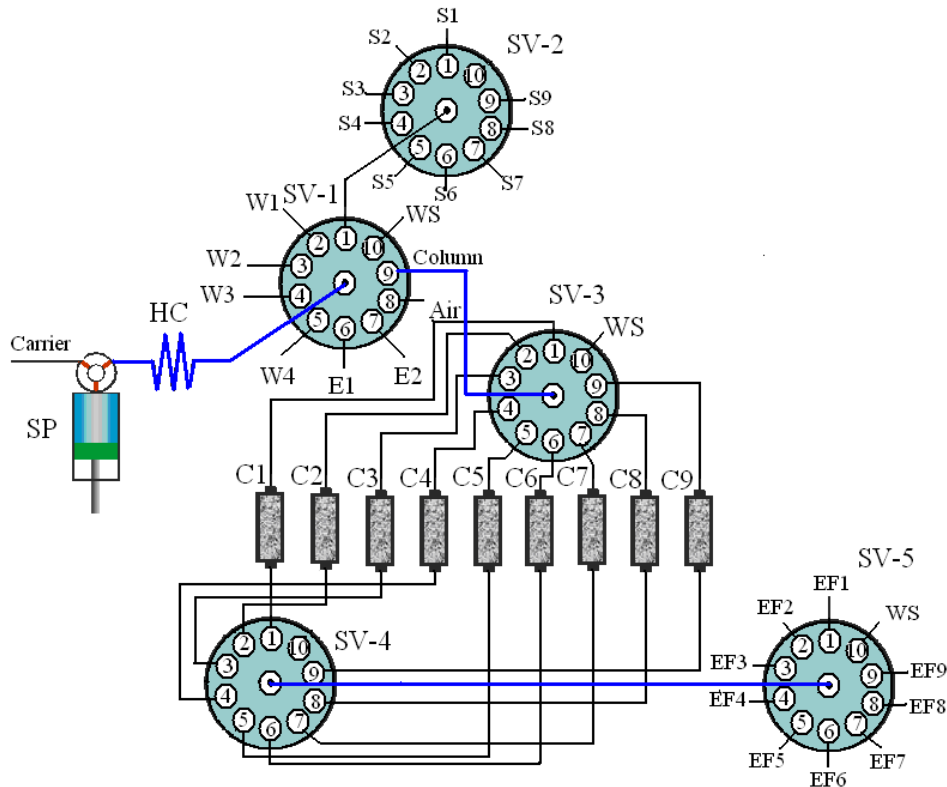


Fig.3. Scheme of the experimental setup

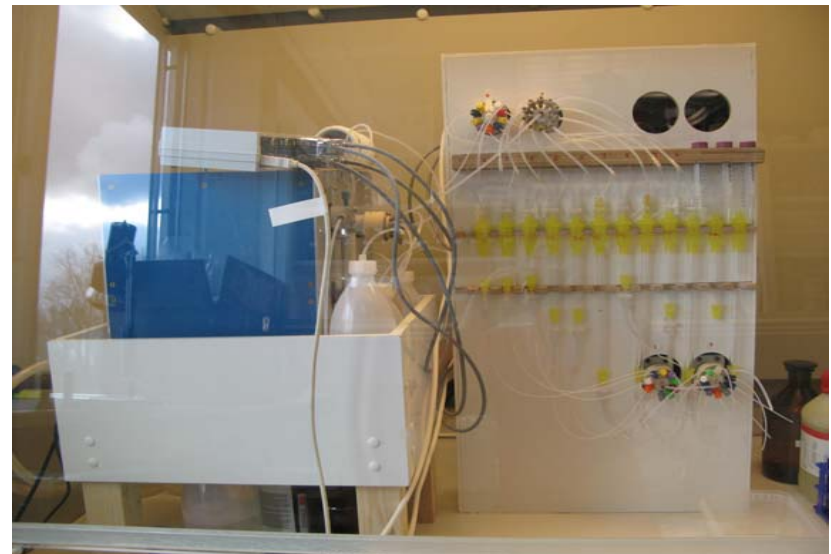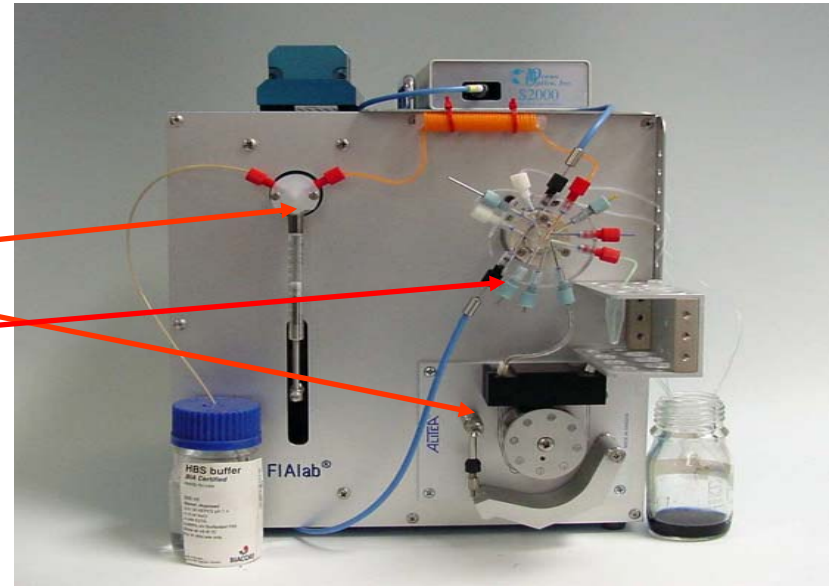


Instrumentation

2) ICP-MS


Thermo X-series inductively coupled plasma mass spectrometry (ICP-MS)



Pre-treatment

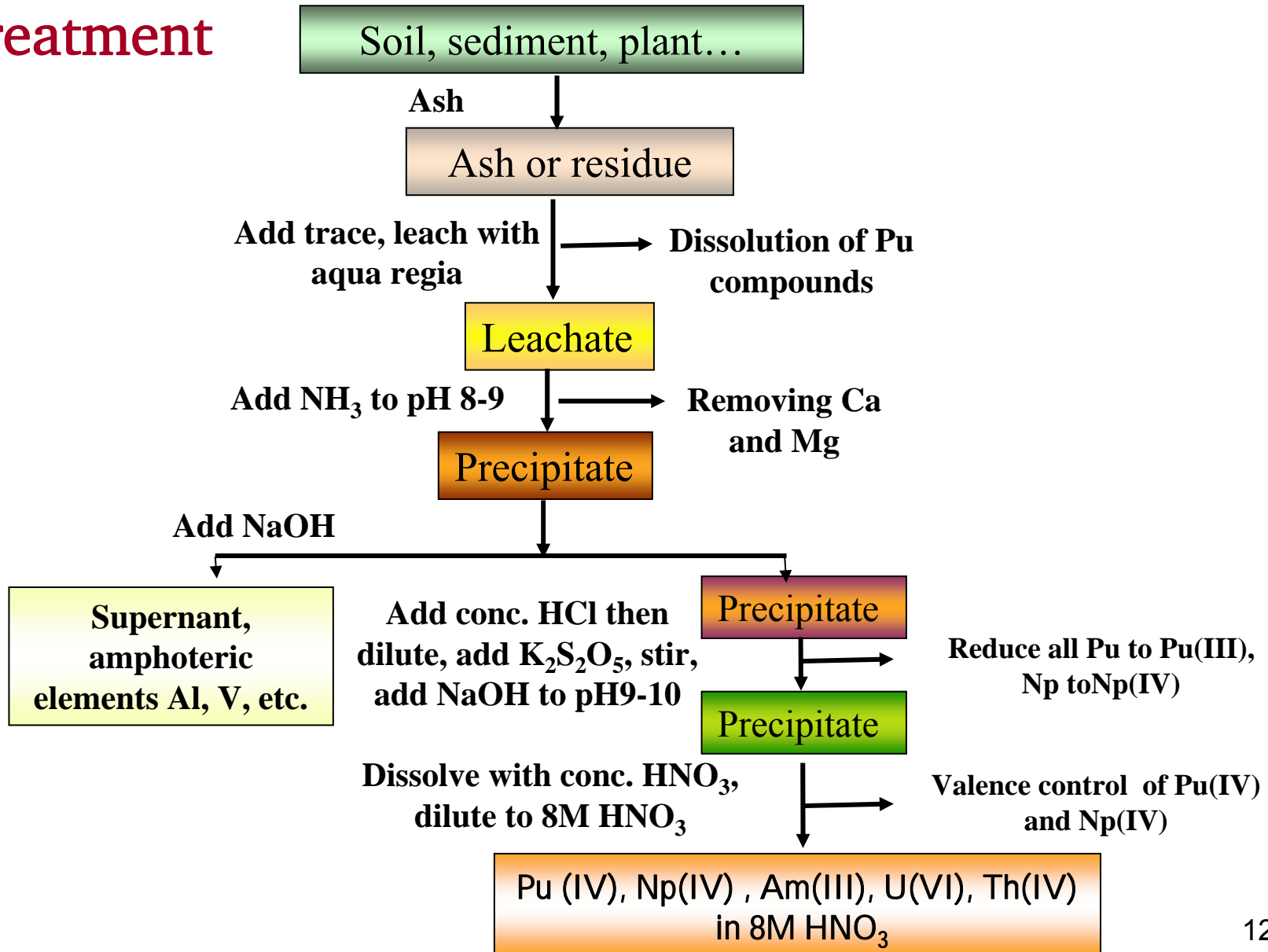

Separation & detection

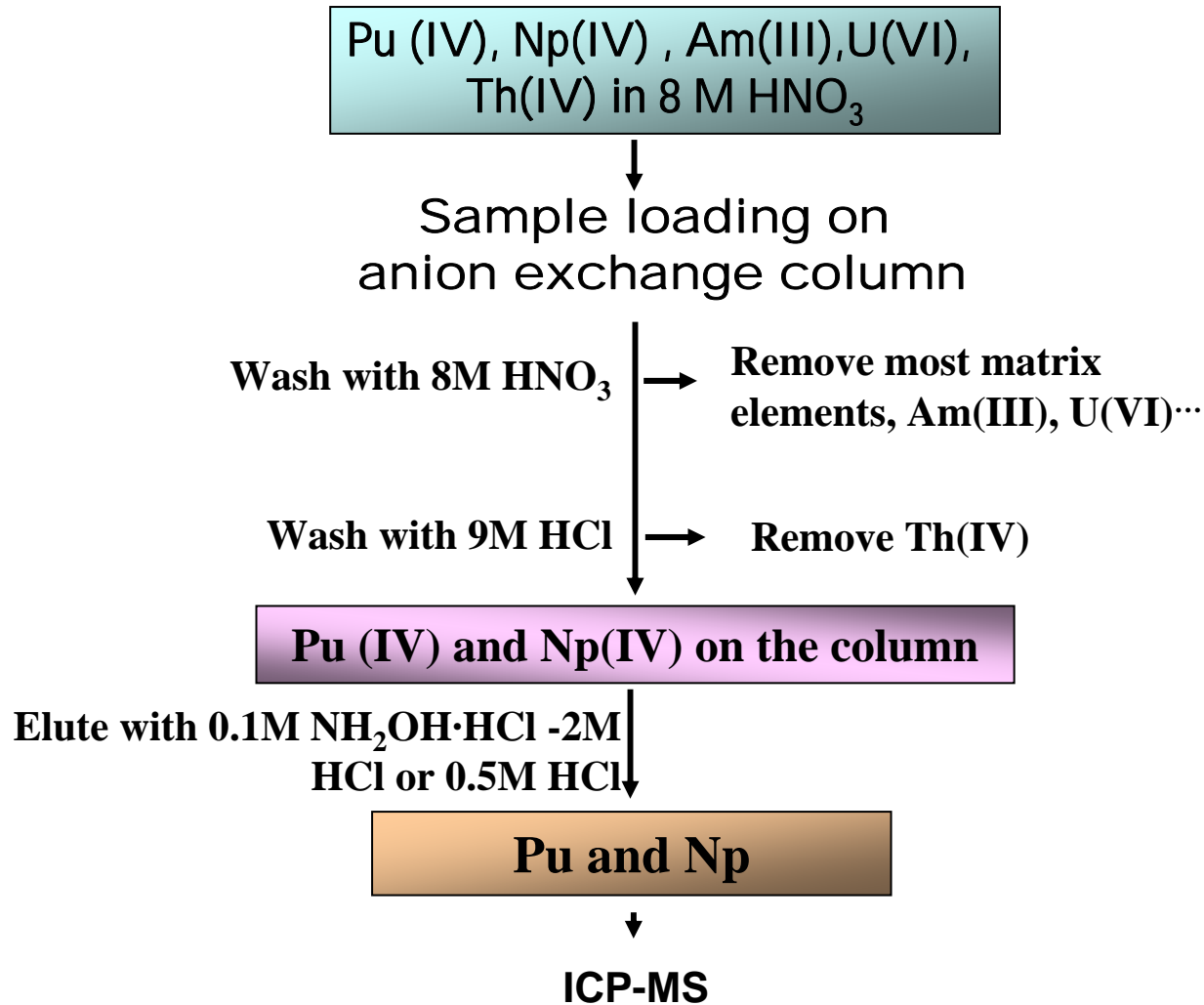

Experimental parameters for comparison

- **Column size**
- **Resin type**
- **Washing solution (1.0-8.0 mol/L HNO₃)**
- **Elution solution (NH₂OHHCl-HCl, 0.1-1.0mol/L HCl)**
- **Flow rate (1.0-5.0 mL/min)**

Key factors for evaluation of experimental results

- Chemical yields of Pu and Np

😊 >85% 😞 <85%

- Ratio between the chemical yield of ^{237}Np and ^{242}Pu

😊 0.9-1.1 😞 otherwise

- Measured values of ^{239}Pu and ^{240}Pu

😊 agree well with the reference values

😞 otherwise

- Decontamination factors for U, Th and Pb

😊 > 10^3 😞 < 10^3

Table 3. comparison of different experimental conditions for the separation of Pu and Np (1)

Column size	Resin	Separation condition#	Chemical yield of ^{242}Pu , Y_{Pu} (%)	Chemical yield of ^{237}Np , Y_{Np} (%)	Ration of $Y_{\text{Np}}/Y_{\text{Pu}}$	^{239}Pu measured (Bq/kg) *	^{240}Pu measured (Bq/kg) *	Decontamination factor **		
								^{238}U	^{232}Th	Pb
16mL (1.0 × 20cm)	AG1 × 2	W-1, 2.5 E-1, 1.0	102.3 ± 5.1	95.8 ± 4.8	0.9	0.14 ± 0.0	0.10 ± 0.01	3.3 × 10 ³	1.9 × 10 ⁴	1.4 × 10 ⁴
	AG1 × 4		99.9 ± 5.0	94.8 ± 4.7	0.9	0.23 ± 0.02	0.08 ± 0.01	3.0 × 10 ²	2.9 × 10 ³	5.6 × 10 ³
	AG1 × 8		96.4 ± 4.8	90.9 ± 4.5	0.9	1.39 ± 0.14	0.11 ± 0.01	4.2 × 10 ¹	3.8 × 10 ²	1.2 × 10 ⁴
8mL (0.7 × 20cm)	AG1 × 2	W-1, 2.5 E- 2, 1.0	71.5 ± 3.6	67.4 ± 3.4	0.9	0.17 ± 0.02	0.12 ± 0.01	1.1 × 10 ³	5.2 × 10 ³	3.1 × 10 ⁴
	AG1 × 4		100.0 ± 5.4	100.0 ± 5.3	1.0	0.16 ± 0.02	0.10 ± 0.01	1.6 × 10 ³	6.4 × 10 ³	3.9 × 10 ³
	AG1 × 8		94.2 ± 4.7	87.9 ± 4.4	0.9	0.16 ± 0.02	0.10 ± 0.01	3.2 × 10 ²	7.8 × 10 ²	8.1 × 10 ³
	AG1 × 8	W-1, 5.0 E-2, 2.5	91.9 ± 4.6	80.1 ± 4.0	0.9	0.18 ± 0.02	0.12 ± 0.01	1.1 × 10 ²	9.1 × 10 ¹	8.8 × 10 ³
4mL (0.7 × 10cm)	AG1 × 2	W-1, 2.5 E-2, 1.0	71.2 ± 3.6	48.4 ± 2.4	0.7	0.12 ± 0.01	0.06 ± 0.01	2.1 × 10 ³	4.4 × 10 ³	1.1 × 10 ⁴
	AG1 × 4		100.0 ± 5.0	98.2 ± 4.9	1.0	0.12 ± 0.01	0.10 ± 0.01	1.3 × 10 ³	2.4 × 10 ³	2.2 × 10 ⁴
	AG1 × 8		98.7 ± 4.9	97.2 ± 4.9	1.0	0.16 ± 0.02	0.10 ± 0.01	1.0 × 10 ³	8.9 × 10 ²	8.6 × 10 ³
	AG1 × 8	W-1, 5.0 E-2, 2.5	92.6 ± 4.6	86.3 ± 4.3	0.9	0.17 ± 0.02	0.12 ± 0.01	2.2 × 10 ²	9.8 × 10 ¹	2.7 × 10 ³

Table 3. comparison of different experimental conditions for the separation of Pu and Np (2)

Column size	Resin	Separation condition#	Chemical yield of ²⁴² Pu, Y _{Pu} (%)	Chemical yield of ²³⁷ Np, Y _{Np} (%)	Ration of Y _{Np} /Y _{Pu}	²³⁹ Pu measured (Bq/kg) *	²⁴⁰ Pu measured (Bq/kg) *	Decontamination factor **		
								²³⁸ U	²³² Th	Pb
2mL (0.5 × 10cm)	AG1 × 2	W-2, 2.5 E-2, 2.5	75.0 ± 3.8	19.0 ± 1.0	0.3	0.26 ± 0.03	0.08 ± 0.01	3.2 × 10 ²	2.0 × 10 ²	2.6 × 10 ³
		W-2, 5.0 E-2, 2.5	48.6 ± 2.4	35.7 ± 1.8	0.7	0.15 ± 0.02	0.10 ± 0.01	8.9 × 10 ²	6.8 × 10 ²	9.3 × 10 ³
	AG1 × 4	W-2, 2.5 E-2, 2.5	103.0 ± 5.2	106.0 ± 5.3	1.0	0.14 ± 0.01	0.09 ± 0.01	3.9 × 10 ³	2.4 × 10 ⁴	2.7 × 10 ⁴
		W-2, 5.0 E-2, 2.5	94.0 ± 4.7	89.7 ± 4.5	1.0	0.25 ± 0.03	0.09 ± 0.01	3.9 × 10 ²	6.7 × 10 ³	1.6 × 10 ⁴
	AG1 × 8	W-2, 2.5 E-2, 2.5	90.5 ± 5.0	88.7 ± 4.9	1.0	0.25 ± 0.03	0.09 ± 0.01	3.4 × 10 ²	2.7 × 10 ²	9.0 × 10 ⁵
		W-2, 5.0 E-2, 2.5	100.5 ± 5.0	98.7 ± 4.9	1.0	0.29 ± 0.03	0.07 ± 0.01	4.8 × 10 ¹	2.1 × 10 ²	2.9 × 10 ⁵
	AG1 × 4	W-2, 2.5 E-4, 2.5	72.9 ± 3.6	64.1 ± 3.2	0.9	0.18 ± 0.02	0.12 ± 0.01	1.6 × 10 ³	1.7 × 10 ⁴	6.1 × 10 ³
		W-3, 2.5 E-4, 2.5	81.8 ± 4.1	69.2 ± 4.1	0.8	0.38 ± 0.04	0.10 ± 0.01	2.1 × 10 ³	1.3 × 10 ⁴	7.5 × 10 ³
		W-4, 2.5 E-4, 2.5	80.2 ± 4.0	63.8 ± 4.1	0.8	0.19 ± 0.02	0.09 ± 0.01	2.4 × 10 ³	1.4 × 10 ⁴	1.5 × 10 ⁴
		W-5, 2.5 E-4, 2.5	39.6 ± 2.0	20.6 ± 4.1	0.5	0.20 ± 0.02	0.15 ± 0.02	2.6 × 10 ³	1.1 × 10 ⁴	9.4 × 10 ³
		W-6, 2.5 E-4, 2.5	31.0 ± 1.6	11.7 ± 4.1	0.4	0.19 ± 0.02	0.20 ± 0.02	3.1 × 10 ³	5.0 × 10 ⁴	1.2 × 10 ⁵

Table 3. comparison of different experimental conditions for the separation of Pu and Np (3)

Column size	Resin	Separation condition #	Chemical yield of ²⁴² Pu, Y _{Pu} (%)	Chemical yield of ²³⁷ Np, Y _{Np} (%)	Ration of Y _{Np} /Y _{Pu}	²³⁹ Pu measured (Bq/kg) *	²⁴⁰ Pu measured (Bq/kg) *	Decontamination factor **		
								²³⁸ U	²³² Th	Pb
2mL (0.5 × 10cm)	AG1 × 4, 100-200 mesh	W-2, 2.5 § E-3, 2.5	61.0 ± 3.1	64.9 ± 3.2	1.1	0.18 ± 0.02	0.12 ± 0.01	3.1 × 10 ³	2.8 × 10 ⁴	1.2 × 10 ⁴
		W-2, 2.5 E-4, 2.5	91.6 ± 4.6	91.0 ± 4.6	1.0	0.14 ± 0.01	0.10 ± 0.01	6.9 × 10 ³	1.7 × 10 ⁴	1.0 × 10 ³
		W-2, 2.5 E- 5, 2.5	66.8 ± 3.3	74.6 ± 3.7	1.1	0.14 ± 0.01	0.09 ± 0.01	8.6 × 10 ³	1.2 × 10 ⁴	1.0 × 10 ³
		W-2, 2.5 E-6, 2.5	78.5 ± 3.9	81.6 ± 4.1	1.0	0.14 ± 0.01	0.07 ± 0.01	6.3 × 10 ³	1.9 × 10 ⁴	1.5 × 10 ³
		W-3, 2.5 E-4, 2.5	35.9 ± 1.8	61.2 ± 3.1	1.7	0.22 ± 0.02	0.19 ± 0.02	3.4 × 10 ³	3.1 × 10 ⁴	2.9 × 10 ³
		W-4, 2.5 E-4, 2.5	82.3 ± 4.1	80.9 ± 4.0	1.0	0.21 ± 0.02	0.12 ± 0.01	3.8 × 10 ³	1.4 × 10 ⁴	4.3 × 10 ⁴
		W-5, 2.5 E-4, 2.5	63.6 ± 3.2	30.1 ± 1.5	0.5	0.19 ± 0.02	0.10 ± 0.01	4.5 × 10 ³	3.1 × 10 ⁴	1.8 × 10 ⁴
2mL (0.7 × 5cm)	AG1 × 4, 100-200 mesh	W-2, 2.5 E-4, 2.5	40.4 ± 2.0	37.1 ± 1.9	0.9	0.18 ± 0.02	0.07 ± 0.01	1.1 × 10 ⁴	1.6 × 10 ³	7.2 × 10 ³
1mL (0.5 × 5cm)	AG1 × 4, 100-200 mesh	W-2, 2.5 E-4, 2.5	50.7 ± 2.5	44.8 ± 2.2	0.9	0.14 ± 0.02	0.10 ± 0.01	4.1 × 10 ³	4.0 × 10 ³	1.5 × 10 ⁴

The reference values of ²³⁹Pu and ²⁴⁰Pu concentration in the Danish soil were reported to be 0.140 ± 0.008 Bq/kg and 0.098 ± 0.006 Bq/kg.* Experimental results are given as the average of three replicates ± standard deviation. ** The relative standard deviations were in all instances better than 10%. § flow rate, mL/min.# W-1: washing sequence 200mL of 8 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-2: 100mL of 8 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-3: 100mL of 6 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-4: 100mL of 4 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-5: 100mL of 2 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-6: 100mL of 1 mol/L HNO₃ + 100mL of 9 mol/L HCl; Pu eluting solution: E-1:Pu elution solution 200mL of 0.1 mol/L NH₂OH-HCl-2 mol/L HCl; E-2: 100mL of 0.1 mol/L NH₂OH-HCl-2 mol/L HCl; E-3: 40mL of 0.1 mol/L NH₂OH-HCl-2 mol/L HCl; E-4: 40mL of 0.5 mol/L HCl; E-5: 40mL of 0.1 mol/L HCl; E-6: 40mL of 1.0 mol/L HCl.

Table 3. Comparison of different experimental conditions for the separation of Pu and Np (1)

Column size	Resin	Separation condition#	Chemical yield of ²⁴² Pu, Y _{Pu} (%)	Chemical yield of ²³⁷ Np, Y _{Np} (%)	Ration of Y _{Np} /Y _{Pu}	²³⁹ Pu measured (Bq/kg) *	²⁴⁰ Pu measured (Bq/kg) *	Decontamination factor **		
								²³⁸ U	²³² Th	Pb
16mL (1.0 × 20cm)	AG1 × 2	W-1, 2.5 E-1, 1.0	😊	😊	😊	😊	😊	😊	😊	😊
	AG1 × 4		😊	😊	😊	😞	😞	😞	😊	😊
	AG1 × 8		😊	😊	😊	😞	😊	😞	😞	😊
8mL (0.7 × 20cm)	AG1 × 2	W-1, 2.5 E- 2, 1.0	😞	😞	😊	😊	😊	😊	😊	😊
	AG1 × 4		😊	😊	😊	😊	😊	😊	😊	😊
	AG1 × 8		😊	😊	😊	😊	😊	😞	😞	😊
	AG1 × 8	W-1, 5.0 E-2, 2.5	😊	😞	😊	😞	😞	😞	😞	😊
4mL (0.7 × 10cm)	AG1 × 2	W-1, 2.5 E-2, 1.0	😞	😞	😞	😊	😞	😊	😊	😊
	AG1 × 4		😊	😊	😊	😊	😊	😊	😊	😊
	AG1 × 8		😊	😊	😊	😞	😊	😊	😞	😊
	AG1 × 8	W-1, 5.0 E-2, 2.5	😊	😊	😊	😞	😞	😞	😞	😊

Table 3. Comparison of different experimental conditions for the separation of Pu and Np (2)

Column size	Resin	Separation condition#	Chemical yield of ²⁴² Pu, Y _{Pu} (%)	Chemical yield of ²³⁷ Np, Y _{Np} (%)	Ration of Y _{Np} /Y _{Pu}	²³⁹ Pu measured (Bq/kg) *	²⁴⁰ Pu measured (Bq/kg) *	Decontamination factor **		
								²³⁸ U	²³² Th	Pb
2mL (0.5 × 10cm)	AG1 × 2	W-2, 2.5 E-3, 2.5	☹	☹	☹	☹	☹	☹	☹	☺
		W-2, 5.0 E-3, 2.5	☹	☹	☹	☺	☺	☹	☹	☺
	AG1 × 4	W-2, 2.5 E-3, 2.5	☺	☺	☺	☺	☺	☺	☺	☺
		W-2, 5.0 E-3, 2.5	☺	☺	☺	☹	☺	☹	☺	☺
	AG1 × 8	W-2, 2.5 E-3, 2.5	☺	☺	☺	☹	☺	☹	☹	☺
		W-2, 5.0 E-3, 2.5	☺	☺	☺	☹	☹	☹	☹	☺
	AG1 × 4	W-2, 2.5 E-4, 2.5	☹	☹	☺	☹	☹	☺	☺	☺
		W-3, 2.5 E-4, 2.5	☹	☹	☹	☹	☺	☺	☺	☺
		W-4, 2.5 E-4, 2.5	☹	☹	☹	☹	☺	☺	☺	☺
		W-5, 2.5 E-4, 2.5	☹	☹	☹	☹	☹	☺	☺	☺
W-6, 2.5 E-4, 2.5		☹	☹	☹	☹	☹	☺	☺	☺	

Table 3. Comparison of different experimental conditions for the separation of Pu and Np (3)

Column size	Resin	Separation condition #	Chemical yield of ²⁴² Pu, Y _{Pu} (%)	Chemical yield of ²³⁷ Np, Y _{Np} (%)	Ration of Y _{Np} /Y _{Pu}	²³⁹ Pu measured (Bq/kg) *	²⁴⁰ Pu measured (Bq/kg) *	Decontamination factor **		
								²³⁸ U	²³² Th	Pb
2mL (0.5 × 10cm)	AG1 × 4, 100-200 mesh	W-2, 2.5 § E-3, 2.5	☹	☹	☺	☹	☹	☺	☺	☺
		W-2, 2.5 E-4, 2.5	☺	☺	☺	☺	☺	☺	☺	☺
		W-2, 2.5 E- 5, 2.5	☹	☹	☺	☺	☺	☺	☺	☺
		W-2, 2.5 E-6, 2.5	☹	☹	☺	☺	☹	☺	☺	☺
		W-3, 2.5 E-4, 2.5	☹	☹	☹	☹	☹	☺	☺	☺
		W-4, 2.5 E-4, 2.5	☹	☹	☺	☹	☹	☺	☺	☺
		W-5, 2.5 E-4, 2.5	☹	☹	☹	☹	☺	☺	☺	☺
		W-2, 2.5 E-4, 2.5	☹	☹	☺	☹	☹	☺	☺	☺
2mL (0.7 × 5cm)	AG1 × 4, 100-200 mesh	W-2, 2.5 E-4, 2.5	☹	☹	☺	☹	☹	☺	☺	☺
1mL (0.5 × 5cm)	AG1 × 4, 100-200 mesh	W-2, 2.5 E-4, 2.5	☹	☹	☺	☺	☺	☺	☺	☺

The reference values of ²³⁹Pu and ²⁴⁰Pu concentration in the Danish soil were reported to be 0.140 ± 0.008 Bq/kg and 0.098 ± 0.006 Bq/kg.* Experimental results are given as the average of three replicates \pm standard deviation. ** The relative standard deviations were in all instances better than 10%. § flow rate, mL/min.# W-1: washing sequence 200mL of 8 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-2: 100mL of 8 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-3: 100mL of 6 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-4: 100mL of 4 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-5: 100mL of 2 mol/L HNO₃ + 100mL of 9 mol/L HCl; W-6: 100mL of 1 mol/L HNO₃ + 100mL of 9 mol/L HCl; Pu eluting solution: E-1:Pu elution solution 200mL of 0.1 mol/L NH₂OH·HCl-2 mol/L HCl; E-2: 100mL of 0.1 mol/L NH₂OH·HCl-2 mol/L HCl; E-3: 40mL of 0.1 mol/L NH₂OH·HCl-2 mol/L HCl; E-4: 40mL of 0.5 mol/L HCl; E-5: 40mL of 0.1 mol/L HCl; E-6: 40mL of 1.0 mol/L HCl.

Main Results

- 1) ^{242}Pu performs well as a tracer for both Pu isotope and ^{237}Np .
- 2) Cross-link of the resins has significant influence on the separation efficiency. **Finally, AG 1x4 resin was chosen as the optimum.**
- 3) **Small-sized column packed with 2mL resin suffices up to 10g of soil.**

Table 4. Selected results from the experiment (10g of soil)

Column size	Resin	Chemical yield of ^{242}Pu , Y_{Pu} (%)	Chemical yield of ^{237}Np , Y_{Np} (%)	Ration of $Y_{\text{Np}}/Y_{\text{Pu}}$	^{239}Pu measured (Bq/kg) *	^{240}Pu measured (Bq/kg)**	Decontamination factor ***		
							^{238}U	^{232}Th	^{208}Pb
2mL (0.5 × 10cm)	AG1 × 4, 50-100 mesh	103.0 ± 5.2	106.0 ± 5.3	1.0	0.14 ± 0.01	0.09 ± 0.01	3.9 × 10 ³	2.4 × 10 ⁴	2.7 × 10 ⁴
	AG1 × 4, 100-200 mesh	91.6 ± 4.6	91.0 ± 4.6	1.0	0.14 ± 0.01	0.10 ± 0.01	6.9 × 10 ³	1.7 × 10 ⁴	1.0 × 10 ³

*The reference value is 0.140 ± 0.008 Bq/kg. **The reference value is 0.098 ± 0.006 Bq/kg.

** The relative standard deviations were in all instances better than 10%.

Main Results

- 4) The total time of on-line separation for a single sample is ~ **2.5h**.
For comparison: 2-3days is need for off-line separation.
- 5) Chemical yields of Pu and Np equally range from 90% to 100%.
- 6) Decontamination factor for ^{238}U , ^{232}Th and ^{208}Pb are in the range of 10^3 to 10^4 .

Table 4. Selected results from the experiment (10g of soil)

Column size	Resin	Chemical yield of ^{242}Pu , Y_{Pu} (%)	Chemical yield of ^{237}Np , Y_{Np} (%)	Ration of $Y_{\text{Np}}/Y_{\text{Pu}}$	^{239}Pu measured (Bq/kg) *	^{240}Pu measured (Bg/kg)**	Decontamination factor ***		
							^{238}U	^{232}Th	^{208}Pb
2mL (0.5 × 10cm)	AG1 × 4, 50-100 mesh	103.0 ± 5.2	106.0 ± 5.3	1.0	0.14 ± 0.01	0.09 ± 0.01	3.9 × 10 ³	2.4 × 10 ⁴	2.7 × 10 ⁴
	AG1 × 4, 100-200 mesh	91.6 ± 4.6	91.0 ± 4.6	1.0	0.14 ± 0.01	0.10 ± 0.01	6.9 × 10 ³	1.7 × 10 ⁴	1.0 × 10 ³

*The reference value is 0.140 ± 0.008 Bg/kg. **The reference value is 0.098 ± 0.006 Bg/kg.

*** The relative standard deviations were in all instances better than 10%.

- 1) Innovation:** Automatic
Rapid
Simultaneous
Low consumption of resins
Low generation of wastes
- 2) Next step:** Stability of Np(IV) and Pu(IV)
Capacity of the SI system
Reusability of the resin

- **Xiaolin Hou**
- **Per Roos**
- **Manuel Miró**
- **Radioecology and Tracers Programme (headed by Sven P. Nielsen), Radiation Research Division, Risø-DTU, Denmark.**

THANK YOU !

Personal e-mail: qiaojixin2004@gmail.com