A promising new ELISA diagnostic test for cattle babesiosis based on Babesia bigemina Apical Membrane Antigen-1

Alessandra Torina¹, Antonio Cordaro¹, Valeria Blanda^{1*}, Rosalia D'Agostino¹, Salvatore Scimeca¹, Maria E. Scariano¹, Guido Sireci² & Rossella C. Lelli¹

¹ Istituto Zooprofilattico Sperimentale della Sicilia, Via G. Marinuzzi 3, 90129 Palermo, Italy. ² Dipartimento di Biopatologia e Biotecnologie Mediche e Forensi, Università di Palermo, Corso Tukory 211, Palermo, Italy.

* Corresponding author at: Istituto Zooprofilattico Sperimentale della Sicilia, Via G. Marinuzzi 3, 90129 Palermo, Italy. Tel.: +39 091 6565364, e-mail: valeria.blanda@gmail.com.

> Veterinaria Italiana 2016, **52** (1), 63-69. doi: 10.12834/Vetlt.74.237.2 Accepted: 21.08.2015 | Available on line: 23.02.2016

LXVII Meeting of the Italian Society for Veterinary Sciences (SISVet) - Selected papers

Keywords

Antibodies, Apical Membrane Antigen-1, Babesia bigemina, Babesiosis, Diagnostic test, ELISA, Piroplasmosis, Tick borne pathogens.

Summary

Babesiosis due to *Babesia bigemina* is a relevant tick-borne disease, affecting cattle worldwide. Many surface proteins of the pathogen including the Apical Membrane Antigen 1 (AMA-1) have been analysed for vaccine and diagnostic purposes. This study focused on *B. bigemina* AMA-1 and on its use for the assessment of diagnostic tests. After bioinformatic analyses, AMA-1 codifying region was amplified and cloned into an expression vector used to induce protein synthesis in *Escherichia coli* cells. AMA-1 was purified by affinity chromatography and used to set up the best condition for an ELISA protocol. Bovine field sera positive to *B. bigemina* were used to evaluate the presence of anti-AMA-1 antibodies. In order to verify the assay specificity, sera positive to *Babesia bovis* or to the piroplasm *Theileria annulata* were also included. Significant differences were obtained between sera negative to both *B. bigemina* and *B. bovis* and samples positive to *B. bigemina*, to *B. bovis* or to both pathogens. No significant reaction was observed with *T. annulata* positive sera. The results showed that AMA-1 protein is suitable to be used as antigen in diagnostic assays for babesiosis diagnosis in cattle, as it does not show any cross reaction with anti-*T. annulata* antibodies.

Un nuovo test ELISA per la diagnosi della babesiosi bovina basato sull'Apical Membrane Antigen-1 di Babesia bigemina

Parole chiave

Anticorpi, Apical Membrane Antigen-1, *Babesia bigemina*, Babesiosi, ELISA, Patogeni trasmessi da zecche, Piroplasmosi, Test diagnostici.

Riassunto

La Babesiosi causata da Babesia bigemina è una malattia trasmessa da zecche che colpisce gli allevamenti bovini in diverse aree del mondo. Molte proteine di superficie del patogeno - tra cui anche l'Apical Membrane Antigen 1 (AMA-1) - sono state studiate al fine di realizzare nuovi vaccini e ottenere nuovi metodi diagnostici. Questo studio è stato condotto sulla proteina AMA-1 di B. bigemina e sul suo possibile utilizzo nella messa a punto di nuovi test diagnostici. In seguito a studi bioinformatici, la regione codificante la proteina AMA-1 è stata amplificata e clonata all'interno di un vettore di espressione. Il vettore ricombinante è stato utilizzato per indurre la sintesi proteica in cellule di Escherichia coli. La proteina AMA-1 è stata purificata mediante cromatografia d'affinità e utilizzata per ottimizzare le condizioni per il protocollo ELISA. Sieri di bovini di campo, risultati positivi a B. bigemina, sono stati utilizzati per rilevare la presenza di anticorpi anti-AMA-1. Al fine di verificare la specificità del saggio, sono stati inclusi nell'analisi sieri positivi a Babesia bovis o al piroplasma Theileria annulata. I sieri positivi a B. bigemina, B. bovis o ad entrambi i patogeni hanno mostrato risposte significativamente differenti rispetto ai sieri negativi sia a B. bigemina che a B. bovis. Non è stata osservata alcuna reazione significativa in presenza di sieri positivi a T. annulata. Lo studio ha mostrato che la proteina AMA-1 di B. bigemina può essere impiegata come antigene nei test per la diagnosi della babesiosi nei bovini, in quanto non mostra alcuna cross-reattività nei confronti degli anticorpi anti-T. annulata.

Introduction

Cattle babesiosis is a tick-borne disease caused by apicomplexan protozoan such as *Babesia bovis* and *Babesia bigemina*. The pathology affects cattle in many areas of the world, strongly reducing meat and milk production (Bock *et al.* 2004, Torina and Caracappa 2007). Understanding the basic molecular mechanism of merozoite invasion of erythrocytes may improve the development of an effective control strategy. Host-cell invasion by apicomplexan parasites is a complex, multistep process (Black and Boothroyd 2000), as a number of proteins are implicated in it, although in most cases their precise functions remain unknown (Zhou *et al.* 2006).

One of the molecules potentially involved in the erythrocyte invasion by *B. bigemina* is the Apical Membrane Antigen-1 (AMA-1), an apically located protein that is shared by many Apicomplexa and has a role in the host red blood cells invasion process. While the biological function of AMA-1 is not completely understood, its cellular localization, a phase-specific expression, and the secretion during the host cell invasion suggest that this protein could have a key role in the red blood cell penetration by the parasite.

The most characterized AMA-1 protein is the one codified by Plasmodium falciparum (PfAMA) (Peterson et al. 1989). This protein has been extensively studied as it is one of the most promising erythrocyte stage candidates among the few anti-malarial vaccine candidates under development (Latitha et al. 2008). Babesia bovis AMA-1 protein (BbAMA-1) was identified and characterized in 2004 (Gaffar et al. 2004) and it showed structural features similar to those of other AMA-1 proteins (Yokoyama et al. 2006). It is supposed that also Babesia bigemina invades erythrocytes using an Apical Membrane Antigen-1. Babesia bigemina AMA-1 is a 65.9KDa protein showing several common features with other AMA-1 family members. It also exhibits a high degree of similarity with the orthologue gene of B. bovis. Structural features of B. bigemina AMA-1 protein were recently described (Torina et al. 2010).

This study was devoted to the synthesis of *B. bigemina* AMA-1 recombinant protein into a prokaryotic expression system and its use to asses a new ELISA test useful for diagnosis of cattle babesiosis.

Material and methods

B. bigemina AMA-1 production in E. coli

Bioinformatic analyses were performed using ClustalW2.0.10 (Larkin *et al.* 2007) software for sequence alignment with the purpose to select a *B. bigemina* surface antigen conserved enough among different strains of the pathogen. DAMBE (Xia and Xie 2001) and MEGA (Kumar *et al.* 2008, Tamura *et al.* 2007) software were used to obtain similarity percentage among the analysed sequences.

Apical Membrane Antigen-1 codifying region was amplified from DNA extracted from blood of a bovine naturally infected by B. bigemina and the amplification product was cloned into the expression vector pET160/GW/D-TOPO (Life Technologies Corporation[™], Carlsbad, CA, United States), following the manufacturer's instruction. The construct was initially used to transform One Shot® TOP10 Chemically Competent Escherichia coli cells (Life Technologies Corporation™, Carlsbad, CA, United States). The plasmid was thus extracted, purified, and sent to Macrogen Europe (Amsterdam, The Netherlands) for sequencing using T7F and T7R universal primers in order to test if the PCR product was cloned in the correct direction. The recombinant plasmid was then used to transform E. coli BL21 Star[™] (DE3) cells (Life Technologies Corporation[™], Carlsbad, CA, United States). Bacterial cells were cultured in LB-Broth at 37°C. Bacterial cell growth was monitored by reading the Optical Density (O.D.) at 600 nm up to the value of approximately 0.6 and then the protein production was induced by 0.75 mM IPTG (isopropil-tio-β-Dgalattoside). After 3 hours of induction, bacteria were collected and treated as previously described (Gallo et al. 2012). Briefly, pellets were cleaned 3 times with a washing solution (10 mM Tris-HCl pH 7.5, 5 mM EDTA, 1 mM DTT, 0.5 mM PMSF, 4 mg/ ml leucopeptin, 0.7 mg/ml pepstatin, 5 mg/ml benzamidin). Following re-suspension in washing solution containing 0.3% SDS, bacterial cells were disrupted by sonication on ice (output control 4, 4 x 15 s, Vibra Cell, Sonics & Materials, Inc., Newtown, CT, USA). Samples were then boiled (5 minutes) and rapidly cooled down on ice (15 minutes). DNase (100 µg/ml) and RNase (50 µg/ml) were added in ice for 20 minutes. Cell debris and non-broken cells were separated by centrifugation at 15,000 x g, for 15 minutes, at 4°C. Protein solutions were precipitated by using 3 volumes of acetone, at -20°C, overnight, and re-suspended in mQ H₂O.

Recombinant AMA-1 protein was thus purified using HisGraviTrapTM kit (GE HealthCare, Chalfont, niUted Little Kingdom) following manufacturer's instruction, and then digested with AcTEV Protease (Life Technologies Corporation™, Carlsbad, CA, United States) for 3 hours at 30°C to remove the His Tag from the fusion protein. After the digestion, proteins were newly purified with HisGraviTrapTM kit (GE HealthCare, Little Chalfont, Buckinghamshire, United Kingdom) to eliminate the protease, the undigested recombinant proteins and the cleaved tags. Protein extracts were analysed by SDS-gel electrophoresis after staining with Lumio[™] Green Detection Kit (Life Technologies Corporation[™], Carlsbad, CA, United States). The Lumio[™] Green Detection Reagent is non-fluorescent and becomes fluorescent upon binding to the tetracysteine Lumio[™] tag codified by the expression vector, allowing the specific visualization of the recombinant protein at UV light. As an alternative to the Lumio[™] Green Detection, the protein was stained with 0.005% Blue Coomassie (Coomassie Brillant Blue R-250) (w/v). The protein was quantified by Bradford protein assay.

B. bigemina AMA-1 based ELISA optimization

The obtained protein was used to set up an ELISA protocol. ELISA tests were performed in 96 well plates (MaxiSorp, NUNC, Roskilde, Denmark). Apical Membrane Antigen-1 protein in carbonate/ bicarbonate buffer (0.1 M NaHCO, in distilled water pH 8.2 Sigma-Aldreich, St. Louis, Missouri, United States) was used to coat plate wells at 4°C over nigth. Then, 5% BSA solution to saturate free sites of the wells, serum samples, and Anti-bovine lgG secondary antibody conjugated with Peroxidase (Sigma- Aldreich, St. Louis, Missouri, United States) were added in this order. After the addition of the specific substrate (SigmaFastTM OPD, Sigma-Aldreich, St. Louis, Missouri, United States), O.D. values were determined at 405 nm at spectrophotometer (Mod 680, Biorad, Hercules, CA, United States) or, if the plate could not be read immediately, 3M H₂SO₄ were added and the stopped reactions were read at 490 nm.

Different concentrations and incubation times for all of these reagents were tested to find the optimal conditions for the assay. To validate our results, serum samples were tested in duplicate and in all plates, foetal bovine serum (FBS) was used as negative control, while blank samples were included replacing sera with the same volume of Phosphate-buffered saline (PBS). Field serum samples were previously selected and tested by serological and molecular techniques, as described in the section 'Field samples tested in the study', in order to verify their positivity/ negativity to *B. bigemina* and their closely related organisms *B. bovis* and *Theileria annulata*.

One sample, which resulted positive at both serological and molecular tests, was then used as positive control. Another sample, which was negative at serological and molecular tests, was subsequently used as negative control. These 2 samples were included in all the plates to verify result reproducibility. Antibody titers were considered positive when they yielded an O.D. 405 nm value at least twice as high as the negative field sera.

Field samples tested in the study

In order to investigate the specificity of the new B. bigemina AMA-1 based ELISA, the assay was used to analyse field serum samples that had been previously tested by serological and molecular techniques and to verify their positivity/negativity to B. bigemina and its closely related organisms B. bovis and T. annulata. In particular, bovine serum samples were analysed using a commercial ELISA kit for the presence of antibodies against B. bigemina (Svanova Biotech AB, Uppsala, Sweden) and an IFA kit for the detection of antibodies against B. bovis (Fuller Laboratories, Fullerton, CA, United States) following the manufacturer's instructions. The presence of antibodies against T. annulata was assessed by immunofluorescence using a protocol modified from Burridge and Kimber (Burridge and Kimber 1972) so to find out the serum antibody titres to Theileria infection. For molecular analysis, DNA was extracted from EDTA-treated blood samples using the PureLink Genomic Mini kit (Life Technologies Corporation[™], Carlsbad, CA, United States) following the manufacturer's instructions. DNA samples were analysed by polymerase chain reaction (PCR) to detect the presence of DNA from Babesia bovis (Figueroa et al. 1993), B. bigemina (Figueroa et al. 1993), and T. annulata (D'Oliveira et al. 1995). PCRs were performed in a reaction buffer containing 1.5 mM MgCl₂, 0.2 mM dNTPs, forward and reverse primers at a concentration of 0.4 mM, and 0.025 U/µl of Tag DNA polymerase (5 U/µl) (Promega, Madison, WI, United States). For each reaction, a positive control consisting of pathogen DNA and a negative control in which DNA was replaced by water were used. Polymerase chain reaction products were visualized after agarose gel electrophoresis containing 10 μ g/ml ethidium bromide.

Following serological and molecular analyses, 136 sera were divided into 6 groups according to their positivity/negativity for *B. bigemina*, *B. bovis*, and *T. annulata* pathogens. In detail, Group 1 included samples positive for *B. bigemina*, *B. bovis* and *T. annulata* (26 samples). Group 3 comprised 16 samples positive for *B. bovis* and *B. bigemina* but negative for *T. annulata*. Group 2 comprehended samples that were positive only for *B. bovis* (26 samples), Group 4 comprehended samples positive only for *B. bigemina* (28 samples). Samples that were positive (22 samples) only for *T. annulata* were enclosed in Group 5 while Group 6 included 18 samples negative to all three pathogens.

Statistical analysis of the results was performed by ANOVA test. Threshold O.D. values were determined to discriminate *B. bigemina/B. bovis* positive from negative sera. The threshold O.D. values were obtained by analysing the distribution frequencies of the O.D. values normalized according to Gauss' curves.

Results

Recombinant *B. bigemina* AMA-1 in *E. coli* production

The region codifying AMA-1 protein was amplified by PCR and the amplification product was inserted into the expression vector. The correct orientation of the cloned fragment was confirmed by sequence analysis. Obtained recombinant plasmid was used to transform BL21 Star™ (DE3) E. coli cells (Life Technologies Corporation™, Carlsbad, CA, United States). Protein expression was induced by IPTG and, after 3 hours, pellets were collected. Obtained cell suspensions were lysed and treated to extract the total protein pattern that was separated by SDS-PAGE and visualized by staining with Lumio Green[™] (Life Technologies Corporation[™], Carlsbad, CA, United States) as shown in Figure 1A. Protein purification by affinity chromatography was verified by SDS-PAGE (Figure 1B).

B. bigemina AMA-1 ELISA results

Recombinant protein was used for the assessment of an ELISA test. The optimum conditions of the test included the coating of individual wells of a microtiter plate (Nunc) with 0.1 μ g/well of AMA-1 antigen overnight at 4°C in an antigen

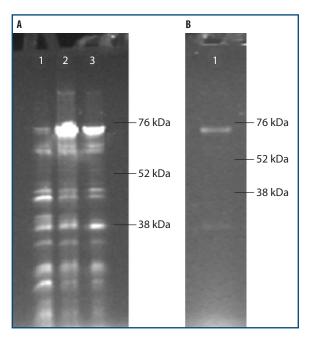


Figure 1. A. SDS-PAGE of total protein extracts from E.coli cells allowed to grow for 3 hours in absence (lane 1) or in presence (lane 2 and 3) of IPTG, stained with LumioTM Green Detection Kit (Life Technologies CorporationTM, Carlsbad, CA, United States) and observed at UV light. **B.** SDS-PAGE of the recombinant protein after purification by affinity chromatography, stained with LumioTM Green Detection Kit (Life Technologies CorporationTM, Carlsbad, CA, United States) and observed at UV light.

coating buffer (bicarbonate/carbonate at pH 8.2, Sigma-Aldreich, St. Louis, Missouri, United States). After 5 washes with PBS Tween, plated were saturated with 5% BSA in PBS (100 µl/well) for 2 hours and 30 minutes at 37°C. Plates were washed 5 times with PBS Tween, sera diluted (1:100) in PBS Tween were added and incubated for 1 hour at 37°C (100 µl/well). After 5 washes with PBS Tween, secondary antibody (Anti-bovine IgG-Peroxidase antibody, Sigma-Aldreich, St. Louis, Missouri, United States) diluted 1:10,000 in PBS-Tween was added and incubated for 1 hour at 37°C (100 µl/well). Plates were washed with PBS Tween 5 times and, then, fast OPD substrate (100 µl/well) was added (Sigma-Aldreich, St. Louis, Missouri, United States) and plates were read at a spectrophotometer at 405 nm after 20 minutes of incubation at room temperature. All the plates showed O.D. values of 0.1 for FBS and between 0.2 and 0.1 for blank samples.

The average values of each individual group are shown in Table I.

The statistical analysis performed by ANOVA test among the different groups confirmed a significant difference between samples that were positive for *B. bigemina* and/or *B. bovis* and samples that were negative for these 2 pathogens. No significant reaction was observed in presence of *T. annulata* positive sera.

Moreover, we did not find significant differences between the samples of Groups 1 and 3 positive for both *B. bovis* and *B. bigemina*, which were then grouped in Group A, so it was the case for Group 2 and Group 4 (*B. bovis* or *B. bigemina* positive, respectively), grouped in Group B, and samples of Group 5 and Group 6 (negative for both *B. bovis* and *B. bigemina*, positive or negative for *T. annulata*) unified in group C. It was not possible to unify the samples of Groups A and B into a single group because the values of the samples of Group A showed O.D. values significantly higher (1.15 \pm 0.12) than those of Group B (0.94 \pm 0.15).

Values obtained by the ELISA test for groups A, B, and C with their means, standard deviations, and number of samples are summarized in Table II.

Cut-off value to discriminate *B. bigemina/B. bovis* positive sera from the negative ones was recognised at 0.55 O.D., in such a way as to favour the sensitivity with respect to the specificity of the assay. In particular, it can be established that, following the specific optimised protocol, samples with O.D. lower than 0.55 can be considered reasonably negative. In contrast, values of O.D. greater than 0.69 indicate positivity to at least 1 of the 2 considered *Babesia*. The intermediate value of O.D. (between 0.55 and 0.69) could be considered doubtful, so for these

Table I. Results of B. bigemina Apical Membrane Antigen-1 ELISA assay performed on serum samples divided into 6 groups according to their positivity/ negativity to B. bigemina, B. bovis and T. annulata pathogens. The field serum samples included in the study were taken from a panel of sera and DNA samples collected by the lstituto Zooprofilattico Sperimentale della Sicilia from 2011 to 2012 for diagnostic purposes. The samples had been previously tested by serological and molecular techniques.

Group number	B. bigemina	B. bovis	T. annulata	N° of samples	Samples tested by serological analyses	Samples tested by molecular analyses	Mean O.D. value	Standard Deviation
1	POS	POS	POS	26	26	1	1.15	0.11
2	NEG	POS	NEG	16	9	7	0.99	0.18
3	POS	POS	NEG	26	10	16	1.15	0.14
4	POS	NEG	NEG	28	21	7	0.9	0.13
5	NEG	NEG	POS	22	10	12	0.4	0.13
6	NEG	NEG	NEG	18	12	6	0.42	0.11
			Total	136	87	49		

samples it may be appropriate to repeat the test on a withdrawal made after 2 weeks.

Discussion and conclusions

This work describes the optimisation of a prokaryotic system for AMA-1 recombinant protein production and the use of the produced antigen for a new diagnostic ELISA test. Obtained results showed a significant signal difference among the analysed serum samples, confirming the presence of B. bigemina antibodies directed against AMA-1 protein in infected bovine sera and suggesting that B. bigemina AMA-1 protein is sufficiently immunogenic and conserved to be recognised by antibodies against B. bovis, as it can be observed also in the amino acidic alignment showed in Figure 2. The higher O.D. values of the samples of Group A are probably due to presence in these sera of antibodies to both B. bigemina and B. bovis able to determine a stronger response to the antigen.

The use of the negative controls (SBF) and blank samples showed that the values obtained were consistent and very low in all the plates, suggesting the absence of non-specific signals in the plate wells.

In this paper a new ELISA test reliable, reproducible, and specific for *B. bovis* and *B. bigemina* is described based on *B. bigemina* AMA-1 protein, which is an apically located protein shared by many Apicomplexa. These proteins show common features and are enclosed in the AMA-1 protein family whose members are characterised by the presence of a Signal Peptide for membrane targeting, a Trans-membrane helix (TM), and disulphure bonds dividing the protein in functional domains. Apical Membrane Antigen-1 members have a key role in the host red cell invasion process and are extensively studied for their possible inclusion in a subunit vaccine. Sequence analyses showed the presence of regions highly conserved **Table II.** Results of B. bigemina Apical Membrane Antigen-1 ELISA assay performed on serum samples divided into 3 groups obtained by the association of the samples of Group 1 and Group 3 (Group A), Group 2 and Group 4 (Group B), Group 5 and Group 6 (Group C), after ANOVA analysis. The field serum samples included in the study were taken from a panel of sera and DNA samples collected by the Istituto Zooprofilattico Sperimentale della Sicilia from 2011 to 2012 for diagnostic purposes. The samples had been previously tested by serological and molecular techniques.

Group	B. bigemina	B. bovis	T. annulata	N° of samples	Mean O.D. value	Standard Deviation
٨	POS	POS	POS	52	1.15	0.12
A	POS	POS	NEG	52		
В	NEG	POS	NEG	44	0.94	0.15
D	POS	NEG	NEG			
c	NEG	NEG	POS	40	0.41	0.12
C	NEG	NEG	NEG	40		0.12

among the different examined organisms (Torina *et al.* 2010), so the choice to include samples positive for *B. bovis* and *T. annulata* in the analysis is justified.

Obtained results suggest the use of this assay for screening purposes in order to enable a rapid diagnosis and lower costs of the reaction. In case of a negative reaction, in fact, it should not be necessary to carry out species-specific serological tests that would be performed only on those serum samples positive at the first screening.

Acknowledgements

Authors would like to thank Mrs. Rosa Filippi and Mr. Pippo Bono for their technical support.

Grant support

Research supported by the Italian Ministry of Health (IZSSI 11/10).

CLUSTAL 2.1 multip	le sequence alignment	
BbigeminaAMA1 BbovisAMA1	MQCILNRITLLATPVIFFLWLSTEISPTGCAFVTFQNEPTST MQLHYKMQSTSPKYNYKRMLCMVFVPVILSSFFAEDALASNSTLFAFHREPTNR **. :. :: : : :: : : : : : : : : : : : :	
BbigeminaAMA1 BbovisAMA1	RGTRRSSRSSRHQQATSSTSEAGAGEATERAGGRTAGSKLIPQTPWTRYMIKYDIARCHG RLTRRASR-GQLLNSRRGSDDTSESSDRYSGRSGGSKNSGQSPWIKYMQKFDIPRNHG * ***:** . * .*:.*:::* .**:.*** *:** *:** *:**	
BbigeminaAMA1 BbovisAMA1	SGIFVDIGGYEAVGNKYYRMPTGKCPVMGKVISLASGADLLEPISADNPRYRGLAFPETV SGIYVDLGGYESVGSKSYRMPVGKCPVVGKIIDLGNGADFLDPISSEDPSYRGLAFPETA ***:**:***:**.* ****.****:**:**.**.********	
BbigeminaAMA1 BbovisAMA1	IKHTGASAGALTNAGNI-HGNLSPVSAADLRKWGYKGNAVTNCAEYASNIVPGSDQRTKY VDSNIPTQPKTRGSSSVTAAKLSPVSAKDLRRWGYEGNDVANCSEYASNLIPASDKTTKY :: .:: .:: .:****** ***:************	
BbigeminaAMA1 BbovisAMA1	RYPFVYDGKEEMCYILYSPMQYNQGTRYCDEDGSAKEGPSSLLCMKPYKSEADAHLYYGS RYPFVFDSDNQMCYILYSAIQYNQGNRYCDNDGSSEEGTSSLLCMKPYKSAEDAHLYYGS *****:*:****************************	
BbigeminaAMA1 BbovisAMA1	ARIDPKWDQNCPMKPIKDAIFGTWVSGACVALESAFEEYVDSAEECAAILFENSAADVDI AKVDPDWEENCPMHPVRDAIFGKWSGGSCVAIAPAFQEYANSTEDCAAILFDNSATDLDI *::**.*:****:*************************	
BbigeminaAMA1 BbovisAMA1	DIDSERYNEISELYNGLKNLQLQQIAFSLFAPMAKSAASATLSKGVGMNWANYESETGVC EVVNEEFNELKELTSGLKRLNLSKVANAIFSPLSNVAGTSRISRGVGMNWATYDKDSGMC :: .*.:**:.** .***.*::* ::*:*:: *.:: *:::*:*******	
BbigeminaAMA1 BbovisAMA1	RILNATPTCLIINAGSLAMTALGSPLESDAINYPCHIDTLGYVEPRKRNSREDGDRNSGI ALINETPNCLILNAGSIALTAIGSPLEYDAVNYPCHIDTNGYVEPRAKTTNKYLDVPFEV ::* **.***:****:*:***** **:************	
BbigeminaAMA1 BbovisAMA1	TTALNMKTLKCTKYVHSKYSESCGTYYYCSEEKSGYLSRLYQFMCSHNVKKAAV-ISTAL TTALSTKTLKCNAYVHTKYSDSCGTYFLCSDVKPNWFIRFLHMIGLYNTKRIVIFVCCTT ****. *****. ***:***:***** **: *: *: ::: :*.*: :::	
BbigeminaAMA1 BbovisAMA1	VLLCLAIYWIYQRLWSTKKGRQHDDYDRLMSKYEYD-DVSHDNIEPEHQLRTDAYIWGEA TAIVLTI-WIWKRFIKAKKEPAPPSFDKYLSNYDYDTTLDADN-ETEQRLDSSAYSWGEA . : *:* **::*: .:** .:** :: *:**** :. ** *.*::* :.** ****	
BbigeminaAMA1 BbovisAMA1	AARPSDITFVHMAKLN 595 VQRPSDVTFVKLSKIN 605 . ****:***:::*:*	

Figure 2. Amino acidic alignment of the B. bigemina and B. bovis Apical Membrane Antigen 1 (AMA-1).

References

- Black M.W. & Boothroyd J.C. 2000. Lytic cycle of *Toxoplasma* gondii. Microbiol Mol Biol Rev, **64**, 607-623.
- Bock R., Jackson L., de Vos A. & Jorgensen W. 2004. Babesiosis of cattle. *Parassitologia*, **129**, Suppl. S2, 47-69.
- Burridge M.J. & Kimber C.D. 1972. The indirect fluorescent antibody test for experimental East Coast fever (*Theileria parva* infection of cattle). Evaluation of a cell culture schizont antigen. *Res Vet Sci*, **13**, 451-455.
- D'Oliveira C., van der Weide M., Habela M.A., Jacquiet P. & Jongejan F. 1995. Detection of *Theileria annulata* in blood samples of carriers by PCR. *J Clin Microbiol*, **33**, 1665-2669.
- Figueroa J.V., Chieves L.P., Johnson G.S. & Buening G.M. 1993. Multiplex polimerase chain reaction based assay for the detection of *Babesia bigemina*, *Babesia bovis* and *Anaplasma marginale* DNA in bovine blood. *Vet Parasitol*, **50**, 69-81.
- Gaffar F.R., Yatsuda A.P., Franssen F.F. & de Vries E. 2004. Erythrocyte invasion by *Babesia bovis* merozoites is inhibited by polyclonal antisera directed against peptides derived from a homologue of *Plasmodium falciparum* apical membrane antigen 1. *Infect Immun*, **72**, 2947-2955.
- Gallo G., Baldi F., Renzone G., Gallo M., Cordaro A., Scaloni A. & Puglia A.M. 2012 Adaptative biochemical pathways and regulatory networks in *Klebsiella oxytoca* BAS-10 producing a biotechnologically relevant exopolysaccharide during Fe(III)-citrate fermentation. *Microb Cell Fact*, **23**, 152.
- Kumar S., Dudley J., Nei M. & Tamura K. 2008. MEGA: a biologist-centric software for evolutionary analysis of DNA and protein sequences. *Brief Bioinform*, 9, 299-306.
- Larkin M.A., Blackshields G., Brown N.P., Chenna R.,

McGettigan P.A., McWilliam H., Valentin F., Wallace I.M., Wilm A., Lopez R., Thompson J.D., Gibson T.J. & Higgins D.G. 2007. ClustalW and ClustalX version. *Bioinformatics*, **23**, 2947-2948.

- Latitha P.V., Biswas S., Pillai C.R. & Saxena R.K. 2008. Immunogenicity of a recombinant malaria vaccine candidate, domain I + II of AMA-1 ectodomain, from Indian *P. falciparum* alleles. *Vaccine*, **26**, 4526-4535.
- Peterson M.G., Marshall V.M., Smythe J.A., Crewther P.E., Lew A., Silva A., Anders R.F. & Kemp D.J. 1989. Integral membrane protein located in the apical complex of *Plasmodium falciparum*. *Mol Cell Biol*, **9**, 3151-3154.
- Tamura K., Dudley J., Nei M. & Kumar S. 2007. MEGA4: molecular Evolutionary Genetics Analysis software version 4.0. *Mol Biol Evol*, 24, 1596-1599.
- Torina A. & Caracappa S. 2007. Babesiosis in Italy: an overview. *Parassitologia*, **49**, 23-28.
- Torina A., Agnone A., Sireci G., Mosqueda J.J., Blanda V., Albanese I., La Farina M., Cerrone A., Cusumano F. & Caracappa S. 2010. Characterization of the Apical Membrane Antigen-1 in Italian strains of *Babesia bigemina*. *Transbound Emerg Dis*, **57**, 52-56.
- Xia X. & Xie Z. 2001. DAMBE: data analysis in molecular biology and evolution. *J Hered*, **92**, 371-373.
- Yokoyama N., Okamura M. & Igarashi I. 2006. Erythrocyte invasion by *Babesia* parasites: current advances in the elucidation of the molecular interactions between the protozoan ligands and host receptors in the invasion stage. *Vet Parasitol*, **138**, 22-32.
- Zhou J., Yang J., Zhang G., Nishikawa Y., Fujisaki K. & Xuan X. 2006. *Babesia gibsoni*: an apical membrane antigen-1 homologue and its antibody response in the infected dogs. *Exp Parasitol*, **114**, 329-333.