



The University of  
**Nottingham**

UNITED KINGDOM · CHINA · MALAYSIA

Llibre, Alba and Shimakawa, Yusuke and Mottez, Estelle and Ainsworth, Shaun and Buivan, Tan-Phuc and Firth, Rick and Harrison, Elliott and Rosenberg, Arielle R. and Meritet, Jean-Francois and Fontanet, Arnaud and Castan, Pablo and Madejon, Antonio and Laverick, Mark and Glass, Allison and Viana, Raquel and Pol, Stanislas and McClure, C. Patrick and Irving, William Lucien and Miele, Gino and Albert, Matthew L. and Duffy, Darragh (2018) Development and clinical validation of the Genedrive point-of-care test for qualitative detection of hepatitis C virus. Gut . ISSN 1468-3288

**Access from the University of Nottingham repository:**

<http://eprints.nottingham.ac.uk/51091/1/gutjnl-2017-315783.pdf>

**Copyright and reuse:**

The Nottingham ePrints service makes this work by researchers of the University of Nottingham available open access under the following conditions.

This article is made available under the Creative Commons Attribution Non-commercial licence and may be reused according to the conditions of the licence. For more details see: <http://creativecommons.org/licenses/by-nc/2.5/>

**A note on versions:**

The version presented here may differ from the published version or from the version of record. If you wish to cite this item you are advised to consult the publisher's version. Please see the repository url above for details on accessing the published version and note that access may require a subscription.

For more information, please contact [eprints@nottingham.ac.uk](mailto:eprints@nottingham.ac.uk)



OPEN ACCESS

ORIGINAL ARTICLE

# Development and clinical validation of the Genedrive point-of-care test for qualitative detection of hepatitis C virus

Alba Llibre,<sup>1,2</sup> Yusuke Shimakawa,<sup>3</sup> Estelle Mottez,<sup>4,5</sup> Shaun Ainswroth,<sup>6</sup> Tan-Phuc Buivan,<sup>4,7</sup> Rick Firth,<sup>6</sup> Elliot Harrison,<sup>6</sup> Arielle R Rosenberg,<sup>8</sup> Jean-François Meritet,<sup>8</sup> Arnaud Fontanet,<sup>3,9</sup> Pablo Castan,<sup>10</sup> Antonio Madejón,<sup>11</sup> Mark Laverick,<sup>6</sup> Allison Glass,<sup>12</sup> Raquel Viana,<sup>12</sup> Stanislas Pol,<sup>1,2,4,7</sup> C Patrick McClure,<sup>13</sup> William Lucien Irving,<sup>13</sup> Gino Miele,<sup>6</sup> Matthew L Albert,<sup>1,2,4,14</sup> Darragh Duffy<sup>1,2,4</sup>

► Additional material is published online only. To view please visit the journal online (<http://dx.doi.org/10.1136/gutjnl-2017-315783>).

For numbered affiliations see end of article.

## Correspondence to

Darragh Duffy, Immunobiology of Dendritic Cells, Institut Pasteur, Paris 75015, France; [darragh.duffy@pasteur.fr](mailto:darragh.duffy@pasteur.fr)

AL, YS, MLA and DD contributed equally.

Received 30 November 2017

Revised 9 March 2018

Accepted 14 March 2018

## ABSTRACT

**Objective** Recently approved direct acting antivirals provide transformative therapies for chronic hepatitis C virus (HCV) infection. The major clinical challenge remains to identify the undiagnosed patients worldwide, many of whom live in low-income and middle-income countries, where access to nucleic acid testing remains limited. The aim of this study was to develop and validate a point-of-care (PoC) assay for the qualitative detection of HCV RNA.

**Design** We developed a PoC assay for the qualitative detection of HCV RNA on the PCR Genedrive instrument. We validated the Genedrive HCV assay through a case-control study comparing results with those obtained with the Abbott RealTime HCV test.

**Results** The PoC assay identified all major HCV genotypes, with a limit of detection of 2362 IU/mL (95% CI 1966 to 2788). Using 422 chronically infected patients with HCV and 503 controls negative for anti-HCV and HCV RNA, the Genedrive HCV assay showed 98.6% sensitivity (95% CI 96.9% to 99.5%) and 100% specificity (95% CI 99.3% to 100%) to detect HCV. In addition, melting peak ratiometric analysis demonstrated proof-of-principle for semiquantification of HCV. The test was further validated in a real clinical setting in a resource-limited country.

**Conclusion** We report a rapid, simple, portable and accurate PoC molecular test for HCV, with sensitivity and specificity that fulfils the recent FIND/WHO Target Product Profile for HCV decentralised testing in low-income and middle-income countries. This Genedrive HCV assay may positively impact the continuum of HCV care from screening to cure by supporting real-time treatment decisions.

## INTRODUCTION

Chronic infection with hepatitis C virus (HCV) is a major public health problem, estimated to infect 1.0% of the world's population (71 million people)<sup>1</sup> and to be responsible for 400 000 annual deaths as a result of cirrhosis and liver cancer.<sup>2</sup> With the advent of highly potent direct-acting antiviral (DAA) therapy, cure is now possible in >95% of treated individuals achieving a

## Significance of this study

### What is already known on this subject?

- Recent direct-acting antiviral therapies have provided transformative therapeutic options for chronic hepatitis C virus (HCV) infection, yet the challenge remains to identify the 1% of the world's population that is chronically infected with HCV.
- Current WHO guidelines highlight the lack of field-based platforms for point-of-care (PoC) HCV RNA testing.
- The Cepheid Xpert HCV Viral Load assay (Cepheid) is currently the only CE-IVD-certified assay for decentralised HCV viral load determination, yet according to the WHO, it presents many limitations.

### What are the new findings?

- We developed a PoC assay for the qualitative detection of HCV RNA on the PCR Genedrive instrument. The PoC assay identified all six major HCV genotypes, with a limit of detection of 2362 IU/mL, 98.6% sensitivity (95% CI 96.9% to 99.5%) and 100% specificity (95% CI 99.3% to 100%) to detect HCV viraemia.
- The test was also validated in a real-life clinical setting in a resource-limited country.
- The Genedrive instrument is a 600g portable device that can be battery operated, thus making it highly suitable for decentralised testing in field settings, and requires only 30 µL of sample.
- We demonstrated a proof of concept for semiquantification of HCV viral load using melting peak ratiometric analysis.

sustained virological response,<sup>3–7</sup> thus changing the landscape of clinical HCV management in recent years. Consequently, the WHO published in 2016 a global strategy to eliminate HCV as a public health threat by 2030.<sup>8</sup>

To achieve this ambitious goal, it is essential to scale up diagnosis and treatment capacities in low-income and middle-income countries,



**To cite:** Llibre A, Shimakawa Y, Mottez E, et al. Gut Epub ahead of print: [please include Day Month Year]. doi:10.1136/gutjnl-2017-315783

AUTHOR PROOF

## Significance of this study

**How might it impact on clinical practice in the foreseeable future?**

► The Genedrive instrument provides a rapid, simple, portable and accurate PoC molecular test for HCV and demonstrated sensitivity and specificity fulfilling the requirements published by FIND, the Foundation for Innovation in New Diagnostics, and WHO for decentralised testing in low- and middle-income countries. Thus, the Genedrive HCV assay could help to support the decentralisation of HCV testing through rapid diagnosis, resulting in improved patient outcome and public health. In particular its portability, rapid results, and ease of use make it ideal for confirming infection cases in population-based screening approaches. This is crucial for ensuring that diagnosis translates into treatment and that patients in low- and middle-income settings undergo screening and treatment initiation during a single medical intervention. This highly sensitive and specific test has recently obtained CE-IVD certification and is therefore positioned to impact the continuum of HCV care from screening to cure through real-time treatment decisions.

particularly in Africa and Central Asia where HCV is highly endemic.<sup>9</sup> However, in these regions, access to HCV diagnostic tools is severely limited.<sup>10</sup> The current diagnostic algorithm is based on serological screening of HCV-specific antibodies (two different positive serological tests may be mandatory in some countries), followed by nucleic acid testing to detect presence of HCV RNA.<sup>11 12</sup> Although rapid serological diagnostic tests exist for anti-HCV antibodies,<sup>10 13</sup> nucleic acid testing for HCV RNA often involves dedicated facilities and highly qualified personnel, even if platforms such as Cepheid, Panther or VERIS/MDx are changing such requirements. In low-income and middle-income countries, HCV RNA tests are currently only available in centralised laboratories, resulting in less than 1% of infected people in these regions being aware of their infection.<sup>14</sup> In certain at-risk patient subgroups (eg, migrants and prisoners) with current testing practices, it requires multiple visits to make a diagnosis, and as a result, many positives cases are missed. Therefore, to improve the 'suboptimal' cascade from screening to linkage to care and efficient therapy, there is an urgent unmet clinical need for decentralised HCV RNA testing.<sup>11</sup> Point-of-care (PoC) molecular technologies have emerged as a viable strategy for improving clinical management of HCV, as they facilitate diagnosis and treatment programmes with the possibility to work in rural, resource-limited settings. In addition to decentralised nucleic acid amplification testing (NAAT), serological rapid diagnostic tests are also contributing to the improvement of diagnosis strategies.

We developed the Genedrive HCV assay for the detection of HCV RNA and conducted analytical and clinical validation studies to assess its diagnostic accuracy. We used plasma and serum samples from a large cohort of treatment-naïve patients with chronic HC and HCV-negative controls at two centres (Institut Pasteur, Paris, France, and National Health Service (NHS), Nottingham, UK). Samples were assayed using the newly developed PoC assay, and results were compared with the commercial Abbott RealTime HCV Assay. We further validated this new PoC test in a resource-limited setting (Lancet Laboratories, Johannesburg, South Africa).

## METHODS

**Genedrive HCV assay development**

The Genedrive HCV assay is a two-step procedure requiring a plasma or serum preparation step, followed by a reverse transcription (RT) reaction to generate the complementary DNA (cDNA) from the target HCV RNA. This cDNA undergoes asymmetric PCR to generate linear amplification of single stranded products, followed by detection using a secondary hybridisation probe and dissociation curve analysis. The Genedrive HCV assay amplifies and detects a 91 base-pair from the 5' untranslated region of the HCV genome.<sup>15</sup> The pan-genotypic HCV target selection was performed through an iterative process comprising 339 sequences from the Los Alamos National Security Database, representing all major HCV genotypes (1a and 1b, 2, 3, 4, 5 and 6) (online supplementary table S1 A). Diagnostic accuracy was shown for all six major genotypes (online supplementary table S1 B).

To enable PoC testing, the Plasma Preparation Cartridge contains lyophilised reagents for plasma processing (protease and reagents to prevent thermal-induced coagulation of plasma proteins), as well as the internal positive control (IPC) oligonucleotides and template (pUCV57 plasmid). The IPC reaction targets a DNA sequence located in the backbone of the pUCV57, which shows non-cross reactivity for all HCV genotypes.

**Genedrive HCV assay workflow**

Following dilution with nuclease free water (1:2 ratio) of the plasma, 15 µL was added to each of the three channels of the plasma preparation cartridge. The protease incubation step was then performed within the Genedrive instrument, consisting of a 5 min incubation at 37°C followed by 5 min at 95°C, thus rendering the sample non-infectious. During this incubation time, 100 µL of nuclease free water was added to resuspend the RT-PCR reagents. After removal of the cartridge from the device, 30 µL of the RT-PCR suspension was dispensed into each channel, which was resealed with a new lid. The cartridge was placed back in the Genedrive instrument, and the programme resumed. The target RNA was converted to DNA (via RT) and subsequently amplified by asymmetric PCR. Following amplification, a fluorescent probe was used to detect target-specific sequences by monitoring changes in fluorescent signal intensity that occurred during dissociation of the fluorescent probe from its hybridised target sequence (if present) as temperature increased. The Genedrive instrument detects this fluorescence at defined melt temperature positions (T<sub>m</sub>) for both the HCV target and for the IPC. A schematic workflow of the Genedrive HCV assay and a representative outcome for a positive and a negative sample are illustrated in online supplementary figures S1 and S2. The melting temperatures (T<sub>m</sub>) and peak height values for both the IPC and the HCV peak are also illustrated (online supplementary figure S2). Each of the three channels can result in either one or two separate melt peaks: the IPC peak and the HCV peak (if positive). The calling logic for determining a result and the four possible outcomes from the sample test replicates and IPC are depicted in online supplementary figures S3 and S4. In all cases, a *positive/negative* outcome on the first test run was taken as a result. *Control fail-retest* and *indeterminate-retest* outcomes were retested, and the second outcome was taken as the final result. If that second run returned an *indeterminate/control fail*, no further retests were performed, and no result was obtained. The instrument runs for 88 min (from 1 hour 40 min to 2 hours including experimental manipulation).

The runs were monitored, and the data were saved using the Genedrive Engineering software.

### Analytical validation

Analytical specificity was assessed using a group of HCV RNA negative plasma samples, some of which were positive for hepatitis B virus (HBV), HIV-1 or flaviviral samples (dengue, yellow fever, Zika and GBV-C) as determined by serology, PCR or both. These specimens were obtained at either NHS, Nottingham or Instituto de Salud Carlos III & Hospital la Paz, Madrid, as part of routine diagnostic procedures, with clearance for samples surplus to diagnostic use being permitted for assay development and validation.

Analytical sensitivity was assessed for each of six major HCV genotypes using 16 commercial standard plasma panels (Seracare (USA) and Tissue Solutions (UK)). Each sample was diluted to 5000, 3000 and 1000 IU/mL, and at each concentration, 24 replicates were tested to obtain the limit of detection (LoD) using probit analysis. To examine the effect of anticoagulant on the LoD, analytical sensitivity was further assessed using synthetic, stable, homogenous and non-infectious armoured RNA in EDTA and Na-heparin plasma (Assuragen, Quant HCV GT 2b).<sup>16</sup>

### Case-control study

We defined cases as patients with treatment-naïve chronic HCV infection (anti-HCV positive and HCV RNA positive). At the Institut Pasteur (Paris, France), pretreatment heparinised plasma samples were obtained from three approved cohort studies in which consent was obtained for secondary use of samples: Inserm C10-08, Inserm C10-54 and ANRS Cupic CO20. The negative control group consisted of samples from healthy donors, negative for both anti-HCV antibody and HCV RNA. Control heparinised plasma samples were obtained from the CoSImmGEn cohort of the *Investigation Clinique et Accès aux Ressources Biologiques* platform (Centre de Recherche Translationnelle, Institut Pasteur, Paris) and the *Etablissement français du sang*. Nottingham samples were derived from those sent to the diagnostic laboratory for routine investigation of HCV or other blood-borne virus infection that were surplus to requirement and could be used for assessment of an in vitro assay. The daily positive controls were obtained from Seracare (USA) and used at a viral load (VL) of 10 000 IU/mL. The negative controls were from a healthy donor plasma (K2EDTA) internally sourced by genedrive plc. Samples were stored at  $-80^{\circ}\text{C}$ , thawed, aliquoted (operators were blinded to infection status) and kept frozen at  $-80^{\circ}\text{C}$  prior to use for both Genedrive and Abbott RealTime HCV Assay testing. Additionally, to test a potential freeze cycle effect, matched fresh and frozen plasma from HCV RNA positive and negative individuals were assessed. The study was sponsored and approved by the ANRS and registered with clinicaltrials.gov (NCT02992184: PoC-HCV Genedrive Viral Detection Assay Validation Study).

### Reference standard test for HCV RNA

VLs were measured by RT real-time PCR using the Abbott RealTime HCV assay, according to manufacturer's instructions. The Abbott m2000sp and m2000rt instruments were used for automated sample preparation and real-time amplification and detection, respectively. The lower limit of quantification of the assay is 12 IU/mL. HCV genotyping was performed based on NS5b region sequence phylogenetic analysis using the Abbott RealTime HCV Genotype II on the Abbott m2000, following manufacturer's instructions.

### Robustness, repeatability and stability

Daily controls (positive and negative samples) were performed at both sites, and results were as expected (online supplementary table S2). Nine hundred and twenty-five case and control samples were tested across 12 Genedrive instruments, five operators and three different assay batches. Although our study was not specifically designed and powered for this analysis, no effect of the Genedrive instruments, the operator and the assay batch were observed on assay performance. To determine the effect of interfering substances on assay performance, 37 different substances including a wide range of antiviral drugs and several endogenous substances such as bilirubin or haemoglobin were tested using HCV positive plasma at  $3 \times \text{LoD}$  (online supplementary table S3).

### Real-life clinical setting validation

A total of 130 plasma and serum samples were tested across three Genedrive instruments and four operators. They were collected as part of a routine HCV RNA diagnostic testing using Abbott RealTime HCV Genotype II on the Abbott m2000, in different African countries: Ghana (9 cases, 3 controls), Kenya (37 cases, 1 control), Mauritius (18 cases), Mozambique (1 case), Nigeria (3 cases), South Africa (33 cases, 6 controls), Tanzania (10 cases, 2 controls), Uganda (1 case), Zambia (1 case) and Zimbabwe (5 cases). Samples were stored at  $-70^{\circ}\text{C}$  prior to testing.

### Data analysis

The accuracy of the Genedrive HCV assay to detect HCV viraemia was assessed by diagnostic sensitivity and specificity, and positive and negative likelihood ratios. Factors associated with the lack of result or false-negative results were identified using Wilcoxon rank-sum test for continuous and Fisher's exact test for categorical variables. All the variables found to be significantly associated with the lack of result in the crude analysis ( $P < 0.05$ ) were further assessed in the multivariable logistic regression. The sample size of the case-control study was determined based on the FIND/WHO Target Product Profile indicating the minimally acceptable sensitivity and specificity as 95% and 98%, respectively.<sup>17,18</sup> Four hundred and nine cases were required to demonstrate that sensitivity is at least higher than 95% at a two-sided significance level of 5% with a power of 90% when the true sensitivity was 98%. Similarly, 495 controls were needed to show that the specificity is at least higher than 98% with a power of 80%, when the true specificity was 99.5%.

To assess the performance of the Genedrive HCV assay to quantify/semiquantify HCV RNA levels, we split the HCV RNA-positive cases into a derivation (Paris site) and validation set (Nottingham site). The correlation between HCV/IPC peak ratio (calculated as the mean value of ratios of HCV peak and IPC peak at each channel) and VL was evaluated using Pearson's correlation coefficient. To predict VL, a conversion factor from the HCV/IPC peak ratio to HCV RNA levels ( $\log_{10}$  IU/mL) was developed by fitting a linear regression model in the Paris cohort. The conversion factor was then applied to the HCV/IPC peak ratio of the Nottingham cohort, and the agreement between the observed HCV RNA levels using RT-PCR and the estimated HCV RNA levels from the Genedrive was assessed using Bland and Altman plots. All analyses were performed using STATA, V.13.0. The study was reported in accordance with the standards for reporting of diagnostic accuracy checklists.<sup>19</sup>

**Table 1** Study sample characteristics

Study site	Institut Pasteur (Paris, France)	Queen's Medical Centre (Nottingham, UK)	Lancet Laboratories (Johannesburg, South Africa)
Gold standard assay	Abbott RealTime	Abbott RealTime	Abbott RealTime
HCV-positive cases (n)	221	201	118
Controls (n)	190	313	12
Condition of sample storage	Frozen	Frozen (201 cases and 313 controls) Fresh (51 cases and 49 controls)	Frozen
Type of sample	Plasma	Plasma (124 cases and 234 controls) Serum (77 cases and 79 controls)	Plasma/serum
Anticoagulant	Na-heparin	EDTA and SST tubes	Information not recorded
Genotype distribution	g1a and g1b (88%) and g4 (12%)	g1a and g1b (62%), g3 (33%), g2 (2.2%), g4 (2.2%), g5 (1.1%)	Data available for 98% of the cases: g1a (40.5%), g4 (27.5%), g1b (15.5%), g3 (6%), g2 (5%), g5 (5%).
No of daily controls	1 positive and 3 negative	4 positive and 4 negative	None
No of operators	2	3	4
No of Genedrive units	4	8	3

HCV, hepatitis C virus.

## RESULTS

### Analytical sensitivity

The Genedrive HCV assay detected all major HCV genotypes, with a LoD ranging from 1406 IU/mL (genotype 6) to 3203 IU/mL (genotype 5) (online supplementary table S4). The LoD of the 16 standard plasma samples was 2362 IU/mL (95% CI 1966 to 2758). Using synthetic armoured RNA, the LoD in plasma derived from K2 EDTA collected blood (1918 IU/mL, 95% CI 1789 to 2205) was significantly lower ( $P < 0.0001$ ) than the LoD in Na-heparin derived plasma (2359 IU/mL, 95% CI 2147 to 2826).

### Analytical specificity

Analytical specificity was examined for 17 different pathogens including samples infected with HIV-1, HBV and various *Flaviviridae* (including: dengue, Zika, yellow fever and GBV-C) (online supplementary table S5). The Genedrive HCV assay was

shown to be non-cross-reactive for any of the aforementioned samples with the sole exception of one positive out of two replicate runs on a GBV-C sample.


### Clinical validation study (case-control study)

A summary of all samples tested in all study sites is illustrated in table 1. The characteristics of case and control subjects including age, sex, median HCV RNA values and HCV genotype are summarised in table 2. For the first part of the study performed in Western laboratories (Institut Pasteur, Paris, and Queen's Medical Centre, Nottingham), a total of 422 cases and 503 negative controls were included, and the flow chart of the study samples is illustrated in figure 1. Results were obtained at the first attempt for 97.2% of the samples (899/925). Consequently, 26 samples were retested, of which 61.5% (16/26) returned a result with the second test. The 10 samples (4 cases/6 controls) that failed such retesting were excluded from subsequent analysis

**Table 2** Characteristics of cases and control samples

	European study			African study		
	Cases (n=422)	Control (n=503)	P values	Cases (n=118)	Control (n=12)	P values
Mean age (SD)	50.4 (14.3)	47.9 (15.7)	0.02	ND	ND	N/A
Male sex (%)	254/380 (66.8)	272/433 (62.8)	0.2	ND	ND	N/A
Study site						
Paris	221 (52.4)	190 (37.8)	<0.001	N/A	N/A	N/A
Nottingham	201 (47.6)	313 (62.2)		N/A	N/A	N/A
Positive anti-HCV at baseline (%)	422/422 (100)	0/503 (0)	1.0	ND	ND	N/A
Mean HCV RNA (log <sub>10</sub> IU/mL) (SD)	5.7 (1.0)	N/A	N/A	5.3 (0.8)	N/A	N/A
HCV RNA (log <sub>10</sub> IU/mL) (%)						
<3.0	13 (3.1)	N/A	N/A	0	N/A	N/A
3.0–6.0	227 (53.8)	N/A		99 (83.9)	N/A	
≥6.0	182 (43.1)	N/A		19 (16.1)	N/A	
HCV genotype (%)						
1	301 (71.3)	N/A	N/A	65 (55.1)	N/A	N/A
2	4 (1.0)	N/A		6 (5.1)	N/A	
3	59 (14.0)	N/A		7 (5.9)	N/A	
4	31 (7.3)	N/A		32 (27.1)	N/A	
5	2 (0.5)	N/A		6 (5.1)	N/A	
Not determined	25 (5.9)	N/A		2 (1.7)	N/A	
HIV infection	0/422	29/503		ND	ND	N/A
HBV infection	0/422	25/503		ND	ND	N/A

HCV, hepatitis C virus; N/A, not applicable; ND, no data.

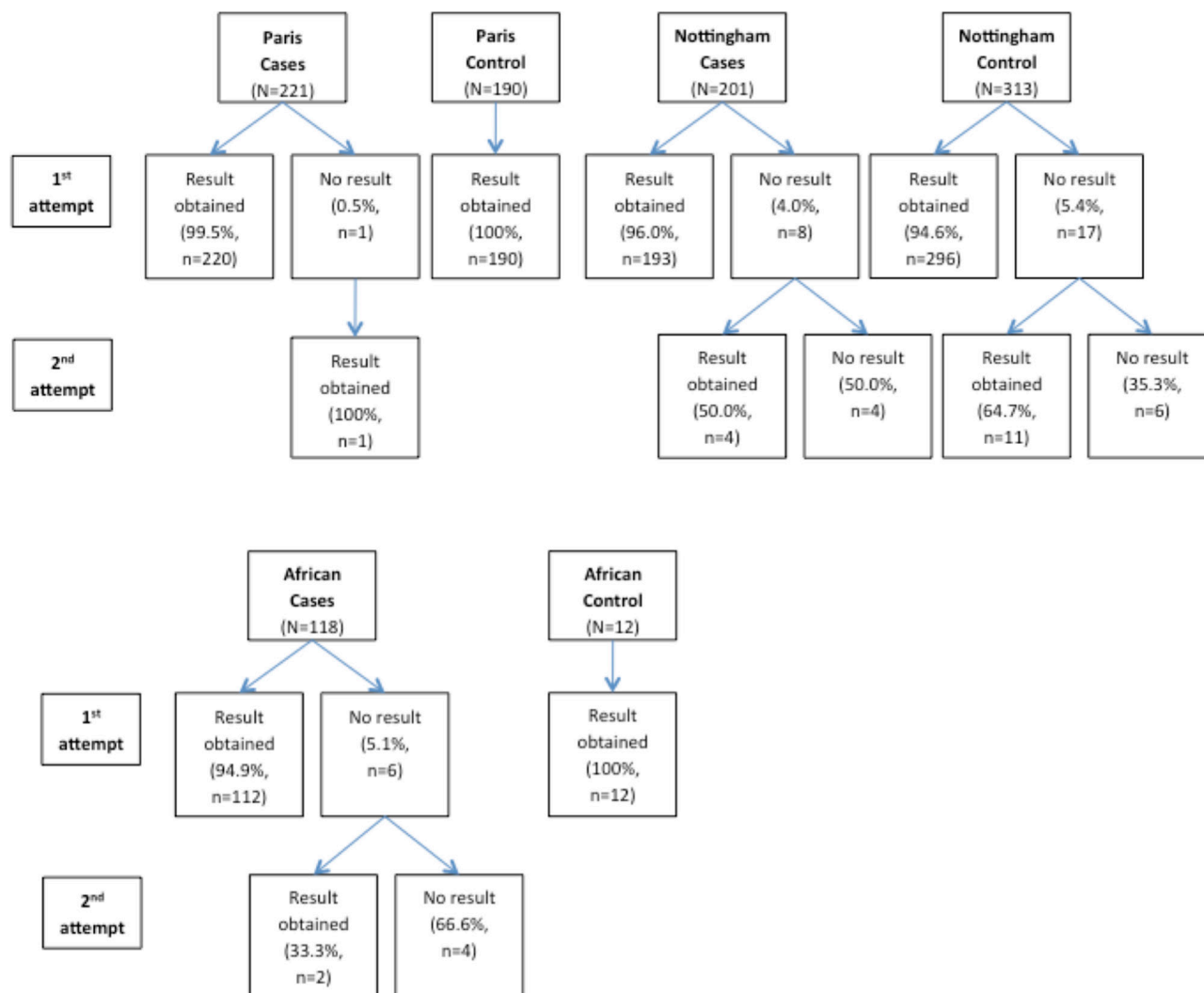

**Figure 1** Flow chart of the study participants.

for sensitivity and specificity. Multivariable analysis identified low HCV RNA levels as the sole independent factor significantly associated ( $P < 0.001$ ) with the lack of a result at the first attempt in the HCV RNA positive case samples (online supplementary table S6). The median HCV RNA level (IQR) was  $5.9 \log_{10}$  IU/mL ( $5.4$ – $6.3$ ) in cases with results and  $2.9 \log_{10}$  IU/mL ( $2.6$ – $3.5$ ) in cases without results. No factor was identified as associated with a lack of result in the control samples.

#### Diagnostic performance of the HCV Genedrive assay

The diagnostic sensitivity and specificity of the Genedrive HCV assay to detect HCV RNA in Western laboratories was 98.6%

(412/418; 95% CI 96.9% to 99.5%) and 100% (497/497; 95% CI 99.3% to 100%), respectively (table 3). Apart from one sample with a VL of 2508 IU/mL, all of the false-negative samples were below 1000 IU/mL. Low HCV RNA level was the only factor significantly associated ( $P < 0.0001$ ) with false-negative results; age, sex, study site, sample type or viral genotype were not found to influence the accuracy of Genedrive HCV assay (online supplementary table S7). No differences in sensitivity were observed between the two sites, which tested samples with difference anticoagulants. To examine for potential effects of a freeze/thaw cycle on the assay, the diagnostic performance was further assessed on 100 freshly collected (non-frozen)


**Table 3** Diagnostic performance of Genedrive using frozen and fresh samples

	Frozen samples from Europe (n=915)*		Fresh samples from Europe (n=96)†		Frozen samples from Africa (n=126)‡	
	Point estimate	95% CI	Point estimate	95% CI	Point estimate	95% CI
Sensitivity	98.6% (412/418)	96.9% to 99.5%	98.0% (48/49)	89.1% to 99.9%	100% (114/114)	96.8% to 100%
Specificity	100% (497/497)	99.3% to 100%	100% (47/47)	92.5% to 100%	100% (12/12)	73.5% to 100%
Positive likelihood ratio	N/A	N/A	N/A	N/A	N/A	N/A
Negative likelihood ratio	0.014	0.007 to 0.032	0.020	0.003 to 0.142	0.000	N/A

\*Excluding four cases and six controls who did not give any result after two attempts.

†Excluding two cases and two controls who did not give any result after two attempts.

‡Excluding four cases who did not give any result after two attempts.

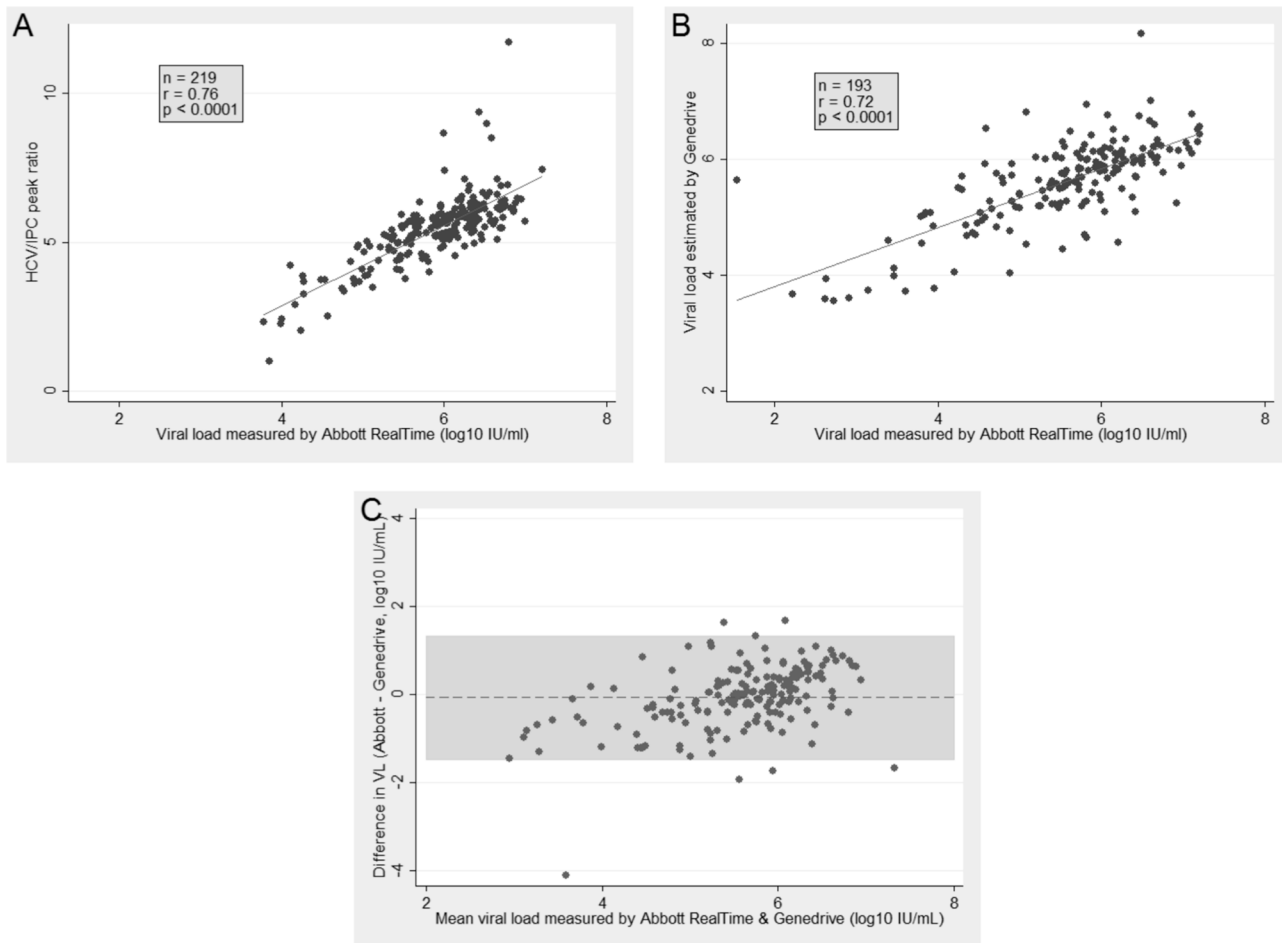

**Figure 2** Comparison between the HCV/IPC peak ratio (Genedrive) and HCV RNA levels (Abbott RealTime). (A) Correlation between the HCV/IPC peak ratio and the viral load measured by Abbott RealTime in the Paris cohort (derivation set). (B) Correlation between the viral load estimated by Genedrive and viral load measured by Abbott RealTime in Nottingham cohort (validation set). (C) Agreement between the viral loads estimated by Genedrive and viral loads observed by Abbott RealTime in Nottingham cohort (validation set).

plasma samples. After excluding four samples (two cases/two controls) that did not give a result, the sensitivity and specificity was 98.0% (48/49; 95% CI 89.1% to 99.9%) and 100% (47/47; 95% CI 92.5% to 100%), respectively (table 3), with one false negative sample that had a VL of 230 IU/mL. This demonstrates that the accuracy of the assay was not affected by a freezing step.

### HCV VL semiquantification

In the derivation data set (Paris cohort), the mean of the HCV/IPC peak ratio was closely correlated with VL as measured by the Abbott platform (Pearson's correlation coefficient:  $r=0.76$ ,  $P<0.0001$ ) (figure 2A). By fitting a linear regression, we obtained a conversion factor from the mean of the HCV/IPC peak ratio to the HCV RNA levels. The estimated  $\log_{10}$  VL was equal to  $3.5+0.4 \times$  HCV/IPC peak ratio. This conversion factor was then applied to the mean HCV/IPC peak ratio in the validation set, for which historical VL measurements were used. There was a linear correlation between the VL measured by the Abbott RealTime platform and the value estimated using the Genedrive results (Pearson's correlation coefficient:  $r=0.72$ ,  $P<0.0001$ ) (figure 2B). The agreement between the observed and estimated VL was assessed using Bland and Altman plots; mean difference was  $-0.06$   $\log_{10}$  IU/mL (SD:  $\pm 0.72$ ; 95% limits of agreement:  $-1.47, 1.33$ ) (figure 2C). Low HCV VL was associated with non-agreement (supplementary table S8). These

results support the potential use of the Genedrive HCV assay to semiquantify VL.

### Real-life clinical setting validation

In order to further assess the diagnostic accuracy of the Genedrive HCV assay outside Western laboratories, we extended our study to a real-life clinical setting. Samples tested in Johannesburg (South Africa) are summarised in tables 1 and 2. A total of 118 cases and 12 negative controls were included, and the flow chart of the study samples is illustrated in figure 1. Results were obtained at the first attempt for 95.4% of the samples (124/130). Consequently, six samples were retested, of which 33.3% (2/6) returned a result with the second test. The four samples (four cases) that failed such retesting were excluded from subsequent analysis for sensitivity and specificity. All the samples that required retesting were positive for HCV RNA. As summarised in table 3, the Genedrive HCV assay demonstrated 100% sensitivity and 100% specificity when compared with the Abbott m2000 in real-life conditions.

### DISCUSSION

We developed and validated the Genedrive HCV assay for the qualitative detection of HCV RNA. The Genedrive instrument is a handheld rapid thermocycler that permits NAAT on plasma

or serum samples, which fulfils the target product profile for HCV diagnosis drafted by FIND/WHO (diagnostic specificity >98% and sensitivity >95%).<sup>17,18</sup> Diagnostic accuracy was shown for all six major genotypes, although genotypes 2, 5 and 6 had a limited presence in the samples tested (online supplementary table S1 B).

Although DAAs for HCV cure are now available, the challenge for HCV infection, and especially for HCV elimination programmes, is the screening of undiagnosed HCV-infected subjects and the linkage to care, namely the availability of DAA at reasonable costs. The majority of infected individuals are unaware of their infection and as such do not benefit from this unprecedented medical advance. HCV infection is often first identified by serological testing, followed by diagnostic confirmation of HCV nucleic acid in circulation.<sup>20–22</sup> At present, HCV RNA testing is only available in centralised medical settings resulting in less than 1% of infected people in low-income and middle-income countries being aware of their infection.<sup>14</sup> Therefore, there is an urgent unmet clinical need for decentralised HCV testing that can provide rapid diagnosis that may result in improved patient outcome and public health.<sup>23</sup>

Initial assay validation was conducted to establish specificity for HCV with clinical validation performed in a two-site retrospective clinical study to determine sensitivity and specificity. In our treatment-naïve case subjects, the lack of a result and false-negative results were only observed in samples with RNA viral levels below the LoD of the Genedrive HCV assay. It is unlikely that this analytical sensitivity would impact diagnostic sensitivity, as the vast majority of people with pretreatment chronic HCV infection have high VL (>5 log IU/mL). This is supported by a study of 2400 treatment-naïve patients with chronic HCV, where it was shown that <2% of patients had viraemia below 4 log<sub>10</sub> IU/mL.<sup>24</sup> Furthermore, recent WHO guidelines explicitly acknowledged that ‘a limit of detection of 3000 IU/mL or lower would be acceptable and would identify 95% of those with viremic infection’.<sup>10</sup>

Linking diagnosis to treatment decisions is now the major clinical challenge: in Europe, the percentage of cured individuals from initial viraemic infections is only 4.1%.<sup>25</sup> In the best cases, the treatment rate is around 5% (eg, France and Germany) and does not exceed 1% in resource-limited countries.<sup>26,27</sup> Furthermore, it has been well described how HCV prevalence negatively correlates with treatment rate,<sup>28,29</sup> as illustrated by countries such as Egypt and Mongolia, which have very high HCV prevalence and treatment rates between 0.1% and 1.2%. These data show how both initial screening and linkage to care after diagnosis are still major challenges, one of the main consequences being the low number of individuals that are prescribed and adhere to treatment.<sup>30</sup> This is partially due to a lack of PoC diagnostic tests, as well as the requirement for multiple clinical visits (eg, diagnosis by serology followed by NAAT). Devices like Genedrive could have a major impact in improving the chronic HCV care continuum, as diagnosis through HCV RNA detection could be performed in a wide range of field settings and no intermediate visits would be required. Thus, screening in resource-limited settings could be easy, effective and rapid and, consequently, treatment could be started as soon as a positive result has been obtained.

While several NAAT PoC assays are in development, such as Alere q (Alere), EOSCAPE (Wave 80 Biosciences), PanNAT platform (Micronics), Truelab PCR (Molbio Diagnostics) and RT CPA (Ustar Biotechnologies), the Cepheid Xpert HCV Viral Load assay (Cepheid) is currently the only CE-IVD certified assay for decentralised HCV VL determination. The Xpert HCV test has been reported to have good performance but also presents important

limitations.<sup>31–33</sup> The conventional Xpert HCV assay requires a large volume of sample (1 mL plasma) and an electrical power supply, meaning that basic laboratory infrastructure is needed,<sup>34,35</sup> although a recent study carried out by Grebeley and colleagues validated the use of 100 µL of capillary blood.<sup>25</sup> Furthermore, the Cepheid Xpert cartridge contains guanidinium thiocyanate as the lysis reagent, which is highly toxic and therefore needs special precautions when handling and should be disposed of by incineration.<sup>36</sup> In contrast, the Genedrive assay does not contain any highly toxic chemicals and only requires 30 µL of sample providing the potential for blood pin prick testing, although a requirement for serum/plasma remains an important limitation. Genedrive is a 600 g portable device that can be battery operated, thus making it highly suitable for decentralised testing and use in field settings. While the Genedrive HCV assay is not fully automated, it is a user-friendly device and requires less than a day’s training to operate. While the Xpert HCV assay is quantitative, the Genedrive HCV assay has potential semiquantitative properties as demonstrated herein. Nevertheless, the clinical significance of quantifying HCV VL is less important in the era of highly potent DAAs at least for diagnosing active HCV infection. The aim of this study was to develop and assess the diagnostic accuracy of the Genedrive HCV assay for the qualitative detection of HCV RNA. Thus, we compared Genedrive directly to the Abbott RealTime HCV test, a well-established PCR test that offers HCV VL and genotype testing and is considered the current gold standard. However, it would also be of interest to perform a direct comparison with Cepheid Xpert, to test for equivalence or superiority between the two platforms.

Finally, a crucial factor when assessing the feasibility of PoC technologies is their cost, both for platform and assays. However, direct comparisons of costs to other platforms are challenging as there are geographic differences and unknowns such as subsidised or negotiated pricing, distributor margins and import duty. Furthermore, additional costs related to the platform should be considered, including installation, maintenance, training and lifespan of the equipment and tests. While the final Genedrive costs remain to be defined, the instrument is approximately \$5000, which is considerably lower than other molecular systems, and each Genedrive HCV test will be available for \$30–\$40, depending on the country or region in question, which is equivalent to other products on the market.

The results obtained in Western laboratory-based settings are highly encouraging; however, the accuracy of PoC tests can be often reduced when used in real clinical settings.<sup>34,35</sup> To address this point, we performed a validation study in South Africa with samples obtained from various African countries. In this real-life clinical setting, where there is often a large time between sample collection and testing—as illustrated by the high levels of hemolyses (27/130 samples)—the performance of the Genedrive HCV assay was equally specific and sensitive. Although we demonstrated the efficacy of the Genedrive test in real-life clinical conditions, including different temperatures, different time delays from sample withdrawal and processing and different operators, further validation of the whole system for infield use in decentralised settings is necessary to achieve its full potential of supporting wide spread screening campaigns.

Here, we provide proof of concept in a real-life clinical setting that the Genedrive HCV assay has great potential to provide an affordable and robust instrument for decentralised HCV NAAT testing. This highly sensitive and specific test has recently obtained CE-IVD certification and is positioned to enable real-time treatment management of patients with chronic HCV in any clinical setting. Therefore, the next step with the Genedrive



HCV assay requires prospective validation in real-life decentralised settings in low-income and middle-income countries.

#### Author affiliations

- <sup>1</sup>Immunobiology of Dendritic Cells, Institut Pasteur, Paris, France
- <sup>2</sup>InsERM U1223, Institut Pasteur, Paris, France
- <sup>3</sup>Unité d'Epidémiologie des Maladies Emergentes, Institut Pasteur, Paris, France
- <sup>4</sup>Centre for Translational Research, Institut Pasteur, Paris, France
- <sup>5</sup>INSERM UMS20, Institut Pasteur, Paris, France
- <sup>6</sup>genedrive plc, Manchester, UK
- <sup>7</sup>Groupe Hospitalier Cochin Hôtel-Dieu, Université Paris Descartes et Département d'hépatologie, Paris, France
- <sup>8</sup>Service de Virologie, Université Paris Descartes, EA4474 "Hepatitis C Virology", AP-HP, Hôpital Cochin, Paris, France
- <sup>9</sup>PACRI Unit, Conservatoire National des Arts et Métiers, Paris, France
- <sup>10</sup>Hospital Universitario La Paz, Madrid, Spain
- <sup>11</sup>Liver Unit, Hospital Universitario La Paz, CIBERehd, IdiPAZ, Madrid, Spain
- <sup>12</sup>Lancet Laboratories, Johannesburg, South Africa
- <sup>13</sup>Gastrointestinal and Liver Disorders Theme, NIHR Nottingham Biomedical Research Centre, Nottingham University Hospitals NHS Trust, University of Nottingham, Nottingham, UK
- <sup>14</sup>Department of Cancer Immunology, Genentech Inc, San Francisco, California, USA

**Contributors** AL, T-PB, ARR, J-FM, PC, AM, CPM, WLI, AG and RV analysed clinical samples. AL, YS and SA analysed the data. SA, RF, EH, ML and GM developed the assay. YS, EM, AF, SP, GM, MLA and DD designed the study. AF, SP, GM, MLA and DD secured funding. AL, YS and DD wrote the manuscript. All authors approved the final version and revised it critically before submission.

**Funding** The overall study was funded by the European Commission FP7 (n° 601851). Work performed in Nottingham and South Africa received additional financial support from genedrive plc.

**Competing interests** SA, RF, EH, ML and GM are full time employees of genedrive plc. The remaining authors declare that they have nothing to disclose regarding funding or conflict of interest with respect to this manuscript.

**Patient consent** Detail has been removed from this case description/these case descriptions to ensure anonymity. The editors and reviewers have seen the detailed information available and are satisfied that the information backs up the case the authors are making.

**Ethics approval** ANRS.

**Provenance and peer review** Not commissioned; externally peer reviewed.

**Open Access** This is an Open Access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 4.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/4.0/>

© Article author(s) (or their employer(s) unless otherwise stated in the text of the article) 2018. All rights reserved. No commercial use is permitted unless otherwise expressly granted.

#### REFERENCES

- 1 Polaris Observatory HCV Collaborators. Global prevalence and genotype distribution of hepatitis C virus infection in 2015: a modelling study. *Lancet Gastroenterol Hepatol* 2017;2:161–76.
- 2 World Health Organization. Hepatitis C Fact sheet. Updated Oct 2017 <http://www.who.int/mediacentre/factsheets/fs164/en>
- 3 Lawitz E, Sulkowski MS, Ghalib R, et al. Simeprevir plus sofosbuvir, with or without ribavirin, to treat chronic infection with hepatitis C virus genotype 1 in non-responders to pegylated interferon and ribavirin and treatment-naïve patients: the COSMOS randomised study. *Lancet* 2014;384:1756–65.
- 4 Afdhal N, Zeuzem S, Kwo P, et al. Ledipasvir and sofosbuvir for untreated HCV genotype 1 infection. *N Engl J Med* 2014;370:1889–98.
- 5 Kowdley KV, Lawitz E, Poordad F, et al. Phase 2b trial of interferon-free therapy for hepatitis C virus genotype 1. *N Engl J Med* 2014;370:222–32.
- 6 Ferenci P, Bernstein D, Lalezari J, et al. ABT-450/r-ombitasvir and dasabuvir with or without ribavirin for HCV. *N Engl J Med* 2014;370:1983–92.
- 7 Zeuzem S, Dusheiko GM, Salupere R, et al. Sofosbuvir and ribavirin in HCV genotypes 2 and 3. *N Engl J Med* 2014;370:1993–2001.
- 8 World Health Organization. *Global health sector strategy on viral hepatitis 2016–2021. Towards ending viral hepatitis. Section 01: Towards eliminating viral hepatitis.* Geneva: WHO, 2016:10–16.
- 9 Gower E, Estes C, Blach S, et al. Global epidemiology and genotype distribution of the hepatitis C virus infection. *J Hepatol* 2014;61:S45–S57.
- 10 World Health Organization. *Guidelines on Hepatitis B and C testing. Section 5: Background - Diagnostics for Testing for Hepatitis B and C Infection.* Geneva: WHO, 2017:30–4.
- 11 Di Bisceglie AM. Natural history of hepatitis C: its impact on clinical management. *Hepatology* 2000;31:1014–8.
- 12 Micallef JM, Kaldor JM, Dore GJ. Spontaneous viral clearance following acute hepatitis C infection: a systematic review of longitudinal studies. *J Viral Hepat* 2006;13:34–41.
- 13 Shivkumar S, Peeling R, Jafari Y, et al. Accuracy of rapid and point-of-care screening tests for hepatitis C. *Ann Intern Med* 2012;157:558–66.
- 14 Foundation for Innovative New Diagnostics. *Strategy for Hepatitis C 2015–2020. Section 5: Global Diagnostic Needs in Hepatitis C: The Current Landscape.* Geneva Foundation for Innovative New Diagnostics, 2014:6–12.
- 15 Genedrive. *Performance evaluation study: use of genedrive HCV ID kit for the detection of hepatitis C virus in plasma.* Manchester: Genedrive plc, 2016.
- 16 Pasloske BL, Walkerpeach CR, Obermoeller RD, et al. Armored RNA technology for production of ribonuclease-resistant viral RNA controls and standards. *J Clin Microbiol* 1998;36:3590–4.
- 17 Foundation for Innovative New Diagnostics. *High-priority target product profile for hepatitis C diagnosis in decentralized settings: report of a consensus meeting. Section 7: revised target product profile for a test for diagnosis of active HCV infection and test of cure.* Geneva: Foundation for Innovative New Diagnostics, 2015:11–20.
- 18 Ivanova Reipold E, Easterbrook P, Trianni A, et al. Optimising diagnosis of viraemic hepatitis C infection: the development of a target product profile. *BMC Infect Dis* 2017;17:707.
- 19 Bossuyt PM, Reitsma JB, Bruns DE, et al. STARD 2015: an updated list of essential items for reporting diagnostic accuracy studies. *Radiology* 2015;277:826–32.
- 20 World Health Organization. *Guidelines for the screening, care and treatment of persons with chronic hepatitis c infection. Section 5: Recommendations on Screening.* Geneva: World Health Organization, 2016:43–9.
- 21 European Association for the Study of the Liver. *EASL recommendations on treatment of hepatitis C 2016.* *J Hepatol* 2017;66.
- 22 AASLD & IDSA. *Recommendations for testing, managing, and treating hepatitis C. Section 6: HCV Testing and linkage to care.* American association for the study of liver diseases and infectious diseases society of America. [www.hcvguidelines.org/full-report-view](http://www.hcvguidelines.org/full-report-view) (accessed 22 May 2017).
- 23 Schito M, Peter TF, Cavanaugh S, et al. Opportunities and challenges for cost-efficient implementation of new point-of-care diagnostics for HIV and tuberculosis. *J Infect Dis* 2012;205(Suppl 2):S169–S180.
- 24 Ticehurst JR, Hamzeh FM, Thomas DL. Factors affecting serum concentrations of hepatitis C virus (HCV) RNA in HCV genotype 1-infected patients with chronic hepatitis. *J Clin Microbiol* 2007;45:2426–33.
- 25 European Union HCV Collaborators. Hepatitis C virus prevalence and level of intervention required to achieve the WHO targets for elimination in the European Union by 2030: a modelling study. *Lancet Gastroenterol Hepatol* 2017;2:325–36.
- 26 Razavi H, Waked I, Sarrazin C, et al. The present and future disease burden of hepatitis C virus (HCV) infection with today's treatment paradigm. *J Viral Hepat* 2014;21(Suppl 1):34–59.
- 27 Hatzakis A, Chulanov V, Gadano AC, et al. The present and future disease burden of hepatitis C virus (HCV) infections with today's treatment paradigm - volume 2. *J Viral Hepat* 2015;22(Suppl 1):26–45.
- 28 Mark E, Dore G, Dore GJ, et al. Hepatitis C disease burden and strategies to manage the burden. *J. Viral Hepat* 2014;21:1–4.
- 29 Wedemeyer H, Dore GJ, Ward JW. Estimates on HCV disease burden worldwide - filling the gaps. *J Viral Hepat* 2015;22(Suppl 1):1–5.
- 30 Zhou K, Fitzpatrick T, Walsh N, et al. Interventions to optimise the care continuum for chronic viral hepatitis: a systematic review and meta-analyses. *Lancet Infect Dis* 2016;16:1409–22.
- 31 Gupta E, Agarwala P, Kumar G, et al. Point-of-care testing (POCT) in molecular diagnostics: Performance evaluation of GeneXpert HCV RNA test in diagnosing and monitoring of HCV infection. *J Clin Virol* 2017;88:46–51.
- 32 McHugh MP, Wu AHB, Chevaliez S, et al. Multicenter Evaluation of the Cepheid Xpert Hepatitis C Virus Viral Load Assay. *J Clin Microbiol* 2017;55:1550–6.
- 33 Grebely J, Lamoury FMJ, Hajarizadeh B, et al. Evaluation of the Xpert HCV Viral Load point-of-care assay from venepuncture-collected and finger-stick capillary whole-blood samples: a cohort study. *Lancet Gastroenterol Hepatol* 2017;2:514–20.
- 34 Pai NP, Vadnais C, Denking C, et al. Point-of-care testing for infectious diseases: diversity, complexity, and barriers in low- and middle-income countries. *PLoS Med* 2012;9:e1001306.
- 35 Hänscheid T, Rebelo M, Grobusch MP. Point-of-care tests: where is the point? *Lancet Infect Dis* 2014;14:922.
- 36 World Health Organization. *WHO Prequalification of In Vitro Diagnostics PUBLIC REPORT.* 2017:1–34. Product: Xpert® HCV Viral Load with GeneXpert® Dx, GeneXpert® Infinity-48s, and GeneXpert® Infinity-80 WHO reference number: PQDx 0260-070-00.