


Risk stratification of symptomatic patients suspected of colorectal cancer using faecal and urinary markers

M. M. Widlak*[†], M. Neal[‡], E. Daulton[§], C. L. Thomas[¶], C. Tomkins[¶], B. Singh^{**},
C. Harmston^{††}, A. Wicaksono[§], C. Evans^{††}, S. Smith^{¶‡‡}, R. S. Savage[‡], J. A. Covington[§] and
R. P. Arasaradnam*^{†§§} 

*Department of Gastroenterology, University Hospitals Coventry and Warwickshire NHS Trust, Coventry, UK, [†]Medical School, University of Warwick, Coventry, UK, [‡]Department of Statistics, University of Warwick, Coventry, UK, [§]School of Engineering, University of Warwick, Coventry, UK, [¶]Department of Biochemistry, University Hospitals Coventry and Warwickshire NHS Trust, Coventry, UK, ^{**}Department of Colorectal Surgery, Leicester General Hospital, Leicester, UK, ^{††}Department of Colorectal Surgery, University Hospitals Coventry and Warwickshire NHS Trust, Coventry, UK, ^{‡‡}Midlands and North West Bowel Cancer Screening Hub, University Hospitals Coventry and Warwickshire, Coventry, UK, and ^{§§}Applied Biological and Experimental Sciences, University of Coventry, Coventry, UK

Received 20 June 2018; accepted 17 September 2018; Accepted Article online 24 September 2018

Abstract

Aim Faecal markers, such as the faecal immunochemical test for haemoglobin (FIT) and faecal calprotectin (FCP), have been increasingly used to exclude colorectal cancer (CRC) and colonic inflammation. However, in those with lower gastrointestinal symptoms there are considerable numbers who have cancer but have a negative FIT test (i.e. false negative), which has impeded its use in clinical practice. We undertook a study of diagnostic accuracy CRC using FIT, FCP and urinary volatile organic compounds (VOCs) in patients with lower gastrointestinal symptoms.

Method One thousand and sixteen symptomatic patients with suspected CRC referred by family physicians were recruited prospectively in accordance with national referring protocol. A total of 562 patients who completed colonic investigations, in addition to providing stool for FIT and FCP as well as urine samples for urinary VOC measurements, were included in the final outcome measures.

Results The sensitivity and specificity for CRC using FIT was 0.80 [95% confidence interval (CI) 0.66–0.93]

and 0.93 (CI 0.91–0.95), respectively. For urinary VOCs, the sensitivity and specificity for CRC was 0.63 (CI 0.46–0.79) and 0.63 (CI 0.59–0.67), respectively. However, for those who were FIT-negative CRC (i.e. false negatives), the addition of urinary VOCs resulted in a sensitivity of 0.97 (CI 0.90–1.0) and specificity of 0.72 (CI 0.68–0.76).

Conclusions When applied to the FIT-negative group, urinary VOCs improve CRC detection (sensitivity rises from 0.80 to 0.97), thus showing promise as a second-stage test to complement FIT in the detection of CRC.

Keywords Faecal biomarker, faecal immunochemical test for haemoglobin, urinary volatile organic compounds, colorectal cancer

What does this paper add to the literature?

This novel study has applied both faecal and urinary markers to improve the diagnostic accuracy for CRC from 80% to 97%, with a negative predictive value of 100% in those with lower gastrointestinal symptoms. Application of low-cost (< £50) diagnostic tests could result in changes to the current clinical pathway.

Introduction

Colorectal cancer accounts for 12% of all cancers in the UK with as strong age relationship. The clinical presentation of CRC is varied; hence the dilemma for the clinician is to distinguish those with significant versus

nonsignificant pathology without recourse to invasive and costly investigations [1].

Symptoms alone are not sufficiently sensitive to diagnose CRC, and up to a third of patients who undergo invasive investigations have normal outcomes (i.e. colonoscopy reported as macroscopically and microscopically normal) [2]. For those with lower gastrointestinal symptoms there is increasing evidence for the use of the faecal immunochemical test for haemoglobin (FIT) as

Correspondence to: Professor Ramesh P. Arasaradnam, Department of Gastroenterology, University Hospitals Coventry and Warwickshire NHS Trust, Clifford Bridge Road, Coventry CV2 2DX, UK.
E-mail: r.arasaradnam@warwick.ac.uk

first-line testing; this has a high negative predictive value (NPV) of 0.99, but sensitivity is relatively low (0.80–0.90) [2–4].

In symptomatic patients it remains uncertain as to how to interpret a FIT-negative result – in other words, can such patients be sufficiently reassured that they do not have cancer? Alternative or additional tests are therefore required. For example, a stool marker for inflammation – faecal calprotectin (FCP) – or urinary volatile metabolic markers could be utilized to risk stratify those with suspected CRC. Metabolic markers such as urinary, faecal or breath volatile organic compounds (VOCs) have previously been shown to aid in CRC detection [5–8]. To our knowledge, this is the first study to use a combination of faecal tests (FIT and FCP) as well as urinary VOC testing in a symptomatic population suspected of CRC.

Method

Design and setting

This was a single-centre, prospective, blinded study of patients with lower gastrointestinal symptoms referred by family physicians to tertiary care with suspected CRC. Ethical approval was granted by Coventry and Warwickshire Research Ethics Committee, UK as part of the FAMISHED (Food and Fermentation using Metagenomics in Health and Disease) multicentre study (09/H1211/38). The study protocol conforms to ethical guidelines of the 1975 Declaration of Helsinki as reflected by the institution's human research committee.

Participants

A total of 1850 patients were approached with criteria for inclusion based on national referral criteria (Appendix S1 in the online Supporting Information). Of these, 834 were excluded for a combination of reasons including physical frailty, illness, language barriers, etc. One thousand and sixteen patients provided consent and underwent colonic investigations (endoscopic or radiological colonic cross-sectional imaging). Sixteen patients withdrew from the study and 310 failed to provide stool samples (69% return rate). A further 78 patients were excluded as only one stool sample was provided (both the FIT sample collection device and stool for FCP were required). Spot urine samples were received from 762 patients (76% return rate), but 39 were excluded due to insufficient sample volume or failed VOC urine analysis. This left a group of 562 patients with matching urine and stool samples (FIT and FCP) who were included for final statistical analysis

(Fig. 1). Those who were under the age of 18, pregnant, did not meet the referral criteria for urgent review for lower gastrointestinal symptoms or had incomplete colonic examinations were excluded from the study.

All study participants were given a pack containing a FIT sample collection device (Extel Hemo-auto MC A device; Kyowa Medex, Tokyo, Japan via Alpha Laboratories Ltd., Eastleigh, UK), which holds 2 mg of faeces in 2 ml of buffer, and a Universal Sterilin 30 ml stool pot for the FCP sample. Written and pictorial instructions for collections were provided, with the Fe-Col® sample collection aid (Alpha Laboratories Ltd.). Time of collection and time of receipt of samples at the laboratory were recorded and all samples were stored refrigerated at 2–8°C until analysis. Patients were asked to return the sample prior to colonic investigations. Samples returned more than 4 days after collection were excluded from analysis.

All study participants were also provided with a Universal Sterilin 30 ml pot to collect a spot urine sample at the time of the clinic visit. Timing of collection was recorded and urine samples were stored with sealed caps in a freezer at –80°C. For patients who were unable to provide the sample in the clinic, a urine container was provided as well as a return envelope, and the sample was sent via courier from their primary physician and stored at –80°C.

Intervention

Quantitative FIT was performed on automated HM-JACKarc analyser (Kyowa Medex) by the Midlands and North West Bowel Cancer Screening Hub, Rugby, UK on a weekly basis. Stool samples for FCP analysis were extracted manually by trained laboratory staff alongside the routine calprotectin service. Extracted calprotectin was measured using the EliA Calprotectin fluoroenzyme immunoassay on the automated Thermo Fisher ImmunoCap 250 analyser (Thermo Fisher Scientific, Waltham, Massachusetts, USA).

In the analyses for CRC, thresholds were determined from the data (*a priori*) to maximize sensitivity under the constraint that the NPV was ≥ 0.99 . For FCP, no threshold achieved an NPV ≥ 0.99 . The lowest detection limit of this assay for FIT is 3 µg/g faeces.

Samples were initially stored at –80°C. Prior to urinary VOC analysis, they underwent a graded defrost process (based on our unpublished established protocol). Ten-millilitre glass vial aliquots (Thermo-Fisher Scientific, Loughborough, UK) suitable for use with an autosampler (MPS, Gerstel GmbH, Mülheim an der Ruhr, Germany) were used. Crimp caps (Chromacol Ltd., Merck, UK) with silicone polytetrafluoroethylene

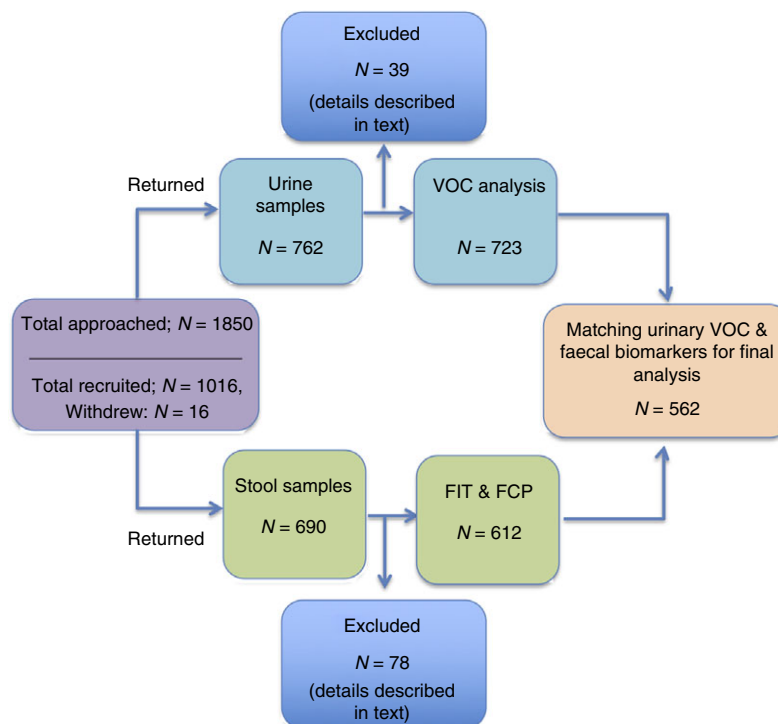


Figure 1 Study flow diagram of the total number of patients recruited via the urgent colorectal lower gastrointestinal pathway who met the inclusion criteria.

septa were used to seal each sample. Each septum and crimp cap was baked for 6 h at 200°C prior to use to remove any potential interfering molecules. Control blanks of air were prepared using the same method. A commercial gas analysis instrument [Lonestar (FAIMS), Owlstone, Cambridge, UK], based on ion mobility spectroscopy (IMS), was used to analyse VOCs emanating from the urine samples. Details of the Lonestar and its application in medical diagnostics have been reviewed by Covington *et al.* [9]. The setup was bespoke for this application, to detect unique VOC chemical ‘fingerprints’. Our previous study has determined the optimal sample capture and storage methods to minimize diurnal and day to day variation [10] (see Appendix S1 for further methodological details on the use of FAIMS).

Outcomes

Participating clinicians, endoscopists and radiologists were blinded to the results of the urinary VOCs and faecal tests (FIT/FCP). Diagnosis of CRC and adenomas was confirmed histologically. High-risk adenomas were defined as lesions with high-grade dysplasia and/or serrated, villous histology, ≥ 10 mm in size or the presence of three or more adenomas. Hyperplastic polyps were excluded. Faecal and urinary VOC results were compared with the outcome of the colonic investigations and divided into clinical groups: CRC, high-risk adenoma, all adenomas, others and normal.

All analysis was carried out using R 3.4.2 (R Foundation for Statistical Computing, Vienna, Austria). To combine FIT and FCP measurements, 10-fold cross validation and a Bayesian robust logistic regression model were used to generate predictive scores for each patient using the above criteria.

The VOC data variables were assessed for predictive value with a Mann–Whitney *U*-test at a significance level of 0.8, with multi-test correction performed using sequential goodness of fit (SGoF) [11]. This yielded a set of about 10^4 candidate variables. The dimensionality of the selected variables was then reduced using principal component analysis (PCA) [12], with the number of principal components chosen using cross-validation (typically around 35) [13]. For testing using VOCs alone, a logistic regression model was trained on the results. For the two-phase test incorporating FIT as a first phase a support vector machine, radial kernel, was used on the top 128 variables, without using SGoF for variable selection. Ten-fold cross-validation was used to generate predictive scores for each patient, and a threshold was selected to maximize sensitivity.

Results

The main patient characteristics are shown in Table 1. A total of 562 patients, who completed colonic investigations, in addition to providing stool for FIT and FCP

Table 1 Patient demographic to include clinical characteristics and main diagnostic outcomes.

Characteristic	<i>n</i> (%)	Median FIT ($\mu\text{g/g}$ faeces)	Median FCP ($\mu\text{g/g}$ faeces)
Total	562		
Male gender	286 (51)		
Median age (years) (range)	68 (29–89)		
Presenting symptoms*			
Altered bowel habit	369 (66)		
Weight loss	87 (15)		
Rectal bleeding	232 (41)		
Anaemia	121 (22)		
Iron-deficiency anaemia	91 (16)		
Abdominal pain	164 (29)		
CRC	35 (6.2)	270	86
High-risk adenoma	27 (4.8)	14	23
All adenomas (low and high risk)	94 (17)	3.5	19
Others (DD, IBD, MC etc.)	173 (31)	1.8	21
Normal	233 (41)	1.4	16

Normal – no colonic pathology identified.

DD, diverticular disease; IBD, inflammatory bowel disease; MC, microscopic colitis.

*Note that some patients may present with a constellation of symptoms.

as well as urine samples for urinary VOC measurements, were included in the final statistical analysis.

Faecal biomarkers (FIT and FCP)

The use of FIT alone for CRC detection revealed a sensitivity of 0.80 [95% confidence interval (CI) 0.66–0.93] and specificity of 0.93 (95% CI 0.91–0.95). The area under the curve (AUC) was 0.90 (95% CI 0.85–0.96). The NPV was 0.99 (95% CI 0.98–1.0).

For high-risk adenomas, applying threshold of 3 $\mu\text{g/g}$ faeces which is the lowest detection limit of the assay, the sensitivity was 0.70 (CI 0.52–0.87), specificity of 0.66 (CI 0.623–0.703) with a NPV of 0.98

(CI 0.96–0.99) (see Table 2 for the diagnostic performance of FIT).

FCP alone performed less well than FIT in the CRC and adenoma groups. Test performances combining FIT and FCP are shown in Table 2. The respective receiver operator curve characteristics (ROCs) for CRC, high-risk adenomas and all adenomas using FIT are shown in Fig. 2.

Urinary VOCs

For urinary VOC analysis, logistic regression using SGoF selection and PCA was applied to form a unique ‘chemical fingerprint’. The sensitivity of urinary VOCs

Table 2 Diagnostic performance of FIT alone and in combination with FCP for CRC, high-risk adenoma and all adenomas.

Disease group	Sensitivity (CI)	Specificity (CI)	NPV (CI)	PPV (CI)	AUC (CI)
FIT					
CRC (<i>n</i> = 35)	0.80 (0.66–0.93)	0.93 (0.91–0.95)	0.99 (0.98–1.0)	0.44 (0.32–0.56)	0.90 (0.85–0.96)
High-risk adenoma (<i>n</i> = 27)	0.63 (0.44–0.81)	0.76 (0.72–0.79)	0.98 (0.96–0.99)	0.11 (0.07–0.17)	0.73 (0.63–0.83)
All adenomas (<i>n</i> = 94)	0.53 (0.43–0.63)	0.68 (0.63–0.72)	0.88 (0.84–0.91)	0.25 (0.19–0.31)	0.63 (0.58–0.69)
FIT + FCP					
CRC (<i>n</i> = 35)	0.80 (0.66–0.93)	0.93 (0.91–0.95)	0.99 (0.97–1.0)	0.43 (0.31–0.55)	0.91 (0.86–0.96)
High-risk adenoma (<i>n</i> = 27)	0.93 (0.81–1.0)	0.25 (0.21–0.29)	0.99 (0.96–1.0)	0.06 (0.04–0.08)	0.69 (0.59–0.79)
All adenomas (<i>n</i> = 94)	0.86 (0.79–0.93)	0.26 (0.22–0.30)	0.90 (0.85–0.95)	0.19 (0.15–0.23)	0.60 (0.54–0.67)

High-risk adenoma: adenoma with high-grade dysplasia, villous histology, ≥ 10 mm or ≥ 3 adenomas.

Thresholds applied to achieve the highest sensitivity under the constraint of keeping the negative predictive value > 0.99 .

Note that the lowest limit of detection for the FIT assay is 3 $\mu\text{g/g}$ feces.

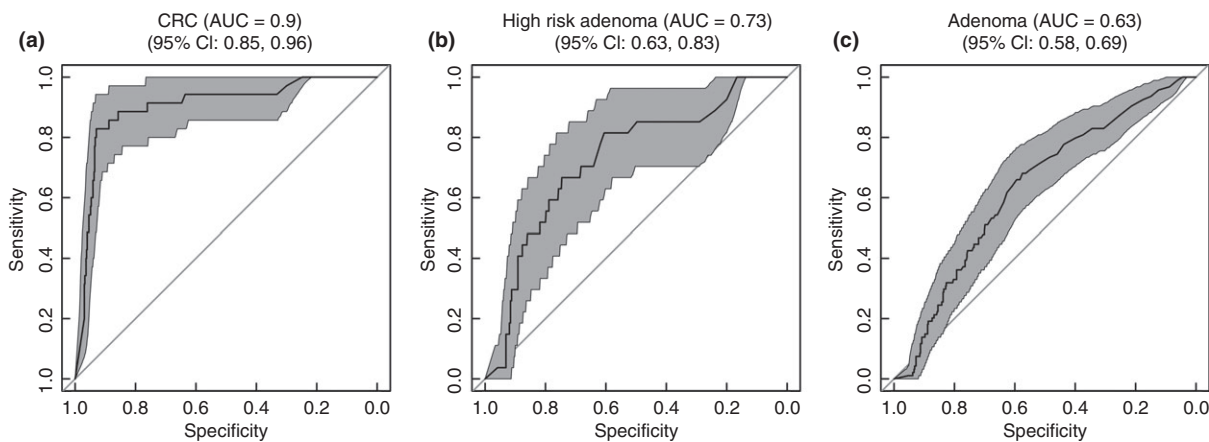


Figure 2 Respective ROC curves for (a) CRC, (b) high-risk adenoma and (c) all adenomas using FIT.

for the detection of CRC was 0.63 (95% CI 0.46–0.79) and specificity was 0.63 (95% CI 0.59–0.67). The AUC was 0.67 (95% CI 0.57–0.77) and the NPV was 0.96 (95% CI 0.94–0.98).

For high-risk adenoma and all adenomas, using urinary VOCs the sensitivity was 0.93 (95% CI 0.81–1.0) and 0.91 (95% CI 0.85–0.97), respectively, with specificity of 0.16 (95% CI 0.13–0.20) and 0.15 (95% CI 0.12–0.19), respectively (Table 3).

Combining FIT and urinary VOCs

When used in combination (FIT and urinary VOCs) for CRC, the sensitivity was 0.80 (95% CI 0.66–0.93) and specificity was 0.89 (95% CI 0.87–0.93). The AUC was 0.86 (95% CI: 0.77–0.94) and the NPV 0.99 (95% CI: 0.97–1.0).

Further analysis of urinary VOCs in a setting of FIT-negative CRC (missed cancer) improved the sensitivity to 0.97 (95% CI 0.90–1.0) and specificity to 0.72 (95% CI 0.68–0.76) with a NPV of 1.0 (95% CI 0.99–1.0). Figure 3 shows a box plot applying the two-stage filter process to depict improvement in CRC detection when urinary VOCs are used in FIT-negative CRC patients compared with FIT alone. The decision threshold line in each case is the value (of either FIT or predicted

probability) that divides prediction of cancer. Thus, for any patient above the line the test indicates the likelihood of having cancer, and vice versa. Overall, patients are classed as negative for the two-stage test if they are negative for both FIT and VOC screening, and positive if they are positive for either test – FIT (above the threshold) or VOC screening. We did not observe any differences even after stratifying for age or gender. For all adenomas, urinary VOC did not improve detection in those with false-negative FIT.

Discussion

Whilst studies evaluating various faecal markers (methylated genes, microRNA and protein markers) have shown promise for the detection of CRC and adenomas their application within a clinical setting has been limited due to high cost and poor sensitivity, especially when applied in areas with a low disease prevalence [14]. Thus, the emphasis has been on low-cost, non-invasive testing such as FIT for detection of CRC and adenomas. Experience from using FIT in the screening population reveals that it has a relatively high specificity at the expense of sensitivity. Various FIT devices have been trialled and it has become evident that there is considerable heterogeneity in FIT devices for the

Table 3 Diagnostic performance of urinary VOCs for CRC, high-risk adenoma and all adenomas.

Disease group	Sensitivity (CI)	Specificity (CI)	NPV (CI)	PPV (CI)	AUC (CI)
CRC ($n = 35$)	0.63 (0.46–0.79)	0.63 (0.59–0.67)	0.96 (0.94–0.98)	0.10 (0.06–0.14)	0.67 (0.57–0.77)
High-risk adenoma ($n = 27$)	0.93 (0.81–1.0)	0.16 (0.13–0.20)	0.98 (0.94–1.0)	0.05 (0.03–0.07)	0.56 (0.45–0.68)
All adenomas ($n = 94$)	0.91 (0.85–0.97)	0.15 (0.12–0.19)	0.90 (0.83–0.96)	0.18 (0.15–0.21)	0.55 (0.49–0.61)

High-risk adenoma (HRA) - adenoma with high-grade dysplasia, villous histology, ≥ 10 mm or ≥ 3 adenomas.

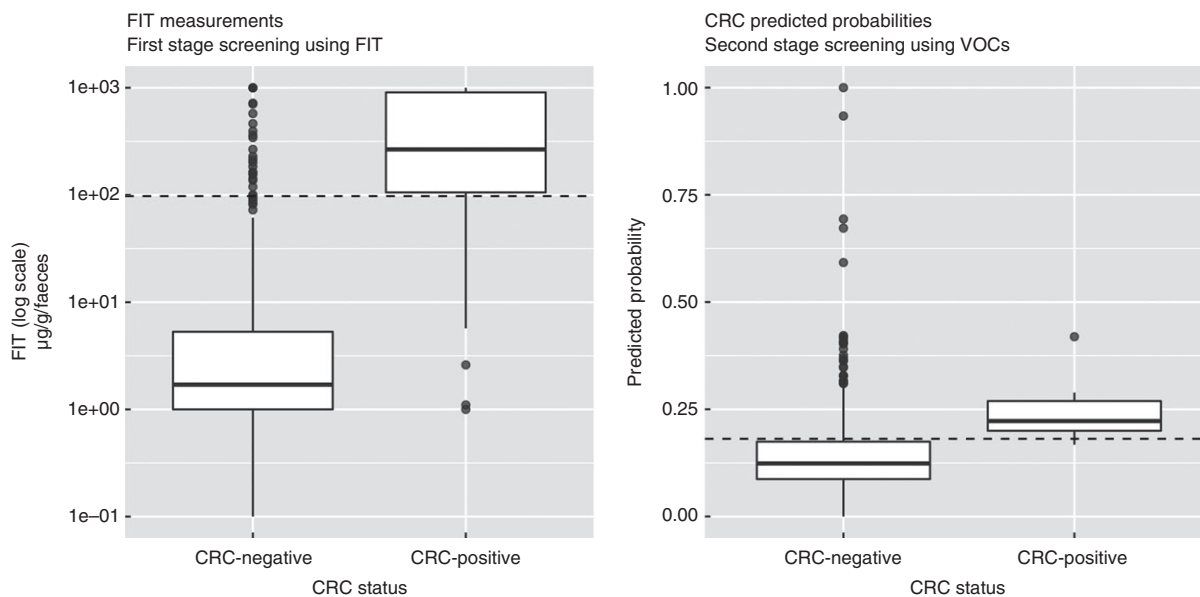


Figure 3 Two plots showing the stages of the FIT/VOC two-stage test. The left-hand plot shows the distribution of FIT measurements for CRC-positive and CRC-negative patients. The dotted line shows the threshold for the first stage test. Patients above this threshold are classed as positive for the test while patients below this threshold go on to the second stage of VOC screening, shown in the right-hand plot. For patients who are screened using VOCs, the distribution of predicted probabilities of CRC is shown in the right-hand plot, with the decision threshold again shown by a dotted line. Patients above the line are positive for the test, while patients below the line are negative.

detection of CRC. This is further compounded by the fact that various threshold levels are applied and there is uncertainty surrounding those who test negative with FIT [15].

In the UK, CRC detection based on symptoms alone is low; ranging from 4% to 8% [16,17]. This suggests that over 90% of patients undergo negative tests for exclusion of CRC. Thus, it is imperative that alternative noninvasive prescreening markers such as FIT and other faecal or urinary markers are used in those with gastrointestinal symptoms so as to minimize unnecessary investigations.

In this study we have shown for the first time the value of dual-modality testing using noninvasive markers (FIT and urinary VOCs) as a two-stage process to exclude CRC. In those with lower gastrointestinal symptoms suspicious of CRC, the use of FIT alone revealed a sensitivity of 80% [meta-analysis suggests 90% (CI 87–92%); unpublished]. Thus, there is the potential to miss one or two cancers out of every ten, which is not sufficiently robust for everyday clinical use. Urinary VOC on its own was less sensitive (0.63; CI 0.46–0.79) and in combination with FIT did not show any improvement in sensitivity (0.80; CI 0.60–0.93). However, the use of a two-stage test, namely the addition of urinary VOC testing in those who test negative for FIT (i.e. false-negative CRC), increases the combined

sensitivity to 97%; this is more acceptable for clinical use and comparable to the performance of colonoscopy at a fraction of the cost.

Urinary VOC analysis on its own (VOC ‘positive’) only provides 63% sensitivity and specificity for CRC detection, and performance does not improve when it is used in combination with FIT or FCP. This may be due to the heterogeneity of the chemical fingerprint that is produced by CRC. However, following preselection by FIT it performs well, as there is reduced background ‘volatile noise’ making it more specific to detect either haemoglobin moiety/breakdown products or glycation end-products. Whilst urinary VOCs demonstrate a high sensitivity for adenoma detection, the lack of specificity and high false-positive rate suggests that this marker may not perform as well for adenoma detection.

The VOCs that are detected using our preanalytical method and FAIMS are unique, as detection is based on volatiles existing in the gaseous phase rather than in the liquid phase. The disease separation is characterized by the mobility of individual ions (i.e. physical rather than chemical properties of the ion), which have low molecular weights (20–200 kDa). Unlike conventional gas chromatography and mass spectroscopy (GC-MS), the specific chemicals are not identified but a ‘chemical fingerprint’ is formed. GC-MS is limited by its high running and labour cost as well as run time. This

impedes its use within routine clinical practice where rapid, low cost and simple operation (nonskilled operators) would be preferable. The composition of key volatile compounds gives rise to the unique chemical fingerprint identified in this study, allowing classification into 'VOC positive or negative' outcomes.

VOCs reflect metabolic cellular changes within the host; for example, detection of advanced glycation end-products which have been implicated in colon carcinogenesis [18]. Our previous work [19] [undertaken using a Bruker Scion GCMS, fitted with dynamic head space sampling and solid-phase micro-extraction (SPME) pre-concentration system] has identified three chemicals which are modulated in CRC – the 'VOC positive' signature. In particular we noted a high incidence of 1,3,5,7-cyclooctatetraene and a low incidence of 1,3-propanediamine and 4-methylbenzoic acid (dietary metabolite). Methylbenzene has also recently been reported to provide a unique chemical signature in those with CRC when a different technology (non-GCMS) is used [20]. A higher incidence of acetone was noted in those with colorectal adenomas; this has been shown by members of our group (unpublished) to be produced by *C. difficile* and in other Clostridiales. Allyl isothiocyanate was also detected in those with CRC but not at elevated levels; the latter is produced by certain *E. coli* strains as dietary substrates with can affect the integrity of the gut mucosa [21].

The use of dual-modality testing, initially with FIT followed by urinary VOCs, enables 97% sensitivity with 100% NPV if both tests are negative for the detection of CRC and high-grade adenomas. Furthermore, findings from this study suggest that the combination of FIT and VOC offers the option for personalized strategies for CRC detection in those with symptoms and avoids the need for repeat FIT testing (if FIT is negative the test probability is unlikely to improve unless there are preanalytical errors).

It is envisaged that both these noninvasive tests (FIT and urinary VOCs) can be undertaken within primary care and analysed within a central laboratory (as FIT currently is) at low cost to guide secondary care referral patterns. The FAIMS unit is commercially available and urine VOCs are deemed stable up to 12 months when stored frozen [10]. It has been purported within a simulation model that for an equivalent biomarker to compare with FIT (£18/test) it should not exceed seven-fold the unit cost of FIT [22] – this is fulfilled in urinary VOCs (£28/test), which cost less than twice the unit cost of FIT. A proposed clinical algorithm is outlined, highlighting the use of a dual noninvasive diagnostic approach in those with lower gastrointestinal symptoms (Figure S1).

Acknowledgements

Patients attending University Hospitals Coventry and Warwickshire NHS Trust for their altruistic contributions. Matthew Davies at Alpha Laboratories Ltd for provision of the reagents and FIT sample collection devices. The NHS Bowel Cancer Screening Hub at Rugby for allowing us to use their HM-JACKarc and Angela Ryder for analysing the FIT sample collection devices. Outpatient and surgical departments (Mr P Baragwanath, H Kashi, K Leong, A Majeed, P McCullough, C Wallis, N Williams and L Wong) for helping to recruit patients. Research nurses Subiatu Wurie, Nicola O'Connell, Bei Bei Liu and research technician Leighanne Burns. Tissue bank and biochemistry staff at the University Hospitals Coventry and Warwickshire NHS Trust. All authors approved the final version of the article, including the authorship list

Author contributions

RPA initiated the project, led the grant funding application, ethics approval and manuscript writing. MMW recruited patients, performed sample handling, produced documentation and contributed to manuscript writing. ED, AW and JAC developed and performed the VOC analysis and contributed to manuscript writing. MN and RSS performed statistical analysis and contributed to manuscript writing. RPA, JAC, CH, CE, SS and BS reviewed the manuscript for intellectual content. CT and CLT were involved in the study design, performed biochemical analysis and contributed to manuscript writing. SS supervised the biochemical analysis. All authors reviewed the final manuscript.

Funding

Bowel Disease Research Foundation, UK, Alphas Labs Ltd.

Disclosure and conflict of interest statement

RPA has provided educational lectures on behalf of Alpha Labs Ltd & Thermo Fisher Scientific Ltd. CT has provided educational lectures on behalf of Thermo Fisher Scientific Ltd. All remaining authors disclose no conflicts of interest.

References

- 1 Arasradnam RP, Brown S, Forbes A *et al.* Guidelines for the investigation of chronic diarrhoea in adults: British Society of Gastroenterology, 3rd edition. *Gut* 2018; **67**: 1380–99

- 2 Widlak MM, Thomas CL, Thomas MG *et al.* Diagnostic accuracy of faecal biomarkers in detecting colorectal cancer and adenoma in symptomatic patients. *Aliment Pharmacol Ther* 2017; **45**: 354–63.
- 3 Mowat C, Digby J, Strachan JA *et al.* Faecal haemoglobin and faecal calprotectin as indicators of bowel disease in patients presenting to primary care with bowel symptoms. *Gut* 2016; **65**: 1463–9.
- 4 Westwood M, Corro Ramos I, Lang S *et al.* Faecal immunochemical tests to triage patients with lower abdominal symptoms for suspected colorectal cancer referrals in primary care: a systematic review and cost-effectiveness analysis. *Health Technol Assess* 2017; **21**: 1–234.
- 5 Arasaradnam RP, McFarlane MJ, Ryan-Fisher C *et al.* Detection of colorectal cancer (CRC) by urinary volatile organic compound analysis. *PLoS One* 2014; **9**: e108750.
- 6 Westenbrink E, Arasaradnam RP, O'Connell N *et al.* Development and application of a new electronic nose instrument for the detection of colorectal cancer. *Biosens Bioelectron* 2015; **67**: 733–8.
- 7 Altomare DF, Di Lena M, Porcelli F *et al.* Exhaled volatile organic compounds identify patients with colorectal cancer. *Br J Surg* 2013; **100**: 144–50.
- 8 de Meij TG, Larbi IB, van der Schee MP *et al.* Electronic nose can discriminate colorectal carcinoma and advanced adenomas by fecal volatile biomarker analysis: proof of principle study. *Int J Cancer* 2014; **134**: 1132–8.
- 9 Covington JA, van der Schee MP, Edge AS *et al.* The application of FAIMS gas analysis in medical diagnostics. *Analyst* 2015; **140**: 6775–81.
- 10 Esfahani S, Sagar NM, Kyrou I *et al.* Variation in gas and volatile compound emissions from human urine as it ages, measured by an electronic nose. *Biosensors* 2016; **6**: 1–11.
- 11 Carvajal-Rodriguez A, de Una-Alvarez J. Assessing significance in high-throughput experiments by sequential goodness of fit and q-value estimation. *PLoS ONE* 2011; **6**: e24700.
- 12 Minka P. Automatic choice of dimensionality for PCA. *MIT Media Lab Percept Comput Sect Tech Rep* 2000; **514**: 1–16.
- 13 Arasaradnam RP, Wicaksono A, O'Brien H *et al.* Noninvasive diagnosis of pancreatic cancer through detection of volatile organic compounds in urine. *Gastroenterology* 2018; **154**: 485–7. e1.
- 14 Iannone A, Losurdo G, Pricci M *et al.* Stool investigations for colorectal cancer screening: from occult blood test to DNA analysis. *J Gastrointest Cancer* 2016; **47**: 143–51.
- 15 Gies A, Cuk K, Schrotz-King P *et al.* Direct comparison of diagnostic performance of 9 quantitative fecal immunochemical tests for colorectal cancer screening. *Gastroenterology* 2018; **154**: 93–104.
- 16 Pande R, Leung E, McCullough P *et al.* Impact of the United Kingdom national bowel cancer awareness campaign on colorectal services. *Dis Colon Rectum* 2014; **57**: 70–5.
- 17 Vulliamy P, McCluney S, Raouf S *et al.* Trends in urgent referrals for suspected colorectal cancer: an increase in quantity, but not in quality. *Ann R Coll Surg Engl* 2016; **98**: 564–7.
- 18 Turovskaya O, Foell D, Sinha P *et al.* RAGE, carboxylated glycans and S100A8/A9 play essential roles in colitis-associated carcinogenesis. *Carcinogenesis* 2008; **29**: 2035–43.
- 19 Westenbrink E. *Electronic Nose Technology as a Point-of-Care Solution for the Detection of Colorectal Cancer*. School of Engineering, University of Warwick, 2016.
- 20 Bhattacharyya D, Kumar P, Mohanty SK *et al.* Detection of four distinct volatile indicators of colorectal cancer using functionalized titania nanotubular arrays. *Sensors* 2017; **17**: 1795.
- 21 Nyangale EP, Mottram DS, Gibson GR. Gut microbial activity, implications for health and disease: the potential role of metabolite analysis. *J Proteome Res* 2012; **11**: 5573–85.
- 22 Lansdorp-Vogelaar I, Goede SL, Bosch LJW *et al.* Cost-effectiveness of High-performance Biomarker Tests vs Fecal Immunochemical Test for Noninvasive Colorectal Cancer Screening. *Clin Gastroenterol Hepatol* 2018; **16**: 504–12.

Supporting Information

Additional Supporting Information may be found in the online version of this article:

Figure S1. Clinical pathway – lower gastrointestinal symptoms suspicious of colorectal cancer.

Appendix S1. Methods.