

Engineering Conferences International ECI Digital Archives

Integrated Continuous Biomanufacturing II

Proceedings

Fall 11-2-2015

Continuous purification of hepatitis C virus-like particles by multi-column chromatography

Ricardo Silva
iBET, rsilva@ibet.pt

Alex Xenopoulos
EMD Millipore

Jose Mota
FCT-UNL

Cristina Peixoto
iTQB/iBET

Manuel Carrondo
FCT-UNL/iBET

See next page for additional authors

Follow this and additional works at: http://dc.engconfintl.org/biomanufact_ii

 Part of the [Biomedical Engineering and Bioengineering Commons](#)

Recommended Citation

Ricardo Silva, Alex Xenopoulos, Jose Mota, Cristina Peixoto, Manuel Carrondo, and Paula Alves, "Continuous purification of hepatitis C virus-like particles by multi-column chromatography" in "Integrated Continuous Biomanufacturing II", Chetan Goudar, Amgen Inc. Suzanne Farid, University College London Christopher Hwang, Genzyme-Sanofi Karol Lacki, Novo Nordisk Eds, ECI Symposium Series, (2015). http://dc.engconfintl.org/biomanufact_ii/131

This Conference Proceeding is brought to you for free and open access by the Proceedings at ECI Digital Archives. It has been accepted for inclusion in Integrated Continuous Biomanufacturing II by an authorized administrator of ECI Digital Archives. For more information, please contact franco@bepress.com.

Authors

Ricardo Silva, Alex Xenopoulos, Jose Mota, Cristina Peixoto, Manuel Carrondo, and Paula Alves

Introduction

- VLPs present challenges for DSP, thus alternative strategies are desired;
- We evaluated semi-continuous, multicolumn chromatography techniques, that have the potential to improve process economics;
- A simple serial connection of two chromatographic columns directs the effluent of one column to an adjacent one, capturing the breakthrough from the previous column and avoiding product loss;
- This simple setup modification overcomes the limitations of single-column processes by achieving saturation of the first column;
- It also benefits from the countercurrent flow between the mobile and the stationary phases, which optimizes the driving force for mass transfer throughout the overall trajectory of the two phases;

AIM: Evaluate continuous multicolumn chromatography for hepatitis C VLP purification

Materials & Methods

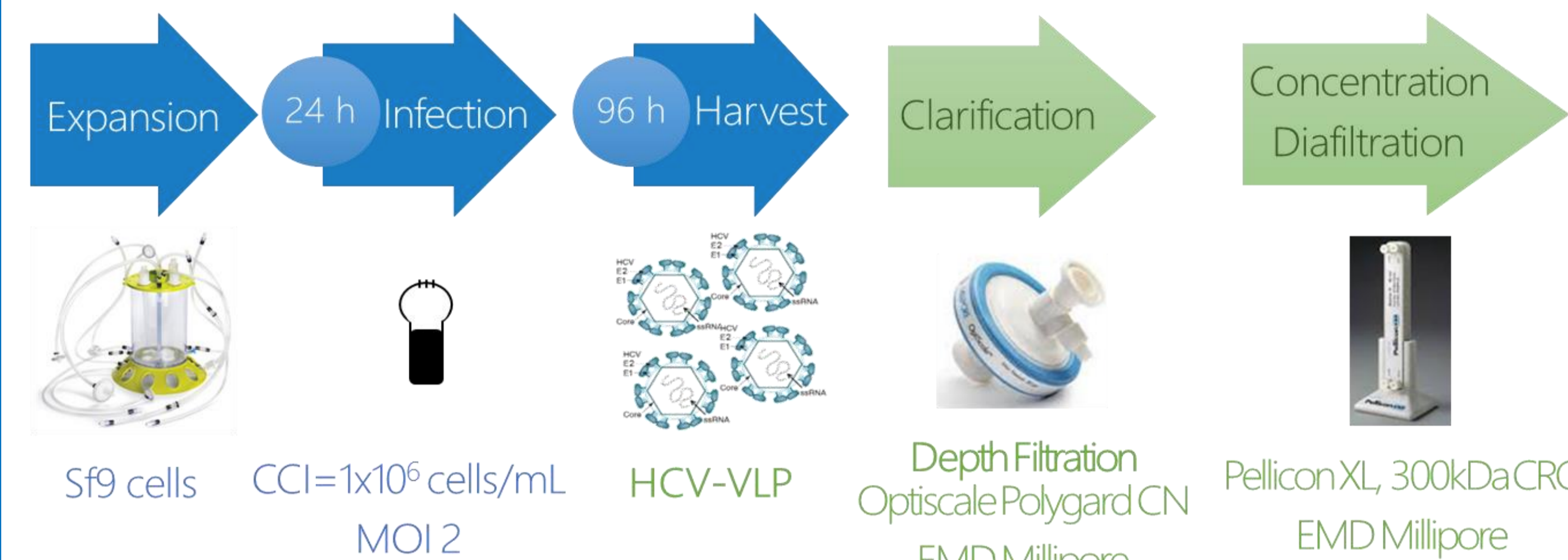


Figure 1: Steps and materials used for preparation of the HCV-VLP bulk

Anion exchange resins: two commercially available matrices (Fractogel® TMAE and DMAE) and three Eshmuno® Q prototypes (1, 2, 3) with different ligand densities (Merck Millipore)

Buffers: 50 mM HEPES at a pH of 7.4 with 0-1 M NaCl

DBC determination for single-column experiments using 1 mL scout columns on ÄKTA Explorer 10s system. This system was modified to support the multicolumn process operation

Columns used in the multicolumn process were slurry packed with Eshmuno® prototype 1 to a final bed volume of 4 mL

HCV-VLP quantification: ELISA kit targeting the GAG p30 protein of the VLPs capsid (QuickTiter™ MuLV, Cell Biolabs, Inc.)

Total Protein quantification: BCA protein assay kit (Pierce Biotechnology)

Total double stranded DNA quantification: PicoGreen™ dsDNA Assay kit (Life Technologies)

Total host cell protein quantification: two-site immunoenzymetric assay (Cygnus Technologies, Inc.)

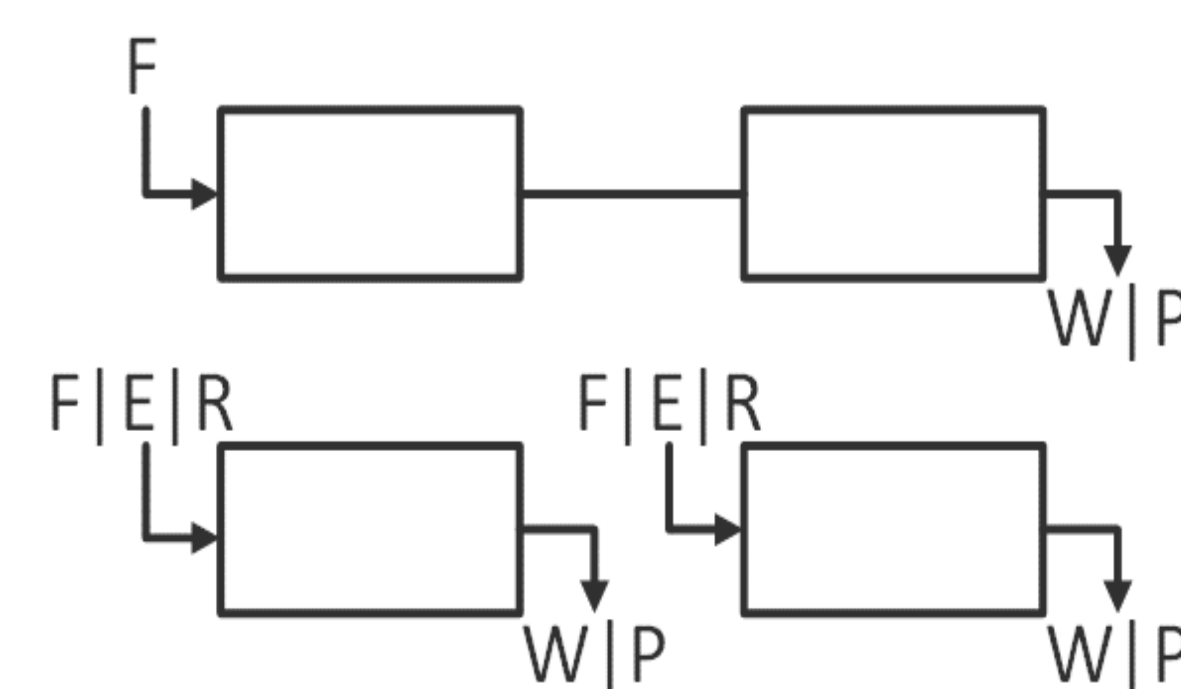
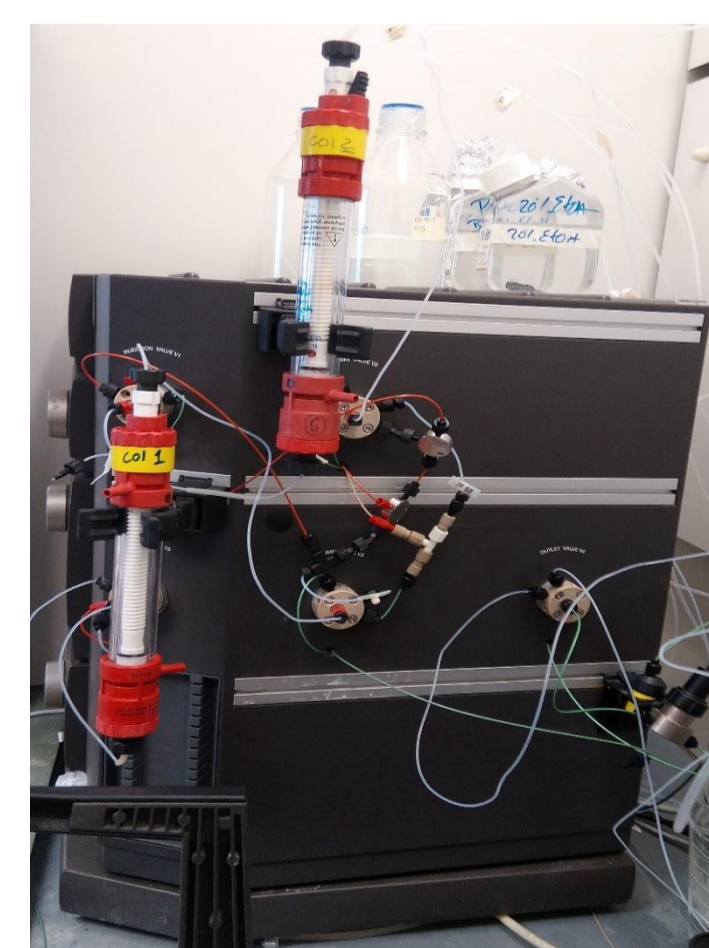


Figure 2: Picture of modified ÄKTA Explorer 10s (left) for implementation of the multicolumn flow schemes (top)

Resin screening

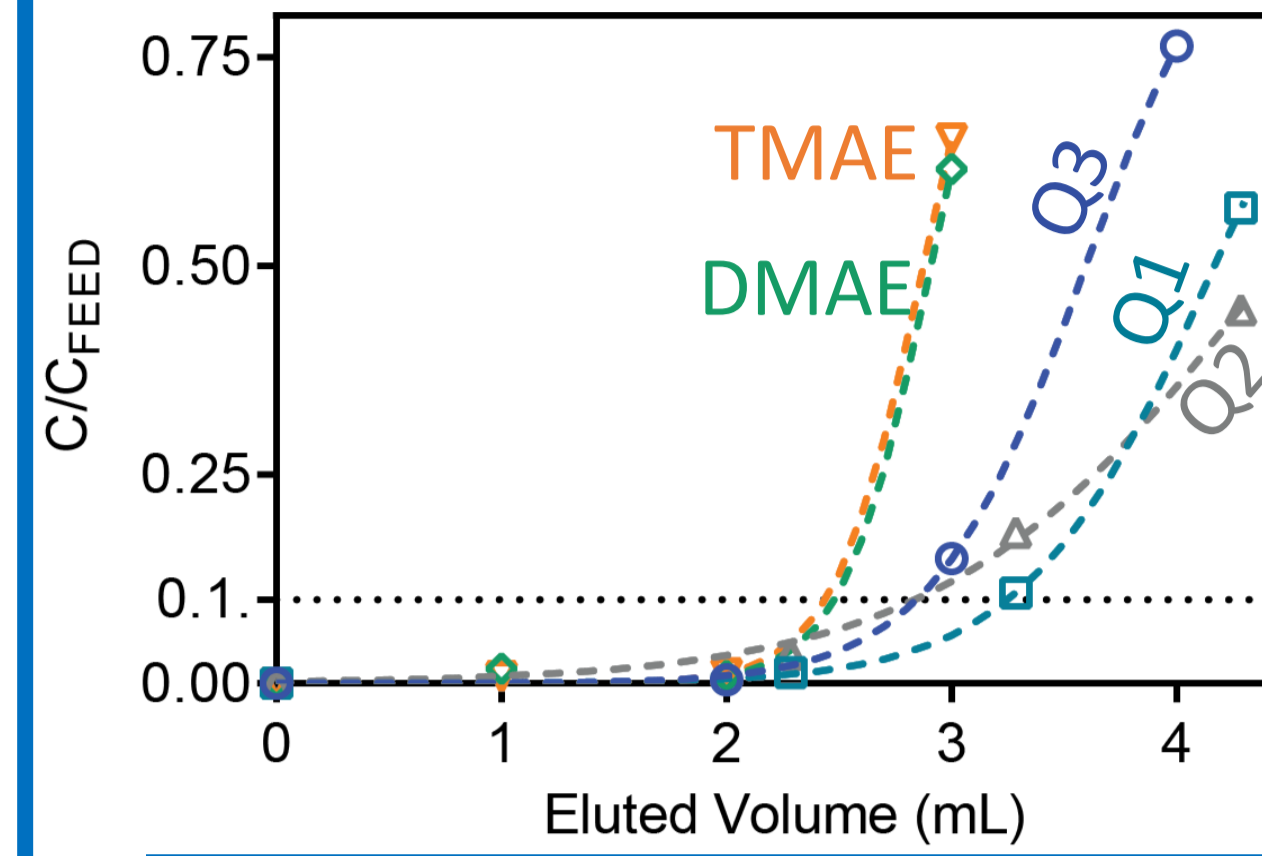


Table 1: DBC (ng p30/mL packed resin) for evaluated AEX resins and relation to ligand density

AEX Resin	Q1	Q2	Q3	TMAE	DMAE
DBC _{10%}	1310	1168	1175	923	917
LD	+++++	++	+++	+	+

Figure 3: Frontal experiments for DBC determination using scout columns (1mL)

Eshmuno 1 was selected due to its higher capacity

Single Column Batch Experiments

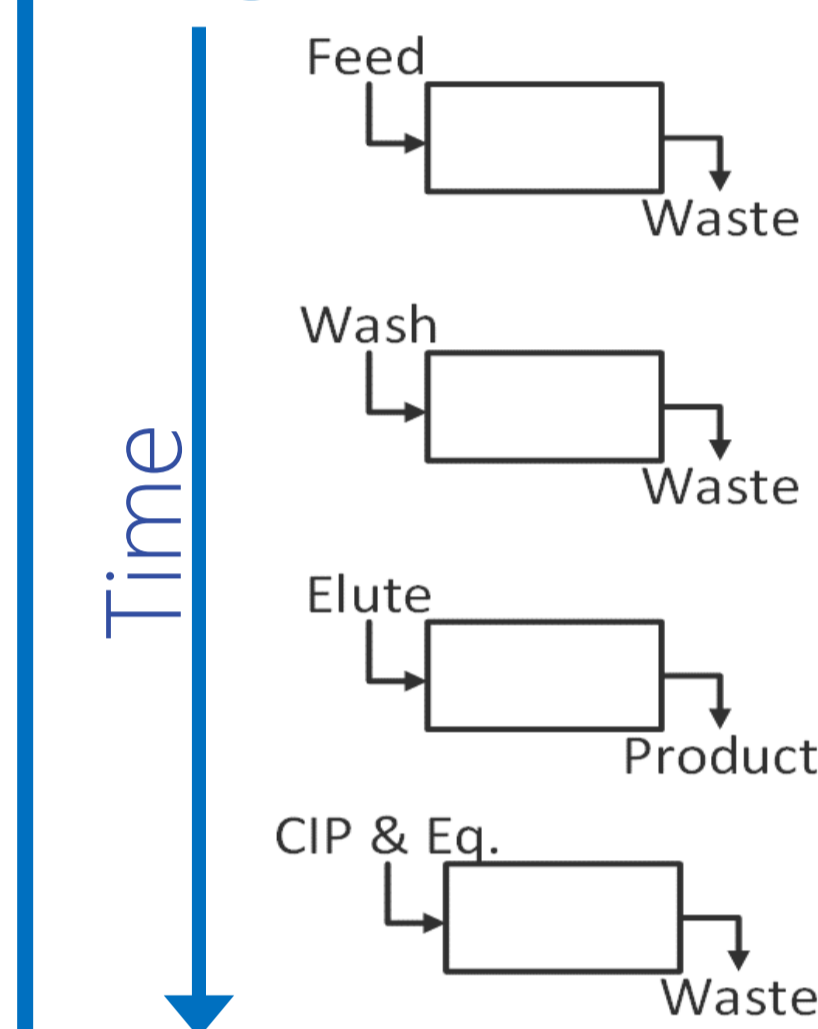


Table 2: Summary of experiments to determine optimal buffer elution conditions; Load up to 90% of column breakthrough

Elution step (mM NaCl)	400	550	700
VLPs Recov.	19 %	35 %	37 %
1 M (NaCl) Elution step Recov.	11 %	5 %	5 %
DNA Remov.	67 %	60 %	53 %
Host cell Prot. Remov.	42 %	33 %	27 %
Total Prot. Remov.	92 %	63 %	46 %

Figure 3: Sequence of steps in Bind & Elute chromatography

- A step elution with 550 mM of NaCl displayed best compromise between recovery and impurity removal;
- Loading below 10% of column breakthrough resulted in lower VLP recoveries (< 19%, data not shown);
- Loads of 90% of column breakthrough lead to a loss of VLPs in the flowthrough fraction, since the mass transfer zone (MTZ) exited the column.

Multicolumn Process Shortcut Design Approach


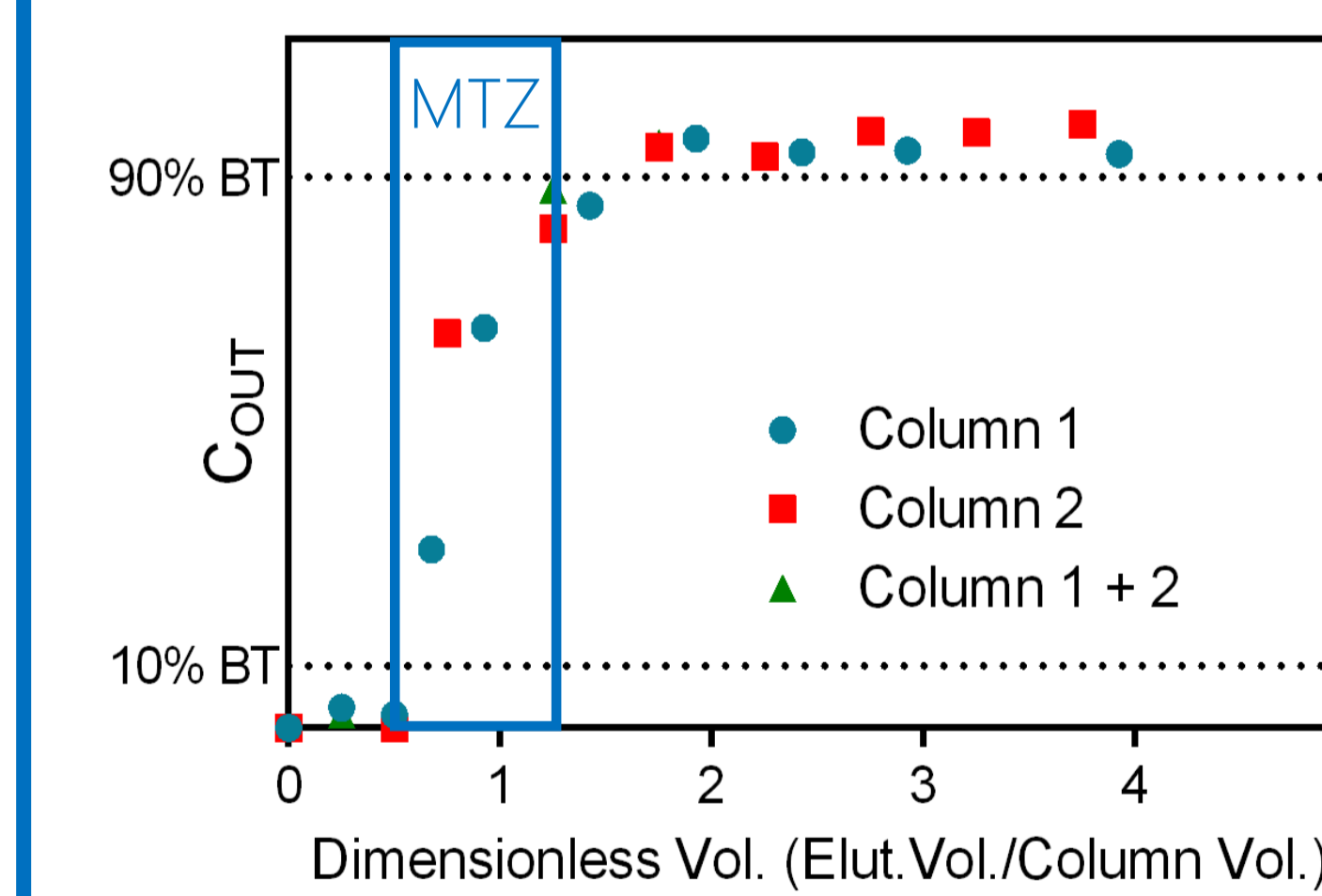


Figure 4: Frontal experiments performed in the chromatographic beds used for the multicolumn process (left); Proposed step sequence for the first half cycle of the multicolumn process (right); In the second half-cycle columns exchange position.

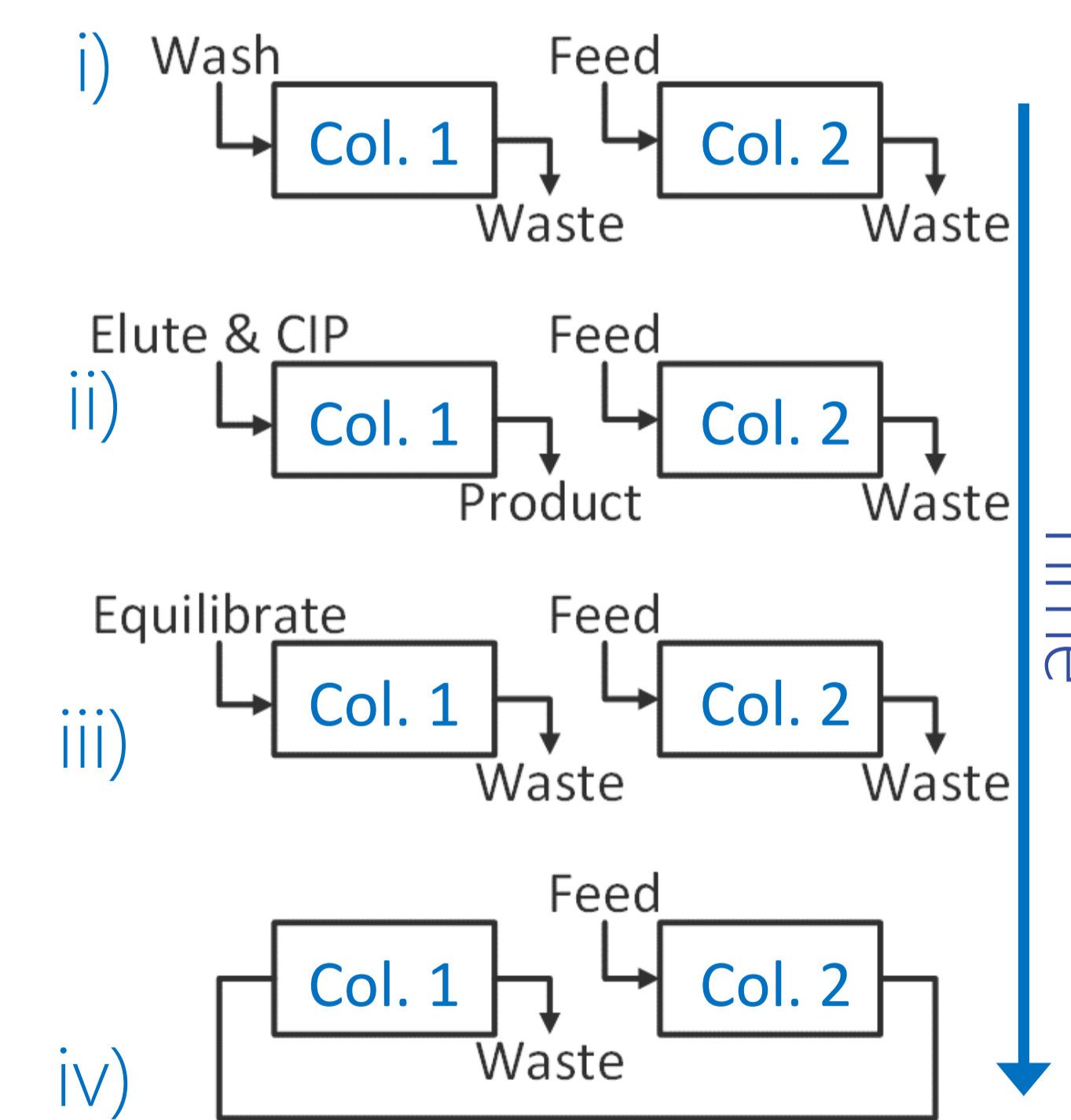

- In steps i) to iii) the downstream column is loaded up to 10% of breakthrough, while the upstream column, loaded in the previous operation cycle follows a sequence of steps similar to a batch operation;
- The effluent of the downstream column is connected to the previous column while loading is maintained, thus avoiding the exit of the MTZ and allowing a 90% of breakthrough column load;

Experimental Run

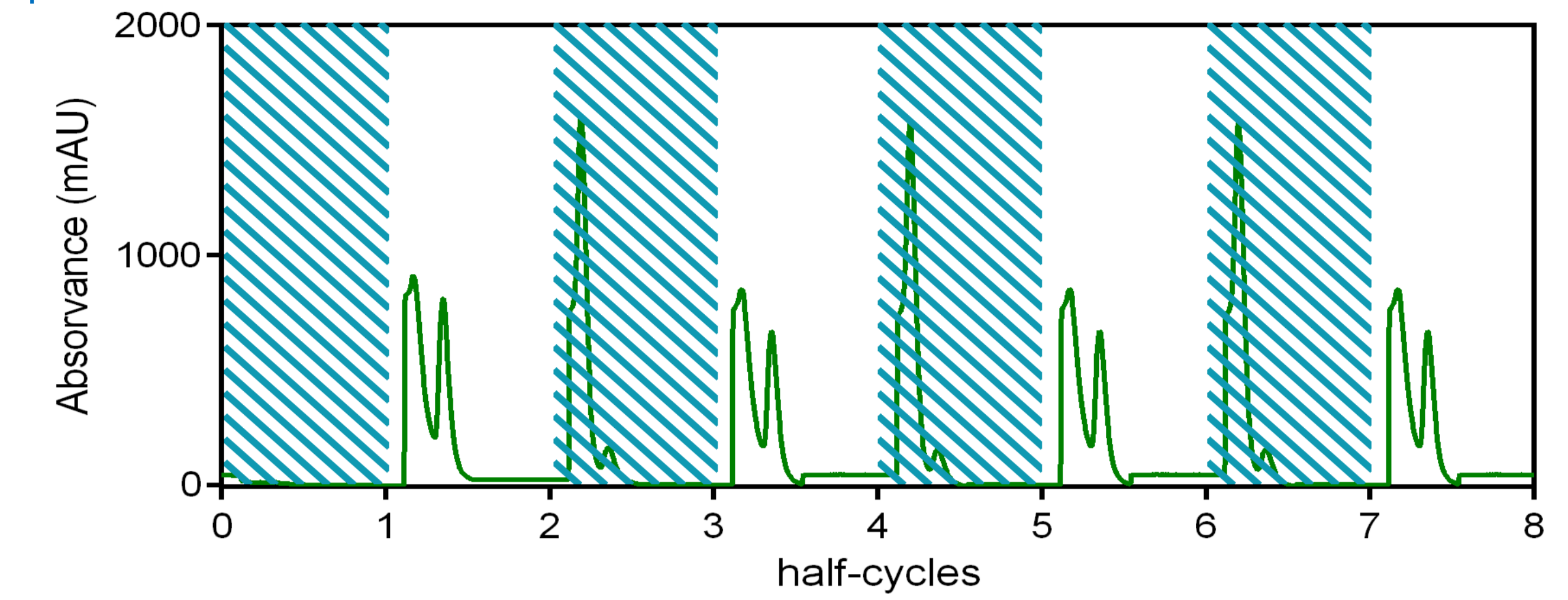


Figure 5: Experimental run of the proposed multicolumn process for continuous purification of HCV-VLPs; Shaded areas correspond to product collection in column 1; non-shaded areas correspond to the product collection in column 2;

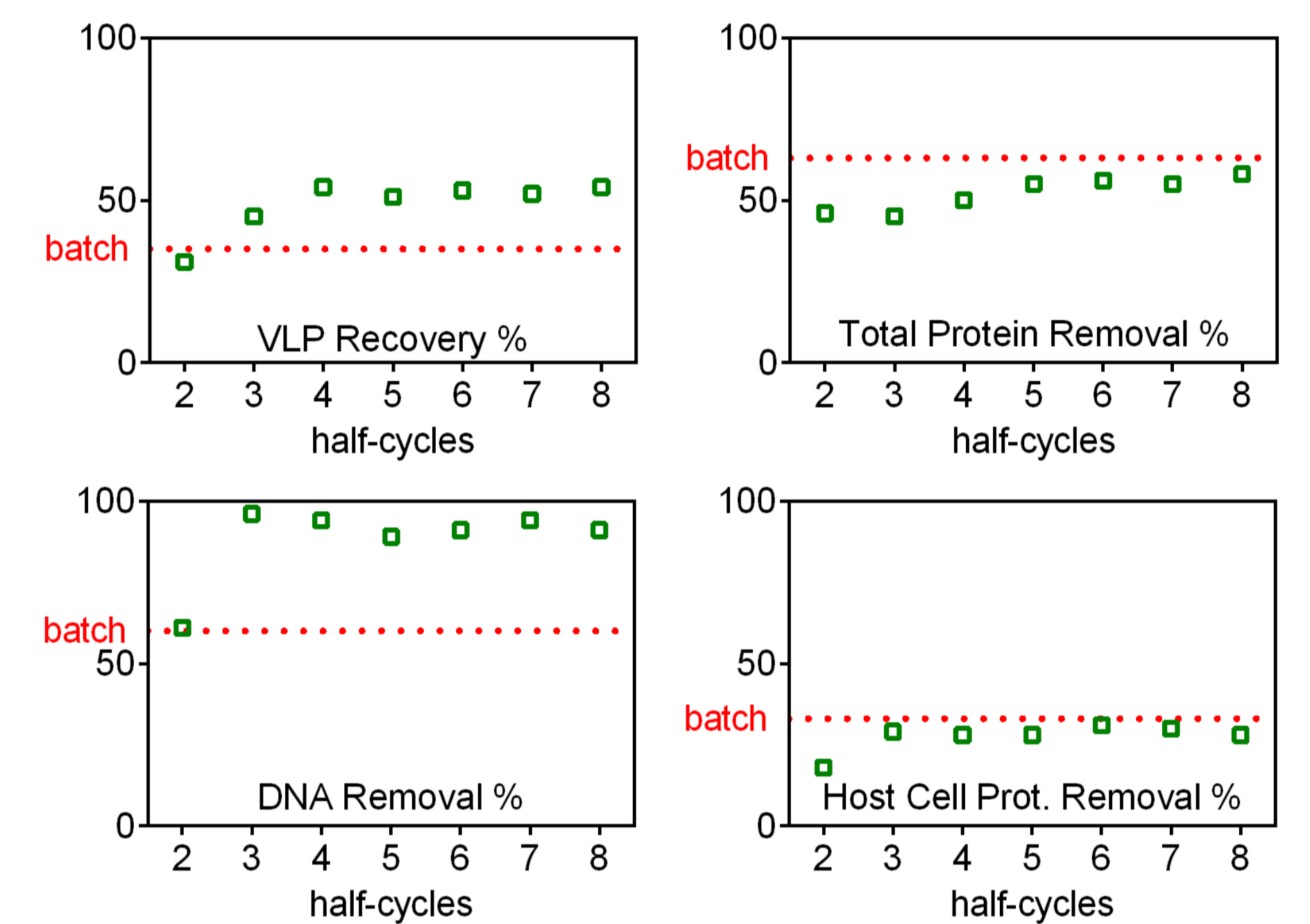


Figure 6: Comparison between batch and multicolumn runs

- Similar HCP and total protein clearance was achieved;
- Increased clearance of DNA;
- Increased recovery of HCV-VLPs

Final Remarks & Future Work

- Successful implementation of a fully continuous multicolumn chromatographic step was achieved;
- The used shortcut design method allowed to increase HCV-VLPs recovery from 35% to 55% without cycle optimization;
- By maintaining the MTZ inside the column loop it was possible to improve the resin capacity utilization;
- Due to the cyclic exchange of the columns position, the multicolumn process benefits from the advantageous simulated countercurrent movement of fluid and solid phases
- Mathematical modelling can be used to further improve cycle design and optimize both step sequence and duration;

Acknowledgements

The authors acknowledge the financial support received from the European Commission through the project H2020-BIOTEC-2014-2015/H2020-LEIT-BIO-2014-1-635770. Also, the authors acknowledge R. Castro and J. Clemente for the production of the HCV-VLPs bulk material.