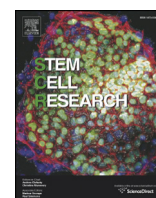


Contents lists available at [ScienceDirect](http://ScienceDirect.com)

Stem Cell Research

journal homepage: www.elsevier.com/locate/scr

Lab Resource: Stem Cell Line

Establishment of an induced pluripotent stem cell (iPSC) line from a patient with Clozapine-responder Schizophrenia

Matteo Marcatili ^{a,b,c}, Fabio Marsoner ^{a,1}, Armando D'Agostino ^{b,c}, Thodoris Karnavas ^{d,2}, Daniele Bottai ^b, Silvio Scarone ^{b,c,*}, Luciano Conti ^{a,*}^a Laboratory of Stem Cell Biology, Centre for Integrative Biology - CIBIO, Università degli Studi di Trento, Trento, Italy^b Department of Health Sciences, Università Degli Studi di Milano, Milan, Italy^c Department of Mental Health, San Paolo Hospital, Milan, Italy^d Chromatin Dynamics Unit, San Raffaele University and Research Institute, Milan, Italy

ARTICLE INFO

Article history:

Received 11 October 2016

Accepted 1 November 2016

Available online 09 November 2016

ABSTRACT

Peripheral blood mononuclear cells (PBMCs) were collected from a patient with treatment-refractory Schizophrenia who presented an exceptional clinical response to Clozapine. iPSC lines were established with a non-integrating reprogramming system based on Sendai virus. A footprint-free hiPSC line was characterized to confirm the expression of the main endogenous pluripotency markers and have a regular karyotype. Pluripotency was confirmed by differentiation into cells belonging to the three germ layers. This hiPSC line represents a valuable tool to study the molecular, biochemical and electrophysiological properties of mature neuronal populations belonging to Clozapine responder patients with a severe form of Schizophrenia.

© 2016 Published by Elsevier B.V. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Resource table.

Name of stem cell line	SCZ#4-9 iPSC
Institution	University of Trento
Person who created resource	Silvio Scarone, Luciano Conti
Contact person and email	luciano.conti@unitn.it ; silvio.scarone@unimi.it
Date archived/stock date	December 2015
Origin	Peripheral Blood Mononuclear Cells (PBMCs)
Type of resource	Induced Pluripotent Stem Cells (iPSCs) derived from a schizophrenic patient (confirmed with SCID-I) Clozapine Non-Responder
Sub-type	Induced Pluripotent Stem Cells (iPSCs)
Key transcription factors	hOCT4, hSOX2, hC-MYC, hKLF4 (CytoTune™-iPS 2.0 Sendai Reprogramming Kit - Thermo Fisher Scientific)
Authentication	Identity and purity of the cell lines was confirmed by SeV specific Polymerase chain reaction (PCR), pluripotent markers detection (Western Blot and immunocytochemistry), karyotyping, expression of specific markers of the three germ layers by means of <i>in vitro</i> differentiation

Link to related literature /
Information in public /
databases

Ethics Patient informed consent obtained; Ethics Review Board-competent authority approval was obtained from the San Paolo Hospital Ethical Board

Resource details

Blood sample was collected by a 40-year old male patient with a diagnosis of treatment-refractory Schizophrenia at the Department of Mental Health of the San Paolo Hospital, Milan (Italy). The diagnosis of Schizophrenia was confirmed by the assessment of two independent psychiatrists with the Structured Clinical Interview for DSM-IV Axis I Disorders (SCID-I).

According to specific criteria (Caspi et al., 2004), the treatment-resistance to conventional antipsychotics justified the introduction of Clozapine which resulted in an extraordinary clinical response. Indeed, the patient underwent multiple hospitalizations for approximately a decade after onset of symptoms associated with a marked psychosocial impairment. After stabilization on Clozapine, the patient was never again hospitalized and satisfied specific criteria for full recovery from Schizophrenia (Jaaskelainen et al., 2013).

To generate the SCZ#4-9 iPSC line, an integration-free gene delivery method (CytoTune™-iPS 2.0 Sendai Reprogramming Kit; Thermo Fisher Scientific) of the four Yamanaka reprogramming factors OCT4, SOX2, KLF4, and C-MYC (Takahashi et al., 2007) based on Sendai viral

* Corresponding authors.

E-mail addresses: fmar@uni-bonn.de (F. Marsoner), tk2708@cumc.columbia.edu (T. Karnavas), silvio.scarone@unimi.it (S. Scarone), luciano.conti@unitn.it (L. Conti).¹ Current address: Institute of Reconstructive Neurobiology, LIFE and BRAIN Center, University of Bonn, Bonn, Germany.² Current address: Department of Genetics & Development, Columbia University, Medical Center, New York, USA.

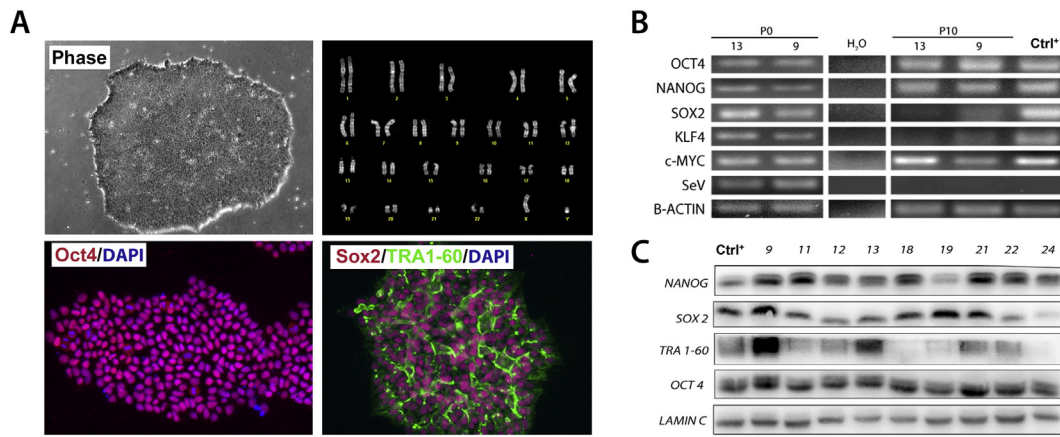


Fig. 1. Characterization of SCZ#4-9 iPSC line. A: Representative image of a SCZ#4-9 hiPSC colony (5 \times) and its karyogram showing a normal diploid 46, XY karyotype with no noticeable cytogenetic aberrations. Immunocytochemical analysis presenting the expression of the pluripotency markers OCT4, TRA-1-60, SOX2 (20 \times). B: RT-PCR presenting the expression of the pluripotency-associated transcripts in SCZ#4-9 iPSCs cultures and in another hiPSC clone (SCZ#4-13) derived from the same patient at passages 0 and 10. Absence of Sendai virus genome conservation is shown in passage 10 cultures. A commercial hiPSC line was used as positive control for pluripotency-associated transcripts. C: Western blot analysis showing expression levels of pluripotency-associated markers (NANOG, SOX2, TRA1-60, OCT4) in SCZ#4-9 iPSCs and in nine other iPSC clones obtained from the same patient.

particles (Fusaki et al., 2009; Yang et al., 2008–2012) was used. iPSC-like colonies appeared after 10–12 days and were picked one week later (Fig. 1A). One of the selected clones gave rise to the established SCZ#4-9 iPSC line with a characteristic iPSC-like morphology (Fig. 1A) and specific immunoreactivity for OCT4, SOX2 and TRA1-60 pluripotency-associated markers (Fig. 1A). PCR (Fig. 1B) and Western Blot (Fig. 1C) analyses confirmed the expression of pluripotency transcripts. SCZ#4-9 iPSC line displayed a normal diploid 46, XY karyotype, without noticeable abnormalities (Fig. 1B). The expression of Sendai virus genome in cultures at passage 0 and passage 10 was analyzed by PCR and the loss of the viral genome was confirmed in passage 10 SCZ#4-9 iPSCs (Fig. 1B).

Embryoid body assay was performed to assess the pluripotent competence SCZ#4-9 iPSC line. Cells were cultured for one week in EB suspension and for additional 7 days in adhesion to promote the *in vitro* maturation towards the three germ layer derivatives (Carpenter et al., 2003). EBs cultures at 14 days displayed the presence of differentiated cells immunoreactive for ectodermal (β 3-Tubulin), mesodermal (α -SMA) and endodermal (AFP) markers (Fig. 2A). The differentiation competency of SCZ#4-9 iPSCs was comparable to that observed for a commercial hiPSC line, as shown by the similar expression levels of transcripts for FGF5 (ectoderm marker), Nestin (neuro-ectoderm

marker), T-Brachyury (mesoderm marker), SOX-17 (endoderm marker) assessed by qRT-PCR (Fig. 2B).

Materials and methods

PBMCs collection and freezing

Peripheral Blood Mononuclear Cells (PBMCs) from patients were isolated in BD Vacutainer CPT Cell Preparation tubes with sodium citrate, after 30 min centrifugation (1800 \times g at room temperature). PBMCs were collected in PBS for a total volume of 35 ml and centrifuged at 300 \times g for 15 min RT and resuspended in fetal bovine serum (FBS) with 10% DMSO. 2×10^6 cells were aliquoted and frozen.

PBMCs thawing and reprogramming with Sendai virus particles

PBMCs were thawed at 37 $^{\circ}$ C and centrifuged at 200 \times g for 10 min in expansion medium (EM) made of StemPro-34 Serum Free Medium (SFM, Thermo Fisher Scientific) Basal Medium, StemPro-34 Nutrient Supplement, 200 mM GlutaMAX, 1% Penicillin/Streptomycin, 100 ng/ml Stem Cell Factor (SCF, Prepotech), 100 ng/ml FLT-3 (Thermo Fisher Scientific), 20 ng/ml Interleukin-6 (IL-6) (Thermo Fisher Scientific), 20 ng/ml

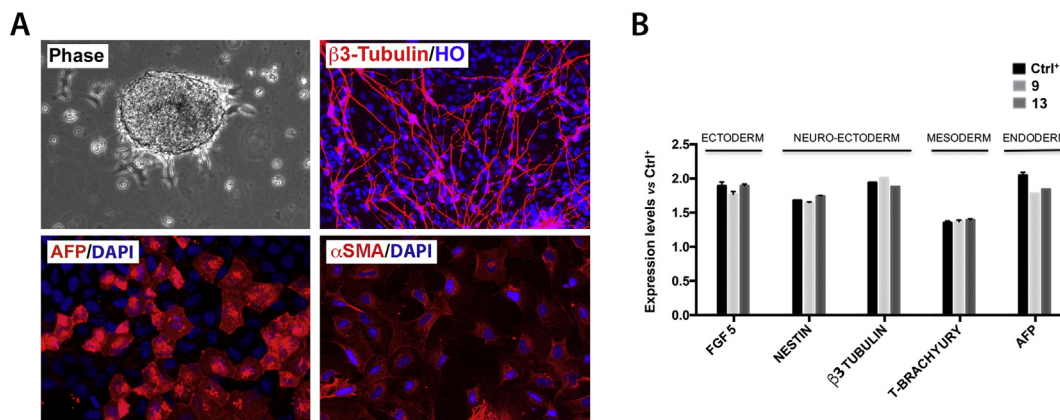


Fig. 2. *In vitro* differentiation SCZ#4-9 iPSC line. A: Embryoid Bodies obtained after 4 days of suspension culture (5 \times). D14 cultures exhibit cells immunoreactive for ectodermal (β 3-Tubulin), mesodermal (α -SMA) and endodermal (AFP) germ layer markers (20 \times). B: qRT-PCR presenting comparable expression levels of transcripts for the 3 germ-layers, FGF5 (ectoderm), Nestin and β 3-Tubulin (neuro-ectoderm), T-Brachyury (mesoderm) and AFP (endoderm) between 14 days differentiated SCZ#4-9 iPSCs, SCZ#4-13 iPSCs and differentiated reference commercial hiPSCs.

Table 1

List of primers sequences, amplicons size and number of PCR cycles.

Gene	Primer sequence	Amplicon size	Cycles (PCR)
BETA ACTIN	F: GACAGGATGCAGAAGGAGATTACTG R: CTCAGGAGGAGCAATGATCTTGAT	72 bp	25
OCT4	F: GGAAGGAATTGGGAACACAAGG R: AACTTCACCTTCCCTCCAACCA	71 bp	30
SOX2	F: GCTACAGCATGATGCAGGACCA R: TCTGCGAGCTGGTCATGGAGTT	135 bp	30
c-myc	F: CCTGGTGCTCCATGAGGAGAC R: CAGACTCTGACCTTTTGCCAGG	128 bp	30
KLF4	F: CATCTCAAGGCACACCTGCGAA R: TCGGTGCGATTTTGGCACTGG	156 bp	30
NANOG	F: CCTGTGATTTGTGGGCTG R: GACAGTCTCCGTGTGAGGCAT	78 bp	30
SEV	F: GGATCACTAGGTGATATCGAGC R: ACCAGACAAGAGTTTAAGAGATATGTATC	181 bp	35
NESTIN	F: GGAGAAGGACCAAGAAGT R: ACCTCTCTGTGGCATT	153 bp	qRT-PCR
T-Brachyury	F: CCTTCAGCAAAGTCAAGCTCACC R: TGAAGTGGGTCTCAGGGAAGCA	153 bp	qRT-PCR
FGF5	F: GGAATACGAGGAGTTTTCAGCAAC R: CTCCTGAACCTGCAGTCATCTG	99 bp	qRT-PCR
AFP	F: GCAGAGGAGATGTCTGGATTG R: CGTGGTCAGTTTGACAGCATCTG	113 bp	qRT-PCR
β3-TUBULIN	F: TCAGCGTCTACTACAACGAGGC R: GCCTGAAGAGATGTCCAAAGGC	120 bp	qRT-PCR

Interleukin-3 (IL-3) (Thermo Fisher Scientific). The medium was replaced daily for the following 3 days.

In order to deliver reprogramming genes in PBMCs, viral particles provided with the CytoTune-iPS 2.0 Sendai Reprogramming Kit (Thermo Fisher Scientific) were used following the manufacturer's protocol. Three weeks post-transduction colonies with iPSCs morphology appeared and were picked, transferred onto a new well and cultured on Geltrex-coated plastic dish in E8 according to the manufacturer's protocol.

In vitro differentiation

Embryoid body (EB) formation assay was performed by gently resuspending iPSCs clumps in 100-mm non-tissue culture-treated dish in Essential 6 medium (E6 medium; Thermo Fisher Scientific). Medium was changed daily. At day 7, EBs were collected and plated on Geltrex-coated dishes in E6 medium to allow growth in adhesion for further 7 days. Medium was changed every other day.

RNA isolation, polymerase chain reaction (PCR) and quantitative-PCR (qPCR)

RNA was isolated with the TRIzol Reagent (Thermo Fisher Scientific) following the manufacturer's protocol and reverse transcribed using iScript cDNA Synthesis Kit (BioRad). Transcripts of interest were amplified using EURO TAQ Thermostable DNA polymerase (EUROCLONE) and detection of genes of interest was confirmed with specific primers (Table 1). Quantitative Reverse Transcription-Polymerase Chain Reaction (qRT-PCR) was performed using the SsoAdvanced Universal SYBR Green Supermix Kit following the manufacturer's instructions. Beta-

actin was used as housekeeping gene to normalize data. Amplification was performed on a CFX96 BioRad RealTime PCR machine. Results were analyzed with BioRad CFX Manager dedicated software.

Immunofluorescence assay

Cells were fixed with PFA 4% for 15 min RT, permeabilized with Triton 0.5% for 15 min RT and blocked with blocking solution (10% FBS in PBS) for 1 h RT. Cultures were then incubated with specific primary antibodies overnight at 4 °C (Table 2) and stained for 45 min at RT with secondary antibody and Hoechst 33,258 1 µg/ml (Thermo Fischer Scientific). Images were detected with the microscope Leica DM IL Led Fluo with Leica DFC450 C camera (Leica Microsystem).

Western blot assay

Cultures were lysated in SDS Sample Buffer (62.5 mM Tris-HCl pH 6.8; 2% SDS; 10% Glycerol; 50 mM DTT; Bromophenol Blue). Samples were boiled at 95 °C for 5 min and loaded in the 8% polyacrylamide gel and proteins blotted on a PVDF membrane by means of Trans Blot Turbo apparatus (BioRad). Primary antibodies (Table 2) were incubated overnight at 4 °C in agitation and secondary antibody for 45 min at RT. Signal was detected with the ECL Clarity system (BioRad) in dark chamber UVITECH Cambridge (Uvitech) and Uvitech software was used to acquire and analyze the data.

Karyotyping

Cell cultures were treated with colcemid (Gibco KaryoMAX Colcemid solution in PBS, Thermo Fischer Scientific) at a final

Table 2

List of the antibodies used in for immunocytochemistry (IC) and Western Blot (WB) assays, working dilution and species in which they are produced.

Antibody	Company	Dilution	Species
TRA 1-60	Santa Cruz Biotech	1:1000 (WB) 1:200 (IC)	Mouse IgM
NANOG	Santa Cruz Biotech	1:1000	Mouse
SOX2	Millipore	1:2000 (WB) 1:300 (IC)	Rabbit
OCT4 (WB)	Santa Cruz Biotech	1:1000	Mouse
LAMIN A/C	Santa Cruz Biotech	1:1000	Rabbit
OCT4 (IC)	Santa Cruz Biotech	1:100	Rabbit
α-SMA	Sigma	1:100	Mouse
TROMA-1	Iowa DHB	1:100	Mouse
AFP	Abnova	1:50	Rabbit
β3-TUBULIN	Promega	1:1000	Mouse
Anti-rabbit HRP	BioRad	1:3000	Goat
Anti-mouse HRP	BioRad	1:3000	Goat
Anti-mouse IgM FITC	Santa Cruz Biotech	1:200	Goat
Alexa Fluor IgG anti-rabbit 568	Life Technologies	1:300	Goat
Alexa Fluor IgG anti-rabbit 488	Life Technologies	1:300	Goat
Anti-mouse IgM FITC	Santa Cruz Biotech	1:200	Goat
Alexa Fluor IgG Anti-mouse 568	Life Technologies	1:400	Donkey

concentration of 10 ng/ml for 16 h (overnight) at 37 °C and metaphases harvest was carried out according to standard protocols. Briefly, PBS washed cells were treated with hypotonic solution (0.075 M KCl for 15 min at RT) and fixed in acetic acid/methanol (1:3 v/v). Air-dried metaphase spreads slides were analyzed by QFQ banding following standard procedures. Microscope observation was performed using a Nikon Eclipse 90i (Nikon Instruments, Japan) equipped with the acquisition and analysis Genikon software (Nikon Instruments S.p.a. Italy).

Acknowledgements

We are grateful to Riccardo Ghidoni and Michele Samaja for the use of laboratory equipment at the Department of Mental Health. This study was supported by a CIBIO start-up grant from the University of Trento (40201803) (L.C.).

References

- Carpenter, M.K., Rosler, E., Rao, M.S., 2003. Characterization and differentiation of human embryonic stem cells. *Cloning Stem Cells* 5 (1), 79–88.
- Caspi, A., Davidson, M., Tammimga, C.A., 2004. Treatment-refractory schizophrenia. *Dialogues Clin. Neurosci.* 6 (1), 61–70.
- Fusaki, N., Ban, H., Nishiyama, A., Saeki, K., Hasegawa, M., 2009. Efficient induction of transgene-free human pluripotent stem cells using a vector based on Sendai virus, an RNA virus that does not integrate into the host genome. *Proc. Jpn. Acad. Ser. B Phys. Biol. Sci.* 85 (8), 348–362.
- Jaaskelainen, E., Juola, P., Hirvonen, N., McGrath, J.J., Saha, S., Isohanni, M., Veijola, J., Miettinen, J., 2013. A systematic review and meta-analysis of recovery in schizophrenia. *Schizophr. Bull.* 39 (6), 1296–1306.
- Takahashi, K., Tanabe, K., Ohnuki, M., Narita, M., Ichisaka, T., Tomoda, K., Yamanaka, S., 2007. Induction of pluripotent stem cells from adult human fibroblasts by defined factors. *Cell* 131 (5), 861–872.
- Yang, W., Mills, J.A., Sullivan, S., Liu, Y., French, D.L., Gadue, P., 2008–2012. iPSC Reprogramming from Human Peripheral Blood Using Sendai Virus Mediated Gene Transfer. *StemBook* [Internet]. Harvard Stem Cell Institute, Cambridge (MA) .