

Software

Open Access

## SimArray: a user-friendly and user-configurable microarray design tool

Richard P Auburn\*, Roslin R Russell, Bettina Fischer, Lisa A Meadows, Santiago Sevillano Matilla and Steven Russell

Address: Department of Genetics, University of Cambridge, Cambridge, UK

Email: Richard P Auburn\* - [r.auburn@gen.cam.ac.uk](mailto:r.auburn@gen.cam.ac.uk); Roslin R Russell - [r.russell@gen.cam.ac.uk](mailto:r.russell@gen.cam.ac.uk); Bettina Fischer - [b.fischer@gen.cam.ac.uk](mailto:b.fischer@gen.cam.ac.uk); Lisa A Meadows - [l.meadows@gen.cam.ac.uk](mailto:l.meadows@gen.cam.ac.uk); Santiago Sevillano Matilla - [s.sevillano@gen.cam.ac.uk](mailto:s.sevillano@gen.cam.ac.uk); Steven Russell - [s.russell@gen.cam.ac.uk](mailto:s.russell@gen.cam.ac.uk)

\* Corresponding author

Published: 01 March 2006

Received: 11 October 2005

BMC Bioinformatics 2006, 7:102 doi:10.1186/1471-2105-7-102

Accepted: 01 March 2006

This article is available from: <http://www.biomedcentral.com/1471-2105/7/102>

© 2006 Auburn et al; licensee BioMed Central Ltd.

This is an Open Access article distributed under the terms of the Creative Commons Attribution License (<http://creativecommons.org/licenses/by/2.0>), which permits unrestricted use, distribution, and reproduction in any medium, provided the original work is properly cited.

### Abstract

**Background:** Microarrays were first developed to assess gene expression but are now also used to map protein-binding sites and to assess allelic variation between individuals. Regardless of the intended application, efficient production and appropriate array design are key determinants of experimental success. Inefficient production can make larger-scale studies prohibitively expensive, whereas poor array design makes normalisation and data analysis problematic.

**Results:** We have developed a user-friendly tool, SimArray, which generates a randomised spot layout, computes a maximum meta-grid area, and estimates the print time, in response to user-specified design decisions. Selected parameters include: the number of probes to be printed; the microtitre plate format; the printing pin configuration, and the achievable spot density. SimArray is compatible with all current robotic spotters that employ 96-, 384- or 1536-well microtitre plates, and can be configured to reflect most production environments. Print time and maximum meta-grid area estimates facilitate evaluation of each array design for its suitability. Randomisation of the spot layout facilitates correction of systematic biases by normalisation.

**Conclusion:** SimArray is intended to help both established researchers and those new to the microarray field to develop microarray designs with randomised spot layouts that are compatible with their specific production environment. SimArray is an open-source program and is available from <http://www.flychip.org.uk/SimArray/>.

### Background

The full utility of the spotted microarray format is clearly reflected in the range of its applications. Transcriptome arrays, containing cDNA, gDNA, or oligonucleotide probes, are used to measure differential gene expression [1-5]. Whole-genome arrays, typically composed of tiled gDNA or oligonucleotides [6], have been used to identify *in vivo* sites of protein-DNA interactions [7,8] or allelic

variation [9,10]. Whilst these applications dominate, other formats, for example antibody arrays, facilitate analysis of protein and small-molecule analytes [11,12]. Thus, spotted microarrays enable high-throughput, cost-effective, and large-scale analysis of molecular interactions.

Robotic spotters deposit probes as an ordered array by repetition of a simple multi-step procedure [13-15]. First,

the print tool is positioned over the first batch of probes to be printed, arranged in 96-, 384- or 1536-well microtitre plates. Second, the spotting pins are filled by capillary action with probe material. This step is often called a source visit. Third, the probe material is deposited on chemically-modified glass slides [14,16]. Finally, the pins are cleaned, to prevent cross-contamination between subsequent spot depositions, before re-filling and printing the next batch of probes. The diversity of instrumentation, spotting pins, and reagents available, mean that procedures may be refined for optimal throughput, spot density, morphology, and consistency [16]. Whilst this facilitates production of high-quality arrays, it can also lead to significant differences between facilities with regard to the instrumentation, protocols, and reagents employed.

Robotic spotters are supplied with sophisticated software to convert operator inputs to the precise list of instructions needed by the arrayer, *e.g.*, how often each source visit is to be printed, and at which spot location [13-15]. Most spotters, however, are not supplied with adequate array design tools. Operators are instead left to develop suitable spot layouts in an *ad hoc* fashion. This often leads to sub-optimal designs with spots positioned according to print order, thus juxtaposing replicates, when a non-sequential or randomised spot layout can help to control for confounding spatial effects [17,18]. For example, biases caused by inconsistent probe concentrations in the microtitre plates [18,19] and local variations in hybridisation or washing efficiency [20-23], also see Additional file 1. Whilst random noise can be overcome with simple replication and averaging, systematic biases must be specifically addressed by randomisation and normalisation [24-27].

There is therefore a need for a microarray design tool that can generate randomised source visits and, in the case of some instruments, permit the use of variable numbers of replicates, since current spot density constraints and the need for genome-wide coverage mean that replication is often limited to the normalisation controls. Exogenous or 'spike' controls, *i.e.*, probes that are complementary to targets not present in the genome of interest, can be employed for this purpose [28,29]. Print time and maximum meta-grid area estimates enable users to evaluate the suitability of the array design.

We have addressed the current lack of such microarray design tools by developing SimArray, a user-friendly and user-configurable program that generates a randomised spot layout, computes the maximum meta-grid area, and estimates printing time, in response to user-defined design decisions. The user enters these parameters by running SimArray twice. The first run produces the source

visit list that can be edited to include variable numbers of replicates, or for specific source visits to be omitted when plates are partially filled. The second run processes this source list to create the spot layout, maximum meta-grid area and estimated print time. User-configurable files mean that SimArray can be adapted to most production environments.

### Implementation

SimArray was developed in Perl version 5.6.1 and 5.8.3, under both Windows and UNIX operating systems. SimArray can be run under Windows (after installing Perl; for example, ActivePerl [30]) and UNIX.

#### Before the first run

Download the three configuration files and an 'index.sa' file from the SimArray web site [31]. Configuration files that describe instrument-specific pin configurations and achievable spot densities have already been created. If a suitable file is not available, an existing one should be downloaded and edited. An example file for the user to record their specific print cycle times is also available for editing. The 'index.sa' file should then be updated, to record the locations and names of the configuration files. The configuration files will then only need to be re-edited if the printing environment is altered.

#### First run

**Probe number:** enter the number of microarray probes (or wells) to be printed.

**Plate format:** select an appropriate plate format.

**Tools available:** select an appropriate pin configuration.

**Source visits:** the source visit list is generated for editing.

#### Second run


**Required spot density:** SimArray counts the number of spots to be printed.

**Pin type:** select the spotting pin to be used.

**Evaluate pin selection:** SimArray evaluates whether the selected pin is compatible with the required spot density.

**Compute spots\_x and spots\_y:** SimArray computes and displays the sub-grid dimensions that fall between the target spot number and a user-specified upper limit, for the user to select.

**Compute print time:** select an appropriate print set-up, SimArray then calculates the estimated print time.

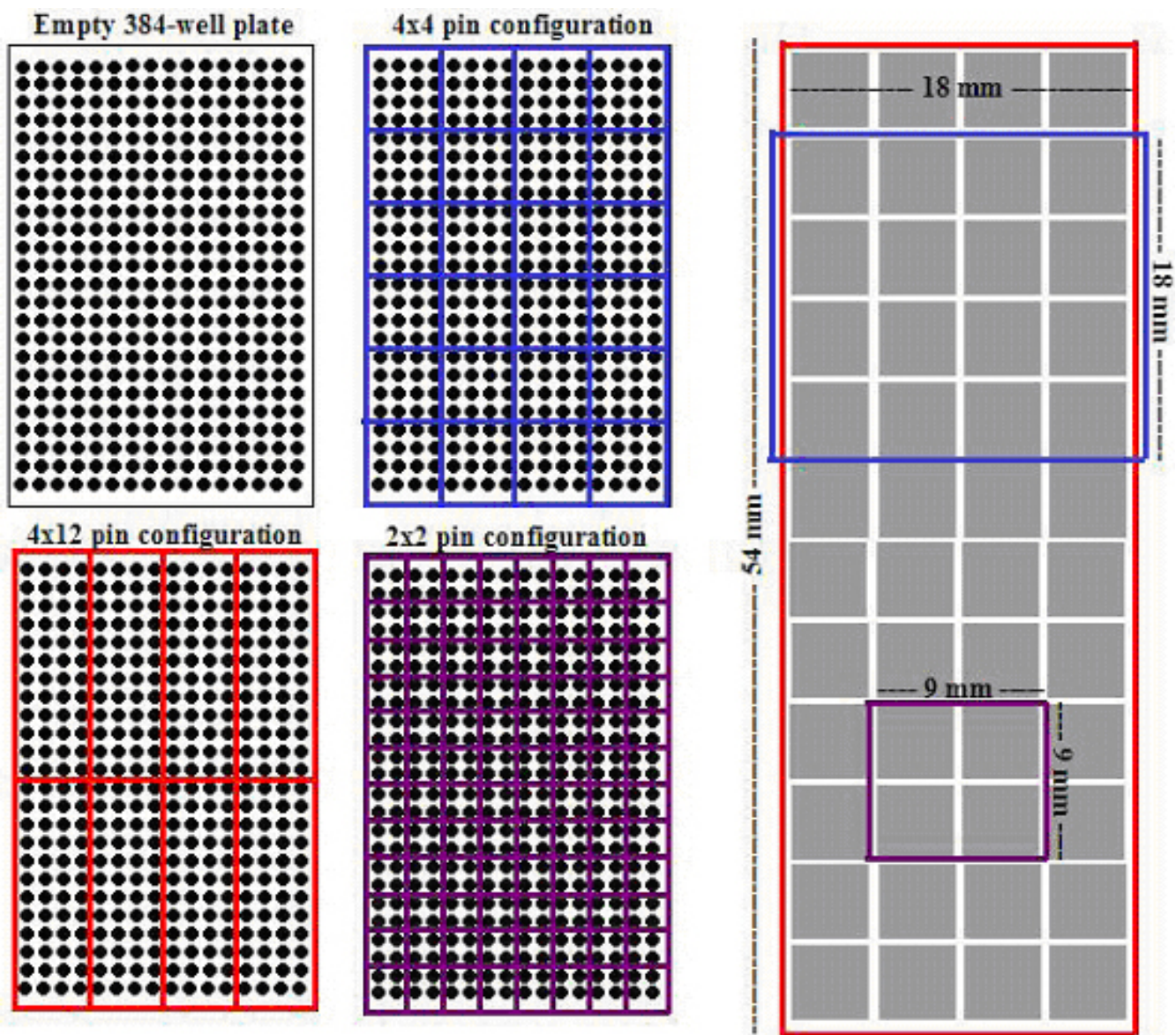

**Figure 1**  
**Pin configuration directly affects array area and source visit number.** Each pin in the print tool will print a single patch or sub-grid of spots. The sum total of sub-grids printed by one print tool is called a meta-grid. Pin configuration directly affects maximum meta-grid area because the spotting pins have a fixed plate-specific pitch to enable adjacent pins in the tool to enter adjacent wells of the microtitre plates. 96-well print tools have a pitch of 9 mm, 384-well print tools have a pitch of 4.5 mm, and 1536-well print tools have a 2.25 mm pitch. SimArray estimates maximum meta-grid area by simply multiplying the number of pins in each axis by the pin pitch. SimArray will thus overestimate the meta-grid area, when a reduced spot pitch is employed because spot pitch is not taken into consideration. Pin configuration also defines how many spots need to be printed per sub-grid, as more pins means that fewer pin loadings or source visits need to be performed. For example, a 4 × 12 print tool will print all probes in a 384-well plate with just 8 source visits, meaning that each pin will print just 8 spots per sub-grid. Whereas, a 4 × 4 tool would require 24 source visits, and a 2 × 2 tool requires 96 source visits.

**Summary report:** SimArray generates a report containing a summary of the user's responses, the randomised spot layout, an estimated print time, and the maximum meta-grid area.

**After the second run**  
 The randomised source visit map can be directly uploaded to instruments that accept either comma-, tab-, or space-separated values source files, or manually entered. Micro-

```
#Version 1.0 (Standard)
#
#Plate : Microtitre plate format
#Tool : Pin configuration name
#Pins(X) : Number of pins (x-axis)
#Pins(Y) : Number of pins (y-axis)
#Grids : Number of sub-grids, per meta-grid
#
#Plate Tool Pins(X) Pins(Y) Grids
96 1x1 1 1 1
96 2x2 2 2 4
96 2x3 2 3 6
96 2x4 2 4 8
96 2x6 2 6 12
384 1x1 1 1 1
384 2x2 2 2 4
384 4x4 4 4 16
384 4x6 4 6 24
384 4x8 4 8 32
384 4x12 4 12 48
1536 1x1 1 1 1
1536 2x2 2 2 4
1536 4x4 4 4 16
1536 8x8 8 8 64
1536 8x16 8 16 64
1536 8x24 8 24 192
```

**Figure 2**

**Example 'tool' user-configurable file.** Print tools are composed of spotting pins that are arranged in a pre-defined plate-format-specific configuration (Fig. 1.). Each user records the pin configurations that are available to them by adding their robotic spotter's list of options to the tool file. Pin configuration directly affects source visit number and maximum meta-grid area (Fig. 1). Tool files for the Qarray2 (Genetix), MicroGrid, and OmniGrid (Genomic Solutions) robotic spotters are available from the SimArray web site [31].

arrays can then be manufactured, with spots no longer positioned according to print order. Standard robotic spotter data tracking software can be used to record which probe is present at each spot location.

**Results and Discussion**

Pin configuration affects the number of source visits that must be performed and the maximum meta-grid area (Fig. 1). For these reasons, the user is required to enter the probe number (an integer) and to select a microtitre plate format (e.g., 96, 384 or 1536), before selecting a compatible pin configuration (Fig. 2). SimArray then prints a source visit list for the user to edit (Fig. 3). If the number of probes to be printed is not compatible with the selected pin configuration, SimArray will round up the source visit number to the nearest whole number, as robotic spotters can only print with a full complement of spotting pins. In such instances, the last source visit to be printed would include some empty wells. Users are, however, able to specify any number of replicates, for any number of source visits. Additionally, source visits can be omitted by setting the replicate number to zero. These features enable array designs with odd numbers of probes, variable num-

```
#Probes: 15360
#Plate format: 384
#Print tool: 4x12
#Max area (mm): 18x54
#Plate Source Replicates
0 0 0
1 1 4
1 2 4
1 3 4
1 4 4
1 5 4
1 6 4
1 7 4
1 8 4
2 9 4
2 10 4
2 11 4
2 12 4
2 13 4
2 14 4
2 15 0
2 16 0
3 17 1
3 18 1
3 19 1
3 20 1
3 21 1
3 22 1
3 23 1
3 24 1
4 25 1
4 26 1
4 27 1
4 28 1
4 29 1
4 30 1
4 31 1
4 32 1
Etc.
```

**Figure 3**

**Example SimArray source visit list.** The source list is generated after the first run of SimArray. Headers are marked with a hash at the beginning of the line and provide a summary of the specified design decisions. The replicate column can be edited to record how many times each source visit is to be printed. Source visits with zero replicates are ignored and omitted from the final array design. Other columns and the header should not be modified. An example file is available from the SimArray web site [31].

```
#Version 1.0 (Genetix, aQu)
#
#Pin      : Spotting pin (type, model, etc.)
#Arrayer  : Spotting instrument (type, model, etc.)
#SpotsX   : Maximum number of spots per sub-grid (x-axis)
#SpotsY   : Maximum number of spots per sub-grid (y-axis)
#Total    : Total number of spots per sub-grid
#Comment  : Free text 'comment' for the user
#
#Pin      Arrayer  SpotsX   SpotsY   Total    Comment
aQu65    Qarray2   44        44        1936    Genetix
aQu75    Qarray2   34        34        1156    Genetix
aQu150   Qarray2   18        18        325     Genetix
```

**Figure 4**  
**Example 'pins' user-configurable file.** Achievable spot densities vary according to the production environment, *i.e.*, spotting pin, spotting buffer, substrate, temperature, relative humidity, *etc.* The available pins and their achievable spot densities are recorded by editing the pins file. These data will directly affect selection of the sub-grid dimensions, *i.e.*, the number of spots that the user would like to print in the x- and y-axis. Pin files for Telechem, Genetix, Matrix, and Genomic Solutions spotting pins are available from the SimArray web site [31].

bers of replicates, and non-sequential source visits to be processed.

Maximum meta-grid area is calculated by simply multiplying the number of pins in each axis by the pin tool's pre-defined pin pitch (Fig. 1). Consequently, SimArray does not take spot pitch into consideration and can overestimate meta-grid areas, especially for low-density arrays that are printed with reduced spot pitches. Since high-density arrays limit the scope for reducing spot pitch, we believe this is a reasonable approach because SimArray will be of most use when designing higher-density arrays. Additionally, most operators print microarrays with the spot pitch set to the near-maximum distance permissible to reduce the probability that neighbouring spots printed by the same pin will be merged together. Prediction of the maximum meta-grid area will at least allow users to decide whether it is possible for them to hybridise the array, *e.g.*, when the hybridisation area is constrained by automated hybridisation stations.


The number of spots to be printed per sub-grid is calculated by counting the number of spots that are specified in the source visit list. This total is displayed and users are asked to select a suitable spotting pin (Fig. 4). The selected pin is evaluated and the script exits if the pin's achievable spot number per sub-grid is incompatible with the required target spot number per sub-grid. Exiting the script at this stage, if a problem is found, removes the risk of downstream errors and provides an opportunity for the array design to be modified, or for a different spotting pin to be selected.

```
#Version 1.0 (BioRobotics, MicroGrid II)
#
#Set_up   : Print setting employed
#Arrayer  : Spotting instrument (type, model, etc.)
#Time     : Time for one print cycle, in minutes
#Comment  : Free text 'comment' for the user
#
#Set_up   Arrayer  Mins     Comment
Salt      MGII610   1.8      90 slides
Salt      MGII610   1.4      40 slides
DMSO     MGII610   1.7      90 slides
DMSO     MGII610   1.3      40 slides
```

**Figure 5**  
**Example 'time' user-configurable file.** Print time is dependent on a range of instrument and user-defined parameters. To estimate the print time, users record the time it takes to perform a single print cycle for a given print set-up by editing the time file. Single print cycle duration is best measured by recording how long it takes to perform a full print-run (under defined conditions) and then dividing by the number of the print cycles. Estimated print time is limited by the user-specified print cycle times. An additional example file is available from the SimArray web site [31].

Sub-grid dimensions, *i.e.*, the number of spots in the x and y-axis, which are compatible with the target spot number per sub-grid are then calculated, and users select an option from a list of compatible choices. To limit the length of this list, SimArray will only display sub-grid dimensions that are equal to or greater than the target spot number per sub-grid and less than a user specified limit. The upper limit is the target spot number per sub-grid, plus the user-specified 'spot number margin'. SimArray prevents users from selecting grid dimensions that are incompatible with the spotting pins' maximum achievable spot density. If, however, the selected sub-grid dimensions permit more than the required number of spots to be printed, additional spot locations are flagged as blanks by assigning them a source visit number of zero, *i.e.*, not printed. SimArray will fail at this stage if there are no viable sub-grid dimension between the minimum and maximum target. We therefore recommend using a 'spot number margin' of at least ten.


Print time is dependent on the number of source visits (Fig. 1), the number of slides to be printed, a range of (perhaps) user-defined options, *e.g.*, pre-blotting, contact speed, *etc.*, and the hardware itself, *e.g.*, microtitre plate handling, x-y-z-axis motor speeds, pin (or tool) travel distances, *etc.* Additionally, wash conditions vary according to the production environment, *i.e.*, spotting pin, spotting buffer, *etc.* Print time is therefore calculated after the user has selected the intended print setting from a list of available options (Fig. 5). The list of options includes the user defined 'single print cycle duration', *i.e.*, the time it takes to perform a single print cycle from source visit to wash/dry cycle. This has to be determined empirically because it

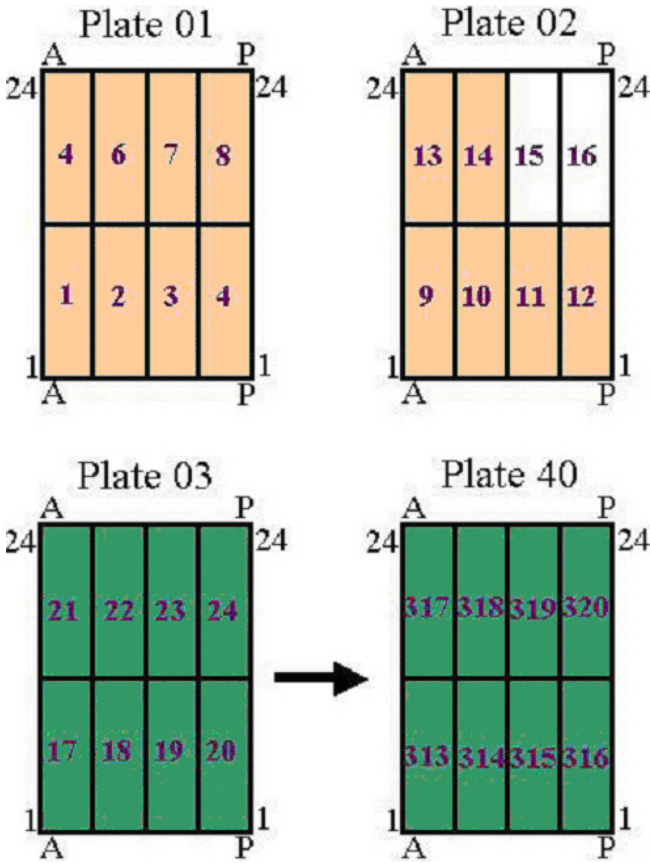

**Figure 6**  
**Example library.** These probes are to be printed by a 48-pin print-head, with the pins arranged in a 4 × 12 configuration. The library consists of 14 independent exogenous controls (orange) and 14,592 transcript-specific probes (green). Since the exogenous controls must be printed by all pins, these must be arranged in 4 × 12 blocks of 48 wells. Conversely, the 14,592 target probes only need to be randomly distributed, to facilitate printing once per meta-grid and to ameliorate systematic biases [18]. Source visit number and hence print order are indicated in the figure.

is production environment and instrument dependent. Print time is calculated by simply multiplying the number of print cycles that need to be performed, by the time taken to perform each print cycle. Print time accuracy is therefore dependent on the user-specified single print cycle estimates.

Finally, a report containing the randomised source visit map along with a summary of the user's responses, an estimated print time, and the estimated maximum meta-grid area is generated (Fig. 6). The user can specify a comma-, tab-, or space-separated source visit map with the command line keys -C, -T, or -S (default), respectively. The source visit map can then either be directly uploaded to instruments that accept source files in these formats, or

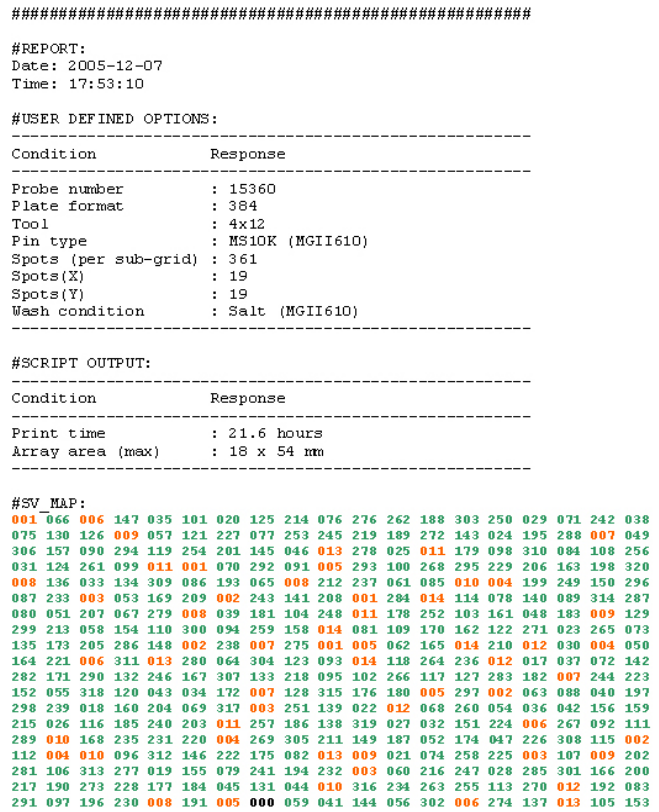

**Figure 7**  
**Example SimArray summary report.** The report is generated after the second run of SimArray and includes a summary of the user-specified design decisions, estimated print run time, estimated maximum meta-grid area, and the randomised source visit map. The source visit map can be either comma-, tab- or space-delimited (default). In this example, and in Figure 8, the design requirement was to print 4 × copies of the exogenous spike control source visits (orange), 1 copy of each transcript-specific probe source visit (green), and omit the empty wells. This report was generated whilst simulating printing by a MicroGrid II 610 robotic spotter (Genomic Solutions). This figure has been colour-coded to aid comparison with Figures 6 and 8.

manually entered. Microarrays can then be manufactured, with the spots and replicates positioned randomly, rather than according to their print order. Standard instrument data tracking software can be used to document what probe is present in each spot. New array designs are appended to the existing report to provide a full record of all array designs. Users are therefore able to perform multiple 'simulated print runs', with different configurations to compare the results, *i.e.*, the estimated print times and maximum meta-grid areas. Each array design includes a date and time stamp.

The user-configurable files described above and in Figures 2, 3 and 4, maximise the utility of this tool because they

```
#####
#REPORT:
Date: 2005-12-07
Time: 17:14:18

#USER DEFINED OPTIONS:
-----
Condition      Response
-----
Probe number   : 15360
Plate format   : 384
Tool           : 4x12
Pin type       : aQu75 (Qarray2)
Spots (per sub-grid) : 361
Spots(X)       : 19
Spots(Y)       : 19
Wash condition : Salt (Qarray2)
-----

#SCRIPT OUTPUT:
-----
Condition      Response
-----
Print time     : 13.8 hours
Array area (max) : 18 x 54 mm
-----

#SV_MAP:
095 009 138 230 311 035 314 025 099 268 102 212 214 001 128 006 120 100 043
266 167 122 113 000 269 030 150 198 301 097 209 060 006 007 134 036 260 073
009 315 172 287 005 277 027 013 164 126 204 158 180 010 285 318 003 175 020
105 174 238 171 202 286 279 006 185 270 297 093 293 163 137 123 151 264 092
010 034 101 190 003 159 013 005 098 220 298 176 111 005 208 139 216 039 283
045 114 169 274 153 064 021 155 278 165 072 316 284 080 001 081 067 041 047
011 058 003 179 031 210 295 002 129 267 004 259 065 199 147 008 087 182 023
024 312 273 289 233 244 201 094 104 162 306 242 275 001 004 253 052 272 003
002 046 265 079 183 319 022 091 082 149 118 002 207 280 228 135 181 044 004
117 140 048 085 011 090 005 006 254 308 227 009 192 032 271 157 109 217 282
205 076 281 130 249 014 168 276 004 017 303 061 074 071 056 258 222 059 014
026 194 075 195 189 057 068 237 300 196 232 263 213 001 156 112 302 245 215
110 040 257 241 224 243 294 290 211 166 193 012 009 145 028 033 262 012 136
218 116 115 223 197 161 007 089 148 225 125 008 051 146 014 310 008 011 231
246 007 133 240 234 236 141 010 107 251 062 119 292 186 055 042 142 178 086
078 131 007 077 200 177 296 054 038 034 037 203 083 252 012 247 221 170 132
154 239 011 313 256 317 261 248 188 010 309 029 320 053 160 307 088 173 066
049 013 019 050 144 106 127 152 084 124 012 013 018 063 288 206 299 108 096
008 226 219 235 070 069 304 291 143 002 184 250 103 229 255 187 305 191 121
```

**Figure 8**  
**Alternative SimArray array report.** The parameters described in Figure 7 were used to simulate printing the library described in Figure 6 with a Qarray2 (Genetix). In this case, the print time has reduced from almost 22 hours to 14 hours, showing that either the Qarray2 is the faster instrument or the MicroGrid II print conditions are not optimal. This example demonstrates that SimArray rapidly facilitates refinement of microarray designs, whilst generating randomised source visit maps that can be used to program robotic spotters. This figure has been colour-coded to aid comparison with Figures 6 and 7.

enable a range of production environments to be explicitly modelled. All SimArray configuration files contain a header, which includes a key to the file's contents to help with this task. If required, additional comments can also be added because all lines marked with a hash at the start are automatically ignored when SimArray reads these files. However, column meaning and order is fixed and must be preserved. We aim to develop an on-line library of configuration files at the SimArray web site [31].

**A fully worked simple example**

For this worked example, we compare the performance of a MicroGrid II (Genomic Solutions) and a Qarray2 (Genetix) instrument, printing a 15 K probe library. The

configuration files were edited to reflect the specific set of printing conditions for each robotic spotter. The example library consists of two probe types: transcript-specific probes and exogenous controls, along with some empty wells (Figure 6). The design requirement is to print single copies of the transcript-specific probes, and for every pin to print quadruplet spots for the exogenous control probes, randomising the distribution of elements on the array, whilst omitting the empty wells. The probes were arranged in the spotting plates, according to these design criterion (Figure 6).

SimArray was used to generate a randomised spot layout for each instrument, to assess which would be better suited to printing this library. The SimArray simulated print-runs indicated that the MicroGrid II spotter would take 56% longer to print microarrays according to the specified criterion (Figs 7 and 8). The estimated print times agreed with how long it would take to print the arrays, provided no manual intervention, e.g., refilling of wash solutions, etc., was required. The user can now either enter the randomised spot layout and print this microarray design with the Qarray2, or re-evaluate whether the print settings for the MicroGrid II were optimal.

SimArray also permits further simulations, allowing an evaluation of alternative pin configurations, replicate numbers, slide numbers and instrument configurations, subject to the availability of appropriately annotated configuration files. This further ensures that microarrays of an ideal design can be generated, whilst permitting each to be evaluated for its compatibility to the local production environment. Printing with different pin configurations, however, requires the spotting library itself to be redesigned, as spotting pins enter adjacent wells of the microtitre plate and the probes must be arranged accordingly (Figs 1 and 6). This suggests that spot layouts should be defined before the spotting probes are transferred to microtitre plates for printing.

**Conclusion**

We have developed a user-friendly microarray design tool, SimArray, which generates a randomised spot layout, computes a maximum meta-grid area, and an estimated print time, in response to user-specified design decisions. SimArray is of general utility for all users of robotic spotters and can be configured to suit individual production environments.

**Availability and requirements**

Project name: SimArray

Project home page: <http://www.flychip.org.uk/SimArray/>

Operating system: Windows and UNIX

**Programming language:** Perl version 5.6.1 (and more recent versions)

**Other requirements:** tool (Fig. 2), pins (Fig. 4), and time (Fig. 5) configuration files

**Licence:** free

**Any restrictions to use by non-academics:** none

### List of abbreviations used

cDNA Complementary DNA

gDNA Genomic DNA

**Meta-grid** The sub-grids (and spots) that are printed by one print tool

**Sub-grid** The patch of spots that is printed by a single pin

### Authors' contributions

RRR suggested that array design could be automated. RPA wrote the source code, with technical guidance from RRR. RPA, BF, LAM, and SSM tested the code and validated performance. The manuscript was prepared by RPA. SR is the group leader, providing funding, critical assessment and general guidance. All authors read and approved the final manuscript.

### Additional material

#### Additional File 1

This dataset includes heat maps for six different microarrays that were printed and hybridised by the UK Drosophila microarray facility. The heat maps demonstrate that spot location has a direct impact on the measured differential gene expression ratios and hence supports the argument for randomisation of the spot layout.

Click here for file

[<http://www.biomedcentral.com/content/supplementary/1471-2105-7-102-S1.doc>]

### Acknowledgements

The authors would like to thank Gos Micklem, David Kreil, and the four anonymous referees for constructive criticism of the manuscript, and François Guillier for help with the web site. This work was supported by research grants from the BBSRC.

### References

- Schena M, Shalon D, Davis RW, Brown PO: **Quantitative monitoring of gene expression patterns with a complementary DNA microarray.** *Science* 1995, **270**:467-470.
- Schena M, Shalon D, Heller R, Chai A, Brown PO, Davis RW: **Parallel human genome analysis: microarray-based expression monitoring of 1000 genes.** *Proc Natl Acad Sci U S A* 1996, **93**:10614-10619.
- Hayward RE, Derisi JL, Alfidhli S, Kaslow DC, Brown PO, Rathod PK: **Shotgun DNA microarrays and stage-specific gene expression in Plasmodium falciparum malaria.** *Mol Microbiol* 2000, **35**:6-14.
- Zaigler A, Schuster SC, Soppa J: **Construction and usage of a one-fold-coverage shotgun DNA microarray to characterize the metabolism of the archaeon Haloferax volcanii.** *Mol Microbiol* 2003, **48**:1089-1105.
- Religio A, Schwager C, Richter A, Ansoorge W, Valcarcel J: **Optimization of oligonucleotide-based DNA microarrays.** *Nucleic Acids Res* 2002, **30**:e51.
- Mockler TC, Ecker JR: **Applications of DNA tiling arrays for whole-genome analysis.** *Genomics* 2005, **85**:1-15.
- Ren B, Robert F, Wyrick JJ, Aparicio O, Jennings EG, Simon I, Zeitlinger J, Schreiber J, Hannett N, Kanin E, Volkert TL, Wilson CJ, Bell SP, Young RA: **Genome-wide location and function of DNA binding proteins.** *Science* 2000, **290**:2306-2309.
- Mukherjee S, Berger MF, Jona G, Wang XS, Muzzev D, Snyder M, Young RA, Bulky ML: **Rapid analysis of the DNA-binding specificities of transcription factors with DNA microarrays.** *Nat Genet* 2004, **36**:1331-1339.
- Pinkel D, Segraves R, Sudar D, Clark S, Poole I, Kowbel D, Collins C, Kuo WL, Chen C, Zhai Y, Dairkee SH, Ljung BM, Gray JW, Albertson DG: **High resolution analysis of DNA copy number variation using comparative genomic hybridization to microarrays.** *Nat Genet* 1998, **20**:207-211.
- Pinkel D, Albertson DG: **Comparative genomic hybridization.** *Annu Rev Genomics Hum Genet* 2005, **6**:331-354.
- Lueking A, Cahill DJ, Mullner S: **Protein biochips: A new and versatile platform technology for molecular medicine.** *Drug Discov Today* 2005, **10**:789-794.
- Chiosis G, Brodsky JL: **Small molecule microarrays: from proteins to mammalian cells - are we there yet?** *Trends Biotechnol* 2005, **23**:271-274.
- Cheung VG, Morley M, Aguilar F, Massimi A, Kucherlapati R, Childs G: **Making and reading microarrays.** *Nat Genet* 1999, **21**:15-19.
- Hegde P, Qi R, Abernathy K, Gay C, Dharap S, Gaspard R, Hughes JE, Snesrud E, Lee N, Quackenbush J: **A concise guide to cDNA microarray analysis.** *Biotechniques* 2000, **29**:548-50, 552-4, 556 passim.
- Affara NA: **Resource and hardware options for microarray-based experimentation.** *Brief Funct Genomic Proteomic* 2003, **2**:7-20.
- Auburn RP, Kreil DP, Meadows LA, Fischer B, Matilla SS, Russell S: **Robotic spotting of cDNA and oligonucleotide microarrays.** *Trends Biotechnol* 2005, **23**:374-379.
- Wernisch L, Kendall SL, Soneji S, Wietzorrek A, Parish T, Hinds J, Butcher PD, Stoker NG: **Analysis of whole-genome microarray replicates using mixed models.** *Bioinformatics* 2003, **19**:53-61.
- Qian J, Kluger Y, Yu H, Gerstein M: **Identification and correction of spurious spatial correlations in microarray data.** *Biotechniques* 2003, **35**:42-4, 46, 48.
- Spruill SE, Lu J, Hardy S, Weir B: **Assessing sources of variability in microarray gene expression data.** *Biotechniques* 2002, **33**:916-20, 922-3.
- Yang YH, Buckley MJ, Dudoit S, Speed TP: **Comparison of methods for image analysis on cDNA microarray data.** *J Comp Graph Stat* 2002, **11**:108-136.
- Yang YH, Dudoit S, Luu P, Lin DM, Peng V, Ngai J, Speed TP: **Normalization for cDNA microarray data: a robust composite method addressing single and multiple slide systematic variation.** *Nucleic Acids Res* 2002, **30**:e15.
- Kreil DP, Russell RR: **There is no silver bullet--a guide to low-level data transforms and normalisation methods for microarray data.** *Brief Bioinform* 2005, **6**:86-97.
- Futschik ME, Crompton T: **OLIN: optimized normalization, visualization and quality testing of two-channel microarray data.** *Bioinformatics* 2005, **21**:1724-1726.
- Kerr MK, Churchill GA: **Statistical design and the analysis of gene expression microarray data.** *Genet Res* 2001, **77**:123-128.
- Churchill GA: **Fundamentals of experimental design for cDNA microarrays.** *Nat Genet* 2002, **32** Suppl:490-495.
- Brodsky L, Leontovich A, Shtutman M, Feinstein E: **Identification and handling of artifactual gene expression profiles emerging in microarray hybridization experiments.** *Nucleic Acids Res* 2004, **32**:e46.
- Le Meur N, Lamirault G, Bihouee A, Steenman M, Bedrine-Ferran H, Teusan R, Ramstein G, Leger JJ: **A dynamic, web-accessible**



**resource to process raw microarray scan data into consolidated gene expression values: importance of replication.** *Nucleic Acids Res* 2004, **32**:5349-5358.

28. van Bakel H, Holstege FC: **In control: systematic assessment of microarray performance.** *EMBO Rep* 2004, **5**:964-969.
29. Baker SC, Bauer SR, Beyer RP, Brenton JD, Bromley B, Burrill J, Causton HC, Conley MP, Elespuru R, Fero M, Foy C, Fuscoe J, Gao X, Gerhold DL, Gilles P, Goodsaid F, Guo X, Hackett J, Hockett RD, Ikononi P, Irizarry RA, Kawasaki ES, Kaysser-Kranich T, Kerr K, Kiser G, Koch WH, Lee KY, Liu C, Liu ZL, Lucas A, Manohar CF, Miyada G, Modrusan Z, Parkes H, Puri RK, Reid L, Ryder TB, Salit M, Samaha RR, Scherf U, Sendera TJ, Setterquist RA, Shi L, Shippy R, Soriano JV, Wagar EA, Warrington JA, Williams M, Wilmer F, Wilson M, Wolber PK, Wu X, Zadrow R: **The External RNA Controls Consortium: a progress report.** *Nat Methods* 2005, **2**:731-734.
30. **ActivePerl.** : [ <http://www.activestate.com/ActivePerl/> ].
31. **SimArray.** : [ <http://www.flychip.org.uk/SimArray> ].

Publish with **BioMed Central** and every scientist can read your work free of charge

*"BioMed Central will be the most significant development for disseminating the results of biomedical research in our lifetime."*

Sir Paul Nurse, Cancer Research UK

Your research papers will be:

- available free of charge to the entire biomedical community
- peer reviewed and published immediately upon acceptance
- cited in PubMed and archived on PubMed Central
- yours — you keep the copyright

Submit your manuscript here:  
[http://www.biomedcentral.com/info/publishing\\_adv.asp](http://www.biomedcentral.com/info/publishing_adv.asp)

