

*Review*

Luminescence- and Fluorescence-Based Complementation Assays to Screen for GPCR Oligomerization: Current State of the Art

Elise Wouters ^{1,†}, Lakshmi Vasudevan ^{1,2,†}, René A.J. Crans ¹, Deepak K. Saini ²
and Christophe P. Stove ^{1,*}

¹ Laboratory of Toxicology, Department of Bioanalysis, Faculty of Pharmaceutical Sciences, Ghent University, Ottergemsesteenweg 460, 9000 Ghent, Belgium; elise.wouters@ugent.be (E.W.); lakshmi.vasudevan@ugent.be (L.V.); rene.crans@ugent.be (R.A.J.C.)

² Department of Molecular Reproduction, Development and Genetics, Indian Institute of Science, Bangalore 560012, India; deepaksaini@iisc.ac.in

* Correspondence: christophe.Stove@ugent.be; Tel.: +32-9-264-81-35

† These authors contributed equally to this work.

Received: 30 April 2019; Accepted: 12 June 2019; Published: 17 June 2019

Abstract: G protein-coupled receptors (GPCRs) have the propensity to form homo- and heterodimers. Dysfunction of these dimers has been associated with multiple diseases, e.g., pre-eclampsia, schizophrenia, and depression, among others. Over the past two decades, considerable efforts have been made towards the development of screening assays for studying these GPCR dimer complexes in living cells. As a first step, a robust in vitro assay in an overexpression system is essential to identify and characterize specific GPCR–GPCR interactions, followed by methodologies to demonstrate association at endogenous levels and eventually in vivo. This review focuses on protein complementation assays (PCAs) which have been utilized to study GPCR oligomerization. These approaches are typically fluorescence- and luminescence-based, making identification and localization of protein–protein interactions feasible. The GPCRs of interest are fused to complementary fluorescent or luminescent fragments that, upon GPCR di- or oligomerization, may reconstitute to a functional reporter, of which the activity can be measured. Various protein complementation assays have the disadvantage that the interaction between the reconstituted split fragments is irreversible, which can lead to false positive read-outs. Reversible systems offer several advantages, as they do not only allow to follow the kinetics of GPCR–GPCR interactions, but also allow evaluation of receptor complex modulation by ligands (either agonists or antagonists). Protein complementation assays may be used for high throughput screenings as well, which is highly relevant given the growing interest and effort to identify small molecule drugs that could potentially target disease-relevant dimers. In addition to providing an overview on how PCAs have allowed to gain better insights into GPCR–GPCR interactions, this review also aims at providing practical guidance on how to perform PCA-based assays.

Keywords: G protein-coupled receptor (GPCR), dimerization; oligomerization; protein complementation assay; bimolecular fluorescence complementation (BiFC) assay; bimolecular luminescence complementation (BiLC) assay

1. Introduction

Membrane receptors are the key players in mediating communication between the cell and the extracellular space. The G protein-coupled receptor (GPCR) protein family represents one of the

largest group of cell membrane signaling proteins in the human genome [1]. Currently, approximately 800 members of this superfamily are classified according to the International Union of Pharmacology, Committee on Receptor Nomenclature and Classification (NC-IUPHAR). Herein, GPCRs are grouped into three classes, namely Class A (rhodopsin-like), Class B (secretin receptor family), and Class C (metabotropic glutamate-like) [2]. The remaining GPCRs—fungal mating pheromone receptors, cyclic AMP receptors, and frizzled/smoothed—are also included in alternative classification systems [3].

G protein-coupled receptors are activated by a panel of different triggers, including photons, nucleotides, ions, hormones, peptides, chemokines, and others. Notwithstanding the broad collection of ligands to activate specific GPCRs, they all share a common structure, namely, an N-terminal domain, seven transmembrane α -helices, extra- (ECL1–3) and intracellular loops (ICL1–3) and a C-terminal domain. Specific domains within these GPCRs show plasticity, since conformational changes are clearly distinguishable between their activated and inactive states, as observed from the crystal structures [4]. Upon activation, the motile domains, such as the ICL3 and the C-terminus, are crucial for interaction with intracellular signaling partners, i.e., G-proteins, GPCR kinases (GRKs) and β -arrestins. Given their crucial role in cellular signaling, altered activity or deviating expression levels of GPCRs have frequently been correlated with diseases, including neurodegenerative disorders, depression, and cancer [5–10].

Overall, their abundance, regulation in pathophysiology of diverse disease areas, accessibility at the cell surface, and availability of druggable binding sites have made GPCRs the largest drug target family. It is, hence, not surprising that 34% of the current FDA-approved drugs on the market target GPCRs [11]. Interestingly, virtually all pharmaceutical companies have primarily focused on GPCR monomers as a target for drug development. However, the recent knowledge that GPCRs form homo-oligomeric or hetero-oligomeric complexes has opened up unexplored avenues in the development of new drug therapies (e.g., small ligands, peptides, and bivalent ligands) [12–15].

1.1. GPCR Dimerization

There are approximately 150,000 to 300,000 protein–protein interactions in the human interactome [16,17] and constant efforts are ever-expanding this number. This also holds true for the interest in GPCR–GPCR interactions, as evidenced by the remarkable rise in publications in this field of research, with over 300 reports (Figure 1) covering GPCR multimers during the last 5 years.

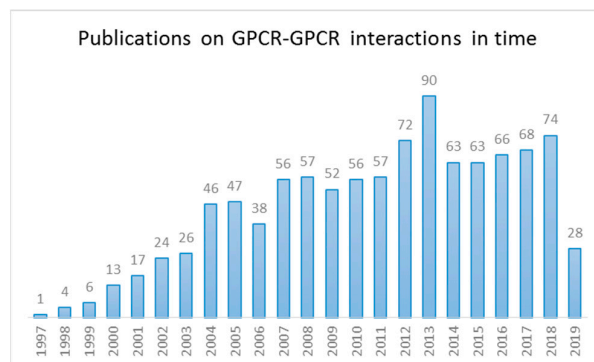


Figure 1. Number of publications on G protein-coupled receptor (GPCR) dimer and oligomers from 1997–2019 (data incomplete for 2019). (Source: NCBI—a search for the terms “GPCR* AND DIMER*” or “GPCR* AND OLIGOMER*” in the PubMed database).

While the traditional view on GPCRs was a “one ligand-one monomeric GPCR unit-one signaling protein” model, an emerging body of evidence shows that GPCRs can form homo- or hetero-oligomeric complexes. This phenomenon adds an additional layer to the complexity of GPCR signaling. Dysfunction of these dimers has been associated with multiple diseases, e.g., pre-eclampsia (angiotensin II type 1 and bradykinin receptor B2 heterodimer complexes) [18],

schizophrenia (dopamine D2 receptor (D₂R) heterocomplexes) [19], depression (serotonin 1A receptor 5-HT_{1A} and serotonin 7 receptor 5-HT₇ heterocomplexes, galanin receptor type 1 (GalR1), type 2 (GalR2), and 5-HT_{1A} heterocomplexes) [20–26].

Although the concept of GPCR dimerization has been controversial over the past two decades, several sophisticated studies have led scientists to the validation of the GPCR dimerization theory. An indisputable and widely accepted example is the obligate dimer of Class C GABA receptors [27]. These receptors are obligatory dimers, as GABAB1 is necessary for ligand binding and GABAB2 ensures efficient cell trafficking and downstream signaling. Both protomers are not functional when expressed alone. In addition to obligate GPCR dimers observed in Class C receptors, an increasing amount of convincing evidence suggests there is also -transient dimer formation between Class A GPCRs [28–30]. Conflicting observations both in favor of and against oligomerization have been reviewed by Bouvier and Herbert (2014) [31] and by Lambert and Javitch (2014) [32]. Recent insights into Class A GPCR dimerization have shed new light on this concept, i.e., a large-scale comparative study of 60 members of the Rhodopsin family revealed that only a small fraction, i.e., 20%, forms dimers [33]. The dimer formation of Class A GPCRs seems to have a clustered distribution, and furthermore, a striking correlation has been suggested between receptor organization and size of the GPCR family. A restricted family size, as for the small Glutamate family, correlates with predominantly dimeric behavior. This also seems to be true for the receptors exhibiting less diversity like the frizzled GPCRs [33].

Besides the fact that Class A GPCRs seem to dimerize to a lesser extent as compared to Class C receptors, the dimerization event itself is also much more complicated, i.e., Class A GPCRs often seem to interact at multiple interfaces, thus causing a broad impact on functionality [34]. This is noted in the case for the muscarinic acetylcholine M₁ or M₂ dimerization, wherein reports have suggested that the sites of contact are transient and could involve multiple regions of the receptor [35–37]. For the D₂R, this transient dimerization phenomenon has been studied in more detail, with a lifetime of 68 ms being assigned to the dimer. This dimer lifetime could be prolonged with a factor of approximately 1.5 by agonist addition (i.e., dopamine and quinpirole) [38,39], whereas incubation with the antagonist spiperone reduced the level of D₂R homodimerization by more than 40% [40], indicating a potential role of the dimer in signaling cascades. Dijkman and colleagues introduced the “rolling dimer” interface model for the neurotensin receptor 1 homodimer. This concept unites earlier seemingly conflicting data since multiple dimer conformations could co-exist that interconvert during the dimer lifetime, by rotation of the monomers relative to one another.

In general, the phenomenon of oligomerization can exert a significant impact on receptor–ligand binding, downstream signaling, crosstalk, internalization, and trafficking. For instance, the interaction of D₂R with dopamine D3 receptor (D₃R) results in a higher potency of certain anti-parkinsonian agents, like pramipexole, compared to their monomers [28]. Furthermore, selective stimulation of D₁R or D₂R or both by the neurotransmitter dopamine triggered co-internalization of the D₁R–D₂R heterodimer [29].

1.2. Studying GPCR–GPCR Interactions: Biochemical Methods

In order to target and study GPCR–GPCR interactions, multiple *in vitro* (and *in vivo*) biochemical and physical, and functional studies have been reported in the literature, which described approaches such as co-immunoprecipitation (co-IP), proximity ligation assays (PLA), and fluorescence- or bioluminescence-based techniques [41].

The presence of multiple protein–protein interactions (PPIs) has been demonstrated by co-IP or pull-down assays [42–44]. On a non-reducing SDS-PAGE (sodium dodecyl sulfate-polyacrylamide gel electrophoresis), it has been reported that GPCRs often migrate in a way that corresponds to twice their molecular weight. Co-IP of endogenous GPCRs entirely relies on the selectivity of the antibodies used. Furthermore, this technique necessitates the disruption of biological samples and the solubilization of membrane proteins and generally provides no information about the physical interaction.

On the other hand, fluorescence- or bioluminescence-based energy transfer techniques can be applied in living cells, thereby overcoming some of the limitations of the classical biophysical techniques. For GPCR dimerization, this has been reviewed by Ciruela et al. (2010) [45]. Fluorescence and bioluminescence resonance energy transfer (FRET and BRET) involve the non-radiative energy transfer between a donor and acceptor, which allows the examination of PPIs in their native cell context by monitoring emissions from the acceptor species. Fluorescence resonance energy transfer microscopy offers the advantage that it allows visualization of the interaction inside living cells. Bioluminescence resonance energy transfer avoids some of the troubleshooting involved with FRET such as photobleaching, spectral bleed-through, etc., but does not provide spatial information on the interaction. Below, we will focus primarily on recent advances in yet another biochemical method, namely, bioluminescence and fluorescence-based complementation techniques that have been applied to study GPCR–GPCR interactions.

1.3. Protein Complementation Assays

Protein complementation assays (PCAs), also referred to as split systems, have been implemented over the past 2 decades to study PPIs. In these assays, a reporter protein with enzymatic or fluorescent properties is engineered or “split” into non-active or non-fluorescent fragments. These moieties are fused to potential interacting proteins. Upon interaction, the fragments will be brought into close proximity and re-assemble spontaneously into a functional biosensor. Although these assays do not formally prove direct protein–protein interactions, they do suggest co-localization of the GPCRs of interest since PCAs rely on the fusion partners “interacting”. Complementation-based assays comprise bimolecular fluorescence and luminescence complementation (BiFC/BiLC) as well as β -Lactamase and β -Galactosidase complementation assays. Since 1997, several “split” fluorescent and luminescent proteins have been developed, of which the most commonly used ones are shown in the timeline in Figure 2.

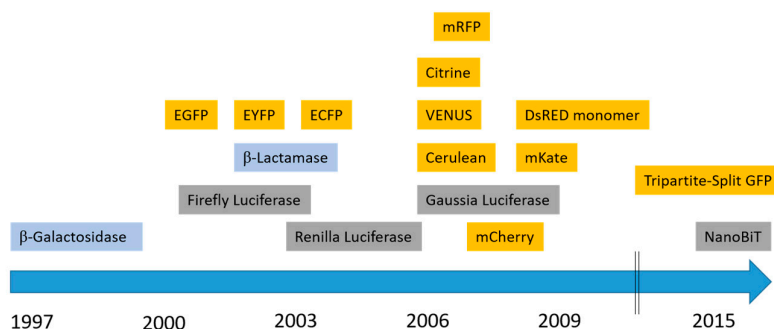


Figure 2. History of the development of split protein biosensors. Commonly used fluorescent split proteins (yellow), luminescent split reporters (grey), and enzymes with colorigenic substrates (blue) are shown.

These assays have the potential to be implemented as screening assays to identify GPCR–GPCR interactions, which later on should be confirmed by a variety of complementary approaches, for example, Spatial Intensity Distribution Analysis (SpIDA), advanced microscopy techniques, and molecular dynamics, among others [35,39,46–49]. A benefit of a robust complementation assay is not only its operational simplicity but even more the possibility to test if dimers occur constitutively or whether ligands can alter the oligomeric state of GPCRs. Moreover, these cell-based assays lend themselves for early-stage drug discovery, since molecules that potentially exert an impact on the level of dimerization, whether desired or undesired, can be rapidly identified. In addition, off-target effects can be revealed, such as altered receptor expression, localization or cellular signaling. An overview of the available PCAs as of today and their characteristics are presented in Table 1.

Table 1. Overview of proteins, already implemented in fluorescence- and luminescence-based complementation assays. N/A = not applicable, nd = not determined.

Reporter Protein	Source Species	Readout	Excitation Wavelength (nm)	Emission Wavelength (nm)	Substrate	Cofactor	Stability (h)	Maturation Time (t1/2) (min)	MW (kDa)
VENUS	<i>Aequorea victoria</i>	Fluorescence	515	528	-	N/A	-	40 (in vitro)	27
GFP	<i>Aequorea victoria</i>	Fluorescence	488	510	-	N/A	-	53 (in vitro)	27
mCherry	<i>Discosoma</i>	Fluorescence	587	610	-	N/A	-	17 + 30 (<i>S. cerevisiae</i>)	29
Cerulean	<i>Aequorea victoria</i>	Fluorescence	433	475	-	N/A	-	nd	27
Tripartite-Split GFP	<i>Aequorea victoria</i>	Fluorescence	488	530	-	N/A	-	nd	23
EYFP	<i>Aequorea victoria</i>	Fluorescence	514	527	-	N/A	-	23 (in vitro)	26.4
ECFP	<i>Aequorea macrodactyla</i>	Fluorescence	405	485	-	N/A	-	49 (<i>S. cerevisiae</i>)	26.8
Citrine	<i>Aequorea victoria</i>	Fluorescence	516	529	-	N/A	-	nd	27
mRFP	<i>Discosoma striata</i>	Fluorescence	584	607	-	N/A	-	<60	25.9
mKate	<i>Discosoma striata</i>	Fluorescence	588	635	-	N/A	-	20	26
DsRed monomer	<i>Discosoma striata</i>	Fluorescence	558	583	-	N/A	-	600	28
Renilla luciferase	<i>Renilla reniformis</i>	Luminescence	-	480	Coelenterazine	N/A	4.5 h (cell)	-	36
Firefly Luciferase	<i>Photinus pyralis</i>	Luminescence	-	550–570	D-luciferin	ATP, O ₂	4.0 h (cell)	-	62
Gaussia Luciferase	<i>Gaussia princeps</i>	Luminescence	-	485	Coelenterazine	N/A	60 h (cell media)	-	20
NanoBiT	<i>Oplophorus gracilirostris</i>	Luminescence	-	460	Furimazine	N/A	6.0 h (cell)	-	19
β-lactamase	<i>Bacillus licheniformis</i>	Luminescence	-	492	Nitrocefin	N/A	nd	-	29
β-Galactosidase	<i>Escherichia coli</i>	Fluorescence	Reliant on the substrate	Reliant on the substrate	FDG, MUG a.o.	Mg ²⁺	1.1 h (yeast cells)	-	464
Click Beetle luciferase	<i>Pyrophorus plagiophthalmus</i>	Luminescence	-	-	D-luciferin	Mg ²⁺ , ATP	-	-	64

Many questions remain unanswered in the field of GPCR di- and oligomerization. Therefore, it is a highly warranted first step to fully characterize the *in vitro* assays featuring robust sensitivity in detecting PPIs of interest. This review aims at offering insight into the progress that has been made in the field of protein complementation assays to study GPCR dimerization. This progress concerns, for example, optimization of complementation assays in stability and light output by a consensus sequence driven strategy. Several mutations applied in fluorescent and luminescent complementation assays in the field of GPCR dimerization will be discussed throughout the following sections of this review, for additional information concerning mutational optimization of split reporter genes in general, we refer to the excellent review by Wehr and Rossner [50].

2. Fluorescence-Based Complementation Assays

One of the well-studied PCAs is BiFC, which is based on the structural and functional complementation of two non-fluorescent protein fragments brought into close proximity by their interacting fusion partners (Figure 3). Upon complementation of the fluorescent protein, the fluorescence can be read via a plate-reader, imaged using a fluorescence microscope or even analyzed by flow cytometry without the need for any treatment of the cells.

In 1999, Regan and colleagues reported BiFC for the first time in *Escherichia coli* (*E. coli*) by using a strategy based on the non-covalent association of the split fragments of green fluorescent protein (GFP), fused to two antiparallel leucine zippers [51]. Subsequently, Hu et al. [52] reported on the development of a BiFC assay based on the split fragments of yellow fluorescent protein (YFP), tagged to two interacting transcription factors in living mammalian cells. Shyu et al. [53] developed Venus, one of the brightest fluorescent proteins, for studying PPIs at physiological conditions using BiFC. Since then, numerous types of fluorescence complementation-based assays have been developed to visualize signaling events involving two or even more interacting proteins, ranging from bacteria to mammalian cells.

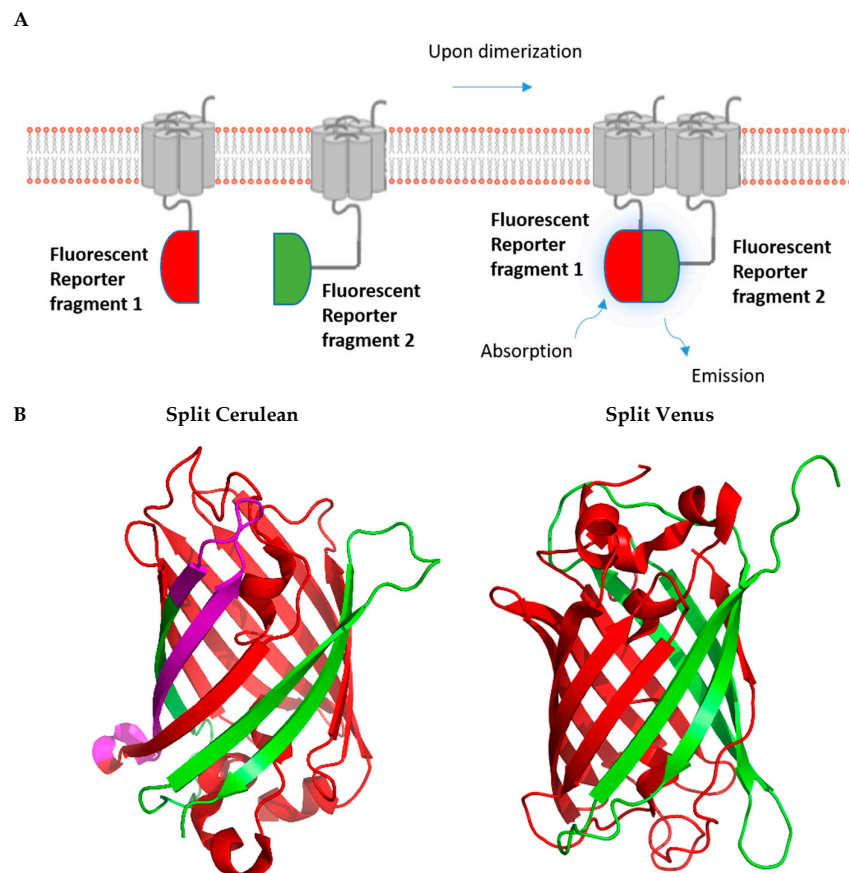


Figure 3. Fluorescence-based complementation assay. (A). The principle of the fluorescent complementation assay is shown schematically. (B). Split versions of the fluorescent reporters Cerulean and Venus are shown in green and red. Purple refers to overlapping sections. (PDB; Accession no. 3AKO for Venus; PDB: Accession no. 5OXB for Cerulean).

2.1. Fluorescent Proteins

2.1.1. Green Fluorescent Protein (GFP)

Regan and colleagues used GFP, the first fluorescent protein, for characterizing BiFC by tagging the split fragments of GFP, i.e., NGFP and CGFP, split at residue 157–158, to strongly interacting antiparallel leucine zippers. The split fragments did not self-assemble in bacteria, whereas upon fusion to antiparallel leucine zippers, the fragments could properly fold to form a functional GFP protein, and hence, gain fluorescence. A recent study performed by Son and colleagues [54] made use of enhanced GFP (EGFP), an improved version of GFP that is brighter and more photostable. The EGFP molecule was split to yield N-EGFP (1–158) and C-EGFP (159–238), respectively. These split parts were fused to α -factor receptor (Ste2p), to study homodimerization in *S. cerevisiae*. Both full length and C-terminally truncated forms of the receptors led to fluorescence, indicating dimerization.

A recent advance in BiFC based on GFP is Tripartite GFP, introduced by Cabantous et al. [55] in 2013 in an attempt to cope with the problem of poor folding and self-assembly seen with split GFP. Tripartite GFP is based on the tripartite association of two twenty amino-acid residue fragments of GFP, coined as GFP 10 and GFP 11, that can be fused to the proteins under investigation, and the complementary GFP 1–9 detector. Upon interaction between the protein partners, GFP 10 and GFP 11 associate with GFP 1–9 to form the functional GFP molecule. This system was initially characterized in *E. coli* and in mammalian cells, where the rapamycin induced association of FRB (FKBP-rapamycin binding domain of mTOR) and FKBP12 (FK506- and rapamycin-binding protein) was studied. This system has extensively been exploited for the validation of multiple reporter systems. Although use of this Tripartite GFP system has not been reported for studying GPCR dimerization, it could be an advantageous platform to be used over the conventional split GFP.

2.1.2. Yellow Fluorescent Protein (YFP) Variants

In 2001, Nagai et al. [56] first demonstrated the use of the yellow fluorescent protein (YFP) in BiFC, wherein they used a variant of YFP to study the effect of Ca^{2+} on PPI in living cells. Shortly after, Griesbeck and colleagues [52,57] tried several permutations and combinations of split regions in EYFP (S65G, S72A, T203Y) to improve the fluorescence signal along with minimal self-assembly of the split fragments. Several variants of YFP have been developed such as Citrine (a variant having a Q69M mutation) and Venus (a variant carrying a F46L mutation), that folds even at a physiological temperature of 37 °C [58]. This led to the usage of these fluorescent protein fragments to study GPCR–GPCR interactions.

Vidi et al. (2008) [59,60] used split fragments of Venus fused to the adenosine $\text{A}_{2\text{A}}$ receptor ($\text{A}_{2\text{A}}$) to address the question about $\text{A}_{2\text{A}}$ oligomerization in the CAD neuronal cell line. Co-expression of $\text{A}_{2\text{A}}$ -VN and $\text{A}_{2\text{A}}$ -VC in these cells, as well as in HEK293 cells, resulted in fluorescence, indicating $\text{A}_{2\text{A}}$ dimerization, in contrast to negative controls wherein cells transfected with $\text{A}_{2\text{A}}$ -VN/D1R-VC/M4R-VC did not show any fluorescence, pointing at the specificity of the interaction.

Similarly, Kilpatrick et al. (2014) [61] assessed dimerization of a receptor for NPY, Neuropeptide Y Y1/Y5. This neuropeptide is a widely expressed modulator of the central nervous system, mainly known for its role in the regulation of appetite upon its release from hypothalamic arcuate neurons. The receptors, tagged with split YFP, showed heterodimerization and the irreversible nature of these fragment tags also helped to elucidate the cellular localization of these dimers. Evaluation of the effect of ligands on the internalization of dimers, as opposed to individual

protomers, demonstrated that Y1/Y5 receptor dimers displayed altered ligand pharmacology, indicative of allosteric interaction.

Przybyla and Watts [62] in 2010 used Venus-based BiFC to demonstrate heterodimerization between the cannabinoid receptor 1 (CB₁) and the dopamine D_{2L} receptor (D_{2L}R) in CAD cells, wherein CB₁ and D_{2L}R receptors were fused at their C-terminus to VN and the C-terminal fragment of Cerulean (CC), respectively (see Section 2.2.1).

Ang et al. [63] in 2018 assessed homo- and heterodimerization of free fatty acid receptors, GPCRs that are expressed on mammalian cells to sense the short-chain fatty acids derived from the microbiota. To study dimerization, HEK293T cells were transfected with VN- or VC-tagged (free fatty acid receptor 2/3) FFAR2/3 receptor. Homodimerization of FFAR3 and FFAR2, as well as heterodimerization between FFAR2-VN and FFAR3-VC were observed, as compared to negative control FFAR2/3 tagged to VN or VC plus P2RY1-VN or P2RY1-VC (for heterodimerization) and P2RY1-VN plus P2RY1-VC (for homodimerization).

Xue et al. [64] in 2018 used a BiFC assay based on split Venus to detect heterodimerization between growth hormone secretagogue receptor 1 α (GHSR1 α) and Orexin 1 receptor (OX1R). Recently, Navarro et al. [65] (2018) used interfering peptides corresponding to regions in TM1–7 domain for A_{2A} and TM5–7 domain for A₁, each fused to TAT peptides, to disrupt the interaction between the A₁–A_{2A} heteromer. They therefore expressed receptors tagged with split fragments of YFP (nYFP and cYFP) in HEK293T cells. In the absence of the TM peptides, fluorescence complementation was detected for A_{2A}-nYFP and A_{2A}-cYFP in the presence of untagged A₁, suggesting homodimerization of A_{2A}. The interaction sites were mapped to TM4 and 5 of A_{2A} with the help of interfering peptides. In the same way, a fluorescence signal corresponding to the heterodimer (A₁-nYFP and A_{2A}-cYFP) was also detected. Upon addition of TM peptides corresponding to regions 5/6, a reduction in fluorescence was observed, suggesting the involvement of these regions in the formation of the heterodimer. In the same year (2018), BiFC using split fragments of YFP was utilized by Hinz et al. [66] for studying homodimerization between A_{2A} and heterodimerization with A_{2B} receptor. For this purpose, CHO cells were transfected with a constant amount of A_{2A}-NYFP and increasing amounts of HA-A_{2A}-CYFP and a strong fluorescence signal was observed, indicative of the A_{2A}-A_{2A} homodimer. For the purpose of studying A_{2B}-A_{2A} heterodimerization, the GPCRs were tagged with NYFP and CYFP, respectively. CHO cells transfected with constant amounts of A_{2B}-NYFP and increasing amounts of HA-A_{2A}-CYFP showed significantly higher fluorescence when compared to the negative control comprising of GABAB2-NYFP and increasing amounts of HA-A_{2A}-CYFP.

An interesting study was performed by Song et al. [67] in 2019, where a homo-molecular fluorescence complementation (homo-FC) probe was designed by splitting a single fluorescent protein (superfolder GFP, sfGFP) into two fragments and linking the C-terminal fragment (strands 8–11) with a short linker N-terminally to the N-terminal fragment (strands 1–7). This “flopped fusion” construct, which does not show intramolecular self-complementation, was used to tag the β_2 adrenergic receptor (β_2 -AR), to study its homodimerization. A fluorescent signal was generated by complementation of strands 1–7 of one “flopped” sfGFP molecule, fused to one β_2 -AR, with strands 8–11 of another “flopped” sfGFP molecule, fused to a second β_2 -AR. The advantage of this approach lies in the fact that only one receptor construct needs to be generated, to allow the assessment of homodimerization.

2.1.3. Cyan Fluorescent Protein (CFP)

Apart from using the aforementioned split Venus, Vidi et al. [59,68] also used split fragments of Cerulean to show that A_{2A} homodimerizes in CAD cells (see Section 2.2.1). The same fluorescent protein fragments were also used in 2010 by Przybyla and Watts [62] to illustrate the homodimerization of D_{2L}R.

2.1.4. Red, Far-Red and Near-Infrared Fluorescent Proteins

Red fluorescent proteins encompass a couple of fluorescent proteins that emit in the red region of the visible spectrum. They have evolved through time to overcome issues such as slow maturation or the tendency to oligomerize (as exhibited by Discosoma Red (DsRed)) [69]. Red variants such as monomer red fluorescent protein 1 (mRFP1) are not appropriate for BiFC due to their low extinction coefficient, quantum yield, and low photostability, despite faster maturation than DsRed and better tissue penetration. Although several variants of RFP have been used for FRET, BiFC assays with the split parts have not been reported for demonstrating GPCR–GPCR interactions.

For imaging deep tissues in animals, fluorescent proteins with emission in the far-red region of the spectrum are chosen. RFP, mcherry and DsRed, have limited application in BiFC due to the requirement of low temperatures for their maturation [70,71]. This led to the development of mKate, a far-red, monomeric form of Katushka fluorescent protein, that matures faster and is photostable [72]. BiFC based on split fragments of mKate fused to transcription factors was successfully applied in COS-7 cells [73], but has not been applied yet for studying GPCR oligomerization.

2.2. BiFC Assays

2.2.1. MBiFC

Multicolor BiFC (MBiFC) has gained a lot of attention in recent years as it essentially gives the freedom to explore multiple protein–protein interactions inside the same living cell [74]. It has been studied with split parts of Venus and Cerulean, wherein the C terminal part of Cerulean is tagged to protein A (CC), the N-terminal part of Venus is tagged to protein B (VN), and finally, the N-terminal fragment of Cerulean is fused to protein C (CN). The C-terminal fragment of Cerulean (CC), comprising residues 155–238 (C155), can functionally complement with the N-terminal fragment of VN and CN, comprising residues 1–172 (N173), to produce a Venus (VN+CC) signal or Cerulean (CN+CC) signal. Since these parts are non-fluorescent per se, simple employment of different excitation and emission wavelengths allows to study simultaneously the interaction between A–B and A–C within the same cell [59,75,76].

Multicolor BiFC has aided in the understanding of GPCR oligomerization amongst all its other applications. Vidi et al. [59] (2008) used MBiFC to visualize the hetero- and homodimerization capacity of A_{2A} in the CAD neuronal cell line. To this end, D_{2L}R-VN was co-transfected with A_{2A}-CN and A_{2A}-CC in CAD cells. The detection of Venus (D_{2L}R-VN/A_{2A}-CC) and Cerulean (A_{2A}-CN/A_{2A}-CC) fluorescence signals revealed the presence of the A_{2A}–D_{2L}R heteromer and A_{2A} homomer, respectively. In addition, fluorescence microscopy also revealed that the dimers were localized at the plasma membrane as well as in intracellular compartments.

Przybyla and Watts [62] in 2010 showed oligomerization of CB₁ and D_{2L}R. To achieve this, they used CAD cells that expressed CB₁-VN and D_{2L}R fused to the split fragments of Cerulean (D_{2L}R-CN, D_{2L}R-CC) and these cells were imaged using fluorescence microscopy. The heterodimers in the cell formed by CB₁-VN–D_{2L}R-CC produced a Venus signal that was mainly localized in the intracellular compartment, whereas the homodimers of D_{2L}R produced a Cerulean signal that was localized both in the intracellular compartment and on the plasma membrane.

2.2.2. BiFC-RET

The fluorescent proteins used in BiFC could also be coupled to FRET or BRET to study GPCR oligomerization, e.g., functionally complemented YFP could act as a FRET/BRET acceptor. BiFC-BRET [77] employs the energy transfer between GPCR-A fused to a luciferase and the complemented fragments of YFP that are fused to GPCR-B and GPCR-C that are as close as 10 nm apart, indicating a trivalent complex of GPCRs.

Navarro et al. [78] in 2008 used BiFC-BRET to demonstrate the interaction between D₂R, A_{2A}, and CB₁ receptors in HEK cells. The split fragments of YFP fused to CB₁ and A_{2A} acted as the

acceptor fluorophore in BRET, and *Renilla* luciferase (RLuc) fused to D₂R (D₂R-RLuc) acting as the donor.

BiFC-FRET [79] is similar, offering the additional benefit that the complex can be visualized within the cell. BiFC-FRET was used to understand the existence of higher order oligomers in the case of A_{2A} by using fusion constructs with Cerulean and split parts of Venus. Prior to this, to test if the system could record FRET, A_{2A} dimerization was tested by fusing A_{2A} to Cerulean and Venus. Based on this, CAD neuronal cells were co-transfected with A_{2A}-Cerulean and A_{2A} fused to split parts of Venus. FRET signals were detected, suggesting the existence of A_{2A} oligomers, as opposed to a negative control wherein cells transfected with M₄ as the acceptor only showed very low FRET [60].

In 2017, Bagher et al. [80] combined BRET with BiFC which forms the basis of “sequential resonance energy transfer 2” (SRET²). Herein, the donor for BRET was the *Renilla* luciferase (RLuc) fused to D_{2L}R (D_{2L}R-RLuc), which upon addition of the substrate coelenterazine 400a, excited the acceptor GFP² fused to D_{2L}R (D_{2L}R-GFP²), thereby confirming the homodimerization of D_{2L}R. Sequentially, the GFP² acted as FRET donor for the CB₁-VN/CB₁-VC homodimer, thus providing evidence for heterotetramerization. With increasing concentrations of CB₁-VN/CB₁-VC, there was a hyperbolic increase in net SRET², as opposed to negative control cells that were transfected with D_{2L}R-RLuc and mGLuR6-GFP² and increasing concentrations of CB₁-VN/CB₁-VC, confirming that the homodimers of D_{2L}R and CB₁ associate to form oligomers.

2.3. Ligand-Dependent Modulation of Dimerization

The level of GPCR dimer formation could possibly be altered by ligands interacting with the receptor(s). The potential of split-protein sensors to monitor the dynamic changes in dimerization, provoked by chemical inducers or inhibitors, was demonstrated using the rapamycin-dependent FKBP/FRB Chemically Induced Dimerization (CID) system [81]. The first demonstration for the use of MBiFC to monitor drug-modulated GPCR oligomerization was published by Vidi et al. [59] in 2008 for the A_{2A}-D₂R interaction. Stimulation of D₂R with the D₂R agonist quinpirole led to internalization of D₂R homodimers as well as of A_{2A}-D₂R oligomers, which was blocked by the D₂R antagonist, sulpiride [59]. Furthermore, treatment with quinpirole decreased the formation of A_{2A}-D₂R heterodimers, as compared to A_{2A} homodimers. Also, treatment with 5-N-methylcarboxamidoadenosine (MECA), an agonist for the adenosine receptor, increased the proportion of A_{2A}/D_{2L}R heterodimers, as compared to homodimers of A_{2A}, while treatment with its antagonist (CGS15943) had no effect on BiFC. Interestingly, these drug-induced alterations on the formation of oligomers could not be supported by changes in the receptor density.

As discussed above, Przybyla and Watts [62] demonstrated heterodimerization of CB₁-D_{2L}R using MBiFC. Also, the influence of ligands on the balance between hetero- and homodimers was analyzed. Quinpirole and CP55,940 shifted the balance towards the formation of CB₁-D_{2L}R heterodimers, as opposed to D_{2L}R-D_{2L}R homodimers. The D₂R antagonist sulpiride favored D_{2L}R homodimerization. Consequently, it was shown that sustained treatment with quinpirole influenced the expression of D_{2L}R, which in turn might affect the dimerization status inside the cell as: D_{2L}R-CB₁ > D_{2L}R-D_{2L}R > A_{2A}-D_{2L}R. In a similar way, studies also indicated that a CB₁ antagonist (1-(7-(2-chlorophenyl)-8-(4-chlorophenyl)-2-methylpyrazolo (1,5-a)-(1,3,5)triazin-4-yl)-3-ethylaminoazetidino-3-carboxylic acid amide benzenesulfonate) could have a potential therapeutic role in Parkinsonism by enhancing the activity of L-DOPA [82].

In 2018, Navarro et al. [65] evaluated if activation of A₁ or A_{2A} alone or by a combination of agonists could modulate the TM interface of the heteromer formed by A_{2A}-nYFP and A₁-cYFP. Selective agonists for A₁, N⁶-cyclopentyladenosine (CPA), and A_{2A}, 4-(2-((6-Amino-9-(N-ethyl-β-D-ribofuranuronamidosyl)-9H-purin-2-yl) amino) ethyl) benzenepropanoic acid (CGS-21680), were used alone or in combination. None of the agonists, either used in combination or alone, could modify the TM interface of the A₁-A_{2A} heterodimer.

Various ligands were tested for their capacity to affect β₂AR oligomerization by Song et al. [67] in 2019. A panel of compounds consisting of five agonists, two antagonists, and two inverse agonists

were screened using the abovementioned homo-FC assay. They observed that agonists induced oligomerization, while inverse agonists reduced the signals, and a combination of both also resulted in reduced intensity of signal, thus confirming the same.

3. Luminescence-Based Complementation Assays

Luciferases, which have seen an expansive growth in use as reporter proteins in biological research, are attractive due to the high signal-to-background ratio associated with their usage, as no excitation light is required to generate a signal [83]. Similar to fluorescent proteins, these enzymes can also be used in applications where the luminescent protein itself is split into two fragments, which are conjugated to proteins of interest (Figure 4). Several luminescent proteins have served this purpose, with Firefly luciferase (FLuc) [84,85] and *Renilla* luciferase (RLuc) [86,87] being the two most commonly used bioluminescent proteins. However, many other novel luciferases have also been developed with favorable characteristics in terms of stability, substrate requirement, brightness, and emission spectrum, e.g., *Gaussia princeps* luciferase (GLuc) [88] and NanoLuciferase [89–91]. All mentioned split-protein reporters were initially validated using the aforementioned FKBP/FRB CID system [84,88,92–94].

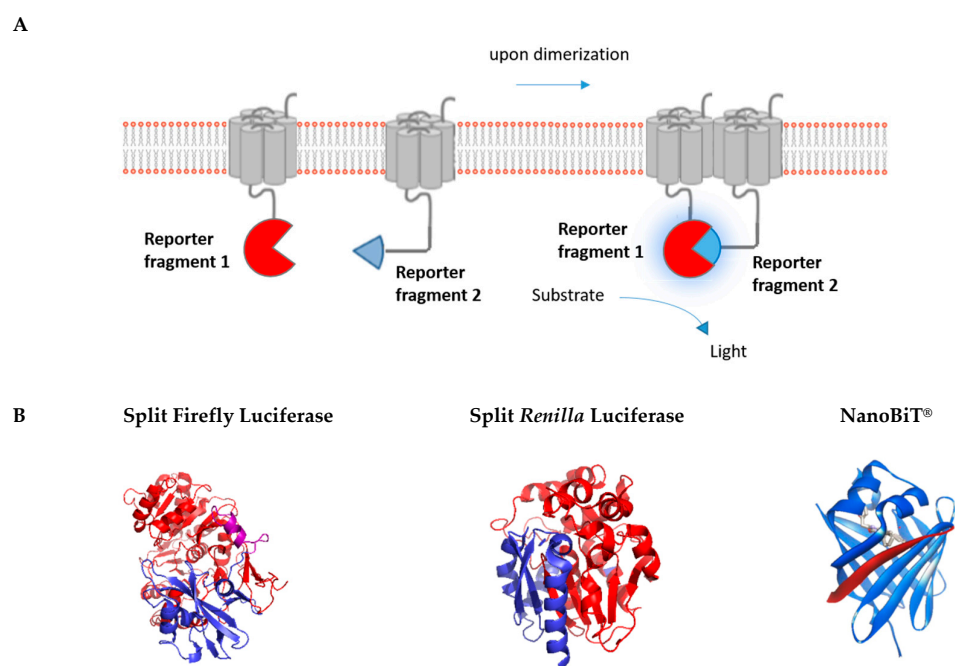


Figure 4. Luminescence-based protein complementation assay. **(A)** The principle of the luminescent complementation assay is shown schematically. **(B)** Split versions of the luminescent reporters Firefly, *Renilla*, and Nanoluciferase are shown in blue and red. Purple refers to overlapping sections. Split luminescent biosensors are depicted in proportion to their size (Fluc: 62 kDa, RLuc: 36 kDa, and NanoBiT®: 19 kDa) (PDB Accession no. 1LCI for Firefly luciferase, PDB Accession no. 2PSD for *Renilla* luciferase) (NanoBiT®: source Promega).

3.1. *Renilla*/Firefly Luciferase

FLuc is the least optimal reporter for employment as a bioluminescent tag, due to its size (± 60 kDa) and its dependence on ATP, molecular oxygen, and magnesium for activity. Luciferases that use coelenterazine as a substrate, such as the luciferase from the sea pansy *Renilla reniformis* RLuc [95,96], have an advantage over FLuc in that they are not ATP dependent and only require the presence of molecular oxygen for the enzyme-catalyzed conversion of a substrate to a luminescent reaction product. As the bioluminescent activity of formerly used non-truncated luciferases is a

limiting factor to permit their application in high throughput screenings, specific mutations of certain luciferases, like in *RLuc*, were selected using a consensus sequence-driven strategy and screened for their ability to confer stability in activity as well as for their light output [97].

For GPCR dimerization purposes, the split-luciferase reporter approaches have not been applied to the same extent as their fluorescence counterparts. Nevertheless, these reporter systems also allow the rapid detection of macromolecular GPCR–GPCR interactions. For instance, the split *FLuc* methodology was used in cell culture models as well as in tumor xenograft models of breast cancer (see Section 8) to measure changes in chemokine receptor *CXCR4* (and *CXCR7*) homodimerization in response to pharmacological agents [98]. Likewise, the *CXCR4* homodimer has also been demonstrated with a split *RLuc* assay [99]. For this *Renilla* luciferase complementation assay, *RlucII* was implemented, which is derived from *RLuc*, wherein two mutations (C124A and M185V) were introduced to make it brighter.

For studying dopamine *D₂R* oligomers, an alternative *RLuc* construct was implemented, namely, *Rluc8* [100]. Owing to eight favorable mutations, this *RLuc8* has a 4-fold improved light output, compared with the parental enzyme [97]. More recently, the same complementation assay was also used by Casado–Anguera et al. [101] (2016) to demonstrate the existence of a therapeutically relevant GPCR dimer, namely, the *A_{2A}–D₂R* dimer. This dimer and its inter-protomer allosteric mechanisms have been proposed as a new model that could contribute to our knowledge concerning drug dosage for the treatment of Parkinson’s disease.

3.2. NanoLuciferase

One of the most recent and undoubtedly the most optimized luminescent PCA so far is “NanoLuciferase Binary Technology” or the NanoBiT® system, developed by Promega. The assay is based on complementation of the split fragments of NanoLuciferase or NanoLuc. This engineered luciferase reporter is a small (19 kDa), ATP-independent luminescent protein, originating from a luminous deep-sea shrimp, *Oplophorus gracilirostris* [102]. In combination with the development of a novel cell-permeable imidazopyrazinone substrate, furimazine, this bioluminescence system generates a glow-type luminescent signal that is over 150-fold greater compared to that of the former *Renilla* and Firefly luciferases. NanoLuc exhibits high physical stability, in a wide range of environmental conditions such as temperature, pH, urea, and ionic strength [103]. As opposed to luminescent protein fragments described earlier, the NanoBiT® subunits do not consist of two fragments similar in size but correspond to a small 1.3 kDa subunit (Small BiT; SmBiT) and an 18 kDa subunit (Large BiT; LgBiT). Both subunits have been thoroughly characterized and have a low affinity ($K_D = 190 \mu\text{M}$) for one another, thus providing the ability to follow kinetics of PPIs.

Although the NanoBiT® PCA has only recently been developed [92], its broad applicability has proven successful in numerous research fields. NanoBiT® has, for instance, been applied to develop bio-assays based on the recruitment of β -arrestins [104–110] and G-proteins [111] for the detection or activity profiling of certain compounds or to elucidate the molecular interaction between the transducer and a GPCR. For GPCR dimerization purposes, NanoBiT® has been implemented to detect ligand-dependent modulation of *D_{2L}R–D_{2L}R* homodimers [40,112]. The methodology has also been used to perform live-cell monitoring of the dynamics of the interaction between the Melanocortin 4 receptor (*MC4R*), and the Melanocortin 2 receptor accessory protein 2 (*MRAP2*), as well as of *MC4R* homodimerization [113].

The NanoBiT® system has been further optimized into a tri-part protein fragment complementation assay by two independent groups [114,115]. For this purpose, the LgBiT was dissected into a smaller C-terminal part of 11 amino acids (LcBiT) and an N-terminal derivative (LnBiT), serving as a “detector” protein of 16.5 kDa, similar to the Tripartite GFP mentioned in Section 2.1.1. Initially, this assay was designed to facilitate the purification of fusion proteins to procure significant quantities and to avoid the lapse of detection of the PPI of interest due to the steric hindrance of LgBiT, due to its higher molecular weight. In addition, another novel split-luciferase reporter based on NanoBiT® was developed that implements GFP- and mCherry-recognizing nanobodies fused to LgBiT and SmBiT [116]. Using this experimental set-up,

GFP- or mCherry-tagged dimers or higher order oligomers can be detected. As fluorescently tagged proteins are already often in use, this luminescent PCA can be directly used and consequently bypasses the tedious protein re-cloning to explore the different possible configurations for setting up a PCA (also see Section 6). Furthermore, due to the strong affinity of the GFP- and mCherry-recognizing nanobodies, a low limit of detection of PPIs or protein aggregations, from micromolar up to low nanomolar, can be achieved.

3.3. BiLC-RET

The application of BiLC has already shown its efficiency in combination with a BRET assay, allowing to detect ternary protein complexes. Sahlholm et al. [117] (2018) demonstrated, by combining BiLC of D₂R–A_{2a} heterodimers with a YFP-tagged β -arrestin, a BRET signal could be obtained. Furthermore, it was postulated that D₂R agonists quinpirole or UNC9994 require the formation of D₂R–A_{2a}R heterodimers, to promote β -arrestin2 recruitment.

4. Combinatorial Assays: BiFC and BiLC

For the detection of higher order oligomerization of GPCRs, a combinatorial assay with both BiFC and BiLC can be implemented (Table 2). Rebois et al. [77] (2008) introduced the concept of detection of a tetravalent complex using a combination of BiFC and BiLC, wherein four β_2 ARs were fused to split parts of Venus or GLuc and the homodimerized β_2 AR-GLuc acted as the donor for the homodimerized β_2 AR-Venus acceptor, the presence of a BRET signal indicating the presence of the tetramer. In the same year (2008), Guo and colleagues [100] showed the existence of a D₂R homo-oligomer, using the same technique.

A complemented donor–acceptor resonance energy transfer (CODA-RET) assay showed that the A_{2A}–D₂R dimer not only forms dimers but can also assemble to form a heterotetramer, composed of two receptor homodimers [118,101]. For this experimental setup, a complemented YFP was implemented and combined with Rluc8.

Table 2. Combinatorial assays with BiFC and BiLC for GPCR oligomerization purposes.

GPCR Dimer	Oligomeric Type	PCA Type	Split Biosensor	Fragments	Negative Control	Cell-Type	Year	Ref.
CXCR4–CXCR4/CC2	Hetero-oligomer with CC2, homotetramer	BiLC and BiFC	RLucII, vYFP	RLucII: 1–330, 331–936 vYFP: 1–465, 466–720	D ₂ R	HEK293	2014	[99]
A _{2a} –D ₂ R	heterotetramer	BiLC and BiFC	RLuc8, YFP	RLuc8: 1–229, 230–311 YFP: 1–155, 156–238	A ₁ R, D ₁ R	HEK293	2015/2016	[118,101]
D _{2s} R–D _{2s} R	Homo-oligomer	BiLC and BiFC	RLuc8, mVenus	RLuc8: 1–229, 230–311 mVenus: 1–155, 156–240	CD8, TSHr	HEK293T	2008	[100]
β ₂ AR–β ₂ AR	homotetramer	BiLC and BiFC	GLuc, Venus	GLuc: 1–63, 64–185	GLucN, VN, VC	HEK293	2008	[77]

5. Comparison of Split Protein Approaches

Protein complementation techniques have shown their applicability in many fields to unravel PPIs. However, every technique has its advantages and limitations, as summarized in Table 3 for BiFC and BiLC.

5.1. Advantages of PCA

A key advantage of BiFC and BiLC lies in the fact that complementation assays offer a robust and straightforward approach to evaluate PPIs in living cells. Consequently, their applicability in high-throughput experiments is a valuable quality. The development of MBiFC, BiFC-RET, or CODA-RET has made the detection of multiple protein–protein interactions feasible, such as GPCR oligomers, without the need for cell lysis or fixation. Due to the simplicity of performing these assays, they can be used as a screening platform for drugs [59]. When used with a plate reader, BiFC and BiLC also offer quantitative and rapid results on relatively large cell populations.

The most valuable aspect of BiLC is that kinetic measurements can be performed due to its sensitivity and the limited propensity of self-association events, at least for certain BiLC reporters (e.g., NanoBiT®). For the purpose of studying GPCR dimerization, the NanoBiT® assay has been compared to related bioluminescence and fluorescence PCAs, including split Venus and RLuc [112]. As a benchmark, the generally accepted D₂R–D₂R dimer was compared in the different PCAs. The NanoBiT® assay presented the best signal-to-noise ratio and was considered the most optimal candidate assay for targeting GPCR dimers. In addition, this method can also be implemented to analyze the kinetics of ligand-dependent modulation of dimerization, broadening its application potential [40]. Even high-throughput screenings can be performed, which is highly relevant, given the growing interest and effort to identify small molecule drugs that can target disease-relevant dimers (or even selectively alter GPCR dimer function). Moreover, BiLC has also proven its applicability in an *in vivo* setup [85,87].

On the other hand, BiFC has the capacity to visualize and localize the interaction, also offering insight into potential protein aggregation artefacts inside the cell. BiFC does not require any exogenous stains or substrate. Thus, depending on the intrinsic fluorescence, it allows the direct measurement of the interaction between proteins.

5.2. Limitations of PCA

BiFC and BiLC also come with their own set of drawbacks. First of all, both techniques require the fusion of the target proteins with split fragments of a fluorescent or luminescent protein, which could affect the original proteins' interaction dynamics. Moreover, the amount of proteins expressed by the cell may potentially cause false positive readouts. Therefore, a good set of negative controls is obligatory (see Section 6.3). In addition, one should bear in mind that these techniques, though indicating the close proximity of the proteins being studied, are not an actual assessment of the physical contact between these proteins.

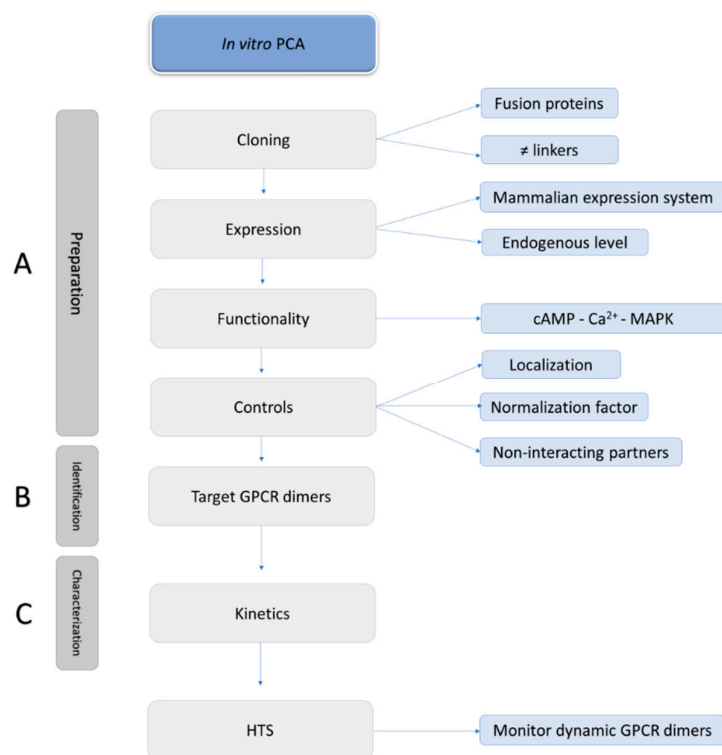
BiFC needs a longer maturation time, has an irreversible nature of complementation, shows a high degree of self-assembly of the split fragments of the fluorescent protein, and requires molecular oxygen for the maturation of the fluorophore, thereby making it unsuitable for obligate anaerobes. Since maturation is a time-dependent process, the data procured from BiFC are not in real time. The irreversibility of the technique also makes it more difficult to study the dynamics of interactions inside the cell and the influence of drugs on the interaction. Despite these limitations, BiFC is still the only protein complementation technique that gives information about the location of PPIs in living cells. For BiLC, its limitations are the need of a substrate and the rather limited possibility to detect the localization of the interaction.

Table 3. Advantages and disadvantages of BiFC and BiLC techniques.

BiFC	
Advantages	Disadvantages
Straightforward technique	Need for tagged proteins (GPCRs)
High-throughput experiments	Autofluorescence
Imaging microscopy: localization of the interaction	Photobleaching
Study intact cells	Measuring dynamics: limited (maturation time)
	Not applicable for studying inhibition of interactions
BiLC	
Advantages	Disadvantages
Straightforward technique	Need for tagged proteins (GPCRs)
High-throughput experiments	Requires a substrate
Kinetic measurements	Detection of localization: limited
Study intact cells	
In vivo application	

6. Guidelines to Perform Accurate PCA-Based Assays

A potential but inherent drawback of PCA-based assays is the requirement to fuse fragments derived from fluorescent or luminescent proteins to the GPCR, which inherently may alter the GPCR's function or cellular localization. Therefore, a list of recommendations for proper implementation of PCA-based assays can be put forward: 1) examine all possible combinations of fusion proteins; 2) verify the functionality and localization of the fused proteins; 3) include proper controls, i.e., non-interacting partners, to control for self-assembly; 4) include a normalization reporter, to allow compensation for differences in transfection efficiencies; 5) preferentially transfect an amount of DNA near to endogenous expression levels of GPCRs, to avoid false positives due to the random collisions or potentially generate stable cell lines; and 6) follow the kinetics of the GPCR interaction (Figure 5).

**Figure 5.** A guideline for the set-up of an optimized PCA.

6.1. Possible Fusions

The development of PCAs requires multiple rounds of cloning, as one needs to design the most optimal “split” system for each GPCR of interest. Several configurations should be explored, being (i) different reporter systems, (ii) (N- or) C-terminally tagged proteins, and (iii) diverse linkers (i.e., length and composition). While this is a time-consuming step, having a certain workload, it is of paramount importance for the development of a robust, sensitive, and successful complementation assay.

6.1.1. Selection of the Reporter System

A broad variety of fluorescent and luminescent split reporter fragments is available to unravel PPIs. For studying GPCR dimerization, different split reporter systems might be advised, depending on the characteristics of the GPCR dimer of interest (if known). For obligate GPCR dimers, the rules could be less stringent for the choice of the reporter system to be used, whereas for transient GPCR dimers, a dynamic reporter assay with kinetic measurements close to real-time is required. When BiFC is preferred, the truncated YFP fragments YN155 and YC155 are recommended, due to the high signal-to-noise ratio [119]. Venus, a brighter version of GFP, split at 155 or 173, is also commonly used for studying protein interactions. The benefit of using Venus fragments over YFP is that the complementation can be read at 37 °C, thus avoiding the incubation at 30 °C which is necessary when dealing with split-YFP.

When BiLC is preferred, the NanoBiT® system offers advantages in flexibility and in the low affinity of the SmBiT–LgBiT fragments for each other, so dynamic interactions can be monitored in living cells in real time [92].

6.1.2. (N- or) C-Terminally Tagged GPCRs

The site at which the tag is fused to a GPCR plays an important role in the functionality of the GPCR under study and the purpose of the study (Figure 6). All the published research on PCAs with GPCR dimers involves fusion proteins where the split biosensors are tagged to the C-terminal end of GPCRs, so transport and expression at the membrane surface are not hampered. Especially for monitoring interactions of Class C GPCRs, like the metabotropic glutamate receptors, C-terminal tagging of biosensors is recommended, as modification of the long N-terminus of the receptor could lead to deviating outcomes.

On the other hand, split fragments tagged to the N-termini of GPCRs, though not yet implemented for research on GPCR dimerization, might be an interesting tool as well, as these fusion proteins may enable the simultaneous monitoring of GPCR-GPCR interactions and the signaling properties of the dimer. For example, when a split biosensor is fused to the N-terminus of the GPCRs to monitor the GPCR-GPCR interaction, another split biosensor fused to the C-terminus of the GPCR can be applied to monitor G-protein or β -arrestin recruitment.

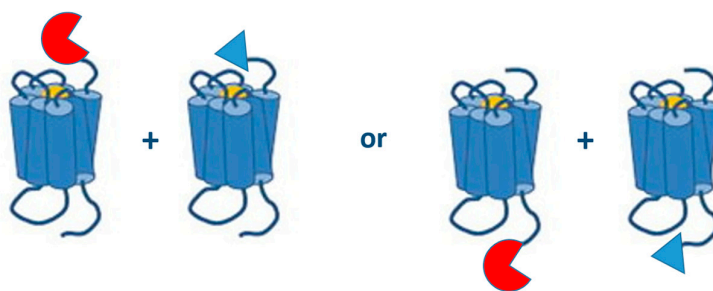


Figure 6. N- or C-terminally split biosensor GPCRs.

6.1.3. Linkers

GPCRs of interest and the fragments of the split luminescent or fluorescent partners are usually separated by a linker sequence to assure the flexibility required for the proper folding and “settling” of the biosensor. Multiple linkers with alternating sequences and conformations have been designed, the design of these linkers often being empirical but based on linker sequences derived from natural multi-domain proteins [120].

The most frequently applied linkers in PCAs are composed of non-polar glycine and polar serine amino acids (“GS” linker), developed by Argos [121]. The flexibility of these linkers originates from the incorporation of these small residues. Moreover, the use of serine (or threonine) residues in these linkers contributes to the stability in aqueous conditions by the formation of hydrogen bonds with water molecules [122]. A model linker design involves the (GGGS)_n template, where n represents the copy number. By adjusting this copy number, one might obtain the most optimal “GS” linker length to achieve appropriate separation of the GPCR and the split biosensor. For PCA purposes, often utilized GS linkers comprise roughly eight to fifteen amino acids, such as DGGSGGG [123], GGGSGGG [99,124] or GSSGGGGSGGGSSG [105]. Other -more rigid linkers have also been implemented in PCA, for example ATGLDLELKASNSAVDGTAGPVAT [117]. The proline residues in linkers might increase stiffness to keep fusion moieties at a distance, whereas lysine residues are often added to improve solubility [122]. Furthermore, proper codon usage (i.e., avoiding rare ones) is advised to not compromise expression levels. Also, depending on the expression system, sequences which could represent potential protease recognition sites should be avoided.

Overall, the selection of the linker length is crucial to allow proper folding and accomplish optimal biological activity of the fusion proteins. For fusion proteins, the linker length is an essential feature of the PCA setup since linkers that are too short often result in impaired biological activity due to the inability to accommodate the complementation of the two split proteins, whereas linkers that are too long can create false positives. Therefore, the linker length will vary on a case-by-case basis. With this design in mind, the NanoBiT® system offers different restriction sites in the provided biosensor templates, resulting in various GS linkers, ranging from 15 to 21 amino acids [112]. Some examples of linkers that have been applied in PCA assays to detect GPCR dimers are given in Table 4.

Table 4. Linker sequences applied in PCA assays for studying GPCR dimerization.

GPCR Dimer	PCA Type	Linker	Ref
AT ₁ -AT ₂	BiFC, Venus	GGGSGGGG	[125]
CXCR4-CXCR4	BiLC, Rluc	(GGGS) ₂	[99]
D ₂₁ R-D ₂₁ R	BiFC, Venus	LG	[100]
D ₂₁ R-D ₂₁ R A _{2a} -D ₂ R	BiLC, Rluc	ATGLDLELKASNSAVDGTAGPVAT	[117,112]
D ₂₁ R-D ₂₁ R	BiLC, NanoBiT	GNS-GSSGGGGSGGGSSG	[112]
MOP-NPFF2	BiFC, Venus	DGGSGGGG	[123]

6.2. Functionality and Localization of the Fusion Proteins

Occasionally, the fusion of biosensors to a protein of interest (POI) may interfere with the expression level, activity or function of the latter [126,127]. This may, for example, highly depend on the sequence and length of the linker between the POI and the biosensor (as discussed in Section 6.1.). To verify whether the functional integrity of the fusion POI is not affected, several techniques should ideally be implemented, such as ligand binding as well as functional assays (e.g., calcium signaling, cAMP, β -arrestin recruitment, MAPK activation, etc.). An immunofluorescence assay can also be used to confirm correct localization of the recombinant proteins at the plasma membrane.

6.3. Non-Interacting Partners

Studying GPCR-GPCR interactions by complementation assays is relatively straightforward. However, the importance of including appropriate controls (i.e., GPCRs that do not bind to the GPCR of interest) should not be underestimated. It is also essential that the split reporter proteins do

not spontaneously associate in the absence of binding partners (that normally drive the complementation). If so, a high number of false positives would be generated [112,128,129]. Ideally, when the site of interaction is known, a mutation at the specific interaction site of the receptor, is performed to evaluate whether this indeed disrupts the interaction. In those cases where the interaction site is unknown, a screening could be performed to elucidate the perfect non-interacting GPCR partner. Negative controls implemented in already published research articles are shown in Table 5 for the BiFC and in Table 6 for BiLC.

Table 5. An overview of GPCR–GPCR interactions, shown by fluorescence complementation assays.

GPCR Dimer	Oligomeric Type	PCA	Split Biosensor	Fragments	Negative Control	Cell-Type	Year	Ref.	In Vivo or Native Tissue Evidence	Ref.
mGluR5–D2R	Heterodimer	BiFC	YFP	1–155, 155–231	GABA _{B2}	HEK	2009	[130]	Yes	[130]
D2R–D2R	Homodimer	BiFC	YFP	1–155, 156–238	D1R	HEK	2015	[118]	Yes	[131,132]
A2A–D2L R	Heterodimer	MBiFC	Venus/Cerulean	1–172, 155–238	D1	CAD	2008	[59]	Yes	[130,133–137]
D2L R–CB1	Heterodimer	MBiFC	Venus/Cerulean	1–172, 155–238	M4	CAD	2010	[62]	Yes	[133]
D2L R–D2L R	Oligomer	MBiFC	Venus/Cerulean	1–172, 155–238	-	CAD	2010	[62]	-	-
D25R–D25R	Homodimer	BiFC	Venus	1–155, 156–240	CD8	HEK293T	2008	[100]	Yes	[131,132]
AT1–AT2	Homo- and heterodimer	BiFC	Venus	1–158,159–239	ATIP	HEK293FT	2011	[125]	Yes	[138–140]
CXCR4–CXCR4	Homodimer	BiFC	vYFP	1–465, 466–720	D2R	HEK293	2014	[99]	-	-
A2A–A2A	Homodimer	MBiFC	Venus/Cerulean	1–172, 155–238	-	CAD	2008	[59]	-	-
A2A–A2A	Homodimer	BiFC	YFP	1–155, 155–238	Non-fused A1 (competition)	HEK293T	2018	[65]	-	-
A2A–A1	Heterodimer	BiFC	YFP	1–155, 155–238	-	HEK293T	2018	[65]	Yes	[141]
GHSR1a–OX1R	Heterodimer	BiFC	Venus	1–172, 156–239	-	HEK293T	2018	[64]	-	-
β2AR–β2AR	Oligomer	BiFC	-15sfGFP	-	-	HeLa	2019	[67]	-	-
A2A–A2A	Homodimer	BiFC	YFP	1–155,156–239	GABA _{B2}	CHO	2018	[66]	-	-
A2B–A2A	Heterodimer	BiFC	YFP	1–155,156–239	GABA _{B2}	CHO	2018	[66]	-	-
FFAR3–FFAR3	Homodimer	BiFC	Venus	1–155 (I152L), 155–239	P2RY1	HEK293T	2018	[63]	-	-
FFAR2–FFAR3	Heterodimer	BiFC	Venus	1–155 (I152L), 155–238	P2RY1	HEK293T	2018	[63]	-	-
mGluR2–mGluR2	Homodimer	BiFC	mCitrine	1–172, 155–238	-	HEK293T	2016	[142]	-	-
α1b–α1b	Homodimer	BiFC	eYFP	1–172, 155–238	-	HEK293T	2007	[143]	-	-

Table 6. An overview of GPCR–GPCR interactions, shown by luminescence complementation assays.

GPCR Dimer	Oligomeric Type	PCA Type	Split Biosensor	Fragments	Negative Control	Cell-Type	Year	Ref.	In vivo or Native Tissue Evidence	Ref.
CXCR4–CXCR4/CC2	Homodimer									
CXCR7–CXCR7	Homodimer	BiLC	FLuc	NLuc-416 and CLuc-398	β2-AR	HEK293T	2009	[98]	Yes	[144–146]
CXCR4–CXCR7	Heterodimer									
A2a–D2R	heterodimer	BiLC	RLuc8	1–229, 230–311	-	HEK293T	2018	[117]	Yes	[118,130,134–

										137]
A _{2a} -A _{2a}	homodimer	BiLC	RLuc8	1-229, 230-311	A ₁	HEK293	2016	[118,101]	-	-
D _{2L} R-D _{2L} R	homodimer	NanoBiT	NanoLuc	1-11, 12-167	CB ₂	HEK293T	2018/2019	[40,112]	Yes	[131,132]

6.4. Normalization Factor

To allow better comparison between different experimental setups, it is advisable to include a normalization factor, i.e., a reporter gene such as a fluorescent marker, as previously described [112,147]. In a transient expression system, co-transfection of this fluorescent marker, which should not cause cross-excitation, is recommended. A list of PCAs with corresponding fluorescent markers is given in Table 7. For BiLC in principle, all fluorescent reporters can be used, although it is recommended to avoid cross-excitation by choosing an excitation wavelength more than 20 nm before or after the emission wavelength of the luminescent protein. A potential strategy to cope with alterations in expression levels of the fusion proteins, is the implementation of ratiometric expression systems, consisting of two GPCRs of interest and a fluorescent marker for optical expression control. The so-called ratiometric BiFC (rBiFC) makes it feasible to control the expression through FACS or flow cytometry analysis [147].

Table 7. The application of normalization factors. * Given as examples.

Fluorescence	PCA (Excitation/Emission)	Fluorescent Marker (Excitation/Emission)
	Venus (515/528)	CFP (433/475)
	mCherry (587/610)	CFP (433/475)/GFP (488/510)/YFP (514/527)
	GFP (488/510)	CFP (433/475)
	Cerulean (452/478)	mTagBFP (402/457)
Luminescence	PCA	Fluorescence Marker
	Rluc	CFP (433/475)/mCherry (587/610)*
	Fluc	CFP (433/475)/mCherry (587/610)*
	Nluc	CFP (433/475)/mCherry (587/610)*

6.5. Endogenous Expression Levels

In an arbitrary system, PCAs often involve overexpression of split biosensors since these reporter genes are typically introduced in the cells by transfection. To circumvent excessively high levels of membrane expression, which may ultimately result in random collisions rather than a real interaction, low levels of biosensors should ideally be expressed, preferably close to endogenous expression levels. White et al. [148] (2017) have demonstrated that PPIs can be endogenously monitored by BRET when a Nanoluciferase reporter is genetically fused to a natively expressed GPCR of interest by CRISPR/Cas9-mediated homology-directed repair. This methodology also offers the potential to study GPCR–GPCR interactions at endogenous levels.

6.6. Kinetics

To follow the dynamics of GPCR interactions in living cells in real time or to monitor the influence of ligand-dependent modulation on the level of dimerization, one can use a real-time PCA-based method for a certain amount of time. For example, to evaluate the effect of the α -melanocyte-stimulating hormone (α -MSH) on the interaction between MC4R and MRAP2, Habara et al. [113] (2018) measured this interaction with a NanoBiT®-based PCA for up to 120 s before and after the addition of α -MSH. By doing so, the authors could demonstrate that the stimulation of MC4R with α -MSH slightly decreased the NanoBiT® signal, which led to the postulation that the activated structural change of MC4R negatively impacts the interaction with MRAP2. Overall, depending on the half-life of the substrate and the PPI of interest, different timeframes of real-time measurement can be implemented.

7. HTS with Cell-Based PCAs

High-throughput screening (HTS) assays are powerful ways to assess the influence of ligands that affect the protein complexes, thereby providing a method to measure spatial and temporal changes in the protein association in response to drugs (Figure 7). This section will focus on the potential use of HTS assays based on BiLC and BiFC for GPCRs.

7.1. GPCR Oligomerization Screening

Given that the protein interaction partners are known, PCAs based on BiLC and BiFC could aid in studying the strength of the interaction compared to a negative control. To find new (unknown) interactions, a screening assay would be handy, in which case the bait protein is fused to a split luminescent/fluorescent fragment and screened against a cDNA library fused to the other half of the reporter fragment. This way, one can identify which GPCRs interact with each other. In conjunction to this is the functional validation of the detected interactions. False readouts can be sidelined upon employment of fusion proteins with mutations at the site of interaction, which consequently inhibit the interaction and decrease the level of complementation of the luminescent or fluorescent fragments. Clearly, this is only feasible if the interaction sites are known, which may be unraveled by alanine screening. For such screening assays, a robust plate reader format would be the best choice.

A non-GPCR example of such a workflow was published by Zych et al. in 2013 [149]. These authors established a high-throughput imaging based-screening that was based on the principles of BiFC, with as the aim the study of the dynamics of Vpr (a nonstructural protein encoded by HIV-1) dimerization.

Screening based on BiFC coupled to flow cytometry is a fast way to capture weak interactions, and hence, this has been used to screen for mutations that modulate the affinity and specificity of PPIs. For example, in 2013, Morell et al. [150] coupled flow cytometry to BiFC to study the weak interaction between the SH3 domain of C-Abl kinase and its natural or mutated binding partners. This combination helps to select for good interacting partners, even if they were deficient in the whole cell population. Hence, this combination has proven to be fast, specific, and sensitive and can capture even weak, transient interactions, thus opening up a new dimension in the field of proteomics.

7.2. GPCR Drug Discovery

Since GPCR–GPCR interactions have been implied in several disease patterns, e.g., neurodegenerative disorders, a strong interest in compounds which could interfere with or alter these GPCR–GPCR interactions has arisen. Monitoring these GPCR–GPCR interactions has the potential of unveiling differences in interaction specifics, such as dynamics, as well as identifying potential therapeutic agents.

For an *in vitro* HTS assay, PCAs can be used to screen for compounds that inhibit (BiLC) or enhance (BiLC and BiFC) the PPIs. Upon addition of the modulating compound, the recorded signal changes, depending on its potential to influence the GPCR–GPCR interaction. Counter screens should be implemented to negate the possibility of an artefact, such as compounds which possibly interfere with reporter complementation, show autofluorescence, scatter light, quench luminescence or have toxic side-effects [151–153]. A key advantage of PCA-based approaches is that these include formats using living cells, so GPCR receptors are in their native environment.

Overall, this type of HTS strategy would be of great help to further enhance the interest in GPCR drug-discovery research.

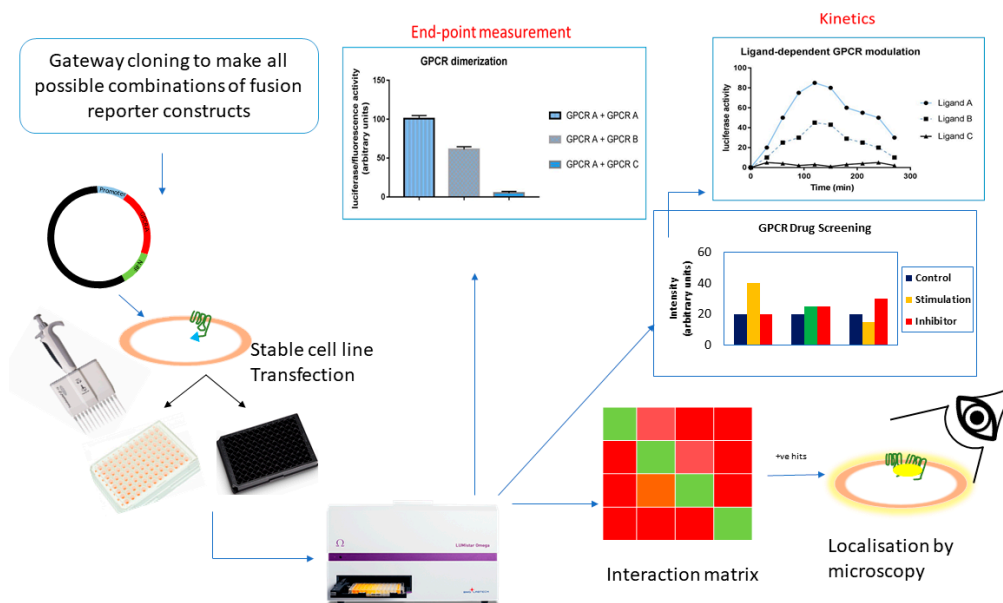


Figure 7. High-throughput screening.

8. In Vivo Application

Besides the potential as HTS assays, in recent years, BiFC and BiLC complementation techniques have also been applied in vivo, aiming at a better understanding of the dynamics of homo- or heterodimers in living animals. Several studies established an imaging reporter strategy, which can monitor the specific pharmacological regulation of oligomer complexes.

For GPCR dimerization purposes, Luker et al. [98] (2009) applied the FLuc-based PCA (NFLuc-416 and CFLuc-398) to investigate the homo- and heterodimerization of CXCR4 and CXCR7 in a tumor xenograft model of breast cancer. Type 231 cells stably expressing CXCR4-NFLuc and CXCR4-CFLuc or CXCR7-NFLuc and CXCR7-CFLuc were injected bilaterally into 6 to 8 weeks-old female nude mice. To test whether treatment with chemokine ligands produced time- and dose-dependent changes in reporter signal due to the alterations in the level of dimerization, fibroblasts stably expressing the chemokine CXCL12 were co-implanted. This chemokine reduced the level of both CXCR4 homodimerization and CXCR4–CXCR7 heterodimerization. Blocking the binding of CXCL12 to CXCR4 resulted in a time-dependent increase in the level of CXCR4 homodimer formation, but no increase in the CXCR4–CXCR7 heterodimerization was detected. In another setup, mice were injected with 231 cells stably expressing CXCR4 or CXCR7 homodimers and with fibroblasts only transduced with GFP. After 1 h of treatment with the CXCR7 modulator CCX754 (100 mg/kg subcutaneous), an increase in CXCR7 homodimer formation, but not CXCR4 homodimer formation, was detected. Overall, it was observed that the chemokine CXCL12 decreased the level of CXCR4 homodimerization and CXCR4–CXCR7 heterodimerization and the CXCR7 modulator CCX754 increased the level of CXCR7 homodimerization. The application of such newly designed drugs that specifically modify the levels of certain dimers could be beneficial in the treatment of not only breast cancer, but also in rheumatoid arthritis, HIV disease, lung and prostate cancers, where these GPCRs play an important role as well.

Rather than trying to bridge the gap between cell culture studies and *in vivo* physiology by the injection of cells stably expressing GPCRs fused to the split biosensors, use of the CRISPR/Cas9 methodology (previously described in Section 6.5) to fuse native receptors to split protein fragments would allow the detection of dimerization at endogenous protein levels. This has already been successfully implemented for the β -arrestin2 recruitment to CXCR4 [148], thus, it could likewise be tested for GPCR dimerization.

Most importantly, it has been shown that the bioluminescent reporter is sensitive at all the depths and locations within a live mouse [154], thus making this one of the best platforms to study PPIs as it better mimics a “real-life scenario”.

In the field of fluorescence, Han et al. [155] in 2014 applied BiFC based on split fragments of mNeptune, a far-red monomeric Neptune variant (600 nm/650 nm), for in vivo applications. The main obstacle of imaging in live animals, being the tissue opacity to excitation light below 600 nm, was overcome by this fluorescent protein. This system was tested by imaging the interaction between bFos and bJun in live mice and could also serve to unravel GPCR–GPCR interactions.

Finally, luminescence- or fluorescence-based systems can be applied to study GPCR oligomerization in vivo as these assays allow real-time detection of the interaction and subsequently would help in the screening of pharmacological compounds acting on these oligomers. These non-invasive animal imaging assays, for the quantification of GPCR–GPCR interactions, might aid in the transformation of treatments targeting specific oligomers (e.g., small molecules or bivalent ligands) from cell-based assays to clinical trials.

9. Conclusion

GPCRs are highly dynamic proteins that are subject to multiple spatiotemporal interactions inside the cell to trigger appropriate cellular responses. These interactions involve transducers, kinases, as well as other GPCRs, among others [156–158]. While great progress has been made during the past decade in targeting and understanding GPCR–GPCR interactions, some key questions remain unanswered. Although increasing evidence supports the existence of GPCR–GPCR interactions in vitro and in vivo [159], the sophisticated characteristics of these interactions, such as their variation in dynamics, complexity (i.e., dimerization or higher-order oligomerization), binding affinities, and the lack of a consensus sequence in the interfaces, make GPCR–GPCR interactions a challenging research field. Multiple studies based on PCA techniques have shed new light on the composition, localization, and even ligand-induced modulation of these GPCR complexes. Moreover, the potential of these PCA-based methodologies has been further explored by their combination in a variety of ways to investigate not only binary but also ternary or even quaternary partnerships. The current review describes the journey taken by these BiLC- and BiFC-based PCAs and the GPCR–GPCR interactions that have been identified by applying these techniques. PCAs are straightforward and flexible, and provide on the one hand a convenient approach to visualize the oligomerization processes in living cells and on the other hand an ideal tool to monitor the interaction dynamics. Nevertheless, certain milestones could still be achieved, such as the development of reversible split-fluorescent proteins or the discovery of split-luminescent proteins and their respective ligands that allow enhanced imaging of the GPCR–GPCR interactions. Finally, we anticipate a growing interest in the application of PCA-based high-throughput drug screenings, which makes it possible to identify compounds that selectively modulate the composition of GPCR oligomers in living cells.

Funding: This research was funded by IWT/SBO, grant number 140028 and BOF.DCV.20—BOF15/DOS/021 of Ghent University.

Conflicts of Interest: The authors declare no conflict of interest. The funders had no role in the design of the study; in the collection, analyses, or interpretation of data; in the writing of the manuscript, or in the decision to publish the results.

Abbreviations

A ₁ /A _{2A}	Adenosine receptor type 1/type 2A
α _{1b}	α _{1b} adrenergic receptor
AT ₁ /AT ₂	Angiotensin II receptor type 1/type 2
β ₂ AR	β ₂ Adrenergic receptor
BiFC	Bimolecular fluorescence complementation assay
BiLC	Bimolecular luminescence complementation assay

BRET	Bioluminescence resonance energy transfer
CAD	Cath.-a-differentiated
cAMP	Cyclic adenosine monophosphate
CB ₁	Cannabinoid receptor 1
CC	C-terminal fragment of split Cerulean
CFP	Cerulean fluorescent protein
CC2	C-C chemokine receptor type 2
CFLuc	C-terminal fragment of split Firefly luciferase
CGFP	C-terminal fragment of split green fluorescent protein
CID	Chemically induced dimerization
CN	N-terminal fragment of split Cerulean
CODA-RET	Complemented donor–acceptor resonance energy transfer
co-IP	Co-immunoprecipitation
CPA	N6-cyclopentyladenosine
CXCR4/CXCR7	C-X-C chemokine receptor type 4/type 7
D ₁ R/D ₂ R/D ₃ R	Dopamine receptor type 1/type 2/type 3
D _{2L} R	Dopamine receptor type 2 long isoform
D _{2S} R	Dopamine receptor type 2 short isoform
DsRed	Discosoma Red
ECL	Extracellular loop
EGFP	Enhanced green fluorescent protein
FLuc	Firefly Luciferase
FFAR2/3	Free fatty acid receptor type 2/type 3
FRB	FKBP-rapamycin binding domain of mTOR
FKBP12	FK506- and rapamycin-binding protein
FRET	Förster resonance energy transfer
GHSR1a	Growth hormone secretagogue receptor 1α
GLuc	<i>Gaussia princeps</i> luciferase
GPCR	G protein-coupled receptor
GRK	G protein-coupled receptor kinase
HEK	Human embryonic kidney
Homo-FC	Homo-molecular fluorescence complementation
5-HTR	5-hydroxytryptamine receptor
HTS	High-throughput screening
L-DOPA	L-3,4-dihydroxyphenylalanine
LcBiT	C-terminal fragment of split large BiT
LgBiT	Large BiT, large subunit of Nanoluciferase
LnBiT	N-terminal fragment of split large BiT
M ₁ /M ₂ /M ₄	Muscarinic acetylcholine receptor type 1/type 2/type 4
MAPK	Mitogen activated protein kinase
MBiFC	Multicolor Bimolecular fluorescence complementation assay
MC4R	Melanocortin 4 receptor
MOP	Mu-opioid receptor
MRAP2	Melanocortin 2 receptor accessory protein 2
α-MSH	α-Melanocyte-stimulating hormone
mGluR ₂ /mGluR ₅ /mGluR ₆	Metabotropic glutamate receptor type 2/type 5/type 6
mRFP	Monomer red fluorescent protein
NanoBiT®	NanoLuc Binary Technology®
NFLuc	N-terminal fragment of split Firefly luciferase
NGFP	N-terminal fragment of split green fluorescent protein
NPF2	Neuropeptide FF receptor 2
NPY	Neuropeptide Y
OX1R	Orexin 1 receptor
PCA	Protein complementation assay

PLA	Proximity ligation assay
PPI	Protein–protein interaction
P2RY1	P2Y purinoceptor 1
RLuc	<i>Renilla</i> Luciferase
SDS-PAGE	Sodium dodecyl sulfate–polyacrylamide gel electrophoresis
SH3	SRC Homology 3 Domain
SmBiT	Small BiT, small subunit of Nanoluciferase
SpIDA	Spatial intensity distribution analysis
SRET	Sequential resonance energy transfer
TAT	Transactivator of transcription
TM	Transmembrane
TSHr	Thyrotropin receptor
VC	C-terminal fragment of split Venus
VN	N-terminal fragment of split Venus
YC	C-terminal fragment of split yellow fluorescent protein
YFP	Yellow fluorescent protein
YN	N-terminal fragment of split yellow fluorescent protein

References

- Venter, J.C.; Adams, M.D.; Myers, E.W.; Li, P.W.; Mural, R.J.; Sutton, G.G.; Smith, H.O.; Yandell, M.; Evans, C.A.; Holt, R.A.; et al. The sequence of the human genome. *Science* **2001**, *291*, 1304–1351.
- Foord, S.M.; Bonner, T.I.; Neubig, R.R.; Rosser, E.M.; Pin, J.P.; Davenport, A.P.; Spedding, M.; Harmar, A.J. International Union of Pharmacology. XLVI. G protein-coupled receptor list. *Pharmacol. Rev.* **2005**, *57*, 279–288.
- Bjarnadottir, T.K.; Gloriam, D.E.; Hellstrand, S.H.; Kristiansson, H.; Fredriksson, R.; Schioth, H.B. Comprehensive repertoire and phylogenetic analysis of the G protein-coupled receptors in human and mouse. *Genomics* **2006**, *88*, 263–273.
- Rasmussen, S.G.; Choi, H.J.; Fung, J.J.; Pardon, E.; Casarosa, P.; Chae, P.S.; Devree, B.T.; Rosenbaum, D.M.; Thian, F.S.; Kobilka, T.S.; et al. Structure of a nanobody-stabilized active state of the beta(2) adrenoceptor. *Nature* **2011**, *469*, 175–180.
- Reyes-Resina, I.; Aguinaga, D.; Labandeira-Garcia, J.L.; Lanciego, J.L.; Navarro, G.; Franco, R. Usefulness of identifying G-protein-coupled receptor dimers for diagnosis and therapy of neurodegenerative diseases and of gliomas. *Histol. Histopathol.* **2018**, *10.14670/HH-11-963*, 11963.
- Zhao, J.; Deng, Y.; Jiang, Z.; Qing, H. G Protein-Coupled Receptors (GPCRs) in Alzheimer's Disease: A Focus on BACE1 Related GPCRs. *Front. Aging Neurosci.* **2016**, *8*, 58.
- Tomita, H.; Ziegler, M.E.; Kim, H.B.; Evans, S.J.; Choudary, P.V.; Li, J.Z.; Meng, F.; Dai, M.; Myers, R.M.; Neal, C.R.; et al. G protein-linked signaling pathways in bipolar and major depressive disorders. *Front. Genet.* **2013**, *4*, 297.
- Bar-Shavit, R.; Maoz, M.; Kancharla, A.; Nag, J.K.; Agranovich, D.; Grisaru-Granovsky, S.; Uziely, B. G Protein-Coupled Receptors in Cancer. *Int. J. Mol. Sci.* **2016**, *17*, 1320.
- Nieto Gutierrez, A.; McDonald, P.H. GPCRs: Emerging anti-cancer drug targets. *Cell. Signal.* **2018**, *41*, 65–74.
- Huang, Y.H.; Todd, N.; Thathiah, A. The role of GPCRs in neurodegenerative diseases: avenues for therapeutic intervention. *Curr. Opin. Pharmacol.* **2017**, *32*, 96–110.
- Hauser, A.S.; Chavali, S.; Masuho, I.; Jahn, L.J.; Martemyanov, K.A.; Gloriam, D.E.; Babu, M.M. Pharmacogenomics of GPCR Drug Targets. *Cell* **2018**, *172*, 41–54.e19.
- Waldhoer, M.; Fong, J.; Jones, R.M.; Lunzer, M.M.; Sharma, S.K.; Kostenis, E.; Portoghese, P.S.; Whistler, J.L. A heterodimer-selective agonist shows in vivo relevance of G protein-coupled receptor dimers. *Proc. Natl. Acad. Sci. USA* **2005**, *102*, 9050–9055.
- Vinals, X.; Moreno, E.; Lanfumey, L.; Cordini, A.; Pastor, A.; de La Torre, R.; Gasperini, P.; Navarro, G.; Howell, L.A.; Pardo, L.; et al. Cognitive Impairment Induced by Delta9-tetrahydrocannabinol Occurs through Heteromers between Cannabinoid CB1 and Serotonin 5-HT2A Receptors. *PLoS Biol.* **2015**, *13*, e1002194.

14. Qian, M.C.; Vasudevan, L.; Huysentruyt, J.; Risseeuw, M.D.P.; Stove, C.; Vanderheyden, P.M.L.; Van Craenenbroeck, K.; Van Calenbergh, S. Design, Synthesis, and Biological Evaluation of Bivalent Ligands Targeting Dopamine D-2-Like Receptors and the -Opioid Receptor. *Chemmedchem* **2018**, *13*, 944–956.
15. Qian, M.; Wouters, E.; Dalton, J.A.R.; Risseeuw, M.D.P.; Crans, R.A.J.; Stove, C.; Giraldo, J.; Van Craenenbroeck, K.; Van Calenbergh, S. Synthesis toward Bivalent Ligands for the Dopamine D2 and Metabotropic Glutamate 5 Receptors. *J. Med. Chem.* **2018**, *61*, 8212–8225.
16. Venkatesan, K.; Rual, J.F.; Vazquez, A.; Stelzl, U.; Lemmens, I.; Hirozane-Kishikawa, T.; Hao, T.; Zenkner, M.; Xin, X.F.; Goh, K.I.; et al. An empirical framework for binary interactome mapping. *Nat. Methods* **2009**, *6*, 83–90.
17. Zhang, Q.C.; Petrey, D.; Deng, L.; Qiang, L.; Shi, Y.; Thu, C.A.; Bisikirska, B.; Lefebvre, C.; Accili, D.; Hunter, T.; et al. Structure-based prediction of protein-protein interactions on a genome-wide scale. *Nature* **2012**, *490*, 556–560.
18. AbdAlla, S.; Lothar, H.; el Massiery, A.; Quitterer, U. Increased AT(1) receptor heterodimers in preeclampsia mediate enhanced angiotensin II responsiveness. *Nat. Med.* **2001**, *7*, 1003–1009.
19. Seeman, P. Dopamine D2 receptors as treatment targets in schizophrenia. *Clin. Schizophr. Relat. Psychoses* **2010**, *4*, 56–73.
20. Renner, U.; Zeug, A.; Woehler, A.; Niebert, M.; Dityatev, A.; Dityateva, G.; Gorinski, N.; Guseva, D.; Abdel-Galil, D.; Frohlich, M.; et al. Heterodimerization of serotonin receptors 5-HT1A and 5-HT7 differentially regulates receptor signalling and trafficking. *J. Cell Sci.* **2012**, *125*, 2486–2499.
21. Borroto-Escuela, D.O.; Romero-Fernandez, W.; Perez-Alea, M.; Narvaez, M.; Tarakanov, A.O.; Mudo, G.; Agnati, L.F.; Ciruela, F.; Belluardo, N.; Fuxe, K. The existence of FGFR1-5-HT1A receptor heterocomplexes in midbrain 5-HT neurons of the rat: relevance for neuroplasticity. *J. Neurosci.* **2012**, *32*, 6295–6303.
22. Borroto-Escuela, D.O.; Corrales, F.; Narvaez, M.; Oflijan, J.; Agnati, L.F.; Palkovits, M.; Fuxe, K. Dynamic modulation of FGFR1-5-HT1A heteroreceptor complexes. Agonist treatment enhances participation of FGFR1 and 5-HT1A homodimers and recruitment of beta-arrestin2. *Biochem. Biophys. Res. Commun.* **2013**, *441*, 387–392.
23. Borroto-Escuela, D.O.; Perez-Alea, M.; Narvaez, M.; Tarakanov, A.O.; Mudo, G.; Jimenez-Beristain, A.; Agnati, L.F.; Ciruela, F.; Belluardo, N.; Fuxe, K. Enhancement of the FGFR1 signaling in the FGFR1-5-HT1A heteroreceptor complex in midbrain raphe 5-HT neuron systems. Relevance for neuroplasticity and depression. *Biochem. Biophys. Res. Commun.* **2015**, *463*, 180–186.
24. Borroto-Escuela, D.O.; Tarakanov, A.O.; Fuxe, K. FGFR1-5-HT1A Heteroreceptor Complexes: Implications for Understanding and Treating Major Depression. *Trends Neurosci.* **2016**, *39*, 5–15.
25. Millon, C.; Flores-Burgess, A.; Narvaez, M.; Borroto-Escuela, D.O.; Santin, L.; Parrado, C.; Narvaez, J.A.; Fuxe, K.; Diaz-Cabiale, Z. A role for galanin N-terminal fragment (1-15) in anxiety- and depression-related behaviors in rats. *Int. J. Neuropsychopharmacol.* **2014**, *18*, pyu064.
26. Millon, C.; Flores-Burgess, A.; Narvaez, M.; Borroto-Escuela, D.O.; Santin, L.; Gago, B.; Narvaez, J.A.; Fuxe, K.; Diaz-Cabiale, Z. Galanin (1-15) enhances the antidepressant effects of the 5-HT1A receptor agonist 8-OH-DPAT: involvement of the raphe-hippocampal 5-HT neuron system. *Brain Struct. Funct.* **2016**, *221*, 4491–4504.
27. Rondard, P.; Huang, S.; Monnier, C.; Tu, H.; Blanchard, B.; Oueslati, N.; Malhaire, F.; Li, Y.; Trinquet, E.; Labesse, G.; et al. Functioning of the dimeric GABA(B) receptor extracellular domain revealed by glycan wedge scanning. *EMBO J.* **2008**, *27*, 1321–1332.
28. Maggio, R.; Millan, M.J. Dopamine D2-D3 receptor heteromers: Pharmacological properties and therapeutic significance. *Curr. Opin. Pharmacol.* **2010**, *10*, 100–107.
29. So, C.H.; Varghese, G.; Curley, K.J.; Kong, M.M.; Alijaniam, M.; Ji, X.; Nguyen, T.; O'Dowd, B. F.; George, S.R. D1 and D2 dopamine receptors form heterooligomers and cointernalize after selective activation of either receptor. *Mol. Pharmacol.* **2005**, *68*, 568–578.
30. Petrin, D.; Hebert, T.E. The functional size of GPCRs - monomers, dimers or tetramers? *Subcell. Biochem.* **2012**, *63*, 67–81.
31. Bouvier, M.; Hébert, T.E. CrossTalk proposal: Weighing the evidence for Class A GPCR dimers, the evidence favours dimers. *J. Physiol.* **2014**, *592*, 2439–2441.
32. Lambert, N.A.; Javitch, J.A. CrossTalk opposing view: Weighing the evidence for class A GPCR dimers, the jury is still out. *J. Physiol.* **2014**, *592*, 2443–2445.

33. Felce, J.H.; Latty, S.L.; Knox, R.G.; Mattick, S.R.; Lui, Y.; Lee, S.F.; Klenerman, D.; Davis, S.J. Receptor Quaternary Organization Explains G Protein-Coupled Receptor Family Structure. *Cell Rep.* **2017**, *20*, 2654–2665.
34. Herrick-Davis, K.; Milligan, G.; Di Giovanni, G. *G-Protein-Coupled Receptor Dimers*; Springer: Berlin/Heidelberg, Germany, 2017.
35. Marsango, S.; Caltabiano, G.; Jimenez-Roses, M.; Millan, M.J.; Pediani, J.D.; Ward, R.J.; Milligan, G. A Molecular Basis for Selective Antagonist Destabilization of Dopamine D3 Receptor Quaternary Organization. *Sci. Rep.* **2017**, *7*, 2134.
36. Nenasheva, T.A.; Neary, M.; Mashanov, G.I.; Birdsall, N.J.; Breckenridge, R.A.; Molloy, J.E. Abundance, distribution, mobility and oligomeric state of M(2) muscarinic acetylcholine receptors in live cardiac muscle. *J. Mol. Cell. Cardiol.* **2013**, *57*, 129–136.
37. Hern, J.A.; Baig, A.H.; Mashanov, G.I.; Birdsall, B.; Corrie, J.E.; Lazareno, S.; Molloy, J.E.; Birdsall, N.J. Formation and dissociation of M1 muscarinic receptor dimers seen by total internal reflection fluorescence imaging of single molecules. *Proc. Natl. Acad. Sci. USA* **2010**, *107*, 2693–2698.
38. Kasai, R.S.; Kusumi, A. Single-molecule imaging revealed dynamic GPCR dimerization. *Curr. Opin. Cell Biol.* **2014**, *27*, 78–86.
39. Kasai, R.S.; Ito, S.V.; Awane, R.M.; Fujiwara, T.K.; Kusumi, A. The Class-A GPCR Dopamine D2 Receptor Forms Transient Dimers Stabilized by Agonists: Detection by Single-Molecule Tracking. *Cell Biochem. Biophys.* **2018**, *76*, 29–37.
40. Wouters, E.; Marin, A.R.; Dalton, J.A.R.; Giraldo, J.; Stove, C. Distinct Dopamine D(2) Receptor Antagonists Differentially Impact D(2) Receptor Oligomerization. *Int. J. Mol. Sci.* **2019**, *20*, 1686.
41. Van Craenenbroeck, K. GPCR oligomerization: contribution to receptor biogenesis. *Sub-Cell. Biochem.* **2012**, *63*, 43–65.
42. Baragli, A.; Alturaihi, H.; Watt, H.L.; Abdallah, A.; Kumar, U. Heterooligomerization of human dopamine receptor 2 and somatostatin receptor 2 Co-immunoprecipitation and fluorescence resonance energy transfer analysis. *Cell. Signal.* **2007**, *19*, 2304–2316.
43. Mizuno, N.; Suzuki, T.; Kishimoto, Y.; Hirasawa, N. Biochemical assay of G protein-coupled receptor oligomerization: Adenosine A1 and thromboxane A2 receptors form the novel functional hetero-oligomer. *Methods Cell Biol.* **2013**, *117*, 213–227.
44. Skieterska, K.; Duchou, J.; Lintermans, B.; Van Craenenbroeck, K. Detection of G protein-coupled receptor (GPCR) dimerization by coimmunoprecipitation. *Methods Cell Biol.* **2013**, *117*, 323–340.
45. Ciruela, F.; Vilardaga, J.P.; Fernandez-Duenas, V. Lighting up multiprotein complexes: lessons from GPCR oligomerization. *Trends Biotechnol.* **2010**, *28*, 407–415.
46. Ward, R.J.; Pediani, J.D.; Godin, A.G.; Milligan, G. Regulation of oligomeric organization of the serotonin 5-hydroxytryptamine 2C (5-HT2C) receptor observed by spatial intensity distribution analysis. *J. Biol. Chem.* **2015**, *290*, 12844–12857.
47. Ward, R.J.; Pediani, J.D.; Harikumar, K.G.; Miller, L.J.; Milligan, G. Spatial intensity distribution analysis quantifies the extent and regulation of homodimerization of the secretin receptor. *Biochem. J.* **2017**, *474*, 1879–1895.
48. Pediani, J.D.; Ward, R.J.; Godin, A.G.; Marsango, S.; Milligan, G. Dynamic Regulation of Quaternary Organization of the M1 Muscarinic Receptor by Subtype-selective Antagonist Drugs. *J. Biol. Chem.* **2016**, *291*, 13132–13146.
49. Pediani, J.D.; Ward, R.J.; Marsango, S.; Milligan, G. Spatial Intensity Distribution Analysis: Studies of G Protein-Coupled Receptor Oligomerisation. *Trends Pharmacol. Sci.* **2018**, *39*, 175–186.
50. Wehr, M.C.; Rossner, M.J. Split protein biosensor assays in molecular pharmacological studies. *Drug Discov. Today* **2016**, *21*, 415–429.
51. Ghosh, I.; Hamilton, A.D.; Regan, L. Antiparallel leucine zipper-directed protein reassembly: Application to the green fluorescent protein. *J. Am. Chem. Soc.* **2000**, *122*, 5658–5659.
52. Hu, C.D.; Chinenov, Y.; Kerppola, T.K. Visualization of interactions among bZip and Rel family proteins in living cells using bimolecular fluorescence complementation. *Mol. Cell* **2002**, *9*, 789–798.
53. Shyu, Y.J.; Liu, H.; Deng, X.H.; Hu, C.D. Identification of new fluorescent protein fragments for bimolecular fluorescence complementation analysis under physiological conditions. *Biotechniques* **2006**, *40*, 61–66.

54. Cevheroglu, O.; Kumas, G.; Hauser, M.; Becker, J.M.; Son, C.D. The yeast Ste2p G protein-coupled receptor dimerizes on the cell plasma membrane. *Bba-Biomembranes* **2017**, *1859*, 698–711.
55. Cabantous, S.; Nguyen, H.B.; Pedelacq, J.D.; Koraichi, F.; Chaudhary, A.; Ganguly, K.; Lockard, M.A.; Favre, G.; Terwilliger, T.C.; Waldo, G.S. A new protein-protein interaction sensor based on tripartite split-GFP association. *Sci. Rep.* **2013**, *3*, 2854.
56. Nagai, T.; Sawano, A.; Park, E.S.; Miyawaki, A. Circularly permuted green fluorescent proteins engineered to sense Ca²⁺. *Proc. Natl. Acad. Sci. USA* **2001**, *98*, 3197–3202.
57. Griesbeck, O.; Baird, G.S.; Campbell, R.E.; Zacharias, D.A.; Tsien, R.Y. Reducing the environmental sensitivity of yellow fluorescent protein. Mechanism and applications. *J. Biol. Chem.* **2001**, *276*, 29188–29194.
58. Nagai, T.; Ibata, K.; Park, E.S.; Kubota, M.; Mikoshiba, K.; Miyawaki, A. A variant of yellow fluorescent protein with fast and efficient maturation for cell-biological applications. *Nat. Biotechnol.* **2002**, *20*, 87–90.
59. Vidi, P.A.; Chemel, B.R.; Hu, C.D.; Watts, V.J. Ligand-dependent oligomerization of dopamine D(2) and adenosine A(2A) receptors in living neuronal cells. *Mol. Pharmacol.* **2008**, *74*, 544–551.
60. Vidi, P.-A.; Chen, J.; Irudayaraj, J.M.K.; Watts, V.J. Adenosine A2A receptors assemble into higher-order oligomers at the plasma membrane. *FEBS Lett.* **2008**, *582*, 3985–3990.
61. Kilpatrick, L.E.; Humphrys, L.J.; Holliday, N.D. A G protein-coupled receptor dimer imaging assay reveals selectively modified pharmacology of neuropeptide Y Y1/Y5 receptor heterodimers. *Mol. Pharmacol.* **2015**, *87*, 718–732.
62. Przybyla, J.A.; Watts, V.J. Ligand-induced regulation and localization of cannabinoid CB1 and dopamine D2L receptor heterodimers. *J. Pharmacol. Exp. Ther.* **2010**, *332*, 710–719.
63. Ang, Z.; Xiong, D.; Wu, M.; Ding, J.L. FFAR2-FFAR3 receptor heteromerization modulates short-chain fatty acid sensing. *FASEB J.* **2018**, *32*, 289–303.
64. Xue, Q.; Bai, B.; Ji, B.; Chen, X.; Wang, C.; Wang, P.; Yang, C.; Zhang, R.; Jiang, Y.; Pan, Y.; et al. Ghrelin Through GHSR1a and OX1R Heterodimers Reveals a Galphas-cAMP-cAMP Response Element Binding Protein Signaling Pathway in Vitro. *Front. Mol. Neurosci.* **2018**, *11*, 245.
65. Navarro, G.; Cordomi, A.; Brugarolas, M.; Moreno, E.; Aguinaga, D.; Perez-Benito, L.; Ferre, S.; Cortes, A.; Casado, V.; Mallol, J.; et al. Cross-communication between Gi and Gs in a G-protein-coupled receptor heterotetramer guided by a receptor C-terminal domain. *BMC Biol.* **2018**, *16*, 24.
66. Hinz, S.; Navarro, G.; Borroto-Escuela, D.; Seibt, B.F.; Ammon, Y.C.; de Filippo, E.; Danish, A.; Lacher, S.K.; Cervinkova, B.; Rafehi, M.; et al. Adenosine A2A receptor ligand recognition and signaling is blocked by A2B receptors. *Oncotarget* **2018**, *9*, 13593–13611.
67. Song, D.; Jung, Y. Homo-molecular Fluorescence Complementation for Direct Visualization of Receptor Oligomerization in Living Cells. *Angew. Chem. (Int. Ed. Engl.)* **2019**, *58*, 2045–2049.
68. Vidi, P.A.; Chemel, B.R.; Watts, V.J. Direct visualization of adenosine A(2A) and dopamine D-2L receptor oligomers in a neuronal cell model. *FASEB J.* **2008**, *22*.
69. Baird, G.S.; Zacharias, D.A.; Tsien, R.Y. Biochemistry, mutagenesis, and oligomerization of DsRed, a red fluorescent protein from coral. *Proc. Natl. Acad. Sci. USA* **2000**, *97*, 11984–11989.
70. Fan, J.Y.; Cui, Z.Q.; Wei, H.P.; Zhang, Z.P.; Zhou, Y.F.; Wang, Y.P.; Zhang, X.E. Split mCherry as a new red bimolecular fluorescence complementation system for visualizing protein-protein interactions in living cells. *Biochem. Biophys. Res. Commun.* **2008**, *367*, 47–53.
71. Jach, G.; Pesch, M.; Richter, K.; Frings, S.; Uhrig, J.F. An improved mRFP1 adds red to bimolecular fluorescence complementation. *Nat. Methods* **2006**, *3*, 597–600.
72. Shcherbo, D.; Merzlyak, E.M.; Chepurnykh, T.V.; Fradkov, A.F.; Ermakova, G.V.; Solovieva, E.A.; Lukyanov, K.A.; Bogdanova, E.A.; Zaraisky, A.G.; Lukyanov, S.; et al. Bright far-red fluorescent protein for whole-body imaging. *Nat. Methods* **2007**, *4*, 741–746.
73. Chu, J.; Zhang, Z.; Zheng, Y.; Yang, J.; Qin, L.; Lu, J.; Huang, Z.L.; Zeng, S.; Luo, Q. A novel far-red bimolecular fluorescence complementation system that allows for efficient visualization of protein interactions under physiological conditions. *Biosens. Bioelectron.* **2009**, *25*, 234–239.
74. Hu, C.D.; Kerppola, T.K. Simultaneous visualization of multiple protein interactions in living cells using multicolor fluorescence complementation analysis. *Nat. Biotechnol.* **2003**, *21*, 539–545.
75. Vidi, P.A.; Przybyla, J.A.; Hu, C.D.; Watts, V.J. Visualization of G protein-coupled receptor (GPCR) interactions in living cells using bimolecular fluorescence complementation (BiFC). *Curr. Protoc. Neurosci.* **2010**, *Chapter 5*, Unit 5.29, doi:10.1002/0471142301.ns0529s51.

76. Vidi, P.A.; Watts, V.J. Fluorescent and bioluminescent protein-fragment complementation assays in the study of G protein-coupled receptor oligomerization and signaling. *Mol. Pharmacol.* **2009**, *75*, 733–739.
77. Rebois, R.V.; Robitaille, M.; Petrin, D.; Zylbergold, P.; Trieu, P.; Hebert, T.E. Combining protein complementation assays with resonance energy transfer to detect multipartner protein complexes in living cells. *Methods (San Diego, Calif.)* **2008**, *45*, 214–218.
78. Navarro, G.; Carriba, P.; Gandia, J.; Ciruela, F.; Casado, V.; Cortes, A.; Mallol, J.; Canela, E.I.; Lluís, C.; Franco, R. Detection of heteromers formed by cannabinoid CB1, dopamine D2, and adenosine A2A G-protein-coupled receptors by combining bimolecular fluorescence complementation and bioluminescence energy transfer. *Sci. World J.* **2008**, *8*, 1088–1097.
79. Shyu, Y.J.; Suarez, C.D.; Hu, C.D. Visualization of AP-1 NF-kappaB ternary complexes in living cells by using a BiFC-based FRET. *Proc. Natl. Acad. Sci. USA* **2008**, *105*, 151–156.
80. Bagher, A.M.; Laprairie, R.B.; Toguri, J.T.; Kelly, M.E.M.; Denovan-Wright, E.M. Bidirectional allosteric interactions between cannabinoid receptor 1 (CB1) and dopamine receptor 2 long (D2L) heterotetramers. *Eur. J. Pharmacol.* **2017**, *813*, 66–83.
81. Spencer, D.M.; Wandless, T.J.; Schreiber, S.L.; Crabtree, G.R. Controlling Signal-Transduction with Synthetic Ligands. *Science* **1993**, *262*, 1019–1024.
82. Cao, X.; Liang, L.; Hadcock, J.R.; Iredale, P.A.; Griffith, D.A.; Menniti, F.S.; Factor, S.; Greenamyre, J.T.; Papa, S.M. Blockade of cannabinoid type 1 receptors augments the antiparkinsonian action of levodopa without affecting dyskinesias in 1-methyl-4-phenyl-1,2,3,6-tetrahydropyridine-treated rhesus monkeys. *J. Pharmacol. Exp. Ther.* **2007**, *323*, 318–326.
83. Fan, F.; Wood, K.V. Bioluminescent assays for high-throughput screening. *Assay Drug Dev. Technol.* **2007**, *5*, 127–136.
84. Luker, K.E.; Smith, M.C.; Luker, G.D.; Gammon, S.T.; Piwnica-Worms, H.; Piwnica-Worms, D. Kinetics of regulated protein-protein interactions revealed with firefly luciferase complementation imaging in cells and living animals. *Proc. Natl. Acad. Sci. USA* **2004**, *101*, 12288–12293.
85. Paulmurugan, R.; Umezawa, Y.; Gambhir, S.S. Noninvasive imaging of protein-protein interactions in living subjects by using reporter protein complementation and reconstitution strategies. *Proc. Natl. Acad. Sci. USA* **2002**, *99*, 15608–15613.
86. Paulmurugan, R.; Gambhir, S.S. Monitoring protein-protein interactions using split synthetic renilla luciferase protein-fragment-assisted complementation. *Anal. Chem.* **2003**, *75*, 1584–1589.
87. Paulmurugan, R.; Massoud, T.F.; Huang, J.; Gambhir, S.S. Molecular imaging of drug-modulated protein-protein interactions in living subjects. *Cancer Res.* **2004**, *64*, 2113–2119.
88. Remy, I.; Michnick, S.W. A highly sensitive protein-protein interaction assay based on Gaussia luciferase. *Nat. Methods* **2006**, *3*, 977–979.
89. Bodle, C.R.; Hayes, M.P.; O'Brien, J.B.; Roman, D.L. Development of a bimolecular luminescence complementation assay for RGS: G protein interactions in cells. *Anal. Biochem.* **2017**, *522*, 10–17.
90. Oh-Hashi, K.; Hirata, Y.; Kiuchi, K. SOD1 dimerization monitoring using a novel split NanoLuc, NanoBit. *Cell Biochem. Funct.* **2016**, *34*, 497–504.
91. Verhoef, L.G.; Mattioli, M.; Ricci, F.; Li, Y.C.; Wade, M. Multiplex detection of protein-protein interactions using a next generation luciferase reporter. *Biochim. Biophys. Acta* **2016**, *1863*, 284–292.
92. Dixon, A.S.; Schwinn, M.K.; Hall, M.P.; Zimmerman, K.; Otto, P.; Lubben, T.H.; Butler, B.L.; Binkowski, B.F.; Machleidt, T.; Kirkland, T.A.; et al. NanoLuc Complementation Reporter Optimized for Accurate Measurement of Protein Interactions in Cells. *ACS Chem. Biol.* **2016**, *11*, 400–408.
93. Galarneau, A.; Primeau, M.; Trudeau, L.E.; Michnick, S.W. beta-Lactamase protein fragment complementation assays as in vivo and in vitro sensors of protein-protein interactions. *Nat. Biotechnol.* **2002**, *20*, 619–622.
94. Paulmurugan, R.; Gambhir, S.S. Novel fusion protein approach for efficient high-throughput screening of small molecule-mediating protein-protein interactions in cells and living animals. *Cancer Res.* **2005**, *65*, 7413–7420.
95. Matthews, J.C.; Hori, K.; Cormier, M.J. Substrate and substrate analogue binding properties of Renilla luciferase. *Biochemistry-US* **1977**, *16*, 5217–5220.
96. Matthews, J.C.; Hori, K.; Cormier, M.J. Purification and properties of Renilla reniformis luciferase. *Biochemistry-US* **1977**, *16*, 85–91.

97. Loening, A.M.; Fenn, T.D.; Wu, A.M.; Gambhir, S.S. Consensus guided mutagenesis of Renilla luciferase yields enhanced stability and light output. *Protein Eng. Des. Sel.* **2006**, *19*, 391–400.
98. Luker, K.E.; Gupta, M.; Luker, G.D. Imaging chemokine receptor dimerization with firefly luciferase complementation. *FASEB J.* **2009**, *23*, 823–834.
99. Armando, S.; Quoyer, J.; Lukashova, V.; Maiga, A.; Percherancier, Y.; Heveker, N.; Pin, J.P.; Prezeau, L.; Bouvier, M. The chemokine CXC4 and CC2 receptors form homo- and heterooligomers that can engage their signaling G-protein effectors and betaarrestin. *FASEB J.* **2014**, *28*, 4509–4523.
100. Guo, W.; Urizar, E.; Kralikova, M.; Mobarec, J.C.; Shi, L.; Filizola, M.; Javitch, J.A. Dopamine D2 receptors form higher order oligomers at physiological expression levels. *EMBO J.* **2008**, *27*, 2293–2304.
101. Casado-Anguera, V.; Bonaventura, J.; Moreno, E.; Navarro, G.; Cortes, A.; Ferre, S.; Casado, V. Evidence for the heterotetrameric structure of the adenosine A2A-dopamine D2 receptor complex. *Biochem. Soc. Trans.* **2016**, *44*, 595–600.
102. Hall, M.P.; Unch, J.; Binkowski, B.F.; Valley, M.P.; Butler, B.L.; Wood, M.G.; Otto, P.; Zimmerman, K.; Vidugiris, G.; Machleidt, T.; et al. Engineered luciferase reporter from a deep sea shrimp utilizing a novel imidazopyrazinone substrate. *ACS Chem. Biol.* **2012**, *7*, 1848–1857.
103. England, C.G.; Ehlerding, E.B.; Cai, W. NanoLuc: A Small Luciferase Is Brightening Up the Field of Bioluminescence. *Bioconjug. Chem.* **2016**, *27*, 1175–1187.
104. Cannaert, A.; Franz, F.; Auwarter, V.; Stove, C.P. Activity-Based Detection of Consumption of Synthetic Cannabinoids in Authentic Urine Samples Using a Stable Cannabinoid Reporter System. *Anal. Chem.* **2017**, *89*, 9527–9536.
105. Cannaert, A.; Storme, J.; Franz, F.; Auwarter, V.; Stove, C.P. Detection and Activity Profiling of Synthetic Cannabinoids and Their Metabolites with a Newly Developed Bioassay. *Anal. Chem.* **2016**, *88*, 11476–11485.
106. Dupuis, N.; Laschet, C.; Franssen, D.; Szpakowska, M.; Gilissen, J.; Geubelle, P.; Soni, A.; Parent, A.S.; Pirotte, B.; Chevigne, A.; et al. Activation of the Orphan G Protein-Coupled Receptor GPR27 by Surrogate Ligands Promotes beta-Arrestin 2 Recruitment. *Mol. Pharmacol.* **2017**, *91*, 595–608.
107. Noble, C.; Cannaert, A.; Linnet, K.; Stove, C.P. Application of an activity-based receptor bioassay to investigate the in vitro activity of selected indole- and indazole-3-carboxamide-based synthetic cannabinoids at CB1 and CB2 receptors. *Drug Test. Anal.* **2019**, *11*, 501–511.
108. Storme, J.; Cannaert, A.; Van Craenenbroeck, K.; Stove, C.P. Molecular dissection of the human A3 adenosine receptor coupling with beta-arrestin2. *Biochem. Pharmacol.* **2018**, *148*, 298–307.
109. Wouters, E.; Mogler, L.; Cannaert, A.; Auwärter, V.; Stove, C. Functional evaluation of carboxy metabolites of synthetic cannabinoid receptor agonists featuring scaffolds based on L-valine or L-tert-leucine. *Drug Test. Anal.* **2019**, doi:10.1002/dta.2607.
110. Cannaert, A.; Vasudevan, L.; Friscia, M.; Mohr, A.L.A.; Wille, S.M.R.; Stove, C.P. Activity-Based Concept to Screen Biological Matrices for Opiates and (Synthetic) Opioids. *Clin. Chem.* **2018**, *64*, 1221–1229.
111. Laschet, C.; Dupuis, N.; Hanson, J. A dynamic and screening-compatible nanoluciferase-based complementation assay enables profiling of individual GPCR-G protein interactions. *J. Biol. Chem.* **2018**, doi:10.1074/jbc.RA118.006231.
112. Wouters, E.; Vasudevan, L.; Ciruela, F.; Saini, D.K.; Stove, C.; Van Craenenbroeck, K. Assessing GPCR Dimerization in Living Cells: Comparison of the NanoBiT Assay with Related Bioluminescence- and Fluorescence-Based Approaches. In *Receptor-Receptor Interactions in the Central Nervous System*; Fuxe, K., Borroto-Escuela, D.O., Eds.; Springer: New York, NY, USA, 2018; pp. 239–250, doi:10.1007/978-1-4939-8576-0_15.
113. Habara, M.; Mori, N.; Okada, Y.; Kawasumi, K.; Nakao, N.; Tanaka, Y.; Arai, T.; Yamamoto, I. Molecular characterization of feline melanocortin 4 receptor and melanocortin 2 receptor accessory protein 2. *Gen. Comp. Endocrinol.* **2018**, *261*, 31–39.
114. Dixon, A.S.; Kim, S.J.; Baumgartner, B.K.; Krippner, S.; Owen, S.C. A Tri-part Protein Complementation System Using AntibodySmall Peptide Fusions Enables Homogeneous Immunoassays. *Sci. Rep.-UK* **2017**, *7*, 8186.
115. Ohmuro-Matsuyama, Y.; Ueda, H. Homogeneous Noncompetitive Luminescent Immunodetection of Small Molecules by Ternary Protein Fragment Complementation. *Anal. Chem.* **2018**, *90*, 3001–3004.
116. Moustaqil, M.; Bhumkar, A.; Gonzalez, L.; Raoul, L.; Hunter, D.J.B.; Carrive, P.; Sierrecki, E.; Gambin, Y. A Split-Luciferase Reporter Recognizing GFP and mCherry Tags to Facilitate Studies of Protein-Protein Interactions. *Int. J. Mol. Sci.* **2017**, *18*, 2681.

117. Sahlholm, K.; Gomez-Soler, M.; Valle-Leon, M.; Lopez-Cano, M.; Taura, J.J.; Ciruela, F.; Fernandez-Duenas, V. Antipsychotic-Like Efficacy of Dopamine D2 Receptor-Biased Ligands is Dependent on Adenosine A2A Receptor Expression. *Mol. Neurobiol.* **2018**, *55*, 4952–4958.
118. Bonaventura, J.; Navarro, G.; Casado-Anguera, V.; Azdad, K.; Rea, W.; Moreno, E.; Brugarolas, M.; Mallol, J.; Canela, E.I.; Lluís, C.; et al. Allosteric interactions between agonists and antagonists within the adenosine A2A receptor-dopamine D2 receptor heterotetramer. *Proc. Natl. Acad. Sci. USA* **2015**, *112*, E3609–E3618.
119. Kerppola, T.K. Design and implementation of bimolecular fluorescence complementation (BiFC) assays for the visualization of protein interactions in living cells. *Nat. Protoc.* **2006**, *1*, 1278–1286.
120. Gokhale, R.S.; Khosla, C. Role of linkers in communication between protein modules. *Curr. Opin. Chem. Biol.* **2000**, *4*, 22–27.
121. Argos, P. An investigation of oligopeptides linking domains in protein tertiary structures and possible candidates for general gene fusion. *J. Mol. Biol.* **1990**, *211*, 943–958.
122. Chen, X.; Zaro, J.L.; Shen, W.C. Fusion protein linkers: Property, design and functionality. *Adv. Drug Deliv. Rev.* **2013**, *65*, 1357–1369.
123. Carayon, K.; Mouledous, L.; Combedazou, A.; Mazeret, S.; Haanappel, E.; Salome, L.; Mollereau, C. Heterologous regulation of Mu-opioid (MOP) receptor mobility in the membrane of SH-SY5Y cells. *J. Biol. Chem.* **2014**, *289*, 28697–28706.
124. Stefan, E.; Aquin, S.; Berger, N.; Landry, C.R.; Nyfeler, B.; Bouvier, M.; Michnick, S.W. Quantification of dynamic protein complexes using Renilla luciferase fragment complementation applied to protein kinase A activities in vivo. *Proc. Natl. Acad. Sci. USA* **2007**, *104*, 16916–16921.
125. Porrello, E.R.; Pflieger, K.D.; Seeber, R.M.; Qian, H.; Oro, C.; Abogadie, F.; Delbridge, L.M.; Thomas, W.G. Heteromerization of angiotensin receptors changes trafficking and arrestin recruitment profiles. *Cell. Signal.* **2011**, *23*, 1767–1776.
126. Lalonde, S.; Weise, A.; Walsh, R.P.; Ward, J.M.; Frommer, W.B. Fusion to GFP blocks intercellular trafficking of the sucrose transporter SUT1 leading to accumulation in companion cells. *BMC Plant Biol.* **2003**, *3*, 8.
127. Meyer, E.; Fromherz, P. Ca²⁺ activation of hSlo K⁺ channel is suppressed by N-terminal GFP tag. *The European journal of neuroscience* **1999**, *11*, 1105–1108.
128. Cabantous, S.; Terwilliger, T.C.; Waldo, G.S. Protein tagging and detection with engineered self-assembling fragments of green fluorescent protein. *Nat. Biotechnol.* **2005**, *23*, 102–107.
129. Morell, M.; Ventura, S.; Aviles, F.X. Protein complementation assays: approaches for the in vivo analysis of protein interactions. *FEBS Lett.* **2009**, *583*, 1684–1691.
130. Cabello, N.; Gandia, J.; Bertarelli, D.C.; Watanabe, M.; Lluís, C.; Franco, R.; Ferre, S.; Lujan, R.; Ciruela, F. Metabotropic glutamate type 5, dopamine D2 and adenosine A2a receptors form higher-order oligomers in living cells. *J. Neurochem.* **2009**, *109*, 1497–1507.
131. Wang, M.; Pei, L.; Fletcher, P.J.; Kapur, S.; Seeman, P.; Liu, F. Schizophrenia, amphetamine-induced sensitized state and acute amphetamine exposure all show a common alteration: increased dopamine D2 receptor dimerization. *Mol. Brain* **2010**, *3*, 25.
132. Zawarynski, P.; Tallerico, T.; Seeman, P.; Lee, S.P.; O'Dowd, B.F.; George, S.R. Dopamine D2 receptor dimers in human and rat brain. *FEBS Lett.* **1998**, *441*, 383–386.
133. Bonaventura, J.; Rico, A.J.; Moreno, E.; Sierra, S.; Sanchez, M.; Luquin, N.; Farre, D.; Muller, C.E.; Martinez-Pinilla, E.; Cortes, A.; et al. L-DOPA-treatment in primates disrupts the expression of A(2A) adenosine-CB1 cannabinoid-D-2 dopamine receptor heteromers in the caudate nucleus. *Neuropharmacology* **2014**, *79*, 90–100.
134. Fink, J.S.; Weaver, D.R.; Rivkees, S.A.; Peterfreund, R.A.; Pollack, A.E.; Adler, E.M.; Reppert, S.M. Molecular-Cloning of the Rat Adenosine-A2 Receptor—Selective Coexpression with D2-Dopamine Receptors in Rat Striatum. *Mol. Brain Res.* **1992**, *14*, 186–195.
135. Fuxe, K.; Ferre, S.; Canals, M.; Torvinen, M.; Terasmaa, A.; Marcellino, D.; Goldberg, S.R.; Staines, W.; Jacobsen, K.X.; Lluís, C.; et al. Adenosine A(2A) and dopamine D-2 heteromeric receptor complexes and their function. *J. Mol. Neurosci.* **2005**, *26*, 209–219.
136. Hillion, J.; Canals, M.; Torvinen, M.; Casado, V.; Scott, R.; Terasmaa, A.; Hansson, A.; Watson, S.; Olah, M.E.; Mallol, J.; et al. Coaggregation, cointernalization, and codesensitization of adenosine A(2A) receptors and dopamine D-2 receptors. *J. Biol. Chem.* **2002**, *277*, 18091–18097.

137. Soriano, A.; Ventura, R.; Molero, A.; Hoen, R.; Casado, V.; Cortes, A.; Fanelli, F.; Albericio, F.; Lluís, C.; Franco, R.; et al. Adenosine A(2A) Receptor-Antagonist/Dopamine D-2 Receptor-Agonist Bivalent Ligands as Pharmacological Tools to Detect A(2A)-D-2 Receptor Heteromers. *J. Med. Chem.* **2009**, *52*, 5590–5602.
138. Riggelman, A.; Harvey, J.; Baylis, C. Endothelin mediates some of the renal actions of acutely administered angiotensin II. *Hypertension* **2001**, *38*, 105–109.
139. Zeng, C.Y.; Hopfer, U.; Asico, L.D.; Eisner, G.M.; Felder, R.A.; Jose, P.A. Altered AT(1) receptor regulation of ETB receptors in renal proximal tubule cells of spontaneously hypertensive rats. *Hypertension* **2005**, *46*, 926–931.
140. Zeng, C.Y.; Wang, Z.; Asico, L.D.; Hopfer, U.; Eisner, G.M.; Felder, R.A.; Jose, P.A. Aberrant ETB receptor regulation of AT(1) receptors in immortalized renal proximal tubule cells of spontaneously hypertensive rats. *Kidney Int.* **2005**, *68*, 623–631.
141. Ciruela, F.; Casado, V.; Rodrigues, R.J.; Lujan, R.; Burgueno, J.; Canals, M.; Borycz, J.; Rebola, N.; Goldberg, S.R.; Mallol, J.; et al. Presynaptic control of striatal glutamatergic neurotransmission by adenosine A(1)-A(2A) receptor heteromers. *J. Neurosci.* **2006**, *26*, 2080–2087.
142. Moreno, J.L.; Miranda-Azpiazu, P.; Garcia-Bea, A.; Younkin, J.; Cui, M.; Kozlenkov, A.; Ben-Ezra, A.; Voloudakis, G.; Fakira, A.K.; Baki, L.; et al. Allosteric signaling through an mGlu2 and 5-HT2A heteromeric receptor complex and its potential contribution to schizophrenia. *Sci. Signal.* **2016**, *9*, ra5, doi:10.1126/scisignal.aab0467.
143. Lopez-Gimenez, J.F.; Canals, M.; Padiani, J.D.; Milligan, G. The alpha(1b)-adrenoceptor exists as a higher-order oligomer: Effective oligomerization is required for receptor maturation, surface delivery, and function. *Mol. Pharmacol.* **2007**, *71*, 1015–1029.
144. Rodriguez-Frade, J.M.; del Real, G.; Serrano, A.; Hernanz-Falcon, P.; Soriano, S.F.; Vila-Coro, A.J.; de Ana, A.M.; Lucas, P.; Prieto, I.; Martinez-A, C.; et al. Blocking HIV-1 infection via CCR5 and CXCR4 receptors by acting in trans on the CCR2 chemokine receptor. *EMBO J.* **2004**, *23*, 66–76.
145. Sohy, D.; Parmentier, M.; Springael, J.Y. Allosteric transinhibition by specific antagonists in CCR2/CXCR4 heterodimers. *J. Biol. Chem.* **2007**, *282*, 30062–30069.
146. Sohy, D.; Yano, H.; de Nadai, P.; Urizar, E.; Guillabert, A.; Javitch, J.A.; Parmentier, M.; Springael, J.Y. Hetero-oligomerization of CCR2, CCR5, and CXCR4 and the Protean Effects of "Selective" Antagonists. *J. Biol. Chem.* **2009**, *284*, 31270–31279.
147. Grefen, C.; Blatt, M.R. A 2in1 cloning system enables ratiometric bimolecular fluorescence complementation (rBiFC). *BioTechniques* **2012**, *53*, 311–314.
148. White, C.W.; Vanyai, H.K.; See, H.B.; Johnstone, E.K.M.; Pflieger, K.D.G. Using nanoBRET and CRISPR/Cas9 to monitor proximity to a genome-edited protein in real-time. *Sci. Rep.* **2017**, *7*, 3187.
149. Zych, C.; Domling, A.; Ayyavoo, V. Development of a robust cell-based high-throughput screening assay to identify targets of HIV-1 viral protein R dimerization. *Drug Des. Dev. Ther.* **2013**, *7*, 403–412.
150. Morell, M.; Espargaro, A.; Aviles, F.X.; Ventura, S. Detection of transient protein-protein interactions by bimolecular fluorescence complementation: The Abl-SH3 case. *Proteomics* **2007**, *7*, 1023–1036.
151. Auld, D.S.; Inglese, J. Interferences with Luciferase Reporter Enzymes. In *Assay Guidance Manual*; Sittampalam, G.S., Coussens, N.P., Brimacombe, K., Grossman, A., Arkin, M., Auld, D., Austin, C., Baell, J., Bejcek, B., Chung, T.D.Y., et al., Eds.; Eli Lilly & Company and the National Center for Advancing Translational Sciences: Bethesda, MD, USA, 2004.
152. Simeonov, A.; Davis, M.I. Interference with Fluorescence and Absorbance. In *Assay Guidance Manual*; Sittampalam, G.S., Coussens, N.P., Brimacombe, K., Grossman, A., Arkin, M., Auld, D., Austin, C., Baell, J., Bejcek, B., Caaveiro, J.M.M., et al., Eds.; Eli Lilly & Company and the National Center for Advancing Translational Sciences: Bethesda, MD, USA, 2004.
153. Wade, M.; Mendez, J.; Coussens, N.P.; Arkin, M.R.; Glicksman, M.A. Inhibition of Protein-Protein Interactions: Cell-Based Assays. In *Assay Guidance Manual*; Sittampalam, G.S., Coussens, N.P., Brimacombe, K., Grossman, A., Arkin, M., Auld, D., Austin, C., Baell, J., Bejcek, B., Caaveiro, J.M.M., et al., Eds.; Eli Lilly & Company and the National Center for Advancing Translational Sciences: Bethesda, MD, USA, 2017.
154. Wu, J.C.; Sundaresan, G.; Iyer, M.; Gambhir, S.S. Noninvasive optical imaging of firefly luciferase reporter gene expression in skeletal muscles of living mice. *Mol. Ther. J. Am. Soc. Gene Ther.* **2001**, *4*, 297–306.

155. Han, Y.; Wang, S.F.; Zhang, Z.P.; Ma, X.H.; Li, W.; Zhang, X.W.; Deng, J.Y.; Wei, H.P.; Li, Z.Y.; Zhang, X.E.; et al. In vivo imaging of protein-protein and RNA-protein interactions using novel far-red fluorescence complementation systems. *Nucleic Acids Res.* **2014**, *42*, e103.
156. Faron-Gorecka, A.; Szlachta, M.; Kolasa, M.; Solich, J.; Gorecki, A.; Kusmider, M.; Zurawek, D.; Dziedzicka-Wasylewska, M. Understanding GPCR dimerization. *Methods Cell Biol.* **2019**, *149*, 155–178.
157. Maurice, P.; Guillaume, J.L.; Benleulmi-Chaachoua, A.; Daulat, A.M.; Kamal, M.; Jockers, R. GPCR-interacting proteins, major players of GPCR function. *Adv. Pharmacol. (San Diego, Calif.)* **2011**, *62*, 349–380.
158. Milligan, G.; Ward, R.J.; Marsango, S. GPCR homo-oligomerization. *Curr. Opin. Cell Biol.* **2018**, *57*, 40–47.
159. Gomes, I.; Ayoub, M.A.; Fujita, W.; Jaeger, W.C.; Pflieger, K.D.; Devi, L.A. G Protein-Coupled Receptor Heteromers. *Annu. Rev. Pharmacol. Toxicol.* **2016**, *56*, 403–425.



© 2019 by the authors. Licensee MDPI, Basel, Switzerland. This article is an open access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (<http://creativecommons.org/licenses/by/4.0/>).