

## Original Article

### Diagnostic Immunology



Ann Lab Med 2017;37:248-253  
<https://doi.org/10.3343/alm.2017.37.3.248>

ISSN 2234-3806 · eISSN 2234-3814

# ANNALS OF LABORATORY MEDICINE

## Comparison of the QIAGEN *artus* HBV QS-RGQ Assay With the Roche COBAS AmpliPrep/COBAS TaqMan HBV Assay for Quantifying Viral DNA in Sera of Chronic Hepatitis B Patients

Mi-Soon Han, M.D.<sup>1,3</sup>, Yongjung Park, M.D.<sup>2</sup>, Hyunjin Nah, M.D.<sup>3</sup>, and Hyon-Suk Kim, M.D.<sup>3</sup>

Medical Clinic Laboratory Department of U2Bio Co. Ltd.<sup>1</sup>, Seoul; Department of Laboratory Medicine<sup>2</sup>, National Health Insurance Service Ilsan Hospital, Goyang; Department of Laboratory Medicine<sup>3</sup>, Severance Hospital, Yonsei University College of Medicine, Seoul, Korea

**Background:** Hepatitis B virus DNA quantification is essential for managing chronic hepatitis B (CHB). We compared the performance of *artus* HBV QS-RGQ (QIAGEN GmbH, Germany) and CAP/CTM v2.0 HBV assays (Roche Molecular Diagnostics, USA) in CHB patients.

**Methods:** A comparative evaluation between two assays was performed with 508 clinical serum samples. Precision, linearity, and the limit of detection (LOD) of QS-RGQ assay was evaluated by using the WHO standard 97/750 and clinical samples.

**Results:** Detection rates and viral loads as determined QS-RGQ assay were significantly lower than those from the CAP/CTM v2.0 assay (52.8% vs 60.6%;  $3.55 \pm 1.77$  IU/mL vs  $4.18 \pm 1.89$  IU/mL,  $P < 0.0001$ ). The kappa coefficient between qualitative results was 0.79 (95% confidence interval, 0.74 to 0.85). Bland-Altman plot found a mean difference of (QS-RGQ - CAP/CTM v2.0) =  $-0.63 \log_{10}$  IU/mL (95% limit of agreement,  $-1.48$  to  $0.22$ ). Repeatability and total imprecision (% CV) of the QS-RGQ assay were 1.0% and 1.1% at 2,000 IU/mL, and 0.7% and 1.4% at 20,000 IU/mL, respectively. Linearity of this assay ranged from 31.6 to  $1.0 \pm 10^7$  IU/mL, and the LOD was 2.95 IU/mL.

**Conclusions:** The *artus* HBV QS-RGQ assay showed good performance but significantly decreased detection rate and viral load compared with CAP/CTM v2.0 assays. This assay recommends using plasma; however, we used stored serum because of the retrospective study design. Usually HBV DNA quantification is performed in plasma or serum, but sample type and clinical relevance of quantitative values should be considered when determining the clinical application of this reagent.

**Key Words:** Hepatitis B virus, HBV DNA quantification, *artus* HBV QS-RGQ assay, CAP/CTM v2.0 assay, Performance

**Received:** June 16, 2016

**Revision received:** September 27, 2016

**Accepted:** January 9, 2017

**Corresponding author:** Hyon-Suk Kim  
Department of Laboratory Medicine,  
Severance Hospital, Yonsei University  
College of Medicine, 50-1 Yonsei-ro,  
Seodaemun-gu, Seoul 03722, Korea  
Tel: +82-2-2228-2443  
Fax: +82-2-364-1583  
E-mail: kimhs54@yuhs.ac

### © Korean Society for Laboratory Medicine

This is an Open Access article distributed under the terms of the Creative Commons Attribution Non-Commercial License (<http://creativecommons.org/licenses/by-nc/4.0>) which permits unrestricted non-commercial use, distribution, and reproduction in any medium, provided the original work is properly cited.

## INTRODUCTION

According to the WHO's global estimates, more than 240 million people are chronically infected with Hepatitis B virus (HBV), and 780,000 people die every year as a consequence [1]. The

prevalence of chronic HBV infection varies geographically, but Southeast Asia is an area with high prevalence (8–15%) [2]. Detection of HBV DNA is essential for diagnosis, and HBV DNA quantification is a key determinant of treatment for both HBV envelope antigen (HBeAg)-negative or HBeAg-positive chronic

Hepatitis B (CHB) [3, 4]. Recent studies suggested that after completing successful treatment, remaining HBV DNA correlates with relapse and recurrence of HBV infection, even at a low concentration [5, 6].

The application of real-time PCR for viral diagnostics has been reported to exhibit high sensitivity, a broad dynamic range, and short turnaround time, and therefore is considered as the standard method for quantification [7, 8]. To date, many HBV DNA assays are commercially available. Although two assays may be comparable, assays may report discrepant viral load levels [9, 10]. Therefore, it is important to adopt a highly reliable PCR-based assay to quantify HBV DNA in order to enable appropriate clinical management of CHB. The QIAGEN *artus* HBV QS-RGQ assay (QIAGEN GmbH, Hilden, Germany) was newly launched but so far, only two abstracts have been published about its performance [11, 12].

Our study aimed to verify the performance of the *artus* HBV QS-RGQ assay compared with the Roche COBAS AmpliPrep/COBAS TaqMan HBV assay (Roche Molecular Diagnostics, Pleasanton, CA, USA) in clinical samples.

## METHODS

### 1. Samples

A total of 508 serum samples where HBV DNA quantification had been requested by clinicians and determined by CAP/CTM v2.0 assay were randomly collected from 2008 to 2014. The remaining samples were stored at  $-70^{\circ}\text{C}$  until examined by the *artus* QS-RGQ assay. This study was performed according to the World Medical Association's Declaration of Helsinki guidelines and waived from the institutional review board of the Severance Hospital.

### 2. HBV DNA quantification assays

1) The COBAS AmpliPrep/COBAS TaqMan HBV assay (CAP/CTM v2.0) The CAP/CTM v2.0 is an automated real-time PCR assay that targets the precore and core regions of the HBV genome. The HBV DNA preparation with the COBAS AmpliPrep instrument requires 650  $\mu\text{L}$  of serum or plasma. DNA is extracted and eluted in a volume of 65  $\mu\text{L}$ , 50  $\mu\text{L}$  of which is then analyzed by PCR. Real-time PCR is performed by using the COBAS TaqMan 96 analyzer with a multiplex TaqMan assay. Two targets are amplified, HBV DNA and the internal quantitation standard. All procedures were performed according to the manufacturer's protocol. HBV DNA levels were expressed in international units per milliliter (IU/mL) with a conversion factor of 5.82 copies of HBV

DNA per IU. The lower limit of detection (LOD) provided by the manufacturer is 20.0 IU/mL, and the dynamic range of quantification is from 20.0 to  $1.7 \times 10^8$  IU/mL ( $1.3\text{--}8.2 \log_{10}$  IU/mL). This assay is capable of analyzing 24 specimens in about two hours [13-15].

### 2) The *artus* HBV QS-RGQ assay (*artus* QS-RGQ)

The *artus* QS-RGQ kit requires an input of 1,000  $\mu\text{L}$  of plasma sample. In this study, clinical serum samples were used instead of plasma. Nucleic acids are extracted and eluted in a final volume of 60  $\mu\text{L}$  in elution buffer. Total nucleic acids were extracted with the QIASymphony SP using the QIASymphony DSP Virus/Pathogen Midi kits. Samples processed on the QIASymphony SP can be transferred automatically to the QIASymphony AS module (integrated operation) for assay setup. The QIASymphony AS sets up the PCR reaction by mixing 30  $\mu\text{L}$  of master mix and 20  $\mu\text{L}$  of DNA template. The master mix contains reagents and enzymes for specific amplification of a 134-bp region of the HBV core gene and for direct detection of the amplicon. The real-time PCR was performed on Rotor-Gene Q instruments under conditions described in the QIAGEN *artus* HBV QS-RGQ kit handbook. Data were analyzed with the Rotor-Gene Q software version 2.02 using thresholds of 0.04 and 0.03 to detect signals from HBV and the internal control, respectively. Quantification of HBV DNA was determined by using 5-point external standards. The 1 IU/mL corresponds to 8.21 copies/mL for HBV DNA detection on the Rotor-Gene Q. The LOD claimed by the manufacturer is 10.0 IU/mL, and the linear range of quantification is from 31.6 to  $2.0 \times 10^7$  IU/mL ( $1.5\text{--}7.3 \log_{10}$  IU/mL) (Table 1).

**Table 1.** Comparison between the characteristics of two Hepatitis B virus (HBV) DNA quantification assays

Variables	CAP/CTM v2.0	<i>artus</i> QS-RGQ
DNA extraction principle	Magnetic particle	Magnetic particle
Sample volume ( $\mu\text{L}$ )	650	1,200
Processed volume ( $\mu\text{L}$ )	500	1,000
Elution volume ( $\mu\text{L}$ )	65	60
Sample capacity per batch	24	24
Extraction runtime (min)*	120	90
Target HBV genome region	Precore and Core	Core
Claimed low limit of detection (IU/mL)	20 <sup>†</sup>	10
Claimed dynamic range (IU/mL)	$2.00 \times 10^1$ to $1.70 \times 10^8$	$3.16 \times 10^1$ to $2.00 \times 10^7$

\*Including the lysis step but not hands-on time; <sup>†</sup>In plasma (vs 12 IU/mL in serum).

### 3. Interpretative criteria of HBV DNA viral load

The qualitative results were interpreted as “Detected” or “Not detected” on the basis of each claimed LOD. The quantitative results were included for further analysis, only if exceeding the limit of quantification (LOQ), i.e. above 20.0 IU/mL (CAP/CTM v2.0) and 31.6 IU/mL (*artus* QS-RGQ).

### 4. Serologic CHB assay

The HBeAg status was simultaneously evaluated by the Architect i2000SR analyzer (Abbott Diagnostics, Abbott Park, IL, USA) using a chemiluminescent microparticle immunoassay in clinical samples before frozen storage. According to HBeAg-negative or -positive status, the differences in viral loads detected by the assays were determined.

### 5. Performance characteristics of the *artus* QS-RGQ assay

The second WHO international standard for HBV DNA (NIBSC code 97/750,  $1 \times 10^6$  IU/mL) and clinical samples with high viral load were serially diluted with human serum matrix. The precision performances near the two clinically important levels (low: 2,000 IU/mL; high: 20,000 IU/mL) were evaluated in four replicates over five days using clinical samples. To evaluate linearity, serial dilutions of the WHO standard (20 to  $10^5$  IU/mL) and clinical samples ( $10$  to  $10^7$  IU/mL) were simultaneously tested in two to four replicate measurements. The diluted WHO standard materials of 1, 2.5, 5, 10, and 20 IU/mL were tested in 12 replicates to determine the LOD.

### 6. Statistical analysis

Data was analyzed by using Analyse-it Method Evaluation Edition, version 3.76 software (Analyse-it Software Ltd., Leeds, UK) and IBM SPSS Statistics 22 (International Business Machines Corp., Armonk, NY, USA). The method comparisons were performed by the Kappa test and Pearson Chi-squared test for qualitative results, and the Spearman's test and Passing-Bablok regression for quantitative results. The differences between two methods were presented in a Bland-Altman plot. Linear regression analysis and Probit analysis were used for verification of the *artus* QS-RGQ assay performance. The *P* values  $<0.05$  were considered significant.

## RESULTS

### 1. Agreement and correlation between the two assays

Among the 508 serum samples tested, HBV DNA was quantified in 308 samples (60.6%) by the CAP/CTM v2.0 assay, and

268 (52.8%) were quantified by the *artus* QS-RGQ assay. The detection rate by the *artus* QS-RGQ assay was significantly lower than that by the CAP/CTM v2.0 assay ( $P < 0.0001$ ). The kappa coefficient between the qualitative results of the two assays was 0.79 (95% confidence interval [CI]=0.74 to 0.85). By both assays (Table 2), 229 (45.1%) samples were detected and 196 (38.6%) samples were not detected. However, the CAP/CTM v2.0 assay detected HBV viral load in 46 (9.1%) samples whereas the *artus* QS-RGQ did not detect HBV viral load. Conversely, the *artus* QS-RGQ assay detected six (1.2%) positive samples but not CAP/CTM v2.0. Overall, 12 samples with  $\geq 20,000$  IU/mL and 29 with 2,000 to 20,000 IU/mL, based on the result from CAP/CTM v2.0 assay, were quantified as lower viral loads by the *artus* QS-RGQ assay. Inversely, one sample (4,230 IU/mL) among the results from *artus* QS-RGQ assay was measured as a lower level by the CAP/CTM v2.0 assay (755 IU/mL).

The mean  $\pm$  SD of HBV DNA was  $4.18 \pm 1.89$  IU/mL for the CAP/CTM v2.0 assay and  $3.55 \pm 1.77$  IU/mL for the *artus* QS-RGQ assay. The viral loads quantified by *artus* QS-RGQ assay were significantly lower than those quantified by the CAP/CTM v2.0 assay ( $P < 0.0001$ ). Further, 38 (7.5%) samples showed significant difference of  $\geq 1 \log_{10}$  IU/mL.

Both assays quantified 105 (20.7%) samples with an HBV DNA level  $< 2,000$  IU/mL ( $3.30 \log_{10}$  IU/mL) and showed fair correlation ( $r = 0.49$ ; 95% CI, 0.33 to 0.63;  $P < 0.0001$ ). Of samples detected by both assays, Passing-Bablok regression analysis and Bland-Altman plot are shown in Fig. 1A and Fig. 1B, respectively.

Of the 508 study samples, 177 samples were evaluated by the HBeAg assay, which found 106 HBeAg-negative and 71 HBeAg-

**Table 2.** Comparison of viral load results in 508 clinical samples using the CAP/CTM v2.0 and *artus* QS-RGQ assays (N)

<i>artus</i> QS-RGQ (IU/mL)	CAP/CTM v2.0 (IU/mL)				Total
	$\geq 20,000$	$\geq 2,000$ to $< 20,000$	$\geq 20$ to $< 2,000$	$< 20$	
$\geq 20,000$	69	0	0	0	69
$\geq 2,000$ to $< 20,000$	9*	15	1 <sup>§</sup>	0	25
$\geq 32$ to $< 2,000$	3 <sup>†</sup>	28 <sup>‡</sup>	104	4	139
$< 32$ , Detected	0	0	33	2	35
$< 10$ , Not Detected	0	1	45	196	240
Total	81	44	183	200	508

Hepatitis B virus DNA levels described below are shown as median (1st to 3rd quartiles) (IU/mL) obtained by the *artus* QS-RGQ and CAP/CTM v2.0, respectively: \*7,280 (3,324–10,500) vs 92,189 (34,467–145,000), <sup>†</sup>1,602 (1,482–1,700) vs 41,200 (29,133–57,717), <sup>‡</sup>812 (437–1,124) vs 6,162 (2,752–6,848), <sup>§</sup>4,230 vs 755.

positive samples. The CAP/CTM assay detected 63 samples (59.4%) and 53 samples (74.6%), meanwhile the *artus* QS-RGQ assay detected 53 samples (50.0%) and 50 samples (70.4%) that were HBeAg-negative and HBeAg-positive, respectively. The tendency for viral load difference between the two assays was consistent regardless of HBeAg negativity or positivity. However, the viral load (mean±SD) in the HBeAg-negative group was lower than that in the HBeAg-positive group in both assays ( $3.97 \pm 0.25$  vs  $5.25 \pm 0.34$  for CAP/CTM assay,  $P < 0.0001$ ;  $3.35 \pm 0.24$  vs  $4.60 \pm 0.31$  for *artus* QS-RGQ assay,  $P < 0.0001$ ).

## 2. Performance characteristics of the *artus* QS-RGQ assay


The mean levels and total % CV for the *artus* QS-RGQ assay at two clinically important HBV DNA levels (target levels), i.e. 3.30  $\log_{10}$  IU/mL (2,000 IU/mL) and 4.30  $\log_{10}$  IU/mL (20,000 IU/mL),

were 3.29  $\log_{10}$  IU/mL and 4.31  $\log_{10}$  IU/mL, and 1.1% and 1.4%, respectively (Table 3). The linearity of the *artus* QS-RGQ assay was verified in the claimed linear range for each HBV DNA level of the WHO standard (ranged 1.30 to 5.00  $\log_{10}$  IU/mL,  $y = 0.94x + 0.28$ ,  $r = 0.995$ ) and clinical samples (ranged 1.00 to 7.00  $\log_{10}$  IU/mL,  $y = 0.97x + 0.10$ ,  $r = 0.997$ ) (Table 4). Unlike results for the WHO standard, the *artus* QS-RGQ assay showed less than

**Table 3.** Precision performance of the *artus* QS-RGQ assay

Target concentration ( $\log_{10}$ IU/mL)	Mean level ( $\log_{10}$ IU/mL)	SD (% CV)		
		Repeatability	Between-day	Total
3.30*	3.29	0.03 (1.0%)	0.01 (0.3%)	0.04 (1.1%)
4.30†	4.31	0.03 (0.7%)	0.05 (1.2%)	0.06 (1.4%)

\*2,000 IU/mL; †20,000 IU/mL.

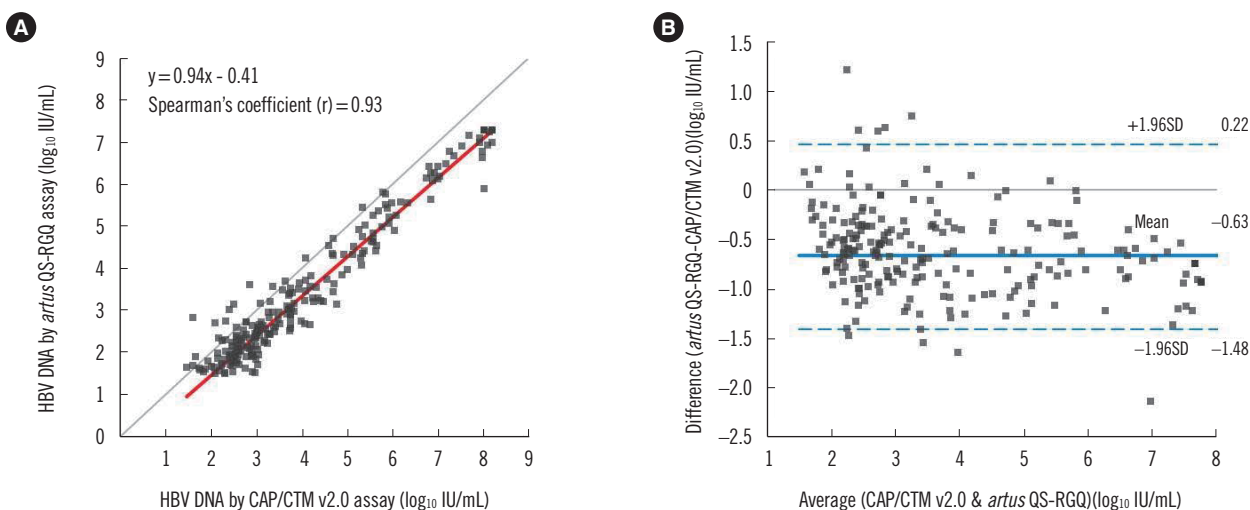

**Fig. 1.** Comparison of Hepatitis B virus (HBV) DNA levels in clinical samples between CAP/CTM v2.0 and *artus* QS-RGQ assays (N=229). (A) The Passing-Bablok equation (95% CI of intercept: -0.52 to -0.27, slope: 0.92 to 0.0.96,  $r = 0.91$  to 0.95). (B) The mean difference (QS/RGQ-CAP/CTM) in Bland-Altman plot was  $-0.63 \pm 0.85 \log_{10}$  IU/mL.

**Table 4.** Linearity of the *artus* QS-RGQ assay for WHO standard and clinical samples

Assigned concentration ( $\log_{10}$ IU/mL)	WHO standard					Clinical samples*				
	Replicates N	Mean	SD	CV (%)	Recovery %	Replicates N	Mean	SD	CV (%)	Recovery %
7.00						3	6.93	0.02	0.32	98.99
6.00						3	5.95	0.01	0.10	99.14
5.00	2	4.87	0.06	1.19	97.50	3	4.98	0.03	0.50	99.52
4.00	4	4.07	0.03	0.64	101.72	2	3.93	0.03	0.87	98.27
3.00	4	3.11	0.02	0.49	103.58	2	2.93	0.03	1.10	97.61
2.00	4	2.07	0.07	3.40	103.69	2	2.06	0.25	12.17	102.97
1.30	4	1.52	0.25	16.52	116.90					
1.00						2	1.15	0.49	43.02	114.88

\*Concentrations measured by CAP/CTM v2.0 assay.

**Table 5.** Lower limit of detection\* of the *artus* QS-RGQ assay

Intended concentration (IU/mL)	Replicate N	Detected N	% Detected	Probability (%)
20.0	12	12	100.0	100.0
10.0	12	12	100.0	100.0
5.0	12	12	100.0	100.0
2.5	12	11	91.7	91.8
1.0	12	7	58.3	58.5

\*The concentration with 95% probability of detection was 2.95 IU/mL according to probit analysis.

100% recovery at  $>10^2$  IU/mL viral loads and increased % CV at  $2.00 \log_{10}$  IU/mL in clinical samples. The probit analysis for the *artus* QS-RGQ assay resulted in an LOD of 2.95 IU/mL (Table 5).

## DISCUSSION

CHB is a major cause of cirrhosis and hepatocellular carcinoma in Southeast Asia, China, and Africa [16]. HBV DNA quantification is essential for monitoring disease status and treatment response in CHB patients. Until now, several commercially available real-time PCR assays for this have been developed. Newly launched *artus* QS-RGQ assay suggested a lower LOD, higher sample volume, and higher LOQ than CAP/CTM v2.0 assay by manufacturer. This new assay has not been researched enough about its performance, except for two abstracts with plasma samples [11, 12]. Brichler *et al* [11] reported good correlation and agreement ( $r^2=0.89$ , mean difference= $0.1 \log_{10}$  IU/mL) between the new assay and CAP/CTM v2.0 assay. However, 23% of the 230 samples showed more than  $\pm 0.5 \log_{10}$  IU/mL difference range, and the HBV DNA levels by the *artus* QS-RGQ assay were lower than the results by the CAP/CTM v2.0 assay in their study (*P* value is not shown). Fielder *et al* [12] suggested that new assay showed detection capability of all eight genotypes, no cross-reactivity with 30 different pathogens and low LOD of 4.1 IU/mL.

In our study, there were significant differences between the CAP/CTM v2.0 and the *artus* QS-RGQ assays in the detection rate and viral load when quantifying HBV DNA levels in clinical serum samples, even though good correlation was observed ( $r=0.93$ ). The HBV DNA levels determined by the *artus* QS-RGQ assay were substantially lower than the results by the CAP/CTM v2.0 assay. In addition, the correlation between the results of the two assays in low viral load samples ( $<2,000$  IU/mL) was not as strong as the correlation for all samples ( $r=0.49$ ). According to the manufacturer's protocol, the LOD of the *artus* QS-RGQ

assay (10 IU/mL) is lower than that of the CAP/CTM v2.0 assay (20 IU/mL). However, our results from the *artus* QS-RGQ assay showed a lower positive rate than those from the CAP/CTM v2.0. Hence, discrepancies between the assays were observed when sample results were plotted against important HBV DNA target levels for clinical management points and each claimed LOD and LOQ (Table 2). These discrepancies may result in different classifications of inactive HBsAg carriers and active CHB patients, and may even change CHB management strategy.

There are several possible factors that may lead to assay discrepancies. Yeh *et al* [10] showed that the B genotype and low HBV viral load were two factors that contribute to significant differences in HBV DNA viral load detection, i.e.  $\geq 1 \log_{10}$  IU/mL, when comparing the CAP/CTM v2.0 assay with the Real-Time HBV assay (Abbott Laboratories, Abbott Park, IL, USA). When Ismail *et al* [9] evaluated the Abbott RealTime HBV assay and *artus* assays on two different platforms, the two *artus* assays showed low quantitative values in comparison with the RealTime assay, especially in samples with low viral load. In our study, the *artus* QS-RGQ assay also determined relatively lower viral loads than other assays, especially in samples with low viral load ( $\leq 2,000$  IU/mL).

Mutations in precore and core promoter regions may occur as CHB progresses [17]. These mutations may influence the result of HBV DNA quantification [18, 19]. Yeh *et al* [19] concluded that the detection difference between the CAP/CTM v2.0 and Abbott Real-Time assays in low viral load samples may be influenced by the tyrosine-methionine-aspartate-aspartate (YMDD) mutation, which confers lamivudine resistance. Because of sample volume limitations, we could not genotype our samples, but almost all HBV genotypes are reported as genotype C in Korea, and mutations in the core promoter are more prevalent than mutations in the precore region. Moreover, two mutations are often associated with HBeAg-negative results or reduced HBeAg production [20, 21]. Thus, detection sensitivity for HBV mutants in the common target region may differ depending on the assay. However, no additional study was performed to identify HBV genotype and mutation of precore and core regions due to insufficient sample volume and the retrospective study design that utilized reserved samples. Finally, the *artus* QS-RGQ assay showed less reliable performance in clinical serum samples than in the WHO standard samples. This discrepancy may be due to the sample type. The QS-RGQ assay recommends the use of plasma samples only, while other manufacturers adopt both sample types, i.e. plasma or serum.

Our study has limitations because this study examined stored

clinical samples to evaluate the newly launched reagents in a clinical laboratory. A well-designed prospective study with freshly drawn blood samples is preferable according to the manufacturer; however, in this case, we retrospectively studied stored samples to determine whether to adopt a clinical test.

In conclusion, the newly launched *artus* QS-RGQ assay displayed good precision at important HBV DNA target levels, linearity over clinically significant ranges, and an acceptable LOD. However, this assay showed a significantly decreased detection rate and viral load compared with results of the CAP/CTM v2.0 assay in serum samples, especially in samples with low viral load. However, this assay recommends plasma only as the sample type. Thus, careful consideration of sample type and sample volume is necessary when evaluating the relevance of these results to determine the clinical application of the QS-RGQ assay.

## Authors' Disclosures of Potential Conflicts of Interest

No potential conflicts of interest relevant to this article were reported.

## Acknowledgments

Miran Kim is acknowledged for their expert technical assistance. The reagents and consumables were provided by QIAGEN Korea. QIAGEN did not influence how the study was conducted or the approval of the manuscript.

## REFERENCES

- World Health Organization. Hepatitis B. <http://www.who.int/mediacentre/factsheets/fs204/en/> (Updated on Jul 2016).
- Oakes K. Hepatitis B: prevalence and pathophysiology. *Nurs Times* 2014; 110:12-6.
- Lok AS and McMahon BJ. Chronic hepatitis B: update 2009. *Hepatology* 2009;50:661-2.
- Martin P, Lau DT, Nguyen MH, Janssen HL, Dieterich DT, Peters MG, et al. A treatment algorithm for the management of chronic hepatitis B virus infection in the United States: 2015 Update. *Clin Gastroenterol Hepatol* 2015;13:2071-87.
- González R, Torres P, Castro E, Barbolla L, Candotti D, Koppelman M, et al. Efficacy of hepatitis B virus (HBV) DNA screening and characterization of acute and occult HBV infections among blood donors from Madrid, Spain. *Transfusion* 2010;50:221-30.
- Jiang JN, Huang ZL, He LX, Huang YH, Su MH, Xie R, et al. Residual amount of HBV DNA in serum is related to relapse in chronic hepatitis B patients after cessation of nucleos(t)ide analogs. *J Clin Gastroenterol* 2015;49:323-8.
- Stelzl E, Muller Z, Marth E, Kessler HH. Rapid quantification of hepatitis B virus DNA by automated sample preparation and real-time PCR. *J Clin Microbiol* 2004;42:2445-9.
- Poljak M, Lepej SZ, Rode OD. Recent developments in serologic and molecular diagnosis of hepatitis B and C. *Acta medica Croatica* 2013; 67:281-90.
- Ismail AM, Sivakumar J, Anantharam R, Dayalan S, Samuel P, Fletcher GJ, et al. Performance characteristics and comparison of Abbott and artus real-time systems for hepatitis B virus DNA quantification. *J Clin Microbiol* 2011;49:3215-21.
- Yeh ML, Huang CF, Hsieh MY, Huang JF, Dai CY, Yu ML, et al. Comparison of the Abbott RealTime HBV assay with the Roche Cobas AmpliPrep/Cobas TaqMan HBV assay for HBV DNA detection and quantification. *J Clin Virol* 2014;60:206-14.
- Brichler S, Dziri S, Duhant A, Faleu G, Gordien E. P0543: Comparative evaluation of the artus HBV QS-RGQ assay and the roche cobas mpli-prep/cobas taqman V.2 HBV test for the quantification of HBV DNA in plasma. *J Hepatol* 2015;62:S518-9.
- Fielder H, Javed S, Marques C, Wall G, Lange C, Tomlinson D. Quantitative detection of hepatitis B virus (HBV) DNA with the novel artus® HBV QS-RGQ Kit version 2. *J Clin Virol* 2015;70:S85.
- Alice T, Cerutti F, Pittaluga F, Varetto S, Gabella S, Marzano A, et al. COBAS AmpliPrep-COBAS TaqMan hepatitis B virus (HBV) test: a novel automated real-time PCR assay for quantification of HBV DNA in plasma. *J Clin Microbiol* 2007;45:828-34.
- Chevaliez S, Bouvier-Alias M, Laperche S, Hézode C, Pawlotsky J.M. Performance of version 2.0 of the Cobas AmpliPrep/Cobas TaqMan real-time PCR assay for hepatitis B virus DNA quantification. *J Clin Microbiol* 2010;48:3641-7.
- Caliendo AM, Valsamakis A, Bremer JW, Ferreira-Gonzalez A, Granger S, Sabatini L, et al. Multilaboratory evaluation of real-time PCR tests for hepatitis B virus DNA quantification. *J Clin Microbiol* 2011;49:2854-8.
- Lok AS. Natural history and control of perinatally acquired hepatitis B virus infection. *Dig Dis* 1992;10:46-52.
- Funk ML, Rosenberg DM, Lok AS. World-wide epidemiology of HBeAg-negative chronic hepatitis B and associated precore and core promoter variants. *J Viral Hepat* 2002;9:52-61.
- Lindh M, Hannoun C, Malmström S, Lindberg J, Norkrans G. Lamivudine resistance of hepatitis B virus masked by coemergence of mutations in probe region of the COBAS AMPLICOR assay. *J Clin Microbiol* 2006;44:2587-9.
- Yeh ML, Huang CF, Huang CI, Liu SF, Yang HL, Hsieh MY, et al. Abbott RealTime HBV assay is more sensitive in detection of low viral load and little impacted by drug resistant mutation in chronic hepatitis B patients under nucleot(s)ide analogues therapy. *PLoS One* 2014;9:e101790.
- Alestig E, Hannoun C, Horal P, Lindh M. Phylogenetic origin of hepatitis B virus strains with precore C-1858 variant. *J Clin Microbiol* 2001;39: 3200-3.
- Alexopoulou A and Karayiannis P. HBeAg negative variants and their role in the natural history of chronic hepatitis B virus infection. *World J Gastroenterol* 2014;20:7644-52.