

PROCTOR, K., PETRIE, B., BARDEN, R., ARNOT, T. and KASPRZYK-HORDERN, B. 2019. Multi-residue ultra-performance liquid chromatography coupled with tandem mass spectrometry method for comprehensive multi-class anthropogenic compounds of emerging concern analysis in a catchment-based exposure-driven study. *Analytical and bioanalytical chemistry* [online], Online First. Available from: <https://doi.org/10.1007/s00216-019-02091-8>.

Multi-residue ultra-performance liquid chromatography coupled with tandem mass spectrometry method for comprehensive multi-class anthropogenic compounds of emerging concern analysis in a catchment-based exposure-driven study.

PROCTOR, K., PETRIE, B., BARDEN, R., ARNOT, T.,
KASPRZYK-HORDERN, B.

2019





Multi-residue ultra-performance liquid chromatography coupled with tandem mass spectrometry method for comprehensive multi-class anthropogenic compounds of emerging concern analysis in a catchment-based exposure-driven study

Kathryn Proctor^{1,2} · Bruce Petrie^{1,3} · Ruth Barden⁴ · Tom Arnot² · Barbara Kasprzyk-Hordern^{1,2}

Received: 21 June 2019 / Revised: 26 July 2019 / Accepted: 13 August 2019
© The Author(s) 2019

Abstract

This paper presents a new multi-residue method for the quantification of more than 142 anthropogenic compounds of emerging concern (CECs) in various environmental matrices. These CECs are from a wide range of major classes including pharmaceuticals, household, industrial and agricultural. This method utilises ultra-performance liquid chromatography coupled with tandem mass spectrometry (UPLC-MS/MS) for analysis of five matrices (three liquid and two solid) from wastewater treatment processes and the surrounding environment. Relative recoveries were predominantly between 80 and 120%; however, due to the complexity of the matrices used in this work, not all compounds were recovered in all matrices, from 138/142 analytes in surface water to 96/142 analytes in digested solids. Method quantification limits (MQLs) ranged from 0.004 ng L⁻¹ (bisoprolol in surface water) to 3118 ng L⁻¹ (creatinine in wastewater treatment work (WwTW) influent). The overall method accuracy was 107.0%, and precision was 13.4%. To test its performance, the method was applied to the range of environmental matrices at WwTWs in South West England. Overall, this method was found to be suitable for application in catchment-based exposure-driven studies, as, of the total number of analytes quantifiable in each matrix, 61% on average was found to be above their corresponding MQL. The results confirm the need for analysing both the liquid and solid compartments within a WwTW to prevent under-reporting of concentrations.

Keywords Chemicals of emerging concern · Analysis · Environment · Mass spectrometry

Electronic supplementary material The online version of this article (<https://doi.org/10.1007/s00216-019-02091-8>) contains supplementary material, which is available to authorized users.

✉ Barbara Kasprzyk-Hordern
B.Kasprzyk-Hordern@bath.ac.uk

¹ Department of Chemistry, University of Bath, Claverton Down, Bath BA2 7AY, UK

² Water Innovation & Research Centre, Department of Chemical Engineering, University of Bath, Claverton Down, Bath BA2 7AY, UK

³ School of Pharmacy and Life Sciences, Robert Gordon University, Aberdeen AB10 7JG, UK

⁴ Wessex Water Services Ltd., Claverton Down, Bath BA2 7WW, UK

Introduction

The use of anthropogenic, household, industrial or agricultural chemicals such as pharmaceuticals, pesticides, plasticisers, UV filters, industrial chemicals and microplastics is ubiquitous, and they have been recognised as a source of environmental contamination. These compounds have been quantified at levels ranging from ng L⁻¹ to µg L⁻¹, and their impact on the environment is not well known. These compounds are often designed to be biologically active and can be persistent in the environment, where they have the potential to bioaccumulate within the tissues of organisms [1–4]. For these reasons, among many others, they are known as compounds of emerging concern (CECs) [5–7].

CECs are primarily introduced to the environment via point sources such as wastewater treatment works (WwTWs), industrial discharge points and landfill leachates [8–11]. Diffuse sources, such as direct application to land in agriculture, have

been identified as a source of select sub-classes of CECs, such as pesticides and veterinary pharmaceuticals. Additionally, the application of digested sludge, from WwTW processes, directly to the land in modern farming practices is a potential source of other classes of CECs [12]. Bisphenol A (BPA) for instance, has been found at high levels within this matrix [13, 14].

CECs have been detected across the world in a multitude of environmental matrices [6, 15, 16]. This is due to their widespread use and to partitioning that can occur from the aqueous phases into suspended solids and sediments, where it can affect terrestrial organisms and fauna. Whilst their presence does not necessarily mean harm, the ecotoxicological effects of many of the CECs have been quantified in laboratory-based studies for a variety of different organisms across trophic levels and toxic effects have been demonstrated [17–19]. These laboratory-based ecotoxicological studies broadly focus on a single compound versus a single organism, but the environment is a ‘cocktail’ of CECs and different microorganisms. Ermler et al. [20] addressed this lack of knowledge for anti-androgens and found that concentration addition is a good model for predicting the effect of mixtures for up to six compounds. Leading on from this work, Orton et al. [21] tested multi-component mixtures of up to 30 compounds with varying mixture ratios. At the point where one mixture caused a 10% inhibition of the cancer cell assay, the concentrations of the components were a factor of 5.8 lower than the concentration that would be needed for them to individually cause this effect [21]. This highlights the need for current analytical methods to have method quantification limits (MQLs) lower than the no-observed effect concentrations (NOECs) of individual CECs, to enable accurate measurements of very low concentrations, for a better understanding of the risk they pose to the environment in combination.

Work is being done to further understand the fate, behaviour and effect of CECs within our environment. However, the sheer number of them and their everchanging usage makes this a challenge. The European Inventory of Existing Commercial Chemical Substances (EINECS) contains over 100,000 substances [22], with further substances being registered across all EEA countries, through Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH). Currently, there are 94,705 registrations containing 22,257 substances, or 22,096 unique substances that are in use [23]. This number is increasing every year, as new substances are developed and registered, and this pattern can be found across all classes of CECs.

The exponential growth of populations across the world, due to increased life expectancy and decreased infant mortality, is increasing the usage of pharmaceuticals. This puts further pressure on agriculture to produce food faster and cheaper, often via the use of pesticides, herbicides and other anthropogenic compounds. With the number of CECs in use, and more being developed each year, it is not feasible to determine the exposure

and effects of all these compounds in a single catchment, let alone across a country or continent. To further complicate matters, many CECs degrade through different processes such as hydrolysis, photolysis and metabolism. These form transformation products and metabolites, which are often more harmful than the parent compounds [24].

Brack et al. [25] conducted a review of the EU Water Framework Directive (WFD) and concluded that there are specific challenges both at European scale and at a local scale. Regulation and national monitoring schemes such as the WFD and the UK Chemical Investigations Programme (CIP) are suitable for furthering the understanding of this problem on a wider scale. Through identification and assessment of the most widespread CECs, high-risk compounds can be identified and managed through harmonised methods across diverse areas. However, they are often limited by the sampling method, i.e. a few samples collected at many locations, using only grab samples, which are not very representative [26]. Sampling at a local scale is crucial to determine catchment-specific substances and mixtures that might be a specific problem in the local environment. This can allow more targeted management of risks and hazards at local catchment level. There is a further need for analytical methods which can detect CECs down to, or below, the predicted no-effect concentration (PNEC), as this will allow a more adequate risk assessment [27].

To gain a better understanding of the fate of CECs within an environmental catchment, analysis of the influent and effluent of contributing WwTWs is required, including analysis of solid particulate matter (SPM) and digested solids, as well as corresponding surface water. Overall, in the literature, there is still a lack of multi-residue methods for quantification of CECs in solid matrices. Even fewer publications consider both the liquid and solid phases in the WwTW or the environment. Many studies focus specifically on pharmaceuticals, personal care products, industrial chemicals, veterinary pharmaceuticals or pesticides. However, a catchment potentially contains compounds from a variety, if not all, classes of CECs, though it is rare that they are all the focus of analysis in a single campaign. One study by Gustavsson et al. [28] covers many different classes of CEC, leading to the analysis of 172 compounds; however, this was only achieved through multiple sample preparation methods. For a comprehensive understanding of the exposure and environmental risk of chemical mixtures, multi-residue quantitative methods covering a large variety of CECs are required.

The aim of this work is to develop and validate a new multi-residue method (< 190 compounds including internal standards) for a wide range of CECs prioritised for risk assessment at a catchment scale, and accounting for the highly urbanised and agricultural areas of one catchment. The classes of CECs covered by this method include the following: UV filters, parabens, plasticisers, steroid estrogens, antibacterials/antibiotics, antifungals, hypertension drugs, non-steroidal anti-inflammatory drugs

(NSAIDs), lipid regulators, antihyperlipidaemics, antihypertensives, antihistamines, drugs for erectile dysfunction, drugs for diabetes, cough suppressants, beta blockers, H₂ receptor agonists, X-ray contrast media, drug precursors, anticancer drugs, anaesthetics, antidepressants, anti-epileptics, calcium channel blockers, hypnotics, antipsychotics, drugs for dementia, human indicators, analgesics, stimulants, opioids, drugs used in veterinary medicine, pesticides, fungicides and herbicides and metabolites. The selection of analyte groups was based not only on prioritisation including existing and proposed legislation, European and national watch lists (UKWIR CIP, EU Watch List) [29–32] and the literature [33, 34] but also on exploring usage statistics (NHS prescriptions [35]), entry into the environment (metabolism, excretion from DrugBank [36]), persistence, bioaccumulation, transport throughout the environment and toxicity of organisms [mammals, aquatic and benthic (log K_{ow} , log K_{oc} , log D_{ow} , water solubility, vapour pressure, Henry's law constant, bioconcentration factor, EPI Suite, ACD/Labs [37, 38]).

Materials and methods

This paper provides an expanded and broader scope method based on the method published by Petrie et al. [39]. Electronic Supplementary Material (ESM) Table S1 contains data on the suppliers of the compounds as well as their physiochemical properties.

The analytes were primarily purchased in solid form, before being accurately weighed and dissolved in HPLC grade methanol (MeOH) (Sigma-Aldrich), or other suitable solvents, at a concentration of 1.0 mg mL⁻¹. These stock solutions were stored in silanised glass vials in the dark at -20 °C, unless otherwise stated. Multi-analyte mixtures were prepared from these stock solutions. The CECs, their corresponding internal standards, data acquisition method and MS/MS detection parameters can be found in ESM Table S2. Chromatograms of all analytes can be found in ESM Fig. S1.

The internal standards 1*S*,2*R*-(+) ephedrine-d3, amphetamine-d5, benzoylecgonine-d8, cocaethylene-d8, cocaine-d3, codeine-d6, cotinine-d3, 2-ethylidene-1,5-dimethyl-3,3-diphenylpyrrolidine-d3 (EDDP-d3), estradiol (2,4,16,16-d4), estrone (2,4,16,16-d4), heroin-d9, ketamine-d4, 3,4-methylenedioxyamphetamine-d5 (MDA-d5), 3,4-methylenedioxymethamphetamine (MDMA-d5), mephedrone-d3, methadone-d9, methamphetamine-d5, methylparaben-13C, morphine-d3, norketamine-d4, quetiapine-d8 and tempazepam-d5 were purchased from LGC Standards (Middlesex, UK). Amitriptyline-d3, amoxicillin-d4, capecitabine-d11, ciprofloxacin-d8, citalopram-d6, diazepam-d5, erythromycin-13C,d3, fluoxetine-d5, gabapentin-d4, imidacloprid-d4, metazachlor-d6, metoprolol-d7, mirtazapine-d3, nortriptyline-d4, nortriptyline-d3, ofloxacin-d3, oxazepam-d5, sildenafil-d8 and verapamil-d7 were purchased from Toronto Research Chemicals (TRC)

(Toronto, Canada). Acetaminophen-d4, atenolol-d5, bisphenol A-d16, carbamazepine-13C6, ibuprofen-d3, ketoprofen-d3, metformin (dimethyl-d6), methiocarb-d3, naproxen-d3, propranolol-d7, sertraline-d3 and tamoxifen-13C2,15N were purchased from Sigma-Aldrich (Gillingham, UK), and bezafibrate-d6 was purchased from QMX Laboratories (Thaxted, UK). These were purchased as solutions at a concentration of 0.1 mg mL⁻¹ or 1.0 mg mL⁻¹ in methanol or other appropriate solvents. If no solutions were available, 1.0 mg powder was purchased, and the entire contents of the vial were dissolved in methanol. The MS/MS detection parameters for the internal standards can be found in ESM Table S3.

All glassware in this paper was silanised to prevent the analytes and internal standards from absorbing to the surface. This was done by coating the internal surfaces of the glassware with 5% dimethylchlorosilane (DMDCS) in toluene (Sigma-Aldrich), rinsing with toluene (Sigma-Aldrich) twice, then rinsing again with MeOH three times, and leaving to dry between each coating or rinse.

Methods for sample collection

Sampling was carried out in a river catchment in South West England. Sampling involved collection of samples from all major WwTWs and receiving environmental waters. Twenty-four-hour composites using a 3700 ISCO sampler (RS Hydro) were collected for both influent and effluent wastewaters in each case. Grab samples were utilised in surface water samples. SPM was collected from influent samples.


SPM per litre was calculated by filtering 30 mL through a pre-dried and pre-weighed GF/F glass microfibre filter. This was then re-dried at 105 °C for 8 h, left to cool and re-weighed to quantify SPM in grams per litre. Data are shown in ESM Table S4. Digested sludge was collected on three consecutive days, both directly after digestion and prior to disposal. Data are shown in ESM Table S5.

All samples were kept on ice during sampling (ice was placed within the composite sampler to maintain a cool temperature of 0–4 °C and promote stability) or placed in a cool box and kept on ice until the samples were transported to the lab. Once at the lab, samples for liquid analysis were transferred to 125-mL PPE bottles (Fisherbrand) and frozen (-20 °C) for further preparation and analysis at a later date. For the influent samples, the remainder of the sample was filtered to collect the SPM, which was then frozen (-20 °C). Most compounds do not adsorb to the PPE bottles, with very few exceptions (ESM Table S6).

Methods for extraction and analysis

The methods used for sample preparation of both liquid and solid matrices, as well as their analysis, can be found

Fig. 1 Flow chart from sample preparation to analysis, for analysis of liquid and solid samples by ESI- and ESI+ methods

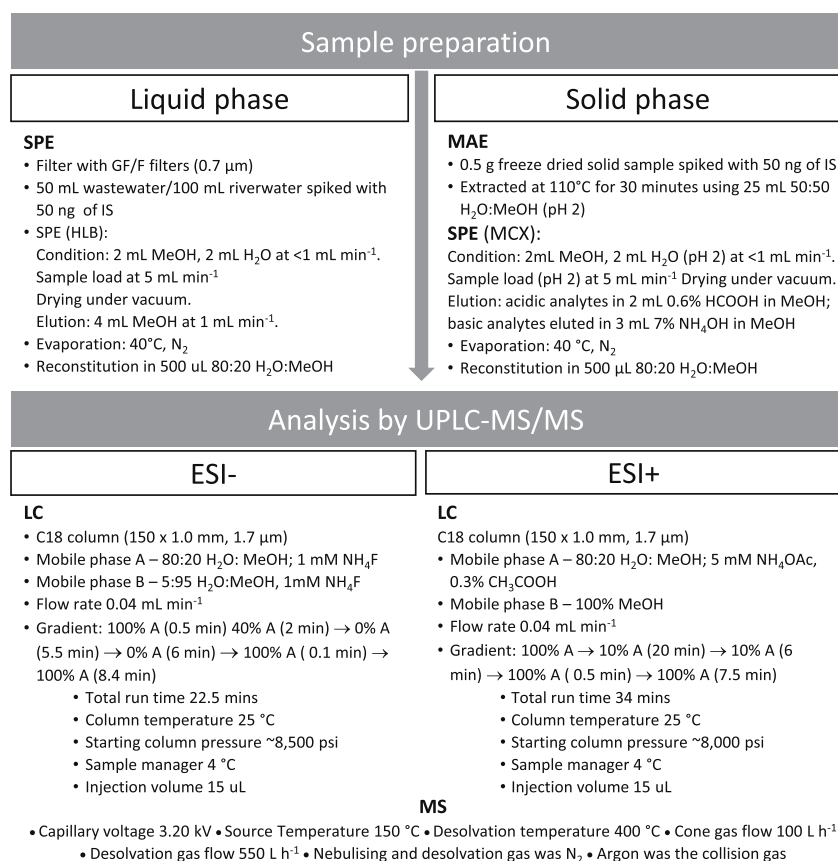

in Fig. 1; this is also discussed in more detail below. Development of an extraction method for liquid matrices with hydrophilic-lipophilic-balanced solid-phase extraction (HLB SPE) was developed based on the method published by Kasprzyk-Hordern et al. [40]. The microwave-assisted extraction (MAE) for solid matrices was developed based on the method published by Petrie et al. [39].

Extraction for liquid matrices

The samples were filtered with a GF/F glass fibre filter (0.7 μm) (Whatman, Sigma-Aldrich) and adjusted to pH 7.5–8.5. 50 mL of influent or effluent (100 mL surface water) was accurately measured and spiked with 50 μL of 1 $\mu\text{g mL}^{-1}$ (50 ng) internal standard mixture in MeOH. The 60-mg HLB SPE cartridges (OASIS, Waters, UK) were conditioned and equilibrated with 2 mL of MeOH, followed by 2 mL of deionised water (H₂O) at a rate of < 1 mL min⁻¹ (under gravity). The deionised H₂O was obtained from a Milli-Q system (18.2 M Ω). The sample was loaded at 5 mL min⁻¹ before the cartridges were dried under vacuum. Analytes were eluted using 4 mL MeOH at 1 mL min⁻¹ (under gravity). Once eluted, the extracts were evaporated to dryness at 40 °C, with a steady flow of nitrogen using a TurboVap LV concentration workstation. Finally, the samples were reconstituted in 500 μL of

80:20 H₂O:MeOH, mixed thoroughly to ensure completed dissolution and transferred to LC vials (polypropylene) (Waters, UK).

Extraction for solid matrices

The solid samples, digested sludge and SPM were initially frozen and freeze-dried (ScanVac, CoolSafe freeze dryer, Lyngø, Denmark). The freeze-dried samples were homogenised, and 0.25–0.5 g was weighed out and spiked with 50 μL of 1 $\mu\text{g mL}^{-1}$ (50 ng) internal standard mixture in MeOH. This was left for 30 min to 1 h, for the methanol to evaporate off. The samples were then transferred to MAE perfluoroalkoxy (PFA) tubes with 25–30 mL of 50:50 acidified H₂O (pH 2):MeOH. The MAE tubes were placed in 800 W MARS 6 microwave (CEM, UK). The temperature was ramped to 110 °C, over 10 min, then held at this temperature for 30 min, before allowing the samples to cool. SPE was then carried out using 60-mg mixed-mode cationic exchange (MCX) cartridges (Oasis, Water, UK). These were conditioned and equilibrated using 2 mL MeOH and 2 mL acidified H₂O (pH 2) at < 1 mL min⁻¹ under gravity. The samples were then loaded at 5 mL min⁻¹ and dried under vacuum. Once dried, the acid analytes were eluted first with 2 mL of 0.6% formic acid (HCOOH) (Sigma-Aldrich) in MeOH. The basic analytes were eluted second with 3 mL of 7% ammonium hydroxide (NH₄OH)

in MeOH (Sigma-Aldrich). These extracts were evaporated to dryness at 40 °C, with a steady flow of nitrogen using a TurboVap. The residue was reconstituted in 500 µL of 80:20 H₂O:MeOH, mixed thoroughly to ensure completed dissolution and transferred to LC vials.

Analysis of samples

The analytes were separated by UPLC performed on the Waters ACQUITY UPLC™ system (Waters, UK). The column used was a reversed-phase C18 column (Waters, UK), 150 mm × 1.0 mm, with a particle size of 1.7 µm. The samples were analysed with Xevo Triple Quadrupole (TQD) Mass Spectrometer (Waters, UK), equipped with an electrospray ionisation (ESI) source in positive and negative modes. To optimise ionisation, two different sets of parameters were used for the ESI-positive and ESI-negative modes. The parameters for these can be seen in Fig. 1. The systems were controlled using MassLynx (Waters, UK). Argon (99.998%) gas, supplied by a BOC cylinder, was used as a collision gas. The nebulising gas was nitrogen, provided by a high-purity nitrogen generator (Waters, Manchester, UK). Two mobile phases were used in the gradient mode ESI+: mobile phase A contained 80:20 H₂O:MeOH with 5 mM ammonium acetate (NH₄OAc) and 3 mM acetic acid (CH₃COOH) (pH 4.7), and mobile phase B contained 100% MeOH. Starting conditions were 100% A, decreasing to 10% A over 20 min, maintained at this level for 6 min, before increasing back to 100% A over 0.5 min and held for 7.5 min to return the column to equilibrium; in ESI-: mobile phase A contained 80:20 H₂O:MeOH with 1 mM ammonium fluoride (NH₄F) and mobile phase B, which was 5:95 H₂O:MeOH containing 1 mM NH₄F. The gradient began at 100% A for 0.5 min and reduced to 40% A over 2 min, before being further reduced to 0% A over the next 5.5 min. It was held at 0% A for 6 min before increasing back up to 100% A over 0.1 min. This was maintained for 8.4 min to re-equilibrate the column. The HPLC grade MeOH, NH₄OAc (Fluka) and CH₃COOH were obtained from Sigma-Aldrich (Gillingham, UK). The NH₄F (Fluka) was obtained from Fisher Scientific (Loughborough, UK) and deionised water (18.2 MΩ) obtained from a Milli-Q system. The temperature of the built-in sample manager was 4.0 °C with an injection volume of 20 µL. The mobile phases were run at a rate of 0.04 mL min⁻¹, as a gradient of a high ratio of aqueous (80%) to 100% MeOH in both positive and negative ionisation modes. The exact gradients and composition of the mobile phases can be also seen in Fig. 1. Data processing was carried out using TargetLynx software, which is an extension to MassLynx (version 4.1, Waters).

Instrument performance

To quantify the analytes, an internal standard approach with an 18-point calibration curve was used. For the majority of the

compounds, the analysed range covers 6 orders of magnitude, from ng L⁻¹ to mg L⁻¹. Each point was repeated every 24 h over 3 days. These calibration samples were prepared in a ratio of 80:20 H₂O:MeOH unbuffered solutions. The signal-to-noise ratios of these samples were used to determine the concentration of the instrument detection limit (IDL) and instrument quantification limit (IQL), where S/N ≥ 3 or S/N ≥ 10, respectively. Determination coefficients (*r*²) were calculated for the full linear range (IDL to ≤ 1000 µg L⁻¹). Inter- and intra-day precision and accuracy were calculated from repeated injections, at regular intervals (*n* = 3) of three concentrations (10 µg L⁻¹, 100 µg L⁻¹, 500 µg L⁻¹) in 80:20 H₂O:MeOH, across 24 h (intra-day) and across 72 h (inter-day).

Method performance

As SPE was used as a preconcentration and clean-up step, the recovery of each analyte must be assessed. Absolute and relative recoveries for SPE of liquid matrices were calculated from matrices spiked in duplicate (*n* = 2) at three concentrations for A-ESI+ (100 ng L⁻¹, 1000 ng L⁻¹ and 5000 ng L⁻¹ for effluent and influent, and 50 ng L⁻¹, 500 ng L⁻¹ and 2500 ng L⁻¹ for surface water). For B-ESI+ and C-ESI-, the matrices were spiked at two concentrations (100 ng L⁻¹ and 1000 ng L⁻¹ for effluent and influent, and 50 ng L⁻¹ and 500 ng L⁻¹ for surface water). For influent suspended particulate matter (SPM) and digested sludge (DS), the absolute and relative recoveries take into account MAE and SPE. The samples were spiked at 50 ng g⁻¹ and 100 ng g⁻¹.

Method detection limits (MDLs) and MQLs were calculated using Eq. 1

$$\text{MDL} = \frac{\text{IDL} \times 100}{\text{Rec} \times \text{Cf}} \quad (1)$$

where IDL is the instrumental detection limit, which is calculated as discussed in the section “Method performance”; 100 is the conversion factor for recovery of a specific matrix (Rec); and Cf is the concentration factor for the specific liquid matrix, e.g. 200 for surface waters or 100 for effluent of effluent. For solid matrices, Cf is replaced with a conversion factor of 2, which converts the volume into grams, based on the 0.25 g of solid matrix being extracted into a 0.5-mL vial for analysis. MQL is calculated with the same equation but by replacing IDL with IQL.

Furthermore, the accuracy and precision of the overall method, including SPE, are also required. These were calculated from samples of 3 matrices taken from 5 different WwTWs in the South West UK. These samples were spiked at 100 ng L⁻¹ and 1000 ng L⁻¹ for influent and effluent samples and at 50 ng L⁻¹ and 500 ng L⁻¹ for surface water for A-ESI+ and at 100 ng L⁻¹ for all matrices for B-ESI+ and C-ESI-. The accuracy of the method was determined from the

percentage deviation from the known concentration of analyte added to the sample. Precision was calculated as the relative standard deviation (RSD) of the replicates.

It has been found that complex matrices such as influent can affect the detection of analytes, especially when these samples have undergone SPE with HLB cartridges, as these cartridges extract a huge range of compounds. Therefore, matrix suppression was determined for the liquid matrices, which were extracted by employing a method using this approach. Samples for calculating matrix suppression were prepared by spiking samples with 50 μL of 1 $\mu\text{g mL}^{-1}$ of internal standards, after the elution step of SPE, and prior to evaporation and reconstitution. Once analysed, matrix suppression for each analyte was calculated using the following equation:

$$\text{Matrix suppression} = 1 - \left(\frac{\text{PA in SS} - \text{PA in US}}{\text{PA in MP QC}} \right) \quad (2)$$

where PA is the peak area of the analyte in spiked sample (SS), unspiked sample (US) and mobile phase quality control (MP QC) sample.

All matrices used were collected via grab sampling and homogenised, and all analyses were carried out on this single sample to ensure consistent results. Due to the use of these environmental matrices and the potential presence of analytes within the matrices prior to spiking, the 'blank' or unspiked (with analytes) portion of the sample was spiked with internal standards and analysed to confirm the concentration of analytes, prior to spiking, for recoveries and matrix effect analysis.

Quality control

Quality control samples were analysed before and after each batch at three concentrations (10 $\mu\text{g L}^{-1}$, 100 $\mu\text{g L}^{-1}$, 500 $\mu\text{g L}^{-1}$) along with procedural blanks, to ensure the method and instrumental performance and to monitor for contamination of the equipment.

All samples were spiked with the internal standards listed in ESM Table S3, for accurate quantification, and to account for loss from the point of spiking until analysis and to limit matrix effects.

Results and discussion

This paper provides an expanded and broader scope method based on a method published by Petrie et al. [39], enabling the analysis of household and agricultural chemicals whilst utilising one sample preparation protocol and comprehensive UPLC-MS/MS methodology.

UPLC-MS/MS method

All analytes were analysed using MRM and ESI⁻ and ESI⁺ modes and conditions previously selected by Petrie et al. [39]. Following EU guidelines, two MRM transitions were used for most of the 195 compounds (142 analytes and 53 ISTDs). This is except for cefalexin, ketoprofen, diclofenac, ibuprofen, 1,7-dimethylxanthine and norfluoxetine, which are to be considered as semi-quantitative. For ESI⁻, the parent ion $[\text{M}-\text{H}]^{-}$ was selected, and for ESI⁺ mode, $[\text{M}-\text{H}]^{+}$ was selected. The most abundant daughter ion was used for quantification and the second most for confirmation. ESM Table S2 includes the MRMs for all analytes, the acquisition method was used to analyse them and the assigned internal standard. The MRMs for the internal standards can be found in ESM Table S3.

The reversed-phase BEH C18 column provided good separation and sensitivity for all compounds. These conditions achieved good separation and peak shape for most analytes. Further information can be found in Fig. 1. A more acidic mobile phase was also trialled in the ESI⁺ method, and although the peak shape and separation for the acidic compounds, particularly the fluoroquinolones, were improved slightly, they only benefitted a small number of compounds and had a detrimental effect against many others. Satisfactory separation of analytes was achieved within a < 40-min retention time.

Quality control criteria as recommended by the EU directive [41] utilised quality control samples, standard tolerances of ion ratio, chromatographic retention time, relative retention time and signal-to-noise ratio.

Instrument performance

The instrument performance was assessed by considering linear response, inter- and intra-day precision and accuracy and, finally, instrumental detection and quantification limits. This data can be found in Table 1.

For the linear response, a linear range covering several orders of magnitude with the r^2 value of ≥ 0.997 was ideal. However, for a few analytes (triclosan, benzoylcegonine and mirtazapine), this could only be achieved through splitting the linear range into two overlapping ranges, each with the r^2 value ≥ 0.997 . This allowed adequate quantification across the entire range. Of the 142 compounds, 120 have r^2 values ≥ 0.997 . Twenty-one of the remaining compounds have r^2 values ≥ 0.992 . Although this is not ideal, it is still adequate for accurate quantification, as the other parameters indicate. Chlorpyrifos has the lowest r^2 value, likely due to the lack of an analogous internal standard; however, it passes further instrumental performance criteria.

The intra-day instrumental performance is high across many compounds. Out of the 142 analytes, the majority of which are very precise, with a deviation of $\leq 5\%$. The

Table 1 Instrumental performance data for ECs of interest in the mobile phase (ordered by class)

Class of analyte	Analyte	RT	RRT	Linearity	r^2	Intra-day instrumental performance		Inter-day instrumental performance		IDL _{SN} ^N ($\mu\text{g L}^{-1}$)	IQL _{SN} ^N ($\mu\text{g L}^{-1}$)
						Precision (deviation) (%)	Accuracy (%)	Precision (%)	Accuracy (%)		
UV filter	Benzophenone-1	9.6	0.9	0.06–684.0	0.996	2.3	106.8	3.3	106.7	0.01	0.06
	Benzophenone-2	7.7	1.0	0.05–583.8	0.997	1.1	99.6	4.2	97.6	0.01	0.05
	Benzophenone-3	21.2	1.2	0.05–404.0	0.995	3.2	84.9	4.5	86.8	0.01	0.05
	Benzophenone-4	6.9	0.9	1.01–502.5	0.997	2.3	103.0	3.8	105.1	0.31	1.01
Parabens	Methylparaben	7.5	1.0	0.06–1122	0.998	1.1	93.3	6.0	97.4	0.01	0.06
	Ethylparaben	8.3	1.0	0.11–663.6	0.997	2.6	112.3	2.1	113.1	0.03	0.11
	Propylparaben	9.2	1.0	0.12–462.0	0.997	5.7	96.4	4.3	98.4	0.04	0.12
Plasticizer	Butylparaben	10.1	1.0	0.06–696.6	0.997	5.0	97.1	3.6	100.3	0.01	0.06
	Bisphenol A	9.0	1.1	0.10–626.4	0.997	2.4	103.6	1.3	104.6	0.03	0.10
Steroid estrogens	E1	9.8	1.0	0.49–989.0	0.998	1.8	96.9	2.1	98.6	0.10	0.49
	E2	9.8	1.0	0.47–949.0	0.997	3.1	96.6	2.6	96.3	0.09	0.47
Antibiotics and antibacterials	EE2	9.7	1.0	0.48–950.0	0.997	2.6	94.6	3.3	93.2	0.10	0.48
	Sulfasalazine	7.1	0.8	0.90–904.0	0.999	3.9	105.2	2.4	104.7	0.27	0.90
	Clarithromycin	18.9	1.1	0.06–561.0	0.999	2.6	99.8	2.4	101.8	0.01	0.06
	Azithromycin	14.0	0.9	0.001–1000	0.998	4.5	108.9	1.5	102.0	0.01	0.05
	Trimethoprim	8.4	1.0	0.10–500.0	0.998	3.0	96.9	2.2	99.5	0.03	0.10
	Sulfamethoxazole	9.6	1.0	0.10–1000	0.999	3.5	95.1	2.4	96.0	0.03	0.10
	Triclosan ^a	12.3	1.2	1.13–225.6 112.8–1128	0.997/0.998	9.4	69.1	6.5	71.4	0.34	1.13
	Amoxicillin	3.1	0.2	0.06–439.5	0.995	5.3	105.7	6.7	94.4	0.02	0.06
	Metronidazole	5.3	1.0	1.00–1000	0.999	2.5	105.0	1.2	102.9	0.06	0.21
	Sulfadiazine	4.8	0.9	0.05–795.2	0.999	2.8	105.3	1.5	104.4	0.01	0.03
Antibiotics and antibacterials	Cefalexin ^b	9.2	0.3	15.9–200	0.995	9.5	111.3	12.3	102.9	4.78	15.94
	Ofloxacin	9.6	1.0	0.23–986.0	0.998	4.2	97.4	2.8	95.9	0.07	0.23
Antibiotics and antibacterials	Ciprofloxacin	9.9	1.0	1.18–902	0.999	8.7	89.0	5.5	90.2	0.35	1.18
	Tetracycline	10.0	1.0	0.06–864.0	0.999	6.8	115.1	8.5	113.1	0.02	0.06
Antibiotics and antibacterials	Danofloxacin	10.2	1.0	1.05–1000	0.998	7.3	106.0	6.0	99.2	0.32	1.05
	Oxytetracycline	10.4	1.1	2.36–800.8	0.997	4.6	93.5	3.0	88.9	0.71	2.36
Antibiotics and antibacterials	Chloramphenicol	12.6	0.6	1.74–400	0.999	3.8	103.5	3.0	100.8	0.52	1.74
	Penicillin G	13.1	0.5	4.68–93.6	0.994	10.3	115.5	4.4	111.7	0.02	0.07
Antibiotics and antibacterials	Penicillin V	14.5	0.8	5.00–200	0.993	4.4	88.5	15.0	96.8	0.15	0.49
	Erythromycin	17.2	1.0	204.4–1022	0.999	2.3	94.4	2.9	95.2	0.20	0.65
Antibiotics and antibacterials	Prulifloxacin	18.0	1.9	100–1000	0.997	4.4	98.7	8.9	86.4	2.44	8.13

Table 1 (continued)

Class of analyte	Analyte	RT	RRT	Linearity Range ($\mu\text{g L}^{-1}$)	r^2	Intra-day instrumental performance		Inter-day instrumental performance		IDL _{IS/N} ($\mu\text{g L}^{-1}$)	IQL _{IS/N} ($\mu\text{g L}^{-1}$)
						Precision (deviation) (%)	Accuracy (%)	Precision (%)	Accuracy (%)		
Antifungal	Norfloxacin	9.7	1.0	0.01–1000	0.996	4.1	85.5	4.4	85.1	0.002	0.01
	Griseofulvin	17.2	0.9	0.26–205.2	0.999	1.6	89.2	3.0	91.6	0.08	0.26
Hypertension	Ketoconazole	21.7	1.2	0.02–800.0	0.999	3.8	94.8	2.5	91.7	0.01	0.02
	Valsartan	7.6	0.9	1.12–1122	0.998	1.9	115.8	3.5	118.6	0.34	1.12
NSAIDs	Irbesartan	8.6	1.0	0.50–603.6	0.998	2.6	96.9	4.1	98.3	0.10	0.50
	Lisinopril	7.1	0.9	0.93–372.5	0.995	2.2	97.2	7.2	95.2	0.09	0.93
	Ketoprofen ^b	7.9	0.9	0.54–1085	0.998	2.2	99.9	2.6	99.4	0.11	0.54
	Ibuprofen ^b	9.8	1.0	0.05–1071	0.998	2.4	93.7	2.3	94.2	0.01	0.05
	Naproxen	8.1	1.0	0.49–989.0	0.998	1.5	97.7	2.5	98.3	0.10	0.49
Lipid regulator	Diclofenac ^b	9.0	1.0	0.10–619.2	0.997	7.9	89.6	4.5	91.8	0.03	0.10
	Acetaminophen	5.1	1.0	0.54–1070	0.998	1.6	97.4	2.6	99.0	0.11	0.54
Antihyperlipidaemic	Bezafibrate	7.9	1.0	0.10–976.0	0.998	2.3	97.8	2.8	97.9	0.03	0.10
	Atorvastatin	9.3	1.1	0.05–500.0	0.997	2.6	98.0	3.5	100.9	0.01	0.05
Antihyperintensive	Gemfibrozil	23.3	1.2	1.01–100.5	0.994	7.8	118.5	6.9	121.1	0.11	0.35
	Candesartan cilexetil	23.0	0.9	226.8–680.4	0.995	5.2	100.5	0.9	106.9	1.58	5.28
Antihistamine	Fexofenadine	8.4	1.0	0.09–937.5	0.998	2.1	106.3	6.5	104.6	0.03	0.09
	Cetirizine	18.7	1.0	0.08–417.7	0.999	1.3	100.5	1.3	100.8	0.02	0.08
Diabetes	Sildenafil	18.3	1.0	0.01–1000	1.000	3.5	99.5	3.0	99.1	0.002	0.01
	Metformin	2.8	1.0	0.43–862.5	0.998	1.5	96.3	1.3	97.0	0.09	0.43
Cough suppressant	Gliclazide	17.8	1.0	0.05–508.0	0.997	2.1	93.2	2.8	95.3	0.01	0.05
	Sitagliptin	11.8	0.7	0.08–646.4	0.998	3.2	111.7	3.0	110.3	0.01	0.02
Beta blocker	Pholcodine	3.7	0.9	1.14–570.0	0.999	4.7	99.5	3.3	99.2	0.35	1.14
	Atenolol	4.3	1.0	0.10–502.5	0.999	2.1	95.3	2.3	96.8	0.03	0.10
H ₂ receptor agonist	Metoprolol	11.2	1.0	0.05–507.5	0.999	1.3	96.8	2.0	96.1	0.01	0.05
	Propranolol	15.1	1.0	0.09–434.9	0.999	2.0	105.4	1.0	106.2	0.03	0.09
X-ray contrast media	Bisoprolol	13.7	0.8	0.10–1004	0.999	4.8	100.4	2.0	96.0	0.0004	0.0012
	Ranitidine	4.6	1.1	5.17–517.0	0.998	2.5	100.1	9.7	97.4	1.03	5.17
Various	Cimetidine	5.3	1.0	0.52–1043	0.999	4.2	104.1	9.0	99.3	0.10	0.52
	Iopromide	4.9	0.9	5.79–1158	0.997	5.0	101.2	12.0	105.4	1.16	5.79
Drug precursor	Buprenorphine	21.8	1.2	0.08–100	0.996	8.9	94.5	11.5	88.2	0.02	0.08
	Ephedrine/pseudoephedrine	7.2	1.0	0.10–500.0	0.997	4.1	94.0	3.4	97.3	0.03	0.10

Table 1 (continued)

Class of analyte	Analyte	RT	RRT	Linearity Range ($\mu\text{g L}^{-1}$)	r^2	Intra-day instrumental performance		Inter-day instrumental performance		IDL _{LS/N} ($\mu\text{g L}^{-1}$)	IQL _{LS/N} ($\mu\text{g L}^{-1}$)	
						Precision (deviation) (%)	Accuracy (%)	Precision (%)	Accuracy (%)			
Anticancer	Norephedrine	6.3	0.9	0.50–1000	0.999	4.3	96.3	5.1	95.2	0.01	0.50	
	Azathioprine	7.8	0.9	0.10–490.0	0.999	7.6	97.5	13.9	97.4	0.03	0.10	
	Methotrexate	7.9	1.0	0.92–458.0	0.997	8.7	108.0	4.1	112.2	0.28	0.92	
	Ifosfamide	12.7	1.1	0.05–509.0	0.999	2.4	93.6	2.7	95.3	0.01	0.05	
	Tamoxifen	22.4	1.0	0.03–668.4	0.998	4.0	96.0	2.4	96.8	0.01	0.03	
	Imatinib	15.4	0.8	0.88–88.4	0.994	2.5	103.8	1.5	101.3	0.08	0.28	
	Capecitabine	16.1	0.9	0.01–594.6	0.999	2.3	89.2	2.8	89.7	0.001	0.004	
	Bicalutamide	18.2	0.9	0.10–784.0	0.995	2.7	90.1	2.9	92.0	0.03	0.10	
	Ketamine	10.6	1.0	0.05–500.0	0.998	1.8	92.5	1.3	93.6	0.01	0.05	
	Norketamine	11.1	1.0	0.10–500.0	0.999	1.8	94.1	3.2	94.0	0.03	0.10	
	Venlafaxine	14.1	1.3	0.04–434.8	0.998	2.5	91.2	1.7	90.5	0.01	0.04	
	Desmethylvenlafaxine	10.8	1.0	0.10–500.0	0.998	2.8	101.3	2.1	102.3	0.03	0.10	
	Fluoxetine	18.4	1.0	0.05–1000	0.999	1.7	96.8	1.8	98.3	0.01	0.05	
	Norfluoxetine ^b	18.4	1.0	0.05–500.0	0.998	1.5	102.7	3.1	103.1	0.01	0.05	
Anaesthetic and metabolite	Sertraline	19.2	1.0	0.05–500.0	1.000	1.6	95.3	1.7	95.7	0.01	0.05	
	Mirtazapine ^a	13.5	1.0	0.05–100.0, 50.0–500.0	0.999/0.997	3.4	94.8	2.7	97.6	0.01	0.05	
	Citalopram	15.1	1.0	0.50–1000	0.999	0.7	101.2	2.6	101.8	0.05	0.50	
	Desmethylcitalopram	15.2	1.0	0.05–500.0	0.998	1.8	103.0	3.0	103.4	0.01	0.05	
	Paroxetine	17.3	0.9	5.00–600	0.998	3.2	103.4	1.3	102.1	0.01	0.03	
	Duloxetine	17.8	1.0	1.00–1000	0.997	3.0	91.2	13.6	78.3	0.003	0.01	
	Amitriptyline	18.2	1.0	0.11–885.0	1.000	4.5	99.6	2.4	96.8	0.03	0.11	
	Nortriptyline	18.4	1.0	0.22–800	0.999	4.0	95.5	3.1	92.9	0.07	0.22	
	Norsertaline	19.8	1.0	0.23–100	0.999	8.7	99.0	11.0	91.8	0.07	0.23	
	Carbamazepine	16.2	1.0	0.05–514.0	1.000	2.0	91.7	1.6	92.7	0.01	0.05	
	Carbamazepine-10,11-epoxide	13.5	0.8	0.10–1000	0.997	1.6	88.9	2.1	89.9	0.03	0.10	
	10,11-Dihydro-10-hydroxycarbamazepine	13.5	0.8	0.50–100.0	0.997	2.8	92.2	5.6	93.8	0.05	0.50	
	Calcium channel blocker	Diltiazem	16.7	1.0	0.10–486.2	0.996	2.3	92.7	2.3	93.6	0.01	0.10
	Hypnotic	Verapamil	16.2	1.0	0.01–600	0.998	2.9	103.1	2.4	101.9	0.001	0.004
Temazepam		18.2	1.0	0.05–500.0	0.998	1.0	97.0	1.6	97.9	0.01	0.05	
Oxazepam		17.8	1.0	0.10–800	0.999	3.3	94.8	3.4	94.3	0.02	0.08	
Diazepam		19.5	1.0	0.01–1000	1.000	1.6	100.7	4.5	99.6	0.003	0.01	

Table 1 (continued)

Class of analyte	Analyte	RT	RRT	Linearity Range ($\mu\text{g L}^{-1}$)	r^2	Intra-day instrumental performance		Inter-day instrumental performance		IDL _{LS/N} ($\mu\text{g L}^{-1}$)	IQL _{LS/N} ($\mu\text{g L}^{-1}$)
						Precision (deviation) (%)	Accuracy (%)	Precision (%)	Accuracy (%)		
Antipsychotic	Quetiapine	17.9	1.0	0.05–1000	0.997	1.4	95.3	1.2	96.4	0.01	0.05
	Risperidone	13.7	0.8	0.01–200	0.997	3.2	101.6	1.2	96.8	0.002	0.01
Dementia	Donepezil	13.9	0.9	0.01–1000	0.998	2.6	110.8	1.3	107.7	0.17	0.58
	Memantine	15.7	1.0	0.05–506.4	0.998	3.5	106.3	0.9	104.3	0.02	0.05
Human indicators	Creatinine	2.7	1.0	1.00–1000	0.999	1.4	100.5	2.8	100.1	0.30	1.00
	Nicotine	3.3	0.8	1.00–500.0	0.998	1.2	98.3	2.4	98.4	0.30	1.00
	Caffeine	8.3	1.2	0.50–500.0	0.999	1.7	99.6	2.8	100.4	0.10	0.50
	Cotinine	7.2	1.0	0.05–1000	0.999	1.5	98.4	1.5	98.8	0.01	0.05
Analgesics and metabolites	1,7-Dimethylxanthine ^b	6.8	0.9	1.00–500.0	0.999	6.0	94.3	9.9	94.9	0.30	1.00
	Morphine	3.5	1.0	1.00–500.0	0.998	2.9	99.1	2.5	97.5	0.30	1.00
	Dihydromorphine	3.3	1.0	0.05–500.0	0.997	4.4	106.0	2.7	108.5	0.01	0.05
	Normorphine	3.4	1.0	1.00–500.0	0.999	1.5	100.9	2.2	99.8	0.30	1.00
	Methadone	17.6	1.0	0.05–400.0	0.998	1.5	98.7	1.4	100.2	0.01	0.05
	EDDP	14.8	1.0	0.05–500.0	0.999	1.2	96.5	1.1	96.4	0.01	0.05
	Codeine	6.1	1.0	0.50–500.0	0.997	2.0	93.5	4.0	95.1	0.10	0.50
	Norcodeine	6.5	1.1	1.00–500.0	0.998	2.8	98.5	4.8	98.6	0.30	1.00
	Dihydrocodeine	5.5	0.9	0.10–500.0	0.999	1.6	94.2	2.1	94.6	0.03	0.10
	Tramadol	11.0	1.0	1.00–500.0	0.999	1.6	100.1	1.9	98.4	0.01	1.00
	<i>N</i> -Desmethyltramadol	11.9	1.1	0.50–500.0	0.998	2.5	92.5	2.2	94.4	0.01	0.50
	<i>O</i> -Desmethyltramadol	8.3	1.2	1.00–400.0	0.997	3.3	95.3	4.9	98.5	0.01	1.00
Stimulants and metabolites	Amphetamine	8.4	1.0	0.10–500.0	0.999	4.4	100.8	1.6	100.7	0.03	0.10
	Methamphetamine	8.5	1.0	0.10–500.0	0.999	2.2	101.0	1.3	101.1	0.03	0.10
	MDMA	8.6	1.0	0.05–1000	0.999	1.3	99.2	1.7	99.8	0.01	0.05
	MDA	8.6	1.0	0.10–1000	0.998	1.1	98.4	0.7	100.0	0.03	0.10
	Cocaine	11.3	1.0	0.05–500.0	0.999	2.2	97.2	1.5	99.0	0.01	0.05
	Benzoyllecgonine ^a	9.7	1.0	0.05–100.0, 50.0–500.0	0.998/0.999	2.4	103.4	0.9	103.2	0.01	0.05
	Anhydroecgonine methyl ester	3.5	1.3	0.50–500.0	0.999	2.3	101.1	2.4	98.7	0.10	0.50
	Cocaine/lecgonine	12.9	1.0	0.05–500.0	0.999	2.8	95.1	1.7	94.7	0.01	0.05
	Mephedrone	9.8	1.0	0.05–500.0	0.998	1.8	87.1	2.9	85.7	0.01	0.05
	MDPV	12.1	0.9	0.05–500.0	0.999	2.2	99.6	0.7	101.4	0.01	0.05
	Heroin	10.9	1.0	0.50–500.0	0.999	1.9	98.2	1.8	99.3	0.10	0.50

Table 1 (continued)

Class of analyte	Analyte	RT	RRT	Linearity Range ($\mu\text{g L}^{-1}$)	r^2	Intra-day instrumental performance		Inter-day instrumental performance		IDL _{S/N} ($\mu\text{g L}^{-1}$)	IQL _{S/N} ($\mu\text{g L}^{-1}$)
						Precision (deviation) (%)	Accuracy (%)	Precision (%)	Accuracy (%)		
Opioid and metabolite Pesticides, fungicides and herbicides	6-Acetylmorphine	7.7	1.1	0.10–500.0	0.997	6.1	95.3	5.1	100.1	0.03	0.10
	Thiamectoxam	8.3	0.4	1.00–100	0.994	4.7	93.8	5.4	96.9	0.02	0.06
	Imidacloprid	10.1	0.6	0.10–595.2	0.996	2.8	100.5	5.5	103.5	0.01	0.04
	Clothianidin	10.4	0.5	1.00–800	0.999	3.2	97.9	3.3	98.6	0.01	0.04
	Metazachlor	17.1	1.0	0.05–1011	0.999	2.5	106.0	2.6	104.7	0.004	0.01
	Terbutylazine	19.3	1.0	0.05–519	1.000	2.4	99.8	3.3	97.5	0.01	0.02
	Methiocarb	19.4	1.0	0.08–1007	0.999	1.9	101.8	1.8	100.6	0.02	0.08
	Dichlofluanid	20.4	1.1	6.83–1092	0.994	3.8	94.9	4.4	90.9	1.29	4.30
	Flufenacet	20.5	1.2	0.01–986.0	0.997	2.0	104.2	2.9	106.2	0.002	0.01
	Oxadiazon	24.2	1.2	1.00–99.6	0.996	4.0	95.5	2.8	97.1	0.02	0.08
	Chlorpyrifos ^c	24.8	1.5	1.87–98.5	0.985	11.8	80.7	7.8	83.3	0.56	1.87
	Triallate	24.9	1.3	0.03–79.0	0.992	7.6	81.3	13.2	70.6	0.01	0.03
Veterinary pharmaceuticals	Tylosin	17.3	1.0	0.56–560.0	0.999	2.2	99.5	4.0	100.2	0.11	0.56
	Sulfapyridine	6.4	1.2	0.05–800	0.999	2.6	110.7	1.1	109.5	0.01	0.03
	Sarafloxacin	10.9	0.7	0.88–442	0.995	5.2	112.1	2.3	107.1	0.22	0.75
	Ceftiofur	12.1	1.3	0.28–800.0	0.993	3.6	89.5	2.0	86.4	0.08	0.28
	Diazinon	21.9	1.2	0.11–2100	0.998	2.7	98.9	4.1	96.0	0.01	0.02

Where possible, instrumental performance was determined at concentrations of 10 $\mu\text{g L}^{-1}$, 100 $\mu\text{g L}^{-1}$ and 500 $\mu\text{g L}^{-1}$; i.e. those analytes where these concentrations were outside the range of linearity or results were <LOQ were not included

IDL instrumental detection limit, IQL instrumental quantification limit

^a Linearity range was split into two overlapping ranges to ensure the r^2 value ≥ 0.997

^b Semi-quantitative, due to only one MRM transition

^c Semi-quantitative, due to poor r^2 value

precision of 20 compounds was between 5 and 10%, and only one compound was > 10% (chlorpyrifos (11.8%), likely due to poor r^2). However, all deviations are inside the recommended maximum of 20%. For 119/142 compounds, the accuracy lay within the ideal range of 90–110%. Of the remaining compounds, only triclosan has an inaccuracy of > 20%.

Regarding the inter-day instrumental performance, the precision is high across many analytes with 115 analytes with a deviation of $\leq 5\%$. Of the 27 analytes with > 5% deviation, only 8 were > 10%. The inter-day accuracy of these compounds was also high, with only 20 compounds that deviated from the QC by > 10%, 4 of which were only slightly greater than 20%. Data for both intra- and inter-day precision and accuracy can be found ESM Table S7.

IDLs ranged from 0.4 ng L^{-1} (bisoprolol) to 4783 ng L^{-1} (cefalexin), and IQLs ranged from 1.2 ng L^{-1} (bisoprolol) to $15,940 \text{ ng L}^{-1}$ (cefalexin). Whilst many of these IQL concentrations are very low, the samples may still need to be concentrated (utilising SPE), as the concentrations of most compounds in environmental matrices are likely to be even lower.

Method performance

Sample extraction was carried out using SPE following the method shown in Fig. 1 (sections “Extraction for liquid matrices” and “Extraction for solid matrices”) and has shown good extraction performance for many CECs. The Oasis HLB sorbent is essential to multi-residue methods, as it has the ability to retain a large range of polar analytes at neutral pH and is therefore widely used in analysis of environmental matrices. However, two drawbacks have been found in the use of this sorbent. Firstly, HLB might result in low recovery of very polar compounds such as metformin and creatinine. This is easily remedied, as metformin and creatinine along with acetaminophen, nicotine, caffeine and 1,7-dimethylxanthine are present in the environment at such high levels that direct injection are utilised instead. The second drawback of HLB is its lack of selectivity. Whilst enabling the extraction of a large variety of polar analytes, in complex matrices, much of the matrix can be extracted along with the chosen analytes, causing significant signal interference. It is therefore important to assess matrix suppression. Previously, it was found that the use of HLB with digested solids provided poor results. This not only necessitated the use of an alternative sorbent, MCX, but also splitting the eluents into acidic and basic fractions (see Fig. 1) [39].

Figure 2 shows relative recoveries for all matrices: surface water, effluent, influent, SPM and digested solids. These very different matrices exhibit similar ranges of relative recovery. Most of the compounds had recoveries between 80 and 120%. Due to the complexity of environmental matrices, not all compounds were recovered adequately across all matrices; therefore, the number of compounds accurately quantifiable in each matrix varies from 138 analytes in surface water to 96 analytes in

digested solids. The data for relative and absolute recoveries for all matrices can be found in ESM Tables S8 and S9 and Fig. S2.

Matrix suppression was analysed for liquid matrices (see Fig. 3 and ESM Table S10). Proximity to zero shows limited matrix effects. Most analytes experienced matrix suppression, shown in Fig. 3 as a positive percentage. For all matrices, at least 108 analytes were below 70% suppression. However, a few compounds, primarily from ESI⁻, experienced significant (> 20%) signal enhancement. Some compounds experienced signal enhancement in some matrices and suppression in others (naproxen, erythromycin, EE2, bicalutamide, candesartan cilexetil, gemfibrozil, chlorpyrifos). This considerable variation is to be expected for a multi-residue method with this variety of analytes. It especially highlights the importance of using relevant analogous internal standards, as well as thoroughly exploring the different matrices relevant to this work.

The method accuracy of the method was 107.0%, and the precision was 13.4%. MQLs for liquid matrices range from 0.004 ng L^{-1} for bisoprolol, in surface water, to 3118 ng L^{-1} for creatinine, in influent (Table 2). For the remaining analytes, there is not enough environmental data on their presence in the UK to comment on whether the MDLs and MQLs are low enough to measure these levels in the catchment of interest. Therefore, this method was applied to five different matrices at a WWTW in the South West UK, which were influent, effluent, surface water, SPM and digested solids, to underline the utilisation of this method to relevant environmental matrices.

Overall, this method provides a clear benefit when used in a catchment-based study, compared to more complex approaches with multiple sample preparation and analytical methods for different classes of CECs from the same matrices. The main advantage of the multi-residue method presented in this paper is the applicability of this method to liquid and solid matrices. In particular, SPM, which is an often forgotten but critical aspect of wastewater, as it contains many more hydrophobic compounds at significant concentrations within this matrix [28, 42, 43].

Application to environmental matrices

The data is presented as average concentrations with variation across a week (Table 3 and ESM Tables S11–S15). As expected, the liquid influent fraction shows the highest variation in concentrations both across the week and between different CECs. Of the 138 analytes of this method, 70% were quantified in influent at this site. Creatinine had the highest average concentration of any analyte in this study ($1.3 \pm 0.3 \text{ mg L}^{-1}$) and is often used a human indicator due to its correlation to population [44, 45]. Other human indicators were present at high levels throughout the week. The other CECs were present at a range of concentrations. Methiocarb, for example, is present intermittently across the week with an average concentration of $3.7 \pm 0.6 \text{ ng L}^{-1}$, which is very close to its MQL of 0.86 ng L^{-1} . However, acetaminophen reaches concentrations several magnitudes higher

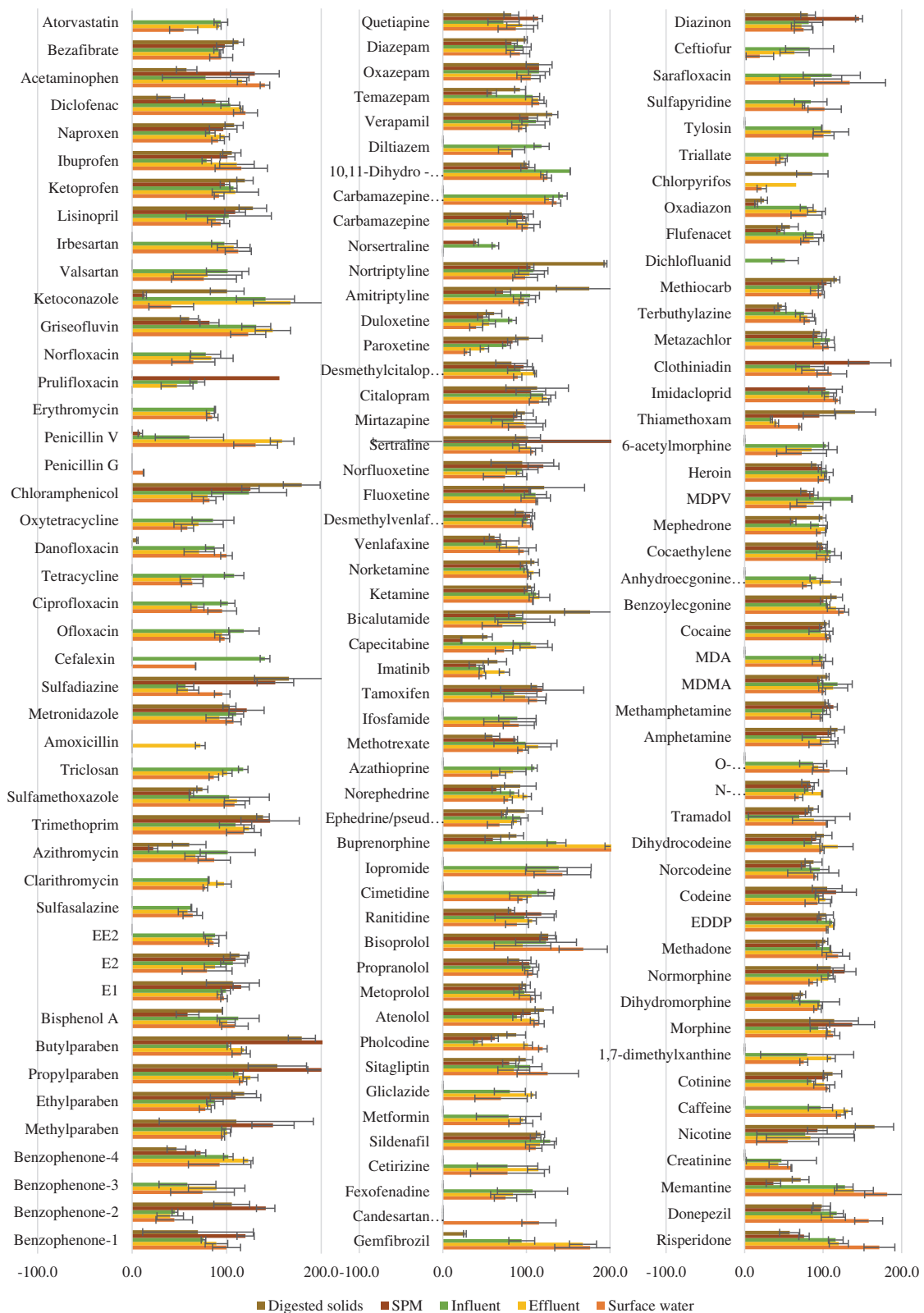


Fig. 2 Relative recoveries for all matrices. Error bars show the range of standard deviation

than methiocarb, at $393.6 \pm 100.0 \mu\text{g L}^{-1}$. This is slightly higher than those published in a review by Verlicchi et al. [11] which

had an absolute highest concentration of $246 \mu\text{g L}^{-1}$. It is interesting to see the presence of the (pesticide) methiocarb in

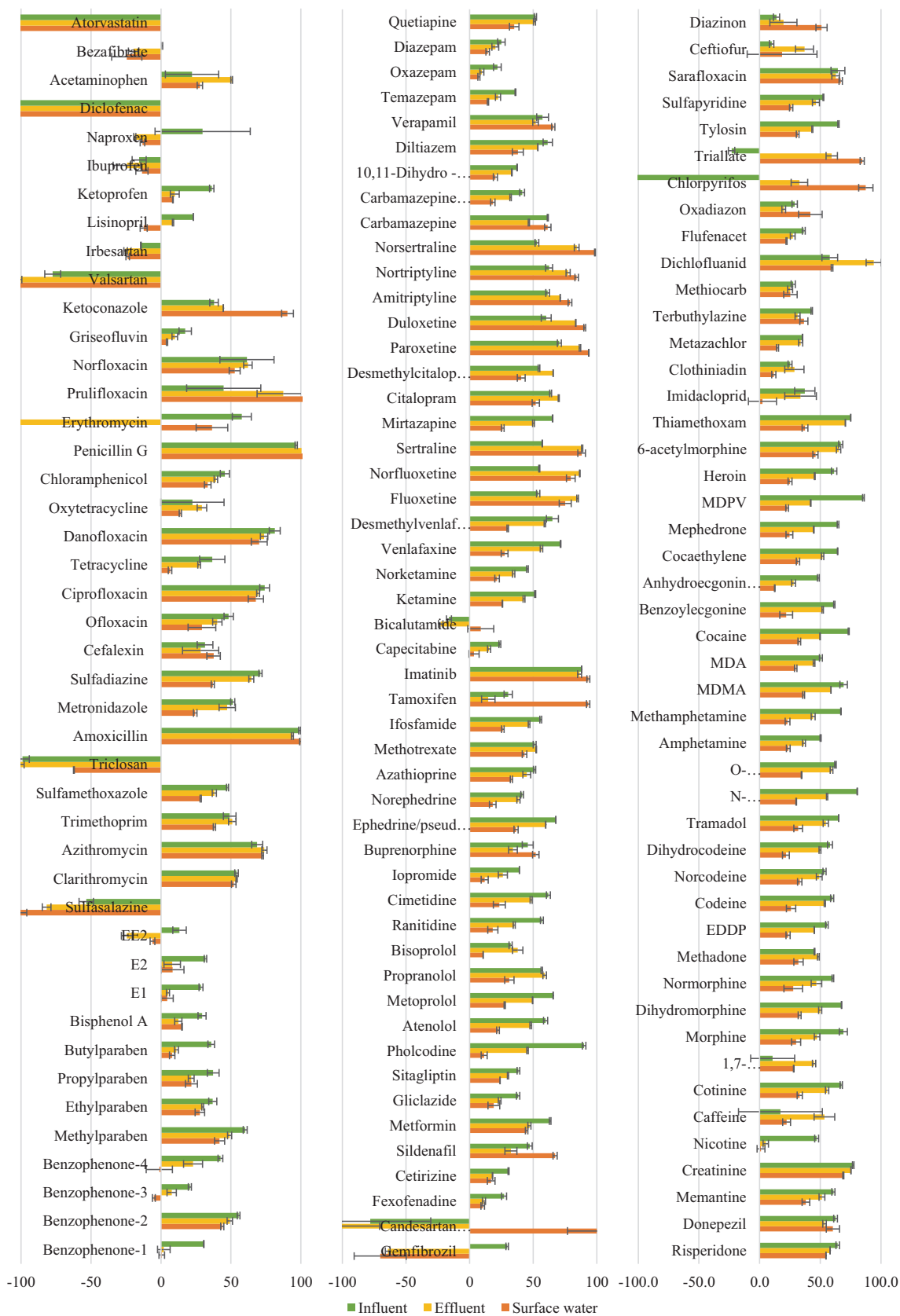


Fig. 3 Matrix suppression for all liquid matrices. Error bars show the range of standard deviation

Table 2 Method performance data for ECs of interest in the mobile phase (ordered by class)

Class of analyte	Analyte	Surface water (ng L ⁻¹)		Effluent (ng L ⁻¹)		Influent (ng L ⁻¹)		Solid particulate matter (ng g ⁻¹)		Digested solids (ng g ⁻¹)	
		MDL	MQL	MDL	MQL	MDL	MQL	MDL	MQL	MDL	MQL
UV filter	Benzophenone-1	0.07	0.35	0.14	0.71	0.23	1.15	0.004	0.02	0.14	0.70
	Benzophenone-2	0.16	0.79	0.34	1.68	0.36	1.82	0.004	0.02	0.09	0.44
	Benzophenone-3	0.15	0.77	0.19	0.97	0.37	1.87	–	–	–	–
	Benzophenone-4	2.09	6.90	5.78	19.1	7.83	25.8	0.21	0.70	4.01	13.2
Parabens	Methylparaben	0.08	0.40	0.19	0.94	0.28	1.41	0.003	0.02	0.06	0.31
	Ethylparaben	0.24	0.79	0.46	1.52	0.49	1.61	0.01	0.05	0.17	0.57
	Propylparaben	0.25	0.83	0.47	1.54	0.63	2.08	0.01	0.03	0.22	0.72
	Butylparaben	0.08	0.38	0.14	0.71	0.24	1.21	0.002	0.01	0.10	0.52
Plasticizer	Bisphenol A	0.26	0.86	0.56	1.84	0.85	2.79	0.03	0.09	0.27	0.88
Steroid estrogens	E1	0.78	3.92	0.15	7.69	1.96	9.78	0.04	0.21	1.68	8.38
	E2	0.90	4.48	1.41	7.03	1.84	9.22	0.04	0.21	1.48	7.41
	EE2	0.98	4.91	1.46	7.32	1.83	9.15	–	–	–	–
Antibiotics and antibacterials	Sulfasalazine	4.31	14.2	9.66	31.9	12.6	41.4	–	–	–	–
	Clarithromycin	0.18	0.90	0.28	1.40	0.34	1.69	–	–	–	–
	Azithromycin	0.08	0.26	0.21	0.68	0.14	0.45	0.03	0.10	0.01	0.04
	Trimethoprim	0.26	0.85	0.51	1.67	0.73	2.41	0.01	0.03	0.07	0.22
	Sulfamethoxazole	0.19	0.63	0.47	1.56	0.72	2.38	0.02	0.08	0.12	0.41
	Triclosan	2.93	9.68	4.55	15.0	4.93	16.3	–	–	–	–
	Amoxicillin	–	–	0.26	0.86	–	–	–	–	–	–
	Metronidazole	0.29	0.98	0.68	2.27	0.57	1.90	0.03	0.09	0.03	0.10
	Sulfadiazine	0.05	0.18	0.18	0.59	0.18	0.62	0.003	0.01	0.003	0.01
	Cefalexin	35.6	118.7	10.2	33.9	18.9	63.1	–	–	–	–
	Ofloxacin	0.35	1.17	0.72	2.40	0.58	1.93	–	–	–	–
	Ciprofloxacin	1.85	6.17	5.10	17.0	3.48	11.6	–	–	–	–
	Tetracycline	0.15	0.50	0.30	1.01	0.18	0.59	–	–	–	–
	Danofloxacin	1.58	5.28	4.45	14.85	3.62	12.08	–	–	2.84	9.45
	Oxytetracycline	6.04	20.1	10.1	33.6	8.26	27.5	–	–	–	–
	Chloramphenicol	3.18	10.6	6.52	21.7	4.21	14.0	0.21	0.69	0.15	0.48
	Penicillin G	0.89	2.98	–	–	–	–	–	–	–	–
	Penicillin V	0.56	1.86	0.92	3.06	2.40	8.00	0.84	2.80	–	–
	Erythromycin	1.15	3.83	2.35	7.85	2.22	7.41	–	–	–	–
Prulifloxacin	–	–	51.3	171.0	35.3	117.6	–	–	–	–	
Norfloxacin	0.01	0.04	0.02	0.06	0.02	0.07	–	–	–	–	
Antifungal	Griseofulvin	0.32	1.06	0.52	1.74	0.59	1.98	0.05	0.16	0.06	0.21
	Ketoconazole	0.06	0.21	0.03	0.10	0.04	0.12	0.02	0.07	0.00	0.01
Hypertension	Valsartan	2.81	9.26	6.40	21.1	7.24	23.9	–	–	–	–
	Irbesartan	0.89	4.47	1.88	9.38	2.50	12.5	–	–	–	–
	Lisinopril	2.17	21.7	4.25	42.5	3.25	32.5	0.04	0.43	0.25	2.47
NSAIDs	Ketoprofen	0.74	3.72	1.60	8.00	2.38	11.9	0.06	0.28	0.47	2.35
	Ibuprofen	0.06	0.31	0.08	0.42	0.19	0.93	0.005	0.02	0.07	0.36
	Naproxen	0.61	3.07	1.17	5.85	6.29	31.5	0.05	0.25	0.60	3.02
	Diclofenac	0.22	0.73	0.44	1.44	0.67	2.22	0.02	0.06	0.75	2.46
	Acetaminophen	1.20	6.02	2.39	12.0	138.0*	1017*	0.04	0.21	2.74	13.7
Lipid regulator	Bezafibrate	0.22	0.66	0.38	1.25	0.64	2.11	0.02	0.05	0.18	0.60
	Atorvastatin	0.14	0.70	0.17	0.84	0.17	0.85	–	–	–	–
Antihyperlipidemic	Gemfibrozil	0.30	1.00	0.63	2.11	1.12	3.75	–	–	0.20	0.67

Table 2 (continued)

Class of analyte	Analyte	Surface water (ng L ⁻¹)		Effluent (ng L ⁻¹)		Influent (ng L ⁻¹)		Solid particulate matter (ng g ⁻¹)		Digested solids (ng g ⁻¹)	
		MDL	MQL	MDL	MQL	MDL	MQL	MDL	MQL	MDL	MQL
Antihyperintensive	Candesartan	6.89	23.0	–	–	–	–	–	–	–	–
Antihistamine	Cilexetil										
	Fexofenadine	0.21	0.69	0.40	1.32	0.56	1.85	–	–	–	–
	Cetirizine	0.26	0.87	0.32	1.06	0.52	1.72	–	–	–	–
GUD/ED	Sildenafil	0.01	0.03	0.02	0.05	0.01	0.05	0.001	0.003	0.001	0.003
Diabetes	Metformin	156.0*	515.0*	163.0*	460.0*	457.0*	1509*	–	–	–	–
	Gliclazide	0.15	0.77	0.16	0.82	0.22	1.09	–	–	–	–
	Sitagliptin	0.03	0.09	0.08	0.27	0.06	0.22	0.004	0.01	0.003	0.01
Cough suppressant	Pholcodine	2.25	7.42	8.02	26.5	25.3	83.3	0.28	0.92	1.52	5.00
Beta blocker	Atenolol	0.20	0.66	0.56	1.84	0.71	2.35	0.01	0.05	0.10	0.33
	Metoprolol	0.07	0.35	0.19	0.96	0.28	1.40	0.01	0.03	0.03	0.14
	Propranolol	0.29	0.96	0.73	2.41	0.68	2.25	0.01	0.04	0.13	0.42
	Bisoprolol	0.001	0.004	0.004	0.01	0.003	0.01	0.0001	0.0005	0.0001	0.0005
H ₂ receptor agonist	Ranitidine	7.96	39.8	22.3	111.4	14.8	73.8	0.44	2.19	4.81	24.1
	Cimetidine	1.60	7.98	3.12	15.6	5.06	25.3	–	–	–	–
X-ray contrast media	Iopromide	5.97	29.9	14.1	70.6	24.5	123.0	–	–	–	–
Various	Buprenorphine	0.06	0.20	0.11	0.36	0.18	0.61	0.02	0.07	0.01	0.05
Drug precursor	Ephedrine/pseudoephedrine	0.60	1.97	1.62	5.36	1.32	4.36	0.02	0.07	0.11	0.35
	Norephedrine	0.18	8.82	0.35	17.3	0.37	18.6	0.01	0.39	0.04	1.85
Anticancer	Azathioprine	0.17	0.55	0.36	1.20	0.41	1.36	–	–	–	–
	Methotrexate	6.13	20.2	9.04	29.8	7.11	23.5	0.16	0.53	1.64	5.42
	Ifosfamide	0.08	0.40	0.24	1.22	0.31	1.53	–	–	–	–
	Tamoxifen	14.5	72.6	0.76	3.82	0.70	3.50	0.004	0.01	2.23	11.14
	Imatinib	0.88	2.93	1.13	3.76	1.78	5.95	0.10	0.35	0.06	0.21
	Capecitabine	0.01	0.02	0.01	0.03	0.01	0.03	0.002	0.01	0.001	0.003
	Bicalutamide	0.22	0.72	0.31	1.02	0.32	1.07	0.02	0.06	0.01	0.03
Anaesthetic and metabolite	Ketamine	0.07	0.37	0.19	0.93	0.24	1.20	0.005	0.02	0.03	0.17
	Norketamine	0.23	0.76	0.56	1.86	0.72	2.37	0.02	0.05	0.10	0.33
	Venlafaxine	0.07	0.37	0.24	1.20	0.37	1.83	0.01	0.03	0.08	0.38
	Desmethylvenlafaxine	0.24	0.80	0.66	2.18	0.85	2.79	0.01	0.05	0.09	0.29
	Fluoxetine	1.14	5.71	1.42	7.08	0.50	2.52	0.005	0.02	0.11	0.53
	Norfluoxetine	1.64	8.21	1.27	6.35	0.42	2.12	0.004	0.02	0.14	0.68
	Sertraline	1.61	8.07	1.21	6.05	0.74	3.72	0.002	0.01	0.17	0.86
	Mirtazapine	0.09	0.44	0.25	1.25	0.39	1.94	0.01	0.03	0.05	0.27
	Citalopram	0.61	6.08	1.41	14.1	1.24	12.4	0.02	0.24	0.16	1.64
	Desmethylcitalopram	0.14	0.69	0.36	1.82	0.31	1.54	0.01	0.03	0.05	0.24
	Paroxetine	0.18	0.59	0.21	0.69	0.13	0.45	0.01	0.02	0.005	0.02
	Duloxetine	0.04	0.13	0.05	0.18	0.04	0.12	0.003	0.01	0.002	0.01
	Amitriptyline	0.16	0.55	0.33	1.09	0.30	1.02	0.02	0.07	0.01	0.03
	Nortriptyline	0.33	1.11	0.63	2.11	0.61	2.03	0.03	0.10	0.02	0.06
Norsertaline	–	–	–	–	1.07	3.58	0.09	0.28	–	–	
Anti-epileptic	Carbamazepine	0.08	0.38	0.19	0.93	0.27	1.37	0.01	0.03	0.10	0.48
	Carbamazepine-10,11-epoxide	0.16	0.53	0.55	1.82	0.53	1.76	–	–	–	–
	10,11-Dihydro-10-hydroxycarbamazepine	0.34	3.37	0.84	8.41	0.99	9.94	0.02	0.25	0.43	4.35
Calcium channel blocker	Diltiazem	0.11	1.11	0.32	3.23	0.27	2.68	–	–	–	–
	Verapamil	0.01	0.02	0.01	0.04	0.01	0.03	0.001	0.002	0.0004	0.001

Table 2 (continued)

Class of analyte	Analyte	Surface water (ng L ⁻¹)		Effluent (ng L ⁻¹)		Influent (ng L ⁻¹)		Solid particulate matter (ng g ⁻¹)		Digested solids (ng g ⁻¹)	
		MDL	MLQ	MDL	MLQ	MDL	MLQ	MDL	MLQ	MDL	MLQ
Hypnotic	Temazepam	0.08	0.38	0.14	0.69	0.18	0.92	0.01	0.04	0.16	0.82
	Oxazepam	0.11	0.36	0.22	0.72	0.20	0.66	–	–	0.01	0.03
	Diazepam	0.02	0.06	0.04	0.13	0.04	0.12	0.002	0.01	0.002	0.01
Antipsychotic	Quetiapine	0.10	0.48	0.21	1.07	0.26	1.32	0.004	0.02	0.05	0.26
	Risperidone	0.01	0.02	0.02	0.06	0.02	0.06	0.001	0.004	0.002	0.01
Dementia	Donepezil	0.55	1.83	1.54	5.12	1.48	4.93	0.09	0.30	0.09	0.29
	Memantine	0.04	0.14	0.11	0.36	0.12	0.39	0.02	0.07	0.01	0.04
Human indicators	Creatinine	511*	1686*	771*	2544*	945*	3118*	–	–	–	–
	Nicotine	3.34	11.0	5.44	18.0	508*	2296*	0.16	–	0.66	2.19
	Caffeine	0.37	1.83	1.11	5.57	121*	581*	–	–	–	–
	Cotinine	0.07	0.35	0.21	1.06	0.27	1.34	0.005	0.02	0.24	1.22
Analgesics and metabolites	1,7-Dimethylxanthine	3.19	10.5	11.4	37.6	560*	2165*	–	–	–	–
	Morphine	2.65	8.75	6.34	20.9	8.85	29.2	0.11	0.37	1.92	6.33
	Dihydromorphine	0.11	0.55	0.32	1.59	0.05	2.51	0.01	0.04	0.09	0.45
	Normorphine	3.54	11.7	7.84	25.9	9.99	33.0	0.12	0.39	1.74	5.75
	Methadone	0.11	0.54	0.21	1.04	0.20	1.01	0.01	0.03	0.03	0.17
	EDDP	0.21	1.05	0.29	1.47	0.23	1.13	0.01	0.03	0.04	0.20
	Codeine	0.74	3.71	1.46	7.31	2.56	12.8	0.04	0.21	0.33	1.66
	Norcodeine	2.88	9.52	8.32	27.4	8.53	28.2	0.19	0.64	1.26	4.17
	Dihydrocodeine	0.23	0.75	0.55	1.83	0.88	2.89	0.02	0.05	0.11	0.36
	Tramadol	0.08	8.20	0.21	21.3	0.30	30.0	0.01	0.62	0.03	3.26
	<i>N</i> -Desmethyltramadol	0.12	5.92	0.30	15.0	0.56	27.9	0.01	0.30	0.04	2.02
	<i>O</i> -Desmethyltramadol	0.09	8.53	0.28	27.8	0.31	31.4	–	–	–	–
Stimulants and metabolites	Amphetamine	0.68	2.23	1.11	3.65	1.23	4.07	0.01	0.05	0.09	0.29
	Methamphetamine	0.32	1.05	0.71	2.35	0.95	3.13	0.01	0.04	0.09	0.30
	MDMA	0.10	0.50	0.27	1.35	0.34	1.70	0.01	0.03	0.04	0.18
	MDA	0.53	1.74	1.00	3.30	0.99	3.26	–	–	–	–
	Cocaine	0.07	0.35	0.22	1.11	0.46	2.31	0.01	0.03	0.03	0.15
	Benzoylcegonine	0.07	0.34	0.18	0.91	0.21	1.07	0.005	0.02	0.03	0.14
	Anhydrocegonine methyl ester	0.93	4.67	1.99	9.96	2.95	14.8	–	–	–	–
	Cocaethylene	0.07	0.35	0.21	1.04	1.31	6.54	0.01	0.03	0.03	0.17
	Mephedrone	0.22	1.09	0.44	2.19	0.55	2.75	0.01	0.04	0.06	0.31
	MDPV	0.04	0.22	0.12	0.59	0.48	2.41	0.01	0.03	0.04	0.20
	Opioid and metabolite	Heroin	0.92	4.62	3.44	17.2	4.18	20.9	0.05	0.25	0.56
6-Acetylmorphine		0.28	0.94	0.76	2.50	0.89	2.95	–	–	–	–
Pesticides, fungicides and herbicides	Thiamethoxam	0.13	0.42	0.44	1.46	0.53	1.76	0.01	0.03	0.01	0.02
	Imidacloprid	0.04	0.15	0.10	0.33	0.10	0.33	0.01	0.02	–	–
	Clothianidin	0.06	0.19	0.14	0.47	0.15	0.50	0.004	0.01	–	–
	Metazachlor	0.02	0.06	0.04	0.14	0.04	0.13	0.002	0.01	0.002	0.01
	Terbutylazine	0.03	0.11	0.07	0.22	0.07	0.23	0.01	0.02	0.01	0.02
	Methiocarb	0.13	0.43	0.27	0.91	0.26	0.86	0.01	0.04	0.01	0.04
	Dichlofluanid	–	–	–	–	25.2	83.8	–	–	–	–
	Flufenacet	0.01	0.04	0.02	0.07	0.02	0.07	0.002	0.01	0.002	0.01
	Oxadiazon	0.15	0.49	0.26	0.85	0.30	0.98	0.08	0.26	0.05	0.16
	Chlorpyrifos	12.9	42.9	8.54	28.5	–	–	–	–	0.33	1.09
	Triallate	0.11	0.37	0.20	0.68	0.09	0.31	–	–	–	–

Table 2 (continued)

Class of analyte	Analyte	Surface water (ng L ⁻¹)		Effluent (ng L ⁻¹)		Influent (ng L ⁻¹)		Solid particulate matter (ng g ⁻¹)		Digested solids (ng g ⁻¹)	
		MDL	QL	MDL	QL	MDL	QL	MDL	QL	MDL	QL
Veterinary pharmaceuticals	Tylosin	1.28	6.39	2.23	11.1	3.27	16.3	–	–	–	–
	Sulfapyridine	0.04	0.14	0.11	0.37	0.10	0.33	–	–	–	–
	Sarafloxacin	0.83	2.78	2.66	8.86	2.01	6.72	–	–	–	–
	Ceftiofur	2.17	7.23	1.32	4.41	1.02	3.39	–	–	–	–
	Diazinon	0.03	0.11	0.07	0.23	0.06	0.21	0.00	0.01	0.003	0.01

*Calculated for direct injection

influent, as this is generally used for plant protection, particularly against slugs. It is thought to be an unlikely compound to make its way into influent and thought to enter surface waters via direct application to the environment. Its presence in the sewage treatment works, although in a low concentration, is notable as it may indicate incorrect disposal, although a higher concentration would be expected from this. It was recently banned for use as a molluscicide in 2014, the grace period of which ended the month before these samples were collected; however, it could still be used as an insect repellent and seed treatment [46]. Other pesticides found in influent at this site include imidacloprid and flufenacet, which are more widely used for vegetable and fruit crops and may be due to rinsing of food prior to consumption. However, this needs further detailed investigation.

Of further interest is the presence of veterinary pharmaceuticals in wastewater, and the presence of these can be justified by considering the number of household pets and the possible disposal route of pet waste down the lavatory. However, it is interesting to note that sulfapyridine (1339.5 ± 330.3 ng L⁻¹) and ceftiofur (451.8 ± 129.2 ng L⁻¹) are usually reserved for the use with individual pigs and cows. These levels suggest potential herd applications, incorrect disposal or unknown contribution of livestock wastewater to this WwTW. Further investigation is needed to determine the source and persistence of these levels.

For the solid portion of influent (SPM), 64% of the 98 analytes quantifiable with this method were found. Only a small fraction of the total concentration can be found in the SPM, as most CECs, particularly pharmaceuticals, prefer to partition to the aqueous phase. For example, only a fraction of the total concentration (6.4%) of bisphenol A ($3.64 \log K_{ow}$) appears in SPM. However, in this case, the concentrations are so high in influent that this results in bisphenol A contributing a large portion of the total measurable concentration of CECs. Concentrations range from an average concentration of 0.1 ng L⁻¹ to 1383 ng L⁻¹ in SPM (converted to ng L⁻¹ for simple comparison to influent concentrations). Ketoconazole prefers to partition to the solid phase with 31% higher average concentration present in SPM. This is not surprising when considering its $\log K_{ow}$ of 4.45. The antidepressants fluoxetine ($4.65 \log K_{ow}$) and amitriptyline ($4.95 \log K_{ow}$) partition

partially with SPM concentrations at 82% and 73% to that in influent. The results for fluoxetine are in line with the results published by Baker and Kasprzyk-Hordern [47]. However, sorption of amitriptyline is far higher in this study. This to be expected, as influent is highly variable, and many factors can affect sorption to solids. Therefore, it is important to analyse both the liquid and solid compartments of this matrix.

There are fewer CECs in effluent than in influent. Of the 137 analytes quantifiable in effluent, 62% were found at the WwTW. Generally, these results show lower concentrations after treatment. However, imatinib, *O*-desmethytramadol, carbamazepine and its epoxide metabolite and venlafaxine and its metabolite increase in effluent. For metabolites, this may be due to degradation/metabolism of the parent compound during treatment. For the parent compounds, this may be due to the undetected presence of conjugated metabolites in influent transforming back to the parent compound. This phenomenon requires more detailed investigation. Comparison of the influent and effluent concentrations can provide data on treatment efficiencies. Trickling filters are employed at the WwTW for the treatment of influent. The results show poor removal of imidacloprid, tramadol, *N*-desmethytramadol, bicalutamide, ranitidine, cetirizine and fexofenadine. Acetaminophen, on the other hand, shows high removal of 99%, similar to what is often seen in the literature for conventional activated sludge treatment [11, 48]. Caffeine, its metabolite 1,7-dimethylxanthine and metformin are highly removed but still at high concentrations of 2.9 ± 0.4 µg L⁻¹, 8.4 ± 1.1 µg L⁻¹ and 14.1 ± 0.9 µg L⁻¹, respectively.

In the surface waters downstream of the WwTW, 51% of the 138 analytes that can be quantified are in this matrix. Of particular noteworthiness is the antidiabetic metformin, which is found at levels > 3000 ng L⁻¹. Metformin has very high usage, and after administration, 100% of the dose is excreted unchanged [49]. Caldwell et al. [50] calculated the PNEC of metformin to be 1 mg L⁻¹, based on critical evaluation of previously published work. This suggests that despite the high environmental levels, metformin is currently of low risk to this catchment. Within the literature, many compounds can be found to have different PNECs due to a lack of consistency in (a) assessment factors

Table 3 Weekly average concentrations of ECs found in several matrices at site A (mean concentration \pm variation from mean across the week)

Class of analyte	Compound	Surface water (ng L ⁻¹)	Effluent (ng L ⁻¹)	Influent (ng L ⁻¹)	SPM (ng L ⁻¹)	Digested solids (ng g ⁻¹)
UV filter	Benzophenone-1	<MQL	<MQL	949.2 \pm 1134.5	13.8 \pm 12.0	N/A
	Benzophenone-2	<MQL	62.0 \pm 36.0	1701.0 \pm 2176.6	22.0 \pm 25.9	10.2 \pm 5.2
	Benzophenone-3	19.0 \pm 1.9	96.7 \pm 36.7	2051.5 \pm 897.3	N/A	N/A
	Benzophenone-4	607.9 \pm 344.8	4224.7 \pm 3455.7	15,817.4 \pm 13,404.6	9.9 \pm 4.0	<MQL
Parabens	Methylparaben	7.0 \pm 2.5	42.3 \pm 49.4	23,401.4 \pm 36,282.5	122.1 \pm 99.8	360.2 \pm 132.5
	Ethylparaben	<MQL	<MQL	2166.0 \pm 2427.8	21.0 \pm 17.8	<MQL
	Propylparaben	3.7 \pm 0.7	32.1 \pm 11.4	5456.2 \pm 3190.2	33.2 \pm 32.6	<MQL
Plasticizer	Butylparaben	<MQL	<MQL	337.4 \pm 434.5	40.3 \pm 66.7	<MQL
	Bisphenol-A	38.1 \pm 15.1	473.6 \pm 256.2	20,395.5 \pm 15,740.7	1383.9 \pm 1204.0	4365.6 \pm 309.6
	E1	<MQL	6.3 \pm 1.6	47.0 \pm 10.7	<MQL	41.5 \pm 2.9
Steroid estrogens	E2	<MQL	<MQL	<MQL	<MQL	<MQL
	EE2	<MQL	<MQL	<MQL	N/A	N/A
	Sulfasalazine	53.6 \pm 6.8	241.4 \pm 27.5	612.9 \pm 248.7	N/A	N/A
Antibiotics and antibacterials	Clarithromycin	102.0 \pm 26.8	776.2 \pm 121.4	1321.1 \pm 293.4	N/A	N/A
	Azithromycin	3.9 \pm 3.2	108.7 \pm 48.4	793.1 \pm 256.1	<MQL	5.7 \pm 5.5
	Trimethoprim	25.9 \pm 5.7	175.6 \pm 30.2	626.2 \pm 151.9	16.7 \pm 6.3	11.1 \pm 2.8
	Sulfamethoxazole	33.2 \pm 3.1	54.3 \pm 6.2	126.7 \pm 58.6	<MQL	<MQL
	Triclosan	<MQL	220.9 \pm 52.2	2480.5 \pm 801.6	N/A	N/A
	Amoxicillin	N/A	<MQL	N/A	N/A	N/A
	Metronidazole	10.1 \pm N/A	37.0 \pm 6.5	81.6 \pm 23.6	<MQL	<MQL
Antifungal	Sulfadiazine	<MQL	<MQL	<MQL	8.5 \pm 3.3	<MQL
	Cefalexin	N/A	N/A	N/A	N/A	N/A
	Ofloxacin	71.5 \pm 29.5	108.7 \pm 52.7	<MQL	N/A	N/A
	Ciprofloxacin	N/A	<MQL	N/A	N/A	N/A
	Tetracycline	<MQL	<MQL	251.5 \pm 70.5	N/A	N/A
	Danofloxacin	<MQL	<MQL	372.3 \pm 414.1	N/A	<MQL
	Oxytetracycline	<MQL	<MQL	N/A	N/A	N/A
	Chloramphenicol	<MQL	157.9 \pm 23.3	<MQL	<MQL	377.6 \pm 46.3
	Penicillin G	N/A	<MQL	<MQL	N/A	N/A
	Penicillin V	<MQL	<MQL	<MQL	<MQL	N/A
	Erythromycin	2148.2 \pm 168.5	6503.8 \pm 619.1	13,591.5 \pm 2445.9	N/A	N/A
	Prulifloxacin	N/A	<MQL	N/A	N/A	N/A
	Norfloxacin	N/A	N/A	N/A	N/A	N/A
Grisofulvin	<MQL	<MQL	<MQL	1.4 \pm N/A	<MQL	

Table 3 (continued)

Class of analyte	Compound	Surface water (ng L ⁻¹)	Affluent (ng L ⁻¹)	Influent (ng L ⁻¹)	SPM (ng L ⁻¹)	Digested solids (ng g ⁻¹)
Hypertension	Ketoconazole	46.5 ± 6.1	64.3 ± 5.2	131.0 ± 40.9	171.6 ± 48.3	929.3 ± 212.1
	Valsartan	<MQL	236.9 ± 23.0	854.9 ± 306.9	N/A	N/A
NSAIDs	Ibuprofen	111.4 ± 14.6	287.9 ± 30.0	380.7 ± 73.3	N/A	N/A
	Lisinopril	<MQL	93.9 ± 16.1	890.5 ± 273.9	12.3 ± 2.7	<MQL
	Ketoprofen	<MQL	<MQL	<MQL	<MQL	<MQL
	Ibuprofen	60.5 ± 21.3	1029.2 ± 357.7	18,840.6 ± 5449.2	42.8 ± 15.9	231.8 ± 39.9
Lipid regulator	Naproxen	234.9 ± 47.7	1452.5 ± 213.7	15,238.5 ± 3111.2	51.7 ± 17.9	108.4 ± 5.3
	Diclofenac	88.2 ± 7.3	345.6 ± 52.7	916.2 ± 367.3	<MQL	26.2 ± 3.4
Antihypertlipidaemic	Acetaminophen	193.0 ± 49.7	1840.7 ± 449.3	393,559.0 ± 100,020.4	28.1 ± 23.7	<MQL
	Bezafibrate	103.6 ± 15.8	591.7 ± 76.7	2038.3 ± 468.2	2.8 ± 1.2	7.5 ± 0.7
Antihyperintensive	Atorvastatin	44.4 ± 8.1	172.9 ± 27.0	1338.2 ± 456.4	N/A	N/A
	Gemfibrozil	N/A	<MQL	N/A	N/A	849.6 ± 183.8
Antihistamine	Candesartan	N/A	N/A	N/A	N/A	N/A
	Cilexetil	N/A	N/A	N/A	N/A	N/A
GUD/ED	Fexofenadine	209.1 ± 31.9	732.6 ± 103.2	1012.1 ± 689.0	N/A	N/A
	Cetirizine	239.6 ± 35.8	1051.1 ± 154.5	1088.4 ± 106.6	N/A	N/A
Diabetes	Sildenafil	2.2 ± 1.3	25.5 ± 12.4	29.0 ± 41.6	2.6 ± 2.9	20.8 ± 1.0
	Metformin	3607.1 ± 413.6	14,050.8 ± 897.2	117,927.4 ± 24,688.7	N/A	N/A
Cough suppressant	Gliclazide	47.3 ± 7.6	84.3 ± 7.5	136.4 ± 41.6	N/A	N/A
	Sitagliptin	142.5 ± 13.5	441.4 ± 67.2	524.7 ± 76.8	11.2 ± 2.9	27.6 ± 2.1
Beta blocker	Pholcodine	<MQL	<MQL	<MQL	<MQL	<MQL
	Atenolol	53.9 ± 10.6	328.9 ± 31.4	2581.7 ± 584.4	23.3 ± 9.6	<MQL
H ₂ receptor agonist	Metoprolol	<MQL	13.3 ± 2.9	32.1 ± 9.3	<MQL	<MQL
	Propranolol	20.6 ± 3.5	106.6 ± 6.1	182.3 ± 30.9	21.5 ± 6.4	198.5 ± 22.9
X-ray contrast media	Bisoprolol	<MQL	14.7 ± N/A	30.6 ± 29.3	<MQL	<MQL
	Ranitidine	148.2 ± 31.8	1163.0 ± 98.1	1545.0 ± 277.6	<MQL	<MQL
Various	Cimetidine	<MQL	<MQL	170.3 ± 64.4	N/A	N/A
	Iopromide	<MQL	<MQL	<MQL	N/A	N/A
Drug precursor	Buprenorphine	6.3 ± 0.2	5.7 ± N/A	60.8 ± 37.7	0.6 ± 0.2	13.2 ± 4.5
	Ephedrine/pseudoephedrine	25.6 ± 2.0	115.5 ± 18.6	655.5 ± 248.6	4.5 ± 1.3	<MQL
Anticancer	Norephedrine	<MQL	<MQL	<MQL	<MQL	<MQL
	Azathioprine	<MQL	<MQL	<MQL	N/A	N/A
Anticancer	Methotrexate	<MQL	<MQL	<MQL	<MQL	<MQL
	Ifosfamide	<MQL	<MQL	<MQL	N/A	N/A
	Tamoxifen	<MQL	<MQL	<MQL	<MQL	<MQL

Table 3 (continued)

Class of analyte	Compound	Surface water (ng L ⁻¹)	Effluent (ng L ⁻¹)	Influent (ng L ⁻¹)	SPM (ng L ⁻¹)	Digested solids (ng g ⁻¹)
Anaesthetic and metabolite	Imatinib	38.3 ± 1.3	143.3 ± 39.6	88.8 ± 14.4	47.6 ± 17.1	123.0 ± 40.3
	Capecitabine	<MQL	<MQL	7.3 ± 2.8	<MQL	<MQL
	Bicalutamide	59.4 ± 1.1	115.6 ± 6.7	154.3 ± 10.0	7.5 ± 2.0	50.0 ± 11.6
	Ketamine	11.8 ± 2.7	75.4 ± 11.3	166.2 ± 65.6	0.4 ± 0.2	2.8 ± 0.5
	Norketamine	<MQL	5.7 ± 2.8	12.2 ± 7.0	<MQL	0.8 ± 0.2
	Venlafaxine	85.5 ± 13.0	509.5 ± 135.8	462.4 ± 137.3	8.6 ± 2.0	126.3 ± 26.0
	Desmethylvenlafaxine	229.2 ± 23.9	981.0 ± 128.6	733.9 ± 129.3	2.9 ± 1.1	22.4 ± 1.8
	Fluoxetine	1.2 ± 0.1	46.2 ± 4.0	53.2 ± 7.1	43.7 ± 32.1	193.7 ± 33.2
	Norfluoxetine	<MQL	<MQL	31.8 ± 6.5	9.6 ± 2.2	91.4 ± 26.5
	Sertraline	<MQL	13.1 ± 1.7	46.3 ± 7.0	114.1 ± 34.8	565.2 ± 71.3
	Mirtazapine	4.2 ± 1.6	41.3 ± 2.4	78.2 ± 13.1	6.9 ± 2.3	69.1 ± 8.7
	Citalopram	<MQL	321.9 ± 19.9	557.5 ± 109.9	87.5 ± 25.7	782.9 ± 75.2
	Desmethylcitalopram	12.0 ± 4.5	119.7 ± 10.7	216.2 ± 59.1	31.7 ± 8.8	295.7 ± 62.6
	Paroxetine	<MQL	<MQL	N/A	1.8 ± N/A	3.2 ± 2.9
Anti-epileptic	Duloxetine	<MQL	<MQL	N/A	<MQL	17.5 ± 5.8
	Amitriptyline	10.3 ± 3.2	53.3 ± 6.9	159.4 ± 25.2	115.6 ± 26.3	471.9 ± 44.1
	Nortriptyline	7.8 ± 4.0	30.6 ± 3.0	9.6 ± 4.2	9.9 ± 3.0	63.2 ± 15.2
	Norsertaline	N/A	N/A	157.8 ± N/A	133.8 ± 37.8	N/A
	Carbamazepine	160.4 ± 16.8	626.8 ± 60.5	521.2 ± 115.3	7.2 ± 2.6	118.7 ± 11.5
	Carbamazepine-10, 11-epoxide	30.9 ± 4.2	138.2 ± 33.3	98.6 ± 29.9	N/A	N/A
	10,11-Dihydro-10-hydroxycarbamazepine	5.7 ± 2.1	43.2 ± 13.1	145.1 ± 55.1	<MQL	<MQL
	Diltiazem	4.4 ± 2.0	59.5 ± 10.3	246.2 ± 55.5	N/A	N/A
	Verapamil	<MQL	<MQL	<MQL	5.8 ± 2.5	51.1 ± 6.1
	Temazepam	<MQL	21.5 ± 13.9	11.1 ± 7.4	<MQL	<MQL
Hypnotic	Oxazepam	2.0 ± 0.6	17.3 ± 4.1	16.6 ± 3.9	N/A	<MQL
	Diazepam	6.2 ± 0.8	11.6 ± 1.2	<MQL	0.9 ± 0.2	4.2 ± 0.4
	Quetiapine	<MQL	<MQL	63.0 ± 31.4	2.9 ± 1.6	17.3 ± 4.7
	Risperidone	1.2 ± 1.4	2.6 ± 0.8	1.0 ± N/A	0.1 ± 0.1	<MQL
Dementia	Donepezil	<MQL	<MQL	N/A	0.9 ± 0.3	8.3 ± 0.2
	Memantine	<MQL	<MQL	<MQL	<MQL	<MQL
	Creatinine	<MQL	<MQL	1,341,104.8 ± 314,821.7	N/A	N/A
Human indicators	Nicotine	35.9 ± 16.1	183.9 ± 36.8	3461.2 ± 688.7	49.4 ± 26.3	227.1 ± 42.7
	Caffeine	426.5 ± 91.2	2938.9 ± 401.5	118,596.4 ± 24,610.9	N/A	N/A

Table 3 (continued)

Class of analyte	Compound	Surface water (ng L ⁻¹)	Effluent (ng L ⁻¹)	Influent (ng L ⁻¹)	SPM (ng L ⁻¹)	Digested solids (ng g ⁻¹)
Analgesics and metabolites	Cotinine	37.8 ± 5.7	214.5 ± 28.6	3028.3 ± 729.6	9.9 ± 2.9	40.1 ± 4.9
	1,7-Dimethylxanthine	923.5 ± 208.7	8352.1 ± 1091.1	153,776.3 ± 36,949.3	N/A	N/A
	Morphine	<MQL	141.7 ± 18.0	1366.1 ± 276.4	10.1 ± 3.8	63.2 ± 43.3
	Dihydromorphine	<MQL	<MQL	133.6 ± 32.7	<MQL	<MQL
	Normorphine	<MQL	<MQL	186.3 ± 50.7	<MQL	<MQL
	Methadone	3.0 ± 0.8	22.1 ± 1.5	41.6 ± 7.1	1.6 ± 0.5	12.9 ± 1.8
	EDDP	10.0 ± 1.5	50.3 ± 3.0	64.8 ± 12.9	11.0 ± 5.1	38.3 ± 22.3
	Codeine	85.6 ± 18.3	577.9 ± 58.6	2570.4 ± 555.0	23.3 ± 5.9	45.0 ± 42.5
	Norcodeine	<MQL	<MQL	181.0 ± 30.9	<MQL	<MQL
	Dihydrocodeine	22.1 ± 4.6	154.0 ± 19.3	445.2 ± 119.4	3.4 ± 0.9	24.8 ± 3.8
	Tramadol	321.7 ± 27.7	1273.7 ± 166.3	1247.5 ± 308.1	3.7 ± 1.0	34.5 ± 7.7
	N-Desmethyltramadol	255.3 ± 24.7	975.5 ± 118.3	1076.5 ± 355.9	1.5 ± 0.5	17.1 ± 5.6
	O-Desmethyltramadol	297.0 ± 32.9	1300.4 ± 200.0	861.6 ± 133.9	N/A	N/A
	Amphetamine	<MQL	<MQL	478.8 ± 135.7	<MQL	<MQL
Stimulants and metabolites	Methamphetamine	<MQL	9.1 ± 0.8	13.9 ± 3.7	<MQL	<MQL
	MDMA	4.3 ± 1.9	59.4 ± 33.3	183.2 ± 127.7	1.4 ± 0.9	6.1 ± 1.0
	MDA	10.0 ± 0.9	30.1 ± 7.6	45.5 ± N/A	N/A	N/A
	Cocaine	4.9 ± 1.6	47.1 ± 11.7	754.6 ± 295.3	6.9 ± 3.2	<MQL
	Benzoyllecgonine	45.5 ± 9.9	291.2 ± 95.7	2078.4 ± 925.3	1.9 ± 1.1	<MQL
	Anhydroecgonine methyl ester	<MQL	<MQL	<MQL	N/A	N/A
	Cocacethylene	<MQL	3.6 ± 1.1	29.3 ± 18.5	0.4 ± 0.3	<MQL
	Mephedrone	<MQL	<MQL	<MQL	<MQL	<MQL
	MDPV	<MQL	<MQL	<MQL	<MQL	<MQL
	Heroin	<MQL	<MQL	<MQL	<MQL	<MQL
	6-Acetyl morphine	<MQL	<MQL	<MQL	N/A	N/A
	Thiamethoxam	<MQL	<MQL	<MQL	<MQL	<MQL
	Imidacloprid	53.0 ± 18.5	346.6 ± 133.7	339.2 ± 178.8	<MQL	N/A
	Clofianidin	<MQL	<MQL	<MQL	1.8 ± 0.8	N/A
Metazachlor	4.5 ± N/A	<MQL	<MQL	4.2 ± 4.1	<MQL	
Terbutylazine	<MQL	<MQL	<MQL	<MQL	<MQL	
Methiocarb	<MQL	3.0 ± N/A	3.7 ± 0.6	0.8 ± 0.6	1.7 ± 0.0	
Dichlofluanid	N/A	N/A	N/A	N/A	N/A	
Flufenacet	24.0 ± 2.4	61.8 ± 3.5	57.1 ± 12.9	4.1 ± 0.7	14.7 ± 0.5	

Table 3 (continued)

Class of analyte	Compound	Surface water (ng L ⁻¹)	Effluent (ng L ⁻¹)	Influent (ng L ⁻¹)	SPM (ng L ⁻¹)	Digested solids (ng g ⁻¹)
Veterinary pharmaceuticals	Oxadiazon	16.9 ± 3.4	29.8 ± 3.4	<MQL	9.3 ± 0.4	<MQL
	Chlorpyrifos	<MQL	<MQL	N/A	N/A	93.8 ± N/A
	Triallate	<MQL	<MQL	<MQL	N/A	N/A
	Tylosin	<MQL	<MQL	<MQL	N/A	N/A
	Sulfapyridine	128.7 ± 19.6	576.0 ± 121.2	1339.5 ± 330.3	N/A	N/A
	Sarafloxacin	<MQL	<MQL	N/A	N/A	N/A
	Ceftiofur	31.2 ± N/A	216.1 ± 77.2	451.8 ± 129.2	N/A	N/A
	Diazinon	6.1 ± 0.3	12.0 ± 0.3	17.5 ± 0.4	11.9 ± 3.9	15.2 ± 3.6

used, (b) limited sources or databases used in studying ecotoxicity data and (c) varying criteria in accepting ecotoxicity study results. This has highlighted the need to harmonise these methods to ensure PNECs are calculated consistently, to provide comparable comparisons between studies and to be clearer of the risks CECs pose to the environment.

Regarding digested solids, the concentrations are in ng g⁻¹ of solid material, which cannot be directly compared to the concentrations of the other matrices, as it is a combination of sludge from various parts of the wastewater treatment process. However, it is an important consideration as a source of CECs in the environment, due to subsequent direct application to the land in agricultural practices. Of the 96 analytes that can be quantified with this method, 55% were found at this site. Antidepressants are high in concentration, with average concentrations of five of the analytes between 126.3 and 782.9 ng g⁻¹. Of the more industrial CECs, methylparaben and bisphenol A are present in solids at high levels. In particular, bisphenol A has been quantified at levels exceeding 4000 ng g⁻¹ at this site. Gemfibrozil, an antihyperlipidaemic, although not found in other matrices at this site, was quantified in digested solids at a concentration of 849.6 ± 183.8 ng g⁻¹. This suggests accumulation of gemfibrozil in other sections of the WwTW that were not analysed. Due to the low MQL, the lack of incoming concentration during the sampling period suggests occasional loads high in gemfibrozil before the study.

Conclusions

This work presents a validated multi-residue method for the analysis of 195 compounds in five matrices (3 liquid and 2 solid). These CECs cover a variety classes, both veterinary and human pharmaceuticals, industrial chemical, personal care products and pesticides. Application of the method to environmental matrices has shown that the method is appropriate for assessing treatment efficiency, partitioning to solids, and environmental concentrations. Discussion of the results has identified several key areas regarding environmental risk assessment, e.g. PNECs that need to be addressed; however, that is outside the scope of this paper. The achieved MDL and MQL concentrations appear low enough to be used to assess the environmental risk of these CECs. The results show a need for analysing both the liquid and solid phases within a WwTW; however, it also indicates a need to monitor all outgoing treated waste materials, e.g. effluent and digested sludge. This was due to the appearance of gemfibrozil, which was not present at quantifiable levels in any other matrix at this site. Overall, this method is suitable to be used in catchment-based exposure-driven studies to further increase knowledge of the contribution of CECs by WwTWs to the environment and the risk they pose.

Acknowledgements The support of Wessex Water Services Ltd. and the University of Bath EPSRC Impact Acceleration Account (Project number: EP/K503897/1) is greatly appreciated. All data supporting this study are provided as supporting information accompanying this paper.

Funding information Kathryn Proctor was funded by NERC and Wessex Water for her CASE studentship (NE/L009579/1).

Compliance with ethical standards

Conflict of interest The authors declare that they have no conflict of interest.

Open Access This article is distributed under the terms of the Creative Commons Attribution 4.0 International License (<http://creativecommons.org/licenses/by/4.0/>), which permits unrestricted use, distribution, and reproduction in any medium, provided you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons license, and indicate if changes were made.

References

- Vorkamp K, Rigét FF. A review of new and current-use contaminants in the Arctic environment: evidence of long-range transport and indications of bioaccumulation. *Chemosphere*. 2014;111:379–95.
- Zenker A, Cicero MR, Prestinaci F, et al. Bioaccumulation and biomagnification potential of pharmaceuticals with a focus to the aquatic environment. *J Environ Manag*. 2014;133:378–87.
- Puckowski A, Mioduszewska K, Łukaszewicz P, et al. Bioaccumulation and analytics of pharmaceutical residues in the environment: a review. *J Pharm Biomed Anal*. 2016;127:232–55.
- Lagesson A, Fahlman J, Brodin T, et al. Bioaccumulation of five pharmaceuticals at multiple trophic levels in an aquatic food web—insights from a field experiment. *Sci Total Environ*. 2016;568:208–15.
- Godoy AA, Kummrow F, Augusto P, et al. Chemosphere occurrence, ecotoxicological effects and risk assessment of antihypertensive pharmaceutical residues in the aquatic environment—a review. *Chemosphere*. 2015;138:281–91.
- Sousa JCG, Ribeiro AR, Barbosa MO, et al. A review on environmental monitoring of water organic pollutants identified by EU guidelines. *J Hazard Mater*. 2018;344:146–62.
- Miller TH, Bury NR, Owen SF, et al. A review of the pharmaceutical exposome in aquatic fauna. *Environ Pollut*. 2018;239:129–46.
- Stuart M, Lapworth D, Crane E, et al. Review of risk from potential emerging contaminants in UK groundwater. *Sci Total Environ*. 2012;416:1–21.
- Ramakrishnan A, Blaney L, Kao J, et al. Emerging contaminants in landfill leachate and their sustainable management. *Environ Earth Sci*. 2014;73:1357–68.
- Castiglioni S, Valsecchi S, Polesello S, et al. Sources and fate of perfluorinated compounds in the aqueous environment and in drinking water of a highly urbanized and industrialized area in Italy. *J Hazard Mater*. 2015;282:51–60.
- Verlicchi P, Al Aukidy M, Zambello E. Occurrence of pharmaceutical compounds in urban wastewater: removal, mass load and environmental risk after a secondary treatment—a review. *Sci Total Environ*. 2012;429:123–55.
- Petrie B, Barden R, Kasprzyk-Hordern B. A review on emerging contaminants in wastewaters and the environment: current knowledge, understudied areas and recommendations for future monitoring. *Water Res*. 2014;72:3–27.
- Petrie B, Lopardo L, Proctor K, et al. Assessment of bisphenol-A in the urban water cycle. *Sci Total Environ*. 2019;650:900–7.
- Lopardo L, Petrie B, Proctor K, et al. Estimation of community-wide exposure to bisphenol A via water fingerprinting. *Environ Int*. 2019;125:1–8.
- aus der Beek T, Weber FA, Bergmann A, et al. Pharmaceuticals in the environment—global occurrences and perspectives. *Environ Toxicol Chem*. 2016;35:823–35.
- Gogoi A, Mazumder P, Tyagi VK, et al. Occurrence and fate of emerging contaminants in water environment: a review. *Groundw Sustain Dev*. 2018;6:169–80.
- US EPA National Health and Environmental Effects Research Laboratory's Mid-Continent Ecology Division (MED). EPA ECOTOX Knowledgebase, <https://cfpub.epa.gov/ecotox/index.cfm> (2000, accessed 13 April 2019).
- Brausch JM, Connors KA, Brooks BW, et al. Human pharmaceuticals in the aquatic environment: a review of recent toxicological studies and considerations for toxicity testing. *Rev Environ Contam Toxicol*. 2012;218:1–99.
- Brausch JM, Rand GM. A review of personal care products in the aquatic environment: environmental concentrations and toxicity. *Chemosphere*. 2011;82:1518–32.
- Ermler S, Scholze M, Kortenkamp A. The suitability of concentration addition for predicting the effects of multi-component mixtures of up to 17 anti-androgens with varied structural features in an in vitro AR antagonist assay. *Toxicol Appl Pharmacol*. 2011;257:189–97.
- Orton F, Ermler S, Kugathas S, et al. Mixture effects at very low doses with combinations of anti-androgenic pesticides, antioxidants, industrial pollutant and chemicals used in personal care products. *Toxicol Appl Pharmacol*. 2014;278:201–8.
- Loos R, Gawlik BM, Locoro G, et al. EU-wide survey of polar organic persistent pollutants in European river waters. *Environ Pollut*. 2015;157:561–8.
- ECHA EU. Registration statistics infograph—ECHA, <https://echa.europa.eu/registration-statistics-infograph#> (2009, accessed 13 April 2019).
- Ccancapa A, Masiá A, Navarro-Ortega A, et al. Pesticides in the Ebro River basin: occurrence and risk assessment. *Environ Pollut*. 2016;211:414–24.
- Brack W, Dulio V, Ågerstrand M, et al. Towards the review of the European Union Water Framework Directive: recommendations for more efficient assessment and management of chemical contamination in European surface water resources. *Sci Total Environ*. 2017;576:720–37.
- Ort C, Lawrence MG, Reungoat J, et al. Sampling for PPCPs in wastewater systems: comparison of different sampling modes and optimization strategies. *Environ Sci Technol*. 2010;44:6289–96.
- Geissen V, Mol H, Klump E, et al. Emerging pollutants in the environment: a challenge for water resource management. *Int Soil Water Conserv Res*. 2015;3:57–65.
- Gustavsson BM, Magnér J, Carney Almroth B, et al. Chemical monitoring of Swedish coastal waters indicates common exceedances of environmental thresholds, both for individual substances as well as their mixtures. *Mar Pollut Bull*. 2017;122:409–19.
- Gardner M. The National Chemical Investigations Programme 2015-2020—volume 2: monitoring of substances of emerging concern (18/EQ/01/13): UKWIR; 2018.
- Comber S, Gardner M, Sörme P, et al. Active pharmaceutical ingredients entering the aquatic environment from wastewater treatment works: a cause for concern? *Sci Total Environ*. 2018;613–614:538–47.
- Gardner M, Jones V, Comber S, et al. Performance of UK wastewater treatment works with respect to trace contaminants. *Sci Total Environ*. 2013;456:359–69.
- European Commission. Directive 2000/60/EC of the European Parliament and of the Council of 23 October 2000 establishing a

- framework for community action in the field of water policy. Off J Eur Parliam. 2000;L327:1–82.
33. Roos V, Gunnarsson L, Fick J, et al. Prioritising pharmaceuticals for environmental risk assessment: towards adequate and feasible first-tier selection. *Sci Total Environ*. 2012;421–422:102–10.
 34. Guo J, Sinclair CJ, ABA B. Risk based prioritisation of pharmaceuticals: UKWIR; 2014.
 35. NHSBSA. NHSBSA practice prescribing dispensing data. NHSBSA Copyright, <https://www.nhsbsa.nhs.uk/prescription-data/dispensing-data/dispensing-contractors-data> (2019, accessed 21 May 2019).
 36. Wishart DS, Feunang YD, Guo AC, et al. DrugBank 5.0: a major update to the DrugBank database for 2018. *Nucleic Acids Research*. 2018;46(D1):D1074–D1082.
 37. US EPA. Estimation Programs Interface (EPI) Suite (TM) for Microsoft(R) Windows, v4.11, <http://www.epa.gov/oppt/exposure/pubs/episuite.htm> (2015).
 38. ACD/Labs (C) 1994–2015. Advanced Chemistry Development (ACD/Labs) software V11.02.
 39. Petrie B, Youdan J, Barden R, et al. Multi-residue analysis of 90 emerging contaminants in liquid and solid environmental matrices by ultra-high-performance liquid chromatography tandem mass spectrometry. *J Chromatogr A*. 2016;1431:64–78.
 40. Kasprzyk-Hordern B, Kondakal VVR, Baker DR. Enantiomeric analysis of drugs of abuse in wastewater by chiral liquid chromatography coupled with tandem mass spectrometry. *J Chromatogr A*. 2010;1217:4575–86.
 41. Official Journal of the European Communities. Commission Decision implementing Council Directive 96/23/EC concerning the performance of analytical methods and the interpretation of results (2002/657/EC). Off J Eur Communities 2002; 8–36.
 42. Campos-Mañas MC, Plaza-Bolaños P, Sánchez-Pérez JA, et al. Fast determination of pesticides and other contaminants of emerging concern in treated wastewater using direct injection coupled to highly sensitive ultra-high performance liquid chromatography-tandem mass spectrometry. *J Chromatogr A*. 2017;1507:84–94.
 43. Hermes N, Jewell KS, Wick A, et al. Quantification of more than 150 micropollutants including transformation products in aqueous samples by liquid chromatography-tandem mass spectrometry using scheduled multiple reaction monitoring. *J Chromatogr A*. 2018;1531:64–73.
 44. Daughton CG. Monitoring wastewater for assessing community health: sewage chemical-information mining (SCIM). *Sci Total Environ*. 2018;619–620:748–64.
 45. Brewer AJ, Ort C, Banta-Green CJ, et al. Normalized diurnal and between-day trends in illicit and legal drug loads that account for changes in population. *Environ Sci Technol*. 2012;46:8305–14.
 46. European Commission. Commission Implementing Regulation (EU) No 187/2014 of 26 February 2014 amending Implementing Regulation (EU) No 540/2011 as regards the conditions of approval of the active substance methiocarb. 2014. <https://doi.org/10.2903/j.efsa.2012.2758>.
 47. Baker DR, Kasprzyk-Hordern B. Multi-residue analysis of drugs of abuse in wastewater and surface water by solid-phase extraction and liquid chromatography-positive electrospray ionisation tandem mass spectrometry. *J Chromatogr A*. 2011;1218:1620–31.
 48. Wilkinson J, Hooda PS, Barker J, et al. Occurrence, fate and transformation of emerging contaminants in water: an overarching review of the field. *Environ Pollut*. 2017;231:954–70.
 49. Scheurer M, Sacher F, Brauch H-J. Occurrence of the antidiabetic drug metformin in sewage and surface waters in Germany. *J Environ Monit*. 2009;11:1608–13.
 50. Caldwell DJ, D'Aco V, Davidson T, et al. Environmental risk assessment of metformin and its transformation product guanylurea:

II. Occurrence in surface waters of Europe and the United States and derivation of predicted no-effect concentrations. *Chemosphere*. 2019;216:855–65.

Publisher's note Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.



Behaviour of Toxic Organic Chemicals in the Waste Stream).

Kathryn Proctor is currently finishing her PhD, joint funded by NERC and Wessex Water, at the University of Bath. She has been working for several years on developing and applying methods for exposure and effects-driven analysis (UPLC-MS/MS, toxicity assays) of chemicals of emerging concern in the environment, both in Bath and abroad in Argentina as part of INTERWASTE (Synergising International Research Studies into the Environmental Fate and



Bruce Petrie is Lecturer in Analytical Chemistry at Robert Gordon University, Aberdeen, with research interests in the sources, fate and impact of trace environmental pollutants.



Ruth Barden is Director of Environmental Strategy at Wessex Water and a Visiting Industrial Fellow at the University of Bath. Her key areas of interest are environmental regulation, emerging issues and innovative wastewater treatment options.



Tom Arnot is Senior Lecturer in Chemical Engineering at the University of Bath. He is Co-Director of the Water Innovation & Research Centre at Bath and Co-Director of the EPSRC-funded Water Informatics: Science & Engineering Centre for Doctoral Training. His research interests focus on a circular economy approach to wastewater treatment, including MBR and membrane technology for water recycling, nutrient recovery, anaerobic digestion for production

of biogas and feedstock chemicals.



Barbara Kasprzyk-Hordern is Professor in Environmental and Analytical Chemistry in the Department of Chemistry, University of Bath. Her principal research interests fall into the three interrelated research areas of environmental, analytical and water sciences.

Analytical and Bioanalytical Chemistry

Electronic Supplementary Material

Multi-residue ultra-performance liquid chromatography coupled with tandem mass spectrometry method for comprehensive multi-class anthropogenic compounds of emerging concern analysis in a catchment-based exposure-driven study

Kathryn Proctor, Bruce Petrie, Ruth Barden, Tom Arnot, Barbara Kasprzyk-Hordern

Contents	Page
<i>Table S1</i> General and chemical information (Class, Chemical, CAS No. Mol Formula, Molecular Weight, Water sol. Log Kow, Log Koc, Log Dow.)	2
<i>Table S2</i> MS/MS detection parameters for analytes of interest sorted by class, includes analytes from ESI –ve and ESI +ve methods	8
<i>Table S3</i> MS/MS detection parameters for internal standards sorted by analogue analyte class, includes internal standards from ESI –ve and ESI +ve methods	13
<i>Figure S1</i> Chromatograms of all analytes in spiked influent sample from the WwTW using methods: A-ESI+, B-ESI+ and C-ESI-	14
<i>Table S4</i> Solid sample weights	15
<i>Table S5</i> Suspended solids weights (g L^{-1})	15
<i>Table S6</i> Adsorbance	16
<i>Table S7</i> Instrument parameters (LC MSMS): Intra- and Inter-day precision and accuracy	20
<i>Table S8</i> Liquid matrix absolute and corrected recoveries	24
<i>Table S9</i> Solid matrix absolute and corrected recoveries	30
<i>Figure S2</i> Absolute recoveries for all matrices, error bars indicate the range of standard deviation	34
<i>Table S10</i> Matrix suppression	35
<i>Table S11</i> Surface water data	39
<i>Table S12</i> Effluent data	45
<i>Table S13</i> Influent data	51
<i>Table S14</i> Solid particulate matter data	57
<i>Table S15</i> Digested solids data	63

Table S1 General and chemical information (Class, Chemical, CAS No. Mol Formula, Molecular Weight, Water sol. Log Kow, Log Koc, Log Dow.)

Compound	CAS no.	Molecular Formula	Molecular Weight	Water Solubility (mg/L) ^j	Log Kow ^h	Log Koc ⁱⁱⁱ	Log Dow ^{iv}	Henry's Law Constant (atm m ³ mol ⁻¹) ^v	Vapour Pressure (Torr) ^{vi}	pKa (Most acidic) ^{vii}	pKa (Most basic) ^{viii}	Supplier
<i>Benzophenone-1</i>	131-56-6	C ₁₃ H ₁₀ O ₃	214.22	413	2.96	3.092	2.9	2.65x10 ⁻¹¹	2.84x10 ⁻⁷	7.72±0.35	-	Sigma Aldrich
<i>Benzophenone-2</i>	131-55-5	C ₁₃ H ₁₀ O ₅	246.22	399	2.78	3.326	2.35	3.61x10 ⁻¹⁶	6.69x10 ⁻¹²	6.98±0.35	-	Sigma Aldrich
<i>Benzophenone-3</i>	131-57-7	C ₁₄ H ₁₂ O ₃	228.25	68.6	3.52	3.44	3.06	1.50x10 ⁻⁸	5.26x10 ⁻⁶	7.56±0.35	-	Sigma Aldrich
<i>Benzophenone-4</i>	4065-45-6	C ₁₄ H ₁₂ O ₆ S	308.31	2.03x10 ⁴	0.37	1.548	-0.53	7.06x10 ⁻¹⁵	-	-0.70±0.50	-	Sigma Aldrich
<i>Methylparaben</i>	99-76-3	C ₈ H ₈ O ₃	152.15	5.98x10 ³	2	2.11	1.63	3.61x10 ⁻⁹	5.55x10 ⁻³	8.31±0.13	-	Sigma Aldrich
<i>Ethylparaben</i>	120-47-8	C ₉ H ₁₀ O ₃	166.18	1.89x10 ³	2.49	2.393	1.99	4.79x10 ⁻⁹	7.59x10 ⁻⁴	8.31±0.13	-	Sigma Aldrich
<i>Propylparaben</i>	94-13-3	C ₁₀ H ₁₂ O ₃	180.21	529	2.98	2.708	2.51	4.25x10 ⁻¹²	9.30x10 ⁻⁴	8.23±0.15	-	Sigma Aldrich
<i>Butylparaben</i>	94-26-8	C ₁₁ H ₁₄ O ₃	194.23	159	3.47	3.001	2.96	6.00x10 ⁻⁹	3.56x10 ⁻⁴	8.22±0.15	-	Sigma Aldrich
<i>Bisphenol A</i>	80-05-7	C ₁₅ H ₁₆ O ₂	228.29	173	3.64	3.095	4.04	9.16x10 ⁻¹²	5.34x10 ⁻⁷	10.29±0.10	-	Sigma Aldrich
<i>E1</i>	53-16-7	C ₁₈ H ₂₂ O ₂	270.37	147	3.43	3.019	4.31	3.80x10 ⁻¹⁰	1.54x10 ⁻⁸	10.25±0.40	-	Sigma Aldrich
<i>E2</i>	50-28-2	C ₁₈ H ₂₄ O ₂	272.39	82.0	3.94	2.899	3.74	1.41x10 ⁻¹²	9.82x10 ⁻⁹	10.27±0.60	-	Sigma Aldrich
<i>EE2</i>	57-63-6	C ₂₀ H ₂₄ O ₂	296.41	116	4.12	2.710	3.9	7.94x10 ⁻¹²	3.74x10 ⁻⁹	10.24±0.60	-	Sigma Aldrich
<i>Sulfasalazine</i>	599-79-1	C ₁₈ H ₁₄ N ₄ O ₅ S	398.4	2.40	3.81	3.177	-1.76	2.19x10 ⁻¹⁸	5.95x10 ⁻²⁰	2.70±0.10	0.90±0.10	Sigma Aldrich
<i>Clarithromycin</i>	81103-11-9	C ₃₈ H ₆₀ NO ₁₃	747.97	0.30	3.18	1.371	2.31	1.73x10 ⁻²⁹	5.06x10 ⁻³⁰	13.08±0.70	8.16±0.70	Sigma Aldrich
<i>Azithromycin</i>	83905-01-5	C ₃₈ H ₇₂ N ₂ O ₁₂	749	0.06	3.24	1.676	-2.48	5.30x10 ⁻²⁹	2.51x10 ⁻³¹	13.28±0.70	8.59±0.70	Sigma Aldrich
<i>Trimethoprim</i>	738-70-5	C ₁₄ H ₁₈ N ₄ O ₃	290.32	2.33x10 ³	0.73	1.896	1.13	2.39x10 ⁻¹⁴	3.74x10 ⁻¹¹	-	7.04±0.10	Sigma Aldrich
<i>Sulfamethoxazole</i>	723-46-6	C ₁₀ H ₁₁ N ₃ O ₃ S	253.28	3.94 x10 ³	0.48	1.536	-0.03	9.56x10 ⁻¹³	1.87x10 ⁻⁹	5.81±0.50	1.39±0.10	Sigma Aldrich
<i>Triclosan</i>	3380-34-5	C ₁₂ H ₇ Cl ₃ O ₂	289.55	4.60	4.66	3.925	4.76	2.13x10 ⁻⁸	3.26x10 ⁻⁵	7.80±0.35	-	Sigma Aldrich
<i>Amoxicillin</i>	26787-78-0	C ₁₆ H ₁₉ N ₃ O ₅ S	365.4	3.43 x10 ³	0.97	0.709	-2.36	2.49x10 ⁻²¹	3.39x10 ⁻²¹	2.44±0.50	7.14±0.10	Sigma Aldrich
<i>Metronidazole</i>	443-48-1	C ₆ H ₉ N ₃ O ₃	171.16	2.57x10 ⁴	0	1.12	-0.46	1.69x10 ⁻¹¹	2.67x10 ⁻⁷	14.44±0.10	2.58±0.34	Sigma Aldrich
<i>Sulfadiazine</i>	68-35-9	C ₁₀ H ₁₀ N ₄ O ₂ S	250.28	2.81x10 ⁴	-0.34	1.392	-0.1	1.58x10 ⁻¹⁰	1.28x10 ⁻¹⁰	6.81±0.10	1.64±0.10	Sigma Aldrich (VETANAL)
<i>Cefalexin</i>	15686-71-2	C ₁₆ H ₁₇ N ₃ O ₄ S	347.39	1.79 x10 ³	0.4	0.421	-2.19	2.77x10 ⁻¹⁷	3.27x10 ⁻²²	3.12±0.50	6.84±0.10	Sigma Aldrich
<i>Ofloxacin</i>	82419-36-1	C ₁₈ H ₂₀ FN ₃ O ₄	361.37	2.83x10 ⁴	-0.2	-0.004	-0.39	4.98x10 ⁻²⁰	6.70x10 ⁻¹⁴	5.19±0.40	7.37±0.42	Sigma Aldrich
<i>Ciprofloxacin</i>	85721-33-1	C ₁₇ H ₁₈ FN ₃ O ₃	331.34	1.15 x10 ⁴	0	-0.004	-0.94	5.09x10 ⁻¹⁹	2.24x10 ⁻¹⁴	6.43±0.41	8.68±0.10	Sigma Aldrich (Fluka)
<i>Tetracycline</i>	60-54-8	C ₂₂ H ₂₄ N ₂ O ₈	444.43	3.88x10 ³	-1.33	-0.128	-2.35	4.66x10 ⁻²⁴	2.40x10 ⁻²⁶	4.50±1.00	11.02±0.70	Sigma Aldrich

Compound	CAS no.	Molecular Formula	Molecular Weight	Water Solubility (mg/L) ⁱ	Log Kow ⁱⁱ	Log Koc ⁱⁱⁱ	Log Dow ^{iv}	Henry's Law Constant (atm m ³ mol ⁻¹) ^v		Vapour Pressure (Torr) ^{vi}	pKa (Most acidic) ^{vii}	pKa (Most basic) ^{vii}	Supplier
								Log Dow ^{iv}	Log Dow ^{iv}				
Danofloxacin	112398-08-0	C ₁₉ H ₂₀ FN ₃ O ₃	357.38	-	-	-	-	-	-	8.41x10 ⁻¹⁴	6.43±0.41	9.00±0.20	Sigma Aldrich
Oxytetracycline	79-57-2	C ₂₂ H ₂₄ N ₂ O ₉	460.43	1.40x10 ³	-2.87	0.093	-3.35	1.70x10 ⁻²⁵	6.27x10 ⁻³⁰	6.27x10 ⁻³⁰	4.50±1.00	10.80±0.70	Sigma Aldrich
Chloramphenicol	56-75-7	C ₁₁ H ₁₂ Cl ₂ N ₂ O ₅	323.13	389	1.14	0.926	-2.19	2.29x10 ⁻¹⁸	1.63x10 ⁻¹⁷	1.63x10 ⁻¹⁷	11.03±0.46	-1.73±0.70	Sigma Aldrich
Penicillin G	61-33-6	C ₁₆ H ₁₈ N ₂ O ₄ S	334.39	210	1.85	1.095	-2.45	1.16x10 ⁻¹⁴	1.69x10 ⁻¹⁸	1.69x10 ⁻¹⁸	2.45±0.50	-1.32±0.60	Sigma Aldrich
Penicillin V	87-08-1	C ₁₆ H ₁₈ N ₂ O ₅ S	350.39	101	1.87	1.295	-2.77	4.42x10 ⁻¹⁵	1.69x10 ⁻¹⁹	1.69x10 ⁻¹⁹	2.44±0.50	-1.68±0.60	Sigma Aldrich
Erythromycin	114-07-8	C ₃₇ H ₆₇ NO ₁₃	733.93	0.52	2.48	1.406	1.66	5.42x10 ⁻²⁹	4.94x10 ⁻³¹	4.94x10 ⁻³¹	13.09±0.70	8.16±0.70	Sigma Aldrich
Prulifloxacin	123447-62-1	C ₂₁ H ₂₀ FN ₃ O ₆ S	461.46	-	-	-	-	-	6.62x10 ⁻¹⁷	6.62x10 ⁻¹⁷	5.85±0.40	7.66±0.10	Sigma Aldrich
Norfloxacin	70458-96-7	C ₁₆ H ₁₈ FN ₃ O ₃	319.33	1.78x10 ⁵	-0.31	-0.392	-1.05	8.70x10 ⁻¹⁹	3.45x10 ⁻¹³	3.45x10 ⁻¹³	0.16±0.20	8.68±0.10	Sigma Aldrich
Griseofulvin	126-07-8	C ₁₇ H ₁₇ ClO ₆	352.77	82.0	1.92	2.543	2.17	2.35x10 ⁻¹¹	5.07x10 ⁻¹³	5.07x10 ⁻¹³	-	-	Sigma Aldrich
Ketoconazole	65277-42-1	C ₂₆ H ₂₈ Cl ₂ N ₄ O ₄	531.43	0.09	4.45	3.492	4.14	5.59x10 ⁻²⁰	1.39x10 ⁻²²	1.39x10 ⁻²²	-	6.88±0.12	Sigma Aldrich
Valsartan	137862-53-4	C ₂₄ H ₂₉ N ₅ O ₃	435.53	1.41	3.65	2.2	1.15	1.82x10 ⁻¹⁸	1.06x10 ⁻¹⁹	1.06x10 ⁻¹⁹	3.56±0.10	0.60±0.10	Sigma Aldrich
Irbesartan	138402-11-6	C ₂₅ H ₂₈ N ₆ O	428.53	-	-	-	4.97	-	1.05x10 ⁻¹⁶	1.05x10 ⁻¹⁶	4.16±0.10	2.60±0.20	LGC
Lisinopril	76547-98-3	C ₂₁ H ₃₁ N ₃ O ₅	405.5	8.60	-0.94	-0.516	-3.23	1.89x10 ⁻²²	1.14x10 ⁻¹⁸	1.14x10 ⁻¹⁸	2.18±0.10	10.50±0.10	LGC
Ketoprofen	22071-15-4	C ₁₆ H ₁₄ O ₃	254.29	120	3	2.077	0.08	2.12x10 ⁻¹¹	3.32x10 ⁻⁸	3.32x10 ⁻⁸	4.23±0.10	-	Sigma Aldrich
Ibuprofen	15687-27-1	C ₁₃ H ₁₈ O ₂	206.29	41.1	3.79	2.352	1.25	1.52x10 ⁻⁷	1.39x10 ⁻⁴	1.39x10 ⁻⁴	4.41±0.10	-	Sigma Aldrich
Naproxen	22204-53-1	C ₁₄ H ₁₄ O ₃	230.27	145	3.1	1.971	-0.54	3.39x10 ⁻¹⁰	3.01x10 ⁻⁷	3.01x10 ⁻⁷	4.84±0.30	-	TRC
Diclofenac	15307-86-5	C ₁₄ H ₁₁ Cl ₂ NO ₂	296.15	4.50	4.02	2.607	0.73	4.73x10 ⁻¹²	1.59x10 ⁻⁷	1.59x10 ⁻⁷	4.18±0.10	-2.26±0.50	Sigma Aldrich
Acetaminophen	103-90-2	C ₈ H ₉ NO ₂	151.17	3.04x10 ⁴	0.27	1.321	0.9	6.42x10 ⁻¹³	1.43x10 ⁻⁶	1.43x10 ⁻⁶	9.86±0.13	1.72±0.50	Sigma Aldrich
Bezafibrate	41859-67-0	C ₁₉ H ₂₀ ClNO ₄	361.83	1.20	4.25	2.311	0.46	2.12x10 ⁻¹⁵	6.29x10 ⁻¹⁴	6.29x10 ⁻¹⁴	3.29±0.10	-2.06±0.70	Sigma Aldrich
Atorvastatin	134523-00-5	C ₃₃ H ₃₅ FN ₂ O ₅	558.66	1.12x10 ⁻³	6.36	2.6	1.86	2.41x10 ⁻²³	6.84x10 ⁻²²	6.84x10 ⁻²²	4.29±0.10	0.38±0.50	LGC
Gemfibrozil	25812-30-0	C ₁₅ H ₂₂ O ₃	250.33	4.96	4.77	2.85	0.86	1.19x10 ⁻⁸	6.13x10 ⁻⁷	6.13x10 ⁻⁷	4.75±0.45	-	Sigma Aldrich
Candesartan Cilexetil	145040-37-5	C ₃₃ H ₃₄ N ₆ O ₆	610.66	-	-	-	1.19	-	1.11x10 ⁻²⁸	1.11x10 ⁻²⁸	4.16±0.10	3.94±0.10	Sigma Aldrich
Fexofenadine	83799-24-0	C ₃₂ H ₃₉ NO ₄	501.67	0.02	2.81	0.822	2.93	1.19x10 ⁻¹⁸	2.08x10 ⁻²⁰	2.08x10 ⁻²⁰	4.43±0.10	9.42±0.10	LGC
Cetirizine	83881-51-0	C ₂₁ H ₂₅ ClN ₂ O ₃	388.9	1.10	-0.61	0.875	0.69	4.19x10 ⁻¹⁷	1.39x10 ⁻¹²	1.39x10 ⁻¹²	3.46±0.10	6.71±0.10	LGC
Sildenafil	139755-83-2	C ₂₂ H ₃₀ N ₆ O ₄ S	474.58	11.3	2.3	2.661	1.13	7.22x10 ⁻²¹	-	-	10.05±0.20	6.03±0.42	Sigma Aldrich (Cerilliant)
Meformin	657-24-9	C ₄ H ₁₁ N ₃	129.17	1.00x10 ⁶	-2.64	-0.666	-6.45	7.64x10 ⁻¹⁶	1.33	1.33	-	12.27±0.10	Sigma Aldrich
Gliclazide	21187-98-4	C ₁₃ H ₂₁ N ₃ O ₃ S	323.41	138	2.12	2.205	0.79	7.95x10 ⁻¹³	-	-	6.07±0.10	3.89±0.20	LGC

Compound	CAS no.	Molecular Formula	Molecular Weight	Water Solubility (mg/L) ⁱ	Log Kow ⁱⁱ	Log Koc ⁱⁱⁱ	Log Dow ^{iv}	Henry's Law Constant (atm m ³ mol ⁻¹) ^v	Vapour Pressure (Torr) ^{vi}	pKa (Most acidic) ^{vii}	pKa (Most basic) ^{vii}	Supplier
Staglipitin	486460-32-6	C ₁₆ H ₁₅ F ₆ N ₅ O	407.31	179	1.39	1.698	-0.51	2.35x10 ⁻¹⁵	2.59x10 ⁻¹¹	7.2±0.10	-	TRC
Pholcodine	509-67-1	C ₂₃ H ₃₀ N ₂ O ₄	398.51	1.01x10 ⁴	0.59	0.731	-0.48	3.42x10 ⁻¹⁹	3.44x10 ⁻¹⁴	13.40±0.20	8.22±0.40	Sigma Aldrich
Atenolol	29122-68-7	C ₁₄ H ₂₂ N ₂ O ₃	266.34	685	-0.03	0.611	-1.71	1.37x10 ⁻¹⁸	3.82x10 ⁻¹¹	13.88±0.20	9.43±0.10	Sigma Aldrich
Metoprolol	51384-51-1	C ₁₅ H ₂₅ N ₃ O ₃	267.37	4.77x10 ³	1.69	1.475	-0.38	1.40x10 ⁻¹³	4.52x10 ⁻⁷	13.89±0.20	9.43±0.10	Sigma Aldrich
Propranolol	525-66-6	C ₁₆ H ₂₁ N ₂ O ₂	259.35	228	2.6	2.451	0.45	7.98x10 ⁻¹³	2.48x10 ⁻⁸	13.84±0.20	9.50±0.30	Sigma Aldrich
Bisoprolol	66722-44-9	C ₁₈ H ₃₁ N ₃ O ₄	325.44	2.24 x10 ³	1.84	1.379	0.06	2.89x10 ⁻¹⁵	1.06x10 ⁻⁸	13.86±0.20	9.42±0.10	Sigma Aldrich
Ranitidine	66357-35-5	C ₁₃ H ₂₂ N ₄ O ₃ S	314.41	2.47x10 ⁴	0.29	1.141	-1.56	3.42x10 ⁻¹⁵	7.66x10 ⁻⁸	-	8.35±0.28	Sigma Aldrich
Cimetidine	51481-61-9	C ₁₀ H ₁₆ N ₆ S	252.34	1.05x10 ⁴	0.57	1.85	-0.3	9.55x10 ⁻¹⁶	3.13x10 ⁻⁹	14.13±0.10	7.07±0.61	Sigma Aldrich
Topromide	73334-07-3	C ₁₈ H ₂₄ N ₃ O ₈	791.12	23.8	-2.49	-1.672	-0.44	1.00x10 ⁻²⁸	5.00x10 ⁻³⁰	10.62±0.70	-2.60±0.70	LGC
Buprenorphine	52485-79-7	C ₂₉ H ₄₁ N ₃ O ₄	467.64	0.66	4.9	3.335	1.37	7.78x10 ⁻¹⁸	-	9.47±0.60	8.31±0.60	Sigma Aldrich
Ephedrine/pseudoephedrine	299-42-3	C ₁₀ H ₁₅ NO	165.24	7.15x10 ⁴	0.68	1.095	-1.13	8.65x10 ⁻¹¹	8.65x10 ⁻³	13.96±0.20	9.38±0.10	Sigma Aldrich
Norephedrine	492-39-7	C ₉ H ₁₃ NO	151.21	1.49x10 ⁵	0.22	0.951	-1.4	3.94x10 ⁻¹¹	1.10x10 ⁻³	12.07±0.45	8.47±0.10	Sigma Aldrich
Azathioprine	446-86-6	C ₉ H ₇ N ₃ O ₂ S	277.26	272	-0.09	2.395	1.21	2.64x10 ⁻¹⁵	5.94x10 ⁻¹¹	-	7.47±0.20	Sigma Aldrich
Methotrexate	59-05-2	C ₂₀ H ₂₃ N ₅ O ₅	454.45	2.60x10 ³	-1.28	-0.387	-7.06	1.54x10 ⁻³¹	-	3.47±0.10	5.56±0.10	LGC
Ifosfamide	3778-73-2	C ₇ H ₁₅ C ₃ N ₂ O ₂ P	261.09	3.78x10 ³	0.97	1.439	0.1	1.36x10 ⁻¹¹	1.15x10 ⁻⁴	-	1.44±0.20	Sigma Aldrich
Tamoxifen	10540-29-1	C ₂₆ H ₂₉ NO	371.53	0.20	6.3	4.4	5.07	4.49x10 ⁻¹⁰	1.85x10 ⁻⁹	-	8.69±0.28	LGC
Imatinib	152459-95-5	C ₂₉ H ₃₁ N ₇ O	493.6	-	-	-	2.04	-	-	13.28±0.70	7.55±0.10	Sigma Aldrich
Capecitabine	154361-50-9	C ₁₅ H ₂₂ FN ₃ O ₆	359.35	1.82x10 ³	0.56	0.173	0.74	2.92x10 ⁻¹⁹	-	5.41±0.40	1.75±0.20	Sigma Aldrich
Bicalutamide	90357-06-5	C ₁₈ H ₁₄ F ₄ N ₂ O ₄ S	430.37	11.8	2.3	2.177	2.71	2.82x10 ⁻¹⁵	8.47x10 ⁻¹⁸	11.49±0.29	-3.75±0.50	Sigma Aldrich
Ketamine	6740-88-1	C ₁₃ H ₁₆ CINO	237.73	3.87x10 ³	3.12	2.283	3.18	1.38x10 ⁻⁸	1.76x10 ⁻⁵	-	6.46±0.20	Sigma Aldrich
Norketamine	35211-10-0	C ₁₅ H ₁₆ CIN ₂ O	233.7	-	-	-	2.71	1.78x10 ⁻¹⁰	1.26x10 ⁻⁵	-	6.25±0.20	Sigma Aldrich
Yenlafaxine	93413-69-5	C ₁₇ H ₂₇ N ₁ O ₂	277.41	267	3.28	2.318	1.32	2.87x10 ⁻¹¹	4.92x10 ⁻⁷	14.84±0.20	9.26±0.28	Sigma Aldrich
Desmethylvenlafaxine	93413-62-8	C ₁₆ H ₂₅ N ₂ O ₂	263.38	-	-	-	1.17	-	3.03x10 ⁻⁷	10.04±0.26	9.33±0.28	Sigma Aldrich
Fluoxetine	54910-89-3	C ₁₇ H ₁₈ F ₃ NO	309.33	60.3	4.65	3.05	1.92	8.90x10 ⁻⁸	1.88x10 ⁻⁶	-	10.05±0.10	LGC (Cerilliant)
Norfluoxetine	83891-03-6	C ₁₆ H ₁₆ F ₃ NO	295.3	-	-	-	1.54	-	5.21x10 ⁻⁶	-	9.05±0.13	LGC (Cerilliant)
Serrraline	79617-96-2	C ₁₇ H ₁₇ Cl ₂ N	306.24	3.50	5.29	3.808	3.11	5.10x10 ⁻⁸	3.85x10 ⁻⁷	-	9.47±0.40	LGC
Mirtazapine	85650-52-8	C ₁₇ H ₁₉ N ₃	265.35	-	-	-	3.15	-	1.11x10 ⁻⁷	-	8.10±0.20	Sigma Aldrich

Compound	CAS no.	Molecular Formula	Molecular Weight	Water Solubility (mg/L) ⁱ	Log Kow ⁱⁱ	Log Koc ⁱⁱⁱ	Log Dow ^{iv}	Henry's Law Constant (atm m ³ mol ⁻¹) ^v	Vapour Pressure (Torr) ^{vi}	pKa (Most acidic) ^{vii}	pKa (Most basic) ^{vii}	Supplier
Citalopram	59729-33-8	C ₂₀ H ₂₁ FN ₃ O	324.4	31.1	3.74	3.23	1.5	2.69x10 ⁻¹¹	1.53x10 ⁻⁷	-	9.57±0.28	Sigma Aldrich
Desmethylcitalopram	62498-67-3	C ₁₉ H ₁₉ FN ₃ O	310.37	-	-	-	0.14	-	1.40x10 ⁻⁷	-	10.50±0.10	TRC
Paroxetine	61869-08-7	C ₁₉ H ₂₀ FN ₃ O ₃	329.37	35.3	3.95	3.088	0.66	1.78x10 ⁻¹²	2.39x10 ⁻⁸	9.68±0.10	-	Sigma Aldrich
Duloxetine	116539-59-4	C ₁₈ H ₁₉ NOS	297.41	13.0	4.68	3.526	1.77	5.43x10 ⁻¹⁰	7.23x10 ⁻⁹	10.02±0.10	-	LGC
Amitriptyline	50-48-6	C ₂₀ H ₂₃ N	277.4	0.82	4.95	3.581	2.57	6.85x10 ⁻⁸	1.50x10 ⁻⁶	9.18±0.28	-	Sigma Aldrich
Nortriptyline	72-69-5	C ₁₉ H ₂₁ N	263.38	2.22	4.74	3.376	1.64	3.12x10 ⁻⁸	1.02x10 ⁻⁶	10.00±0.10	-	Sigma Aldrich
Norserrtraline	87857-41-8	C ₁₆ H ₁₅ Cl ₂ N	292.2	-	-	-	1.68	-	6.41x10 ⁻⁷	9.13±0.40	-	LGC
Carbamazepine	298-46-4	C ₁₅ H ₁₂ N ₂ O	236.28	17.7	2.25	2.227	2.77	1.08x10 ⁻¹⁰	5.78x10 ⁻⁷	13.94±0.20	-0.49±0.20	Sigma Aldrich
Carbamazepine 10,11-epoxide	36507-30-9	C ₁₅ H ₁₂ N ₂ O ₂	252.27	-	-	-	1.97	-	2.69x10 ⁻⁶	13.91±0.20	-0.50±0.20	LGC
10,11-Dihydro-10-hydroxycarbamazepine	29331-92-8	C ₁₅ H ₁₄ N ₂ O ₂	254.28	-	-	-	1.73	-	3.33x10 ⁻⁸	13.75±0.20	-0.53±0.40	LGC
Diltiazem	42399-41-7	C ₂₂ H ₂₆ N ₂ O ₄ S	414.52	12.3	2.79	2.296	1.97	8.61x10 ⁻¹⁷	4.27x10 ⁻¹⁴	-	8.94±0.28	Sigma Aldrich
Verapamil	52-53-9	C ₂₇ H ₃₈ N ₂ O ₄	454.6	4.47	4.8	3.46	2.42	8.79x10 ⁻¹⁵	1.01x10 ⁻¹³	8.97±0.50	-	LGC
Temazepam	846-50-4	C ₁₆ H ₁₃ ClN ₂ O ₂	300.75	164	2.15	1.678	2.79	1.13x10 ⁻⁸	6.33x10 ⁻¹³	11.66±0.40	1.58±0.50	Sigma Aldrich (Cerilliant)
Oxazepam	604-75-1	C ₁₅ H ₁₁ ClN ₂ O ₂	286.71	179	2.32	1.727	2.92	5.53x10 ⁻¹⁰	4.40x10 ⁻¹¹	10.94±0.70	1.17±0.50	Sigma Aldrich (Cerilliant)
Diazepam	439-14-5	C ₁₆ H ₁₃ ClN ₂ O	284.74	58.8	2.7	2.438	3.08	3.64x10 ⁻⁹	4.98x10 ⁻¹⁰	3.40±0.10	-	Sigma Aldrich (Cerilliant)
Quetiapine	111974-69-7	C ₂₁ H ₂₅ N ₃ O ₂ S	384.52	0.60	3.17	2.111	2.67	7.45x10 ⁻¹⁸	3.22x10 ⁻¹³	14.41±0.10	6.74±0.10	LGC
Risperidone	106266-06-2	C ₂₃ H ₂₇ FN ₄ O ₂	410.48	2.76	3.49	2.743	0.86	2.17x10 ⁻¹⁶	4.14x10 ⁻¹³	8.07±0.10	-	LGC
Donepezil	120014-06-4	C ₂₀ H ₂₉ NO ₃	379.49	2.93	4.86	3.855	0.71	1.22x10 ⁻¹²	3.11x10 ⁻¹¹	8.84±0.10	-	LGC
Memantine	19982-08-2	C ₁₂ H ₂₁ N	179.3	894	3.34	2.382	-0.97	1.47x10 ⁻⁵	3.93x10 ⁻²	10.79±0.60	-	Sigma Aldrich
Creatinine	60-27-5	C ₄ H ₇ N ₃ O	113.12	1.66x10 ⁵	-1.21	-0.39	-1.07	2.42x10 ⁻¹²	0.31	-	6.89±0.20	Sigma Aldrich
Nicotine	54-11-5	C ₁₀ H ₁₄ N ₂	162.24	1.00x10 ⁶	1	1.683	-2.34	3.00x10 ⁻⁹	3.03x10 ⁻²	-	8.00±0.50	Sigma Aldrich
Caffeine	58-08-2	C ₈ H ₁₀ N ₄ O ₂	194.19	2.63x10 ³	0.16	0.98	-0.55	3.58x10 ⁻¹¹	3.72x10 ⁻⁷	-	0.52±0.70	Sigma Aldrich
Cotinine	486-56-6	C ₁₀ H ₁₂ N ₂ O	176.22	9.99x10 ⁵	0.34	1.093	0.21	3.33x10 ⁻¹²	4.21x10 ⁻⁴	-	4.72±0.12	Sigma Aldrich (Cerilliant)
1,7-dimethylxanthine	611-59-6	C ₇ H ₈ N ₄ O ₂	180.17	4.14x10 ³	-0.39	0.919	0.24	1.75x10 ⁻¹²	-	8.50±0.50	0.21±0.70	Sigma Aldrich
Morphine	57-27-2	C ₁₇ H ₁₉ NO ₃	285.35	2.64x10 ⁴	0.72	1.163	-0.37	1.33x10 ⁻¹⁶	7.06x10 ⁻¹⁰	9.48±0.40	8.25±0.40	Sigma Aldrich (Cerilliant)

Compound	CAS no.	Molecular Formula	Molecular Weight	Water Solubility (mg/L) ⁱ	Log Kow ⁱⁱ	Log Koc ⁱⁱⁱ	Log Dow ^{iv}	Henry's Law Constant (atm m ³ mol ⁻¹) ^v	Vapour Pressure (Torr) ^{vi}	pKa (Most acidic) ^{vii}	pKa (Most basic) ^{vii}	Supplier
<i>Dihydromorphine</i>	509-60-4	C ₁₇ H ₂₁ NO ₃	287.36	2.38x10 ⁴	0.93	1.185	-0.26	1.51x10 ⁻¹⁶	7.65x10 ⁻¹⁰	9.56±0.40	8.44±0.40	Sigma Aldrich (Cerilliant)
<i>Normorphine</i>	466-97-7	C ₁₆ H ₁₇ NO ₃	271.32	2.56x10 ⁵	0.5	0.599	-2.05	6.07x10 ⁻¹⁷	2.99x10 ⁻¹⁰	9.17±0.40	9.54±0.40	Sigma Aldrich (Cerilliant)
<i>Methadone</i>	76-99-3	C ₂₁ H ₂₇ NO	309.46	48.5	4.17	3.229	3.39	4.97x10 ⁻¹⁰	2.20x10 ⁻⁷	-	9.50±0.50	Sigma Aldrich (Cerilliant)
<i>EDDP</i>	30223-73-5	C ₂₀ H ₂₃ N	277.4	-	-	-	2.5	-	2.06x10 ⁻⁶	-	7.71±0.60	LGC (Cerilliant)
<i>Codeine</i>	76-57-3	C ₁₈ H ₂₁ NO ₃	299.37	1.22x10 ⁴	1.28	1.218	-0.23	7.58x10 ⁻¹⁴	2.47x10 ⁻⁹	13.40±0.20	8.23±0.40	Sigma Aldrich
<i>Norcodeine</i>	467-15-2	C ₁₇ H ₁₉ NO ₃	285.35	3.92x10 ⁴	1.07	0.964	-2.28	3.45x10 ⁻¹⁴	1.51x10 ⁻⁹	13.34±0.20	9.28±0.40	Sigma Aldrich (Cerilliant)
<i>Dihydrocodeine</i>	125-28-0	C ₁₈ H ₂₃ NO ₃	301.39	6.53x10 ³	1.49	1.384	-0.11	8.61x10 ⁻¹⁴	2.48x10 ⁻⁹	14.22±0.20	8.43±0.40	Sigma Aldrich
<i>Tramadol</i>	27203-92-5	C ₁₆ H ₂₅ NO ₂	263.38	1.15x10 ³	3.01	1.959	0.72	1.54x10 ⁻¹¹	1.02x10 ⁻⁶	14.47±0.40	9.61±0.28	Sigma Aldrich
<i>N-desmethyiltramadol</i>	75377-45-6 144830-15-9(S,S)/ 144830-14-8(R,R)	C ₁₅ H ₂₃ NO ₂	249.35	-	-	-	-0.27	-	9.16x10 ⁻⁶	14.46±0.40	10.56±0.10	LGC
<i>O-desmethyiltramadol</i>		C ₁₅ H ₂₃ NO ₂	249.35	-	-	-	0.38	-	3.14x10 ⁻⁷	10.00±0.10	9.61±0.28	LGC
<i>Amphetamine</i>	300-62-9	C ₉ H ₁₃ N	135.21	2.80x10 ⁴	1.76	1.877	-1.23	1.08x10 ⁻⁶	0.31	-	9.94±0.10	LGC (Cerilliant)
<i>Methamphetamine</i>	537-46-2	C ₁₀ H ₁₅ N	149.24	1.33x10 ⁴	2.22	2.027	-1.01	2.37x10 ⁻⁶	0.15	-	10.38±0.10	LGC (Cerilliant)
<i>MDMA</i>	42542-10-9	C ₁₁ H ₁₅ NO ₂	193.25	7.03x10 ³	2.28	2.183	-1.38	2.75x10 ⁻⁹	3.17x10 ⁻³	-	10.32±0.10	LGC
<i>MDA</i>	101-77-9	C ₁₃ H ₁₄ N ₂	198.27	1.96x10 ³	2.18	1.761	-1.61	1.58x10 ⁻¹¹	1.52x10 ⁻⁶	-	5.32±0.25	LGC (Cerilliant)
<i>Cocaine</i>	50-36-2	C ₁₇ H ₂₁ NO ₄	303.36	1.30x10 ³	2.17	2.001	0.92	4.24x10 ⁻¹¹	1.87x10 ⁻⁶	-	8.97±0.60	LGC (Cerilliant)
<i>Benzoylcegonine</i>	519-09-5	C ₁₆ H ₁₉ NO ₄	289.33	1.61x10 ³	-1.32	-0.705	-0.6	1.03x10 ⁻¹³	1.32x10 ⁻⁸	3.35±0.40	10.83±0.40	Sigma Aldrich
<i>Anhydrocegoninemethyl ester</i>	43021-26-7	C ₁₀ H ₁₅ NO ₂	181.23	-	-	-	0.21	-	2.19x10 ⁻²	-	7.97±0.40	Sigma Aldrich (Cerilliant)
<i>Cocacethylene</i>	529-38-4	C ₁₈ H ₂₃ NO ₄	317.38	-	-	-	1.35	-	6.80x10 ⁻⁷	-	9.04±0.60	Sigma Aldrich (Cerilliant)
<i>Mephedrone</i>	1189805-46-6	C ₁₁ H ₁₅ NO	177.24	-	-	-	1.47	-	3.84x10 ⁻³	-	7.41±0.10	Sigma Aldrich (Cerilliant)
<i>MDPV</i>	687603-66-3	C ₁₆ H ₂₁ NO ₃	275.34	-	-	-	2.77	-	4.09x10 ⁻⁷	-	8.41±0.20	Sigma Aldrich (Cerilliant)
<i>Heroin</i>	561-27-3	C ₂₁ H ₂₃ NO ₅	369.42	2.15x10 ³	1.8	1.658	-0.06	6.15x10 ⁻¹³	7.38x10 ⁻¹⁰	-	7.93±0.40	Sigma Aldrich (Cerilliant)
<i>6-acetylmorphine</i>	2784-73-8	C ₁₉ H ₂₁ NO ₄	327.37	-	-	-	-0.22	-	1.83x10 ⁻⁹	9.41±0.40	8.03±0.40	Sigma Aldrich (Cerilliant)

Compound	CAS no.	Molecular Formula	Molecular Weight	Water Solubility (mg/L) ⁱ	Log Kow ⁱⁱ	Log Koc ⁱⁱⁱ	Log Dow ^{iv}	Henry's Law Constant (atm m ³ mol ⁻¹) ^v	Vapour Pressure (Torr) ^{vi}	pKa (Most acidic) ^{vii}	pKa (Most basic) ^{vii}	Supplier
<i>Thiamethoxam</i>	153719-23-4	C ₈ H ₁₀ ClN ₅ O ₃ S	291.71	2.86x10 ³	0.8	1.365	1.29	6.87x10 ⁻¹⁵	1.36x10 ⁻⁹	0.99±0.10	-	LGC (Ultra)
<i>Imidacloprid</i>	138261-41-3	C ₉ H ₁₀ ClN ₅ O ₂	255.66	7.17x10 ³	0.56	1.527	-3.75	1.04x10 ⁻¹³	5.07x10 ⁻⁸	7.16±0.20	4.80±0.20	Sigma Aldrich (PESTANAL)
<i>Clothianidin</i>	210880-92-5	C ₈ H ₈ ClN ₅ O ₂ S	249.68	6.00x10 ³	0.64	1.444	0.41	9.21x10 ⁻¹⁶	5.16x10 ⁻⁷	2.76±0.50	-0.2±0.10	Sigma Aldrich (PESTANAL)
<i>Metazachlor</i>	67129-08-2	C ₁₄ H ₁₆ ClN ₃ O	277.75	250	2.38	2.455	2.98	5.80x10 ⁻¹¹	6.50x10 ⁻⁸	-	1.54±0.10	LGC (Dr Ehrenstorfer)
<i>Terbutylazine</i>	5915-41-3	C ₉ H ₁₆ ClN ₅	229.71	20.4	3.27	2.49	2.48	5.94x10 ⁻⁹	9.19x10 ⁻⁶	-	2.69±0.10	Sigma Aldrich
<i>Methiocarb</i>	2032-65-7	C ₁₁ H ₁₅ NO ₂ S	225.31	104	2.87	2.444	3.13	1.14x10 ⁻⁹	5.92x10 ⁻⁴	12.16±0.46	-1.52±0.70	Sigma Aldrich (PESTANAL)
<i>Dichlofluanid</i>	1085-98-9	C ₉ H ₁₀ Cl ₂ FN ₂ O ₂ S	333.23	5.41	2.72	3.068	3.22	6.74x10 ⁻⁷	1.10x10 ⁻⁴	-5.37±0.50	-	Sigma Aldrich (PESTANAL)
<i>Flufenacet</i>	142459-58-3	C ₁₄ H ₁₅ F ₂ N ₃ O ₂ S	363.33	9.52	2.39	3.102	3.22	2.00x10 ⁻¹¹	1.18x10 ⁻⁶	-	0.31±0.50	LGC (Dr Ehrenstorfer)
<i>Oxadiazon</i>	19666-30-9	C ₁₅ H ₁₈ Cl ₂ N ₂ O	345.22	0.53	4.81	3.532	5.31	3.22x10 ⁻⁸	3.65x10 ⁻⁷	-2.73±0.40	-	Sigma Aldrich (PESTANAL)
<i>Chlorpyrifos</i>	2921-88-2	C ₉ H ₁₁ Cl ₃ NO ₃ PS	350.59	11.0	5.11	3.898	4.78	2.52x10 ⁻⁶	1.63x10 ⁻⁵	-	-5.28±0.10	LGC (Dr Ehrenstorfer)
<i>Triallate</i>	2303-17-5	C ₁₀ H ₁₆ Cl ₃ NOS	304.66	1.36	4.57	3.343	3.8	5.66x10 ⁻⁶	2.86x10 ⁻⁴	-1.48±0.70	-	Sigma Aldrich (PESTANAL)
<i>Tylosin</i>	1401-69-0	C ₄₆ H ₇₇ NO ₁₇	916.12	0.50	1.05	0.205	2.14	5.77x10 ⁻³⁸	0.00	13.06±0.70	7.39±0.70	LGC
<i>Sulfapyridine</i>	144-83-2	C ₁₁ H ₁₁ N ₃ O ₂ S	249.29	1.20x10 ⁴	0.53	1.413	0.21	1.08x10 ⁻¹³	3.90x10 ⁻⁹	8.54±0.30	2.13±0.10	Sigma Aldrich
<i>Sarafloxacin</i>	98105-99-8	C ₂₀ H ₁₇ F ₂ N ₃ O ₃	385.36	1.14x10 ³	1.07	0.791	0.45	1.92x10 ⁻¹⁹	2.67x10 ⁻¹⁶	6.17±0.41	8.68±0.10	Sigma Aldrich (VETRANAL)
<i>Ceftiofur</i>	80370-57-6	C ₁₉ H ₁₇ N ₅ O 7S ₃	523.56	22.8	1.57	0.7	-2.49	2.18x10 ⁻²⁵	-	2.62±0.50	1.70±0.10	Sigma Aldrich
<i>Diazinon</i>	333-41-5	C ₁₂ H ₂₁ N ₃ O ₃ P S	304.35	6.46	3.86	3.339	4.19	8.73x10 ⁻⁸	7.07x10 ⁻⁵	1.21±0.30	-	Sigma Aldrich (PESTANAL)

[i] As calculated by EPI Suite [1] at 25°C

[ii] As calculated by EPI Suite [1] (KOWWIN v1.68 estimate)

[iii] As calculated by EPI Suite [1] based on Log Kow

[iv] As calculated by Marvin Beans [2] at pH 7.5

[v] As calculated by EPI Suite [1] based on Bond SAR method

[vi] As stated on Scifinder calculated using Advanced Chemistry Development (ACD/Labs) Software v11.02 (© 1994-2015 ACD/Labs)[3]

[vii] As stated on Scifinder calculated using Advanced Chemistry Development (ACD/Labs) Software v11.02 (© 1994-2015 ACD/Labs) at 25°C [3]

Table S2 MS/MS detection parameters for analytes of interest sorted by class, includes analytes from ESI^{-ve} and ESI^{+ve} methods

Class of Analyte	Analyte	Method	Precursor (m/z)	Product		CV (V)	CE (eV)	Product Ion 2 (C)	CV (V)	CE (eV)	Internal Standard
				Ion 1 (Q)	Ion 2 (C)						
UV filters	Benzophenone-1	C-ESI-	213.0	134.8	36	20	90.8	34	34	25	Ibuprofen-d3
	Benzophenone-2	C-ESI-	244.9	134.9	32	13	108.9	32	32	20	Bisphenol A-d16
	Benzophenone-3	B-ESI+	229.0	151.0	35	18	105.0	35	35	20	Methadone-d9
	Benzophenone-4	C-ESI-	307.0	227.1	44	24	211.1	42	42	35	Methylparaben-13C
Parabens	Methylparaben	C-ESI-	150.8	91.8	34	20	135.8	20	20	14	Methylparaben-13C
	Ethylparaben	C-ESI-	164.9	91.9	20	20	136.6	26	26	14	Naproxen-d3
	Propylparaben	C-ESI-	179.0	91.8	34	25	136.0	20	20	16	Bisphenol A-d16
	Butylparaben	C-ESI-	193.1	91.8	34	25	136.0	40	40	16	E1 (2,4,16,16)-d4
Plasticizer	Bisphenol-A	C-ESI-	227.3	212.1	40	22	132.7	40	40	25	Bisphenol A-d16
	E1	C-ESI-	269.0	145.0	55	40	158.0	55	55	40	E1 (2,4,16,16)-d4
	E2	C-ESI-	271.1	183.0	60	40	144.9	60	60	45	E2 (2,4,16,16)-d4
Steroid estrogens	EE2	C-ESI-	295.2	144.9	60	45	158.8	60	60	40	E2 (2,4,16,16)-d4
	Sulfasalazine	C-ESI-	369.8	197.0	45	25	239.9	45	45	25	Naproxen-d3
	Clarithromycin	B-ESI+	748.5	158.1	40	31	590.4	40	40	20	Methadone-d9
Antibacterials/antibiotics	Azithromycin	A-ESI+	749.5	116.1	60	54	83.1	60	60	60	Verapamil-d7
	Trimethoprim	B-ESI+	291.2	320.2	26	26	123.1	26	26	36	Methamphetamine-d5
	Sulfamethoxazole	B-ESI+	254.1	92.2	36	30	156.1	36	36	20	Benzoylcegonine-d8
	Triclosan	C-ESI-	288.8	34.8	18	10	36.8	18	18	10	E1 (2,4,16,16)-d4
	Amoxicillin	A-ESI+	365.8	113.9	22	23	207.8	22	22	13	Methiocarb-d3
	Ciprofloxacin	A-ESI+	331.5	288.0	43	18	313.6	43	43	21	Ciprofloxacin-d8
	Erythromycin	A-ESI+	734.3	576.0	34	21	558.3	34	34	21	Erythromycin-13C-d3
	Ofloxacin	A-ESI+	361.9	312.9	35	19	261.0	35	35	30	Ofloxacin-d3
	Metronidazole	A-ESI+	172.16	82.09	28	26	111.0	28	28	24	Metronidazole-d4
	Sulfadiazine	A-ESI+	250.9	107.9	30	26	157.8	30	30	15	Metronidazole-d4
Fungicide	Norfloxacin	A-ESI+	320.22	233.2	58	26	204.9	58	58	38	Ofloxacin-d3
	Cefalexin	A-ESI+	348.22	158.04	48	6	174.0	48	48	4	Cefalexin-d5
	Penicillin V	A-ESI+	351.16	114.15	54	40	160.1	54	54	40	Metazachlor-d6
	Penicillin G	A-ESI+	335.16	160.06	48	20	217.1	48	48	16	Cefalexin-d5
	Chloramphenicol	A-ESI+	323.0	304.8	20	10	274.8	20	20	10	Methiocarb-d3
	Danofloxacin	A-ESI+	358.1	340.0	65	20	255.2	38	38	38	Ciprofloxacin-d8
	Prulifloxacin	A-ESI+	462.2	444.2	40	30	360.1	40	40	44	Ofloxacin-d3
	Griseofulvin	A-ESI+	353.2	69.0	45	25	165.0	45	45	23	Methiocarb-d3

Class of Analyte	Analyte	Method	Precursor (m/z)	Product Ion 1 (Q)	CV (V)	CE (eV)	Product Ion 2 (C)	CV (V)	CE (eV)	Internal Standard
Hypertension	Ketoconazole	A-ESI+	532.0	82.0	60	50	112.13	60	58	Amitriptyline-d3
	Valsartan	C-ESI-	434.0	350.1	35	20	179.1	35	25	Naproxen-d3
NSAIDs	Irbesartan	C-ESI-	427.1	193.0	50	28	121.0	50	65	Bisphenol A-d16
	Lisinopril	B-ESI+	406.2	84.0	38	27	246.1	38	22	Amphetamine-d5
	Ketoprofen	C-ESI-	253.0	209.0	15	7	-	-	-	Ketoprofen-d3
	Ibuprofen	C-ESI-	204.9	161.5	26	6	-	-	-	Ibuprofen-d3
Lipid regulators	Naproxen	C-ESI-	229.0	185.0	20	7	170.0	20	13	Naproxen-d3
	Diclofenac	C-ESI-	293.8	249.9	22	13	-	-	-	Naproxen-d3
	Acetaminophen	B-ESI+	151.9	110.0	26	16	211.2	22	10	Acetaminophen-d4
	Bezafibrate	C-ESI-	360.1	274.0	30	19	154.0	30	28	Bezafibrate-d6
Anti-hyperlipidemic	Atorvastatin	C-ESI-	557.1	397.1	50	30	278.1	50	45	Naproxen-d3
	Gemfibrozil	A-ESI+	251	205.0	21	9	123.0	21	14	Methiocarb-d3
Anti-hyperintensive	Candesartan	A-ESI+	611.1	567.0	44	7	466.9	44	7	Cefalexin-d5
	Fexofenadine	C-ESI-	500.1	456.1	33	14	378.1	33	19	Naproxen-d3
Antihistamines	Cetirizine	B-ESI+	389.1	201.0	32	21	166.0	32	40	Temazepam-d5
	Sildenafil	A-ESI+	475.3	99.9	60	31	282.9	60	43	Sildenafil-d8
GUD/ED	Metformin	B-ESI+	130.0	60.0	30	15	71.0	30	20	Metformin-d6
	Gliclazide	B-ESI+	324.1	127.0	41	20	110.0	41	20	Quetiapine-d8
Cough suppressant	Sitagliptin	A-ESI+	407.8	234.9	46	19	192.9	46	26	Verapamil-d7
	Pholcodine	B-ESI+	399.2	381.2	55	25	100.1	55	37	Atenolol-d7
Beta-blocker	Atenolol	B-ESI+	267.3	145.1	38	30	190.1	38	16	Atenolol-d7
	Metoprolol	B-ESI+	268.1	116.1	42	20	82.1	40	31	Metoprolol-d7
H ₂ receptor agonists	Propranolol	B-ESI+	260.2	183.1	42	18	116.1	42	16	Propranolol-d7
	Bisoprolol	A-ESI+	326.0	116.0	45	18	204.0	45	19	Sildenafil-d8
X-ray contrast media	Ranitidine	B-ESI+	315.9	176.0	26	17	123.9	26	24	Atenolol-D7
	Cimetidine	B-ESI+	252.9	159.4	22	16	211.2	22	10	Acetaminophen-d4
Various uses	Iopromide	B-ESI+	792.0	573.0	46	25	558.9	46	32	Acetaminophen-d4
	Buprenorphine	A-ESI+	468.4	396.0	66	41	414.0	66	35	Amitriptyline-d3
Drug precursor and metabolite	Ephedrine/pseudoephedrine	B-ESI+	166.1	148.1	23	12	133.1	23	21	1S,2R-(+)-Ephedrine-d3
	Norephedrine	B-ESI+	152.2	134.1	23	10	117.1	23	16	1S,2R-(+)-Ephedrine-d3
Anti-cancer	Azathioprine	B-ESI+	278.0	142.0	28	13	85.0	28	20	Cotinine-d3
	Methotrexate	B-ESI+	455.1	308.1	40	20	175.1	40	35	Amphetamine-d5
Anti-cancer	Ifosfamide	B-ESI+	261.0	92.0	40	28	154.0	40	22	Metoprolol-d7
	Tamoxifen	B-ESI+	372.2	72.0	50	25	129.0	50	28	Tamoxifen-13C2-15N

Class of Analyte	Analyte	Method	Precursor (m/z)	Product Ion 1 (Q)	CV (V)	CE (eV)	Product Ion 2 (C)	CV (V)	CE (eV)	Internal Standard
Anaesthetic and metabolite	Capecitabine	A-ESI+	359.9	243.9	25	11	173.8	25	23	Metazachlor-d6
	Imatinib	A-ESI+	494.2	394	57	27	378.3	57	48	Nortriptyline-d3
	Bicalutamide	A-ESI+	428.8	254.8	23	15	410.7	23	15	Methiocarb-d3
Anti-depressants and metabolites	Ketamine	B-ESI+	238.1	125.0	31	27	220.1	31	15	Ketamine-d4
	Norketamine	B-ESI+	224.0	207.1	23	12	125.0	23	27	Norketamine-d4
	Venlafaxine	B-ESI+	278.2	58.1	27	40	260.1	27	12	Metoprolol-d7
	Desvenlafaxine	B-ESI+	264.0	107.1	25	24	246.3	25	20	Metoprolol-d7
	Fluoxetine	B-ESI+	310.2	44.1	34	10	148.1	34	10	Fluoxetine-d5
Anti-epileptic and metabolites	Norfluoxetine	B-ESI+	296.1	134.1	18	6	-	-	-	Fluoxetine-d5
	Sertraline	B-ESI+	306.0	159.0	23	27	275.0	23	10	Sertraline-d3
	Mirtazapine	B-ESI+	266.1	195.0	44	26	72.0	44	18	Mirtazapine-d3
	Citalopram	B-ESI+	325.1	262.1	46	18	109.9	46	26	Citalopram-d6
	Desmethylcitalopram	B-ESI+	311.4	109.0	46	27	262.0	46	18	Citalopram-d6
	Amitriptyline	A-ESI+	278.2	91.1	37	26	233.2	37	18	Amitriptyline-d3
	Nortriptyline	A-ESI+	264.2	233.1	33	16	91.0	33	23	Nortriptyline-d3
	Duloxetine	A-ESI+	297.9	153.9	16	5	187.9	16	5	Nortriptyline-d3
	Paroxetine	A-ESI+	330.0	192.0	45	21	150.9	45	22	Nortriptyline-d3
	Norsertaline	A-ESI+	291.8	158.6	25	23	128.8	25	22	Norsertaline-d4
Anti-epileptic and metabolites	Carbamazepine	B-ESI+	327.0	194.1	40	20	179.1	40	38	Carbamazepine-13C6
	Carbamazepine 10,11-epoxide	B-ESI+	253.1	180.1	39	25	210.1	39	12	Carbamazepine-13C6
	10,11-Dihydro-10-hydroxycarbamazepine	B-ESI+	255.1	194.1	20	20	179.1	20	40	Carbamazepine-13C6
Calcium channel blocker	Diltiazem	B-ESI+	415.0	178.0	40	25	310.1	40	25	Carbamazepine-13C6
	Verapamil	A-ESI+	455.2	165.0	55	31	303.1	55	26	Verapamil-d7
Hypnotic	Temazepam	B-ESI+	301.1	255.1	37	21	283.1	37	14	Temazepam-d5
	Diazepam	A-ESI+	285.0	153.9	56	29	221.8	56	27	Diazepam-d5
Anti-psychoctic	Oxazepam	A-ESI+	286.9	240.8	38	22	268.8	38	14	Oxazepam-d5
	Quetiapine	B-ESI+	384.1	253.1	50	21	221.1	50	40	Quetiapine-d8
Dementia	Risperidone	A-ESI+	411.0	190.9	49	30	109.9	49	51	Verapamil-d7
	Donepezil	A-ESI+	380.3	362.1	56	23	288.0	56	24	Verapamil-d7
	Memantine	A-ESI+	180.0	107.0	36	24	121.0	36	24	Verapamil-d7
Human indicators and metabolites	Creatinine	B-ESI+	114.0	44.0	30	15	86.1	31	11	Metformin-d6
	Nicotine	B-ESI+	163.1	130.0	37	20	117.0	37	20	Atenolol-d7
	Caffeine	B-ESI+	195.1	138.0	38	15	110.0	38	23	Cotinine-d3

Class of Analyte	Analyte	Method	Precursor (m/z)	Product Ion 1 (Q)	CV (V)	CE (eV)	Product Ion 2 (C)	CV (V)	CE (eV)	Internal Standard
Analgesics and metabolites	Cotinine	B-ESI+	177.1	80.0	34	21	98.1	34	22	Cotinine-d3
	1,7 dimethylxanthine	B-ESI+	181.0	124.1	54	21	-	-	-	Cotinine-d3
	Morphine	B-ESI+	286.2	165.1	53	38	152.1	53	56	Morphine-d3
	Dihydromorphine	B-ESI+	288.2	185.0	28	42	213.0	28	32	Morphine-d3
	Normorphine	B-ESI+	272.1	165.0	45	43	152.1	45	49	Morphine-d3
	Methadone	B-ESI+	310.2	265.1	31	15	105.1	31	28	Methadone-d9
	EDDP	B-ESI+	278.2	234.1	50	29	249.1	50	24	EDDP-d3
	Codeine	B-ESI+	300.2	215.1	49	25	152.1	49	57	Codeine-d6
	Norcodeine	B-ESI+	286.1	165.1	46	40	268.2	46	20	Codeine-d6
	Dihydrocodeine	B-ESI+	302.1	199.1	53	33	128.1	53	60	Codeine-d6
	Tramadol	B-ESI+	264.0	58.0	28	45	120.7	28	46	Metoprolol-d7
	N-desmethyltramadol	B-ESI+	250.1	44.0	25	12	232.1	25	8	Cocaine-d3
	O-desmethyltramadol	B-ESI+	250.2	58.0	30	18	232.1	30	10	Cotinine-d3
	Amphetamine	B-ESI+	136.2	91.1	18	16	119.1	18	8	Amphetamine-d5
Stimulants and metabolites	Methamphetamine	B-ESI+	150.2	91.1	18	16	119.1	24	10	Methamphetamine-d5
	MDMA	B-ESI+	194.1	163.1	24	13	105.1	24	24	MDMA-d5
	MDA	B-ESI+	180.0	163.1	21	11	105.1	21	22	MDA-d5
	Cocaine	B-ESI+	304.2	182.1	40	20	82.1	40	31	Cocaine-d3
	Benzoylcegonine	B-ESI+	290.2	168.1	38	19	105.1	38	30	Benzoylcegonine-d8
	Anhydrocegoninemethylester	B-ESI+	182.1	118.0	39	23	122.1	37	20	Atenolol-d7
	Cocaehtylene	B-ESI+	318.2	196.2	38	20	82.1	38	30	Cocaehtylene-d3
	Mephedrone	B-ESI+	178.1	160.1	10	12	145.0	10	22	Mephedrone-d3
	MDPV	B-ESI+	276.1	126.1	40	28	135.0	40	25	Cocaehtylene-d3
	Heroin	B-ESI+	370.2	165.1	51	50	268.1	51	29	Heroin-d9
	6-acetylmorphine	B-ESI+	328.1	165.1	52	39	211.1	52	26	Cotinine-d3
	Dichlofluanid	A-ESI+	334.7	122.9	29	31	270.8	29	6	Amitriptyline-d3
	Oxadiazon	A-ESI+	348.1	331.0	18	7	313.1	18	12	Methiocarb-d3
	Triallate	A-ESI+	305.8	127.8	34	13	144.7	34	26	Methiocarb-d3
Methiocarb	A-ESI+	225.9	208.8	34	15	174.9	34	19	Methiocarb-d3	
Clothianidin	A-ESI+	249.8	131.8	28	15	112.9	28	25	Methiocarb-d3	
Imidacloprid	A-ESI+	255.9	208.8	34	15	174.9	34	19	Metazachlor-d6	
Thiamethoxam	A-ESI+	291.8	210.9	44	12	131.9	44	22	Methiocarb-d3	
Terbutylazine	A-ESI+	229.9	173.9	35	17	131.9	35	24	Diazepam-d5	
Metazachlor	A-ESI+	277.9	209.8	21	21	133.9	21	10	Metazachlor-d6	

Class of Analyte	Analyte	Method	Precursor (m/z)	Product Ion 1 (Q)	CV (V)	CE (eV)	Product Ion 2 (C)	CV (V)	CE (eV)	Internal Standard
Veterinary Pharmaceuticals	Chlorpyrifos	A-ESI+	349.8	124.8	34	19	197.8	34	16	Metazachlor-d6
	Flufenacet	A-ESI+	363.8	152.0	24	18	194.0	24	10	Metazachlor-d6
	Diazinon	A-ESI+	304.9	168.9	36	22	152.9	36	22	Oxazepam-d5
	Tylosin	B-ESI+	916.5	174.2	80	45	101.0	80	56	Methadone-d9
	Oxytetracycline	A-ESI+	460.9	443.5	32	20	425.8	32	13	Ofloxacin-d3
	Sarafloxacin	A-ESI+	385.9	367.8	49	23	298.8	49	28	Verapamil-d7
	Sulfapyridine	A-ESI+	250.03	156.04	40	16	92.0	42	30	Metronidazole-d4
	Tetracycline	A-ESI+	444.9	409.9	33	20	427.5	33	14	Ciprofloxacin-d8
	Ceftiofur	A-ESI+	524.0	241.0	50	50	241.0	50	20	Ofloxacin-d3

Table S3 MS/MS detection parameters for internal standards sorted by analogue analyte class, includes internal standards from ESI –ve and ESI +ve methods

ISTD	Method	Precursor (m/z)	Product Ion 1 (Q)	CV (V)	CE (eV)	ISTD	Method	Precursor (m/z)	Product Ion 1 (Q)	CV (V)	CE (eV)
Metronidazole-d4	A-ESI+	175.9	127.9	26	15	Norketamine-d4	B-ESI+	228.1	128.9	32	28
Gabapentin-d4	A-ESI+	176	158.1	33	16	Metoprolol-d7	B-ESI+	275.4	123.1	44	20
Methiocarb-d3	A-ESI+	228.9	168.9	25	9	Cocaine-d3	B-ESI+	307.2	185.1	40	20
Nortriptyline-d3	A-ESI+	267.1	191	40	20	Cocaethylene-d3	B-ESI+	321.2	199.1	40	22
Amitriptyline-d3	A-ESI+	281.5	90.9	40	25	Mirtazapine-d3	B-ESI+	269.0	194.9	35	25
Metazachlor-d6	A-ESI+	284	215.9	21	10	EDDP-d3	B-ESI+	281.2	234.1	50	29
Diazepam-d5	A-ESI+	290.1	153.9	56	28	Propranolol-d7	B-ESI+	267.0	188.8	40	18
Oxazepam-d5	A-ESI+	292.0	245.9	38	22	Citalopram-d6	B-ESI+	331.0	109.0	46	28
Norsertaline-d4	A-ESI+	296	278.8	39	8	Carbamazepine-13C6	B-ESI+	243.1	200.1	40	20
Ciprofloxacin-d8	A-ESI+	340.1	321.9	55	21	Quetiapine-d8	B-ESI+	392.1	258.1	50	23
Cefalexin-d5	A-ESI+	353.2	194.5	43	10	Methadone-d9	B-ESI+	319.3	268.2	31	15
Ofloxacin-d3	A-ESI+	364.8	260.8	47	28	Temazepam-d5	B-ESI+	306.7	260.1	37	21
Verapamil-d7	A-ESI+	462.1	165	56	31	Fluoxetine-d5	B-ESI+	315.3	153.2	26	8
Sildenafil-d8	A-ESI+	482.9	108.0	60	31	Sertraline-d3	B-ESI+	309.1	159.0	23	27
Erythromycin-13C,D3	A-ESI+	738.5	162	46	33	Tamoxifen-13C2-15N	B-ESI+	375.1	75.0	50	25
Metformin-d6	B-ESI+	136.1	77.0	30	19	Methylparaben-13C	C-ESI-	156.9	97.9	30	20
Morphine-d3	B-ESI+	289.1	152.1	53	56	Ketoprofen-d3	C-ESI-	265.0	212.0	15	7
Atenolol-d5	B-ESI+	274.3	145.1	44	30	Bezafibrate-d6	C-ESI-	366.0	274.0	30	19
Acetaminophen-d4	B-ESI+	156.0	114.0	26	16	Naproxen-d13	C-ESI-	232.0	188.0	15	8
Codeine-d6	B-ESI+	306.2	218.1	52	28	Bisphenol A-d16	C-ESI-	241.1	223.1	40	20
Cotinine-d3	B-ESI+	180.1	80.0	44	24	Ibuprofen-d3	C-ESI-	208.0	164.0	20	6
1S,2R-(+) Ephedrine-d3	B-ESI+	169.2	151.0	23	18	E2 (2, 4, 16, 16)-d4	C-ESI-	275.1	147.0	60	40
Amphetamine-d5	B-ESI+	141.1	92.8	20	14	E1 (2, 4, 16, 16)-d4	C-ESI-	273.1	147.0	55	40
Methamphetamine-d5	B-ESI+	155.1	91.8	28	18						
MDA-d5	B-ESI+	185.1	168.1	21	11						
MDMA-d5	B-ESI+	199.1	165.1	26	13						
Benzoylcegonine-d8	B-ESI+	298.2	171.1	38	19						
Mephedrone-d3	B-ESI+	181.1	148.0	30	22						
Ketamine-d4	B-ESI+	242.1	129.1	31	27						
Heroin-d9	B-ESI+	379.2	165.8	51	50						

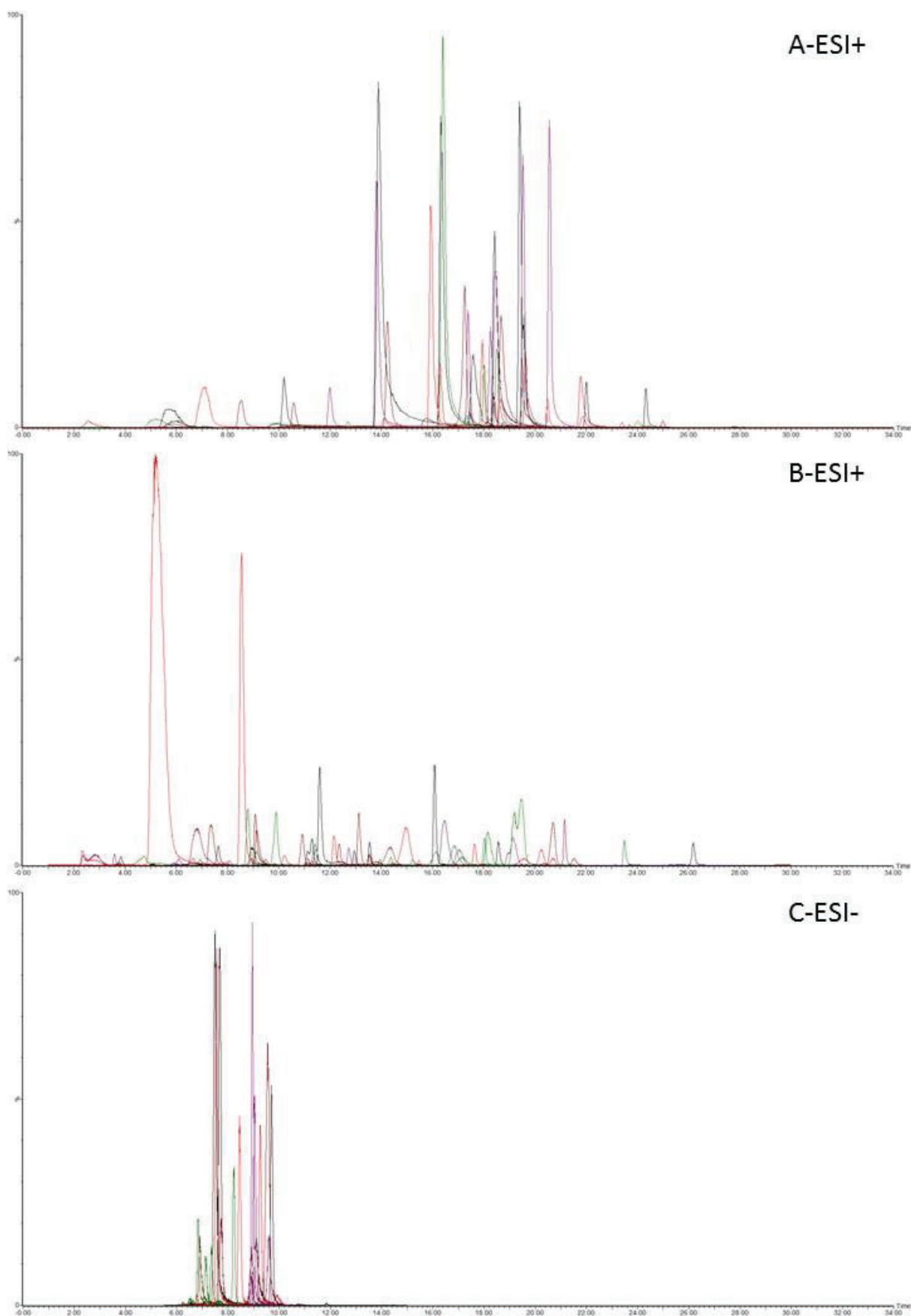


Fig. S1 Chromatograms of all analytes in spiked influent sample from Site A using methods: A-ESI+, B-ESI+ and C-ESI-. The run time for C-ESI- is 22.5 minutes, which is shorter than the 34 minute run time of A-ESI+ and B-ESI+, hence the much more condensed appearance of the peaks.

Table S4 Solid sample weights

Digested solids														
Sample No.	1	2	3	4	5	6								
Day	1	1	2	2	3	3								
Weight (g)	0.24725	0.24989	0.24604	0.24189	0.24855	0.24137								
Solid Particulate matter														
Sample No.	1	2	3	4	5	6	7	8	9	10	11	12	13	14
Day	1	1	2	2	3	3	4	4	5	5	6	6	7	7
Weight (g)	0.23660	0.22825	0.23961	0.23985	0.23428	0.23437	0.23800	0.24990	0.24075	0.23251	0.23590	0.23546	0.24345	0.23930

Table S5 Suspended solids weights (g L⁻¹)

Wed	Thu	Fri	Sat	Sun	Mon	Tue
Day 1	Day 2	Day 3	Day 4	Day 5	Day 6	Day 7
0.2094	0.2072	0.2376	0.2006	0.1634	0.1134	0.147

Table S6 Adsorbance

Compound	Spiked Surface water adsorbance to PPE bottle after 24h storage		Recovery 24h storage on SPE cartridge	
	Ads (%)	SD	Rec (%)	SD
Benzophenone-1	-	-	-	-
Benzophenone-2	-	-	-	-
Benzophenone-3	33.7	18.7	45.3	16.0
Benzophenone-4	-	-	-	-
Methylparaben	-	-	-	-
Ethylparaben	-	-	-	-
Propylparaben	-	-	-	-
Butylparaben	-	-	-	-
Bisphenol A	-	-	-	-
E1	-	-	-	-
E2	-	-	-	-
EE2	-	-	-	-
Sulfasalazine	-	-	-	-
Clarithromycin	-6.4	2.1	106.8	11.8
Trimethoprim	5.2	3.0	94.7	4.3
Sulfamethoxazole	9.3	0.1	98.8	3.9
Triclosan	-	-	-	-
Amoxicillin	-	-	-	-
Metronidazole	-6.8	4.8	107.4	2.1
Sulfadiazine	37.6	3.6	97.9	10.9
Cefalexin	-	-	-	-
Ofloxacin	18.1	16.4	96.8	6.4
Ciprofloxacin	-	-	-	-
Tetracycline	-	-	-	-
Danofloxacin	-	-	-	-
Oxytetracycline	-	-	-	-
Chloramphenicol	-9.7	16.1	126.5	10.7
Penicillin G	-	-	-	-
Penicillin V	-60.9	55.6	149.7	41.9
Erythromycin	-4.0	11.2	146.2	71.1
Azithromycin	16.6	2.3	115.0	21.6
Prulifloxacin	-	-	-	-
Norfloxacin	-	-	-	-
Griseofulvin	-0.3	8.4	98.2	3.2
Ketoconazole	-10.6	2.9	104.6	13.6
Valsartan	-	-	-	-
Irbesartan	-	-	-	-
Lisinopril	73.7	20.8	22.0	1.6
Ketoprofen	-	-	-	-
Ibuprofen	-	-	-	-
Naproxen	-	-	-	-
Diclofenac	-	-	-	-

Compound	Spiked Surface water adsorbance to PPE bottle after 24h storage		Recovery 24h storage on SPE cartridge	
	Ads (%)	SD	Rec (%)	SD
Acetaminophen	1.4	3.1	97.3	2.7
Bezafibrate	-	-	-	-
Atorvastatin	-	-	-	-
Gemfibrozil	-15.5	11.2	107.0	4.2
Candesartan Cilexetil	-141.2	53.2	322.5	208.8
Fexofenadine	-	-	-	-
Cetirizine	1.6	0.7	89.8	6.1
Sildenafil	-9.8	7.6	107.1	3.9
Metformin	-62.3	5.1	71.5	10.4
Gliclazide	6.3	6.6	107.3	3.3
Sitagliptin	15.6	5.0	149.4	34.5
Pholcodine	43.4	1.8	65.2	4.4
Atenolol	-6.8	3.8	109.6	1.5
Metoprolol	7.8	1.0	98.8	1.3
Propranolol	-12.4	7.3	108.2	16.1
Bisoprolol	6.0	10.1	118.5	22.0
Ranitidine	30.7	5.7	118.8	5.8
Cimetidine	52.5	0.7	99.2	2.4
Iopromide	-11.2	2.5	108.1	3.1
Buprenorphine	0.5	1.5	151.6	26.9
Ephedrine/pseudoephedrine	52.9	0.8	93.5	5.1
Norephedrine	18.0	1.8	129.5	23.5
Azathioprine	5.6	3.3	95.6	1.3
Methotrexate	30.2	0.8	57.0	7.2
Ifosfamide	5.5	0.1	99.3	0.3
Tamoxifen	60.2	5.9	73.0	26.7
Imatinib	11.7	2.2	101.9	17.7
Capecitabine	-1.1	2.2	94.2	4.7
Bicalutamide	-19.7	14.4	95.3	14.0
Ketamine	-2.7	0.7	110.2	2.5
Norketamine	5.4	1.3	98.4	1.0
Venlafaxine	3.8	1.3	101.9	2.6
Desvenlafaxine	-1.3	1.8	103.1	3.7
Fluoxetine	-22.2	21.6	121.4	32.9
Norfluoxetine	58.4	9.2	71.3	36.6
Sertraline	1.4	20.6	96.0	25.8
Mirtazapine	6.0	1.6	100.4	6.6
Citalopram	2.5	4.6	92.4	14.5
Desmethylcitalopram	-2.0	7.0	110.1	17.9
Paroxetine	35.8	15.7	85.4	16.4
Duloxetine	-	-	-	-
Amitriptyline	14.0	3.1	95.8	6.2
Nortriptyline	30.0	0.8	110.0	9.8
Norsertaline	-	-	-	-

Compound	Spiked Surface water adsorbance to PPE bottle after 24h storage		Recovery 24h storage on SPE cartridge	
	Ads (%)	SD	Rec (%)	SD
Carbamazepine	6.3	2.1	96.2	0.7
Carbamazepine 10,11 epoxide	-7.8	0.9	109.5	0.8
10,11 dihydro 10 hydroxycarbamazepine	17.4	0.7	82.4	0.1
Diltiazem	-2.0	3.9	91.6	17.3
Verapamil	9.1	3.0	94.5	4.9
Temazepam	-0.5	2.2	100.8	0.7
Oxazepam	14.9	0.4	89.6	0.5
Diazepam	-17.7	0.2	119.4	0.6
Quetiapine	-7.3	5.3	97.0	14.5
Risperidone	14.8	8.3	118.5	19.6
Donepezil	16.4	13.7	118.9	29.4
Memantine	3.2	10.4	174.0	41.0
Creatinine	75.4	0.7	37.8	7.3
Nicotine	-11.8	6.0	147.9	5.6
Caffeine	4.6	0.5	96.8	1.6
Cotinine	0.1	3.0	97.3	1.0
1,7 dimethylxantine	-4.4	3.7	99.1	6.2
Morphine	3.3	10.1	104.7	2.7
Dihydromorphine	24.1	2.6	90.2	0.5
Normorphine	17.0	4.5	107.4	3.3
Methadone	1.5	0.7	101.2	11.4
EDDP	53.7	1.3	98.3	11.8
Codeine	1.2	5.1	107.0	4.9
Norcodeine	6.6	7.0	106.7	5.3
Dihydrocodeine	14.7	0.2	95.1	0.7
Tramadol	1.9	2.3	101.8	0.1
N-desmethyltramadol	16.7	4.3	99.0	3.5
O-desmethyltramadol	4.0	0.4	99.2	2.0
Amphetamine	43.3	1.0	87.1	0.1
Methamphetamine	3.2	0.7	118.2	1.8
MDMA	21.0	2.5	90.2	2.8
MDA	44.4	2.5	82.1	4.1
Cocaine	1.4	0.1	99.7	2.0
Benzoylcegonine	1.0	0.7	95.9	0.2
Anhydrocegonine methylester	10.3	3.7	108.7	3.5
Cocaethylene	6.0	2.0	94.8	2.4
Mephedrone	33.6	4.3	90.1	7.3
MDPV	38.9	3.4	82.8	7.0
Heroin	17.8	0.8	87.7	4.4
O-6-MAM	15.6	2.8	93.8	2.9
Thiamethoxam	-1.7	6.0	106.1	2.2
Imidacloprid	1.1	1.0	95.2	2.3
Clothianidin	-1.8	9.7	107.5	5.9

Compound	Spiked Surface water adsorbance to PPE bottle after 24h storage		Recovery 24h storage on SPE cartridge	
	Ads (%)	SD	Rec (%)	SD
Metazachlor	2.1	6.7	96.9	14.2
Terbutylazine	-13.9	3.7	119.0	0.8
Methiocarb	4.1	5.4	99.3	5.3
Dichlofluanid	-	-	-	-
Flufenacet	-9.7	0.6	103.1	1.3
Oxadiazon	-27.6	28.9	74.1	15.8
Chlorpyrifos	-3.3	0.3	112.0	5.8
Triallate	-44.7	33.0	136.8	9.9
Tylosin	11.4	7.5	86.0	1.4
Sulfapyridine	-2.9	3.9	103.7	1.0
Sarafloxacin	-	-	-	-
Ceftiofur	73.2	2.7	29.9	13.8
Diazinon	-8.9	10.9	106.0	16.8

Table S7 Instrument parameters (LC MSMS): Intra- and Inter-day precision and accuracy

Compound	Inter-day precision (%)			Intra-day precision (%)			Inter-day Accuracy (%)			Intra-day Accuracy (%)		
	Low	Med	High	Low	Med	High	Low	Med	High	Low	Med	High
<i>Benzophenone-1</i>	2.4	6.4	1.0	4.9	1.1	0.9	106.3	116.2	97.6	103.6	119.0	97.9
<i>Benzophenone-2</i>	4.1	6.1	2.6	1.3	1.2	0.8	90.8	106.4	95.6	90.4	110.0	98.2
<i>Benzophenone-3</i>	2.5	3.8	-	4.3	4.6	-	83.4	86.3	-	84.4	89.2	-
<i>Benzophenone-4</i>	7.0	3.3	1.1	2.5	2.7	1.8	117.0	92.3	105.8	115.0	89.0	105.0
<i>Methylparaben</i>	5.2	11.5	1.4	0.7	0.4	2.1	94.4	101.0	96.9	90.7	94.0	95.3
<i>Ethylparaben</i>	2.2	4.1	0.1	4.6	1.2	2.2	120.2	120.9	98.0	117.3	121.6	98.0
<i>Propylparaben</i>	2.8	9.3	-	2.3	12.2	-	99.8	102.9	-	97.1	100.4	-
<i>Butylparaben</i>	3.7	6.2	1.1	5.1	9.2	0.8	97.2	103.9	99.7	93.6	99.2	98.4
<i>Bisphenol A</i>	0.9	2.6	0.3	2.9	2.6	1.8	107.0	108.7	98.0	107.0	106.1	97.7
<i>E1</i>	3.5	1.8	1.0	2.0	1.6	1.7	99.1	97.7	98.9	96.2	96.7	97.9
<i>E2</i>	2.3	5.2	0.4	3.5	4.5	1.1	95.7	93.7	99.5	93.7	96.9	99.2
<i>EE2</i>	7.0	1.3	1.5	2.1	3.9	1.8	84.3	94.3	101.1	86.1	95.2	102.5
<i>Sulfasalazine</i>	3.4	1.4	2.4	4.0	4.3	3.3	112.3	100.4	101.3	115.4	101.3	98.8
<i>Clarithromycin</i>	5.1	1.5	1.3	4.0	2.0	1.2	97.5	99.9	102.0	102.6	99.3	103.6
<i>Azithromycin</i>	-	1.6	1.4	-	3.9	5.2	-	108.1	95.9	-	113.6	104.2
<i>Trimethoprim</i>	1.4	3.5	4.1	1.6	1.8	3.2	104.5	96.8	89.3	104.6	94.2	99.6
<i>Sulfamethoxazole</i>	2.0	3.0	5.6	0.7	2.2	4.3	104.3	91.6	89.3	106.1	93.6	88.3
<i>Triclosan</i>	2.9	13.7	2.8	10.0	16.1	2.0	69.8	72.4	72.1	69.4	66.3	73.2
<i>Amoxicillin</i>	8.1	6.9	5.2	8.8	3.2	3.9	93.5	88.8	100.8	104.2	99.1	113.8
<i>Metronidazole</i>	2.3	0.5	0.7	1.3	2.9	3.3	114.4	99.5	95.0	110.9	104.6	99.5
<i>Sulfadiazine</i>	1.1	1.8	1.6	1.9	3.3	3.0	126.7	90.7	95.9	123.8	93.0	99.2
<i>Cefalexin</i>	-	12.3	-	-	9.5	-	-	102.9	-	-	111.3	-
<i>Ofloxacin</i>	2.8	3.8	1.9	6.8	3.4	2.3	98.0	95.3	94.3	97.6	98.3	96.3
<i>Ciprofloxacin</i>	7.0	6.0	3.3	16.8	6.9	2.4	75.5	97.5	97.7	58.1	104.4	104.4
<i>Tetracycline</i>	13.1	5.0	7.4	12.9	2.9	4.6	151.2	91.3	96.9	134.6	99.6	111.0
<i>Danofloxacin</i>	-	11.1	0.9	-	7.7	6.8	-	95.9	102.5	-	95.0	116.9
<i>Oxytetracycline</i>	-	5.0	1.1	-	3.9	5.2	-	85.7	92.0	-	87.4	99.7
<i>Chloramphenicol</i>	4.4	1.7	-	6.5	1.1	-	102.4	99.2	-	104.8	102.2	-
<i>Penicillin G</i>	6.2	2.7	-	16.1	4.5	-	125.3	98.2	-	128.8	102.2	-
<i>Penicillin V</i>	-	15.0	-	-	4.4	-	-	96.8	-	-	88.5	-
<i>Erythromycin</i>	-	-	2.9	-	-	2.3	-	-	95.2	-	-	94.4
<i>Prulifloxacin</i>	-	3.0	14.7	-	2.2	6.5	-	100.7	72.2	-	106.3	91.0
<i>Norfloxacin</i>	-	1.5	7.3	-	4.1	4.0	-	82.1	88.1	-	84.4	86.6
<i>Griseofulvin</i>	3.4	2.6	-	2.2	1.0	-	88.8	94.4	-	86.1	92.4	-
<i>Ketoconazole</i>	3.5	1.5	2.5	2.4	4.0	4.9	87.2	92.6	95.3	89.7	94.4	100.4

Compound	Inter-day precision (%)			Intra-day precision (%)			Inter-day Accuracy (%)			Intra-day Accuracy (%)		
	Low	Med	High	Low	Med	High	Low	Med	High	Low	Med	High
<i>Valsartan</i>	2.7	3.2	4.8	1.6	0.8	3.4	157.0	98.7	100.0	157.4	95.6	94.5
<i>Irbesartan</i>	6.1	3.8	2.2	4.2	2.9	0.8	85.1	108.6	101.3	83.3	108.7	98.8
<i>Lisinopril</i>	0.9	3.4	-	13.0	1.4	-	103.6	90.7	-	103.0	87.3	-
<i>Ketoprofen</i>	3.9	3.5	0.5	3.6	0.5	2.4	101.4	96.2	100.5	99.6	99.1	101.1
<i>Ibuprofen</i>	3.7	1.3	1.8	3.8	0.9	2.6	83.1	98.4	101.0	81.8	97.8	101.7
<i>Naproxen</i>	3.0	3.6	0.9	1.7	1.0	1.8	91.9	101.8	101.4	90.2	102.5	100.3
<i>Diclofenac</i>	5.5	5.0	2.9	8.6	9.9	5.3	80.4	97.9	97.2	78.0	95.5	95.3
<i>Acetaminophen</i>	1.6	1.6	1.6	4.3	0.9	2.7	96.4	98.0	97.7	98.1	99.3	99.5
<i>Bezafibrate</i>	4.3	1.4	2.7	4.2	2.0	0.7	92.6	100.3	100.8	91.4	99.5	102.3
<i>Atorvastatin</i>	1.6	7.9	1.1	4.7	2.5	0.6	100.9	98.2	103.8	101.9	89.6	102.6
<i>Gemfibrozil</i>	8.0	5.9	-	11.7	3.9	-	144.1	98.1	-	140.5	96.6	-
<i>Candesartan Cilexetil</i>	-	-	0.9	-	-	5.2	-	-	106.9	-	-	100.5
<i>Fexofenadine</i>	4.4	13.6	1.3	1.6	3.4	1.3	101.5	103.9	108.4	98.4	111.4	109.0
<i>Cetirizine</i>	2.0	0.9	0.8	0.8	0.4	2.6	101.8	99.2	100.4	103.8	99.2	99.4
<i>Sildenafil</i>	1.7	4.8	2.6	3.8	2.2	4.6	99.3	106.0	91.9	95.8	110.6	92.2
<i>Metformin</i>	2.6	1.1	0.8	1.7	1.2	0.9	86.2	101.6	101.2	88.7	101.5	100.7
<i>Glitclazide</i>	0.8	4.4	1.1	2.5	5.0	0.8	86.7	92.5	100.3	87.3	97.2	101.5
<i>Sitagliptin</i>	1.8	5.7	1.5	0.9	4.1	4.6	141.1	88.7	101.2	141.5	92.8	100.8
<i>Pholcodine</i>	0.5	6.1	7.3	0.4	6.4	3.0	78.7	105.7	114.1	80.0	100.4	117.3
<i>Atenolol</i>	3.9	1.5	0.9	4.5	2.0	0.5	83.3	99.9	102.8	86.9	100.6	102.8
<i>Metoprolol</i>	2.0	1.2	0.6	0.9	3.2	1.8	91.3	98.1	101.1	90.9	96.7	100.7
<i>Propranolol</i>	2.1	1.5	2.5	0.4	0.7	1.9	122.5	96.8	97.0	125.4	98.4	94.8
<i>Bisoprolol</i>	3.3	0.7	2.1	8.0	2.9	3.4	86.6	109.8	91.5	90.4	115.0	95.7
<i>Ranitidine</i>	3.9	0.5	3.1	3.5	10.9	14.6	100.2	95.0	105.0	95.7	95.4	101.3
<i>Cimetidine</i>	5.2	2.6	4.7	3.0	13.8	10.1	105.3	98.2	108.8	99.0	96.0	102.9
<i>Iopromide</i>	7.0	5.2	2.6	14.7	10.3	11.2	94.0	100.6	109.0	100.1	103.8	112.2
<i>Buprenorphine</i>	13.7	9.4	-	11.4	6.5	-	80.1	96.3	-	92.2	96.7	-
<i>Ephedrine/pseudoephedrine</i>	4.5	3.6	4.1	2.5	4.1	3.5	92.3	95.7	94.1	99.5	95.1	97.4
<i>Norephedrine</i>	5.6	4.1	3.2	7.2	4.8	3.3	101.2	92.3	95.4	97.6	89.5	98.6
<i>Azathioprine</i>	14.8	5.8	2.3	14.6	11.0	16.2	93.9	106.3	92.2	100.1	99.6	92.6
<i>Methotrexate</i>	11.9	7.5	6.8	0.8	5.0	6.5	108.1	103.6	112.3	119.7	110.1	106.8
<i>Ifosfamide</i>	2.9	3.1	1.3	4.0	3.3	0.7	94.5	90.6	95.6	96.4	92.5	96.9
<i>Tamoxifen</i>	4.5	1.8	5.6	2.2	0.7	4.3	88.8	106.1	93.2	93.0	105.2	92.3
<i>Imatinib</i>	-	1.5	-	-	2.5	-	-	101.3	-	-	103.8	-
<i>Capecitabine</i>	3.4	3.7	1.4	2.7	2.1	2.1	79.9	93.9	95.2	79.8	93.4	94.4
<i>Bicalutamide</i>	-	2.9	2.9	-	1.9	3.5	-	80.4	103.6	-	79.9	100.2
<i>Ketamine</i>	2.7	1.7	1.1	0.4	1.8	1.6	75.6	103.6	98.5	77.3	105.3	98.2

Compound	Inter-day precision (%)			Intra-day precision (%)			Inter-day Accuracy (%)			Intra-day Accuracy (%)		
	Low	Med	High	Low	Med	High	Low	Med	High	Low	Med	High
Norketamine	2.9	1.9	0.6	4.8	3.7	0.9	80.3	99.4	102.4	82.2	98.1	101.8
Venlafaxine	3.9	2.0	1.6	1.4	1.3	2.4	74.6	96.2	102.7	72.9	96.7	101.7
Desmethylvenlafaxine	4.0	2.3	2.0	2.4	2.4	1.4	101.4	101.3	101.3	105.9	102.2	98.9
Fluoxetine	1.5	2.1	1.4	2.2	1.3	1.9	91.1	101.1	98.1	92.3	103.0	99.6
Norfluoxetine	2.5	1.6	0.4	4.1	2.9	2.3	107.0	103.3	97.8	106.3	104.8	98.1
Sertraline	3.0	1.0	0.9	0.7	1.7	2.7	84.0	100.8	101.1	84.3	101.9	100.9
Mirtazapine	3.9	4.2	2.1	2.2	3.0	2.9	90.3	99.0	95.0	94.0	101.6	97.3
Citalopram	0.5	0.8	0.8	1.7	4.9	1.2	101.6	100.5	101.7	102.0	101.3	102.1
Desmethylcitalopram	0.1	2.9	2.4	1.7	4.4	3.0	122.1	99.1	97.8	112.2	99.4	98.7
Paroxetine	-	1.5	1.1	-	2.6	3.8	-	101.4	102.7	-	105.4	101.4
Duloxetine	-	13.1	14.0	-	2.9	3.1	-	80.7	75.9	-	97.1	85.4
Amitriptyline	3.6	1.5	2.3	4.3	4.1	4.9	90.8	100.4	99.2	92.2	105.7	100.8
Nortriptyline	4.9	2.8	1.6	5.0	2.6	4.3	78.2	101.4	99.2	80.7	106.7	98.9
Norsertaline	15.2	6.8	-	13.1	4.2	-	76.9	106.7	-	85.8	112.3	-
Carbamazepine	2.6	1.7	1.6	2.4	1.1	1.3	84.1	94.8	96.1	86.3	94.5	97.2
Carbamazepine 10,11-epoxide	0.9	0.1	3.7	2.0	2.4	2.0	94.8	80.4	91.3	94.0	80.5	95.1
10,11-Dihydro-10-												
hydroxycarbamazepine	3.7	2.0	2.8	5.1	3.8	7.9	94.6	96.5	85.4	94.8	98.6	88.1
Diltiazem	0.4	1.5	4.9	3.2	2.2	1.4	93.8	93.4	91.0	94.0	94.6	92.1
Verapamil	-	2.3	2.6	-	3.5	2.3	-	104.9	98.9	-	108.5	97.8
Temazepam	1.7	1.0	0.4	1.1	1.4	2.4	92.3	98.0	100.6	93.8	98.7	101.1
Oxazepam	6.3	2.3	1.8	3.4	3.6	3.0	82.8	100.7	99.4	81.3	103.7	99.5
Diazepam	4.7	5.0	3.7	2.8	0.8	1.3	105.6	99.1	94.2	107.2	101.0	93.9
Quetiapine	1.8	1.1	1.4	1.0	1.0	1.7	83.4	101.3	101.3	85.1	102.5	101.5
Risperidone	-	1.2	-	-	3.2	-	-	96.8	-	-	101.6	-
Donepezil	-	1.8	0.7	-	3.3	1.9	-	111.0	104.3	-	115.5	106.1
Memantine	-	0.6	1.2	-	2.8	4.1	-	105.9	102.7	-	110.4	102.2
Creatinine	1.1	0.9	2.3	2.5	1.9	3.9	99.3	110.8	97.2	93.3	109.7	97.5
Nicotine	1.1	2.1	0.5	2.6	4.1	0.4	102.2	90.0	102.6	103.2	89.6	102.3
Caffeine	2.0	1.8	1.3	1.9	5.4	0.9	113.2	90.1	95.4	115.2	89.5	96.5
Cotinine	0.8	2.2	1.6	1.5	1.2	1.8	94.4	99.6	101.3	95.1	100.9	100.3
1,7-dimethylxanthine	13.0	2.9	2.1	6.8	9.0	13.9	84.2	99.9	98.7	90.0	98.4	96.3
Morphine	6.4	0.2	2.1	3.9	0.1	3.4	94.6	106.3	96.5	87.7	106.3	98.5
Dihydromorphine	3.8	2.5	6.8	0.4	2.5	5.3	110.6	102.0	105.3	115.1	104.0	106.3
Normorphine	2.2	1.5	0.9	3.6	1.1	1.9	101.1	102.7	99.1	98.6	102.3	98.5
Methadone	0.7	2.3	-	1.7	1.2	-	100.4	97.1	-	101.1	99.3	-
EDDP	1.4	0.5	1.8	1.3	1.2	1.0	92.6	102.1	94.5	91.2	102.6	95.5

Compound	Inter-day precision (%)			Intra-day precision (%)			Inter-day Accuracy (%)			Intra-day Accuracy (%)		
	Low	Med	High	Low	Med	High	Low	Med	High	Low	Med	High
<i>Codéine</i>	3.8	1.0	1.1	3.5	4.3	4.1	87.6	98.2	94.6	91.4	99.0	94.8
<i>Norcodéine</i>	3.5	2.4	2.4	6.7	4.2	3.4	105.7	94.7	95.2	106.3	95.1	94.2
<i>Dihydrocodéine</i>	2.7	0.5	1.4	0.8	1.7	3.7	84.8	103.4	94.4	87.4	102.8	93.4
<i>Tramadol</i>	1.7	1.7	1.5	1.2	2.6	1.9	100.6	100.6	99.0	99.0	99.0	97.3
<i>N-desmethyiltramadol</i>	3.4	3.0	1.1	3.2	1.2	2.2	91.1	92.8	93.6	94.7	93.9	94.7
<i>O-desmethyiltramadol</i>	2.0	4.6	-	3.1	6.7	-	92.6	97.9	-	94.2	102.8	-
<i>Amphetamine</i>	0.8	2.1	1.5	1.3	1.3	2.2	100.3	102.9	99.2	100.8	101.6	99.8
<i>Methamphetamine</i>	1.6	0.2	1.3	1.2	0.3	2.4	99.0	101.7	102.3	100.4	101.5	101.3
<i>MDMA</i>	0.6	1.2	2.1	3.3	0.5	1.2	102.5	96.7	98.3	103.2	97.9	98.3
<i>MDA</i>	1.4	1.1	0.7	0.6	0.4	1.1	98.4	99.6	97.3	100.0	100.8	99.3
<i>Cocaine</i>	2.2	2.7	1.7	1.1	1.1	2.2	96.7	98.1	96.8	98.2	100.4	98.4
<i>Benzoylcegonine</i>	1.4	3.3	2.5	0.7	1.1	0.9	104.3	101.6	104.4	105.2	103.2	101.4
<i>Anthydrocegoninemethyl ester</i>	1.6	2.5	2.8	1.2	2.3	3.8	109.1	90.8	103.5	107.1	88.5	100.4
<i>Cocacethylene</i>	5.9	0.8	1.6	2.0	1.2	1.8	88.2	101.8	95.2	88.3	101.2	94.5
<i>Mephedrone</i>	1.8	0.9	2.8	4.9	3.0	0.9	77.3	89.6	94.4	75.9	89.9	91.4
<i>MDPV</i>	4.2	2.3	0.2	0.3	0.6	1.1	105.8	94.1	99.0	108.7	96.5	99.1
<i>Heroin</i>	2.7	2.2	0.8	2.5	1.3	1.7	98.3	97.0	99.4	100.4	98.8	98.6
<i>6-acetyl morphine</i>	1.1	9.6	7.5	3.5	3.4	8.3	104.7	75.0	106.3	105.0	83.0	112.3
<i>Thiamethoxam</i>	8.6	2.3	-	5.6	3.8	-	93.0	100.9	-	89.9	97.7	-
<i>Imidacloprid</i>	-	5.7	5.2	-	3.1	2.5	-	110.2	96.8	-	107.9	93.0
<i>Clothianidin</i>	4.1	3.7	2.0	3.5	2.8	3.3	88.7	106.4	100.8	91.1	105.3	97.4
<i>Metazachlor</i>	-	3.2	2.0	-	1.5	3.5	-	108.2	101.1	-	112.1	99.9
<i>Terbutylazine</i>	1.1	2.5	6.5	1.8	2.3	2.9	91.4	103.7	97.3	94.3	106.1	98.9
<i>Methiocarb</i>	1.4	1.2	2.8	1.5	2.2	2.2	100.6	102.7	98.5	101.1	105.4	98.8
<i>Dichlofluanid</i>	-	3.2	5.5	-	1.3	6.3	-	89.1	92.6	-	93.3	96.5
<i>Flufenacet</i>	1.4	3.1	4.4	1.3	2.3	2.3	125.3	96.6	96.8	121.7	97.5	93.5
<i>Oxadiazon</i>	5.0	0.6	-	3.4	4.6	-	90.2	104.1	-	91.2	99.8	-
<i>Chlorpyrifos</i>	10.4	5.3	-	12.1	11.5	-	71.8	94.9	-	73.1	88.2	-
<i>Triallate</i>	13.2	-	-	7.6	-	-	70.6	-	-	81.3	-	-
<i>Tylosin</i>	2.9	2.2	1.6	5.8	2.8	3.3	104.9	92.2	101.3	108.2	92.7	99.8
<i>Sulfapyridine</i>	1.5	0.5	1.3	2.6	2.3	3.1	145.9	86.6	96.1	143.7	89.4	99.2
<i>Sarafloxacin</i>	3.9	1.8	1.3	5.6	4.8	5.2	138.5	85.0	97.8	134.4	92.9	109.0
<i>Ceftiofur</i>	-	3.4	0.6	-	5.3	2.0	-	81.0	91.7	-	83.4	95.5
<i>Diazinon</i>	3.6	4.2	4.5	3.1	2.3	2.7	101.9	89.1	97.0	102.3	93.2	101.2

Table S8 Liquid matrix absolute and corrected recoveries

Compound	Surface water			Absolute Effluent			Influent			Surface water			Corrected Effluent			Influent		
	Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD	
<i>Benzophenone-1</i>	81.3	0.1		80.2	1.1		49.5	1.4		99.9	14.8		89.2	11.9		75.6	2.4	
<i>Benzophenone-2</i>	31.0	14.4		29.0	1.7		26.7	1.1		44.6	19.3		40.3	14.0		45.4	3.0	
<i>Benzophenone-3</i>	33.0	52.5		51.9	10.0		27.0	32.2		74.4	33.7		89.4	29.9		58.5	30.1	
<i>Benzophenone-4</i>	72.8	3.0		52.7	5.0		38.9	5.0		92.6	33.0		122.8	5.0		101.9	5.0	
<i>Methylparaben</i>	70.2	27.1		59.4	27.4		39.9	37.5		96.3	1.6		100.8	3.3		99.8	5.0	
<i>Ethylparaben</i>	70.1	0.5		73.0	5.1		68.5	0.8		77.5	6.2		83.5	3.8		87.5	8.7	
<i>Propylparaben</i>	69.8	5.8		75.1	0.7		55.6	4.2		117.4	3.8		125.3	7.8		112.6	5.0	
<i>Butylparaben</i>	76.8	0.1		82.3	1.8		47.9	13.9		115.4	9.6		118.3	2.6		102.4	2.0	
<i>Bisphenol A</i>	60.5	3.9		56.8	22.0		37.4	47.1		108.7	14.0		100.7	8.2		112.2	22.3	
<i>E1</i>	63.2	1.7		64.3	3.3		50.6	7.3		97.2	3.6		92.7	2.7		99.1	5.0	
<i>E2</i>	52.9	23.0		67.5	10.9		51.5	15.1		79.3	26.3		87.9	11.4		106.6	13.0	
<i>EE2</i>	48.3	6.1		64.9	11.0		51.9	13.1		86.3	5.4		84.6	7.0		87.9	12.0	
<i>Sulfasalazine</i>	31.8	6.4		28.4	0.2		21.8	3.3		64.2	10.4		58.7	9.8		62.7	0.4	
<i>Clarithromycin</i>	31.0	33.9		40.1	1.2		33.2	20.4		77.6	2.4		97.2	7.7		80.9	0.5	
<i>Azithromycin</i>	21.9	11.4		24.5	6.3		22.2	12.4		86.9	17.0		67.0	11.3		101.1	29.1	
<i>Trimethoprim</i>	59.2	27.5		59.8	0.8		41.4	6.6		118.0	18.3		123.7	5.4		109.6	17.1	
<i>Sulfamethoxazole</i>	78.9	8.5		64.0	14.5		42.1	36.5		108.4	11.0		111.0	13.4		102.8	42.3	
<i>Triclosan</i>	58.3	5.0		75.1	5.0		69.3	5.0		86.5	5.0		100.8	5.0		117.5	5.0	
<i>Amoxicillin</i>	0.4	1.2		30.1	23.1		-	-		-	-		72.3	5.1		-	-	
<i>Metronidazole</i>	63.3	6.9		56.3	10.1		51.8	11.0		107.3	8.0		92.3	14.2		110.2	7.9	
<i>Sulfadiazine</i>	56.0	3.7		28.2	3.2		23.2	4.1		95.4	8.1		59.0	11.4		56.6	8.7	
<i>Cefalexin</i>	30.7	N/A		429.0	413.2		-36.0	135.4		67.2	-		-	-		140.6	5.2	
<i>Ofloxacin</i>	66.3	7.4		50.9	12.1		56.7	7.5		97.9	5.1		95.2	7.9		118.3	16.1	
<i>Ciprofloxacin</i>	28.2	6.6		28.1	3.4		52.3	38.6		95.4	15.2		69.2	6.9		101.3	7.6	
<i>Tetracycline</i>	18.0	14.5		12.3	24.9		14.0	26.9		63.7	11.3		63.0	12.0		108.0	10.2	
<i>Danofloxacin</i>	29.6	6.2		30.4	4.7		30.3	2.9		99.9	6.1		71.0	16.1		87.3	9.7	

Compound	Surface water			Absolute Effluent			Influent			Surface water			Corrected Effluent			Influent		
	Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD	
<i>Oxytetracycline</i>	39.3	6.4	42.3	14.1	101.3	133.7	58.6	6.5	70.3	26.0	85.6	22.1						
<i>Chloramphenicol</i>	131.4	119.6	26.6	35.9	85.9	23.8	82.0	6.8	80.0	16.6	123.8	39.7						
<i>Penicillin G</i>	6.2	3.6	7.3	1.3	3.9	5.1	12.3	0.3	-	-	-	-						
<i>Penicillin V</i>	110.1	22.7	105.5	8.1	25.1	27.8	130.7	23.3	158.7	12.5	60.6	36.1						
<i>Erythromycin</i>	923.8	52.9	858.0	18.4	930.1	20.7	84.8	6.0	82.9	4.3	87.7	0.5						
<i>Prulifloxacin</i>	1.9	3.5	29.2	15.1	44.7	10.7	-	-	47.5	16.9	69.1	7.8						
<i>Norfloxacin</i>	44.7	7.3	48.4	7.4	42.5	5.5	64.8	22.8	83.9	22.9	77.9	15.9						
<i>Griseofulvin</i>	98.2	5.9	89.9	2.3	84.8	1.9	122.9	18.6	149.0	18.8	131.3	15.4						
<i>Ketoconazole</i>	8.9	5.9	67.2	51.7	67.9	40.3	41.4	23.7	167.8	39.5	141.2	30.8						
<i>Valsartan</i>	60.6	4.1	53.1	4.4	46.9	9.1	76.0	34.3	79.9	36.4	101.1	22.4						
<i>Irbesartan</i>	56.2	8.8	53.6	10.1	40.3	12.2	112.3	13.8	107.4	17.3	97.6	13.5						
<i>Lisinopril</i>	21.4	33.2	21.9	4.3	28.6	61.0	93.7	9.5	88.5	8.8	102.3	45.2						
<i>Ketoprofen</i>	72.9	1.1	67.8	2.0	45.6	0.0	92.0	5.2	109.2	24.5	107.2	3.8						
<i>Ibuprofen</i>	85.6	5.6	126.2	18.4	57.9	5.0	115.5	27.6	110.8	18.3	78.9	5.0						
<i>Naproxen</i>	80.5	2.9	84.5	5.0	15.7	5.0	91.0	6.6	98.0	5.0	82.7	5.0						
<i>Diclofenac</i>	70.5	4.2	71.8	5.8	46.5	33.4	120.0	12.7	116.0	1.1	103.9	9.9						
<i>Acetaminophen</i>	44.4	5.0	44.8	5.0	77.8	10.4	140.5	5.0	116.4	5.0	77.8	45.9						
<i>Bezafibrate</i>	73.7	4.3	78.3	18.7	46.3	5.0	94.2	12.1	92.9	0.7	91.4	5.0						
<i>Atorvastatin</i>	35.5	15.9	59.9	2.2	58.6	8.4	54.7	14.9	91.9	2.2	94.0	7.2						
<i>Gemfibrozil</i>	177.3	49.3	129.5	28.1	66.6	12.7	176.0	42.1	167.5	16.2	94.3	15.5						
<i>Candesartan Cilexetil</i>	32.0	7.8	289.8	149.1	401.6	59.5	114.9	20.4	-	-	-	-						
<i>Fexofenadine</i>	67.6	5.5	71.0	22.9	50.7	17.2	74.9	13.4	83.8	26.7	107.3	41.9						
<i>Cetirizine</i>	47.9	50.2	78.8	15.2	48.6	36.0	77.3	44.2	114.3	13.2	77.3	36.2						
<i>Sildenafil</i>	46.9	12.9	79.3	39.3	70.8	26.0	111.1	7.4	116.2	17.2	128.3	7.5						
<i>Metformin</i>	55.2	9.2	53.0	9.2	36.7	23.5	93.3	14.4	96.0	3.3	78.7	38.7						
<i>Gliclazide</i>	32.9	0.8	61.8	7.4	46.4	12.6	69.5	31.3	108.4	2.9	80.1	18.5						
<i>Sitagliptin</i>	38.3	3.1	36.2	3.8	35.1	5.6	125.5	36.8	84.9	19.1	104.3	14.3						
<i>Pholcodine</i>	76.9	5.0	43.1	5.0	13.7	5.0	119.6	5.0	102.0	5.0	42.0	5.0						

Compound	Surface water		Effluent		Influent		Corrected Effluent		Influent			
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD		
<i>Atenolol</i>	76.0	2.0	54.6	5.0	42.8	5.0	115.1	6.1	109.7	5.0	88.6	5.0
<i>Metoprolol</i>	71.8	7.6	53.1	15.0	36.2	16.7	107.8	2.9	104.3	12.9	97.4	12.7
<i>Propranolol</i>	45.4	22.9	36.1	18.9	38.7	12.8	107.2	5.9	102.7	4.8	105.0	9.7
<i>Bisoprolol</i>	73.0	5.9	51.8	7.0	54.3	10.4	167.9	28.8	95.4	33.9	123.4	36.4
<i>Ranitidine</i>	65.0	16.4	46.4	5.0	70.1	72.9	88.6	13.8	106.8	5.0	97.8	35.4
<i>Cimetidine</i>	32.7	7.4	33.4	5.0	20.6	41.0	95.3	5.5	106.2	26.3	123.7	9.7
<i>Iopromide</i>	97.0	11.2	82.1	12.4	47.2	10.1	143.0	33.9	123.7	24.9	138.5	38.7
<i>Buprenorphine</i>	30.6	3.7	39.6	5.7	33.3	10.4	209.6	1.0	228.8	34.9	135.8	11.2
<i>Ephedrine/pseudoephedrine</i>	25.4	15.1	18.7	51.4	23.0	12.3	67.8	14.5	88.6	4.7	92.8	9.2
<i>Norephedrine</i>	28.4	5.0	28.9	5.0	26.9	5.0	78.2	5.0	100.9	5.0	85.4	5.0
<i>Azathioprine</i>	88.9	29.7	81.5	33.9	72.0	25.6	66.4	8.7	83.8	15.7	109.5	3.2
<i>Methotrexate</i>	22.6	11.1	30.7	29.0	39.1	109.1	95.7	6.3	114.1	15.2	98.8	37.8
<i>Isofamide</i>	62.9	24.6	41.9	30.2	33.4	24.6	90.8	18.5	80.2	31.4	89.2	22.4
<i>Tamoxifen</i>	0.2	5.3	8.7	57.5	9.5	43.8	112.8	10.2	98.4	25.3	85.2	27.4
<i>Imatinib</i>	4.8	0.6	44.5	37.7	-12.8	41.0	47.2	4.1	73.6	6.1	46.6	3.2
<i>Capecitabine</i>	62.3	8.8	70.4	6.0	60.4	2.7	73.3	10.5	111.5	19.3	104.8	20.6
<i>Bicalutamide</i>	91.8	32.0	108.2	36.7	104.7	27.6	71.2	24.3	99.8	33.9	95.2	32.9
<i>Ketamine</i>	68.2	7.5	53.8	10.8	41.7	10.9	107.7	0.5	115.6	12.4	111.6	3.2
<i>Norketamine</i>	65.7	17.6	53.9	22.0	42.2	21.2	103.2	1.9	108.6	6.8	98.9	1.8
<i>Venlafaxine</i>	58.8	23.2	36.3	19.4	23.8	3.1	96.8	5.7	89.5	22.1	70.2	5.8
<i>Desmethylvenlafaxine</i>	62.5	24.2	45.9	5.0	35.9	5.0	107.4	0.1	100.0	5.0	101.5	5.0
<i>Fluoxetine</i>	4.4	31.8	7.1	5.5	19.8	6.1	112.3	1.0	111.1	17.6	110.7	13.3
<i>Norfluoxetine</i>	3.0	42.2	7.9	4.3	23.6	1.5	74.4	25.9	91.7	3.6	94.8	19.0
<i>Sertraline</i>	3.1	13.4	8.3	29.5	13.4	6.0	108.0	2.9	105.0	13.8	91.5	8.1
<i>Mirtazapine</i>	57.5	30.1	39.8	26.1	25.8	34.1	99.1	20.6	97.0	25.7	84.7	26.8
<i>Citalopram</i>	41.1	28.8	35.5	5.0	40.3	4.7	115.0	11.6	124.1	5.0	119.9	14.9
<i>Desmethylcitalopram</i>	36.2	5.0	27.5	25.6	32.5	35.1	89.5	5.0	110.9	0.9	93.1	15.6

Compound	Absolute				Corrected							
	Surface water		Effluent		Surface water		Effluent					
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD				
<i>Paroxetine</i>	4.4	1.5	9.8	1.8	19.9	3.3	29.1	3.5	49.7	4.9	76.2	4.7
<i>Duloxetine</i>	5.1	3.3	8.5	4.5	18.7	8.3	39.8	7.7	55.3	7.6	83.2	4.7
<i>Amitriptyline</i>	17.6	5.5	20.7	3.1	25.7	6.7	96.7	5.4	97.5	13.2	104.5	11.1
<i>Nortriptyline</i>	12.6	3.6	16.1	3.1	23.4	5.7	98.4	14.9	103.7	17.8	108.0	17.5
<i>Norsertaline</i>	-	-	-	-	31.8	13.3	-	-	-	-	62.9	4.0
<i>Carbamazepine</i>	68.2	9.7	55.5	13.0	37.4	9.3	101.7	7.0	102.5	13.8	88.3	11.1
<i>Carbamazepine 10,11-epoxide</i>	94.8	5.0	55.0	5.0	56.7	5.0	136.1	5.0	126.9	5.0	143.7	5.0
<i>10,11-Dihydro-10-hydroxycarbamazepine</i>	74.1	13.8	59.5	26.4	50.3	30.2	125.1	4.7	120.6	3.7	152.2	0.1
<i>Diltiazem</i>	43.7	28.2	30.1	24.0	36.3	16.1	82.9	0.0	82.3	15.7	118.1	9.1
<i>Verapamil</i>	25.0	8.2	34.0	4.5	31.8	5.3	95.1	4.4	102.2	20.0	111.2	16.7
<i>Temazepam</i>	66.3	15.5	72.5	0.3	54.5	9.3	115.3	8.5	113.9	7.0	107.3	8.5
<i>Oxazepam</i>	92.2	9.1	82.6	9.8	75.1	3.7	105.1	11.2	105.5	17.0	115.0	12.7
<i>Diazepam</i>	83.5	6.3	70.5	6.3	70.8	2.7	92.4	11.5	85.8	10.2	95.6	9.9
<i>Quetiapine</i>	51.6	28.9	46.8	23.5	37.8	25.4	87.2	21.2	95.1	18.6	72.6	18.7
<i>Risperidone</i>	35.3	6.1	38.1	1.9	30.3	3.8	171.6	19.6	119.9	12.2	116.1	9.4
<i>Donepezil</i>	42.7	8.5	42.1	4.3	41.6	8.6	158.0	17.7	112.9	15.5	117.3	8.4
<i>Memantine</i>	58.6	6.8	42.5	8.0	35.6	8.9	181.6	28.8	139.0	24.8	127.8	10.2
<i>Creatinine</i>	31.0	10.9	24.8	13.0	22.8	17.0	59.6	0.6	42.9	11.8	46.9	44.7
<i>Nicotine</i>	45.3	4.8	55.7	15.1	98.5	11.7	54.9	39.6	83.9	56.1	77.3	61.8
<i>Caffeine</i>	137.0	40.8	89.8	5.0	82.9	11.1	122.8	5.7	131.7	5.0	96.6	15.4
<i>Cotinine</i>	72.0	19.8	47.3	26.8	37.4	5.0	105.9	2.1	101.5	13.4	85.2	5.0
<i>1,7-dimethylxanthine</i>	47.5	5.0	26.6	5.0	89.3	2.1	75.5	5.0	110.6	5.0	79.6	59.3
<i>Morphine</i>	57.1	14.8	47.8	0.9	34.2	6.0	112.5	8.7	109.8	4.7	93.9	10.4
<i>Dihydromorphine</i>	45.4	32.8	31.4	28.1	19.9	42.9	92.8	4.7	96.2	2.8	95.7	25.2
<i>Normorphine</i>	42.8	24.4	38.6	28.6	30.3	30.1	87.3	5.4	106.5	8.6	109.5	3.2
<i>Methadone</i>	46.5	15.0	47.9	6.9	49.6	10.6	119.2	14.9	111.7	13.0	109.0	0.5
<i>EDDP</i>	23.8	32.8	33.9	10.7	44.1	2.0	105.2	1.0	113.6	0.4	111.6	2.9

Compound	Surface water		Effluent		Influent		Corrected Effluent		Influent			
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD		
<i>Cocaine</i>	67.5	10.6	68.4	56.6	39.0	5.0	93.1	16.2	102.8	8.3	94.8	5.0
<i>Norcocaine</i>	52.5	11.1	36.4	21.2	35.5	10.1	90.6	2.4	87.5	32.5	95.7	11.5
<i>Dihydrocodeine</i>	66.4	6.7	54.8	34.0	34.7	13.9	98.5	2.6	118.6	19.7	85.7	9.9
<i>Tramadol</i>	60.9	24.0	47.0	4.1	33.3	10.8	104.8	0.5	88.3	27.3	69.5	64.8
<i>N-desmethyltramadol</i>	42.2	5.0	33.4	15.0	17.9	5.0	69.1	5.0	98.8	0.7	81.3	5.0
<i>O-desmethyltramadol</i>	58.6	19.1	36.0	26.6	31.8	27.6	108.1	21.9	93.9	6.5	87.3	17.3
<i>Amphetamine</i>	22.5	38.1	27.4	23.9	24.5	15.3	98.7	16.7	107.8	11.2	91.7	18.5
<i>Methamphetamine</i>	47.7	4.4	42.5	15.5	31.9	8.2	97.6	1.4	97.1	12.1	102.1	2.4
<i>MDMA</i>	50.0	11.3	37.0	13.2	29.4	10.8	100.9	1.9	112.9	18.5	118.5	18.7
<i>MDA</i>	28.7	18.9	30.3	19.0	30.7	21.1	99.6	1.1	99.1	13.0	99.4	3.8
<i>Cocaine</i>	71.4	5.5	44.9	9.2	21.7	33.0	106.9	2.5	105.2	2.5	97.4	15.3
<i>Benzoyllecgonine</i>	74.2	15.2	55.1	14.6	46.5	3.2	126.2	6.3	116.4	10.8	106.3	2.3
<i>Anhydroecgoninemethylester</i>	53.6	0.6	50.2	10.3	33.9	5.0	79.5	5.8	109.5	13.3	91.3	8.0
<i>Cocacethylene</i>	71.3	7.8	48.1	14.9	7.6	40.9	105.9	2.6	107.4	15.7	109.9	4.8
<i>Mephedrone</i>	22.8	12.9	22.8	11.8	18.2	25.1	97.3	5.6	103.7	0.3	94.8	10.6
<i>MDPV</i>	112.2	17.7	85.5	14.7	20.8	33.8	78.9	20.8	87.9	21.6	136.8	0.2
<i>Heroin</i>	54.1	20.3	29.1	18.3	23.9	16.0	100.9	7.1	103.0	1.2	105.3	7.6
<i>6-acetylmorphine</i>	52.9	27.4	40.0	30.9	33.9	26.6	72.6	31.8	85.4	32.5	103.6	3.1
<i>Thiamethoxam</i>	53.3	11.3	28.3	3.1	23.9	3.5	70.8	1.7	41.0	2.2	34.1	1.3
<i>Imidacloprid</i>	92.7	10.3	62.8	11.8	60.7	15.2	118.6	2.9	107.4	6.0	107.9	9.3
<i>Clothianidin</i>	76.5	7.5	63.1	9.0	64.0	10.5	111.1	18.7	89.5	16.6	83.4	18.4
<i>Metazachlor</i>	85.7	4.4	64.7	6.7	64.7	7.4	107.3	7.7	96.4	8.9	108.5	5.5
<i>Terbutylazine</i>	72.9	3.4	64.1	2.5	56.2	2.8	83.2	7.2	79.8	9.3	75.7	10.8
<i>Methiocarb</i>	79.4	6.9	62.0	4.5	69.1	2.4	97.7	3.8	91.3	7.6	97.1	3.4
<i>Dichlofluanid</i>	23.3	9.6	2.3	2.9	12.6	3.4	-	-	-	-	51.3	16.9
<i>Flufenacet</i>	82.7	10.5	69.7	9.2	62.9	4.3	82.9	11.5	87.9	12.8	87.9	10.2
<i>Oxadiazon</i>	63.9	21.8	60.6	9.8	52.6	6.2	78.9	19.8	91.5	11.3	79.1	8.1

Compound	Surface water			Absolute Effluent			Influent			Surface water			Corrected Effluent			Influent		
	Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD		Rec (%)	SD	
<i>Chlorpyrifos</i>	15.1	5.1		42.2	N/A		138.0	46.3		21.8	5.9		65.7	-		-	-	
<i>Triallate</i>	32.2	9.4		28.7	6.4		73.2	N/A		45.9	6.0		49.6	4.8		106.7	-	
<i>Tylosin</i>	43.8	32.2		50.3	16.6		34.3	13.9		100.4	13.4		110.0	22.6		98.3	0.3	
<i>Sulfapyridine</i>	76.8	4.7		50.0	5.6		47.1	5.9		101.8	21.1		75.0	1.5		84.3	20.8	
<i>Sarafloxacin</i>	37.3	6.6		30.5	8.8		36.3	7.9		134.0	45.5		84.1	39.4		110.9	36.3	
<i>Ceftiofur</i>	11.0	8.9		30.2	7.0		39.1	14.4		19.4	17.5		63.6	18.6		82.8	30.5	
<i>Diazinon</i>	74.8	17.6		67.8	15.0		60.7	14.1		75.2	11.9		72.5	13.0		81.6	18.0	

Table S9 Solid matrix absolute and corrected recoveries

Compound	Absolute				Corrected			
	SPM		Digested Solids		SPM		Digested Solids	
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD
<i>Benzophenone-1</i>	41.4	8.4	8.1	2.2	120.0	9.1	69.7	58.6
<i>Benzophenone-2</i>	34.3	2.8	11.1	0.8	141.3	9.7	105.5	19.0
<i>Benzophenone-3</i>	-	-	-	-	-	-	-	-
<i>Benzophenone-4</i>	17.5	2.8	7.6	1.8	72.6	5.5	47.0	9.9
<i>Methylparaben</i>	44.2	7.7	18.1	11.9	149.1	22.4	110.3	81.5
<i>Ethylparaben</i>	40.2	7.8	19.5	0.7	109.5	26.7	118.7	13.0
<i>Propylparaben</i>	49.5	4.1	16.0	0.6	200.4	14.6	153.6	30.8
<i>Butylparaben</i>	36.8	2.8	11.2	0.7	247.0	10.4	179.6	14.3
<i>Bisphenol A</i>	35.7	10.0	11.9	-	58.8	12.0	95.7	-
<i>E1</i>	16.2	1.9	5.9	2.1	115.3	8.5	106.6	28.0
<i>E2</i>	14.0	1.5	6.4	0.8	109.4	12.6	113.6	9.9
<i>EE2</i>	-	-	-	-	-	-	-	-
<i>Sulfasalazine</i>	-	-	-	-	-	-	-	-
<i>Clarithromycin</i>	-	-	-	-	-	-	-	-
<i>Azithromycin</i>	7.6	6.0	14.0	12.0	22.4	4.8	60.5	17.7
<i>Trimethoprim</i>	69.3	17.2	44.9	2.5	146.0	30.9	138.5	6.6
<i>Sulfamethoxazole</i>	36.1	2.9	24.4	2.4	62.7	2.3	74.4	5.8
<i>Triclosan</i>	-	-	-	-	-	-	-	-
<i>Amoxicillin</i>	22.9	18.7	26.3	25.2	-	-	-	-
<i>Metronidazole</i>	8.0	1.6	5.6	0.3	121.2	18.4	103.1	6.9
<i>Sulfadiazine</i>	9.9	1.7	-	-	151.6	19.3	166.0	36.4
<i>Cefalexin</i>	-	-	-	-	-	-	-	-
<i>Ofloxacin</i>	6.3	3.9	-	-	-	-	-	-
<i>Ciprofloxacin</i>	9.3	3.7	0.9	1.0	-	-	-	-
<i>Tetracycline</i>	-	-	12.5	-	-	-	-	-
<i>Danofloxacin</i>	-	-	0.3	0.4	-	-	5.6	0.7
<i>Oxytetracycline</i>	-	-	-	-	-	-	-	-
<i>Chloramphenicol</i>	26.3	2.3	17.2	3.5	125.3	9.2	179.5	19.5
<i>Penicillin G</i>	-	-	-	-	-	-	-	-
<i>Penicillin V</i>	1.9	0.7	29.0	18.4	8.7	2.3	-	-
<i>Erythromycin</i>	-	-	-	-	-	-	-	-
<i>Prulifloxacin</i>	-0.8	-	-41.8	22.9	156.0	-	-	-
<i>Norfloxacin</i>	-	-	-	-	-	-	-	-
<i>Griseofulvin</i>	17.8	3.1	6.6	1.2	82.0	10.1	60.6	9.5
<i>Ketoconazole</i>	8.5	1.7	-3.2	18.3	13.3	1.8	100.6	17.8
<i>Valsartan</i>	-	-	-	-	-	-	-	-
<i>Irbesartan</i>	-	-	-	-	-	-	-	-
<i>Lisinopril</i>	67.6	6.7	37.7	2.9	109.1	10.9	127.8	14.5
<i>Ketoprofen</i>	28.9	1.0	23.1	3.8	98.2	5.0	119.4	8.6
<i>Ibuprofen</i>	42.3	13.9	14.9	4.0	101.0	7.8	105.5	9.6
<i>Naproxen</i>	44.1	9.9	16.4	2.3	96.5	14.6	107.9	9.6
<i>Diclofenac</i>	41.9	12.3	4.2	1.7	88.3	13.9	40.8	14.6
<i>Acetaminophen</i>	14.8	3.6	3.9	0.5	129.8	25.8	57.7	10.8

Compound	Absolute				Corrected			
	SPM		Digested Solids		SPM		Digested Solids	
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD
<i>Bezafibrate</i>	32.3	2.5	16.4	2.2	98.0	8.7	112.5	5.5
<i>Atorvastatin</i>	-	-	-	-	-	-	-	-
<i>Gemfibrozil</i>	7.5	3.9	2.9	0.3	-	-	26.3	1.6
<i>Candesartan Cilexetil</i>	-	-	-	-	-	-	-	-
<i>Fexofenadine</i>	-	-	-	-	-	-	-	-
<i>Cetirizine</i>	-	-	-	-	-	-	-	-
<i>Sildenafil</i>	35.6	2.2	26.2	1.5	110.9	6.4	116.7	4.9
<i>Metformin</i>	-	-	-	-	-	-	-	-
<i>Gliclazide</i>	-	-	-	-	-	-	-	-
<i>Sitagliptin</i>	42.0	3.7	30.0	3.8	79.4	7.9	99.2	8.1
<i>Pholcodine</i>	29.7	2.6	22.8	2.7	62.0	4.2	87.8	11.0
<i>Atenolol</i>	52.1	4.6	30.8	2.7	105.2	9.4	121.0	10.8
<i>Metoprolol</i>	57.7	2.1	36.6	1.6	95.1	3.9	99.7	4.9
<i>Propranolol</i>	42.7	5.6	20.9	4.8	103.3	8.4	91.5	14.0
<i>Bisoprolol</i>	38.1	5.9	26.6	3.8	125.8	9.9	126.1	8.6
<i>Ranitidine</i>	61.7	8.3	21.5	1.9	117.8	17.8	82.1	3.6
<i>Cimetidine</i>	-	-	-	-	-	-	-	-
<i>Iopromide</i>	-	-	-	-	-	-	-	-
<i>Buprenorphine</i>	18.1	3.6	12.4	2.0	59.9	9.4	88.0	8.5
<i>Ephedrine/pseudoephedrine</i>	40.0	3.1	28.2	5.2	72.9	2.8	97.9	21.0
<i>Norephedrine</i>	37.2	3.9	27.0	5.6	64.0	5.1	92.6	18.9
<i>Azathioprine</i>	-	-	-	-	-	-	-	-
<i>Methotrexate</i>	48.2	1.5	16.9	2.6	87.1	2.2	59.3	8.7
<i>Ifosfamide</i>	-	-	-	-	-	-	-	-
<i>Tamoxifen</i>	15.9	9.5	0.3	0.0	118.7	49.7	112.7	7.0
<i>Imatinib</i>	9.0	2.1	10.5	2.4	39.9	8.4	65.3	10.8
<i>Capecitabine</i>	4.8	2.5	5.0	1.0	22.1	0.5	53.4	5.8
<i>Bicalutamide</i>	23.0	3.5	21.8	4.0	87.0	8.3	176.0	30.6
<i>Ketamine</i>	48.5	3.3	29.0	0.9	101.7	4.3	105.3	4.3
<i>Norketamine</i>	41.8	5.4	30.3	1.6	96.3	3.6	109.7	4.3
<i>Venlafaxine</i>	40.7	10.1	11.3	22.5	70.3	20.8	61.9	6.4
<i>Desmethylvenlafaxine</i>	58.2	5.2	33.9	3.8	104.3	5.6	96.9	10.3
<i>Fluoxetine</i>	28.2	1.6	9.4	2.7	104.1	0.9	121.2	48.4
<i>Norfluoxetine</i>	34.7	6.0	7.4	0.9	120.3	18.5	95.1	37.5
<i>Sertraline</i>	72.1	98.5	5.8	6.2	-	-	101.9	14.8
<i>Mirtazapine</i>	33.0	2.0	18.5	2.3	89.2	4.6	98.2	9.9
<i>Citalopram</i>	46.4	2.5	30.4	5.1	105.3	17.7	112.7	37.3
<i>Desmethylcitalopram</i>	41.5	3.0	20.8	3.0	96.4	10.6	82.2	18.8
<i>Paroxetine</i>	18.3	2.5	15.8	3.4	83.5	6.1	103.0	15.6
<i>Duloxetine</i>	10.2	1.8	9.2	2.1	48.3	6.8	61.3	9.3
<i>Amitriptyline</i>	27.9	3.5	8.9	9.1	72.0	9.2	175.2	39.5
<i>Nortriptyline</i>	27.3	1.7	16.5	3.7	104.8	4.0	194.4	1.8
<i>Norsertaline</i>	22.0	32.0	19.2	12.2	39.7	2.7	-	-
<i>Carbamazepine</i>	30.7	1.9	10.8	0.9	98.7	3.3	94.5	13.9

Compound	Absolute				Corrected			
	SPM		Digested Solids		SPM		Digested Solids	
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD
<i>Carbamazepine 10,11-epoxide</i>	-	-	-	-	-	-	-	-
<i>10,11-Dihydro -10-hydroxycarbamazepine</i>	31.6	3.5	11.5	0.6	101.6	8.3	100.1	2.8
<i>Diltiazem</i>	-	-	-	-	-	-	-	-
<i>Verapamil</i>	46.3	5.9	34.3	2.1	102.5	10.6	130.8	6.7
<i>Temazepam</i>	9.5	2.2	6.1	1.6	58.5	5.5	92.3	6.6
<i>Oxazepam</i>	0.5	N/A	1.4	0.7	115.0		115.1	15.5
<i>Diazepam</i>	24.4	4.0	17.3	2.4	82.2	6.8	99.3	1.7
<i>Quetiapine</i>	51.6	3.5	19.1	3.1	113.9	5.3	81.9	9.0
<i>Risperidone</i>	36.3	7.0	15.7	4.2	75.8	6.6	57.4	12.8
<i>Donepezil</i>	41.6	8.7	26.8	4.2	97.2	11.6	98.3	11.0
<i>Memantine</i>	20.0	1.6	19.6	3.7	36.6	9.2	71.4	10.6
<i>Creatinine</i>	-	-	-	-	-	-	-	-
<i>Nicotine</i>	50.1	7.8	45.6	4.0	92.7	12.5	165.4	24.4
<i>Caffeine</i>	-	-	-	-	-	-	-	-
<i>Cotinine</i>	10.9	1.5	4.1	1.7	102.7	2.8	112.1	11.6
<i>1,7-dimethylxanthine</i>	-	-	-	-	-	-	-	-
<i>Morphine</i>	27.1	3.8	15.8	6.3	137.0	28.6	114.0	30.9
<i>Dihydromorphine</i>	13.2	1.9	11.0	1.2	64.6	3.9	74.3	3.4
<i>Normorphine</i>	25.1	4.3	17.4	4.0	127.2	14.6	109.9	17.2
<i>Methadone</i>	46.3	2.4	29.2	0.7	95.5	3.8	102.8	3.2
<i>EDDP</i>	35.6	3.4	25.5	2.8	98.2	5.7	103.8	9.3
<i>Codeine</i>	58.7	8.3	30.2	1.9	116.5	25.7	105.3	19.0
<i>Norcodeine</i>	39.1	5.3	24.0	1.9	77.9	5.6	87.8	10.8
<i>Dihydrocodeine</i>	44.5	4.4	27.8	2.1	91.1	5.6	101.0	10.4
<i>Tramadol</i>	47.1	3.0	30.7	2.1	80.8	3.6	87.7	6.1
<i>N-desmethyltramadol</i>	44.1	4.5	24.7	1.5	84.3	4.1	83.4	10.7
<i>O-desmethyltramadol</i>	-	-	-	-	-	-	-	-
<i>Amphetamine</i>	60.8	3.7	34.1	2.5	110.7	4.2	118.3	8.7
<i>Methamphetamine</i>	53.4	2.1	33.8	1.8	113.2	5.0	104.4	2.1
<i>MDMA</i>	45.0	2.5	28.1	1.3	96.0	2.6	105.6	1.9
<i>MDA</i>	-	-	-	-	-	-	-	-
<i>Cocaine</i>	52.0	2.3	32.8	1.6	99.2	3.7	104.7	2.9
<i>Benzoylcegonine</i>	54.4	4.4	35.7	1.9	100.4	3.9	117.3	7.3
<i>Anhydroecgoninemethylester</i>	-	-	-	-	-	-	-	-
<i>Cocaethylene</i>	48.5	2.9	30.0	1.3	99.8	3.8	98.9	6.5
<i>Mephedrone</i>	23.3	3.2	15.9	1.5	61.9	3.3	99.4	4.1
<i>MDPV</i>	44.0	4.6	24.6	2.7	88.3	5.2	79.7	8.3
<i>Heroin</i>	27.7	1.8	17.9	2.2	98.7	4.9	91.7	6.5
<i>6-acetylmorphine</i>	-	-	-	-	-	-	-	-
<i>Thiamethoxam</i>	20.1	4.8	15.1	3.5	95.3	20.6	140.7	26.1
<i>Imidacloprid</i>	29.0	7.3	24.2	2.9	103.2	21.4	-	-
<i>Clothianidin</i>	32.6	6.8	29.9	3.0	159.3	27.1	-	-
<i>Metazachlor</i>	26.0	3.4	8.9	1.8	92.7	5.5	96.1	8.0

Compound	Absolute				Corrected			
	SPM		Digested Solids		SPM		Digested Solids	
	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD	Rec (%)	SD
<i>Terbutylazine</i>	13.3	1.6	8.4	1.9	45.3	7.6	47.4	4.9
<i>Methiocarb</i>	22.3	1.6	12.6	1.1	104.3	5.5	117.5	3.7
<i>Dichlofluanid</i>	-	-	0.0	0.2	-	-	-	-
<i>Flufenacet</i>	14.7	2.0	5.6	1.0	45.9	4.6	57.9	10.5
<i>Oxadiazon</i>	3.0	1.1	2.5	0.6	14.8	1.9	25.0	4.0
<i>Chlorpyrifos</i>	4.6	2.4	7.9	1.9	-	-	86.1	20.0
<i>Triallate</i>	2.3	0.9	4.1	1.7	-	-	-	-
<i>Tylosin</i>	-	-	-	-	-	-	-	-
<i>Sulfapyridine</i>	26.5	5.1	23.5	4.5	-	-	-	-
<i>Sarafloxacin</i>	6.3	2.8	-	-	-	-	-	-
<i>Ceftiofur</i>	-	-	-	-	-	-	-	-
<i>Diazinon</i>	0.5	0.3	1.1	0.5	145.8	4.5	80.8	9.3

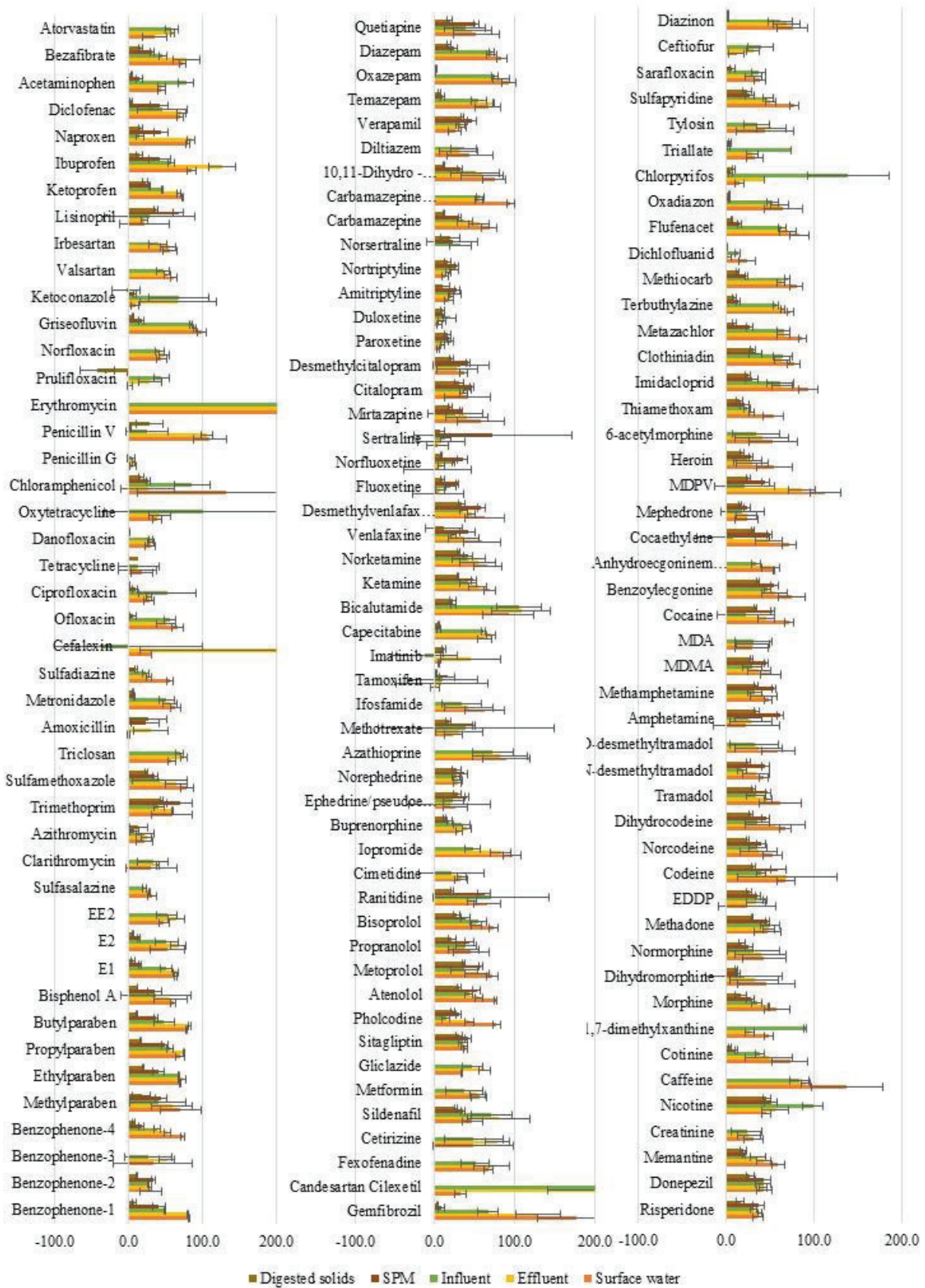


Fig. S2 Absolute recoveries for all matrices, error bars indicate the range of standard deviation

Table S10 Matrix suppression

Compound	Surface water		Effluent		Influent	
	(%)	SD	(%)	SD	(%)	SD
Benzophenone-1	0.6	2	2	4.5	30.6	0.1
Benzophenone-2	43.8	1.1	49.2	1.9	55.5	1
Benzophenone-3	-5.1	1	7.7	3.3	20.7	1.1
Benzophenone-4	-1.2	9.4	23	6.8	42.5	1.7
Methylparaben	41.9	3.6	49.1	1.5	60.2	1.4
Ethylparaben	27.8	3.4	29.8	0.6	37	2.9
Propylparaben	21.7	4.3	21.6	2.1	37.3	4.2
Butylparaben	8	1.9	11	1.4	35.9	2.4
Bisphenol A	15	0.2	12.4	2.6	29.4	2.9
E1	4.5	4.3	5	1.1	28.6	1.4
E2	8.3	8.1	8	5.9	31.8	0.9
EE2	-6.1	1.7	-24.4	3.9	13.2	4.8
Sulfasalazine	-105.8	9.9	-81.7	3.1	-53.4	5.1
Clarithromycin	52	1.6	54.1	0.5	54	1.2
Azithromycin	72.7	0.6	73.9	1.9	68.8	3.9
Trimethoprim	38.2	0.8	51	2.7	49.1	4.5
Sulfamethoxazole	28.4	0.4	38.1	1.5	47.6	0.8
Triclosan	-62.3	0.3	-100.9	3.1	-99	4.8
Amoxicillin	99.8	0.7	93.9	0.7	99.0	0.7
Metronidazole	24.5	1.3	47.4	5.9	51.3	1.5
Sulfadiazine	37.1	1.1	64.6	1.8	71.1	1.1
Cefalexin	37.7	4.8	28.4	12.9	31.5	5.6
Ofloxacin	29.3	9.9	40.2	3.3	48.4	3.3
Ciprofloxacin	67.8	5.4	69.5	1.0	74.1	3.5
Tetracycline	6.5	1.3	27.3	0.9	36.7	9.0
Danofloxacin	70.2	5.7	73.7	2.7	81.3	3.8
Oxytetracycline	13.9	0.8	29.3	3.4	22.5	22.6
Chloramphenicol	33.4	2.4	39.3	1.2	45.7	3.3
Penicillin G	101.1	0.1	100.5	0.2	96.7	0.8
Penicillin V	-15.4	10.8	-13.1	33.7	-2.5	6.9
Erythromycin	36.5	11.2	-195.6	14.7	57.8	6.8
Prulifloxacin	107.9	1.9	87.4	18.7	44.9	26.5
Norfloxacin	52.7	3.9	62.2	3.1	61.5	19.3
Griseofulvin	4.3	0.5	9.9	2.0	17.3	4.4
Ketoconazole	90.4	4.2	44.6	0.1	38.0	3.0
Valsartan	-139.2	2.6	-105.4	5.9	-77.4	5.7
Irbesartan	-23.8	1.2	-25.3	1.3	-14.3	0.2
Lisinopril	-12.3	2.5	8.6	0.5	23.1	0
Ketoprofen	8.5	0.2	9.9	3.2	36.3	1.5
Ibuprofen	-13.6	4.4	-24.3	9.9	-15.6	5
Naproxen	-13.2	1.6	-18.7	0.7	29.8	34.1
Diclofenac	-127.3	0.8	-141.7	5.7	-104.5	2.2
Acetaminophen	27.8	1.9	50.6	0.6	22.2	19
Bezafibrate	-24.5	10.8	-19.6	3.7	0.7	0.5

Compound	Surface water		Effluent		Influent	
	(%)	SD	(%)	SD	(%)	SD
Atorvastatin	-155.7	2.9	-160.3	1.4	-152.8	14.5
Gemfibrozil	-70.3	20.5	-64.4	1.9	29.5	1.0
Candesartan Cilexetil	100.0	23.2	-108.3	37.4	-78.1	47.6
Fexofenadine	10.3	1.4	11.2	1.2	27.1	1.9
Cetirizine	16.9	3	18	0.1	30.5	0.4
Sildenafil	67.2	1.6	32.4	4.7	47.3	1.9
Metformin	44.8	1	47	1.2	63.3	0.6
Gliclazide	19	4.6	23.6	1.1	38.2	1.3
Sitagliptin	23.8	0.0	30.2	0.5	38.4	1.1
Pholcodine	11.4	2.3	45.5	0.6	90	1.4
Atenolol	22.3	1.1	48	0.7	59.7	1.6
Metoprolol	27.6	0.5	49.5	0.1	65.7	0.1
Propranolol	31.2	3.7	59	1.4	56.7	0.6
Bisoprolol	10.5	0.0	37.8	4.0	32.2	1.2
Ranitidine	18.2	4	35	0.8	56.8	1.1
Cimetidine	23.4	4.7	48.4	0.8	61.9	1.6
Iopromide	11.8	2.7	26.1	3.5	39.1	0
Buprenorphine	52.0	2.5	34.1	3.4	45.7	4.3
Ephedrine/pseudoephedrine	36.5	1.6	59.8	0.1	67.7	0.1
Norephedrine	18.2	2.6	38.4	1.2	41.3	1
Azathioprine	32.8	0.9	45	3.1	51.1	0.9
Methotrexate	43.2	1.8	52.3	0.2	51.1	1.1
Ifosfamide	26.1	1	46.7	0.7	55.9	0.8
Tamoxifen	93.2	1.3	14.6	5.3	30.3	3.2
Imatinib	93.3	1.1	86.4	1.5	88.1	0.1
Capecitabine	3.4	4.0	15.3	1.2	23.8	0.7
Bicalutamide	8.8	10.2	-22.8	1.8	-16.3	1.9
Ketamine	25.8	0.2	42.7	0.7	51.5	0.4
Norketamine	21.6	1.8	34.6	0.9	45.4	0.6
Venlafaxine	27.5	2.7	56.5	1	71.5	0.2
Desmethylvenlafaxine	29.9	0.6	59.1	0.6	65.1	4.7
Fluoxetine	75.1	4.7	84.9	0.7	53.9	1.3
Norfluoxetine	79.5	3.4	86.6	0.3	54.8	0.4
Sertraline	88.1	3	88.4	0.5	57.1	0.1
Mirtazapine	26	1.1	50.2	0.7	65.2	0.2
Citalopram	52	2.8	70.1	0.4	63.7	0.8
Desmethylocitalopram	40.6	3.1	65.6	0.1	54.5	0.7
Paroxetine	93.7	0.2	86.7	0.5	70.5	1.6
Duloxetine	90.5	0.7	83.4	0.2	60.2	3.9
Amitriptyline	78.9	1.4	71.2	0.0	61.4	1.5
Nortriptyline	84.1	1.4	77.3	1.6	62.5	2.8
Norsertaline	98.5	0.3	84.2	1.8	52.9	1.5
Carbamazepine	61.4	2.5	46.6	0.5	61.4	0.4
Carbamazepine 10,11-epoxide	18.3	1.7	32.2	0.5	41.3	1.8

Compound	Surface water		Effluent		Influent	
	(%)	SD	(%)	SD	(%)	SD
10,11-Dihydro -10-hydroxycarbamazepine	20.3	1.6	33.2	0.2	37.3	0.2
Diltiazem	37.9	4.3	53.6	0.1	61.4	3.6
Verapamil	65.7	1.3	51.9	2.1	57.3	4.6
Temazepam	14.4	0.5	22.3	2	36	0.3
Oxazepam	7.2	0.8	9.6	1.5	21.8	3.0
Diazepam	14.2	1.4	20.2	2.7	25.0	3.0
Quetiapine	35	3.7	50.9	0.8	51.7	0.9
Risperidone	54.7	0.1	58.1	0.1	64.3	1.4
Donepezil	60.3	5.4	53.5	1.4	62.6	1.5
Memantine	38.1	3.1	51.3	2.4	60.8	1.3
Creatinine	69	0.4	75.2	0.1	77.2	0.6
Nicotine	1.2	3.2	5	2.2	47	1.6
Caffeine	22.2	3.2	53.4	8.5	17.1	34.5
Cotinine	33	2.1	55.6	1.4	67.3	0.9
1,7-dimethylxanthine	28.1	0.2	44.9	1.4	10.7	18.1
Morphine	30.1	3.8	47.3	2.3	69	3.3
Dihydromorphine	32.9	1.2	49.9	1.3	67.5	0.1
Normorphine	27.7	7.6	46.8	4.3	60.3	0.7
Methadone	32	3.8	48.4	0.8	45.2	0.3
EDDP	23.4	1.7	45.1	0.2	55.4	1
Codeine	25.9	3.9	53.7	0.5	59.8	1.3
Norcodeine	33.1	1.8	49.1	2.4	53.5	1.3
Dihydrocodeine	21.7	2.8	49.7	0.8	58	2.1
Tramadol	31.8	3.6	54.5	2.2	65.1	0
N-desmethyltramadol	30.2	0.3	55.6	0.5	80.4	0.2
O-desmethyltramadol	34.6	0.2	59.2	1.1	62.5	0.6
Amphetamine	23.6	1.6	36.5	1.2	50.2	0.4
Methamphetamine	23	1.9	44.1	1.6	67	0.2
MDMA	36.3	0.8	58.8	0	69.2	3.1
MDA	29.7	1	44.9	0.7	50.5	1.2
Cocaine	32.7	1.1	49.8	0.2	73.4	0.4
Benzoylcegonine	21.9	5.2	52	0.6	61.6	0.5
Anhydroecgoninemethylester	12.5	0.3	28	1.6	48.4	0.8
Cocaethylene	31.6	1.4	51.8	1	64.3	0.2
Mephedrone	24.6	2.8	44.6	0.3	64.6	0.7
MDPV	22.7	1.1	42.1	0.2	85.6	0.7
Heroin	25	1.6	45.4	0.4	61.4	2.2
6-acetylmorphine	46	2.1	65.1	1.8	66.9	1.7
Thiamethoxam	37.4	2.3	70.7	0.0	75.0	0.0
Imidacloprid	2.3	11.6	33.7	13.0	37.2	8.3
Clothianidin	11.7	1.7	28.7	7.9	25.1	1.6
Metazachlor	14.8	0.8	33.9	1.7	35.4	0.1
Terbutylazine	36.5	3.4	31.2	2.0	43.0	0.6
Methiocarb	25.4	5.6	25.2	2.0	27.7	1.7

Compound	Surface water		Effluent		Influent	
	(%)	SD	(%)	SD	(%)	SD
Dichlofluanid	59.7	0.7	94.0	6.3	57.8	6.4
Flufenacet	22.1	0.4	27.3	2.0	36.3	1.1
Oxadiazon	41.8	9.7	19.8	1.6	28.9	2.2
Chlorpyrifos	87.5	5.8	32.8	6.7	-222.8	0.5
Triallate	84.5	1.8	59.4	4.7	-23.1	2.6
Tylosin	31.4	1.1	43.4	0.4	64.9	0.5
Sulfapyridine	26.0	1.4	46.5	2.9	52.5	0.4
Sarafloxacin	66.8	1.4	62.7	3.2	64.5	5.7
Ceftiofur	18.5	28.7	37.0	7.4	10.0	1.8
Diazinon	51.0	4.6	19.6	11.2	14.0	2.6

Table S11 Surface water data (ng L⁻¹)

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Benzophenone-1</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Benzophenone-2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Benzophenone-3</i>	21.3	21.5	21	19.2	20.8	21.7	19.7	18.1	17.1	16.4	16.3	15.8	18.7	18.4
<i>Benzophenone-4</i>	476.7	410.3	482	463.4	1339.1	1303.9	878.2	951.3	403	433.7	290.7	280.6	416.1	381.8
<i>Methylparaben</i>	7	6.1	14.3	11.4	4.2	4.3	7.4	7.6	6.4	7	6.2	5.9	5.2	5.5
<i>Ethylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Propylparaben</i>	5.3	4.5	4.9	3.8	3.2	3.1	3.5	4.9	3.6	3.2	2.5	2.7	3.5	3.7
<i>Butylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bisphenol A</i>	75.6	39.6	38.1	27.4	56.5	55.2	22.4	25.9	15.2	15.2	50.4	46.9	39.8	24.7
<i>E1</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>E2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>EE2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfasalazine</i>	49	63.3	58.2	54.9	66.6	57.3	46.9	57.6	54	54.8	47.1	29.6	62.2	48.5
<i>Clarithromycin</i>	115.9	121.3	105.9	98.9	140.1	146.1	105.8	113.5	94.3	93.6	47.4	48.3	98.4	98
<i>Azithromycin</i>	2.5	1	<MQL	<MQL	5.4	5.2	<MQL	0.3	2.5	14.3	<MQL	<MQL	<MQL	<MQL
<i>Trimethoprim</i>	22.8	23.3	23.8	21.8	30.8	25.7	27.8	28.1	25.9	31.2	15.7	15.2	36.3	34.2
<i>Sulfamethoxazole</i>	30	32.2	37.1	36.4	34.6	36	35.8	37.8	34.2	34	29.5	29.6	29	29.1
<i>Triclosan</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Amoxicillin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Metronidazole</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	10.1	<MQL
<i>Sulfadiazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cefalexin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ofloxacin</i>	112.8	<MQL	56.5	<MQL	<MQL	<MQL	45.4	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ciprofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tetracycline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Danofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxytetracycline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Chloramphenicol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Penicillin G</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Penicillin V</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Erythromycin</i>	2285	2380.4	2064.4	2080.1	2435.8	2385.3	2083.6	2048.9	2062.6	2059	1849.3	1912.7	2153.7	2274.1
<i>Prulifloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Norfloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Griseofulvin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ketoconazole</i>	56.9	<MQL	44.4	<MQL	<MQL	42.8	<MQL	<MQL	<MQL	41.9	<MQL	<MQL	<MQL	<MQL
<i>Valsartan</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Irbesartan</i>	161.1	107.8	98	99.9	117.5	122.2	122.3	126.8	108.8	100.8	84.4	93.5	112.3	104.5
<i>Lisinopril</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ketoprofen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ibuprofen</i>	76	75.7	59.9	57.2	65	75.1	52.7	46.5	34.4	43.4	31.8	31.9	96	101.2
<i>Naproxen</i>	243.9	235.3	216.8	193.5	250.2	230	211.9	221.3	249.3	234.8	154.1	179.1	337	332
<i>Diclofenac</i>	82.3	88.3	85.8	77.6	93.4	93.7	98.1	99.7	87.7	96.2	69.1	81.6	92.4	88.5
<i>Acetaminophen</i>	216	214.6	171.7	142.7	223.5	220.3	204.9	173.1	169.3	167.7	127	104.7	260	306
<i>Bezafibrate</i>	114.2	119.8	95.4	96.8	113.8	119.7	111.6	106.7	93.7	96.2	75.2	69.3	124.6	113.7
<i>Atorvastatin</i>	39.3	48.3	40.7	37.7	52.3	50	48.1	46.8	38.8	44	28.9	31.9	60.9	54.1
<i>Gemfibrozil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Candesartan Cilexetil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Fexofenadine</i>	165	240.6	227.1	239.3	225.3	201.6	273.9	253.5	201.8	224.1	135.5	180	225.7	133.7
<i>Cetirizine</i>	188.6	196.3	183.9	205.2	276.8	263.9	268.2	254.1	273.1	297.8	211.1	212.8	269.1	253.1
<i>Sildenafil</i>	2.1	1.1	3.9	3.3	3.7	3.8	1.9	2.1	2.3	4.3	<MQL	0.6	0.6	0.6
<i>Metformin</i>	3750.4	3877.1	3477.7	3526.4	3974.5	3974.5	3594.6	3623.8	3429	3360.8	2727.6	2864	4179.1	4140.1
<i>Gliclazide</i>	51.3	37	36.4	38.6	45	49.9	58.1	54.9	61.2	58.9	44.8	44.7	39.5	42.4
<i>Sitagliptin</i>	170.4	157.2	112.8	153.6	153.3	152.6	150.1	143.8	149.3	142.8	115.3	123.4	128.6	142.2
<i>Pholcodine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Atenolol</i>	61.7	58.2	49.1	50.4	65.3	60.2	53.3	51.3	48	46.7	32.8	37.5	71.7	68.4
<i>Metoprolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Propranolol</i>	23.3	22.4	17.9	19	27.8	25	22.9	21.3	21.3	20.6	15	13.7	19.4	18.8
<i>Bisoprolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ramitidine</i>	135.7	161.6	136.6	140.3	170.9	176.1	174.1	187.7	137.8	179.4	69.2	85.4	150.7	169.3
<i>Cimetidine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Iopromide</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Buprenorphine</i>	<MQL	6.1	<MQL	6.5	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ephedrine/pseudoephedrine</i>	28.3	26.5	24.7	26.2	28.2	26.4	23.6	23.8	24.1	24.5	21.3	23.5	28.2	28.6
<i>Norephedrine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Azathioprine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methotrexate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ifosfamide</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tamoxifen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imatinib</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	37	<MQL	39.6	<MQL	<MQL	<MQL	<MQL
<i>Capecitabine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bicalutamide</i>	58.9	60.1	58.9	59.5	59.9	61.6	59.9	60.7	59	62	56.5	59	59.2	57
<i>Ketamine</i>	12	11	11.2	12	14.9	14.8	14.2	14.6	14.1	14.5	7.1	8.1	8.2	8.3
<i>Nor-ketamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Venlafaxine</i>	87.6	91.7	82.4	83	104.5	100.6	91.7	92.5	93.7	92.1	58.4	57.9	81.8	78.8
<i>Desmethylvenlafaxine</i>	238.9	232.2	214	226.8	266.2	250	246.4	250.4	245.4	248.6	182.7	186.6	218.2	202
<i>Fluoxetine</i>	1.4	1.2	0.8	1.1	1.2	1.2	1.2	1.2	1.2	1.1	1.3	1.4	1.3	1.4
<i>Norfluoxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sertraline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Mirtazapine</i>	5.2	5	3.4	3.9	6.2	6.8	5.3	5.1	3.7	4.3	0.8	1.2	4.5	3.8
<i>Citalopram</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Desmethylcitalopram</i>	11.9	14.7	9.2	9.1	19.9	19.4	15.9	12.5	11	10	<MQL	<MQL	6.5	4.1
<i>Paroxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Duloxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Amitriptyline</i>	16.2	12.7	11.4	12.9	12.2	<MQL	10.4	<MQL	5.5	11.1	<MQL	4.5	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Nortriptyline</i>	<MQL	2.8	<MQL	<MQL	7.9	<MQL	<MQL	<MQL	<MQL	12.6	<MQL	<MQL	<MQL	<MQL
<i>Norsertaline</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Carbamazepine</i>	162.3	164.5	152	153	184.7	185.5	174.6	177.4	168.7	168.1	142.5	144.7	136.7	131.2
<i>Carbamazepine 10,11-epoxide</i>	32.7	33.4	32.8	31.4	36.6	36	32.1	34.7	31.1	33.2	24.3	23.5	25.1	26
<i>10,11-Dihydro-10-hydroxycarbamazepine</i>	4.1	4.5	2.9	3	8.3	8.5	7.8	8.2	6.9	8.2	4.3	3.1	5.2	4.5
<i>Diltiazem</i>	0	5.6	4.6	3.6	8.8	7.7	4	4.9	4.3	5.6	<MQL	<MQL	2.7	1.1
<i>Verapamil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Temazepam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxazepam</i>	3.9	0.9	2.1	1.1	1.2	0.9	2.9	1.6	2.3	3.2	1	1.9	1.7	3.5
<i>Diazepam</i>	6.2	5.5	6	6	5.9	6	5.9	6	5.6	6	5.8	6	8.1	8.5
<i>Quetiapine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Risperidone</i>	5.1	1.3	0.3	0.2	<MQL	<MQL	<MQL	<MQL	<MQL	0.1	<MQL	<MQL	<MQL	<MQL
<i>Donepezil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Memantine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Creatinine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Nicotine</i>	31.1	30.5	24.2	20.5	31.2	32.6	23.4	21.2	23.9	22.8	64.5	65.9	54	56.8
<i>Caffeine</i>	429.2	426.3	409.7	408.2	460.4	467.2	392.8	383.7	375	373.2	309.2	297.2	620.7	617.6
<i>Cotinine</i>	39.4	38.9	36.9	34.8	38.7	42.1	34.2	36.7	35.4	32.9	28.6	31.5	48.9	50
<i>1,7-dimethylxanthine</i>	1051.7	1089.2	902.7	900.6	1038.6	1033.8	883.8	864	750.4	724.2	581.5	581.1	1235.1	1292.8
<i>Morphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dihydromorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Normorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methadone</i>	4.3	3.7	2.8	2.4	4	4.1	3.5	3.5	3.1	3	1.6	1.8	2.5	2.3
<i>EDDP</i>	10.6	10.7	9.7	9.7	12.2	12.6	10.6	11.1	10.4	10	7.4	7.4	8.5	8.5
<i>Codeine</i>	96.7	83.6	80.8	72	100.4	100.2	91.5	87.1	76.3	77.9	51.2	53	119.9	108.2
<i>Norcodeine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dihydrocodeine</i>	30.5	25	20.5	19.4	24.4	28.4	23.9	19.1	19.5	20.5	12.5	14	25.3	25.9

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Tramadol</i>	328.9	323.8	291.6	301	370.9	352.4	344.2	323	352.1	348.9	277	280.5	312.3	296.7
<i>N-desmethyltramadol</i>	251.5	259	254.2	248	291.9	298.4	265.7	272.6	270.9	273.1	218.5	222.1	222.1	226.3
<i>O-desmethyltramadol</i>	310.4	305.2	278.6	273.3	338.8	347.5	327.9	319.2	313.9	316.4	242.8	235.4	277	271.8
<i>Amphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methamphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDMA</i>	3.6	3.6	2.7	2.5	3.7	3.6	4.6	5.1	8.2	8.5	2.3	2.2	5	4.5
<i>MDA</i>	11.7	9.3	11.1	11	9.1	9.1	12.4	9.8	9.1	9.8	8.1	9	9.7	10.3
<i>Cocaine</i>	5	5.1	3.4	3.6	7.1	7.8	5.7	6.2	4.2	4.2	2.3	2.4	5.7	5.8
<i>Benzoylcegonine</i>	54.9	55.9	42.3	41.3	43.3	45	47.4	48.5	44	45.9	25.1	26.3	57.5	59.6
<i>Anhydroecgoninemethylester</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cocathylene</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Mephedrone</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDPV</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Heroin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>6-acetylmorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Thiamethoxam</i>	55.8	55	48.9	55.2	49.7	57.4	60.7	80.2	74	87.5	45.1	33.7	20.3	18.7
<i>Imidacloprid</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Clothianidin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	3.9	5
<i>Metazachlor</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Terbutylazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methiocarb</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dichlofluanid</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Flufenacet</i>	22.4	22.5	22.6	22.3	22.2	22.6	22.4	22.7	22.9	23.4	25.6	26	29.6	28.4
<i>Oxadiazon</i>	<MQL	<MQL	<MQL	<MQL	<MQL	13.6	<MQL	<MQL	<MQL	21.5	<MQL	<MQL	15.5	<MQL
<i>Chlorpyrifos</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Triallate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tylosin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfapyridine</i>	121.9	101	100.5	93.9	136.7	131.5	140.4	121	173.2	148.4	143.8	146.5	119.1	123.2

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Sarafloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ceftiofur</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	31.2	<MQL	<MQL
<i>Diazinon</i>	6	6.1	6.1	6.2	6	6	5.9	5.9	6.2	7.4	6	6	6.5	5.9

Table S12 Effluent data (ng L⁻¹)

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Benzophenone-1</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Benzophenone-2</i>	36.2	27.9	39.2	40.8	104.2	141.9	57.3	48.9	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Benzophenone-3</i>	177.7	185.4	79.4	66.4	114.6	98.3	99.6	87.9	85.4	64.5	76.1	74.5	78.6	64.8
<i>Benzophenone-4</i>	2508.6	2695.1	2816.8	2775.3	6189.3	6422.8	12191.7	11422.4	3158	3357.4	1728.7	1747	1075.4	1058
<i>Methylparaben</i>	14.8	20.2	160.2	164.4	37.8	31.4	24.4	26.5	12.5	19.4	21.6	27.3	16.6	15
<i>Ethylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Propylparaben</i>	23	20.7	54.9	56.8	24.8	32.3	30	25.5	18.1	19.2	36.4	38.2	35.8	33
<i>Butylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bisphenol A</i>	633.8	680.1	493.2	495.8	593.7	611.3	249.9	243	66.7	57.8	857.8	926.1	373.1	347.5
<i>E1</i>	4.1	7	1.4	5.5	8	7.4	5.2	6.8	6.4	5.3	10.2	7.4	4.9	8.8
<i>E2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>EE2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfasalazine</i>	281.6	307.9	219.1	233.4	241.8	237.1	216.2	269.1	241.5	239.2	167.7	222.8	251.3	251.7
<i>Clarithromycin</i>	843.2	870.1	921.1	915.4	854.6	919.7	824.6	837	673.2	708.9	674.3	703.7	563.4	558.4
<i>Azithromycin</i>	120.2	137.1	138.3	103.6	136.8	114.4	161.6	163.9	145.4	150.4	12.2	19.5	59.9	58.5
<i>Trimethoprim</i>	194.2	192.6	165.9	186.7	161.4	172.6	180.7	175.9	162.2	172.5	225	230.3	112.7	125.4
<i>Sulfamethoxazole</i>	55.5	61.8	61.4	60.9	52.3	50.9	51.2	56.2	47.5	54.1	60.3	62.7	43	42.5
<i>Triclosan</i>	279.2	233.6	132.1	154.7	344.9	282.1	267.2	234.5	241.2	157.5	201.5	197.7	175.3	190.9
<i>Amoxicillin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Metronidazole</i>	<MQL	<MQL	<MQL	27.5	40.8	42.5	45.8	47.8	36.9	<MQL	38.4	<MQL	33.7	27.5
<i>Sulfadiazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cefalexin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ofloxacin</i>	241.3	124.1	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	79.3	<MQL	63.9	<MQL
<i>Ciprofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tetracycline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Danofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxytetracycline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Chloramphenicol</i>	<MQL	<MQL	200.3	<MQL	<MQL	138.1	<MQL	<MQL	146.1	<MQL	165	<MQL	<MQL	139.9
<i>Penicillin G</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Penicillin V</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Erythromycin</i>	7087.1	6895.2	7239.8	7380.2	6739.1	6930.4	6278	6391.6	5945.6	5908.5	6752.1	6746.7	5602.2	5156.1
<i>Prulifloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Norfloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Griseofulvin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ketoconazole</i>	82.2	65.2	65.3	67.6	57.5	64	66.3	<MQL	<MQL	<MQL	56.6	60.2	64.2	56
<i>Valsartan</i>	242.9	275.1	208.6	218.2	216.1	229.6	219.6	236.4	237.5	224.4	280.4	285.1	227.7	214.5
<i>Irbesartan</i>	311.5	315.6	307.1	316	280.2	323.7	293.1	291.8	284	298.7	286.8	285.1	223	213.8
<i>Lisinopril</i>	92.3	90.8	62.3	65.6	91.9	110.9	98.2	108.6	91.6	83.7	121.5	119	97.7	80.5
<i>Ketoprofen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ibuprofen</i>	1396	1448.6	786.5	844.8	897.8	925.6	792.3	834.1	648.2	644.8	1769.9	1655.6	828	937.1
<i>Naproxen</i>	1711.4	1853.1	1193.4	1201.8	1292.6	1381.1	1429.2	1440.8	1231.9	1192.9	1661	1787.3	1480.5	1477.3
<i>Diclofenac</i>	335.3	387.8	367.4	395.3	356.8	337.5	410	397.4	344.2	348.3	341	366.8	227	223.8
<i>Acetaminophen</i>	2482.2	2457.3	1634	1634	1771.2	1721.3	1933.4	1521.8	1022.8	935.5	2232.7	2145.4	2157.9	2120.5
<i>Bezafibrate</i>	669	704.1	587.6	619.5	603.8	565.3	625.4	589.6	537.4	547.1	660.2	690.2	442.6	442.2
<i>Atorvastatin</i>	224.3	236.2	171	167.2	169.6	167.3	169.6	195.3	129.3	143.1	159.3	178	155.7	154.5
<i>Gemfibrozil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Candesartan Cilexetil</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Fexofenadine</i>	792.8	874.9	797.7	835.5	702.7	645.9	794.9	740.1	679.1	723.1	803.8	824.3	514.4	526.8
<i>Cetirizine</i>	799.1	839.7	911.2	916.6	1018.7	1055.7	1115.5	1088	1181.3	1176.1	1300.9	1336.5	957.7	1017.8
<i>Sildenafil</i>	11.2	9.8	42.5	40.1	46.6	42.8	27.4	27.2	23.1	22.7	20.3	21.1	12.9	9.9
<i>Metformin</i>	16251.2	14983.1	12922.6	13137.1	14153.3	14367.8	14889.9	14190.6	12829.4	12810.7	15067.1	13845.7	13416.8	13845.7
<i>Gliclazide</i>	88.8	79.6	77.8	81.5	94.3	91.9	93.2	98.1	83.5	83.4	80.4	84.6	72.2	70.7
<i>Sitagliptin</i>	446.5	432	471.8	399.9	422.7	403.5	499.3	475	443.1	519.1	489.1	571.6	324.6	281.4
<i>Pholcodine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Atenolol</i>	346.9	350.7	314	298.9	346.4	306.4	364.6	367.1	289	309.3	371.9	370.5	303.7	265.6
<i>Metoprolol</i>	15	10.2	11	11.4	12.4	15.6	21	13	15.7	15.5	14.3	15.9	8.7	7.1

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Propranolol</i>	103.5	103.5	103.4	105.1	108.2	110.7	121.9	109.8	93.2	98.5	111.7	113.4	111.5	98
<i>Bisoprolol</i>	14.8	14.6	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ramitidine</i>	1277.1	1219.2	1336	1241.8	1227.8	1054.3	1206.6	1161.2	1095.3	980.3	1283.9	1176.8	1053.7	968
<i>Cimetidine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Iopromide</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Buprenorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ephedrine/pseudoephedrine</i>	117.5	109.3	151.5	142.6	135.1	132.7	105.7	108.7	105.4	93.8	122.7	115.5	95.3	81.2
<i>Norephedrine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Azathioprine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methotrexate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ifosfamide</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tamoxifen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imatinib</i>	<MQL	173	188.9	<MQL	150.3	<MQL	143.8	115.1	<MQL	<MQL	<MQL	<MQL	<MQL	75.1
<i>Capecitabine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bicalutamide</i>	126	124.4	117.2	116.5	111.3	120.9	122	117.6	114.3	116.9	114.6	113.8	101.4	101.7
<i>Ketamine</i>	56.7	57.4	81.7	82	74.2	81.6	74.6	74.8	78.3	78.6	94.6	94.1	62.8	64
<i>Nor-ketamine</i>	5	4.1	5.5	3.9	2.9	3.3	3	5.6	7.4	6.8	11.5	12.5	2.9	5.3
<i>Venlafaxine</i>	611.1	622	602	624.3	563.1	574.8	597.4	597.7	567.5	584.1	271.7	276.4	325.5	315.4
<i>Desmethylvenlafaxine</i>	1040.2	1075.7	1063.6	1041.3	977.9	1065.1	1003.3	1072.3	981.5	1017.5	1030.6	1027.1	700.9	637.3
<i>Fluoxetine</i>	37.6	48	45.1	51.9	43.7	44.5	50	53.5	39	47.3	42.8	39.8	56.7	46.1
<i>Norfluoxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sertraline</i>	12.7	16.3	11	11.1	13.4	15	17.7	14.5	11.2	12.9	11.3	12.3	11	12.7
<i>Mirtazapine</i>	41.7	43.9	40	43.4	39.5	42.5	45.6	41.2	39.6	42.7	41.1	45.5	36	35.8
<i>Citalopram</i>	337.1	321.3	301.2	342.3	338.9	340	349.9	324.8	307.1	331.8	314.9	344.7	269.7	282.3
<i>Desmethylcitalopram</i>	121.6	116.2	104.7	112.5	125.7	124.4	127.6	122.2	116.4	114.3	136.3	142.8	104	106.7
<i>Paroxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Duloxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Amitriptyline</i>	61.3	65.3	63.4	55	52	45.8	52.2	46.5	50.5	47.5	57.8	62.3	46.9	39.9

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Nortriptyline</i>	35.6	22.2	35.4	32.5	37.9	30.2	35	25.5	26.6	<MQL	27.6	38.9	33.1	21.1
<i>Norsertaline</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Carbamazepine</i>	649	663.5	653.2	664.8	619.4	633.5	672.6	653.7	658	639.4	652.3	653.9	473	489
<i>Carbamazepine 10,11-epoxide</i>	157.4	161.1	164.3	171.7	163	156.1	152.4	154.5	164.9	145.9	99.9	78.2	81.7	84.2
<i>10,11-Dihydro-10-hydroxycarbamazepine</i>	46.3	45	24.3	22.6	27.6	29	47	48.2	48.4	52.3	66.3	65.2	39.5	43.5
<i>Diltiazem</i>	52.1	54.1	46.6	51.2	53.1	56	60.6	62.5	57.1	58.5	81.2	84.8	57.5	57.9
<i>Verapamil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Temazepam</i>	2.7	3.7	6.4	9.5	16.2	19.3	32.4	27.6	37.4	34.9	42.1	43.6	12.2	13.1
<i>Oxazepam</i>	16.7	22.3	16.6	17.9	17.7	16.6	17.3	18.7	17	22	19.7	24.1	7.4	8.7
<i>Diazepam</i>	14.2	14.2	11.9	<MQL	10.7	10.9	10.4	10.8	11.4	10.9	12.5	11.4	10.9	10.3
<i>Quetiapine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Risperidone</i>	4.2	2.4	1.8	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Donepezil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Memantine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Creatinine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Nicotine</i>	214.1	220.5	137.2	115.3	187.1	181.9	187	198.4	163.6	159.7	242.3	249.3	147.5	170.3
<i>Caffeine</i>	3086.9	3187	2464.7	2453.1	2803.2	2814.7	3062.9	3070.2	2453.1	2540.8	3776.5	3683.5	2852.7	2895.6
<i>Cotinine</i>	261.8	275.2	182.9	197.4	210	204.3	220.4	221.2	189	191.1	241.1	238.9	189.5	180.3
<i>1,7-dimethylxanthine</i>	8959.6	9717.4	7212.3	6632.8	7604.5	7577.8	7827.4	9378.6	7970.1	7675.9	11010.1	9744.2	7488.6	8130.5
<i>Morphine</i>	171.3	178.5	139	139	107.2	138.8	145.9	151.8	139.7	128.4	152.4	155.5	132.2	104.6
<i>Dihydromorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Normorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methadone</i>	21.8	20.4	20.5	23.1	22.5	23.9	24.8	23.8	21.2	23.4	22.5	23.1	19.8	18.8
<i>EDDP</i>	50.2	51.1	50	51.7	51.4	52.7	53.6	53.1	49.1	50.9	50.2	52.8	43.1	43.6
<i>Codeine</i>	597.2	602.1	606.2	584.7	563.6	587.6	623.7	645.3	539.2	509.1	702	594.9	485.4	449.6
<i>Norcodeine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dihydrocodeine</i>	138.2	154.6	167.1	142.4	157.3	173.3	161	170.8	155.7	141.2	182.3	181.3	104.9	126.4

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Tramadol</i>	1370.4	1413.6	1354.8	1377.3	1322.1	1400.5	1400.8	1358.5	1289.4	1316.3	1213.2	1239.8	916.5	858.9
<i>N-desmethyltramadol</i>	1026	1003.8	1054.4	1112.1	1072.1	1045.5	1059	1027	980.8	1010.8	917.8	927.2	719.1	702.2
<i>O-desmethyltramadol</i>	1474	1478.7	1388	1403.6	1370.8	1435.4	1429.8	1473.2	1309	1341.6	1183.8	1198.6	859.6	860.1
<i>Amphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methamphetamine</i>	10.4	10.4	8.5	8.6	8.3	8.8	8.4	7.8	8.8	8.9	9.7	10	9.5	9.3
<i>MDMA</i>	34.6	35.9	25.6	27.3	26.9	27.1	55.4	55	87.1	91.8	124.4	123	56	60.9
<i>MDA</i>	27.8	27.1	18.7	22.9	22.9	21.3	26	27.5	32.4	39.6	44.5	43.2	36.3	31.3
<i>Cocaine</i>	40.8	41.5	35.7	36.6	49.4	50.7	72.5	67.7	54.3	55.4	43.6	45.5	33.3	32.3
<i>Benzoylcegonine</i>	223.7	238.7	206.4	210.5	241.2	242.8	413.8	426.8	406.9	404.8	361	369.9	161.9	168.1
<i>Anhydroecgoninemethylester</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cocaine</i>	2.7	2.6	2.7	2.8	3.2	2.9	4.8	5.5	4.9	5	4.4	4	2.5	2.2
<i>Mephedrone</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDPV</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Heroin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>6-acetyl morphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Thiamethoxam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imidacloprid</i>	437.1	349.5	298.7	289.1	267.5	243	403.8	362	495.4	460.5	501.3	548.3	94	102.9
<i>Clothianidin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Metazachlor</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Terbutylazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methiocarb</i>	<MQL	<MQL	3	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dichlofluanid</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Flufenacet</i>	59.4	<MQL	60.7	<MQL	59.3	59.2	59.9	59.2	59.3	59.4	65.3	65.3	69.1	68.5
<i>Oxadiazon</i>	<MQL	34.5	28.2	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	26.6	<MQL
<i>Chlorpyrifos</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Triallate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tylosin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfapyridine</i>	580.5	540.9	797.7	786.2	506.4	498.5	587.4	497.9	662.1	580.9	722.4	565.2	344.4	393.6

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Sarafloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ceftiofur</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	264.9	123.3	154.4
<i>Diazinon</i>	11.8	13	13.6	11.6	11.9	11.5	12.1	11.6	11.7	11.7	11.5	11.8	12.6	11.8

Table S13 Influent data (ng L⁻¹)

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Benzophenone-1</i>	589.6	557.7	856.8	805.5	2712.7	4675.3	463.7	454.4	230.3	246.9	542.7	518.4	320.9	313.3
<i>Benzophenone-2</i>	906.5	812	2902.1	2825.6	6575.2	6622.6	903.3	818.9	252.2	242.8	155.7	152.8	327.9	316.5
<i>Benzophenone-3</i>	3865.3	3678.5	1422.5	1410.7	3100	2927.2	2007.7	1906.3	1268.2	1216.3	1622	1579.7	1367.7	1349.1
<i>Benzophenone-4</i>	8481.2	8567.9	15251	14741.9	45395.6	40824.6	25292	29646.3	7593	7452.2	3444.5	4625.1	5090.9	5036.7
<i>Methylparaben</i>	6747.9	6546.1	109530. 2	113301. 5	11755.2	11919.2	19411.2	19587.8	5146.1	5827.2	4830.7	4944.2	3846.9	4225.3
<i>Ethylparaben</i>	2267	2564.3	1321.7	1390.4	1382	1278.9	8239.9	7629	1127.9	869.9	834.5	662.6	365	391.6
<i>Propylparaben</i>	4891	4904.7	12625.9	11980.2	3599.6	3324.8	7693.7	7872.3	2898.9	3352.3	3462.2	3379.7	3269.8	3132.4
<i>Butylparaben</i>	162.3	135	69.6	72.7	110.2	117.9	1440.7	1335.7	146.3	149.8	194.1	166.4	311.6	310.9
<i>Bisphenol A</i>	42439.1	40048.6	28356.2	25688.6	29291.6	30071.1	3741.6	3395.1	2650.3	2979.4	34522.9	37814.1	3066	1472.4
<i>E1</i>	47.2	55.1	52.1	70.1	38.5	49.5	55.5	34	64.3	55.2	45.6	34.1	29	27.4
<i>E2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>EE2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfasalazine</i>	1183.7	1093.5	552.8	593	765.3	655.2	601.7	621.4	618.7	416.6	499.5	430.1	296.8	252.4
<i>Clarithromycin</i>	1438	1342	1370.7	1336.8	1904.9	1860.7	1522.5	1461.1	1064.6	1104.9	956.1	900.5	1119.7	1113.5
<i>Azithromycin</i>	1040	985.1	880.7	872.7	506.1	516.5	450.7	363.8	606.1	685.7	1081	1231.7	913	969.9
<i>Trimethoprim</i>	616.9	640.3	615.8	631.6	687.5	745.7	671.2	651.6	866.5	922.1	471.8	459.7	392.4	393.5
<i>Sulfamethoxazole</i>	237	243.7	139.1	139.3	139.3	134.4	122.5	141.9	134.2	128.8	60.8	54.6	43.5	54.7
<i>Triclosan</i>	2161.2	2068.2	2155.5	1862.8	2720.2	2935.1	2527.6	2335.5	2127.8	1935.1	4488.9	4040.8	1811.5	1557.3
<i>Amoxicillin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Metronidazole</i>	103.7	95.3	<MQL	65.3	123.5	126.9	<MQL	<MQL	75.6	63.8	55.9	<MQL	67.6	80.4
<i>Sulfadiazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cefalexin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ciprofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tetracycline</i>	<MQL	<MQL	<MQL	322	<MQL	<MQL	181.1	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Danofloxacin</i>	1468.1	440.3	<MQL	138.7	24.1	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxytetracycline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Chloramphenicol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Penicillin G</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Penicillin V</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Erythromycin</i>	16689.5	17722.8	15790.5	15730.6	15005.6	15297.1	13601.3	13215	11408.6	11267.2	12575.1	12710	9479.8	9788.1
<i>Prulifloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Norfloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Griseofulvin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ketoconazole</i>	140.9	189.5	123	132.8	154.9	163.6	109	112.8	185.7	203.4	87.1	80.1	73.5	78
<i>Valsartan</i>	757.8	709.6	904.6	926.5	935.3	860.6	1422.5	1425.5	1023.4	745.1	929.6	739.6	301.3	287.6
<i>Irbesartan</i>	463.8	421.6	400.9	373.4	425.3	414.1	479.9	452	366.6	378.5	365.2	336.8	235.4	216.3
<i>Lisinopril</i>	740	796.1	926.8	779.5	855.6	928	1497.2	1388.7	1081.7	1067.1	629.1	648.8	607.6	520.5
<i>Ketoprofen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ibuprofen</i>	26399.5	24884.9	23870.1	23627.8	22613	23234	16736.4	17266.5	16448.6	17357.3	17599.7	17175.6	9011.9	7542.7
<i>Naproxen</i>	19304	16694.3	16102.4	16331	17662.8	19573	17582.1	16855.7	16344.5	14528.4	9564.6	10533.1	11555.5	10708
<i>Diclofenac</i>	811.1	671.8	764.6	788.2	1081.2	1041.5	1691.3	1602.2	1092.4	898	991.1	707.3	338.4	347.6
<i>Acetaminophen</i>	497085.7	488975.2	440391.2	429786.3	470886.8	473172.4	453346.9	455736.9	423621.8	420226.2	247334.8	240282.2	236834.3	232145.7
<i>Bezafibrate</i>	2321.8	2286.9	2304.9	2369.4	2137.6	2180	2603.8	2397.9	2302.9	2296.1	1549.7	1528.1	1132.7	1123.8
<i>Atorvastatin</i>	1573.6	1852.7	1431	1556.5	1666.3	1429.8	1927.3	1854.6	1265.8	958.4	1371.6	1041.2	398.2	407.4
<i>Gemfibrozil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Candesartan Cilexetil</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Fexofenadine</i>	802.8	1343.5	891.9	1370.5	167.4	160.4	223.1	219.9	1510.6	1591	2498.4	2025.2	685.1	680.3
<i>Cetirizine</i>	1046.9	1040.6	1099	1144.5	1013.5	1118.4	1204	1159.1	1310.1	1211.9	1084	1005.4	906.1	893.9
<i>Sildenafil</i>	3.7	8	125.2	134.5	18.8	21.4	17.4	18.3	15.1	20.8	6.8	6.5	4.2	5.6
<i>Metformin</i>	145236.8	138418.9	129059.7	122752.1	133968.9	135642.8	137459.5	139755.9	127100	120363.8	87468.7	83620.9	76884.7	73251.2
<i>Gliclazide</i>	118.4	124.4	120.4	125.4	203.5	211.7	195.6	191.2	106.5	113.7	102.9	105.8	94.5	96.3
<i>Sitagliptin</i>	628.8	594.2	547.4	544.3	385.8	418.8	454	442.9	473.5	487.8	607	631.1	549.2	581.1
<i>Pholcodine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Atenolol</i>	3302.6	3260.1	2658.5	2740.7	2768.2	2948.1	3040.6	2944.3	2757.2	2822.4	2196.1	1768.6	1454.6	1482.1

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Metoprolol</i>	29.9	27.3	28.5	32	39.4	40.3	42.9	45.1	43	37	30.5	23.5	16.6	13.2
<i>Propranolol</i>	203.9	205	185.5	209.5	195.7	210.5	187.3	202.3	205.1	202	166.2	143.1	122.4	113.1
<i>Bisoprolol</i>	80	93.3	38.3	37.2	3	<MQL	14.3	6.4	43.6	38.4	<MQL	<MQL	2.7	6.5
<i>Ramitidine</i>	1625.7	1615.8	1940.2	1408.6	1436.7	1626.1	1560.9	1358.7	2033.2	2108.2	1554.5	1163.1	1085.6	1113.1
<i>Cimetidine</i>	277.6	288.4	177.1	231.2	159.9	132.1	202.9	188.1	193.1	177.9	109.4	101.3	77.7	67.8
<i>Iopromide</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Buprenorphine</i>	<MQL	<MQL	41.1	<MQL	102.9	<MQL	<MQL	<MQL	108.7	<MQL	38.1	32.3	16.2	<MQL
<i>Ephedrine/pseudoephedrine</i>	944	920.1	957.4	974.4	879.7	768.3	605.9	608.7	645.5	598.6	347.9	338.5	292.3	296.2
<i>Norephedrine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Azathioprine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methotrexate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Isofamide</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tamoxifen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imatinib</i>	109	89.3	<MQL	62.4	94.4	113.4	89.7	<MQL	63	115.3	<MQL	<MQL	<MQL	<MQL
<i>Capecitabine</i>	7.2	13.3	8.6	8.1	6	3.5	3.9	3.7	10.9	11.8	5.4	5.5	<MQL	<MQL
<i>Bicalutamide</i>	169.6	178.7	161.4	158.4	144.7	152.4	154.6	160	154.6	147.2	142.8	155.2	139.7	141.3
<i>Ketamine</i>	203.6	207.5	249.2	260.6	147.4	157.4	163.4	164.2	246.2	218.6	89.9	85	66.3	67
<i>Norketamine</i>	3	5.3	13.5	13.4	10.9	9	15.3	16	23.5	22.7	18.9	15.9	1.6	1.3
<i>Venlafaxine</i>	532.6	521.4	413.3	397.7	641.1	706.9	588.9	566.8	467.3	510.5	281.3	251	318.7	275.8
<i>Desmethylvenlafaxine</i>	870.3	867.3	821	853.5	700.9	729.6	712.9	744	854.3	770.6	753.1	700.5	439.5	457
<i>Fluoxetine</i>	54.2	53.2	46.8	45.1	57.6	54.3	57.9	46	66	69.4	52.3	52.3	44.5	44.6
<i>Norfluoxetine</i>	31.9	28.1	21.9	24.8	41.3	25.5	30.4	27.8	26.1	34.4	53.8	38.6	30.5	30.4
<i>Sertraline</i>	52.2	48.8	42.2	46.5	54.1	42.8	54.2	56	48.8	54.5	36.6	28.9	41.6	41
<i>Mirtazapine</i>	85.8	78.4	77.5	84.6	87.5	92.9	90.5	87.1	83.5	85.1	75.2	68.9	48.1	50.3
<i>Citalopram</i>	685.9	691.5	566.1	596.9	617.6	654.6	636.7	573.2	594.8	602.5	443.9	443.8	349.8	347.4
<i>Desmethylcitalopram</i>	266.7	268.3	178.6	194	279	280.7	307.5	255.9	216.6	218.1	169.9	156.1	115.3	119.6
<i>Paroxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Duloxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Amiripryline</i>	189	168.7	139.8	161.3	196.9	204.1	143.9	180.3	161.3	179.4	131.6	133	119.9	122.6
<i>Nortriptyline</i>	7.7	15	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	13.3	11.9	10.5	<MQL	2.5
<i>Norsertaline</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	157.8	<MQL	<MQL
<i>Carbamazepine</i>	442.9	439.6	640.2	653.3	580.6	604.8	685	646.7	547.8	540.4	426.2	419.3	334.7	335.8
<i>Carbamazepine 10,11-epoxide</i>	132.1	135.9	102.2	110.7	99.9	102.2	124.7	131.1	108	112.7	41.4	49.8	62.7	67.8
<i>10,11-Dihydro-10-hydroxycarbamazepine</i>	94.4	99.6	91.7	106.4	147.8	152.8	186.3	176.4	257.5	255	143	142.2	90.1	88.7
<i>Diltiazem</i>	278.3	264.2	250.4	261.9	261.5	283.4	308.2	289.5	287.2	269.2	237.4	215.5	119.6	120.1
<i>Verapamil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Temazepam</i>	2.7	1.2	2.4	5.3	12.3	13	25.6	16.8	21.9	17.3	10.6	4.1	<MQL	<MQL
<i>Oxazepam</i>	19.2	25.2	20.1	18.7	19.2	17.9	12.6	11.5	18.6	15.1	16.9	16.9	10.2	10.3
<i>Diazepam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Quetiapine</i>	132	132.8	50.8	59	36.9	40	58.6	55.8	56.2	54.6	76.7	71.7	27.8	29.4
<i>Risperidone</i>	1.2	0.8	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Donepezil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Memantine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Creatinine</i>	134004	149532	141384	157666	157304	165614	162896	164165	143975	290325.	255190.	143501	143201	159747
	9.9	6.2	3.2	9.8	2.9	7.6	3.3	7.4	4.5	8	3	1.7	2.5	2.1
<i>Nicotine</i>	4177.1	3825.2	3496.1	3530.1	3938.7	4074.9	4233.9	4177.1	3496.1	3212.3	3155.5	3076.1	1884.2	2179.4
<i>Caffeine</i>	139854.	138619.	131782.	126261.	132545.	131853	144020.	140777.	126482.	125819.	82621.5	82159.7	79398.9	78154.1
	1	3	7	3	7	7	3	7	1	5				
<i>Cotinine</i>	3561.1	3563.9	3345.8	3426.4	3401.1	3637.9	3601.3	3504.8	3227.8	3262	2507.2	2418.3	1471.8	1467.4
<i>1,7-dimethylxanthine</i>	195989	188964.	168926.	162623.	179266.	182908.	174136	175229.	164709.	173064.	94019.7	101800.	98902.2	92328.3
		1	5	2	7	8	8	8	3	8		2		
<i>Morphine</i>	1588.3	1530.5	1423.7	1483.7	1465.9	1531.7	1651.5	1494.8	1601.6	1548.2	1062.2	1128.6	772	842.6
<i>Dihydromorphine</i>	100.6	118.3	149.7	177.3	156	146.9	200.8	163.2	110.8	131.8	131.2	129.8	72.5	82
<i>Normorphine</i>	198.3	242.8	241.6	230.4	179	226.1	192.2	251.6	209.5	194.2	128.7	139.1	102.8	71.5
<i>Methadone</i>	47.6	45.7	46.6	45.9	49.5	49.5	44.4	43.9	43.8	41.6	34.1	32.4	29.4	28.1
<i>EDDP</i>	75.9	75.5	68.5	75.9	70.4	72.2	78	79	63.7	63.6	53.4	50.7	39.9	40.8

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Codeine</i>	2693.9	2721.7	2939.2	3023.7	3006.8	3053.9	2746.7	2934.8	2908.7	2961.6	2138.8	1991.6	1491.7	1372.6
<i>Norcodeine</i>	213.9	174.6	189.6	228.8	184.5	220.1	201.6	183.5	204.6	197.7	147.7	141.5	129.7	115.7
<i>Dihydrocodeine</i>	558.8	574.8	513.5	555.9	492.5	505.7	480	493.4	452.9	496.7	395.5	332.8	198.8	182
<i>Tramadol</i>	1642.3	1611.5	1434.4	1654.9	1360.6	1413.6	1409.1	1332	1222	1157	889.9	830.8	748.3	758.4
<i>N-desmethylntramadol</i>	1445.5	1423.6	1232.6	1357.8	1210.1	1257.4	1444.8	1322.1	1073.2	1092.2	720.6	679.8	403.1	408.2
<i>O-desmethylntramadol</i>	874.6	816.8	1002.8	890.6	766.1	796.6	1005.8	916.6	1030.2	1065	876.3	814.8	616.8	589.4
<i>Amphetamine</i>	557.9	534.7	475.5	502.1	517.6	593.8	561.4	511.1	631.3	669.5	349	379.3	202.2	217.3
<i>Methamphetamine</i>	20.4	19.8	9.3	11.1	11.2	11.4	14.4	15.3	10.6	11	18.7	17.8	11.7	11.6
<i>MDMA</i>	62.5	65.2	87.5	77.6	126.9	125.8	257.7	249.3	453.5	455	217.5	194.1	95.5	96.5
<i>MDA</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	46.4	44.5	<MQL	<MQL	<MQL	<MQL
<i>Cocaine</i>	748.9	716.5	704	712.1	1067.3	1131.6	1165.9	1098.8	917.8	891.8	367.5	331.5	358.1	352.3
<i>Benzoylcegonine</i>	1931.7	1975.7	1799.7	1761	2424.9	2436.8	3466.5	3491.1	3169.2	2981.7	1205.8	1073.9	687.9	691.2
<i>Anhydroecgoninemethylester</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cocacethylene</i>	17.9	16.3	16.6	17.5	38	40.6	62	55	48.7	53.7	13.2	12.6	8.9	9.4
<i>Mephedrone</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDPV</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Heroin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>6-acetylmorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Thiamethoxam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imidacloprid</i>	433.9	442	295.1	251.9	313.7	344.1	481.6	470.8	624.3	618.6	194.9	204.7	18.8	53.9
<i>Clothianidin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Metazachlor</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Terbutylazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methiocarb</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	3.1	<MQL	4.3	<MQL	<MQL	<MQL	<MQL
<i>Dichlofluanid</i>	29088.8	31767.8	30499.7	32272.8	32626.9	32962.3	25104.3	24874.5	<MQL	24288.3	25926.8	25299.9	21788.1	<MQL
<i>Flufenacet</i>	47.8	47.1	<MQL	<MQL	47.6	<MQL	<MQL	<MQL	49.2	49.1	80.9	81.7	60.5	58.9
<i>Oxadiazon</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Chlorpyrifos</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Triallate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tylosin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfapyridine</i>	1740.3	1871.8	1382.8	1487.1	1695.3	1408.8	1212.3	1358.4	1496	1464.3	1006.5	1255.4	680.2	694.1
<i>Sarafloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ceftiofur</i>	256.9	208.2	309.7	398.4	484.5	471.8	582.5	616.9	285	621.7	661.9	524.2	<MQL	<MQL
<i>Diazinon</i>	17.7	18.2	17.4	17.3	<MQL	<MQL	17	17.1	17.2	19.3	17.3	17.7	17.2	16.9

Table S14 Solid particulate matter data (ng L⁻¹)

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Benzophenone-1</i>	11.9	13.3	8.8	8.5	47.7	36.2	12.4	16.9	8.4	11.1	3.8	4.4	5.5	4
<i>Benzophenone-2</i>	14.4	14.3	33.3	42.4	85.2	72.2	15.1	14.5	1.8	1.8	2.9	3.4	3.1	3.6
<i>Benzophenone-3</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Benzophenone-4</i>	4.7	4.1	11.7	9.5	14.7	16.3	9	8.8	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methylparaben</i>	39.3	43	315.9	345.3	138.7	267.6	84.1	99.4	75.3	62.2	37.7	27.3	76	97.4
<i>Ethylparaben</i>	10.4	8.7	<MQL	<MQL	6.1	91.1	16.7	31.2	<MQL	<MQL	2.5	1.4	<MQL	<MQL
<i>Propylparaben</i>	32	32.8	20.7	27.2	23.4	188.8	36.4	53.1	11.4	14.1	7	5.9	5.9	5.6
<i>Butylparaben</i>	6.6	9.6	2	2.1	16.2	11.6	137.4	262	79.1	18.3	5	4.1	5	5.6
<i>Bisphenol A</i>	3756.1	3209.9	1944.2	2312.4	2416.5	1949.7	264.1	229.7	1211.7	1647	106.7	106.6	101.7	118.8
<i>E1</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>E2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>EE2</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sulfasalazine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Clarithromycin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Azithromycin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Trimethoprim</i>	17.4	16.3	20.5	19.5	20	18.1	40.9	16.9	10.4	12.5	10.5	10	10.9	9.9
<i>Sulfamethoxazole</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Triclosan</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Amoxicillin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Metronidazole</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfadiazine</i>	6	7.5	7.1	7.9	11.2	8.1	23.9	8.2	6.4	8.2	4.9	5.2	5.6	9.2
<i>Cefalexin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Ofloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Ciprofloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Tetracycline</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Danofloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Oxytetracycline</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Chloramphenicol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Penicillin G</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Penicillin V</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Erythromycin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Prulifloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Norfloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Griseofulvin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	1.4	<MQL	<MQL
<i>Ketoconazole</i>	114.4	117	272.1	276.9	172.2	148	169.5	153.9	149.7	122.5	188.3	118.5	202.3	196.8
<i>Valsartan</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Irbesartan</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Lisinopril</i>	14.8	8.2	10.7	19	13.7	16.5	13.7	9.7	12.5	8.7	6.5	8.2	12.3	18.4
<i>Ketoprofen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ibuprofen</i>	50.4	48	50.2	48.3	69	70.3	45.4	44.9	39.1	44.6	29.4	31.1	12.9	15.8
<i>Naproxen</i>	49	55	50	48.6	82.9	88.7	61.7	74.9	32.7	40.4	36	38.6	27.2	38.5
<i>Diclofenac</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Acetaminophen</i>	29.2	27.9	9.7	8.4	50.6	37.2	13.9	22	13	13.8	56.6	99.4	5.3	5.9
<i>Bezafibrate</i>	2.7	2.9	4.6	2.9	4.2	4	4	4.4	1.4	1.8	2.5	2.5	0.9	1
<i>Atorvastatin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Gemfibrozil</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Candesartan Cilexetil</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Fexofenadine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cetirizine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sildenafil</i>	1	0.9	9.8	9.3	1.5	1.9	2.6	2.1	1.5	1.7	1	0.8	1.1	0.9
<i>Metformin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Gliclazide</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sitagliptin</i>	15.5	9.3	10.4	19.7	13.7	15.1	13.2	12.2	8.8	8.9	8.1	6.3	8.1	7.8
<i>Pholcodine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Atenolol</i>	21	11.8	17	14	20.6	25.9	62.4	29.6	17	25.6	11.7	29	25.4	15.6
<i>Metoprolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Propranolol</i>	14.7	13.5	31.9	28.8	23.5	23.5	33.1	27.1	20.3	20.7	14.6	12.9	17.3	19.1
<i>Bisoprolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ranitidine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cimetidine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Iopromide</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Buprenorphine</i>	0.6	0.9	<MQL	0.8	<MQL	<MQL	<MQL	0.5	<MQL	<MQL	<MQL	<MQL	<MQL	0.3
<i>Ephedrine/pseudoephedrine</i>	6.7	5.4	4.5	5	6	6.6	5	4.6	2.8	2.7	3.3	3.5	4	2.7
<i>Norephedrine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Azathioprine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Methotrexate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ifosfamide</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Tamoxifen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imatinib</i>	42.9	40.1	73.8	72.5	71.4	78.3	37.3	28	31.6	31.5	42.2	33.2	39.5	44.6
<i>Capecitabine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bicalutamide</i>	9.4	10.7	8.4	8.6	9.1	<MQL	7.9	7.1	<MQL	<MQL	4	4.1	6.3	5.7
<i>Ketamine</i>	0.3	0.2	0.8	0.6	0.8	0.3	0.6	0.5	0.2	0.2	0.3	0.3	0.2	0.1
<i>Nor-ketamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Venlafaxine</i>	2.7	14.4	8	8.5	5.6	5.9	8.5	8.9	8.9	8.7	6.3	8.1	17.8	7.9
<i>Desmethylvenlafaxine</i>	2	1.4	2.4	2.7	2.4	2.3	4.1	3.7	5.5	4.3	2.6	2.7	2.2	1.8
<i>Fluoxetine</i>	170.7	69.8	46.9	48.8	31.4	28.1	27.8	30.2	27.6	28.5	21.8	21.7	29.2	29.1
<i>Norfluoxetine</i>	12.7	9.7	12.2	12.7	10.7	9.7	9.7	10.1	7	13.6	5.7	5.5	7.1	7.7
<i>Sertraline</i>	93.2	85.8	176.2	157.2	122	95.2	191.9	120.4	87.1	94.2	61.6	63.8	129.8	119.5
<i>Mirtazapine</i>	4.3	4.1	9.1	9.4	8.1	12.4	8.9	7.6	6.7	7.6	4	4	5.3	4.6
<i>Citalopram</i>	64.7	54.2	135.5	142.2	85.8	78.7	98.4	95.7	85.4	92	56.1	53.9	95.4	87.4
<i>Desmethylcitalopram</i>	21.1	18	45.9	47.6	30.4	27.8	35.1	39.2	33.3	34.7	20.9	20.5	37.3	32
<i>Paroxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	1.8	<MQL	<MQL	<MQL
<i>Duloxetine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Amitriptyline</i>	85.2	73.8	156.5	175.4	112.5	110.4	144.9	134.9	104.5	105.6	87.8	120.4	103.9	103.1

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Nortriptyline</i>	6.2	6	15.9	16.7	10.9	10	10.4	8.7	8.1	8.8	6.9	9.3	10.9	10.1
<i>Norserrtraline</i>	104.3	98.7	187	181.1	132.8	99.9	188.5	175.1	91.4	111.7	86	92.7	183.1	140.6
<i>Carbamazepine</i>	5.1	7.6	7	7.8	8.7	8.4	11.1	13.9	4.5	7.2	3.4	3.8	6.6	6.1
<i>Carbamazepine 10,11-epoxide</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>10,11-Dihydro-10-hydroxycarbamazepine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Diltiazem</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Verapamil</i>	3	2.6	8.3	10.9	4.3	3.5	4.2	4.2	5.2	5	5.8	4.6	9.5	9.4
<i>Temazepam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxazepam</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Diazepam</i>	0.7	<MQL	<MQL	<MQL	<MQL	<MQL	1.2	<MQL	<MQL	<MQL	<MQL	0.8	<MQL	<MQL
<i>Quetiapine</i>	4.9	4.5	3.3	3.1	1.2	1.4	2.6	3	7.1	3.8	1.6	1	1.6	1.2
<i>Risperidone</i>	0.1	0.2	0.1	0.1	0.1	0.1	0.1	0.1	0	<MQL	<MQL	<MQL	0.2	<MQL
<i>Donepezil</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	0.6	<MQL	<MQL	<MQL	<MQL	<MQL	1.1	<MQL
<i>Memantine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Creatinine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Nicotine</i>	29.3	31.9	29.7	30.1	80.1	59.5	122	80.8	44.1	63	35.1	46.2	17.4	21.7
<i>Caffeine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cotinine</i>	7.8	7.6	14.1	11.7	14.8	14	12.5	11.4	8.1	9.6	8.1	7.6	5.3	6.3
<i>1,7-dimethylxanthine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Morphine</i>	5.8	6.4	11.1	13.6	7.7	27.8	11.8	8.8	7.1	11.7	5.7	6	7.8	10.5
<i>Dihydromorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Normorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methadone</i>	1.4	1	2.1	2.2	2	2.5	2.2	2.1	1.3	1.4	0.9	0.9	1	1.1
<i>EDDP</i>	5	5	13.6	14.4	10.3	16.2	19.3	22	9	11	4.9	5.2	9.4	8.4
<i>Codeine</i>	19	17	28	27.7	28.1	27.2	34.8	29.4	22.9	25.7	17.2	15.6	16.3	16.7
<i>Norcodeine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dihydrocodeine</i>	2.6	2	4.6	3.6	3.8	3.6	5.1	4.5	3.6	4.9	2.4	2.6	2.4	2.4

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Tramadol</i>	3.8	3	4.6	4.2	4.3	4.6	5.9	5.2	3.1	3.4	2.8	2.8	2.6	2.2
<i>N-desmethyltramadol</i>	1.3	1	2.3	2.1	1.8	1.5	2.3	2	1.2	1.5	1.1	1	0.9	0.9
<i>O-desmethyltramadol</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Amphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methamphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDMA</i>	0.4	0.3	0.8	0.8	1	0.9	1.9	1.8	1.3	5.3	1.8	1.6	0.9	0.9
<i>MDA</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cocaine</i>	5.3	5	6.3	6.1	10.1	10.1	13.9	12.7	4.4	5	5.1	4.8	4.3	3.9
<i>Benzoylcegonine</i>	1.5	1.4	1.1	1.1	2.9	2.7	3.9	3.7	0.9	1	2.1	2.4	0.7	0.5
<i>Anhydroecgoninemethylester</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cocacethylene</i>	<MQL	<MQL	<MQL	<MQL	0.2	0.2	0.8	0.9	0.1	0.4	0.4	0.3	<MQL	<MQL
<i>Mephedrone</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDPV</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Heroin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>6-acetyl/morphine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Thiamethoxam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imidacloprid</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Clothianidin</i>	<MQL	<MQL	0.8	4.6	<MQL	<MQL	<MQL	<MQL	1.3	0.2	<MQL	<MQL	2.6	1.7
<i>Metazachlor</i>	<MQL	1.3	10.7	<MQL	<MQL	4.5	<MQL	<MQL	<MQL	0.1	<MQL	<MQL	<MQL	<MQL
<i>Terbutylazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methiocarb</i>	0.3	0.2	0.9	0.9	0.6	1.4	0.2	0.5	3.2	0.8	<MQL	0.1	0.9	1
<i>Dichlofluanid</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Flufenacet</i>	4.1	4.1	4.5	<MQL	5	5.2	<MQL	<MQL	4.1	3.6	<MQL	<MQL	2.9	2.9
<i>Oxadiazon</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	9.7	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL	8.9
<i>Chlorpyrifos</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Triallate</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Tylosin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sulfapyridine</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A

Compound	Day 1		Day 2		Day 3		Day 4		Day 5		Day 6		Day 7	
	A	B	A	B	A	B	A	B	A	B	A	B	A	B
<i>Sarafloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Ceftiofur</i>	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A	N/A
<i>Diazinon</i>	4.5	14.8	<MQL	17.9	<MQL	<MQL	<MQL	<MQL	<MQL	7.5	12.6	<MQL	<MQL	<MQL

Table S15 Digested solids data (ng g⁻¹)

Compound	Day 1		Day 2		Day 3	
	A	B	A	B	A	B
<i>Benzophenone-1</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Benzophenone-2</i>	21.3	10.9	7.9	8.9	5.2	7.0
<i>Benzophenone-3</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Benzophenone-4</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methylparaben</i>	617.3	431.8	331.8	290.7	215.2	274.3
<i>Ethylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Propylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Butylparaben</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bisphenol A</i>	4984.4	4479.9	4220.1	4200.5	4010.5	4298.3
<i>E1</i>	40.5	47.0	39.0	40.0	43.8	38.9
<i>E2</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>EE2</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sulfasalazine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Clarithromycin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Azithromycin</i>	1.7	16.6	3.1	3.6	8.6	0.3
<i>Trimethoprim</i>	14.4	15.1	8.5	7.9	10.3	10.2
<i>Sulfamethoxazole</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Triclosan</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Amoxicillin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Metronidazole</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Sulfadiazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cefalexin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Ofloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Ciprofloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Tetracycline</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Danofloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxytetracycline</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Chloramphenicol</i>	<MQL	<MQL	359.4	368.8	328.7	453.6
<i>Penicillin G</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Penicillin V</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Erythromycin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Prulifloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Norfloxacin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Griseofulvin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ketoconazole</i>	1164.2	1234.1	602.8	860.5	867.4	847.0
<i>Valsartan</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Irbesartan</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Lisinopril</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ketoprofen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ibuprofen</i>	298.3	266.0	178.4	230.9	213.5	203.8
<i>Naproxen</i>	114.1	110.0	100.6	111.9	101.4	112.2
<i>Diclofenac</i>	33.1	27.7	23.5	24.6	23.1	25.2
<i>Acetaminophen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bezafibrate</i>	7.8	8.8	7.5	7.2	7.2	6.4

Compound	Day 1		Day 2		Day 3	
	A	B	A	B	A	B
<i>Atorvastatin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Gemfibrozil</i>	1115.2	866.8	814.3	1025.1	709.4	566.6
<i>Candesartan Cilexetil</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Fexofenadine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cetirizine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sildenafil</i>	20.8	19.5	22.6	21.3	20.6	20.2
<i>Metformin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Gliclazide</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sitagliptin</i>	29.3	26.9	30.3	29.1	25.9	24.2
<i>Pholcodine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Atenolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Metoprolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Propranolol</i>	236.3	223.6	189.7	183.4	176.0	182.0
<i>Bisoprolol</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ranitidine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Cimetidine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Iopromide</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Buprenorphine</i>	18.7	4.6	13.8	12.1	13.2	16.9
<i>Ephedrine/pseudoephedrine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Norephedrine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Azathioprine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Methotrexate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ifosfamide</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Tamoxifen</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imatinib</i>	164.9	188.9	94.8	113.5	75.5	100.2
<i>Capecitabine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Bicalutamide</i>	71.5	54.1	40.4	53.7	44.5	36.1
<i>Ketamine</i>	3.6	3.4	2.3	2.5	2.8	2.1
<i>Norketamine</i>	1.1	1.2	0.6	0.7	0.7	0.7
<i>Venlafaxine</i>	134.7	180.1	108.5	108.8	105.8	120.0
<i>Desmethylvenlafaxine</i>	25.2	23.6	23.2	21.5	20.0	20.9
<i>Fluoxetine</i>	238.8	242.3	168.2	172.7	169.1	171.3
<i>Norfluoxetine</i>	133.7	123.1	77.2	73.8	68.0	72.4
<i>Sertraline</i>	677.1	653.4	527.3	514.2	506.6	512.8
<i>Mirtazapine</i>	80.5	78.0	60.7	60.6	60.4	74.3
<i>Citalopram</i>	897.2	879.1	746.2	723.6	737.7	714.0
<i>Desmethylcitalopram</i>	392.4	374.1	263.1	257.7	246.9	240.1
<i>Paroxetine</i>	5.4	6.8	0.4	0.4	<MQL	<MQL
<i>Duloxetine</i>	25.9	25.0	16.6	13.8	11.1	12.8
<i>Amitriptyline</i>	546.3	513.1	453.9	458.2	415.5	444.3
<i>Nortriptyline</i>	86.7	82.2	55.8	52.8	48.9	53.0
<i>Norsertaline</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Carbamazepine</i>	137.2	132.2	113.8	110.4	109.1	109.3
<i>Carbamazepine 10,11-epoxide</i>	N/A	N/A	N/A	N/A	N/A	N/A

Compound	Day 1		Day 2		Day 3	
	A	B	A	B	A	B
<i>10,11-Dihydro -10-hydroxycarbamazepine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Diltiazem</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Verapamil</i>	60.6	58.5	46.8	45.8	46.2	48.6
<i>Temazepam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Oxazepam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Diazepam</i>	4.9	4.0	3.9	4.5	3.9	4.3
<i>Quetiapine</i>	18.4	18.0	13.9	13.9	13.1	26.7
<i>Risperidone</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Memantine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Donepezil</i>	8.4	8.1	<MQL	<MQL	<MQL	<MQL
<i>Creatinine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Nicotine</i>	305.7	262.8	208.3	203.6	191.6	190.8
<i>Caffeine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cotinine</i>	49.6	43.0	39.3	37.1	36.7	35.0
<i>1,7-dimethylxanthine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Morphine</i>	104.5	108.5	<MQL	<MQL	22.6	17.2
<i>Dihydromorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Normorphine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methadone</i>	16.4	14.3	12.2	11.7	11.5	11.4
<i>EDDP</i>	75.8	62.4	25.5	24.2	17.3	24.4
<i>Codeine</i>	108.4	101.5	12.6	12.8	19.0	15.8
<i>Norcodeine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Dihydrocodeine</i>	30.8	29.2	23.0	22.1	22.6	21.0
<i>Tramadol</i>	50.1	37.7	31.6	32.2	29.3	26.4
<i>N-desmethyltramadol</i>	27.3	21.8	13.8	13.4	13.7	12.4
<i>O-desmethyltramadol</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Amphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methamphetamine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDMA</i>	7.8	7.1	5.0	5.7	5.6	5.4
<i>MDA</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cocaine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Benzoylcegonine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Anhydroecgoninemethylester</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Cocaethylene</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Mephedrone</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>MDPV</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Heroin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>6-acetylmorphine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Thiamethoxam</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Imidacloprid</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Clothianidin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Metazachlor</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Terbutylazine</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Methiocarb</i>	<MQL	1.7	1.6	<MQL	<MQL	<MQL

Compound	Day 1		Day 2		Day 3	
	A	B	A	B	A	B
<i>Dichlofluanid</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Flufenacet</i>	<MQL	15.4	14.7	<MQL	14.0	14.7
<i>Oxadiazon</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Chlorpyrifos</i>	<MQL	<MQL	93.8	<MQL	<MQL	<MQL
<i>Triallate</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Tylosin</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sulfapyridine</i>	N/A	N/A	N/A	N/A	N/A	N/A
<i>Sarafloxacin</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Ceftiofur</i>	<MQL	<MQL	<MQL	<MQL	<MQL	<MQL
<i>Diazinon</i>	17.4	<MQL	<MQL	10.8	<MQL	<MQL