

## IDENTIFICATION OF A BACTERIOICIN PRODUCING BY *LACTOCOCCUS LACTIS* SUBSP. *LACTIS* PD14

Hoa Thi Minh Tu<sup>1,\*</sup>, Nguyen La Anh<sup>2</sup>, Le Thanh Binh<sup>1</sup>

<sup>1</sup>*Institute of Biotechnology, VAST 18 Hoang Quoc Viet, Cau Giay, Hanoi*

<sup>2</sup>*Food Industrial Research Institute, Ministry of Industry and Trade*

\*Email: [hoathiminhthu@gmail.com](mailto:hoathiminhthu@gmail.com)

Received: 25 November 2012; Accepted for publication: 15 August 2013

### ABSTRACT

The strain *Lactococcus* PD14 isolated from cow's fresh milk was identified as *Lactococcus lactis* subsp. *lactis*. The bacteriocin produced by PD14 was isolated and purified by absorption-desorption method described and then this bacteriocin was purified by solid-phase extraction-SPE and HPLC method with column C18. The result of tricine- SDS-PAGE indicated the molecular weight of bacteriocin was about 3.5 kDa, the same as nisin.

The prenisin encoding gene was cloned and sequenced successfully for PD14. The *nis* gene of the strain PD14 was 99 % homology with *nisA* and *nisZ* while 98 % with *nisF*, 97 % with *nisQ*. The deduced amino acid sequence of prenisin was determined, which had minor difference from published sequence of *nisA* in one amino acid. At the position aminoacid -18, prenisin of PD14 had valine while the known prenisin A had phenylalanine. The matured nisin of PD14 was completely similar to the known amino acid sequence of *nisA*. Thus, the strain *Lactococcus lactis* subsp. *lactis*. PD14 produced a bacteriocin which was identified as nisin A.

**Keywords:** bacteriocin, nisin, identification, amino acid sequence, *Lactococcus lactis* subsp. *Lactis*, homology.

### 1. INTRODUCTION

Nisin is the first lantibiotic was discovered in 1928 but the structure was published by Gross and Morell in 1971, it contains unusual amino acid residues, namely dehydroalanine. Nisin is a polypeptide of size 3.5 kDa, consisted by 34 amino acids produced by *L. lactis* subsp. *lactis* [1, 2], which exhibits antimicrobial activity against a wide range of Gram-positive bacteria, especially that associated with food spoilage, such as *S. aureus*, *B. cereus*, *Clostridium* and *L. monocytogenes* [3, 4]. The unusual biosynthesis of nisin is by post-translational modification of serine and threonine to form dehydro amino acids that react with cysteine to form thioether lanthinine rings makes typical characteristics of this bacteriocin class [4]. Many studies published indicate that nisin is safe and granted GRAS (generally recognized as safe) by the FDA and allowed by the FAO and WHO to use as a food preservative [1, 5, 6]. Nisin is

worldwide used over than 50 countries for different food preservation. The major use of nisin is in dairy industry and canned vegetables and meats. The application of nisin in food preservation has been recently extending to active package field as an antimicrobial packaging material [7, 8], dehydrobutyrine, lantionine, and 3-methylanthionine [1, 4, 9, 10].

The *nis* gene cluster for nisin biosynthesis has been studied in details. It is demonstrated that the cluster is located in a conjugative transposon locates on the chromosome [1, 11, 7]. The transposon also encodes for sucrose fermentation and reduced bacteriophage sensibility. The cluster consists by 11 genes, i.e. *nisA*, *nisB*, *nisT*, *nisC*, *nisI*, *nisP*, *nisR*, *nisK*, *nisF*, *nisE*, *nisG* arranged in three multi-cistronic operon. *NisB* and *nisC* play the role for maturation of the lantibiotic, while *nisT* involved in transport across the cell membrane. *NisI* encodes an immunity protein, *nisP* is responsible for a putative serine protease involving in processing. *NisR*, *nisK* encode a putative regulatory protein and putative histidine kinase, respectively. *NisF*, *nisE*, *nisG* encode ATP-binding cassette transporter, which are together with *nisI* responsible for immunity [12].

Several lantibiotics have been discovered, among them five other natural variants of nisin are A, Z, Q, F and U [3, 13, 14]. The difference between these nisins are from one to few amino acids (aa) in their structure. Nisin A, Z, Q and F have similar activities, but differ in a few aa. Nisin A and nisin Z differ in a single one, while nisin A differs from nisin Q six aa [15, 16, 17] and differs nisin F in two aa. Nisin F differs nisin Q in four amino acids [2].

A bacteriocin was produced by strain *Lactococcus lactis* subsp. *lactis* strain PD14 isolated from cow's fresh milk in the farm at Hanoi outskirts, that was described in the previous paper [18]. The properties of bacteriocin from PD14 showed wide antibacterial range and heat resistance [19]. Thus, it has led to an interest to determine encoding gene sequence for PD14's bacteriocin and its deduced amino acid sequence in comparison with known published nisins.

## 2. MATERIAL AND METHODS

### 2.1. Microbial strains and media

The strain *Lactococcus* PD14 and *Lactobacillus plantarum* JCM1149- indicator strain were received from the collection of biomaterial technology laboratory, Institute of Biotechnology. These strains were cultured on MRS medium those composition is (g/l) peptone 10; meat extract 10; yeast extract 5; glucose 20; Tween 80 1; K<sub>2</sub>HPO<sub>4</sub> 2; CH<sub>3</sub>COONa 5; ammonium citrate 2; MgSO<sub>4</sub> .7 H<sub>2</sub>O 0.2; MnSO<sub>4</sub>.H<sub>2</sub>O 0.05; pH = 6.5 ÷ 6.8.

*Escherichia coli* DH5 $\alpha$  was cultured in Luria-Bertani broth (LB) on a shaker with speed of 220 rpm/min, overnight at 37 °C.

### 2.2. Isolation and purification of bacteriocin

The bacteriocin produced by PD14 was isolated and purified by absorption-desorption method described by Rongguang et al [20] and then this bacteriocin was purified by solid-phase extraction-SPE and HPLC method with column C18.

### 2.3. Determination the size and activity of bacteriocin directly on SDS-polyacrylamide gel

Tricine-SDS-polyacrylamid gel electrophoresis (PAGE) was conducted according to the method of Schagger with some modification for determination of size and activity of bacteriocin [21]. Bacteriocin sample producing by PD14 strain was loaded repeatedly into gel hole. The electrophoresis was stopped running when the dyes was appeared at the position less than 0.5cm from the end of the gel. The gel then cut into two parts as seen in Fig. 1, one half was dyed in coomassie, the other half was washed with 35 % ethanol, 2 % glycerol for 30 - 60 min and rinsed with sterile distilled water 15 - 30 min. The gel finally placed on a sterile pettri dish available a layer of solid agar and then poured over one layer of semi-liquid agar contained *L. plantarum* JCM1149. The dish was kept in the refrigerator for 4 h, and then transferred to 37 °C incubation for further 24 h.

#### **2.4. DNA extraction and Cloning *nis* gene**

Total DNA was prepared by Sambrook et al [22]. Total DNA was isolated by using kit of Promega Wizard® plus SV minipres DNA, A1330.

*Nis* gene was amplified by PCR using a couple of degenerate primers with the following sequences: NisF1:5→3'ATGAGTACAAAAGATTTTNAACTT and NisR1: 3→5' TTATTTNCTTACGTGAA

For the PCR, 2 µl of the total purified DNA and the two primers NisF1, NisR1 were taken. Thermal cycle was performed by Kwaadsteniet [10], denaturation at 94 °C for 4 min, 35 thermal cycles of 94 °C for 1 min, 48 °C for 30 sec, 72 °C for 7 sec and one cycle for 7 min at 72 °C.

The PCR product was purified by QiaQuick PCR Purification (Qiagen, Hilden, Germany). Ligation was conducted by ligation vector pGEM®-T, ligase (Takara) for 1 h at 16 °C. The transformation of ligated product in competent *E. coli* DH5α was performed as described Sambrook et al [22]. The isolated recombinant plasmid from transformed colonies (white colonies) was then purified and cut by restrict enzyme *EcoRI*. The restricted sample was analysed by electphoresis for *nis* gene discovery and it was used for sequence reading. The reading was carried out on Beckam coulter CIQ TM8000 by using GenomeLab Dye Terminator Cycle Sequencing (Quick Start Kit). Nucleotides was treated by SeqEd1.03 and compared by Blast of NCBI.

Sequencing gene *nis* of PD14 strain were compered with the sequences in GenBank by sofwere Clustalw.

### **3. RESULT AND DISCUSION**

#### **3.1. Electrophoresis for determination of bacterioin size**

The result of Tricine - SDS-PAGE was indicated in Fig. 1, where gel A (haft A) showed only one band with the size about 3.5 kDa. On the gel B (haft B), at the same position was appeared an inhibitory zone. It meant that *L. plantarium* JCM1149 was inhibited by the bacteriocin which was located on the gel at the position as 3.5 kDa indicated by the marker. Therefore, it can be concluded the bacteriocin produced by strains PD14 has a molecular weight of about 3.5 kDa as nisin.

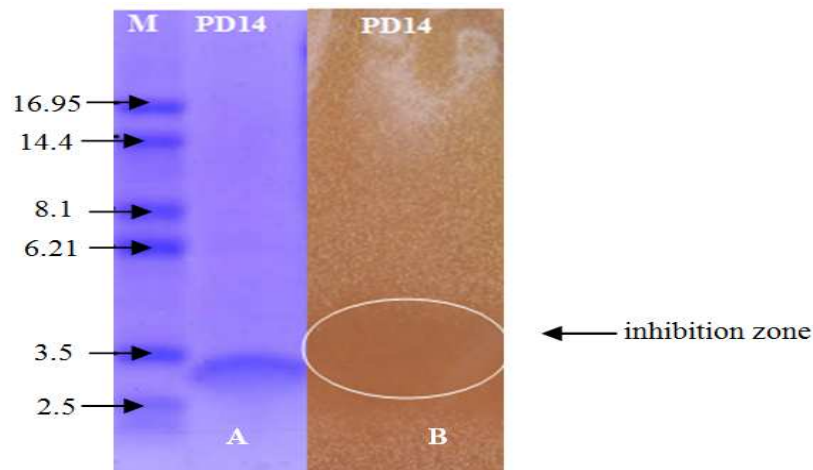


Figure 1. Electrophoresis of bacteriocin producing by *L. lactis* PD14; M: Marker MW-17S (Sigma), A: first part of the gel dyed with comasie, B: second part of the gel with inhibition zone formation (due to *L. plantarum* JCM1149 was inhibited).

### 3.2. Detection of *nis* gene in *Lactococcus lactis* subsp. *lactis* PD14

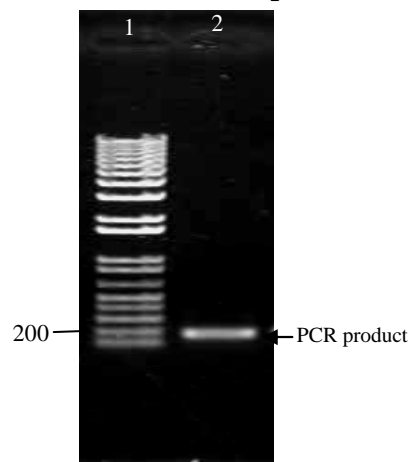


Figure 2. Electrophoresis of *nis* gene amplified on gel agarose 0.8%  
 1. Marker Fermentas GeneRuler 1 kb Ldr. Plus RTU #SM1333  
 2. PCR product.

By PCR method, the use of a couple of degenerate primers as above- mentioned designed on sequences encoded for nisin A, nisin Z, nisin Q and nisin F indicated the existing *nis* gene in the strain PD14 (Fig. 2). The PCR product was a DNA fragment of 171bp (including 69 bp of leader and 102 bp of nisin mature).

The pGEM®-T Easy Vector Systems are convenient systems for the cloning of PCR(c) products and size 3.018 kb. This Vector contains multiple restriction sites within the MCS. These restriction sites allow for the release of the insert by digestion with a single restriction enzyme. In this study, the restriction enzymes *EcoRI* was used for the release of the insert by digestion.

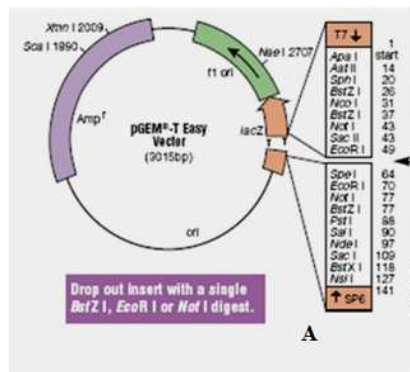


Figure 3. A - pGEM®-T Easy.

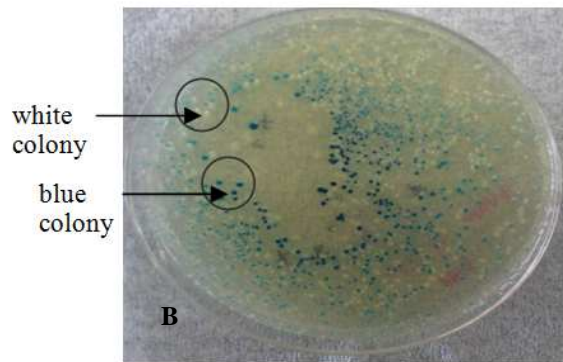


Figure 3. B – Blue and white.

The colonies of *E. coli* containing the recombinant plasmid were selected on solid LB medium containing ampicillin, X-gal and IPTG. The selected eight colonies in which six were white and two were blue colonies. Isolation of recombinant plasmid DNA was described in Promega's Protocols and Applications Guide. The result analyzed by 0.8 % agarose gel electrophoresis illustrated as in Fig. 4.

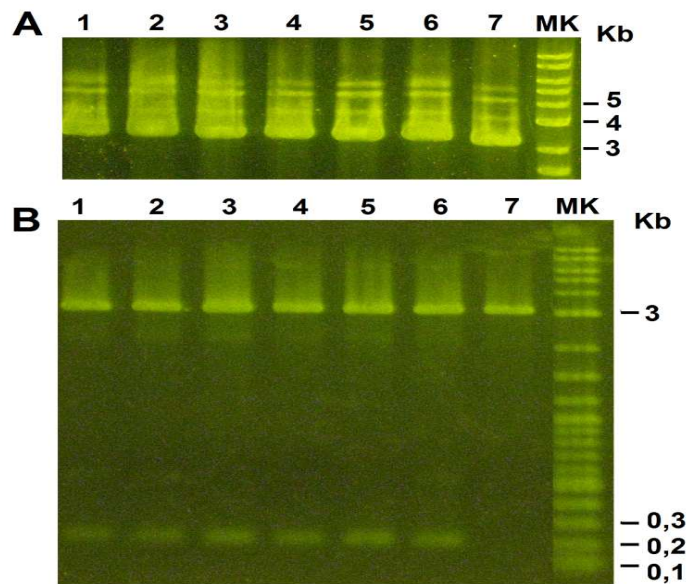


Figure 4. A. Lane: 1, 2, 3, 4, 5, 6: gene *nis* selected from white colonies; lane: 7 from blue colonies): Marker Biolabs (0.1 – 10 kb).

B. Lane: 1, 2, 3, 4, 5, 6, 7 plasmid from colonies containing *nis* gene in which (1, 2, 3, 4, 5, 6 from white colonies while 7 from blue colonies) after digestion with *EcoRI*: Marker Biolabs (0.1 – 10 kb).

To confirm the *nis* gene was inserted into vector, plasmid products were tested with restriction enzymes digestion with *EcoRI*. The results cut with restriction enzymes was checked on 0.8 % agarose gel electrophoresis (Fig. 4A).

Observation of the electrophoresis image as seen in Fig 4B showed that plasmide ADN products after digestion with *EcoRI* which appeared two bands. A band was corresponding to the



subsp. *lactis*. Therefore, it can be confirm that the bacteriocin produced by the strain PD14 is nisin A.

## REFERENCES

1. Chen H., and D. Hoover - Bacteriocins and their Food Applications, Comprehensive reviews in food science and food safety **2** (2003) 82-102.
2. Kwaadsteniet de M., Doeschate K. T. and L. M. T. Dicks - Characterization of the structure gene encoding nisin F, a new lantibiotic produced by a *Lactococcus lactis subsp. lactis* isolate from freshwater catfish (*Clarias gariepinus*), Appl. Environ. Microbiol. **74** (2008) 547-549.
3. Bailey F. J., Hurst A. - Preparation of highly active form of nisin from *Streptococcus lactis*, Canadian J. of Microbiol. **17** (1971) 61-67.
4. Twomey D., Ross R. P., Ryan M., Meaney B., Hill C. - Lantibiotics produced by lactic acid bacteria: structure, function and applications, Antonie van Leeuwenhoek **82** (2002) 165-185.
5. Delves-Broughton J. - Nisin and it's uses as a food preservative, International of Food Technology Volume 43(1990) 73-76.
6. Cleveland J., Montville T. M., Nes I. F., Chikindas M. L. - Bacteriocins: safe, natural antimicrobials for food preservation, Int. J. Food Microbio. **71** (2001) 1-20.
7. McAuliffe O., Ross R. P., Hill C. - Lantibiotics: structure, biosynthesis and mode of action, FEMS Microbiology Rev. **25** (2001) 285-308.
8. Suppakul P., J. Multz, K. Sonneveld, and S. W. Bigger - Active Packaging Technology with emphasis on antimicrobial packaging and its application, J. Food Science **68** (2) (2003) 408-420.
9. Gross E., and J. L. Morell - The structure of nisin, J. Am. Chem. Soc. **93** (1971) 4634-4645.
10. Hurst A. - Nisin. Adv. Appl. Microbiol. **27** (1981) 1285-1231.
11. Chinanoti N. Ishizaki - Nisin Z production by *Lactococcus lactis* IO-1 using xylose as a carbon source, Bioscience biotechnol, Biochem **62** (1998) 1022-1024.
12. Zendo T., Fukao M., Ueda K., Higuchi T., Nakayma J., and Sonomoto K. - Identification of the lantibiotic nisin Q, a new natural nisin variant produced by *Lactococcus lactis* 61-14 isolated from a river in Japan, Biosci. Biotechnol. Biochem. **67** (2003) 1616-1619.
13. Kuipers R. J., O. P. and de W. M. Vos - Comparison of lantibiotic gene clusters and encoded proteins, Antonie von Leeuwenhoek **69** (1996) 171-184.
14. Saris, Sea S. Beasley and Per E. J. - Nisin-producing *Lactococcus lactis* strains isolated from Human Milk, Appl. Environ. Microbiol. **8** (70) (2004) 5051-5053.
15. Buchman G. W., Banerjee S., and Hansen J. N. - Structure, expression, and evolution of a gene encoding the precursor of nisin, a small protein antibiotic, J. Biol. Chem. **263** (1988) 6260-6266.
16. Mulders J. W., I. J. Boerrigter, H. S. Rollema, R. J. Siezen, and W. M. de Vos - Identification and characterization of the lantibiotic nisin Z, a natural nisin variant, Eur. J. Biochem. **201** (1991) 581-584.

17. Zendo T., Fukao M., Ueda K., Higuchi T., Nakayama J., and Sonomoto K. - Identification of the lantibiotic nisin Q, a new natural nisin variant produced by *Lactococcus lactis* 61-14 isolated from a river in Japan, *Biosci. Biotechnol. Biochem* **67** (2003)1616–1619.
18. Hoa T. M. Tu, Nguyen Thi Da, Le Hue Anh and Le Thanh Binh - Screening Lactic Acid Bacteria Producing bacteriocins from cow's raw milk and fecal, *Tạp chí Khoa học và Công nghệ* **47** (3) (2009) 29-38.
19. Hoa Thị Minh Tú, Nguyễn La Anh, Lê Thanh Bình - phân loại chủng vi khuẩn *Lactococcus* PD14 tổng hợp bacteriocin, *Tạp chí Khoa học và Công nghệ* **50**(3) (2012) 633-641.
20. Rongguang Yang, Monty C. Johnson, and Bibek Ray - Novel method to extraction large amounts of bacteriocins from lactic acid bacteria, *Appl. and Environ. Microbiol.* **58** (1992) 3355-3359.
21. Schagger H. - Tricine-SDS-PAGE, *Nat. Protocols* **1** (1) (2006) 16-22.
22. Sambrook J., Russell D. W. - *Molecular Cloning: A Laboratory Manual*, The 3<sup>rd</sup> edition, Cold Spring Harbor Laboratory Press, Cold Spring Harbor, New York, USA, 2001.

## TÓM TẮT

### NGHIÊN CỨU ĐỊNH TÊN BACTERIOCIN DO CHỦNG *LACTOCOCCUS LACTIS* SUBSP. *LACTIS* PD14 TỔNG HỢP

Hoa Thi Minh Tu<sup>1</sup>, Nguyen La Anh<sup>2</sup>, Le Thanh Binh<sup>1</sup>

<sup>1</sup>*Viện Công nghệ sinh học, Viện HLKHCNVN, 18 Hoàng Quốc Việt, Cầu Giấy, Hà Nội*

<sup>2</sup>*Viện Công nghiệp thực phẩm, Thanh Xuân, Hà Nội*

\*Email: [hoathiminhtu@gmail.com](mailto:hoathiminhtu@gmail.com)

Chủng vi khuẩn lactic *Lactococcus* PD14 phân lập từ sữa bò tươi đã được định tên là *Lactococcus lactis* subsp. *lactis* [19]. Chủng này sinh tổng hợp bacteriocin qua phân tích điện di Tricine - SDS- PAGE cho thấy có kích thước khoảng 3,5 kDa, giống như nisin. Bacteriocin này có phổ tác dụng chủ yếu với các vi khuẩn Gram dương, trong đó đặc biệt ức chế mạnh nhóm vi khuẩn gây bệnh, gây ngộ độc thực phẩm như *B. cereus*, *S. aureus*, *L. monocytogenes*. Ngoài ra, các đặc tính chịu nhiệt, bền axit, nhạy cảm với protease rất tương đồng với nisin.

Việc nghiên cứu xác định tên một hoạt chất sinh học luôn là một ưu tiên hàng đầu, làm cơ sở cho các bước tiếp theo, đặc biệt cho những nghiên cứu ứng dụng. Với việc tách và tinh sạch thành công gene *nis* mã hóa tổng hợp nisin đã mở ra cơ hội để xác định trình tự gene *nis*, trình tự axit amin trong việc xác định tên bacteriocin này. Kết quả giải trình tự gene *nis* của chủng PD14 cho thấy 99 % tương đồng với *nisA* và *nisZ*, 98 % với *nisF* và 97 % với *nisQ*. Ngoài ra, phân tích trình tự axit amin của prenisin ở chủng PD14 chứng tỏ nó chỉ khác với nisin A bởi chỉ một axit amin ở vị trí -18 (phần đầu) và có độ tương đồng 98 % với *nisA*, 96 % với *nisZ*, 94 % với *nisF* và 87 % với *nisQ*.



Ngoài ra, chủng PD14 đã được xác định là *Lactococcus lactis* subsp. *lactis*. Vì vậy, với tất cả các kết quả này là những bằng chứng có thể khẳng định bacteriocin do chủng PD14 tổng hợp là nisin A.

*Từ khóa:* amino acid sequence bacteriocin, nisin, identification, *Lactococcus lactis* subsp. *Lactis*, homology.