

Himmelfarb Health Sciences Library, The George Washington University Health Sciences Research Commons

Pediatrics Faculty Publications

Pediatrics

2015

Percutaneous MR guided direct left atrial access to deliver large interventional devices

Toby Rogers

William Schenke

Jonathan R. Mazal

Merdim Sonmez

Ozgur Kocaturk

See next page for additional authors

Follow this and additional works at: http://hsrc.himmelfarb.gwu.edu/smhs_peds_facpubs

 Part of the [Pediatrics Commons](#)

Recommended Citation

Rogers, T., Schenke, W., Mazal, J. R., Sonmez, M., Kocaturk, O., et al. (2015). Percutaneous MR guided direct left atrial access to deliver large interventional devices. *Journal of Cardiovascular Magnetic Resonance*, 17(Suppl 1), O19.

This Journal Article is brought to you for free and open access by the Pediatrics at Health Sciences Research Commons. It has been accepted for inclusion in Pediatrics Faculty Publications by an authorized administrator of Health Sciences Research Commons. For more information, please contact hsrc@gwu.edu.

Authors

Toby Rogers, William Schenke, Jonathan R. Mazal, Merdim Sonmez, Ozgur Kocaturk, Kanishka Ratnayaka, Michael Hansen, Anthony Z. Faranesh, and Robert J. Lederman

ORAL PRESENTATION

Open Access

Percutaneous MR guided direct left atrial access to deliver large interventional devices

Toby Rogers^{1*}, William Schenke¹, Jonathan R Mazal¹, Merdim Sonmez¹, Ozgur Kocaturk^{1,2}, Kanishka Ratnayaka^{1,3}, Michael Hansen¹, Anthony Z Faranesh¹, Robert J Lederman¹

From 18th Annual SCMR Scientific Sessions
Nice, France. 4-7 February 2015

Background

Transcatheter aortic valves have benefitted from device miniaturization and in-situ assembly to reduce delivery system caliber and enable trans-vascular delivery. In contrast, investigational transcatheter mitral prostheses are bulky devices that require large caliber access ports. Trans-apical delivery is undesirable because of increased morbidity associated with rib spreading, closure site bleeding and impact on left ventricle (LV) function. Trans-septal delivery is challenging because of the acute angle required to reach the mitral valve. A 'straight shot' to the mitral valve (Fig 1A) that does not violate the LV myocardium is preferable, both in terms of device engineering and patient outcome. We hypothesized that with realtime MR guidance and by deflating a lung, it is possible to access the left atrium (LA) directly through the posterior chest wall, and close the access port using off-the-shelf nitinol devices.

Methods

LA access was obtained in 8 naïve Yorkshire swine. Animals were positioned on their right side. An active MR needle was used to access the left pleural space and insert a drain (arrow, Fig 1B) to insufflate the pleural space with CO₂. A direct trajectory to the LA was planned (Fig 1C) and the active MR needle was used to puncture through the chest wall, passing through the empty pleural space and enter the LA posteriorly (Fig 1D). Position was confirmed by pressure waveform and by injecting gadolinium. A stiff wire was introduced to the LV apex, over which an 18Fr sheath with a passive MR marker at the tip was advanced into the LA (arrow, Fig 1E-F). Sheath position and relation to the mitral valve were assessed using 3D and cine MRI.

The sheath was withdrawn and the LA puncture was closed with a nitinol closure device (arrow, Fig 1G). Animals were re-imaged 7days later.

Results

Pleural access and left lung deflation was uncomplicated in all 8 animals. Realtime MR guided LA access was successful in all (with a single pass in 7/8, but required a second pass because of rudimentary device failure in 1/8). 3D and cine MRI confirmed that the 18Fr sheath trajectory relative to the plane of the mitral valve was favorable to perform a mitral intervention (Fig 1F). The LA puncture was successfully closed with nitinol closure devices in all animals (under X-ray guidance in 6/8 and under MR guidance in 2/8). The lung was re-inflated by aspirating the pleural CO₂. There were no significant peri-procedural complications or mortality. After 7days, MRI confirmed stable position of the LA closure device (arrow, Fig 1H). Only one significant pericardial effusion was observed in the one animal in which the LA was punctured twice.

Conclusions

Percutaneous (Fig 1I) MR guided direct LA access with large sheaths is feasible in swine to achieve a straight trajectory to the mitral valve without injuring the LV myocardium. The puncture is closed using nitinol closure devices. This technique could provide a simpler and safer access route for transcatheter mitral valve interventions.

Funding

This work was supported by the Division of Intramural Research, National Heart Lung and Blood Institute, National Institutes of Health (Z01-HL005062).

¹National Heart Lung and Blood Institute, National Institutes of Health, Bethesda, MD, USA

Full list of author information is available at the end of the article

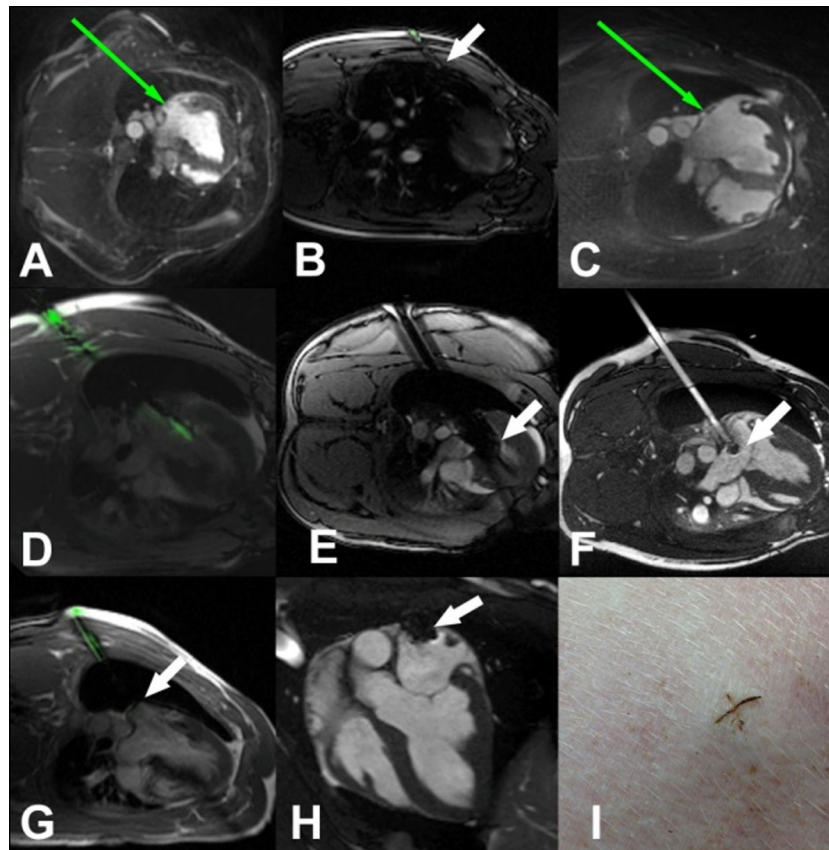


Figure 1

Authors' details

¹National Heart Lung and Blood Institute, National Institutes of Health, Bethesda, MD, USA. ²Institute of Biomedical Engineering, Bogazici University, Istanbul, Turkey. ³Department of Cardiology, Children's National Medical Center, Washington, DC, USA.

Published: 3 February 2015

doi:10.1186/1532-429X-17-S1-O19

Cite this article as: Rogers et al.: Percutaneous MR guided direct left atrial access to deliver large interventional devices. *Journal of Cardiovascular Magnetic Resonance* 2015 **17**(Suppl 1):O19.

Submit your next manuscript to BioMed Central and take full advantage of:

- Convenient online submission
- Thorough peer review
- No space constraints or color figure charges
- Immediate publication on acceptance
- Inclusion in PubMed, CAS, Scopus and Google Scholar
- Research which is freely available for redistribution

Submit your manuscript at
www.biomedcentral.com/submit

