



## A Rapid Bacteriophage DNA Extraction Method

Jakociune, Dziuginta; Moodley, Arshnee

*Published in:*  
Biology Methods and Protocols

*DOI:*  
[10.3390/mps1030027](https://doi.org/10.3390/mps1030027)

*Publication date:*  
2018

*Document version*  
Publisher's PDF, also known as Version of record

*Document license:*  
[CC BY](https://creativecommons.org/licenses/by/4.0/)

*Citation for published version (APA):*  
Jakociune, D., & Moodley, A. (2018). A Rapid Bacteriophage DNA Extraction Method. *Biology Methods and Protocols*, 1(3), [27]. <https://doi.org/10.3390/mps1030027>

Protocol

# A Rapid Bacteriophage DNA Extraction Method

Džiuginta Jakočiūnė  and Arshnee Moodley \*

Department of Veterinary and Animal Sciences, Faculty of Health and Medical Sciences,  
University of Copenhagen, 1870 Frederiksberg C, Denmark; djak@sund.ku.dk

\* Correspondence: asm@sund.ku.dk; Tel.: +45-35-33-27-11

Received: 1 June 2018; Accepted: 24 July 2018; Published: 27 July 2018



**Abstract:** Bacteriophages (phages) are intensely investigated as non-antibiotic alternatives to circumvent antibiotic resistance development as well as last resort therapeutic options against antibiotic resistant bacteria. As part of gaining a better understanding of phages and to determine if phages harbor putative virulence factors, whole genome sequencing is used, for which good quality phage DNA is needed. Traditional phage DNA extraction methods are tedious and time consuming, requiring specialized equipment e.g., an ultra-centrifuge. Here, we describe a quick and simple method (under four hours) to extract DNA from double stranded DNA (dsDNA) phages at titers above  $1.0 \times 10^{10}$  plaque-forming units (PFU)/mL. This DNA was suitable for library preparation using the Nextera XT kit and sequencing on the Illumina MiSeq platform.

**Keywords:** phage; fast; DNA extraction; spin column; genome; sequencing

## 1. Introduction

Bacteriophages (phages) are of great interest because of their ability to specifically attack bacteria. Phages are approved to be used in the food industry to control foodborne pathogens [1]. Moreover, animal and human clinical trials have demonstrated that phages can be efficiently used as alternative to antibiotics [2], including treatment of infections caused by antibiotic resistant bacteria. However, in order to ensure the safety of the phages, i.e., phages do not harbor putative virulence factors or antibiotic resistance genes that could be transferred to bacteria, whole genome analyses are performed [2], requiring extraction of phage DNA. Traditional extraction methods include precipitation of phages, ultra-centrifugation, the use of hazardous solutions, e.g., chloroform, phenol [3], and are time consuming, e.g., the protocol described in Sambrook takes two days [4]. While viral DNA extraction kits based on spin column technology are commercially available, we describe a simple and fast method using the spin column based DNeasy Blood & Tissue Kit (Qiagen, Hilden, Germany), which is used to extract DNA from blood, tissues, bacterial lysates, insects, and yeast. We chose this kit since it is our preferred in-house kit for extraction of bacterial DNA for whole genome sequencing. We applied a modified protocol to extract DNA from streptococcal phages.

## 2. Experimental Design

Required materials: Brain Heart Infusion Broth (BHI, Oxoid, Basingstoke, UK), Agarose (Merck KGaA, Darmstadt, Germany),  $\text{CaCl}_2$  (Merck KGaA, Darmstadt, Germany), SM buffer [5], DNase I (1 U/ $\mu\text{L}$ , Thermo Scientific, Roskilde, Denmark), RNase A (10 mg/mL, Thermo Scientific, Roskilde, Denmark), Proteinase K (3.0–15.0 U/mg solid, lyophilized powder, Merck KGaA, Darmstadt, Germany), 0.5 M ethylenediaminetetraacetic acid (EDTA), DNeasy Blood & Tissue Kit (Qiagen, Hilden, Germany), and 99% ethanol, GeneRuler 1 kb DNA Ladder (250–10,000 bp, Thermo Scientific, Roskilde, Denmark).

Required equipment: Heating block, benchtop centrifuge, shaker, pipettes, filter tips 1  $\mu\text{L}$ –1000  $\mu\text{L}$ , 1.5 mL and 2 mL sterile microcentrifuge tubes, 13 mL centrifuge tubes, 0.22  $\mu\text{m}$  sterile syringe

filters, NanoDrop™ (Thermo Scientific, Roskilde, Denmark), Qubit® fluorometer (Thermo Scientific, Roskilde, Denmark).

### 3. Procedure

#### 3.1. Propagation of Streptococcal Phages by the Plate Method. Time for Completion: 2–3 days

For propagation of our phages, we used a double agarose overlay plaque assay as previously described [5] with some modifications. It should be noted that this method may not apply to all phages or bacteria, hence users should use their preferred method of phage propagation.

(a) To prepare the bottom layer, BHI broth containing a final concentration of 0.8% agarose is heated to melt the agarose, and then cooled to 56 °C before adding CaCl<sub>2</sub> to a final concentration of 5 mM. Twenty-five milliliters are dispensed to sterile petri dishes. Solidified plates are stored at 4 °C until use.

(b) To prepare the top layer, BHI broth containing a final concentration of 0.4% agarose is heated to melt the agarose, and then cooled to 56 °C before adding CaCl<sub>2</sub> to a final concentration of 5 mM.

(c) Prepare SM buffer containing gelatin as described previously [5].

1. Day 1: Determine which dilution of the phage lysate that will produce confluent lysis. Mix 100 µL of this dilution with 4 mL of molten top layer containing 80 µL of an overnight broth culture of the propagation strain. This mix is then poured over the bottom layer. Once the agarose had solidified, plates are incubated upside down at 37 °C overnight.
2. Day 2: Add 5 mL of SM buffer to the plate and incubate at 4 °C overnight with shaking (50 rpm) to elute phages from the top layer.
3. Day 3. Recover the eluted phage by removing the SM buffer to 13 mL centrifuge tubes. Centrifuge at 5 min at 5000× g to remove residual bacteria and filter lysate with a 0.22 µm syringe filter.

**⏸ PAUSE STEP** Lysates can be kept at 4 °C until further use.

#### 3.2. DNA Extraction. Time for Completion: 4 Hours

##### 3.2.1. Removal of Bacterial DNA and RNA, and Digestion of Phage Capsid [4]. Time for Completion: 3 Hours

- To remove any residual bacterial DNA and RNA present in the lysate, 450 µL of the filter-sterilized lysate is incubated with 50 µL DNase I 10x buffer, 1 µL DNase I (1 U/µL), and 1 µL RNase A (10 mg/mL) for 1.5 h at 37 °C without shaking. Thereafter, 20 µL of 0.5 M EDTA (final concentration 20 mM) is added to inactivate DNase I and RNase A.
- To digest the phage protein capsid, 1.25 µL Proteinase K (20 mg/mL) is then added and incubated for 1.5 h at 56 °C without shaking.

##### 3.2.2. DNA Purification with the DNeasy Blood & Tissue Kit (Qiagen, [6]). Time for Completion: 1 Hour for 12 Samples

1. Add 200 µL of lysed phage lysate to 200 µL of AL buffer, vortex to mix thoroughly and incubate for 10 min at 70 °C. Optional: to achieve a higher DNA yield, a larger volume of the lysed phage lysate can be used and add the equivalent volume of AL buffer.
2. Add 200 µL of 99.9% EtOH (or the equivalent volume if using a larger volume), vortex to mix thoroughly.
3. Transfer the mixture to the DNeasy Mini spin column placed in a 2 mL collection tube (provided with the kit) and centrifuge for 1 min at 6000× g. If larger volumes are used, repeat this step several times. Do not load more than 750 µL at a time onto the column. Discard the flow-through and the collection tube.
4. Place the DNeasy Mini spin column in a new 2 mL collection tube (provided), add 500 µL Buffer AW1, and centrifuge for 1 min at 6000× g. Discard the flow-through and the collection tube.

5. Place the DNeasy Mini spin column in a new 2 mL collection tube (provided), add 500  $\mu$ L Buffer AW2, and centrifuge for 3 min at  $20,000\times g$  to dry the DNeasy membrane. Discard the flow-through and the collection tube. Place the DNeasy Mini spin column in a new 2 mL collection tube (provided) and centrifuge for an additional 1 min at  $20,000\times g$  to make sure there is no carryover of ethanol.
6. Place the DNeasy Mini spin column in a sterile 1.5 mL or 2 mL microcentrifuge tube (not provided in the kit), and pipette 30  $\mu$ L AE Buffer directly onto the DNeasy membrane.
7. Incubate at room temperature for 1 min, and then centrifuge for 1 min at  $6000\times g$  to elute the DNA. **OPTIONAL STEP** Repeat the elution step using the eluate from the first elution.

The concentration and quality of the DNA was assessed using the Nanodrop and Qubit according to the manufacturer's instructions. Five microliters of DNA was also run on a 1% agarose gel. Preparation of the sequencing libraries was done according to Kot et al. [7] and sequenced on the Illumina MiSeq platform [8] using MiSeq Reagent Kit v2 (500-cycles) and manufacturer's instructions.

#### 4. Expected Results

We successfully extracted DNA from 16 phages active against *Streptococcus uberis*. The DNA quality and quantity are shown in Table 1 as measured with the Nanodrop, Qubit, and by gel electrophoresis (Figure 1). We have also included the results for the *Salmonella* (S112) and *Escherichia coli* (EC11 and EC80) phage DNA extracted using our method for comparison.


**Table 1.** Overview of DNA quality and quantity compared to starting phage concentration.

| Phage ID | Starting PFU/mL      | Nanodrop (ng/ $\mu$ L) | Qubit (ng/ $\mu$ L) | 260/280 <sup>a</sup> |
|----------|----------------------|------------------------|---------------------|----------------------|
| SU1      | $2.1 \times 10^{11}$ | 225.67                 | 216                 | 1.86                 |
| SU2      | $9.0 \times 10^{10}$ | 81.15                  | 80.8                | 1.82                 |
| SU3      | $9.0 \times 10^{10}$ | 66.87                  | 56.4                | 1.88                 |
| SU4      | $6.0 \times 10^{10}$ | 26.85                  | 17.5                | 1.65                 |
| SU5      | $1.6 \times 10^{11}$ | 268.01                 | 224                 | 1.83                 |
| SU6      | $3.0 \times 10^{11}$ | 69.41                  | 62.4                | 1.82                 |
| SU7      | $6.0 \times 10^{10}$ | 95.73                  | 89.6                | 1.86                 |
| SU8      | $1.0 \times 10^{11}$ | 105.76                 | 78.4                | 1.80                 |
| SU9      | $1.1 \times 10^{11}$ | 59.35                  | 55.2                | 1.84                 |
| SU10     | $2.0 \times 10^{11}$ | 287.35                 | 216                 | 1.84                 |
| SU11     | $1.0 \times 10^{10}$ | 12.67                  | 8.56                | 1.81                 |
| SU12     | $1.4 \times 10^{10}$ | 4.56                   | 1.22                | 1.57                 |
| SU13     | $2.0 \times 10^{11}$ | 85.73                  | 67.6                | 1.82                 |
| SU14     | $1.9 \times 10^{11}$ | 93.16                  | 80.8                | 1.84                 |
| SU15     | $1.3 \times 10^{11}$ | 115.78                 | 91.6                | 1.83                 |
| SU16     | $1.0 \times 10^{12}$ | 73.95                  | 59.6                | 1.88                 |
| *S112    | $1.0 \times 10^{10}$ | 57.60                  | 52                  | 1.60                 |
| *EC11    | $2.8 \times 10^{10}$ | 30.91                  | 35                  | 1.80                 |
| *EC80    | $2.0 \times 10^{11}$ | 33.55                  | 24.8                | 1.84                 |

<sup>a</sup> Ratio =  $\sim 1.8$  is accepted as pure DNA. Lower values might be related to the presence of residues of phenol, guanidine, or other reagents used in DNA extraction or an indication of a very low concentration of DNA. \* *Salmonella* and *Escherichia coli* phages.

The DNA concentration varied from 4.56–287.35 ng/ $\mu$ L. There was no linear correlation observed between the starting concentration of phages and the quantity of DNA extracted. Based on our results, we observed that the optimal concentration of our phages for successful extraction of DNA was above  $1.0 \times 10^{10}$  PFU/mL, as no measurable DNA was obtained when titres lower than this were used (e.g., SU11, starting titer =  $6.0 \times 10^9$  PFU/mL, nanodrop concentration =  $-8.59$  ng/ $\mu$ L). The purity of DNA was determined by the A260/280 ratio. Using this measurement, 16/19 DNA samples were

of high quality. The integrity of the extracted DNA was determined by gel electrophoresis (Figure 1). An intact band suggests that the DNA is not contaminated with RNA nor degraded.

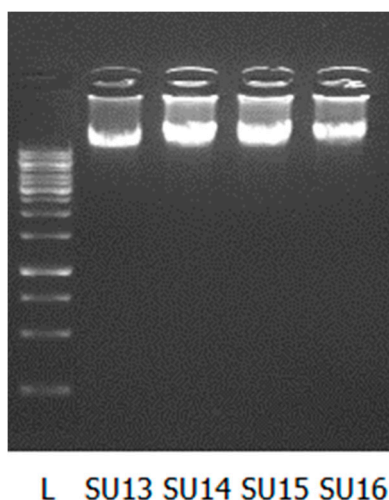

**Figure 1.** Five  $\mu\text{L}$  of extracted DNA was loaded on 1% agarose gel. L—GeneRuler 1 kb DNA Ladder (250–10,000 bp, Thermo Scientific, Roskilde, Denmark); S13–S16 phage DNA.

Despite the low quality scores for a few samples, we were able to prepare libraries and sequence on the Illumina Miseq platform all phage DNA samples. De novo sequence assembly of the trimmed raw reads was performed in the CLC Genomics Workbench v.7 (Qiagen Hilden, Germany) and for each phage,  $\geq 99.85\%$  of reads were assembled into one contig representing the phage genome (Table 2).

**Table 2.** Summary of read and assembly quality.

| Phage ID | Number of Contigs >1000 bp | Max Contig Length | Min Contig Length | N50     | Reads Assembled to Max Contig (%) | Number of Paired Reads | Reads with Phred Score Above 30 (%) |
|----------|----------------------------|-------------------|-------------------|---------|-----------------------------------|------------------------|-------------------------------------|
| SU1      | 1                          | 34,770            | 209               | 34,770  | 99.85                             | 575,584                | 96.93                               |
| SU2      | 3                          | 34,819            | 200               | 34,819  | 99.96                             | 547,242                | 97.58                               |
| SU3      | 1                          | 34,808            | 201               | 34,808  | 99.86                             | 560,430                | 96.76                               |
| SU4      | 1                          | 34,819            | 204               | 34,819  | 99.94                             | 422,742                | 95.65                               |
| SU5      | 1                          | 34,800            | 98                | 34,800  | 99.86                             | 238,880                | 96.95                               |
| SU6      | 1                          | 34,799            | 205               | 34,799  | 99.95                             | 483,526                | 96.61                               |
| SU7      | 1                          | 34,784            | 218               | 34,784  | 99.98                             | 391,000                | 93.43                               |
| SU8      | 1                          | 34,802            | 206               | 34,802  | 99.92                             | 354,100                | 96.86                               |
| SU9      | 1                          | 34,771            | 204               | 34,771  | 99.85                             | 245,804                | 96.83                               |
| SU10     | 1                          | 34,873            | 204               | 34,873  | 99.99                             | 266,872                | 95.73                               |
| SU11     | 2                          | 37,919            | 214               | 37,919  | 99.94                             | 247,626                | 95.79                               |
| SU12     | 1                          | 34,507            | 221               | 34,507  | 99.99                             | 337,214                | 95.70                               |
| SU13     | 1                          | 34,804            | 201               | 34,804  | 99.93                             | 304,082                | 92.69                               |
| SU14     | 1                          | 35,029            | 200               | 35,029  | 99.91                             | 294,802                | 96.27                               |
| SU15     | 1                          | 34,974            | 205               | 34,974  | 99.98                             | 371,720                | 94.62                               |
| SU16     | 1                          | 34,819            | 201               | 34,819  | 99.99                             | 552,896                | 94.31                               |
| S112     | 1                          | 52,664            | 210               | 52,664  | 99.97                             | 100,884                | 93.28                               |
| EC11     | 1                          | 169,478           | 213               | 169,478 | 99.98                             | 206,448                | 97.31                               |
| EC80     | 1                          | 52,706            | 202               | 52,706  | 99.96                             | 5902                   | 95.20                               |

## 5. Conclusions

The DNeasy Blood and Tissue kit from Qiagen is commonly used in research laboratories for extraction of bacterial or eukaryotic DNA. Here, we demonstrate that with some modifications it can also be used for purification of phage DNA that can be sequenced on the Illumina MiSeq platform.

Limitations to the technique. The quality and concentration obtained in this study was sufficient for sequencing on the Illumina MiSeq platform. However, this DNA might not be of sufficient quality for use on other sequencing platforms such as non-amplification, long read technology e.g.,

PacBio [9] or Oxford Nanopore [10]. Furthermore, we only extracted DNA from double stranded DNA (dsDNA) phages and it is unknown if this method would work on single-stranded DNA (ssDNA) or RNA viruses.

Method tried on other types of phages. Our streptococcal phages belonging to the Podoviridae-like family were used in this study. We have also used this method on phages active against *E. coli* and *Salmonella* (Tables 1 and 2) belonging to the Myoviridae and Siphoviridae families.

**Author Contributions:** D.J. and A.M. conceived and designed the experiments; D.J. performed the experiments; D.J. and A.M. analyzed the data and wrote the paper.

**Funding:** This project was funded by Danish Council for Independent Research-Technology and Production Sciences (DFF-FTP).

**Acknowledgments:** We would like to thank Tatjana P. Kristensen (University of Copenhagen) for technical assistance, Michela Gambino and Lone Brøndsted (University of Copenhagen) for allowing us to use S112 DNA information.

**Conflicts of Interest:** The authors declare no conflicts of interest.

## References

1. Sulakvelidze, A. Safety by nature: Potential bacteriophage applications. *Microbe Mag.* **2011**, *6*, 122–126. [CrossRef]
2. Vandenhuevel, D.; Lavigne, R.; Ussow, H.B. Bacteriophage therapy: Advances in formulation strategies and human clinical trials. *Annu. Rev. Virol.* **2015**, *2*, 599–618. [CrossRef] [PubMed]
3. Pickard, D.J.J. Preparation of Bacteriophage Lysates and Pure DNA. In *Bacteriophages*; Clokie, M.R.J., Kropinski, A.M., Eds.; Methods in Molecular Biology; Humana Press: New York, NY, USA, 2009; Volume 2, pp. 3–10, ISBN 978-1-58829-682-5.
4. Sambrook, J.; Russell, D.W. *Molecular Cloning: A Laboratory Manual*; Cold Spring Harbor Laboratory Press: Cold Spring Harbor, NY, USA, 2001; ISBN 9780879695774.
5. Kropinski, A.M.; Mazzocco, A.; Waddell, T.E.; Lingohr, E.; Johnson, R.P. Enumeration of Bacteriophages by Double Agar Overlay Plaque Assay. In *Bacteriophages*; Clokie, M.R.J., Kropinski, A.M., Eds.; Humana Press: New York, NY, USA, 2009; Volume 1, pp. 69–76, ISBN 978-1-58829-682-5.
6. DNeasy®Blood & Tissue Handbook. July 2006. Available online: [http://www.webcitation.org/query?url=http%3A%2F%2Fdiagnostics1.com%2FMANUAL%2FGeneral\\_Qiagen.pdf&date=2018-06-01](http://www.webcitation.org/query?url=http%3A%2F%2Fdiagnostics1.com%2FMANUAL%2FGeneral_Qiagen.pdf&date=2018-06-01) (accessed on 9 September 2017).
7. Kot, W.; Vogensen, F.K.; Sørensen, S.J.; Hansen, L.H. DPS—A rapid method for genome sequencing of DNA-containing bacteriophages directly from a single plaque. *J. Virol. Methods* **2014**, *196*, 152–156. [CrossRef] [PubMed]
8. MiSeq Support—Documentation Literature. Available online: [http://emea.support.illumina.com/sequencing/sequencing\\_instruments/miseq/documentation.html](http://emea.support.illumina.com/sequencing/sequencing_instruments/miseq/documentation.html) (accessed on 26 July 2018).
9. PacBio. Available online: [www.pacb.com](http://www.pacb.com) (accessed on 24 July 2018).
10. Oxford Nanopore. Available online: <https://nanoporetech.com> (accessed on 24 July 2018).



© 2018 by the authors. Licensee MDPI, Basel, Switzerland. This article is an open access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (<http://creativecommons.org/licenses/by/4.0/>).