

RESEARCH ARTICLE

Open Access

Hair-based rapid analyses for multiple drugs in forensics and doping: application of dynamic multiple reaction monitoring with LC-MS/MS

Iltaf Shah¹, Andrea Petroczi¹, Martina Uvacsek², Márta Ránky³ and Declan P Naughton^{1*}**Abstract**

Background: Considerable efforts are being extended to develop more effective methods to detect drugs in forensic science for applications such as preventing doping in sport. The aim of this study was to develop a sensitive and accurate method for analytes of forensic and toxicological nature in human hair at sub-pg levels.

Results: The hair test covers a range of different classes of drugs and metabolites of forensic and toxicological nature including selected anabolic steroids, cocaine, amphetamines, cannabinoids, opiates, bronchodilators, phencyclidine and ketamine. For extraction purposes, the hair samples were decontaminated using dichloromethane, ground and treated with 1 M sodium hydroxide and neutralised with hydrochloric acid and phosphate buffer and the homogenate was later extracted with hexane using liquid-liquid extraction (LLE). Following extraction from hair samples, drug-screening employed liquid chromatography coupled to tandem mass spectrometric (LC-MS/MS) analysis using dynamic multiple reaction monitoring (DYN-MRM) method using proprietary software. The screening method (for > 200 drugs/metabolites) was calibrated with a tailored drug mixture and was validated for 20 selected drugs for this study. Using standard additions to hair sample extracts, validation was in line with FDA guidance. A Zorbax Eclipse plus C18 (2.1 mm internal diameter × 100 mm length × 1.8 μm particle size) column was used for analysis. Total instrument run time was 8 minutes with no noted matrix interferences. The LOD of compounds ranged between 0.05-0.5 pg/mg of hair. 233 human hair samples were screened using this new method and samples were confirmed positive for 20 different drugs, mainly steroids and drugs of abuse.

Conclusions: This is the first report of the application of this proprietary system to investigate the presence of drugs in human hair samples. The method is selective, sensitive and robust for the screening and confirmation of multiple drugs in a single analysis and has potential as a very useful tool for the analysis of large array of controlled substances and drugs of abuse.

Keywords: Dynamic MRM, LC-MS/MS, Human hair, Controlled drugs, Doping, Forensics, Toxicology, Drugs of abuse

Introduction

Hair testing is a convenient, tamper resistant and non-invasive technique for the analysis of many controlled drugs and drugs of abuse as compared to blood tests and urinalysis [1]. Hair testing can be used to screen for the parent drug and for metabolites and could be used to complement urinalysis [1-5]. In addition to urine tests, for the past three decades hair analysis has been employed to detect chronic social drug use and in the

fight against doping in sport. Longer term histories of drug use can be detected as hair grows by approximately 1 cm per month [6-8]. Hair analysis provides a wider window of detection thus hair samples of a few centimetre lengths will provide more accurate information of drug use than blood and urine and it has good potential for out of competition doping [9-12]. Most controlled drugs and drugs of abuse are integrated into the hair matrix in a number of ways: (i) an endogenous-exogenous pathway; transfer or absorption of drug molecules into hair duct in the form of sweat and sebum from transdermal secretion, or (ii) by an endogenous pathway; the drug molecules diffuse into growing hair from circulatory system using

* Correspondence: D.Naughton@kingston.ac.uk

¹School of Life Sciences, Kingston University, Kingston-upon-Thames, Surrey, UK
Full list of author information is available at the end of the article

passive transport [13-17]. Dosage, bioavailability, physiochemical properties and pharmacokinetics also affect the process of drug incorporation into the hair matrix. The concentration of drug detected in hair is also influenced by an individual's metabolic pathway, cosmetic treatments and hair pigmentation. Hair also favours the incorporation of undissociated basic drugs due to carrying no net charge due to an isoelectric pH of about 6. Furthermore, due to the less polar nature of the parent drugs, its incorporation into the keratin matrix is favoured in contrast to their metabolites. While in urine it is the other way around and relatively low concentrations of parent drugs are excreted in urine as compared to their metabolites [18,19].

Doping in sport both in and out of competition is a persistent problem and therefore it is highly desirable to have a high throughput multi-analyte mass spectrometry method for the detection of controlled drugs and drugs of abuse in hair. Gas chromatography–mass spectrometry (GC-MS), requires additional derivatisation and lacks the required sensitivity to simultaneously detect large numbers of drugs in hair samples. Liquid chromatography–mass spectrometry (LC-MS/MS) when used in Multiple Reaction Monitoring (MRM) mode has been a powerful tool for detecting and confirming the presence of drugs in complex biological matrices. Recently, researchers have reported the hybrid triple-quadrupole mass spectrometer (QTrap) to screen and confirm 300, 100 and 88 drugs in human blood samples, respectively [20-26]. While others reported a 3200 Q Trap(R) LC-MS/MS system for analysis of 700 and 301 drugs in serum and urine samples [27,28]. A set of over 500 negative-ion MS–MS spectra was collected from three libraries applied in screening and systematic toxicological analysis is also reported [29]. Moreover, the mass spectra characteristics of more than 2,500 illegal drugs and metabolites have also been measured in urine and plasma using hybrid quadrupole time-of-flight mass spectrometry (LC-QTOF-MS) [30,31]. All of the above techniques are limited to human blood or urine samples. Furthermore, the LOD for most of these drugs were not suitable for doping control purpose as the methods of detections were not sensitive enough. However, there is an analytical method of metabolomic approach used for hair analysis using time of flight detector (TOF) and High resolution mass spectrometry (HRMS), which is limited to metabolite analysis and not the actual compounds and hence not widely used [32,33]. An UPLC-TOF-MS method for simultaneous screening and quantification of 52 drugs in hair was also developed and validated but the analysis of drugs were limited and only 15 autopsy hair samples were tested using this method [34].

The aim of this paper was (i) to develop a hair-based multi-drug/metabolite assay using LC-MS/MS with a dynamic multiple reaction monitoring (DYN-MRM)

[35,36] method using a proprietary software, for both screening and validated confirmatory analyses, and (ii) to apply the method to hair samples for drugs and metabolites relevant to social drug use and doping in sport.

Materials and methods

Reagents and chemicals

LC/MS grade acetonitrile and deionised water together with formic acid, ammonium formate, dichloromethane, sodium hydroxide, hydrochloric acid, pentane, sodium hydrogen phosphate heptahydrate, and sodium phosphate monobasic dihydrate were obtained from Sigma-Aldrich (Poole-Dorset, UK). All drugs including nandrolone, stanozolol, testosterone, boldenone, cocaine, benzoylecgonine, cocaethylene, amphetamine, n-desmethylephedrine, ephedrine, methamphetamine, 3,4-methylenedioxy-N-methylamphetamine (MDMA), 3,4-methylenedioxyamphetamine (MDA), tetrahydrocannabinol (THC), 11-hydroxy-THC (11-OH-THC), THC-carboxylic acid-(THC-COOH), ketamine, norketamine, clenbuterol, propranolol, terbutaline, salbutamol, morphine, codeine and phencyclidine (PCP) together with internal standards stanozolol-d₃, cocaine-d₃, MDMA-d₅, THC-d₃, codeine-d₃, ketamine-d₄ and PCP-d₅ were all purchased from LGC standards (Teddington, UK).

Calibrant and stock solutions

Stock solutions of individual analytes and internal standards were prepared by dissolving or diluting it in methanol obtaining a concentration of 1 mg/mL. Individual stock solutions were diluted to prepare a standard mixture solution by mixing individual drug in methanol arriving at a concentration of 1 µg/mL. Moreover, working solutions of the individual standards were also prepared in methanol at various concentrations. The concentration of the working solution of internal standards mixture used was 10 ng/mL. Calibrants and quality control samples were prepared using blank hair. All standard stock solutions of forensic and toxicological interest were prepared in methanol at a concentration of 1 mg/mL. Furthermore, a bespoke working solution test mixture was prepared from dilution of stock solutions in-house at a concentration of 10 ng/mL. The test mix prepared, consisted of a representative range of forensic analyte classes of 27 components which is used to check the sensitivity and accuracy of the new dynamic MRM assay. All solutions were stored in amber plastic tubes at –80°C until further analysis.

Sample preparation and analysis

Hair samples (50 mg) were collected from 233 volunteers as part of a research project investigating the discriminatory power of drug- and doping-related explicit and implicit social cognitive measurements in behavioural context indexed on the presence or absence of the target

drugs in hair. On the average the hair samples were 3 cm in length (circa 120 mg in weight). The hair samples were cut directly from the vertex posterior of the head where it was possible. The hair samples were stored in labelled, sealed paper envelopes. The study was approved by the Kingston University Research Ethics Committee. Blank hair samples were obtained from healthy volunteers with no previous drugs use.

Extraction of hair samples and screening using dynamic MRM-LC-MS/MS mode

The Agilent Technologies forensic and toxicology dynamic MRM-LC-MS/MS kit was used to screen hair samples for all controlled substances and drugs of abuse of interest by employing a single point calibration curve. For extraction, the hair samples were decontaminated using dichloromethane then milled in a mini ball mill (Fritsch and Gerhardt UK Ltd); 20 mg hair was weighed out and incubated with 1 M sodium hydroxide at 95°C for 5 minutes. The mixture was neutralised with 1 M hydrochloric acid followed by addition of 2 mL of 0.2 M phosphate buffer pH 7 [6]. The sample was further extracted using hexane and 2 µL injected on to the LC-MS/MS system. The kit was found to be a good and convenient starting point for the screening of a large array of toxicological and forensic analytes. As compared to the multiple reactions monitoring (MRM), dynamic MRM has an associated delta retention time window (Del RT) where it dynamically turns on and off without affecting the total cycle time. But some time compromises between, cycle time, dwell time and the total number of MRM transitions are required in a large number of analytes screening assays. For screening purposes we have used mostly 4 MRM transitions for each compound.

Confirmation of extracted hair samples

For quantitative confirmatory analysis of positive hair samples a 6–8 point calibration curve was used along with quality control samples. A further 2 µL of extracted hair sample was injected for quantitative analysis using DYN-MRM-LC-MS/MS. Only the two most abundant MRM transitions were used for confirmatory analysis. Of the two, the most sensitive transition was used as quantifier ions and the other transition as qualifier ions, the ratios of which are indicative of the analyte of interest.

Instrumentation

The LC–MS/MS system comprised of a 1260 infinity LC system (Agilent Technologies UK) coupled to a 6430 triple quadrupole mass spectrometer (Agilent Technologies UK). The LC system consisted of a 1290 infinity thermo-stated autosampler, degasser, binary pump and column heater. Agilent Masshunter Forensics and Toxicology dynamic MRM software kit was used in combination with

Masshunter workstation software. An electrospray ionisation (ESI) source was employed for samples ionisation. The chromatographic separation was achieved using a Zorbax Eclipse plus C18 (2.1 mm × 100 mm, 1.8 µm)

Table 1 Chromatographic and spectrometric characteristics of analytes selected for validation for confirmatory analyses

Compound name	Precursor ion	Product ion	Fragmentor	Collision energy
Nandrolone	275.1	257	105	5
		109.1		21
Stanozolol	329.1	175.2	105	10
		107.1		34
Boldenone	287.2	135.1	105	5
		121.1		20
Cocaine	304.2	182.1	138	17
		77		61
Benzoylcegonine	290.1	168.1	118	17
		77		60
Amphetamine	136.1	119.1	66	5
		91		17
Methamphetamine	150.1	91	92	17
		65.1		45
Ephedrine	166.1	148.1	81	9
		115		25
MDMA	194.1	163	97	9
		105		25
THC	315.2	123.3	150	30
		193.2		20
THC-COOH	345.2	299.2	140	18
		327.2		18
11-Hydroxy THC	331.2	313.2	112	9
		193.1		21
Morphine	286.2	165.1	158	41
		152		60
Codeine	300.2	165.1	158	45
		58.1		29
Clenbuterol	277	203.1	100	15
		240		15
Propranolol	260	115.9	110	22
Salbutamol	240	73.9	90	22
		152.1		20
Terbutaline	226.2	125.1	100	20
		152.1		20
PCP	244.2	91	86	41
		86.1		9
Ketamine	238.1	220.1	105	11
		125		11

column. Column was heated to 65°C for good reproducibility. The analytical column was connected in tandem with a 0.2 µm inline filter to prevent it from blocking. Mobile phase consisted of a Solvent A [10 mM ammonium formate/0.02 M formic acid in water] and solvent B [0.02 M formic acid in acetonitrile].

The flow rate through the column was set at 0.4 mL/min. The gradient flow composition is described as follows. The mobile phase linear gradient was run with 90% solvent A for 0.5 minutes and 10% solvent B, decreasing to 50% solvent A in 3 minutes, decreasing further to 5% solvent A and 95% solvent B, and held from 4 to 6 minutes and then returned to 90% of solvent A until 7 minutes, and stabilised until 8 minutes before next injection. The bypass configurations were set up for the mixer and damper when using this kit with our Agilent HPLC 1260 binary pump. This was to convert the pump to low delay volume mode to achieve better chromatography. The mass spectrometer was operated in positive ionization polarity at a spray capillary voltage of 3000 V. Sheath gas temperature and flow was 350°C, 12 L/min, nozzle voltage was 450 V, drying gas temperature and flow 300°C, 8 L/min, nebulizer gas pressure was 25 psi and fragmentor voltage was set to 150 V, Delta retention time was 1 for all ions in dynamic MRM mode. The total

run time was 8 min together with 2 min post time equilibration. The acquired dynamic MRM transitions, fragmentor voltages and collision energies of quantified analyte are given in Table 1 [36].

Method validation

Given the necessity to adopt tailored test mixtures based on the drug profile of interest it was necessary to validate the methods for confirmation of positive drugs. Repeated injections (n = 6) were performed on a single day to establish the intra-day precision (% CV), while the inter-day precision was evaluated on several days. The accuracy was expressed as the ratio of the compound added to the measured (mean value/nominal value) × 100. Extraction efficiency of the compounds was determined by spiking blank hair samples with 10 pg/mg drugs. Absolute recoveries were determined by comparing the concentrations obtained from peak height ratios of extracted hair samples with added concentrations (standard addition method was used for the analysis of testosterone). Relative recoveries were determined by comparing the % amount of drug recovered from hair matrix against detector response for extracted standards, utilising three different concentrations of quality controls. Stability checks were performed by carrying out three freeze/thaw cycles (between room

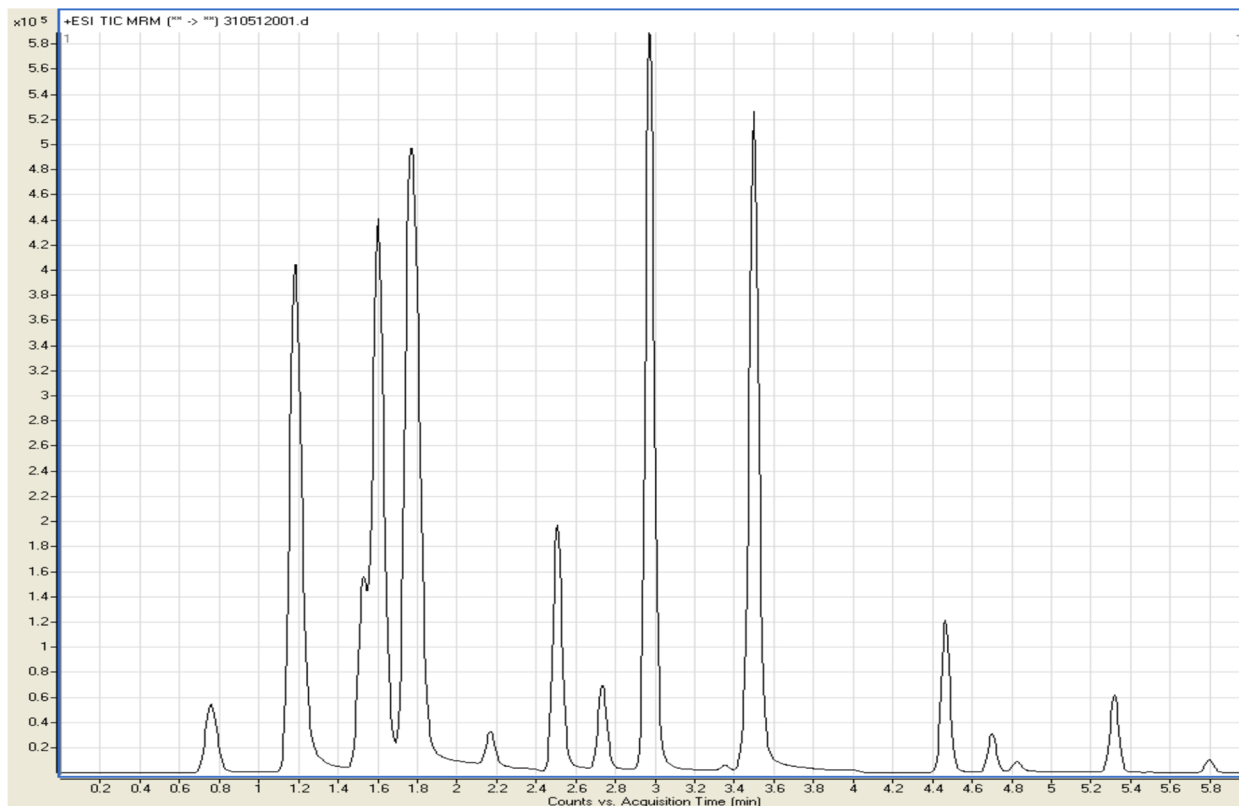


Figure 1 Chromatogram showing peaks for 27 spiked drugs in a single extract of human hair using DYN-MRM-LC-MS/MS. (See text for peak elution times).

temperature to $-20 \pm 5^\circ\text{C}$), in hair samples by spiking blank hair with three different concentrations of quality controls (10, 25 and 50 pg/mg).

Results and discussion

The MassHunter Forensics and Toxicology Dynamic MRM Database Kit was able to screen for 300 analytes of Forensics and Toxicology nature with enhanced sensitivity, all in a single LC/MS analytical run. The Kit also helped to minimize method development for this analysis. The kit came equipped with, MRM transitions for analytes of forensics and toxicological nature along with their optimized fragmentor voltages and collision energy settings. During method development, these MRM transitions were imported from the database to the MassHunter Data Acquisition program where these transitions were used for screening and confirmation purpose. These MRM transitions were also compared against MRM transitions for infused analytes. The Dynamic MRM feature of Agilent Triple Quadrupole mass spectrometer provided us with an adaptive MRM data collection methodology

where the instrument was able to collect data only during a predetermined time window. The advantage is that more than one compound/MRMs can be analyzed in a single run through the Dynamic MRM feature, without losing data quality. The Agilent MassHunter software was used for data generation. Table 1 shows the precursor and product ions of the confirmed analytes along with their fragmentor voltages and collision energies. The protonated molecules $[M + H]^+$ was identified as the precursor ions and the diagnostic product ions were monitored in dynamic multiple reactions monitoring (MRM) mode.

Dynamic MS/MS transitions for the LC/MS Forensics and Toxicology test mixture and their chromatographic evaluation could not be achieved because of the unavailability of this kit in the UK, and due to Home office licensing restrictions in the UK. The commercially available standard Dyn-MRM kit mixture was comprised of about 24 drugs (codeine, oxycodone, amphetamine, methamphetamine, MDA, MDMA, MDEA, hydrocodone, strychnine, cocaine, heroin, meperidine, trazodone, PCP, oxazepam, nitrazepam, verapamil, methadone, lorazepam,

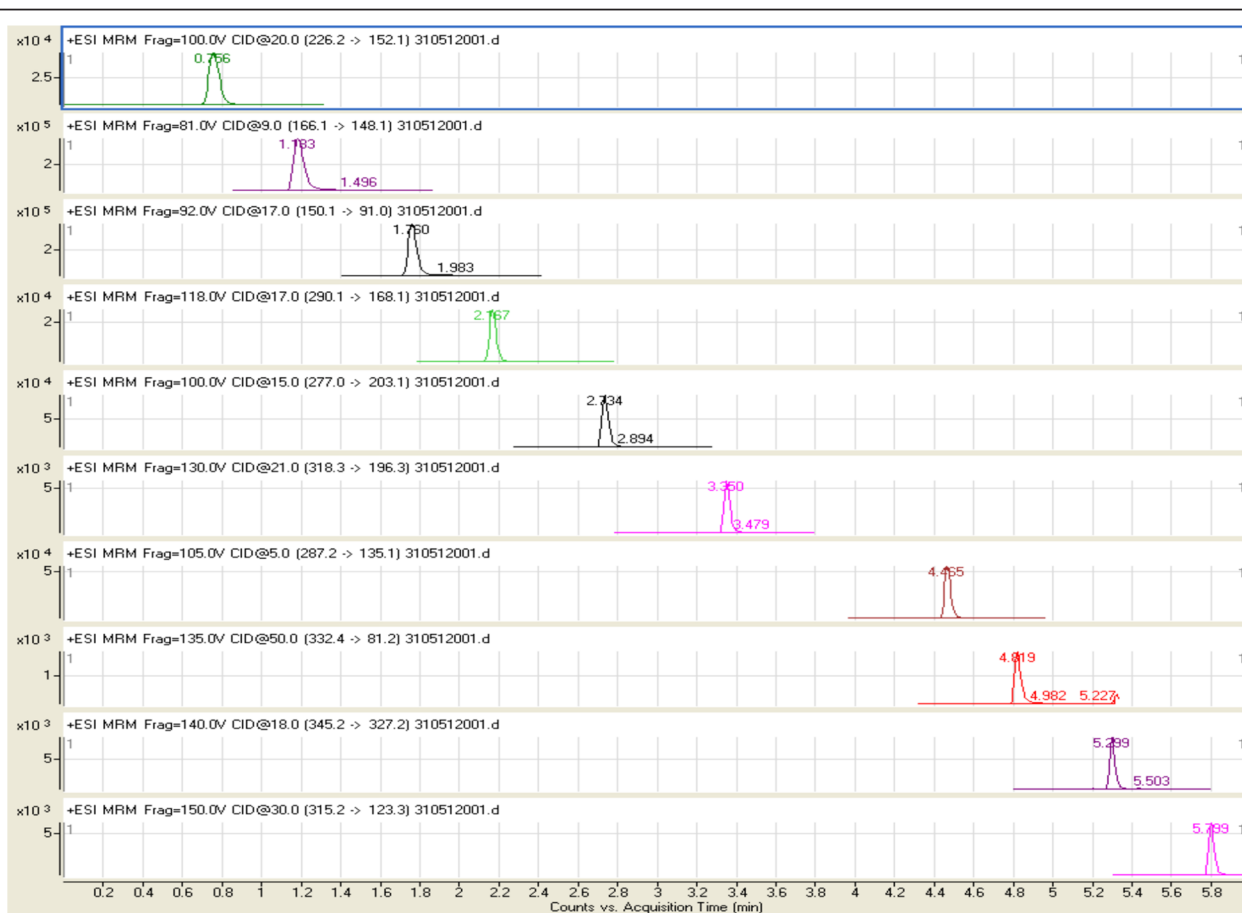


Figure 2 DYN-MRM-LC-MS/MS chromatogram showing delta retention time (Del RT, Δt) window for 10 spiked drugs of abuse in a single extract of human hair sample. From top: terbutaline, ephedrine, methamphetamine, benzoylecgonine, clenbuterol, THC-d₃, boldenone, stanozolol-d₃, THC acid and cannabidiol (THC).

Table 2 Summary of assay validation results

Drugs in (pg/mg)	Concntr. QC's	Linear range	LOD	r ²	Intra-day (n = 6)		Inter-day(n = 6)	
					Precision, % CV	Accuracy,%	Precision, % CV	Accuracy,%
Nandrolone	10	0.1-100	0.05	0.9998	3.23	100.1	4.82	99.8
	25				1.96	100.4	1.95	99.4
	50				1.41	99.9	1.39	100.2
Stanozolol	10	0.1-100	0.05	0.9995	6.3	102.1	9.4	101.2
	25				10.2	98.5	7.1	99.8
	50				8.6	99.4	3.2	101.3
Boldenone	10	0.5-100	0.25	0.9998	4.9	99.9	3.4	101.1
	25				3.3	100.3	2.7	99.8
	50				2.5	99.7	2.3	98.7
Cocaine	10	0.5-100	0.25	0.9998	3.5	98.9	3.2	99.8
	25				4.3	99.4	3.3	98.9
	50				2.9	99.9	2.3	97.8
Benzoylecgonine	10	0.5-100	0.25	0.9987	3.5	107.2	4.7	100.5
	25				2.1	102.3	10.9	102.6
	50				3.3	100.1	8.9	99.9
Amphetamine	10	0.25-100	0.1	0.9979	7.1	98.8	3.2	100.5
	25				5.9	99.6	2.7	100.1
	50				3.2	89.9	3.2	100.3
Methamphetamine	10	0.5-100	0.25	0.9988	5.2	98.3	5.3	100.2
	25				2.3	102.6	2.1	103.3
	50				4.1	100.9	3.6	108.1
Ephedrine	10	1-100	0.5	0.9945	3.8	99.3	4.3	101.2
	25				3.5	98.7	2.1	99.8
	50				3.2	101.2	3.2	99.9
MDMA	10	0.5-100	0.25	0.9968	5.1	99.4	5.1	100.3
	25				2.4	103.2	3.3	110.2
	50				3.7	101.6	6.2	104.2
THC	10	1-100	0.5	0.9996	2.1	107.3	4.3	109.2
	25				3.2	110.2	2.6	105.3
	50				5.1	101.3	3.9	103.2
THC-COOH	10	1-100	0.5	0.9992	4.3	100.1	3.3	101.1
	25				3.6	98.9	3.9	99.9
	50				2.9	98.6	4.3	98.9
11-OH-THC	10	0.5-100	0.25	0.9998	3.6	89.8	5.3	90.5
	25				2.1	89.9	5.9	89.6
	50				5.3	90.9	2.1	99.8
Morphine	10	0.5-100	0.25	0.999	4.2	89.3	3.2	90.9
	25				6.3	88.9	5.3	91.8
	50				3.9	100.9	2.1	99.5
Codeine	10	0.5-100	0.25	0.999	2.9	98.9	1.9	99.5
	25				8.3	99.5	5.6	99.3
	50				2.4	99.3	1.5	9.8

Table 2 Summary of assay validation results (Continued)

Clenbuterol	10	0.5-100	0.25	0.9989	5.1	114.2	4.6	104.3
	25				2.1	89.9	1.3	103.4
	50				1.2	99.9	1.3	105.6
Salbutamol	10	1-100	0.5	0.9998	2.1	99.9	1.9	102.3
	25				1.9	98.4	1.8	101.2
	50				3.2	100.5	2.1	101.1
Propranolol	10	0.1-100	0.05	0.9999	2.2	99.7	1.9	100.5
	25				2.6	89.9	2.5	99.78
	50				1.1	98.9	1.5	98.5
Terbutaline	10	0.5-100	0.25	0.9997	1.7	100.8	1.9	99.9
	25				1.2	99.9	2.1	100.5
	50				2.2	100.5	1.9	102.3
PCP	10	0.5-100	0.25	0.9998	1.4	102.3	2	100.3
	25				2.1	101.3	3.5	100.1
	50				3.2	105.1	1.5	100.5
Ketamine	10	0.5-100	0.25	0.9996	1.6	104.2	1.1	98.9
	25				2.6	102.3	2.1	99.9
	50				2.4	104.2	2.2	99.6

alprazolam, temazepam, proadifen, diazepam and THC). Because of its unavailability in UK, we had to tailor make a test mixture which consisted of the following 21 drugs of interest along with 6 internal standards (nandrolone, stanozolol, boldenone, testosterone, cocaine, benzoylecgonine, amphetamine, methamphetamine, ephedrine, MDMA, THC, THC-COOH, 11-hydroxy THC, morphine, codeine, clenbuterol, salbutamol, propranolol, terbutaline, PCP, ketamine and internal standards like stanozolol-d₃, amphetamine-d₅, methamphetamine-d₅, MDMA-d₅, THC-d₃ and cocaine-d₃). The chromatogram in Figure 1 shows the standard mixture of 27 drugs spiked in human hair extracted and analysed using DYN-MRM-LC-MS/MS method.

The retention times of 21 spiked drugs along with their 6 internal standards in human hair, along with 6 internal standards spiked in human hair where delta retention time was 1, are as follows: Peaks = Morphine (0.71 min), terbutaline (0.75), amphetamine-d₅ (1.67), amphetamine (1.67), methamphetamine (1.76), methamphetamine-d₅ (1.76), MDMA (1.92), 5-MDMA (1.92), ketamine (2.3), benzoylecgonine (2.67), clenbuterol (2.73), cocaine (2.98), cocaine-d₃ (2.98), boldenone (3.22), THC (3.35), THC-d₃ (3.35), nandrolone (3.4), PCP (3.5), codeine (3.54), ephedrine (3.71), salbutamol (3.95), propranolol (4.32), testosterone (4.7), stanozolol (4.81), stanazolol-d₃ (4.81), THC acid (5.29), 11-hydroxy THC (5.3). The Delta retention time for all analytes was 1.

The following representative chromatogram (Figure 2) shows the importance of the short retention time windows

Table 3 The quantitative range and frequency of occurrence of confirmed drugs in human hair

Drugs	Frequency of occurrences	Range of drugs found in hair (pg/mg)
Nandrolone	5	0.12 – 15.28
Stanozolol	14	0.36 – 8.24
Boldenone	1	1.4
Testosterone	233	1.19 – 4.1
Naltrexone	2	8.6 – 16.6
Cocaine	4	1.2 – 10.3
Benzoylecgonine	1	4.4
Amphetamine	6	0.28 – 13.6
Methamphetamine	3	0.884 – 5.6
Ephedrine	1	5.4
MDMA	2	8.6 – 8.76
THC	9	1.2 – 8.4
THC-COOH	1	18.2
11-OH-THC	8	0.61 – 3.1
Morphine	4	1.9 – 10.8
Codeine	5	3.4 – 8.3
Clenbuterol	7	0.95 – 9.2
Salbutamol	7	1 – 12.52
Propranolol	3	1.12 – 12.4
Terbutaline	3	0.79 – 1.8
PCP	4	0.89 – 1.13
Ketamine	4	1.3 – 9.6

functionality (Δt retention time, Δt) for the entire run time. In this example, ten drugs were spiked into a hair sample, extracted and analysed using DYN-MRM-LC-MS/MS method.

Validation results

Intra-day and inter-day precision and accuracy are shown in Table 2. Drug free blank hair was analysed to evaluate the specificity of analytical method. Endogenous drugs (e.g. testosterone) solutions were prepared in human hair samples using the standard additions method. The absolute recoveries of the analytes and internal standards were ranging from 98.2 to 98.7% and 94.3 to 99.6%, respectively. Relative recoveries for analytes and internal standards ranged from 99.8 to 102.2% and 99.8 to 99.9% respectively. In the Table, r^2 represents the linear regression; all the units are represented in pg/mg units. The assay showed excellent sensitivity (LOD = 0.1 pg) and linearity. Where limit of quantitation ranged between 0.1-1 pg/mg. Calibration linearity was not less than ($R^2 = 0.9987$) for the spiked analytes.

Application in hair-analysis

For this study the analysis results of 233 hair samples showed that about 70 volunteers were confirmed positive for some of 20 different of drugs of forensic nature, mainly

steroids and drugs of abuse. Out of these 70 only one volunteer was positive for 5 drugs, 2 volunteers were positive for four different drugs and 11 volunteers were positive for 2 drugs at the same time. The rest of the 56 volunteers were positive for 1 drug each. The range of drugs found in hair and their frequencies are given in Table 3 and Figure 3.

Benefits of dynamic-MRM over MRM technique

Multiple reactions monitoring (MRM) mode has been widely used as a suitable way of conducting quantitative analysis using triple quadrupole mass spectrometry, but the limitation of this method is the use of a series of time segments to monitor MRM transitions. In order to analyse very low quantities of hundreds of compounds of forensic and toxicological nature, in a wide variety of matrices, dynamic MRM is useful. Dynamic MRM method development and modifications are a lot easier and more reliable and provide better quality data compared to traditional time segments in MRM mode. A further benefit of this method is that, rather than performing scans for all the analytes during the entire method, the mass spectrometer only monitors MRM transitions for the analytes that elute in that particular time segment. The result is that there are fewer transitions recorded during each MS scan, allowing the mass spectrometer to use a longer dwell time and more data points per peak.

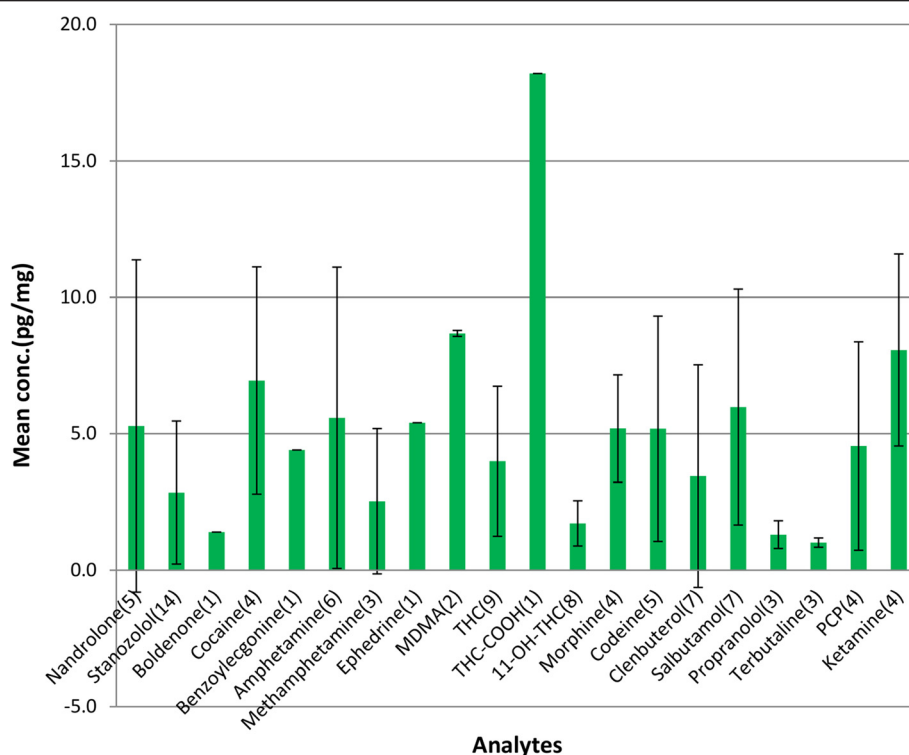


Figure 3 Mean concentrations of analytes along with the frequency of positive occurrences (shown in brackets). The error bars represents standard deviation.

In other words, analytes are only monitored while they are eluting from the LC to focus detection. This approach also allows total exclusion of background and interferences problems. These qualities made it possible to use dynamic MRM mode to screen for about 20 drugs of forensic and toxicological nature in hair samples just over 8 minute time period. Since a quantifier ion and a qualifier ion were necessary for confirmatory purposes, a total of about 40 dynamic MRM transitions were employed during the analysis, covering classes of analytes such as steroids, amphetamines, cannabinoids, cocaine's and bronchodilators etc.

Advantages of hair-dynamic MRM drug analyses method

This study reports the application of a hair-based multi-drug analytical approach which has high sensitivity and accuracy and is capable of analysing some 200 drugs/metabolites in one chromatographic run. The benefits of the novel method over current approaches include: i) long-term such as quarterly detection windows (for 3 cm hair samples) and ease of 'non-invasive hair' collection, reduced infection risk and storage demands, ii) using Dyn-MRM software affords the ability to screen for some 200 compounds of interest simultaneously in one chromatographic run (i.e. 8 min) using a test mixture to calibrate the LC component, iii) develop a full validation for compounds of interest and thus perform full confirmation analysis, iv) adapt the software to add new compounds of interest.

It should be noted that the focus of this study was on the range of drugs and their metabolites listed. For this reason the calibrant test mixture was prepared and validation steps were completed with these selected compounds. The method also allows screening for any of the other 200 plus drugs/metabolites from spectral library and full validation can be added using reference and standard material.

Conclusion

This study reports the development and validation of a method to analyse drugs of forensic and toxicological nature in human hair samples. The method is selective, sensitive and robust for the screening and confirmation of multiple drugs in a single analysis and will be very useful tool for the analysis of large array of controlled substances and drugs of abuse. This is the first report of the application of this proprietary system to drugs of abuse in hair samples. In this example, it employs the Dyn-MRM software to screen and confirm for 20 selected analytes (via full validation) of forensic interest present in the hair samples in a single chromatographic run.

Competing interests

The authors declare that they have no competing interests.

Authors' contributions

DPN and AP initiated the study. All authors contributed to the study design, interpretation of the results and preparation of the manuscript and have read and approved the final version. The method development, validation and sample analyses were conducted by IS with contributions from DPN. MU and MR also contributed by recruitment and in the hair sample collection.

Acknowledgement

The project on drug- and doping-related explicit and implicit social cognition, which provided the hair samples, was funded by the World Anti-Doping Agency Social Science Grant (2010/2011).

Author details

¹School of Life Sciences, Kingston University, Kingston-upon-Thames, Surrey, UK
²Faculty of Physical Education and Sports Sciences, Semmelweis University, Budapest, Hungary.
³Eötvös Lóránd University, Faculty of Education and Psychology, Budapest, Hungary.

Received: 19 September 2014 Accepted: 26 November 2014

Published online: 13 December 2014

References

1. Verstraete AG: **Detection times of drugs of abuse in blood, urine, and oral fluid.** *Ther Drug Monit* 2004, **26**:200–205.
2. Kintz P, Cirimele V, Sachs H, Jeanneau T, Ludes B: **Testing for anabolic steroids in hair from two bodybuilders.** *Forensic Sci Int* 1999, **101**:209–216.
3. Kronstrand R, Nyström I, Strandberg J, Druid H: **Screening for drugs of abuse in hair with ion spray LC–MS–MS.** *Forensic Sci Int* 2004, **145**:183–190.
4. Strahm E, Sottas P-E, Schweizer C, Saugy M, Dvorak J, Saudan C: **Steroid profiles of professional soccer players: an international comparative study.** *B J Sports Med* 2009, **43**:1126–1130.
5. Kintz P: **Testing for anabolic steroids in hair: a review.** *Legal Med* 2003, **5**:S29–S33.
6. Deshmukh N, Hussain I, Barker J, Petroczi A, Naughton DP: **Analysis of anabolic steroids in human hair using LC–MS/MS.** *Steroids* 2010, **75**:710–714.
7. Shen M, Xiang P, Shen B, Bu J, Wang M: **Physiological concentrations of anabolic steroids in human hair.** *Forensic Sci Int* 2009, **184**:32–36.
8. Gosetti F, Mazzucco E, Gennaro MC, Marengo E: **Ultra high performance liquid chromatography tandem mass spectrometry determination and profiling of prohibited steroids in human biological matrices. A review** *J Chromatogr B* 2013, **927**:22–36.
9. Gambelunghe C, Somavilla M, Ferranti C, Rossi R, Aroni K, Manes N, Bacci M: **Analysis of anabolic steroids in hair by GC/MS/MS.** *Biomed Chromatogr* 2007, **21**:369–375.
10. Strano-Rossi S, Castrignanò E, Anzillotti L, Odoardi S, De-Giorgio F, Bermejo A, Pascali VL: **Screening for exogenous androgen anabolic steroids in human hair by liquid chromatography/orbitrap-high resolution mass spectrometry.** *Anal Chim Acta* 2013, **793**:61–71.
11. Vincenti M, Salomone A, Gerace E, Pirro V: **Application of mass spectrometry to hair analysis for forensic toxicological investigations.** *Mass Spectrom Rev* 2013, **32**:312–332.
12. Gomez C, Fabregat A, Pozo ÓJ, Marcos J, Segura J, Ventura R: **Analytical strategies based on mass spectrometric techniques for the study of steroid metabolism.** *TRAC* 2014, **53**:106–116.
13. Pötsch L, Skopp G, Moeller MR: **Biochemical approach on the conservation of drug molecules during hair fiber formation.** *Forensic Sci Int* 1997, **84**:25–35.
14. Pragst F, Balikova MA: **State of the art in hair analysis for detection of drug and alcohol abuse.** *Clin Chim Acta* 2006, **370**:17–49.
15. Nielsen MW, Lasaroms JJ, Mulder PP, Van Hende J, van Rhijn JHA, Groot MJ: **Multi residue screening of intact testosterone esters and boldenone undecylenate in bovine hair using liquid chromatography electrospray tandem mass spectrometry.** *J Chromatogr B* 2006, **830**:126–134.
16. Gaillard Y, Pépin G: **Testing hair for pharmaceuticals.** *J Chromatogr B* 1999, **733**:231–246.
17. Rusevska K, Zdravkovski Z: **Simple extraction method for detecting exogenous substances in scalp hair by GC-MS.** *J Hyg Eng Des* 2013, **2**:59–67.

18. Calvo D, Tort N, Salvador JP, Marco M-P, Centi F, Marco S: Preliminary study for simultaneous detection and quantification of androgenic anabolic steroids using ELISA and pattern recognition techniques. *Analyst* 2011, **136**:4045–4052.
19. Thieme D, Anielski P, Grosse J, Sachs H, Mueller R: Identification of anabolic steroids in serum, urine, sweat and hair: comparison of metabolic patterns. *Anal Chim Acta* 2003, **483**:299–306.
20. Herrin GL, McCurdy HH, Wall WH: Investigation of an LC-MS-MS (QTrap®) method for the rapid screening and identification of drugs in postmortem toxicology whole blood samples. *J Anal Toxicol* 2005, **29**:599–606.
21. Viette V, Guillaume D, Mylonas R, Mauron Y, Fathi M, Rudaz S, Hochstrasser D, Veuthey JL: A multi-target screening analysis in human plasma using fast liquid chromatography–hybrid tandem mass spectrometry (Part I). *Clin Biochem* 2011, **44**:32–44.
22. Leung GN, Leung DK, Wan TS, Wong CH: High throughput screening of sub-ppb levels of basic drugs in equine plasma by liquid chromatography–tandem mass spectrometry. *J Chromatogr A* 2007, **1156**:271–279.
23. Vincenti M, Cavanna D, Gerace E, Pirro V, Petrarulo M, Di Corcia D, Salomone A: Fast screening of 88 pharmaceutical drugs and metabolites in whole blood by ultrahigh-performance liquid chromatography–tandem mass spectrometry. *Anal Bioanal Chem* 2013, **405**:863–879.
24. Maurer HH: What is the future of (ultra) high performance liquid chromatography coupled to low and high resolution mass spectrometry for toxicological drug screening? *J Chromatogr A* 2013, **1292**:19–24.
25. Montesano C, Johansen SS, Nielsen MKK: Validation of a method for the targeted analysis of 96 drugs in hair by UPLC–MS/MS. *J Pharm Biomed Anal* 2014, **88**:295–306.
26. Oberacher H, Schubert B, Libiseller K, Schweissgut A: Detection and identification of drugs and toxicants in human body fluids by liquid chromatography–tandem mass spectrometry under data-dependent acquisition control and automated database search. *Anal Chim Acta* 2013, **770**:121–131.
27. Dresen S, Ferreiros N, Gnann H, Zimmermann R, Weinmann W: Detection and identification of 700 drugs by multi-target screening with a 3200 Q TRAP® LC-MS/MS system and library searching. *Anal Bioanal Chem* 2010, **396**:2425–2434.
28. Mueller C, Weinmann W, Dresen S, Schreiber A, Gergov M: Development of a multi-target screening analysis for 301 drugs using a QTrap liquid chromatography/tandem mass spectrometry system and automated library searching. *Rapid Commun Mass Spectrom* 2005, **19**:1332–1338.
29. Niessen W: Fragmentation of toxicologically relevant drugs in negative ion liquid chromatography–tandem mass spectrometry. *Mass Spectrom Rev* 2012, **31**:626–665.
30. Broecker S, Herre S, Wüst B, Zweigenbaum J, Pragst F: Development and practical application of a library of CID accurate mass spectra of more than 2,500 toxic compounds for systematic toxicological analysis by LC–QTOF-MS with data-dependent acquisition. *Anal Bioanal Chem* 2011, **400**:101–117.
31. Broecker S, Herre S, Pragst F: General unknown screening in hair by liquid chromatography–hybrid quadrupole time-of-flight mass spectrometry (LC–QTOF-MS). *Forensic Sci Int* 2012, **218**:68–81.
32. Polettini A, Gottardo R, Pascali JP, Tagliaro F: Implementation and performance evaluation of a database of chemical formulas for the screening of pharmaco/toxicologically relevant compounds in biological samples using electrospray ionization-time-of-flight mass spectrometry. *Anal Chem* 2008, **80**:3050–3057.
33. Jamin EL, Bonvallot N, Tremblay-Franco M, Cravedi J-P, Chevrier C, Cordier S, Debrauwer L: Untargeted profiling of pesticide metabolites by LC–HRMS: an exposomics tool for human exposure evaluation. *Anal Bioanal Chem* 2014, **406**:1149–1161.
34. Nielsen MKK, Johansen SS, Dalsgaard PW, Linnet K: Simultaneous screening and quantification of 52 common pharmaceuticals and drugs of abuse in hair using UPLC–TOF-MS. *Forensic Sci Int* 2010, **196**:85–92.
35. Agilent Technical Overview (2007). Ion optics innovations for increased sensitivity in hybrid MS systems. Agilent Technologies USA 5989-7408EN. (http://www.chem.agilent.com/Library/technicaloverviews/Public/5989-7408EN_HI.pdf).
36. Peter S, Thomas G, Frank K, Tim S, M K (2009). New dynamic MRM mode improves data quality and triple quad quantification in complex analyses. Agilent Technologies USA 5990-3595EN. (http://www.chem.agilent.com/Library/technicaloverviews/Public/5990-3595en_lo%20CMS.pdf).

doi:10.1186/s13065-014-0073-0

Cite this article as: Shah et al.: Hair-based rapid analyses for multiple drugs in forensics and doping: application of dynamic multiple reaction monitoring with LC-MS/MS. *Chemistry Central Journal* 2014 **8**:73.

Publish with **ChemistryCentral** and every scientist can read your work free of charge

“Open access provides opportunities to our colleagues in other parts of the globe, by allowing anyone to view the content free of charge.”

W. Jeffery Hurst, The Hershey Company.

- available free of charge to the entire scientific community
- peer reviewed and published immediately upon acceptance
- cited in PubMed and archived on PubMed Central
- yours — you keep the copyright



Submit your manuscript here:
<http://www.chemistrycentral.com/manuscript/>

ChemistryCentral