

Characterization of a Nematic Mixture by Reversed-Phase HPLC and UV Spectroscopy: Application to Phase Behavior Studies in Liquid Crystal-CO₂ Systems

Ana R. E. Brás¹, Sara Henriques¹, Teresa Casimiro¹, Ana Aguiar-Ricardo¹,
João Sotomayor¹, Jorge Caldeira^{1,3}, Celina Santos¹, Madalena Dionísio^{1,2}

¹REQUIMTE/CQFB, Departamento de Química, Faculdade de Ciências e Tecnologia
Universidade Nova de Lisboa, 2829-516 Caparica, Portugal

²E-mail address: madalena.dionisio@dq.fct.unl.pt.

³Instituto Superior de Ciências da Saúde – Sul, 2829-516 Caparica, Portugal

In the present work a simple but selective reversed-phase high performance liquid chromatographic method (HPLC) was developed for analysis of the nematic liquid crystal mixture E7 to precisely determine the composition of the liquid crystal mixture used in PDLC films preparation. Ultra-violet absorption spectrophotometry (UV) experiments were carried out to determine the HPLC detection wavelength and to characterize the absorptivity constants of E7 constituents. The developed technique is applied in the case of equilibrium solubility studies for E7 in supercritical carbon dioxide (scCO₂). The results indicate unambiguously that scCO₂ can fractionate the mixture towards the E7 components. The four single components peaks of the E7 mixture were distinctively separated by this method which enabled the determination of the solubility of E7 constituents in the scCO₂.

I. Introduction

Liquid crystals have been widely studied due to their unique electro-optical properties, which resulted in a variety of potential applications ranging from flat panel displays to laser beam steering and optical switches for telecommunications [1,2,3,4].

Among liquid crystals, a multicomponent nematic mixture, E7, has been widely used [5,6,7] due to a high optical anisotropy [8] and a large temperature range where it maintains anisotropic characteristics. Its nematic to isotropic transition occurs at 58 °C [9]. At room temperature it still exhibits a nematic phase, which easily supercools upon cooling, avoiding crystallization. Thus, its liquid crystalline properties are extended down to the glass transition at – 62 °C [9,10].

Most applications rely on the specific nature of electro-optical response functions together with the possibility to control the phase behavior/miscibility and thermophysical properties of these systems.

The present work is part of a research project designed for the development of polymer-dispersed liquid crystals, PDLC, using supercritical fluid technology. Supercritical CO₂ has excellent plasticising properties and can easily penetrate into polymer films and extract or precipitate substances in those matrices, simply by manipulating pressure [11,12]. In order to design these processes and especially for scale-up calculations from laboratory experiments, it is important a detailed knowledge of phase behavior [13,14]. In the present work, a simple

analytical technique was developed to determine the solubility of E7 in scCO₂ and to identify the partitioning of the E7 components into the CO₂ rich phase as a function of pressure and temperature conditions.

The E7 mixture was studied by both ultra-violet absorption spectroscopy and high performance liquid chromatography in order to determine its composition and find a suitable method to quantify its solubility in supercritical CO₂.

The UV spectra also provided the respective molar absorptivities (ϵ) and the choice of adequate wavelength for the HPLC assays.

UV and HPLC analytical techniques allowed the identification of each E7 component and the determination of E7 solubility in supercritical CO₂.

II. Experimental Results

2.1. Reagents

HPLC Milipore water, methanol CHROMASOLV[®] for high-performance Liquid Chromatography grade from Sigma-Aldrich GmbH (Germany), acetonitrile E CHROMASOLV[®] for HPLC from Riedel-de-Haën AG (Germany), and carbon dioxide from Soc. Port. do Ar Líquido, “ARLÍQUIDO”, S.A. (Portugal). The nematic liquid crystal E7 and its pure components (see Figure 1), 4-Cyano-4'-*n*-pentyl-1,1'-biphenyl (5CB), 4-Cyano-4'-*n*-heptyl-1,1'-biphenyl (7CB), 4-Cyano-4'-*n*-octyloxy-1,1'-biphenyl (8OCB) and 4-Cyano-4''-*n*-pentyl-1,1',1''-terphenyl (5CT) were kindly supplied by Merck KGaA (Darmstadt, Germany). All the products were used as received.

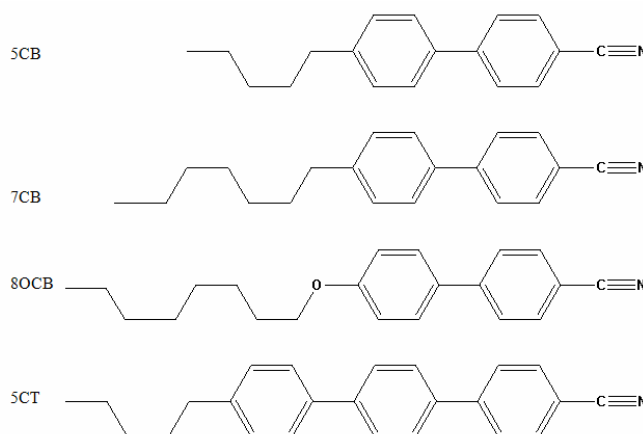


Figure 1: Structures of the components of the nematic liquid crystal mixture E7.

2.2. Sample Preparation

E7 liquid crystal is a highly viscous liquid as well as 5CB and 7CB. 8OCB and 5CT are both solids at room temperature. Solutions of E7 and its pure components within the concentration range 4.0×10^{-6} and 7.5×10^{-4} M were prepared by weighting and dissolved in HPLC grade acetonitrile. The averaged molar mass of E7, \overline{M}_M , was calculated by weighted

sum of the individual molar masses M_{M_i} of the E7 pure components and their molar fractions, x_i , ($\overline{M}_M = \sum_i x_i M_{M_i}$); the following composition (w/w) was assumed: 51 % of 5CB ($M_M=249.36 \text{ g mol}^{-1}$), 25 % of 7CB ($M_M=277.41 \text{ g mol}^{-1}$), 16 % of 8OCB ($M_M=307.44 \text{ g mol}^{-1}$) and 8 % of 5CT ($M_M=325.45 \text{ g mol}^{-1}$).

2.3. UV Spectroscopy

UV Spectrophotometry spectra of samples in HPLC grade acetonitrile were recorded with SHIMADZU UV-2501 PC, UV-VIS spectrophotometer at room temperature.

The absorption spectra of E7 and its components in acetonitrile solutions were measured between 190 and 400 nm. Figure 2 presents the obtained molar absorptivity (ϵ) spectra within 220 and 350 nm. The observed absorption band corresponds to a $\pi-\pi^*$ transition as found in single liquid crystals, with a maximum wavelength located at 281 nm. Two additional bands are also reported in literature, both below 220 nm, due to $\sigma \rightarrow \sigma^*$ and a first $\pi-\pi^*$ transition [15,16]. These bands are not shown in the spectrum since they are located in the cut off region of the solvent.

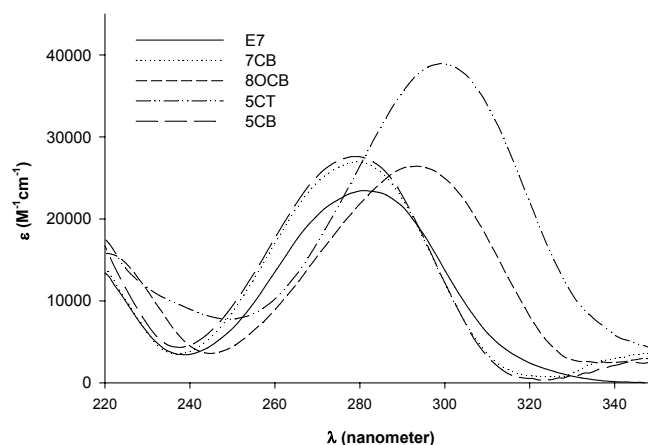


Figure 2: Molar absorptivity spectra of E7 and its components (see legend).

The 5CB spectrum presents a maximum absorption wavelength located at 279 nm, in agreement with λ_{max} reported in literature [15], as well as 7CB. The λ_{max} of both 8OCB and 5CT shifts to higher wavelength (see Table 1) due to the presence of, respectively, the oxygen atom and the third conjugated phenyl ring. The molar absorptivities (ϵ) at the maximum wavelength for the E7 mixture and its constituents were also obtained according to the Beer law within the previous concentration range. The obtained values, ϵ and λ_{max} are listed in Table 1.

Component	λ_{max} (nm)	ϵ ($\text{M}^{-1}\text{cm}^{-1}$)
E7	281	23.4×10^3
5CB	279	27.6×10^3
7CB	279	26.9×10^3
8OCB	294	26.4×10^3
5CT	299	38.9×10^3

Table 1: Molar Absorptivity (ϵ) values of the E7 mixture and its components diluted in acetonitrile.

The UV spectrum presented in Figure 2 shows negligible absorptivity of E7 at wavelengths higher than 365 nm, the wavelength of either the pulsed laser beam used in the liquid crystal deposition or the UV light used for curing sealant in liquid crystal devices [17]. Thus, these results in acetonitrile media seem to confirm the reported stability of pure E7 to photodegradation [18], even in solutions stored for periods as long as 6 months.

2.4. HPLC Analysis

Chromatographic analysis was carried out using a Merck L-7100 HPLC equipped with an L-7400 UV detector and D-7000 computer interface. Compounds were detected at 277 nm. Elution was performed with 79 % methanol (HPLC grade) / 21 % water (Milipore) with a flow rate of 1 mL.min⁻¹, using a water jacket thermostatised at 25 °C. A Purospher[®] 100 STAR RP-18e column (250 x 4 mm, 5µm particle) from Merck KGaA (Darmstadt, Germany) was used. Quantitative determinations were performed with an injection loop of 20 µL, in a time synchronized autosampler (Marathon) from standard compounds dissolved in acetonitrile or from supercritically extracted solutions.

Baseline resolution for all the components including the 7CB and 8OCB critical pair [19], was achieved with the above stated conditions.

The retention times of commercial E7 were 25.57, 55.45, 59.71 and 84.09 min respectively for 5CB, 7CB, 8OCB and 5CT. The attribution of the peaks to the different single components was accomplished by measuring the pure liquid crystals separately under the same conditions.

Through the calibration curves, Figure 3, it was possible to achieve a relation between peak areas and the single components concentrations, which enabled us to determine the exact E7 composition.

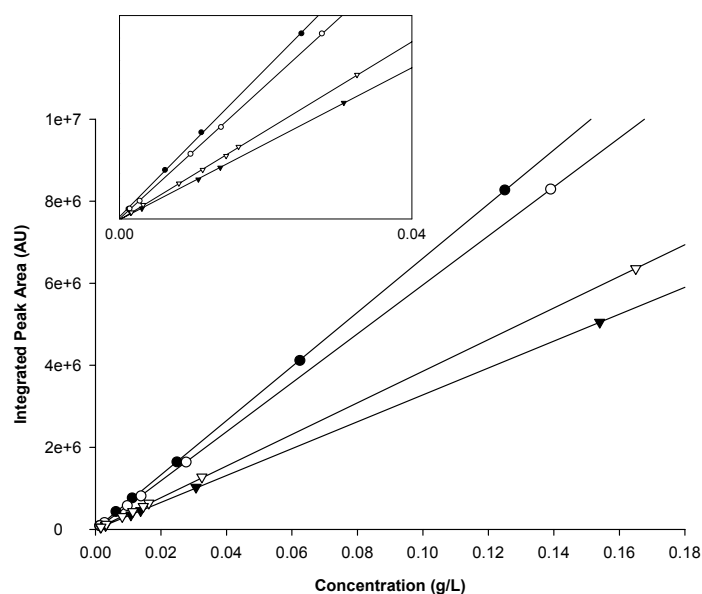


Figure 3: Calibration curves (integrated HPLC peak area vs. concentration) for E7 constituents (5CB-full circles: $y = 6.60 \times 10^7 x + 1.08 \times 10^3$; 7CB-open circles: $y = 5.97 \times 10^7 x - 6.29 \times 10^3$; 8OCB-full triangles: $y = 3.28 \times 10^7 x + 2.57 \times 10^3$ and 5CT-open triangles: $y = 3.85 \times 10^7 x + 3.85 \times 10^3$. All $r^2 = 0.9999$ or better). The inset amplifies the region of lowest concentrations.

2.5 Equilibrium Solubility Measurements

The equilibrium solubility measurements were undertaken in a high-pressure apparatus, which operates according to the static analytical method already described elsewhere [20,21]. This apparatus operates between 290 and 345 K and at pressures between the atmospheric pressure up to 30 MPa [22]. The central piece of the apparatus is a sapphire tube cell of 35 mL useful internal volume, and it allows full observation of the inside. In a typical experiment, the cell is loaded with a known amount of E7 and then purged with a low CO₂ flow to degas. The CO₂ is introduced by means of a HIP manual pressure generator. Pressure is measured by a pressure transducer Newport M3 (Heize HPO and digital panel INFP100 FS). Temperature is measured with a platinum-resistance T (5616 RTD) in contact with the cell, connected to a temperature controller (2100), both from Hart Scientific. Typical temperature stability is ± 0.01 K, at 313.15 K. After the desired temperature had been reached, the pressure is adjusted to the desired value. By stirring the mixture, time for equilibration of the phases is reduced, the pressure reaches a plateau, and the equilibration is continued for at least 40 min. Once the equilibrium is achieved one sample is withdrawn from the rich-CO₂ phase through the six-port valve located on the top of the cell, into the sample loop. The gas in each sample is quantified by expansion into calibrated volumes and the amount of CO₂ calculated from the measurement of the pressure increase at the working temperature. The sampling loop is cleaned with a CO₂ low-flow, after washing the lines with the solvent. The chosen solvent for cleaning the sampling loop was acetonitrile. In this work, the gaseous samples of the vapor-liquid equilibrium were collected at pressures from 5 MPa up to 18 MPa and at 313.15 K.

III. Discussion

The HPLC results allowed the determination of commercial E7 composition: (50.5 \pm 0.2)% 5CB; (24.7 \pm 0.1)% 7CB; (16.8 \pm 0.2)% 8OCB and (8.1 \pm 0.1)% 5CT. These values are in agreement with the nominal values indicated by Merck (51%, 25%, 16% and 8% respectively) [7,23]. As far as we know this is the first time that an accurate analytical determination of E7 composition is reported in literature.

The four single peaks corresponding to the individual E7 components were clearly identified in the HPLC chromatogram (Figure 4 A). This work represents an improvement to other liquid crystal mixtures separation methods reported in literature [19,24] where shoulders and poor peaks were observed, using equipment available in most of laboratories.

The separation efficiency achieved and the accordance of the so determined composition with data furnished by suppliers, allow us to conclude that the proposed reversed phase HPLC method is valid to accurately analyse E7 in different conditions. Thus, the method was applied to determine the equilibrium solubility of E7 and its components in scCO₂.

Samples of liquid crystal mixture collected from the rich-CO₂ phase were diluted in acetonitrile to a convenient known volume. The resulting solutions were then analyzed by HPLC.

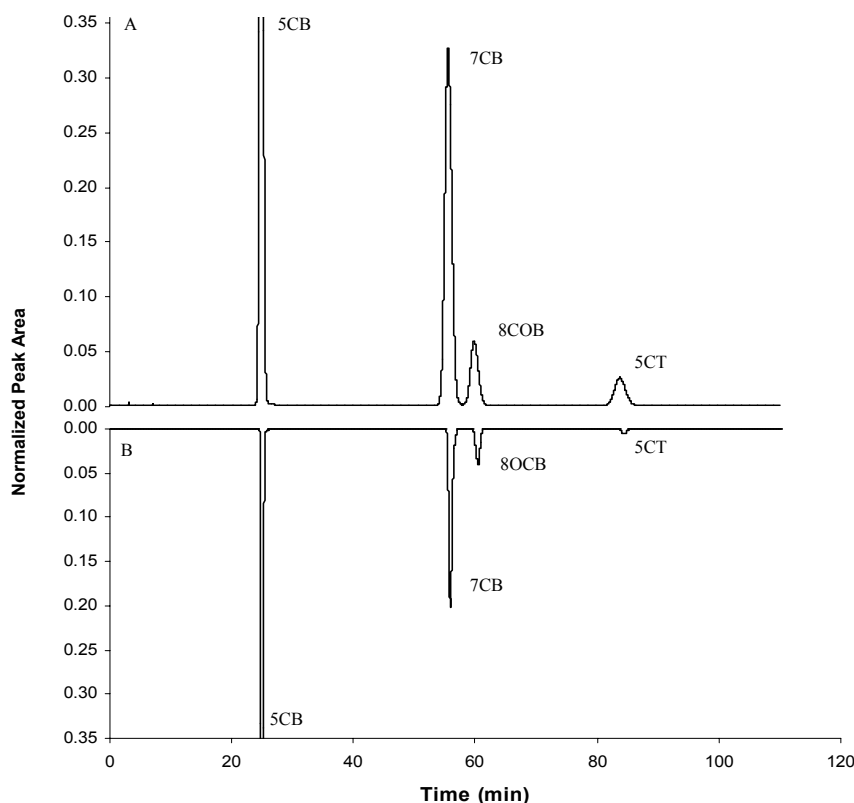


Figure 4: HPLC chromatograms of E7: A- mixture commercially available and B- mixture from VLE sample (experiment at 40 °C and 13.59 MPa). The peaks were normalized by the 5CB peak and were represented in this scale for a better viewing of the smaller peaks.

HPLC chromatograms of commercially available E7 and of the liquid crystal mixture obtained from equilibrium solubility measurements at 313.15 K and 13.59 MPa are compared in Figure 4A and B, respectively. The detailed analysis of the chromatograms shows that the ratio between the sample peaks areas obtained in the solubility measurements is different from the observed ratio in the neat E7.

Using the calibration curves, concentrations of the E7 components were determined for each solubility measurement. In Figure 5, it is represented the relative abundance of each component relative to the total solubilized E7 as a function of pressure, at 313.15 K. Horizontal lines in this figure indicate the composition of 5CB, 7CB, 8OCB and 5CT in parent E7. From this plot it is clear that CO₂ can fractionate E7 components. The 5CB and the 7CB are clearly more soluble in scCO₂ than 8OCB and 5CT in the studied pressure range. Additionally, the composition of the liquid crystal mixture solubilized in the rich-CO₂ phase tends to a constant value at pressures higher than 13.0 MPa: 64.1%5CB, 25.1% 7CB, 9.2% 8OCB and 1.55% 5CT.

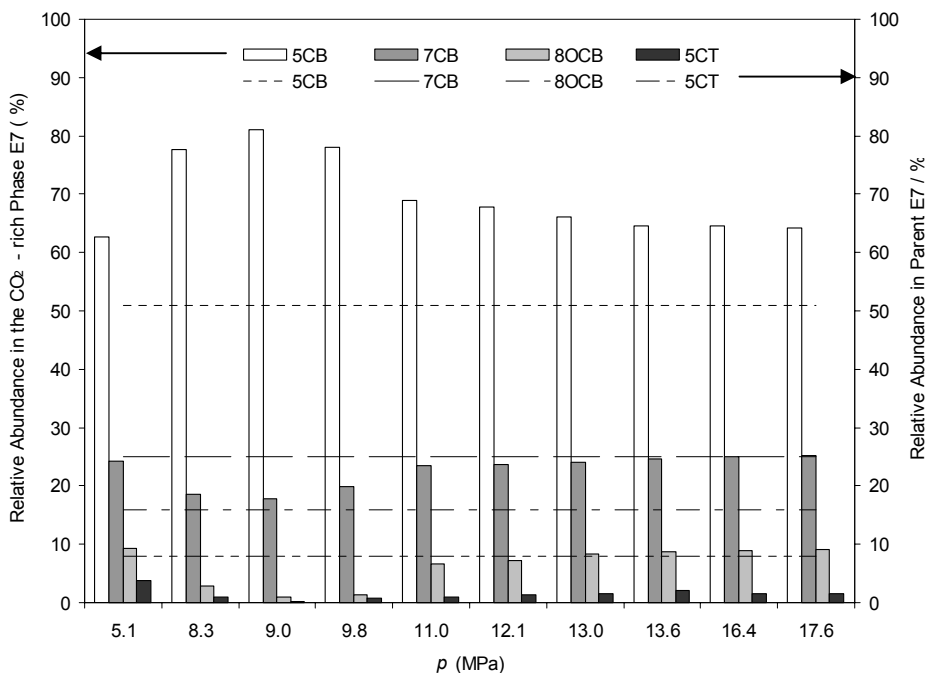


Figure 5: Relative abundance of each E7 component relative to total E7 in the CO₂ rich phase and relative to the initial parent E7, at 313.15 K.

The overall solubility of E7 in CO₂ is presented in Figure 6 as a function of pressure, at 313.15 K. The effect of pressure on the total E7 solubility follows the expected trend as the solvent capacity increases with pressure at constant temperature.

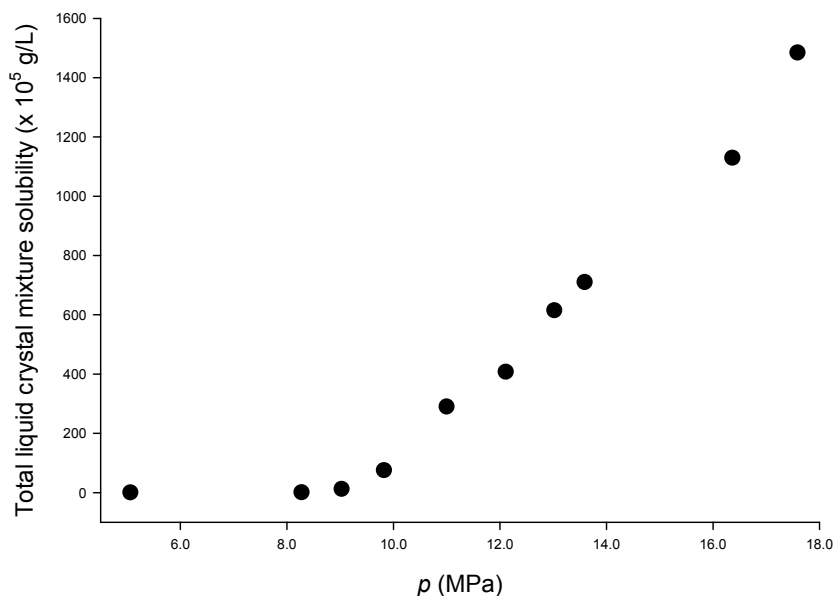


Figure 6: Total liquid crystal mixture solubility as a function of pressure in supercritical CO₂ at 313.15K.

IV. Conclusions

This analytical reversed phase HPLC method enables the quantitative determination and resolution of all E7 components in a reproducible way, using conventional instrumentation. Furthermore, parallel UV studies provided the molar absorptivities and maximum absorption wavelengths of E7 and all its constituents, in diluted acetonitrile solutions.

HPLC procedure meets the demands of analytical determination of E7 mixture components since some of their applications can induce composition variation. This is evident in selective miscibility relatively to the polymer matrix of some PDLC devices or in the studied example of differential extraction with supercritical fluids namely scCO₂. It was demonstrated that selective solubility of E7 occurs in supercritical CO₂, enriching the extracted samples in 5CB and 7CB components and approaching an invariant composition for pressures higher than 13.0 MPa. To our knowledge, this is the first evidence of supercritical CO₂ ability to fractionate E7.

Acknowledgments. This work was supported by Fundação para a Ciência e Tecnologia through the projects POCTI/ QUI/ 35429/ 2000, POCTI/ CTM/ 37435/ 2001, POCTI/ CTM/ 47363/ 2002, grant SFRH/BPD/11665/2002 and by FEDER and FSE.

V. References

-
1. M. Schadt, *Annu. Rev. Mater. Sci.*, **27**, 305, (1997).
 2. J.W. Doane, *Liquid Crystals- Applications and Uses Vol.1* edited by B. Bahadur (Singapore: World Scientific) 361-395, (1990).
 3. H.S. Kitzerow, *Liq. Cryst.*, **16**, 1, (1994).
 4. P.S. Drzaic, *Liquid Crystals Dispersions* (Singapore: World Scientific), (1995).
 5. S. T. Wu, D. K. Yang, *Reflective Liquid Crystal Displays*, Wiley: New York, (2001).
 6. M.H. Godinho, J.L. Figueirinhas, A.F Martins, *Liq. Cryst.*, **20**, 373, (1996).
 7. U. Maschke, M. Benmouna, X. Coqueret, *Macromol. Rapid Commun.*, **23**, 159, (2002).
 8. S. W. Leonard, J. P. Mondia, H. M. van Driel, O. Toader, S. John, K. Busch, A. Birner, U. Gösele, V. Lehmann, *Physical Review B, Condensed Matter and Materials Physics III, series 61*, No. 4, (2000).
 9. F. Roussel, J-M. Buisine, U. Maschke, X. Coqueret, *Mol. Cryst. Liq. Cryst.*, **299**, 321, (1997).
 10. M. T. Viciosa, A. M. Nunes, A. Fernandes, P. L. Almeida, M. H. Godinho, M. D. Dionísio, *Liq. Cryst.*, **29**, 429, (2002).
 11. I. Kikic, F. Vecchione, *Curr. Opin. Solid St. M.*, **7**, 399, (2003).
 12. D.L. Tomasko, H.B. Li, D.H. Liu, X.M. Han, M.J. Wingert, L.J. Lee, K.W. Koelling, *Ind. Eng. Chem. Res.*, **42**, 6431, (2003).
 13. J.A. Lopes, D. Gourguillon, P.J. Pereira, A.M. Ramos, M. Nunes da Ponte, *J. Supercrit. Fluids*, **16**, 261, (2000).
 14. J.C. Gonzáles, M.R. Vieytes, A.M. Botana, J.M. Vieites, L.M. Botana, *J. Chromatogr. A*, **910**, 119, (2001).

-
15. S. T. Wu, *J. Appl. Phys.*, **69**(4), 2080, (1991).
 16. I.C. Khoo, S. T. Wu, *Optics and Nonlinear Optics of Liquid Crystals* (World Scientific, Singapore), (1993).
 17. J. Gonzalo, P.E. Dyer, M. Hird, *Appl. Phys. Lett.*, **71**(19), 2752, (1997).
 18. P.T. Lin, S. T. Wu, C.Y. Shang, C.S. Hsu, *Mol. Cryst. Liq. Cryst.*, **411**, 243, (2004).
 - 19 L Bedjaoui,.; Benmouna, M.; Bouchaour, T.; Coqueret, X.; Maschke, U. *Mol. Cryst. Liq. Cryst.*, **411**, 537-544, (2004).
 20. C. M. M. Duarte, M. Crew, T. Casimiro, A. Aguiar-Ricardo, M. Nunes da Ponte, *J. Supercrit. Fluids*, **22**, 87, (2002).
 21. T. Casimiro, F. Montilla, S. Garcia, T. Avilés, A. Shariati, S. Raessi, C. J. Peters, A. Aguiar-Ricardo, M. Nunes da Ponte, *J. Supercrit. Fluids*, **31**, 1, (2004).
 22. D. Gourgouillon, M. Nunes da Ponte, *Phys. Chem.. Chem. Phys.*, **1**, 5369, (1999).
 23. Information from Merck KGa, Darmstadt, Germany.
 - 24 P.J. Schoenmakers, F.C.C.J.G. Verhoeven, H.M. Van Den Bogaert, *J. Chromatogr.*, **371**, 121, (1986).