Dong et al. Acta Vet Scand (2021) 63:10 https://doi.org/10.1186/s13028-021-00575-1

Acta Veterinaria Scandinavica

RESEARCH Open Access

Prevalence of heterophilic antibodies in serum samples from horses in an equine hospital, and elimination of interference using chicken IgY

Bo Dong^{1,2}, Daniel Bergman¹ and Bodil Ström Holst^{1*}

Abstract

Background: Heterophilic antibodies in serum and plasma can interfere with mammalian antibodies in immunoassays and result in false test results, usually false positive. Although studies screening for heterophilic antibodies as well as elimination studies have been conducted in dogs and cats, knowledge of the presence of heterophilic antibodies in other species in veterinary medicine is limited. In this study, a 2-site sandwich-type interference assay that detects anti-mouse antibodies was used to detect heterophilic antibodies in a population of horses treated in an animal hospital.

Results: A total of 194 serum samples from 127 individual horses were analyzed. There were 11/127 (8.7%) interference-positive horses, and these were analyzed in an assay exchanging the capture mouse IgG with chicken IgY. The positive samples were negative in the chicken IgY assay, indicating elimination of a possible interference, with the chicken-based assay. Four interference-positive samples were from geldings, and anti-Müllerian hormone (AMH) was analyzed from these samples. AMH concentrations were negative in these samples as expected in geldings, indicating that the heterophilic antibodies did not cause interference in the AMH assay.

Conclusion: The present study shows that there are heterophilic antibodies in horse serum samples like in samples from humans, dogs, and cats. The use of chicken-based reagents, such as chicken lgY, which do not cross-react with mammalian lgG, eliminates the effects of interfering antibodies in the samples. Equine heterophilic antibodies do not necessarily cause interference in commercial immunoassays.

Keywords: ELISA, Heterophilic antibodies, Horse, Interference, Serum

Background

Immunoassays are commonly used in veterinary clinical practice, especially for hormone analyses, and provide support for clinical diagnosis and treatment. One assay that is commonly used is the sandwich immunoassay [1]. This assay has the advantage of being very sensitive, but

it is also prone to interference by heterophilic antibodies [1]. Heterophilic antibodies can cross-link capture antibodies with detection antibodies and have been shown to cause false-positive results in human medicine [2–7], for anti-Müllerian hormone (AMH) and B-type natriuretic hormone in dogs [8, 9], and for equine growth hormone (eGH) in horses [10, 11]. In human medicine, heterophilic antibodies can be grouped as true heterophilic antibodies, human antimouse antibodies (HAMA) and rheumatoid factors (RF) [1, 12].

Full list of author information is available at the end of the article



^{*}Correspondence: Bodil.Strom-Holst@slu.se

¹ Department of Clinical Sciences, Swedish University of Agricultural Sciences, Box 7054, 750 07 Uppsala, Sweden

Dong et al. Acta Vet Scand (2021) 63:10 Page 2 of 6

Reported prevalences of heterophilic antibodies vary and depend on methods used. A double-antibody sandwich immunoassay that does not cross-link with any known substance can be used to screen for heterophilic antibodies. In such an assay, signals may be generated by the cross-linking of the assay antibodies by heterophilic antibodies [1]. In veterinary medicine, an interference assay was used to study the prevalence of heterophilic antibodies in the serum of dogs and cats, and the prevalence was reported to be 5–9% [13]. In horses, it has been reported to be 5% [11]. In humans, it has been reported to be as high as 40% [14]. The frequency of interference in human serum samples has been reported to be from 0.5 to 2% to around 4% [15, 16]. It will vary with the assay used but will be lower than the prevalence of heterophilic antibodies [8, 16].

Most previous reports on screening and elimination of interfering antibodies in veterinary clinical laboratories have focused on dogs and cats [9, 13, 15]. In the horse, abnormally high concentrations of eGH analyzed using an in-house enzyme-linked immunosorbent assay (ELISA) have been described to be caused by heterophilic antibodies, and a screening revealed a presence of heterophilic antibodies in 5% of serum samples from healthy horses [10, 11].

Heterophilic antibodies are a heterogeneous group, and multiple strategies are required to eliminate their effect on assay results [17]. One approach is taking advantage of the fact that heterophilic antibodies against mammalian IgG do not cross-react with chicken IgY. The exchange of mouse IgG with chicken IgY has therefore been shown to eliminate the interference of heterophilic antibodies in human samples as well as in samples from dogs and cats [13, 18].

The goals of this study were to use a previously developed species-independent interference assay to screen a population of horses treated in animal hospitals for presence of heterophilic antibodies, to assess whether chicken IgY-based tests eliminate interference and if detected heterophilic antibodies cause interference in a commercial sandwich immunoassay for analysis of AMH.

Methods

Animals

Equine serum that had been analyzed at the Clinical Pathology Laboratory, the University Animal Hospital in Uppsala, Swedish University of Agricultural Sciences, Sweden was used. Exclusion criteria were clearly visible signs of hemolysis or lipemia.

Interference assay

An interference assay was performed as described by Bergman and co-workers [13]. Negative samples from trial runs were pooled and used as negative controls. The chicken anti-mouse IgG was diluted 1:1600 in the pooled negative control sera and this solution represented the positive control. As standard, chicken anti-mouse IgG (Immunsystem AB, Uppsala, Sweden) was used, diluted in phosphate-buffered saline (PBS) 1:1600. The cutoff point was determined by calculating the mean optical density (OD) of the duplicates of a serial dilution seven times of the standard. For each run, the cutoff point was required to be greater than the assay limit of detection (LoD), as determined by the formula LoD = 0-standard + 2 standard deviations (SD) (doing 25 repeat measurements of the 0-standard). The intra-assay coefficients of variation (CVs) for the cutoff point ranged from 0.88 to 6.81%, and the inter-assay CV was 22.1%. A relative OD value for each sample was calculated by dividing the mean OD of the sample by the cutoff value. The cutoff limit for a positive result was > 1. Similarly, the positive control was required to be > 1 and the negative control < 1.

In brief, 100 μL purified mouse IgG (2 μg/mL) was added to each well of a polystyrene microtiter plate and placed overnight at 4 °C. After washing with PBS, a total volume of 50 µL of standards, controls and samples were pipetted in duplicate into the wells followed by incubation at room temperature (RT). After another washing, monoclonal horseradish peroxidase (HRP)conjugated mouse anti-human carcinoembryonic antigen (CEA) IgG was added, selected not to match the specificity of the capture antibody, thereby excluding the possibility that any existing analyte would give rise to a positive result. The plate was then incubated and washed. Finally, 3,3',5,5' tetramethylbenzidinen (TMB) was added and the plate was incubated for 8 min at RT in the dark before adding stop solution (H2SO4) and reading OD at 450 nm.

Interference elimination

Samples that were positive in the interference assay were further tested in a chicken-based assay, as previously described [13]. In brief, every second column of a microtiter plate was coated with either 2 μ g/mL of purified mouse IgG or 2 μ g/mL nonimmunized chicken IgY (Immunsystem AB). 50 μ L of standards, controls and samples were added into each well, and incubated at RT. After washing, 100 μ L of a 1:10,000 dilution of monoclonal HRP-conjugated mouse anti-human CEA IgG was added, followed by incubation at RT. The plate was washed, TMB was added and the plate was incubated at RT in the dark before adding stop solution and reading OD at 450 nm.

Dong et al. Acta Vet Scand (2021) 63:10 Page 3 of 6

Anti-Müllerian hormone assay

Samples from geldings that were positive in the interference assay were analyzed for AMH using an ELISA (AMH Gen II, Beckman coulter), according to the manufacturer. Briefly, 30 µL of standards, controls and samples were mixed with 150 µL assay buffer, and 120 µL of the so achieved solution were pipetted in duplicate into each well and incubated in an anti-AMH antibody coated microtitration plate. After incubation and washing, anti-AMH biotin conjugate was added, and after a second incubation and washing step, streptavidin-HRP was added. After a third incubation and washing step, the substrate, TMB, was added and incubated briefly before adding an acidic stopping solution. The degree of enzymatic turnover of the substrate was determined by dual wavelength absorbance measurement at 450 nm and 620 nm. The intra-assay CV was < 5% and the inter-assay CV was < 15%.

Statistical analysis

Routine descriptive statistical methods were used. The effects of sex, breed, castration status and diagnosis were compared between the positive and the negative group using Student's t-test. In the study of interference elimination, the Wilcoxon signed-rank test was used to compare results using mouse IgG with those using chicken IgY. The level of statistical significance was set at P < 0.05. The GraphPad Prism Software version 8.0 for Windows (San Diego, CA, USA) was used for statistical analyses.

Results

Horses

A total of 194 samples from 127 individual horses were analyzed. The horses were of 27 different breeds. Their median age was 8 years, interquartile range (IQR) 0.5–13.0 years, and included 59 mares, 22 stallions and 46 geldings. Forty-four horses were subjected to multiple

analyses: 27 horses had samples analyzed twice, 10 horses three times, 5 four times and two horses had samples analyzed 5 times. All multiple samples were collected at different times. The horses were classified into 9 different categories according to diagnosis or cause for initial sampling (Table 1).

Interfering antibodies in horses

There were 11/127 (8.7%) interference-positive horses, with a median age of 12.0 years, range 1 month–17 years, (IQR 4.0–15.0 years), seven mares and four geldings of seven different breeds (Fig. 1).

The relative OD for the positive samples ranged between 1.01 and 4.69. A 15-year-old Icelandic horse gelding admitted to the animal hospital for liver disease was sampled three times within two weeks and was only positive for interference the second time (relative OD 1.2). A 14-year-old intact female Icelandic horse with a subcutaneous skin abscess was positive for interference

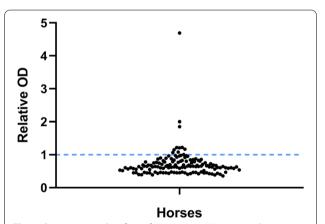


Fig. 1 Screening results of interference assay. Horse samples (n = 127) were screened with a sandwich ELISA using nonimmunized mouse IgG as the capture antibody. Cutoff level is indicated by the dashed horizontal line

Table 1 Characterization of the horses screened for interfering antibodies by disease category

Disease category	Positive samples (n = 11)	%	Negative sampl	es (n = 116) %	Total (n = 127)
Bone, muscle and joint	4	21.1	15	78.9	19
Digestive tract	2	3.9	49	96.1	51
Infection/inflammation	2	9.5	19	90.5	21
Neoplastic disease	1	16.7	5	83.3	6
Neurologic disease	0	0	4	100	4
Reproductive disease	0	0	5	100	5
Respiratory disease	0	0	3	100	3
Skin disease	1	10	9	90	10
Other	1	16.7	5	83.3	6
SUM	11	8.7	116	91.3	127

Dong et al. Acta Vet Scand (2021) 63:10 Page 4 of 6

the first sample (relative OD 1.03), but was negative on three subsequent samplings within three weeks. There was no significant difference of age, sex, breed, neutering status, or diagnostic category between interference-positive and interference-negative horses (P > 0.05).

Interference elimination

When chicken IgY antibody was used as the capture antibody, all the 11 interference-positive samples were negative. When mouse IgG was used, five samples were also detected negative, although they had previously been evaluated positive in screening tests (Fig. 2). The 11 samples produced significantly less signal in the wells containing chicken IgY than in wells containing mouse IgG (Z=-3.059, P<0.01).

Analysis of anti-Müllerian hormone

Enough serum for AMH analysis was available from four geldings (relative OD 1.1; 1.2; 1.2 and 1.8), all of which had AMH concentrations below the lowest standard point.

Discussion

In the present study, equine serum samples were screened for heterophilic antibodies and the prevalence was found to be 8.7% (11/127). This is similar but slightly higher than previously reported for horses [11], and similar to what has been reported for cats and dogs, using the same species-independent ELISA [13]. The horse population studied by Borromeo and co-workers [11] was healthy horses. The samples included in the present study were from horses admitted to an equine hospital, and the

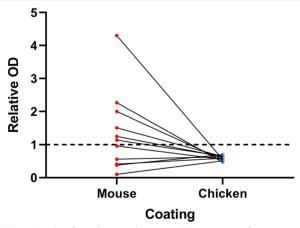


Fig. 2 Results of interference elimination. The reactivity of heterophilic antibodies towards unimmunized mouse IgG and chicken IgY was detected with HRP- conjugated mouse anti-CEA antibody. The positive cutoff level is indicated by a dotted horizontal line

results may therefore better mirror the situation in clinical practice. Variations in prevalence of heterophilic antibodies may also be related to the methods used.

In human medicine, it has been shown that interference in immunoassays is a present and underestimated problem [16, 19]. The result achieved when affected by heterophilic antibodies is usually falsely increased, resulting in overestimation of the analytes in question. Examples of analyses that may be affected in equine medicine are e.g. adrenocorticotrophic hormone (ACTH) and insulin, used in the diagnosis of the two most common endocrine syndromes: pituitary pars intermedia dysfunction (PPID) and equine metabolic syndrome (EMS) [20], and erythropoietin, used for doping [21]. If the test results are clearly above the normal range, interference by heterophilic antibodies may be suspected, as was the case of interference with the in-house eGH ELISA [10]. In other cases, a result that is false high but within the normal range may not raise such a suspicion but lead to an incorrect diagnosis. Evaluating the effect of heterophilic antibodies on immunoassays is challenging, as the true result most often is unknown. Reanalyzing the sample after precipitation with polyethylene glycol (PEG) can result in normalized values in humans [22], but the method does not work well for samples from dogs [8]. One assay with a known expected result for castrated animals is the AMH immunoassay. In horses, this assay can be used to differ geldings, having non-detectable AMH concentrations, from cryptorchid stallions, having high AMH concentrations [23]. It is a sandwich immunoassay using murine antibodies, and thus similar to the interference assay. In the present study, AMH was analyzed in samples from four geldings with heterophilic antibodies. All samples had AMH concentrations below the detection limit, thus no sign of interference was detected. One possible reason for the lack of interference, leading to false positive results, is that the concentrations of heterophilic antibodies were just above the cut-off value, indicating a weak reactivity. Since antibodies affect immunoassay results both by their concentration and their affinity, it is also possible that the affinity of the heterophilic antibodies was too low to cause interference in the AMH assay [1]. In dogs, 2/7 samples with heterophilic antibodies caused a false positive result in the AMH assay, and these were the samples with the strongest anti-mouse reactivity. As commercial assays generally are designed to minimize the effect of heterophilic antibodies, the problem with interference is expected to increase with stronger reactivity.

That the positive results in the interference assay were caused by heterophilic antibodies is supported by the fact that they were eliminated with chicken IgY as capture antibodies, and that samples with visible signs of

Dong et al. Acta Vet Scand (2021) 63:10 Page 5 of 6

hemolysis and lipemia were excluded. In addition, all samples were cryopreserved at $-20\,^{\circ}\text{C}$ to eliminate interference by complement [24]. Clinically relevant interference has previously been described with an in-house ELISA for eGH using mouse and rabbit antibodies [11]. The lack of interference with chicken IgY suggests that exchanging mammalian IgG with chicken IgY in immunoassays may reduce the risk of interference by heterophilic antibodies, as has previously been described [18].

Iatrogenic interference has been reported in human medicine as a result of monoclonal antibody therapy [25], but as no therapeutic antibodies are licensed for use in horses this is not a possible cause for the heterophilic antibodies in the present study. Direct physical contact with other species and indirect exposure through vaccination and food intake has also been suggested to potentially induce heterophilic antibodies [16]. Rheumatoid factors are considered as one of the possible causes of the interference, with age, gender and autoimmunity being main risk factors for interference caused by RF in human medicine [15, 26, 27]. In horses, there have been reports of an age-dependent increase in IgM RF regardless of sex [28]. Anamnestic data from horses with positive samples were reviewed to identify possible risk factors for interference. However, we found no significant relationship between the presence of heterophilic antibodies and age, sex, breed, castration status, or diagnosis.

Of the 11 positive horses, three were sampled and assayed on multiple occasions and tested positive only once, with less than a month between the different samples. A transient character and low affinity of heterophilic antibodies has been described [29, 30], although they have been described to persist for several years in dogs [31]. The reactivity of these samples was rather weak, and a variation between different runs may contribute to the shifting results.

Conclusions

The prevalence of heterophilic antibodies in a hospital population of horses was similar to that previously reported for dogs and cats. Exchanging mammalian IgG with chicken IgY can reduce the problem of interference in immunoassays for equine samples.

Abbreviations

ACTH: Adrenocorticotrophic hormone; AMH: Anti-Müllerian hormone; CEA: Carcinoembryonic antigen; CV: Coefficient of variation; eGH: Equine growth hormone; ELISA: Enzyme-linked immunosorbent assay; EMS: Equine metabolic syndrome; HAMA: Human antimouse antibodies and; HRP: Horseradish peroxidase; IgG: Immunoglobulin G; IgM: Immunoglobulin M; IgY: Immunoglobulin Y; IQR: Interquartile range; LoD: Limit of detection; OD: Optical density; PBS: Phosphate buffered saline; PEG: Polyethylene glycol; PPID: Pituitary pars intermedia dysfunction; RF: Rheumatoid factor; RT: Room temperature; TMB: 3,3′,5,5′-Tetramethylbenzidine.

Acknowledgements

The authors acknowledge the laboratory staff at the University Animal Hospital.

Authors' contributions

BD analyzed the samples, interpreted the results and drafted the manuscript. DB supervised the laboratory work and contributed in writing the manuscript. BSH supervised the study and took a major part in writing the manuscript. All authors read and approved the final manuscript.

Funding

Open access funding provided by Swedish University of Agricultural Sciences. This study was supported by Fujian natural science foundation general project (2019J01804). The funding body had no role in the design of the study and collection, analysis, and interpretation of data or in writing the manuscript.

Availability of data and materials

All data generated or analyzed during this study are included in this published article.

Declarations

Ethics approval and consent to participate

The study involves analyses of left-over serum samples that were collected for diagnostic purposes, and no ethical approval was needed. All horse owners gave their informed written consent. The study was performed in accordance with national animal welfare regulations (SJVFS 2015:38).

Consent for publication

Not applicable.

Competing interests

The authors declare that they have no competing interests.

Author details

¹ Department of Clinical Sciences, Swedish University of Agricultural Sciences, Box 7054, 750 07 Uppsala, Sweden. ² College of Life Science of Longyan University, Longyan 364012, China.

Received: 17 December 2020 Accepted: 3 March 2021 Published online: 12 March 2021

References

- Bolstad N, Warren DJ, Nustad K. Heterophilic antibody interference in immunometric assays. Best Pract Res Clin Endocrinol Metab. 2013:27:647–61.
- Bolstad N, Warren DJ, Bjerner J, Kravdal G, Schwettmann L, Olsen KH, et al. Heterophilic antibody interference in commercial immunoassays; a screening study using paired native and pre-blocked sera. Clin Chem Lab Med. 2011;49:2001–6.
- Censi S, Cavedon E, Fernando SW, Barollo S, Bertazza L, Zambonin L, et al. Calcitonin measurement and immunoassay interference: a case report and literature review. Clin Chem Lab Med. 2016:54:1861–70.
- Mongolu S, Armston AE, Mozley E, Nasruddin A. Heterophilic antibody interference affecting multiple hormone assays: is it due to rheumatoid factor? Scand J Clin Lab Investig. 2016;76:240–2.
- García-González E, Aramendía M, Álvarez-Ballano D, Trincado P, Rello L. Serum sample containing endogenous antibodies interfering with multiple hormone immunoassays laboratory strategies to detect interference. Prac Lab Med. 2016;4:1–0.
- Revet I, Boesten LS, Linthorst J, Yildiz E, Janssen JW, de Rijke YB, et al. Misleading FT4 measurement: assay-dependent antibody interference. Biochem Med. 2016;26:436–43.
- Rotmensch S, Cole LA. False diagnosis and needless therapy of presumed malignant disease in women with false-positive human chorionic gonadotropin concentrations. Lancet. 2000;355:712–5.
- Bergman D, Larsson A, Hansson-Hamlin H, Strom HB. Investigation of interference from canine anti-mouse antibodies in hormone immunoassays. Vet Clin Pathol. 2019;48(Suppl 1):59–69.

Dong et al. Acta Vet Scand (2021) 63:10 Page 6 of 6

- Solter PF, Oyama MA, Sisson DD. Canine heterophilic antibodies as a source of false-positive B-type natriuretic peptide sandwich ELISA results. Vet Clin Pathol. 2008;37:86–95.
- Borromeo V, Berrini A, Gaggioli D, Secchi C. Heterophile antibodies in horse plasma interfering in an enzyme-linked immunosorbent assay. Vet Res Commun. 2007;31(Suppl 1):165–8.
- Borromeo V, Berrini A, Gaggioli D, Secchi C. Heterophile antibody interference in a solid phase sandwich immunoassay for detection of equine growth hormone in plasma. Vet Immunol Immunopathol. 2007;115:1–9.
- Kroll MH, Elin RJ. Interference with clinical laboratory analyses. Clin Chem. 1994;40:1996–2005.
- Bergman D, Larsson A, Hansson-Hamlin H, Svensson A, Holst BS. Prevalence of interfering antibodies in dogs and cats evaluated using a species-independent assay. Vet Clin Pathol. 2018;47:205–12.
- Boscato LM, Stuart MC. Incidence and specificity of interference in twosite immunoassays. Clin Chem. 1986;32:1491–5.
- 15. Bjerner J, Nustad K, Norum LF, Olsen KH, Bormer OP. Immunometric assay interference: incidence and prevention. Clin Chem. 2002;48:613–21.
- Ismail AA, Walker PL, Cawood ML, Barth JH. Interference in immunoassay is an underestimated problem. Ann Clin Biochem. 2002;39:366–73.
- Bergman D, Larsson A, Hansson-Hamlin H, Ahlen E, Holst BS. Characterization of canine anti-mouse antibodies highlights that multiple strategies are needed to combat immunoassay interference. Sci Rep. 2019;9:1–9.
- Larsson A, Mellstedt H. Chicken antibodies: a tool to avoid interference by human anti-mouse antibodies in ELISA after in vivo treatment with murine monoclonal antibodies. Hybridoma. 1992;11:33–9.
- Mongolu S, Armston AE, Mozley E, Nasruddin A. Heterophilic antibody interference affecting multiple hormone assays: is it due to rheumatoid factor? Scand J Clin Lab Invest. 2006;76:240–2.
- 20. McFarlane D. Diagnostic testing for equine endocrine diseases: confirmation versus confusion. Vet Clin N Am Equine Pract. 2019;35:327–38.
- Lönnberg M, Bondesson U, Cormant F, Garcia P, Bonnaire Y, Carlsson J, et al. Detection of recombinant human EPO administered to horses using MAIIA lateral flow isoform test. Anal Bioanal Chem. 2012;403:1619–28.
- Paragliola RM, Corsello A, Papi G, Melfa E, Urbani A, Pontecorvi A, et al. Immunoassay interference on thyroid function tests during treatment with nivolumab. Thyroid. 2020;30:1091–4.

- Claes A, Ball BA, Almeida J, Corbin CJ, Conley AJ. Serum anti-Mullerian hormone concentrations in stallions: developmental changes, seasonal variation, and differences between intact stallions, cryptorchid stallions, and geldings. Theriogenology. 2013;79:1229–35.
- Carlander D, Larsson A. Avian antibodies can eliminate interference due to complement activation in ELISA. Ups J Med Sci. 2001;106:189–95.
- Reinsberg J. Interferences with two-site immunoassays by human anti-mouse antibodies formed by patients treated with monoclonal antibodies: comparison of different blocking reagents. Clin Chem. 1998:44:1742–4
- 26. van Vollenhoven RF. Sex differences in rheumatoid arthritis: more than meets the eye. BMC Med. 2009;7:1–4.
- Rasch EK, Hirsch R, Paulose-Ram R, Hochberg MC. Prevalence of rheumatoid arthritis in persons 60 years of age and older in the United States: effect of different methods of case classification. Arthritis Rheum. 2003;48:917–26.
- Hoikhman R, Kudlackova H, Babak V, Faldyna M, Jahn P. Detection of IgM-rheumatoid factor and anti-citrullinated protein antibodies in healthy horses and their comparison. Vet ImmunolImmunopathol. 2018;202:141–6.
- Kazmierczak SC, Catrou PG, Briley KP. Transient nature of interference effects from heterophile antibodies: examples of interference with cardiac marker measurements. Clin Chem Lab Med. 2000;38:33–9.
- 30. Bjerner J, Børmer OP, Nustad K. The war on heterophilic antibody interference. Clin Chem. 2005;51:9–11.
- 31. Bergman D, Bäckström C, Hansson-Hamlin H, Larsson A, Holst BS. Preexisting canine anti-IgG antibodies: implications for immunotherapy, immunogenicity testing and immunoassay analysis. Sci Rep. 2020;10:1–0.

Publisher's Note

Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Ready to submit your research? Choose BMC and benefit from:

- fast, convenient online submission
- $\bullet\,$ thorough peer review by experienced researchers in your field
- rapid publication on acceptance
- support for research data, including large and complex data types
- gold Open Access which fosters wider collaboration and increased citations
- maximum visibility for your research: over 100M website views per year

At BMC, research is always in progress.

Learn more biomedcentral.com/submissions

