

Baptist Health South Florida

Scholarly Commons @ Baptist Health South Florida

All Publications

5-2021

New Class of Radially Adjustable Stentriever for Acute Ischemic Stroke: Primary Results of the Multicenter Tiger Trial

Guilherme Dabus

Baptist Hospital of Miami; Miami Neuroscience Institute; Miami Cardiac & Vascular Institute,
guilhermed@baptisthealth.net

Italo Linfante

Baptist Hospital of Miami; Miami Neuroscience Institute; Miami Cardiac & Vascular Institute,
italol@baptisthealth.net

Follow this and additional works at: <https://scholarlycommons.baptisthealth.net/se-all-publications>

Citation

Stroke (2021) 1534-1544

This Article -- Open Access is brought to you for free and open access by Scholarly Commons @ Baptist Health South Florida. It has been accepted for inclusion in All Publications by an authorized administrator of Scholarly Commons @ Baptist Health South Florida. For more information, please contact Carrief@baptisthealth.net.



New Class of Radially Adjustable Stentriever for Acute Ischemic Stroke

Primary Results of the Multicenter TIGER Trial

Rishi Gupta¹, MD, MBA*; Jeffrey L. Saver, MD*; Elad Levy, MD; Osama O. Zaidat, MD, MS; Dileep Yavagal², MD; David S. Liebeskind, MD; Ahmad Khaldi, MD, MS; Bradley Gross, MD; Michael Lang, MD; Sandra Narayanan, MD; Brian Jankowitz, MD; Kenneth Snyder, MD; Adnan Siddiqui, MD; Jason Davies, MD; Eugene Lin³, MD; Ameer Hassan, MD; Ricardo Hanel, MD; Amin Aghaebrahim, MD; Ritesh Kaushal, MD; Ali Malek⁴, MD; Nils Mueller-Kronast, MD; Robert Starke, MD; Hormozd Bozorgchami, MD; Gary Nesbit, MD; Masahiro Horikawa, MD; Ryan Priest⁵, MD; Jesse Liu⁶, MD; Ronald F. Budzik, MD; Peter Pema, MD; Nirav Vora, MD; M. Asif Taqi, MD; Edgar Samaniego, MD; Qingliang Tony Wang, MD, PhD; Erez Nossek, MD; Guilherme Dabus⁷, MD; Italo Linfante, MD; Ajit Puri, MD; Eitan Abergel, MD; Sidney Starkman, MD; Satoshi Tateshima, MD; Ashutosh P. Jadhav, MD, PhD, for the TIGER Trial Investigators

BACKGROUND AND PURPOSE: The Tigertriever is a novel, radially adjustable, fully visible, stentriever that permits the operator to align radial expansion with target vessel diameters. This multicenter trial compared the Tigertriever's effectiveness and safety compared with established stent retrievers.

METHODS: Single arm, prospective, multicenter trial comparing the Tigertriever to efficacy and safety performance goals derived from outcomes in 6 recent pivotal studies evaluating the Solitaire and Trevo stent-retriever devices with a lead-in and a main-study phase. Patients were enrolled if they had acute ischemic stroke with National Institutes of Health Stroke Scale score ≥ 8 due to large vessel occlusion within 8 hours of onset. The primary efficacy end point was successful reperfusion, defined as core laboratory-adjudicated modified Thrombolysis in Cerebral Ischemia score 2b-3 within 3 passes of the Tigertriever. The primary safety end point was a composite of 90-day all-cause mortality and symptomatic intracranial hemorrhage. Secondary efficacy end points included 3-month good clinical outcome (modified Rankin Scale score 0-2) and first-pass successful reperfusion.

RESULTS: Between May 2018 and March 2020, 160 patients (43 lead-in, 117 main phase) at 17 centers were enrolled and treated with the Tigertriever. The primary efficacy end point was achieved in 84.6% in the main-study phase group compared with the 63.4% performance goal and the 73.4% historical rate (noninferiority $P < 0.0001$; superiority $P < 0.01$). The first pass successful reperfusion rate was 57.8%. After all interventions, successful reperfusion (modified Thrombolysis in Cerebral Ischemia score $\geq 2b$) was achieved in 95.7% and excellent reperfusion (modified Thrombolysis in Cerebral Ischemia score 2c-3) in 71.8%. The primary safety composite end point rate of mortality and symptomatic intracranial hemorrhage was 18.1% compared with the 30.4% performance goal and the 20.4% historical rate (noninferiority $P = 0.004$; superiority $P = 0.57$). Good clinical outcome was achieved in 58% at 90 days.

CONCLUSIONS: The Tigertriever device was shown to be highly effective and safe compared with Trevo and Solitaire devices to remove thrombus in patients with large-vessel occlusive stroke eligible for mechanical thrombectomy.

REGISTRATION: URL: <https://www.clinicaltrials.gov>; Unique identifier: NCT03474549.

GRAPHIC ABSTRACT: An online [graphic abstract](#) is available for this article.

Key Words: goal ■ intracranial hemorrhage ■ laboratories ■ reperfusion ■ stent

Correspondence to: Rishi Gupta, MD, MBA, Wellstar Medical Group, Department of Neurosurgery, Wellstar Health System Kennestone Hospital Marietta, GA. Email rishi.gupta@wellstar.org

*R. Gupta and J.L. Saver are co-first authors.

This manuscript was sent to Gert Kwakkel, Guest Editor, for review by expert referees, editorial decision, and final disposition.

The Data Supplement is available with this article at <https://www.ahajournals.org/doi/suppl/10.1161/STROKEAHA.121.034436>.

For Sources of Funding and Disclosures, see page 1542.

© 2021 The Authors. *Stroke* is published on behalf of the American Heart Association, Inc., by Wolters Kluwer Health, Inc. This is an open access article under the terms of the [Creative Commons Attribution Non-Commercial-NoDerivs License](#), which permits use, distribution, and reproduction in any medium, provided that the original work is properly cited, the use is noncommercial, and no modifications or adaptations are made.

Stroke is available at www.ahajournals.org/journal/str

Nonstandard Abbreviations and Acronyms

BGC	balloon guide catheter
CT	computed tomography
DSMB	data and safety monitoring board
IQR	interquartile range
IV tPA	intravenous tissue-type plasminogen activator
mRS	modified Rankin Scale
mTICI	modified Thrombolysis in Cerebral Ischemia
NIHSS	National Institutes of Health Stroke Scale
PG	performance goal
sICH	symptomatic intracranial hemorrhage
TIGER	Treatment With Intent to Generate Endovascular Reperfusion

Based on convergent evidence from major clinical trials of reperfusion efficacy and improved clinical outcomes, endovascular thrombectomy with stent retrievers has become the standard of care for patients with large vessel occlusion ischemic strokes.^{1–8} A meta-analysis of 5 randomized trials demonstrated the number needed to treat with thrombectomy is 1 in 2.6 to improve 3-month global disability outcome by 1 or more levels on the modified Rankin Scale (mRS).⁹ However, while quite substantial, the benefits conveyed by first-generation stent retrievers are constrained by less than maximal reperfusion rates achieved with these devices. In pooled, individual participant data meta-analyses of the pivotal trials, failure to achieve successful reperfusion (modified Thrombolysis in Cerebral Ischemia [mTICI] 2b–3) occurred in 29% of patients and failure to achieve complete reperfusion (mTICI 3) in 67%.⁹ Accordingly, developing additional endovascular thrombectomy devices with performance characteristics comparable or better than established stent retrievers is desirable.

The TIGER (Rapid Medical) is a novel operator-adjustable stent retriever that affords the interventionalist incremental control over the radial diameter and radial force of the thrombectomy basket. The design is intended to facilitate alignment of the mesh with the anatomy of the occluded vessel and to increase internalization of the thrombus within the device, thereby facilitating retrieval and reducing downstream embolization (Figure 1A; Figure I in the [Data Supplement](#)). The device has had CE mark, the European Union certification indicating conformity with health, safety, and environmental protection, since 2016 and has shown promising signals of efficacy and safety in preliminary case series reported by European centers.^{10–12}

The TIGER trial (Treatment With Intent to Generate Endovascular Reperfusion) was a single-arm, multicenter, prospective study assessing the efficacy and

safety of the Tiger-21 and Tiger-17 retrievers. The study employed an objective performance criterion, noninferiority design, comparing efficacy and safety of the Tiger-trieter with performance goals derived from adjudicated outcomes in 6 completed prospective trials of 2 FDA-approved predicate stent retrievers, Solitaire (Medtronic) and Trevo (Stryker).^{1–4,13,14}

METHODS

Study Design

The data that support the findings of this study are available from the corresponding author on reasonable request and with approval from the TIGER investigators. The primary TIGER Study hypothesis was that the Tigertriever device would achieve successful reperfusion performance within 3 attempts would be superior to the performance goal established using Bayesian meta-analysis of 6 recent pivotal trials of the Solitaire and Trevo stent-retriever devices—TREVO 2 (Thrombectomy Revascularization of Large Vessel Occlusions in Acute Ischemic Stroke),¹³ SWIFT (SOLITAIRE FR With the Intention for Thrombectomy),¹⁴ MR CLEAN (Multicenter Randomized Clinical Trial of Endovascular Treatment for Acute Ischemic Stroke in the Netherlands),¹ ESCAPE (Endovascular Treatment for Small Core and Anterior Circulation Proximal Occlusion With Emphasis on Minimizing CT to Recanalization Times),⁴ REVASCAT (Randomized Trial of Revascularization With Solitaire FR Device Versus Best Medical Therapy in the Treatment of Acute Stroke due to Anterior Circulation Large Vessel Occlusion Presenting Within Eight Hours of Symptom Onset),³ and SWIFT PRIME (Solitaire With the Intention for Thrombectomy as Primary Endovascular Treatment).²

The study was overseen by an executive committee, composed of academic investigators, and an independent data and safety monitoring board (DSMB) (Table I in the [Data Supplement](#)). The DSMB consisted of an independent stroke neurologist and 2 independent neurointerventionalists, who did not enroll patients, and a biostatistician. At the level of individual events, the DSMB adjudicated all safety end points, including hemorrhages and serious adverse events. At the level of aggregate study data, the DSMB acted in an expert supervisory capacity and monitored subject safety and the conduct of the study. Central readers at a core imaging laboratory assessed qualifying diagnostic catheter angiography to determine location of the target occlusion, and rated reperfusion grades on angiograms obtained after every pass of the Tigertriever device, every pass of rescue therapy, and end of the procedure. The core laboratory imaging readers also assessed baseline and 24-hour post-procedure computed tomography (CT) or magnetic resonance imaging scans to verify study entry criteria and determine the presence and type of any postprocedure intracranial hemorrhage. The sponsor, Rapid Medical, and its contracted Contract Research Organization (Genae Americas) were responsible for logistical operations, data management, and monitoring of the trial. The study protocol was approved by the institutional review board or ethics committee at each participating site.

Population and Participating Centers

All patients or their legally authorized representative provided written informed consent before enrollment. Patients aged 18

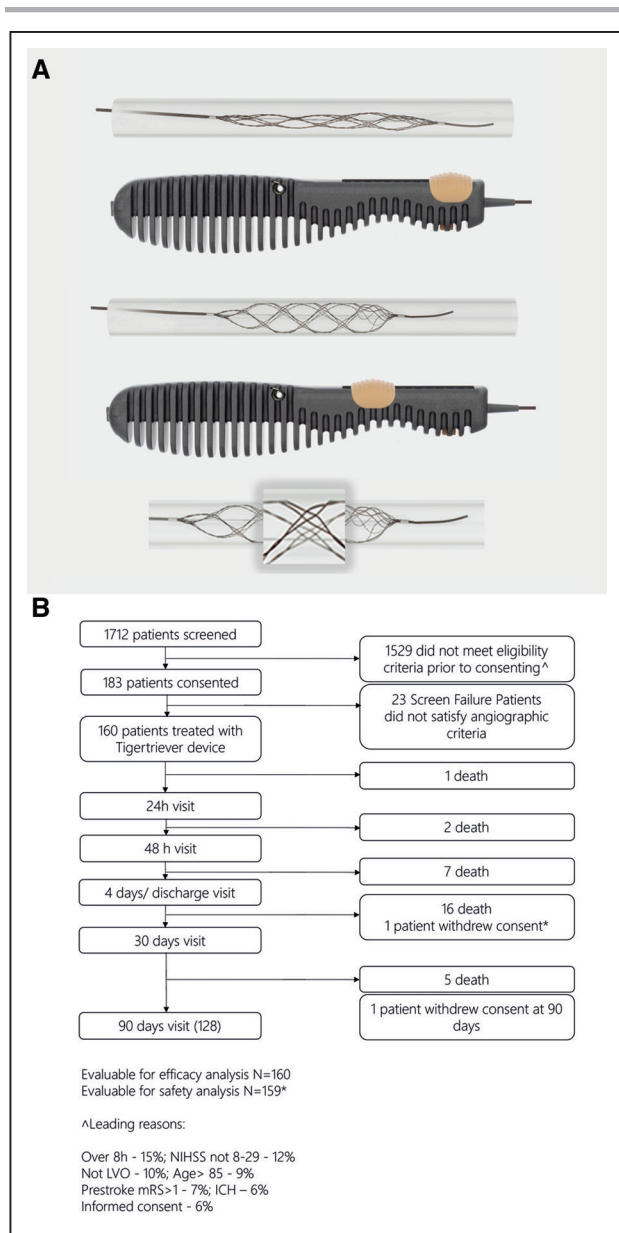


Figure 1. Tigertriever device and chart of patient flow through study.

A, Tigertriever device including external handle for the operator to control the degree of expansion of the stent retriever. Fine wire mesh mounted on a flexible shaft. The design of the wire mesh (magnified view) is optimized to penetrate the clot and encapsulate it during retrieval. **B**, Patient screening, enrollment, treatment, and follow-up. ICH indicates intracranial hemorrhage; LVO, large vessel occlusion; mRS, modified Rankin Scale; and NIHSS, National Institutes of Health Stroke Scale.

to 85 years were eligible for the study if they had acute ischemic stroke with new moderate-to-severe neurological deficits, angiographically confirmed occlusion in the intracranial internal carotid, middle cerebral (M1 or M2 segment), basilar, or vertebral artery, and could undergo endovascular therapy (at least 1 Tigertriever deployment in the target artery) within 8 hours of last known well. Full study entry criteria are shown in the [Data Supplement](#) (Table II in the [Data Supplement](#)). Key inclusion criteria included the following: prestroke mRS score ≤ 1 ; baseline National Institutes of Health Stroke Scale (NIHSS) score 8–29;

Alberta Stroke Program Early CT Score 6–10 for CT-qualifying or diffusion restriction volume ≤ 50 mL for magnetic resonance imaging-qualifying anterior circulation cases; patient being either (1) eligible for, and received, IV tPA (intravenous tissue-type plasminogen activator) within 3 hours of symptom onset or (2) ineligible for IV tPA treatment; and anticipated device deployment within 8 hours of last known well. Key exclusion criteria included stenosis or occlusion in the deployment site or in a proximal vessel that would prevent device access to the thrombus or require angioplasty or stenting to enable device access to the thrombus. The treatment start inclusion criterion of within 8 hours, rather than also 8 to 24 hours with advanced imaging selection, was selected as initial trial design discussions and review with the Food and Drug Administration occurred before the reporting of results of late window trials. Physicians were trained in the use of the device on a bench vascular model before any procedures were done. In addition, participating sites that had no experience with the Tigertriever device were required to first participate in a lead-in phase before participating in the main-study phase for primary end point analysis. During the lead-in phase, up to 4 patients per clinical site were enrolled. The lead-in phase was completed by the site either (1) achieving 2 successive successful reperfusions (mTICI 2b or higher) or (2) performing 4 cases. Patients enrolled in the lead-in phase were consented, treated, and followed according to the clinical protocol, in the same manner as the main-study patients.

Procedure

Site of arterial access and anesthesia mode were based on operator preference at the enrolling center. The selection of guide catheter type (balloon guide catheter (BGC) versus non-BGC) and use of an intermediate catheter was also at the discretion of the operator. After placement of the guide catheter in the target vessel, a roadmap was constructed and the microcatheter was navigated across the clot. The Tigertriever was deployed and the amount of expansion of the device was based on vessel diameter and clot length.

Two versions of the Tigertriever device were available: The standard version (Tigertriever) has a net length of 32 mm (unexpanded form) and can expand up to 6 mm diameter, can be delivered through a microcatheter with an internal diameter of 0.021 inches. And a smaller version (Tigertriever 17) has a net length of 23 mm (unexpanded form) and can be delivered through a microcatheter with an internal diameter of 0.017 inches. It can expand up to 3 mm diameter. The device deployment, dwell time, and number of expansions/deflations were based on the instructions for use. The first pass was required to be with the Tigertriever for inclusion in the trial. Up to 3 passes of the Tigertriever were permitted to achieve successful reperfusion (mTICI, 2b–3). For each pass, a new Tigertriever device was used and the total number of passes and devices used were recorded. After 3 attempts, rescue therapy could be used by the operator. Rescue therapy could consist of use of another mechanical thrombectomy device, angioplasty, intraarterial tPA, or intracranial stenting at the discretion of the treating physician. The Tigertriever could not be used after a rescue attempt and for primary efficacy analysis when revascularization was unsuccessful after Tigertriever attempts it was considered a failure to achieve the end point. Imaging with CT or magnetic resonance imaging was obtained at 24 hours to assess for intracranial hemorrhage.

NIHSS scores were obtained at 24 hours, 48 hours, 4 days (or discharge if sooner), and 90 days; mRS scores were obtained at 4 days (or discharge if sooner), 30 days, and 90 days.

Outcomes

Primary Efficacy End Point

The primary efficacy end point was successful reperfusion, defined as achieving mTICI 2b-3 within 3 passes with the Tigertriever without use of rescue therapy. All angiograms were interpreted by an experienced neurointerventionalist readers at the independent core imaging laboratory, who assigned scores using the mTICI scale.¹⁵ The core imaging laboratory mTICI ratings were used for assessment of the primary efficacy end point.

Primary Safety End Point

The primary safety end point was defined as the composite of all-cause mortality at 90±14 days and symptomatic intracranial hemorrhage (sICH) within 24 (18–36) hours from the study procedure. sICH was defined as any parenchymal hematoma type 2, remote intracerebral hemorrhage, subarachnoid hemorrhage, or intraventricular hemorrhage that is the predominant cause of an NIHSS deterioration of 4 or more points at 24 hours, as adjudicated by the DSMB (modified Heidelberg Bleeding Classification criteria).¹⁶ The independent core imaging laboratory assessed all 24-hour CT/MR images and radiologically classified hemorrhages as hemorrhagic infarction types 1 or 2, parenchymal hematoma types 1 or 2, remote intracerebral hemorrhage, subarachnoid hemorrhage, and intraventricular hemorrhage.

Secondary End Points

Prespecified secondary efficacy end points for the study were (1) good clinical outcome (mRS score 0–2 at 90 days); (2) first pass successful reperfusion (mTICI 2b-3); (3) health-related quality of life at 90 days, assessed with the EQ-5D^{2,17}; and (4) granular level of disability at 90 days, assessed with the Academic Medical Center Linear Disability Score.^{18,19} Prespecified secondary safety end points were (1) asymptomatic intracranial hemorrhage within 24 hours (18–36 hours) of procedure; (2) neurological deterioration within 24 hours after procedure, defined as a NIHSS increase of 4 points or more; and (3) embolization to previously uninvolved vascular territories.

Statistical Analysis

The performance goal for the efficacy primary end point was defined as the incidence in the 6 trials cited below, minus a noninferiority statistical margin of 10%, the same clinically relevant threshold employed in the registration trials of stent-retrievers, SWIFT and TREVO 2,^{9,13} therefore, giving a performance goal (PG) of 73.4%–10%=63.4%. The pooled incidence was derived from the following 6 studies: TREVO 2, SWIFT,¹ MR CLEAN, ESCAPE, REVASCAT, and SWIFT PRIME.^{1–4,9,13} For this efficacy end point, the primary noninferiority hypothesis test was performed at an overall 2-sided alpha level of 0.05 using exact binomial methods, in which the lower confidence bound on the observed incidence of reperfusion is compared with the PG. If noninferiority was demonstrated, superiority was then tested in a post hoc hierarchical manner. The study statistical design is shown graphically in Figure IIA in the [Data Supplement](#).

The PG for safety primary end point was defined as the incidence in the 6 trials cited above, plus a statistical margin of 10%. Thus, the PG was 30.4% (18.2% mortality+2.2% sICH+10% noninferiority statistical margin). For primary safety, the hypothesis test was performed at an overall 2-sided alpha level of 0.05 using exact binomial methods, in which the upper confidence bound on the observed incidence of revascularization is compared with the PG.

The primary study population was the main-study phase group. Sensitivity analysis was performed at regulatory request on the combined lead-in/main-study phase population). To evaluate Tigertriever device performance in vessels with of varied size, the study population was divided into 2 vessel diameter categories, ≥2 and <2 mm, and performance compared across the groups.

RESULTS

Between May 2018 and March 2020, 1712 patients were preliminarily screened, 183 consented before angiography, and 160 met angiographic and device access criteria for full study enrollment. Patients were enrolled at 17 study centers (16 US centers—159 patients, 1 Israeli center—1 patient; Table I in the [Data Supplement](#)) and included 43 lead-in phase and 117 main-study phase patients. Figure 1B shows the screening, enrollment, treatment, and follow-up flow of the study population.

Patient baseline characteristics are shown for the main-study phase and combined lead-in/main-study phase cohorts in Table 1 (Table III in the [Data Supplement](#) shows patient characteristics for the lead-in phase patients alone; Tables IV and V in the [Data Supplement](#) compare patient characteristics and outcomes in the current trial with those of the 6 trials contributing to the performance goal and the more recent ARISE II trial [Analysis of Revascularization in Ischemic Stroke With EmboTrap]). In the main-study phase cohort, patient age was mean 65±15 years, 61.5% were men, and median baseline NIHSS 17 (interquartile range [IQR], 12–21). The median Alberta Stroke Program Early CT score was 9 (IQR, 8–10). Target occlusion location was in the anterior circulation in 97.4% of patients, with the most common sites the M1 middle cerebral artery (57.3%) and intracranial internal carotid artery (20.5%). Time from last known normal to puncture was 172 minutes (IQR, 128–273), and a majority (65.8%) of patients received intravenous t-PA before thrombectomy. A BGC was used in 29.9% of the procedures (21.4% BGC alone and 8.5% BGC+regional aspiration).

Primary Efficacy End Point and Reperfusion Results

In the primary analysis main-study phase population, the Tigertriever achieved the primary end point of successful reperfusion (mTICI 2b-3) within 3 passes and without rescue therapy in 99 of 117 patients, 84.6% (95% CI, 78.1%–91.2%), compared with the 63.4%

Table 1. Patient and Procedure Characteristics for Main-Study and Combined Lead-In/Main-Study Patients

	TIGER	
	Main-study (n=117)	All (n=160)
Age, y; mean (SD)	65 (15)	66 (15)
Male sex, n (%)	72 (61.5%)	92 (58%)
Race, n (%)		
White	99 (84.6%)	134 (83.7%)
Black	17 (14.5%)	24 (15%)
Asian	1 (0.85%)	2 (1.25%)
Hispanic ethnicity, n (%)	16 (13.7%)	24 (15%)
NIHSS score		
Mean (SD)	17.4 (5.6)	17.6 (5.6)
Median (IQR)	17 (12–21)	18 (13–22)
Baseline CT ASPECTS		
Mean (SD)	8.9 (1.1)	8.8 (1.2)
Median (IQR)	9 (8–10)	9 (8–10)
Prestroke mRS, n (%)	(n=117)	(n=156)
0–1	116 (99.2%)	155 (99.4%)
0	93 (79.5%)	127 (81.4%)
1	23 (19.7)	28 (17.8%)
Body mass index, median (IQR)	29.5 (25–35)	30.4 (25.7–36.2)
Medical history, n (%)		
Hypertension	89 (76%)	126 (78.8%)
Diabetes	36 (30.8%)	52 (32.5%)
Atrial fibrillation	47 (40.2)	60 (37.5%)
Dyslipidemia	57 (48.7%)	81 (50.6%)
Previous MI/CAD	23 (19.7%)	34 (21.3%)
Previous ischemic stroke/transient ischemic attack	16 (13.7%)	20 (12.5%)
Intravenous tPA failure	77 (65.8%)	111 (69.4%)
Proximal occlusion location, n (%)		
M1 middle cerebral artery	67 (57.3%)	87 (54.4%)
M2 middle cerebral artery	23 (19.7%)	37 (23.2%)
Internal carotid artery	24 (20.5%)	31 (19.4%)
Basilar	3 (2.6%)	5 (3.2%)
Occlusion side (left)	48 (41%)	80 (50%)
Last known well to arterial puncture, min; median (IQR)	172 (128.3–273)	179 (127.5–293)
Procedure aspects		
General anesthesia	49 (41.9%)	73 (45.6%)
BGC use only	25 (21.4%)	43 (26.9%)
BGC+intermediate catheter use	10 (29.9%)	56 (35%)
Intermediate catheter use only	21 (17.9%)	28 (17.5%)

ASPECTS indicates Alberta Stroke Program Early CT Score; BGC, balloon guide catheter; CAD, coronary artery disease; IQR, interquartile range; MI, myocardial infarction; mRS, modified Rankin Scale; TIGER, Treatment With Intent to Generate Endovascular Reperfusion; and tPA, tissue-type plasminogen activator.

performance goal and the 73.4% historical rate (noninferiority $P<0.0001$; superiority $P<0.01$; Table 2, Figure IIB in the [Data Supplement](#)).

The full distribution of reperfusion outcomes is shown in Figure 2A. The rate of excellent reperfusion (mTICI 2c-3) within 3 passes of the Tigertriever was 63.2% and 71.8% after all interventions. The frequency of first pass excellent reperfusion (mTICI 2c-3) was 41.4% and of first pass substantial reperfusion (mTICI 2b-3) 57.8%. The mean number of passes with the Tigertriever was 1.8 ± 0.9 . The rates of reperfusion achieved with each device size (Tigertriever and Tigertriever 17) and in each target artery are reported in Text I in the in the [Data Supplement](#).

Rescue therapy was employed in 33 of 117 (28.2%) main-study phase patients (frequency of other mechanical thrombectomy devices, angioplasty, intraarterial thrombolysis, and intracranial stenting shown in Table VI in the [Data Supplement](#)). In 5, at the interventionalist's discretion, a different device was used before reaching 3 passes with Tigertriever. Per protocol, these were counted as failures despite not having not made 3 attempts with Tigertriever. In 13, the additional devices were used at interventionalist discretion to further improve reperfusion despite mTICI $\geq 2b$ having been attained within 3 passes of the Tigertriever. Per protocol, these were counted as Tigertriever successes in achieving mTICI $\geq 2b$. Final reperfusion rates after all interventions were 95.7% mTICI 2b-3; 71.8% mTICI 2c-3; and 41.9% mTICI 3. Median time from arterial puncture to successful reperfusion (mTICI 2b-3) was 24 minutes (IQR, 16–38).

Primary efficacy end point and additional reperfusion results in lead-in patients did not statistically differ from main-study phase patient results (Table VII in the [Data Supplement](#)).

Primary Safety End Point

In the primary analysis main-study phase population, the rate for the primary safety composite end point of combined sICH and 90-day all-cause mortality was 18.1% (95% CI, 11.1%–25.1%, Table 3), compared with the PG of 30.4% and historical rate of 20.4% (noninferiority $P=0.004$; superiority $P=0.57$). The primary safety end point in lead-in patients did not statistically differ from main-study phase patient results (Table VIII in the [Data Supplement](#)).

Secondary End Points

Secondary clinical efficacy outcomes are shown in Table 2. In main-study phase patients, good clinical outcome (functional independence, mRS score 0–2) was attained by 58.0% (95% CI, 48.3%–67.3%), a rate superior to the 43.5% in the pooled comparator trials ($P=0.006$). The full distribution of disability levels at 90 days is shown in Figure 2B. Health-related quality of life on the EQ-5D at 90 days was median 80 (IQR, 70–90) and granular disability level on the Academic Medical

Table 2. Angiographic and Clinical Efficacy Outcomes in Main-Study and Combined Lead-In/Main-Study Patients

	TIGER	
	Main-study (n=117)	All (n=160)
Primary efficacy end point		
Successful reperfusion (mTICI 2b-3 within 3 TIGER passes without rescue), n (%)	99 (84.6%) (95% CI, 78.1%–91.2%)	130 (81.3%) (95% CI, 75.2%–87.3%)
Angiographic outcomes within 3 TIGER passes		
Excellent reperfusion (mTICI 2c-3 within 3 TIGER passes without rescue), n (%)	74 (63.2%)	94 (58.8%)
0	9 (7.7%)	17 (10.6%)
1	2 (1.7%)	2 (1.3%)
2a	7 (6%)	11 (6.9%)
2b	25 (21.4%)	36 (22.5%)
2c	28 (23.9%)	30 (18.8%)
3	46 (39.3%)	64 (40%)
Other angiographic and procedural outcomes		
Final successful reperfusion (mTICI 2b-3)	112 (95.7%)	150 (93.8%)
Final excellent reperfusion (mTICI 2-3c)	84 (71.8%)	106 (66.3%)
0	3 (2.6%)	4 (2.5%)
1	0	0
2a	2 (1.7%)	6 (3.8%)
2b	28 (23.9%)	44 (27.5%)
2c	35 (29.9%)	38 (23.8%)
3	49 (41.9%)	68 (42.5%)
First-pass successful reperfusion (mTICI 2b-3)*	67/116 (57.8%)	89/158 (56.3%)
First-pass excellent reperfusion (mTICI 2c-3)*	48/116 (41.4%)	66/158 (41.8%)
Use of rescue therapy†	33 (28.2%)	51 (31.9%)
Time from puncture to mTICI 2b-3, median (IQR)	24 (16–38)	25 (17–43)
	N=103	N=139
Time from puncture to closure, median (IQR)†	43 (28–72.5)	48 (30–84)
Clinical outcomes		
90-d good outcome (mRS score 0–2), n (%) (95% CI)	65/112 (58%)	86/154 (55.8%)
	(95% CI, 48.33%–67.29%)	(95% CI, 47.63%–63.83%)
EQ-5D at 90 d, median (IQR)	0.80 (0.70–0.90)	0.80 (0.68–0.90)
	(N=81)	(N=112)
ALDS at 90 d, median (IQR), mean (SD)	93.3% (57.8–100)	93.3% (60.4–100)
	77.3% (30.2%)	78.1% (30.0%)
	(N=83)	(N=114)

ALDS indicates Academic Medical Center Linear Disability Score; IQR, interquartile range; mRS, modified Rankin Scale; and mTICI, modified Thrombolysis in Cerebral Ischemia.

*Angio after first pass was not available from 1 lead-in and 1 main-study patient.

†Closure time was not available from 1 main-study patient.

‡Stentriever, aspiration retriever, IA tPA, angioplasty or stenting.

Center Linear Disability Score at 90 days was median 93.3 (IQR, 57.8–100).

Secondary safety outcomes are shown in Table 3. Neurological deterioration by 4 or more NIHSS points by 24 hours occurred in 7.7% (95% CI, 3.6%–14.2%). The rate of embolization to a previously uninvolved new territory was 2.6% (95% CI, 0.5%–7.4%), a rate superior to the 7.4% in the pooled comparator trials ($P=0.0403$). Asymptomatic intracranial hemorrhage within 24 hours occurred in 31.0% (95% CI, 22.7%–40.3%). (Additional data on radiological classification of hemorrhages are shown in Tables IX and X in the [Data Supplement](#).)

Secondary clinical efficacy and safety results in lead-in patients did not statistically differ from main-study phase patient results (Tables III and IV in the [Data Supplement](#)). All adverse events are shown in Table XI in the [Data Supplement](#).

Tigertriever Performance in Vessels of Different Size

Tigertriever performance in smaller (<2 mm) and larger (≥ 2 mm) vessels is shown in Table 4. Median vessel diameters were 1.6 mm (IQR, 1.5–1.8) for smaller vessels and 2.5 mm (IQR, 2.2–3.0) for larger vessels. Rates of both the primary efficacy end point (successful reperfusion) and the primary safety end point components (sICH at 24 hours and all-cause mortality at 90 days) were comparable in both vessel size groups.

DISCUSSION

In this prospective, multicenter trial, the Tigertriever yielded high rates of reperfusion among patients with acute ischemic stroke with large vessel occlusion, with substantial reperfusion achieved in >8 of every 10 patients with the Tigertriever alone and >9 of every 10 patients after additional interventions. The rate of successful reperfusion surpassed the predefined performance goal against the predicate Trevo and Solitaire devices derived from regulatory registration and additional pivotal trials. The high successful reperfusion rates with the Tigertriever device translated to favorable functional outcomes, with nearly 6 of every 10 patients achieving functional independence at day 90.

Patients enrolled in the TIGER study were similar to the comparator studies in several key features, including age, presenting deficit severity (NIHSS), extent of ischemic injury on initial imaging (Alberta Stroke Program Early CTs), vascular risk factors, rate of pretreatment with tPA, and frequency of target occlusion in internal carotid artery and M1 middle cerebral artery (Table IV in the [Data Supplement](#)). The frequency of M2 middle cerebral artery occlusions was mildly higher in the present study than several comparator studies but lower than in the most recent comparator investigation, ARISE II.²⁰ The greater proportion of M2

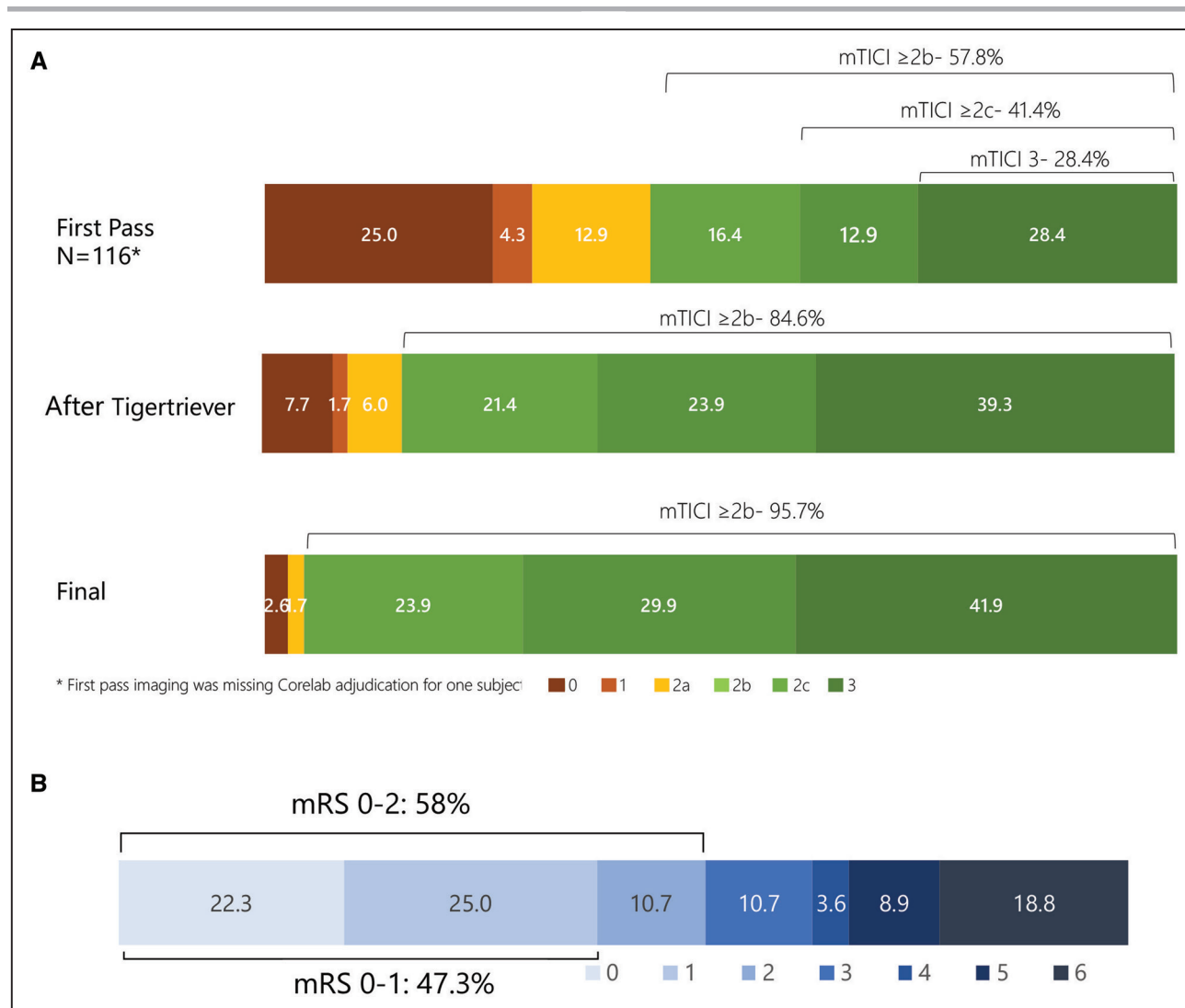


Figure 2. Reperfusion and functional outcomes in main-study phase of TIGER trial (Treatment With Intent to Generate Endovascular Reperfusion).

A, Reperfusion degree after first pass, after up to 3 Tigertriever passes, and after rescue therapy. **B**, Distribution of 90-d global disability outcomes on the modified Rankin Scale (mRS). mRS was available from 114 patients out of the 117 of the main-study group. mTICI indicates modified Thrombolysis in Cerebral Ischemia.

occlusions reflects the increasing treatment of more distal vessel occlusions in clinical practice.²¹ M2 occlusions have features that make them both more difficult targets for EVT, including distal location, smaller size, reduced accessibility, and less difficult, including smaller thrombus burden. As the frequency of reperfusion of the M2 was similar to other occlusion sites in the present study, the distribution of target occlusions did not affect comparison with prior studies. Patients in this study were treated sooner after onset (last known well to puncture 2 hour 59 minutes) than in prior studies, in accord with increasing emphasis in national guidelines and national practice on accelerating speed of endovascular thrombectomy.⁶⁻⁸

Reperfusion rates were better in the current study than in the comparator studies in a variety of parameters, including first pass reperfusion, mean number of passes, reperfusion after up to 3 passes before rescue

therapy, and final reperfusion after rescue therapy. First pass effect, defined as achieving mTICI 2c-3 reperfusion with one pass of the stent retriever, is likely to yield better patient outcomes and is a stringent metric to assess device performance.²² In the current study, a first pass effect was achieved in 41.4% of patients, which compares favorably to the rate of 25% in both the North American Solitaire Acute Stroke registry²³ and Trevo acute ischemic stroke registry.²⁴ In the ARISE II study, the first pass effect was seen in 40% of patients,²⁰ which is comparable to the present study. In the current study, successful reperfusion at procedure end, including rescue therapies, was achieved in nearly every patient (96%), a value that elevates cerebral reperfusion success rates to those of cardiac reperfusion procedures.^{25,26}

Faster time from puncture to achievement of reperfusion has been associated with improved clinical

Table 3. Safety End Points in Main-Study and Combined Lead-In/Main-Study Patients

	TIGER	
	Main-study (n=117)	All (n=160)
Primary safety composite end point, n (%)		
Symptomatic intracranial hemorrhage within 24 h and 90 d all-cause mortality*	21 (18.1%) (95% CI, 11.1%–25.1%)	31/159* (19.5%) (95% CI, 13.3%–25.7%)
Secondary safety end points, n (%)		
Symptomatic intracranial hemorrhage within 24 h†	2 (1.7%) (95% CI, 0.2%–6.03%)	3 (1.9%) (95% CI, 0.39%–5.38%)
90-d all-cause mortality*	21 (18.1%) (95% CI, 11.1%–25.1%)	31/159* (19.5%) (95% CI, 13.3%–25.7%)
Asymptomatic intracranial hemorrhage within 24 h	36 (31.0%) (95% CI, 22.7%–40.29%)	50‡ (31.4%) (95% CI, 24.32%–39.27%)
Neurological deterioration within 24 h	9 (7.7%) (95% CI, 3.61%–14.21%)	14 (8.8%) (95% CI, 4.9%–14.33%)
Embolization to new territory	3 (2.6%) (95% CI, 0.53%–7.37%)	4‡ (2.5%) (95% CI, 0.69%–6.31%)

CT indicates computed tomography; sICH, symptomatic intracranial hemorrhage; and TIGER, Treatment With Intent to Generate Endovascular Reperfusion.

*One patient withdrew consent at the day 4 visit so did not contribute full to 90-d mortality assessment.

†All sICH events occurred in cases treated with Tigertriever plus rescue therapy.

‡N=159. Reduced sample size due to missing 24 h CT or final angiogram.

outcomes.²⁷ The current study results demonstrate the median procedure time (from puncture to achievement of mTICI $\geq 2b$) was 25 (IQR, 17–43) minutes, which compares favorably to the ARISE II study, where the median time was 35 (IQR, 24–58) minutes.²⁰ This time was similar to the contact aspiration first pass thrombectomy time of median 25 minutes as reported in the COMPASS study (Aspiration Thrombectomy Versus Stent Retriever Thrombectomy as First-Line Approach for Large Vessel Occlusion) but shorter than the median 35 minutes reported for that trial's stent retriever group.²⁸ In experienced, high-volume centers performing over 48 thrombectomies annually, the median puncture-to-reperfusion time with the Trevo device was 67 (IQR, 42–105) minutes.²⁹ The reduced procedure times with the Tigertriever may reflect in part the technological advancement of the operator being able to adjust the radial force expansion to interact with the clot in a more productive manner.

The use of a BGC in the proximal vessel of the target occlusion when prior stent retrievers are deployed has been shown to improve the first pass reperfusion effect and clinical outcomes.³⁰ The proximal flow arrest may prevent distal embolization or non-target distal embolization. The unique design of the Tigertriever allows for more effective internalization of the thrombus into the interstices of the device. As the device is collapsed by the operator, it theoretically limits the chance of distal embolization as the clot is completely entwined in the retriever. The current study had a lower rate of BGC use at 35% in the intention to treat population, in comparison to the ARISE II trial in which this technique was used in 73.6% of patients.²⁰ Despite the lower use of BGC in this study, there were only 4 patients (2.5%) who had an embolization to a new territory compared with 6.6% in ARISE II and to 7.4% in the pooled data from of MR CLEAN,¹ ARISE II,²⁰ REVASCAT,³ and TREVO 2.¹³

Table 4. Tigertriever Performance in Vessels of Different Size Among Combined Lead-In/Main-Study Patients

	Vessel diameter <2 mm (N=27)	Vessel diameter ≥ 2 mm (N=133)	P Value
Vessel diameter, mean (SD), mm	1.58 (0.29)	2.71 (0.67)	...
Vessel diameter, median (IQR; range), mm	1.6 (IQR, 1.5–1.8) (Range, 0.7–1.9)	2.5 (IQR, 2.2–3) (Range, 2.0–6.0)	...
Successful reperfusion (mTICI 2b-3) within 3 Tigertriever passes without rescue, n (%)	23 (85.2%)	107 (80.4%)	0.79
90-d all-cause mortality, n (%)	3 (11.1%)	28 (21.0%)	0.29
sICH at 24 h, n (%)	0 (0%)	3 (2.2%)	1.00

IQR indicates interquartile range; mTICI, modified Thrombolysis in Cerebral Ischemia; and sICH, symptomatic intracranial hemorrhage.

Clinical outcomes in the present study were favorable. Functional independence (mRS score 0–2) at 90 days was achieved in 58% of patients, a rate higher than in the comparator studies. Similarly, health-related quality of life showed a utility value of 0.80, higher than the 0.57 in the MR CLEAN trial.³¹ These high rates of good long-term functional outcome likely reflect reduced total brain ischemia time because of high achieved reperfusion rates, early last known well-to-punctures times, and rapid puncture-to-reperfusion times. The faster last known well-to-puncture times versus historical comparator trials likely reflects interval evolution in endovascular workflow and systems of care, while the increased perfusion rates and rapid puncture-to-reperfusion times likely at least in part reflect intrinsic properties of the Tigertriever itself.

Limitations

This study has limitations. The study design was a single-arm trial against objective performance criteria derived from pooled prior studies of predicate devices, rather than a randomized trial with a contemporaneous control group. This approach constrains precision in delineating how well the Tigertriever compares with any particular comparator device. Reflecting the single-arm design, the central imaging and site clinical outcome evaluators and were not blinded to treatment assignment. However, for imaging the use of a central core laboratory was employed to mitigate any resulting bias. Similarly, for clinical outcomes, the primary efficacy outcome was performed by certified raters using a structured assessment system (the Rankin Focused Assessment) and the congruence of rater assessments of global disability and patient-reported health-related quality of life indicates accurate scoring.

Conclusions

In this prospective, multicenter study with a unique operator-controlled stent retriever, the Tigertriever achieved successful reperfusion within 3 passes in >8 of every 10 patients, demonstrating noninferiority over a performance goal and superiority over actual historical rates derived from trials and prospective studies of established devices. In addition, first pass successful reperfusion was achieved in nearly 6 of every 10 patients and final successful reperfusion in >9 of every 10 patients. Rates of embolization to a new territory and sICH were low and rates of good global disability and health-related quality of life outcomes high. Efficacy and safety outcomes were similar in vessels with diameters <2 and ≥2 mm.

ARTICLE INFORMATION

Received January 27, 2021; final revision received March 8, 2021; accepted March 11, 2021.

Presented in part at the International Stroke Conference, virtual, March 17–19, 2021.

Affiliations

Wellstar Medical Group, Department of Neurosurgery, Wellstar Health System Kennestone Hospital Marietta, GA (R.G., A.K.). Department of Neurology and Comprehensive Stroke Center, University of California Los Angeles (J.L.S., D.S.L.). Department of Neurosurgery, State University of New York at Buffalo (K.S., A.S., J.D.). Departments of Endovascular Neurosurgery and Stroke, St. Vincent Mercy Medical Center, Toledo, OH (E.L., O.O.Z.). Department of Neurology (D.Y.) and Department of Neurosurgery (R.S.), University of Miami School of Medicine, FL. Department of Neurosurgery, Stroke Institute, University of Pittsburgh Medical Center, PA (B.G., M.L.). Neurology, University of Pittsburgh Medical Center, PA (S.N.). Department of Neurosurgery, Cooper University Health Care, Camden, NJ (B.J.). Department of Neurology, Valley Baptist Medical Center, Harlingen, TX (A.H.). Stroke and Cerebrovascular Surgery, Lysterly Neurosurgery/Baptist Neurological Institute, Jacksonville, FL (R.H., A.A.). Advanced Neuroscience Network/Tenet South Florida, Delray Beach (R.K., A.M., N.M.-K.). Oregon Health and Science University, Portland (H.B., G.N., M.H., R.P., J.L.). Department of Radiology, Riverside Radiology and Interventional Associates, Columbus, OH (R.F.B., P.P., N.V.). Vascular Neurology of Southern California, Los Robles Hospital, Thousand Oaks (M.A.T.). Departments of Neurology, Neurosurgery and Radiology University of Iowa Hospitals and Clinics, Iowa City (E.S.). Departments of Neurology, Surgery/Neurosurgery, and Comprehensive Stroke Center, Maimonides Medical Center/SUNY Downstate Health Sciences University, Brooklyn, NY (Q.T.W.). Department of Neurosurgery, New York University Medical School (E.N.). Department of Neurointerventional Surgery, Baptist Cardiac and Vascular Institute, Miami, FL (G.D., I.L.). Department of Radiology, University of Massachusetts Medical School, Worcester (A.P.). Department of Neuroradiology, Rambam Health Care, Haifa, Israel (E.A.). Department of Emergency Medicine (S.S.) and Department of Radiology and Neurosurgery (S.T.), University of California Los Angeles. Department of Neurosurgery, Barrow Neurological Institute, Phoenix, AZ (A.P.J.).

Acknowledgments

The members of the Data Safety Monitoring Board were Vineeta Singh, MD (Chair), Colin Derdeyn, MD, M. Shazam Hussain, MD, Scott Hamilton, PhD (Biostatistician).

Sources of Funding

The study was funded by Rapid Medical. The study was designed with the input from an academic steering committee. There was an independent imaging core laboratory, independent statistician, and data safety monitoring board. The sponsor provided operational support with on-site and remote monitoring.

Disclosures

Dr Gupta serves as Principal Investigator (PI) for the TIGER (Treatment With Intent to Generate Endovascular Reperfusion) Study (Rapid Medical), PI for the ASSIST Registry (Stryker Neurovascular), PI for the RECCLAIM II Study (Zoll), CLEAR Study (Vesalio; The Vesalio NeVa Stent Retriever Study for Treatment of Large Vessel Occlusion Strokes), Clinical Events Committee (CEC) for the MIND Trial (Penumbra), consultant for Cerenovus. Dr Saver discloses consultant or advisory boards from Medtronic, Stryker, Cerenovus, and Rapid Medical; institutional conflict of interest from University of California. The University of California has patent rights in retrieval devices for stroke. Dr Levy is consultant or advisory board or Ownership interest for Claret Medical, GLG consulting, Guidepoint Global, Imperative Care, Medtronic, StimMed, Misionix, Mosaic, Clarion, Stryker, NeXtGen Biologics, MEDX, Cognition Medical, Endostream Medical, Rapid Medical, Rebound Therapeutics, Three Rivers Medical, and received Honorarium for training from Medtronic and Penumbra. Dr Zaidat is a consultant for Stryker Neurovascular, Cerenovus, Penumbra, Rapid Medical and Medtronic. He has also received research grants from Penumbra and Stryker Neurovascular. Dr Yavagal is a consultant for Medtronic, Cerenovus, Rapid Medical, Vascular Dynamics, Poseydon, Neurosave, Neural Analytics and Galaxy Therapeutics. Dr Liebeskind is a consultant as the imaging core lab for Cerenovus, Genentech, Medtronic, Stryker and Rapid Medical. Dr Gross is a consultant for Medtronic and Microvention. Dr Jankowitz is a consultant for Medtronic and Stryker. Dr Snyder has reported he is a consultant for Toshiba, Medtronic, EV3, Abbott Vascular, Micrus, Boston Scientific, Codman, Zimmer, Stryker, Vital, Cannon. Significant financial interest in Endo Tex, Micrus, BSC EPI, Access Closure, Inc, Cordis, Primus. He is a major stockholder in Boston Scientific, Access Closure, Inc, Niagara Gorge Medical. Dr Siddiqui is a consultant for Amnis Therapeutics, Apellis Pharmaceuticals, Boston Scientific, Canon Medical Systems, United States, Inc, Cerebrotech Medical Systems, Inc, Cerenovus, Corindus, Inc, Endostream Medical, Ltd, Imperative Care, Inc, Integra Lifesciences Corp, IRRAS, Medtronic, Microvention, Minnetronix Neuro, Inc, Northwest University-DSMB for HEAT trial, Penumbra, Perflow Medical, Ltd, Q'Apel Medical, Inc, Rapid Medical, Rebound Therapeutics Corp, Serenity Medical, Inc, Silk Road Medical, StimMed, Stryker, Three Rivers Medical, Inc, VasSol,

Vizai, Inc. WL. Gore Associates. National PI/Steering Committee for Cerenovous NAPA Trial and ARISE II Trial, Medtronic SWIFT PRIME (Solitaire With the Intention for Thrombectomy as Primary Endovascular Treatment) and SWIFT DIRECT Trials (Bridging Thrombolysis Versus Direct Mechanical Thrombectomy in Acute Ischemic Stroke), Microvention FRED Trial (The Flow Direction Endoluminal Device) and CONFIDENCE Study (Carotid Stent Trial to Evaluate the Safety and Efficacy of the Roadsaver Stent Used in Conjunction with the Nanoparasol Embolic Protection System for Patients at Increased Risk for Adverse Events From Carotid Endarterectomy), Medical University of South Carolina POSITIVE Trial (Perfusion Imaging Selection of Ischemic Stroke Patients for Endovascular Therapy), Penumbra 3D Separator Trial, COMPASS Trial (Aspiration Thrombectomy Versus Stent Retriever Thrombectomy as First-Line Approach for Large Vessel Occlusion), INVEST Trial. Financial interests in Adona Medical, Amnis Therapeutics, Bend IT Technologies, Ltd., BlinkTBI, Inc, Boston Scientific Corp (for purchase of Claret Medical), Buffalo Technology Partners, Inv, Cardinal Consultants, LLC, Cerebrotech Medical Systems, Inv. Cognition Medical, Endostream Medical, Ltd, Imperative Care, Inc, Instylla, Inc, International Medical Distribution Partners, IRRAS, LaunchNY Seed Fund Management, LLC, NeuroRadial Technologies, Inc, Neurovascular Diagnostics, Inc, Perflow Medical, Ltd, QApel, Inc, Radical Catheter Technologies, Inc, Rebound Therapeutics Corp (Purchased 2019 by Integra Lifesciences Corp), Rist Neurovascular, Inc, Sense Diagnostics, Inc, Serenity Medical, Inc, Silk Road Medical, Spinnaker Medical, Inc, StimMed, Synchron, Three Rivers Medical, Inc, Truvic Medical, Inc, Vastrax, LLC, VICIS, Inc, Viseon, Inc, Vizai, Inc. Research grants as co-investigator NIH/NINDS 1R01NS091075 Virtual intervention of aneurysms and Co-Principal Investigator NIH-NINDS R21 NS109575-01 Optimizing approaches to endovascular therapy of acute ischemic stroke. Dr Davies is a consultant for Medtronic, Microvention. Research support NIH RO1. Shares and ownership of Cerebrotech and Rist Neurovascular. Dr Hassan is a consultant for Medtronic, Stryker, Microvention, Penumbra, Cerenovous, Genentech, GE Healthcare, Scientia, Balt, Vizai, Insera therapeutics, Proximie, NovaSignal, Vesalio. PI COMPLETE Study (Penumbra; International Acute Ischemic Stroke Registry With the Penumbra System Aspiration Including the 3D Revascularization Device) and LVO SYNCHRONISE (Vizai; Observational Study of Automated Detection for Identification, Triage, and Timely Intervention in Large Vessel Occlusions). Steering Committee for SELECT (Optimizing Patient Selection for Endovascular Treatment in Acute Ischemic Stroke), DAWN, SELECT 2 (Optimizing Patient Selection for Endovascular Treatment in Acute Ischemic Stroke 2), EXPEDITE II (Evaluate a Next Generation Portable Diagnostic Platform for Determination and Immediate Triage of Emergency Large Vessel Stroke 2), EMBOLISE (Embolization of the Middle Meningeal Artery With ONYX Liquid Embolic System for Subacute and Chronic Subdural Hematoma), CLEAR. Grant support from GE Healthcare and Valley Baptist. Dr Hanel is a consultant for Rapid Medical, Medtronic, Stryker, Cerenovous, Balt, Phenox, Elum, MIVI, ThrombX, Endostream, RIST, REIST, Serenity, BendIT. Dr Malek is a consultant for Stryker, RAPID, InNeuroCo. Dr Mueller-Kronast reports no financial disclosures. Dr Starke has research support from NREG, Joe Niekro Foundation, Brain Aneurysm Foundation, Been Foundation and NIHNS R01 NS111119-01A1 and UL 1TR002736 and KL 2TR002737 through the Miami Clinical and Translational Science Institute, National Center for Advancing Translational Sciences and the National Institute on Minority Health and Health Disparities. Unrestricted research grants from Medtronic. Consultant for Penumbra, Abbott, Medtronic, InNeuroCo and Cerenovous. Dr Bozorgchami is a consultant for Cerenovous. Dr Nesbit reports funding from Oregon Health Sciences. Dr Horikawa is a consultant for Terumo, Inc. Dr Priest is a consultant for Medtronic, Cerenovous and Stryker. Dr Vora is a consultant for Medtronic and Microvention. Dr Taqi is a consultant for Rapid Medical, Stryker and Medtronic. Dr Samaniego is a consultant for Rapid Medical, Medtronic and Microvention. Dr Nossek is a consultant for Rapid Medical. Dr Dabus is a consultant for Microvention, Penumbra, Medtronic, Cerenovous, InNeuroCo. Dr Linfante is a consultant for Medtronic, Stryker, Cerenovous. He is a major shareholder for Three Rivers, Prolong Pharmaceuticals, Prometheus, InNeuroCo. Dr Puri is a consultant for Stryker, Cerenovous, Medtronic, Microvention, QApel, Merit Medical, Arsenal Medical. He has received grant support from SBIR, NIH. Speaker Bureau for Merit Medical, Cerenovous, QApel. Stock options from InNeuroCo, Galaxy Therapeutics, NTI, Agile Medical and Perfuzo. Dr Starkman is a consultant for Medtronic, Microvention. He reports stock ownership with Rist Neurovascular and Cerebrotech. He has research support from the NIH. Dr Tateshima is a consultant for Medtronic, Stryker, Cerenovous, Balt, Phenox, Spartan Micro, and Irvine Neurovascular. He has research support from Biomedical Solutions, Inc. He reports educational support from Microvention. The other authors report no conflicts.

REFERENCES

- Berkhemer OA, Fransen PS, Beumer D, van den Berg LA, Lingsma HF, Yoo AJ, Schonewille WJ, Vos JA, Nederkoorn PJ, Wermer MJ, et al; MR CLEAN Investigators. A randomized trial of intraarterial treatment for acute ischemic stroke. *N Engl J Med*. 2015;372:11–20. doi: 10.1056/NEJMoa1411587
- Saver JL, Goyal M, Bonafe A, Diener HC, Levy EI, Pereira VM, Albers GW, Cognard C, Cohen DJ, Hacke W, et al; SWIFT PRIME Investigators. Stent-retriever thrombectomy after intravenous t-PA vs. t-PA alone in stroke. *N Engl J Med*. 2015;372:2285–2295. doi: 10.1056/NEJMoa1415061
- Jovin TG, Chamorro A, Cobo E, de Miquel MA, Molina CA, Rovira A, San Roman L, Serena J, Abilleira S, Ribó M, et al; REVASCAT Trial Investigators. Thrombectomy within 8 hours after symptom onset in ischemic stroke. *N Engl J Med*. 2015;372:2296–2306. doi: 10.1056/NEJMoa1503780
- Goyal M, Demchuk AM, Menon BK, Eesa M, Rempel JL, Thornton J, Roy D, Jovin TG, Willinsky RA, Sapkota BL, et al; ESCAPE Trial Investigators. Randomized assessment of rapid endovascular treatment of ischemic stroke. *N Engl J Med*. 2015;372:1019–1030. doi: 10.1056/NEJMoa1414905
- Campbell BC, Mitchell PJ, Kleinig TJ, Dewey HM, Churilov L, Yassi N, Yan B, Dowling RJ, Parsons MW, Oxley TJ, et al; EXTEND-IA Investigators. Endovascular therapy for ischemic stroke with perfusion-imaging selection. *N Engl J Med*. 2015;372:1009–1018. doi: 10.1056/NEJMoa1414792
- Powers WJ, Rabinstein AA, Ackerson T, Adeoye OM, Bambakidis NC, Becker K, Biller J, Brown M, Demaerschalk BM, Hoh B, et al. Guidelines for the early management of patients with acute ischemic stroke: 2019 update to the 2018 guidelines for the early management of acute ischemic stroke: a guideline for healthcare professionals from the American Heart Association/American Stroke Association. *Stroke*. 2019;50:e344–e418. doi: 10.1161/STR.0000000000000211
- Turc G, Bhogal P, Fischer U, Khatri P, Lobotesis K, Mazighi M, Schellinger PD, Toni D, de Vries J, White P, et al. European Stroke Organisation (ESO)- European Society for Minimally Invasive Neurological Therapy (ESMINT) guidelines on mechanical thrombectomy in acute ischemic stroke. *J Neurointerv Surg*. 2019;11:535–538. doi: 10.1136/neurintsurg-2018-014568
- Mokin M, Ansari SA, McTaggart RA, Bulsara KR, Goyal M, Chen M, Fraser JF; Society of NeuroInterventional Surgery. Indications for thrombectomy in acute ischemic stroke from emergent large vessel occlusion (ELVO): report of the SNIS Standards and Guidelines Committee. *J Neurointerv Surg*. 2019;11:215–220. doi: 10.1136/neurintsurg-2018-014640
- Goyal M, Menon BK, van Zwam WH, Dippel DW, Mitchell PJ, Demchuk AM, Dávalos A, Majore CB, van der Lugt A, de Miquel MA, et al; HERMES collaborators. Endovascular thrombectomy after large-vessel ischaemic stroke: a meta-analysis of individual patient data from five randomised trials. *Lancet*. 2016;387:1723–1731. doi: 10.1016/S0140-6736(16)00163-X
- Kara B, Selcuk HH, Erbahceci Salik A, Zalov H, Yildiz O, Gul G, Balkan B. Single-center experience with the Tigertriever device for the recanalization of large vessel occlusions in acute ischemic stroke. *J Neurointerv Surg*. 2019;11:455–459. doi: 10.1136/neurintsurg-2018-014196
- Gruber P, Diepers M, von Hesslering A, Weber J, Kahles T, Anon J, Berberat J, Nedeltchev K, Liebeskind DS, Remonda L. Mechanical thrombectomy using the new Tigertriever in acute ischemic stroke patients - A Swiss prospective multicenter study. *Interv Neuroradiol*. 2020;26:598–601. doi: 10.1177/1591019920946499
- Will L, Maus V, Maurer C, Weber A, Weber W, Fischer S. Mechanical thrombectomy in acute ischemic stroke using a maximally expandable stent retriever (Tigertriever): preliminary single center experience [published online June 11, 2020]. *Clin Neuroradiol*. doi: 10.1007/s00062-020-00919-w
- Nogueira RG, Lutsep HL, Gupta R, Jovin TG, Albers GW, Walker GA, Liebeskind DS, Smith WS; TREVO 2 Trialists. Trevo versus Merci retrievers for thrombectomy revascularisation of large vessel occlusions in acute ischaemic stroke (TREVO 2): a randomised trial. *Lancet*. 2012;380:1231–1240. doi: 10.1016/S0140-6736(12)61299-9
- Saver JL, Jahan R, Levy EI, Jovin TG, Baxter B, Nogueira RG, Clark W, Budzik R, Zaidat OO; SWIFT Trialists. Solitaire flow restoration device versus the Merci Retriever in patients with acute ischaemic stroke (SWIFT): a randomised, parallel-group, non-inferiority trial. *Lancet*. 2012;380:1241–1249. doi: 10.1016/S0140-6736(12)61384-1
- Goyal M, Fargen KM, Turk AS, Mocco J, Liebeskind DS, Frei D, Demchuk AM. 2C or not 2C: defining an improved revascularization grading scale and the need for standardization of angiography outcomes in stroke trials. *J Neurointerv Surg*. 2014;6:83–86. doi: 10.1136/neurintsurg-2013-010665
- von Kummer R, Broderick JP, Campbell BC, Demchuk A, Goyal M, Hill MD, Treumiet KM, Majoie CB, Marquering HA, Mazya MV, et al. The Heidelberg bleeding classification: classification of bleeding events after ischemic stroke and reperfusion therapy. *Stroke*. 2015;46:2981–2986. doi: 10.1161/STROKEAHA.115.010049
- Feng YS, Kohlmann T, Janssen MF, Buchholz I. Psychometric properties of the EQ-5D-5L: a systematic review of the literature. *Qual Life Res*. 2021;30:647–673. doi: 10.1007/s11136-020-02688-y

18. Holman R, Weisscher N, Glas CA, Dijkgraaf MG, Vermeulen M, de Haan RJ, Lindeboom R. The Academic Medical Center Linear Disability Score (ALDS) item bank: item response theory analysis in a mixed patient population. *Health Qual Life Outcomes*. 2005;3:83. doi: 10.1186/1477-7525-3-83
19. Weisscher N, Vermeulen M, Roos YB, de Haan RJ. What should be defined as good outcome in stroke trials; a modified Rankin score of 0-1 or 0-2? *J Neurol*. 2008;255:867–874. doi: 10.1007/s00415-008-0796-8
20. Zaidat OO, Bozorgchami H, Ribó M, Saver JL, Mattle HP, Chapot R, Narata AP, Francois O, Jadhav AP, Grossberg JA, et al. Primary Results of the Multicenter ARISE II Study (Analysis of Revascularization in Ischemic Stroke With EmboTrap). *Stroke*. 2018;49:1107–1115. doi: 10.1161/STROKEAHA.117.020125
21. Nogueira RG, Mohammaden MH, Haussen DC, Budzik RF, Gupta R, Krajina A, English JD, Malek AR, Sarraj A, Narata AP, et al. Endovascular therapy in the distal neurovascular territory: results of a large prospective registry [published online December 15, 2020]. *J Neurointerv Surg*. doi: 10.1136/neurintsurg-2020-016851
22. Zaidat OO, Castonguay AC, Linfante I, Gupta R, Martin CO, Holloway WE, Mueller-Kronast N, English JD, Dabus G, Malisch TW, et al. First pass effect: a new measure for stroke thrombectomy devices. *Stroke*. 2018;49:660–666. doi: 10.1161/STROKEAHA.117.020315
23. Zaidat OO, Castonguay AC, Gupta R, Sun CJ, Martin C, Holloway WE, Mueller-Kronast N, English JD, Linfante I, Dabus G, et al. North American Solitaire Stent Retriever Acute Stroke registry: post-marketing revascularization and clinical outcome results. *J Neurointerv Surg*. 2018;10(suppl 1):i45–i49. doi: 10.1136/neurintsurg-2013-010895.rep
24. Zaidat OO, Castonguay AC, Nogueira RG, Haussen DC, English JD, Satti SR, Chen J, Farid H, Borders C, Veznedaroglu E, et al. TREVO stent-retriever mechanical thrombectomy for acute ischemic stroke secondary to large vessel occlusion registry. *J Neurointerv Surg*. 2018;10:516–524. doi: 10.1136/neurintsurg-2017-013328
25. Patel RD, Saver JL. Evolution of reperfusion therapies for acute brain and acute myocardial ischemia: a systematic, comparative analysis. *Stroke*. 2013;44:94–98. doi: 10.1161/STROKEAHA.112.666925
26. Sun CH, Bhatt DL, Nogueira RG, Gupta R. Endovascular therapy for stroke: getting to the "heart" of the matter. *Circulation*. 2014;129:1152–1160. doi: 10.1161/CIRCULATIONAHA.113.003703
27. Jahan R, Saver JL, Schwamm LH, Fonarow GC, Liang L, Matsouaka RA, Xian Y, Holmes DN, Peterson ED, Yavagal D, et al. Association between time to treatment with endovascular reperfusion therapy and outcomes in patients with acute ischemic stroke treated in clinical practice. *JAMA*. 2019;322:252–263. doi: 10.1001/jama.2019.8286
28. Turk AS 3rd, Siddiqui A, Fifi JT, De Leacy RA, Fiorella DJ, Gu E, Levy EI, Snyder KV, Hanel RA, Aghaebrahim A, et al. Aspiration thrombectomy versus stent retriever thrombectomy as first-line approach for large vessel occlusion (COMPASS): a multicentre, randomised, open label, blinded outcome, non-inferiority trial. *Lancet*. 2019;393:998–1008. doi: 10.1016/S0140-6736(19)30297-1
29. Nogueira RG, Haussen DC, Castonguay A, Rebello LC, Abraham M, Puri A, Alsheklee A, Majhoo A, Farid H, Finch I, et al. Site Experience and outcomes in the trevo acute ischemic stroke (TRACK) Multicenter Registry. *Stroke*. 2019;50:2455–2460. doi: 10.1161/STROKEAHA.118.024639
30. Zaidat OO, Mueller-Kronast NH, Hassan AE, Haussen DC, Jadhav AP, Froehler MT, Jahan R, Ali Aziz-Sultan M, Klucznik RP, Saver JL, et al; STRATIS Investigators. Impact of balloon guide catheter use on clinical and angiographic outcomes in the STRATIS Stroke Thrombectomy Registry. *Stroke*. 2019;50:697–704. doi: 10.1161/STROKEAHA.118.021126
31. Schreuders J, van den Berg LA, Franssen PS, Berkhemer OA, Beumer D, Lingsma HF, van Oostenbrugge RJ, van Zwam WH, Majoie CB, van der Lugt A, et al; MR CLEAN investigators. Quality of life after intra-arterial treatment for acute ischemic stroke in the MR CLEAN trial-Update. *Int J Stroke*. 2017;12:708–712. doi: 10.1177/1747493017706244