



Article

# Development of a High-Throughput Calcium Mobilization Assay for CCR6 Receptor Coupled to Hydrolase Activity Readout

Sara Gómez-Melero 1,\*, Fé Isabel García-Maceira 1, Tania García-Maceira 1, Verónica Luna-Guerrero 1, Gracia Montero-Peñalvo 1, Javier Caballero-Villarraso 2,3, Isaac Túnez 2,3 and Elier Paz-Rojas 4

- Canvax Biotech, Parque Científico y Tecnológico Rabanales 21, c/Astrónoma Cecilia Payne s/n, Edificio Canvax, 14014 Córdoba, Spain; fi.garciamaceira@gmail.com (F.I.G.-M.); tagama0802@gmail.com (T.G.-M.); veronicalunaguerrero@gmail.com (V.L.-G.); gracia7432@gmail.com (G.M.-P.)
- <sup>2</sup> Department of Biochemistry and Molecular Biology, Faculty of Medicine and Nursing, University of Córdoba, Avda. Menéndez Pidal s/n, 14004 Córdoba, Spain; bc2cavij@uco.es (J.C.-V.); itunez@uco.es (I.T.)
- Maimonides Biomedical Research Institute of Cordoba, Avda. Menéndez Pidal s/n, 14004 Córdoba, Spain
- <sup>4</sup> Multiplex Biopharma, Parque Científico Rabanales 21, Calle Astronoma Cecilia Payne, 14014 Córdoba, Spain; pazrojaselier@gmail.com
- \* Correspondence: sara.gomez.melero@gmail.com

**Abstract:** CCR6 is a chemokine receptor highly implicated in inflammatory diseases and could be a potential therapeutic target; however, no therapeutic agents targeting CCR6 have progressed into clinical evaluation. Development of a high-throughput screening assay for CCR6 should facilitate the identification of novel compounds against CCR6. To develop a cell-based assay, RBL-2H3 cells were transfected with plasmids encoding β-hexosaminidase and CCR6. Intracellular calcium mobilization of transfected cells was measured with a fluorescent substrate using the activity of released hexosaminidase as readout of the assay. This stable, transfected cell showed a specific signal to the background ratio of 19.1 with low variability of the signal along the time. The assay was validated and optimized for high-throughput screening. The cell-based calcium mobilization assay responded to the specific CCR6 ligand, CCL20, in a dose-dependent manner with an EC50 value of 10.72 nM. Furthermore, the assay was deemed robust and reproducible with a Z′ factor of 0.63 and a signal window of 7.75. We have established a cell-based high-throughput calcium mobilization assay for CCR6 receptor. This assay monitors calcium mobilization, due to CCR6h activation by CCL20, using hexosaminidase activity as readout. This assay was proved to be robust, easy to automate and could be used as method for screening of CCR6 modulators.

**Keywords:** high-throughput screening; HTS; CCR6; hexosaminidase; calcium mobilization; cell-sensor; GPCR

# Citation: Gómez-Melero, S.; García-Maceira, F.I.; García-Maceira, T.; Luna-Guerrero, V.; Montero-Peñalvo, G.; Caballero-Villarraso, J.; Túnez, I.; Paz-Rojas, E. Development of a High-throughput Calcium Mobilization Assay for CCR6 Receptor Coupled to Hydrolase Activity Readout. *Biomedicines* 2022, 10, 422. https://doi.org/ 10.3390/biomedicines10020422

Academic Editors: Yoshitaka Hosokawa

Received: 19 December 2021 Accepted: 7 February 2022 Published: 10 February 2022

**Publisher's Note:** MDPI stays neutral with regard to jurisdictional claims in published maps and institutional affiliations.



Copyright: © 2022 by the authors. Licensee MDPI, Basel, Switzerland. This article is an open access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (http://creativecommons.org/licenses/by/4.0/).

#### 1. Introduction

G-protein-coupled receptors (GPCRs) represent 4% of the human genome and are one of the most attractive therapeutic targets due to their involvement in many important physiological functions [1]. Despite being one of the largest classes of proteins, the majority of GPCRs remains undragged, with only approximately 100 GPCRs targeted with small molecules or peptides [2]. GPCR receptors do not have activity per se, but they are linked to GTPases that participate in signal transduction. Such GTPases are collectively known as G proteins and each G protein comprises three subunits,  $\alpha$ ,  $\beta$ , and  $\gamma$ . There are four families of G proteins that are classified according to their  $\alpha$  subunit:  $G_s$ ,  $G_i$ ,  $G_q$  and  $G_{12/13}$  where  $G_s$  and  $G_i$  participate in cAMP-dependent pathways,  $G_q$  activates phospholipase C-dependent pathways, and  $G_{12/13}$  activates a small family of

Biomedicines **2022**, 10, 422 2 of 12

GTPases [3–5]. In addition, the  $G_q$  family is further divided into four  $G_q$  proteins due to the use of four different  $\alpha$  polypeptides,  $\alpha q$ ,  $\alpha 11$ ,  $\alpha 14$ , and  $\alpha 15/16$  [6]. The  $G_{\alpha 15/16}$  polypeptide has a restricted tissue expression and is predominantly expressed in myeloid and B-cell lineages while  $G_{\alpha q}$  and  $G_{\alpha 11}$  are ubiquitously expressed [6].

There have been extensive efforts to develop new, sensitive and simple-to-use assays to measure second messengers or molecules coupled to the signal transduction as readout for GPCR activation [7]. Based on this signal transduction cascade several assay techniques such as radioligand binding, reporter gene, cAMP detection and calcium mobilization are commonly used [8]. Calcium mobilization is the result of activation of GPCRs through specific subunit  $G_{\alpha q}$  activating classical phospholipase C pathways resulting in the release of calcium from the endoplasmic reticulum into the cytoplasm [7]. In resting cells, the cytosolic calcium concentration is lower (~100–200 nM) than that in the extracellular environment (~2 mM), but when the cells are excited by the activation of GPCRs the concentration of intracellular calcium can increase in seconds to 100  $\mu$ M [9]. Because this change in the concentration of intracellular  $Ca^{2+}$  play a crucial role in GPCR-induced signaling, the develop of measurement assays based on GPCR calcium mobilization has attracted more attention in GPCR-targeted drug discovery and has widely used in the investigation of many GPCRs [10].

The human CC chemokine receptor 6, hCCR6, is a class A of GPCR expressed in a particular diverse range of leukocytes including T cells (specifically Th17 cells and Treg cells), B cell, neutrophils and subsets of dendritic cells [11]. With particular relevance in rheumatoid arthritis, psoriasis, multiple sclerosis, inflammatory bowel disease and cancer; inhibition of the hCCR6 signaling might prove to be a useful therapeutic strategy [11,12]. Some studies indicate that upon CCR6 binding to its ligand, CCL20 is activated a signaling pathway through calcium mobilization [13,14]. Therefore, the development of novel functional assays based on calcium mobilization is of interest for CCR6h-targeted drug discovery [10].

Extensive research has been conducted regarding the molecular determinants of coupling of each GPCR to a particular family of G proteins and today it is known that GPCR and G protein coupling selectivity involves the interaction of GPCR with both the extreme carboxyl terminus and linker I of Gα protein [15] and that development of chimeras of  $G_{\alpha q}$  comprising the C-terminal five amino of  $G_{\alpha i}$  or  $G_{\alpha s}$  and a highly conserved glycine within Linker I is enough for the coupling of non-Gaq coupled GPCR to signaling through phospholipase C-dependent pathways [16]. In addition, it was previously known that  $G_{\alpha 15/16}$  is a promiscuous G protein, that is, it can couple to a wide variety of GPCR to signal through phospholipase C-dependent pathways [17] and this is the molecular basis of the measurement of calcium release upon activation of many nonnaturally  $G_{\alpha q}$  coupled-GPCRs assays that are used in HTS today; for example, the ChemiScreen™ family of HTS products RBL-2H3 is a rat basophilic leukemia cell line that expresses the  $G_{\alpha 15/16}$  of  $G_{\alpha}$  subunit. This promiscuous  $G_{\alpha 15/16}$  protein shifts the receptors coupling to the calcium mobilization pathway [18]. Mast cells release their secretory granules after aggregation of their high affinity receptors, Fc<sub>\varepsilon</sub> receptor I (Fc<sub>\varepsilon</sub>RI), by the antigen-immunoglobulin E (IgE) complexes [13,18]. Aggregation of FcERI by antigens initiates a signaling cascade that induces the elevation of intracellular Ca2+ causing degranulation [18]. In this way, the FcεRI stimulation induces the release of hydrolases, such as  $\beta$ -hexosaminidase (HexB), from intracellular granules [13].

We hypothesized that overexpression of a GPCR receptor on the surface of RBL-HexB cells generates a cell-based sensor useful to determine calcium mobilization through GPCR ligand binding using hydrolase activity as readout of the assay. This assay can be used to measure the effect of compounds on calcium mobilization signaling of CCR6h receptor. When the cell-based sensor is incubated with the specific CCR6h ligand, CCL20, the amount of intracellular calcium increases, and the reporter polypeptide is released. The enzymatic activity of such released non-protease hydrolase reporter polypeptide is detected with the specific substrate. Herein we report the development of a precise, simple and cost-effective assay to measure calcium mobilization for the CCR6h receptor.

Biomedicines **2022**, 10, 422 3 of 12

#### 2. Materials and Methods

#### 2.1. Plasmids

Coding sequence of human  $\beta$ -hexosaminidase (HexB, Genbank: BC017378.2) was synthesized by Integrated DNA Technologies (IDT), flanked with the recognition sites of 5′XhoI and 3′NotI and cloned into the pRV retroviral vector by the restriction with the indicated enzymes. A hygromycin resistance cassette was included in the vector backbone for stable cells selection. The vector also included an IRES-NGFR cassette, cloned downstream of human HexB and, thus, under the control of the same promoter, for selection by flow cytometry of cells expressing HexB stably.

Coding sequence of human CCR6 (Gb: NM\_004367) was obtained by synthesis (IDT), flanked with the recognition sites of 5'XhoI and 3'NotI. It was cloned in frame with the signal peptide IgK, a short sequence (VGS) to improve GPCR expression [19] and a tag sequence (c-myc) into the pcDNA3.1 vector by the restriction with the indicated enzymes. pcDNA3.1 vector included a neomycin resistance cassette for stable cells selection.

The integrity of the constructs was confirmed by sequence analysis (Stabvida, Caparica, Portugal).

#### 2.2. Cell Culture

All cell lines were purchased from the American Type Culture Collection (ATCC, Manassas, VA, USA). RBL-2H3 (rat basophilic leukemia, CRL-2256) cells was cultured in Eagle's minimum essential medium (EMEM)/RPMI-1640 medium. MYC 1-9E10.2 (9E10, CRL-1729) and 200-6-G3-4 (20.4, HB-8737) hybridoma cells were maintained in RPMI-1640 medium. All media (Gibco) were supplemented with 10% fetal bovine serum (FBS, Gibco, New York, NY, USA), 2 mM L-glutamine (Lonza, Basel, Switzerland), 100 U/mL penicillin (Gibco), and 100 µg/mL streptomycin (Gibco). Cells were maintained at 37 °C in a humidified atmosphere of 5% CO<sub>2</sub>.

# 2.3. Antibody Purification

Mouse anti-c-myc and mouse anti-NGFR were obtained from culture supernatant of MYC 1-9E10.2 (9E10, CRL-1729) and 200-6-G3-4 (20.4, HB-8737) hybridoma cells, respectively. Immunoglobulins were affinity-purified from the culture supernatant using Protein G Sepharose columns (HiTrap Protein G HP, GE Healthcare, Madrid, Spain) according to the manufacturer's protocol in AKTA Prime Plus (GE Healthcare). The purity of purified mAbs was checked by SDS-PAGE on 10% acrylamide gel and stained with Coomassie brilliant blue. The concentration of purified immunoglobulins was determined by UV absorbance at 280 nm.

#### 2.4. Flow Cytometry

For positive transfected cells selection, flow cytometric analysis was performed. Mouse anti-c-myc (9E10, in house) or mouse anti-NGFR (HB-8737, in house) antibodies were incubated with RBL-2H3 transfected cells or mock transfectants for 20 min at 4 °C. After washing with PBS, cells were stained with FITC-conjugated goat anti-mouse IgG (1:60, Sigma, Burlington, MA, USA) for 10 min at 4 °C and were analyzed by flow cytometry using a FACSCalibur (Becton Dickinson, Madrid, Spain) cytometer. A total of 10,000 events were acquired and data were analyzed with BD CellQuest Pro software (BD CellQuest Pro, BD Biosciences, Madrid, Spain).

# 2.5. Creation of Stable RBL-2H3 Clone with β-hexosaminidase Activity and CCR6 Expression

RBL-2H3 was transfected, with 2  $\mu g$  of the expression vector pRV-PCMV-HexB-IRES-NGFR-HYG obtained, using Canfast (Canvax Biotech, Córdoba, Spain) according to the manufacturer's instructions. Following transfection, stably transfectants were selected with 2 mg/mL hygromycin B (Duchefa, Torrente, Spain). To confirm the plasmid transfection, flow cytometry against NGFR with mouse anti-NGFR (HB8737, in house)

Biomedicines **2022**, 10, 422 4 of 12

antibody was used. Positive cells were cloned by limiting dilution in 96-well plates and wells with growing colonies were analyzed for NGFR expression by flow cytometry again.

Clones with positive expression for NGFR were expanded to six well plates for determination of HexB activity. The cells were plated in 384-well plates (Corning, 384-well flat-bottom black polystyrene microplates, Corning, New York, USA) in 30 µL per well at 5000 cells per well and cultured for 48 h at 37 °C, 5% CO2. Then, cell media was aspirated and replaced with a mix containing 1 µg/mL mouse anti-TNP antibody (clone IgELb4, Creative Biolabs, Shirley, NY, USA), 1 µg/mL of TNP-BSA (Biosearch Technologies, Hoddesdond, UK) and 1 mM of 4-Methylumbelliferyl N-acetyl-β-D-glucosaminide (4MU-NGlc) substrate (Glycosynth) in BSS (25 mM HEPES/NaOH, 1.2 mM KH2PO4, 65 mM NaCl, 5.65 mM KCl, 0.6 mM MgCl2, 1.8 mM CaCl2, 5.6 mM glucose and 0.1% BSA). Cells were incubated at 37 °C protected from light. The fluorescence was measured at 360 nm excitation and 470 nm emission wavelength in a FLUOStar OPTIMA plate reader (BMG Labtech, Ortenberg, Germany) to measure exocytosis of HexB. A read was taken at 0, 15, 30, 45 and 60 min of incubation. To measure background release, 2 mM of 4MU-NGlc was used without anti-TNP antibody or TNP-BSA. Ionomycin (Sigma, Burlington, MA, USA) at 10 μM was used as a positive control. Wells without cells were used as blank and their fluorescence was subtracted of specific and background release to calculate the specific signal-to-background (S/B) ratio. The stability of HexB expression in 1B7-RBL-2H3-HexB cells' activity in the time was measured for two months.

The selected clone, 1B7-RBL-2H3-HexB, with the best hexosaminidase activity was transfected with 2  $\mu g$  of the expression vector pcDNA3.1-PEF1 $\alpha$ -PSIg $\kappa$ -CMyc-VGS-CCR6-Neo using Canfast (Canvax Biotech, Córdoba, Sapin) according to the manufacturer's instructions. Following transfection, stably transfectants were selected with 0.5 mg/mL G418 (Phytotechnology, Lenexa, KS, USA). Flow cytometry was used to measure surface expression of CCR6 using an anti c-myc antibody (9E10 clone, mouse antibody) followed by anti-mouse-FITC (Sigma).

# 2.6. CCR6h Calcium Mobilization Assay

Mobilization of Ca<sup>2+</sup> by CCL20 ligand was determined using 1B7-RBL-HexB-hCCR6 cells. The cells were plated in 384-well plates (Corning, 384-well flat-bottom black polystyrene microplates) in 30  $\mu$ L per well at 5000 cells per well and cultured for 48 h at 37 °C, 5% CO<sub>2</sub>. Then, cell media was aspirated and replaced with serial dilutions of CCL20 (Peprotech, London, UK) and 1 mM 4-Mthylumbelliferyl N-acetyl- $\beta$ -D-glucosaminide substrate (Sigma) per well in BSS protected from light. Cells were incubated 60 min at 37 °C, 5% CO<sub>2</sub>. The fluorescence was measured (excitation: 360 nm, emission: 470 nm) using a FLUOStar OPTIMA plate reader (BMG Labtech). To generate dose-response curves of CCL20 and calculate EC<sub>50</sub> and Hill slope (H) value non-linear regression analyses were performed.

EC50 and H were used to calculate EC80 using the following equation:

$$EC_{F} = \left(\frac{F}{100 - F}\right)^{\frac{1}{H}} * EC_{50}$$
 (1)

where F is the fraction of maximal response, in this case 80.

# 2.7. Optimization and Performance of the HTS Assay

In order to optimize the incubation time, the assay was run as described above changing the incubation time with the ligand. Furthermore, the method of dispensing cells manually or using a Research Pro 8 Channel Electronic Multichannel Pipette (Eppendorf, Madrid, Spain) was compared. Concentration response curves of CCL20 were run to compare EC50 value and signal-to-background (S/B) ratio.

The assay with optimized conditions was performed to determine the Z' factor and the signal window. The Z' factor [20] was calculated by the following equation: Z' = 1 –

Biomedicines **2022**, 10, 422 5 of 12

((3SDmax + 3SDmin)/|mean max - mean min|) and the signal window (SW) [21] was calculated using the equation: SW = ((mean max - 3SDmax) - (mean min + 3SDmin))/SDmax. In both equations, SDmax is the standard deviation of the positive controls, SDmin is the standard deviation of the negative controls, mean max is the mean value of the positive controls and mean min is the mean value of the negative controls. CCL20 at 60 nM (approximately EC80) was used as a positive control and the assay buffer (BSS) as a negative control.


#### 2.8. Statistical Analysis

Data were represented as mean  $\pm$  SD. The statistical analysis was carried out using GraphPad Prism version 8.0.1 software (GraphPad Software Inc., San Diego, CA, USA). Non-linear regression analyses were performed to generate dose-response curves and calculate EC50 and EC80 values. Linear regression was used to analyze data reproducibility. Statistical significance was evaluated with Student's *t*-test and are indicated with \* p < 0.05; \*\* p < 0.01; \*\*\* p < 0.001; and \*\*\*\* p < 0.0001.

#### 3. Results

#### 3.1. Development of Stable Cell Line Expressing Human HexB

To obtain a cell line with hexosaminidase activity, we transfected RBL-2H3 cells with a plasmid for the expression of human HexB fused to IRES-NGFR cassette (Figure 1a). Transfected cells were kept in antibiotic selection with hygromycin for two weeks, after which, cells were tested in flow cytometry. Positive cells for NGFR expression were selected and cloned to ensure their stability and, growing colonies, were analyzed for NGFR expression again.

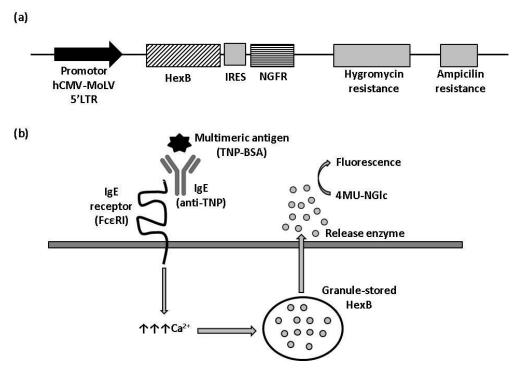

**Figure 1.** β-hexosaminidase cell sensor-based assay. (a) Vector construction. Plasmid vector with hygromycin resistance cassette used to stably express HexB fused to NGFR under control hCMV-MoLV-5'LTR promoter. (b) Illustration of HexB-based calcium mobilization assay. 4MU-NGlc: 4-Methylumbelliferyl N-acetyl- $\beta$ -D-glucosaminide. The arrows represent an increment of calcium.

HexB activity of three clones with positive expression for NGFR, named 1B7-RBL-2H3-HexB, 1C4-RBL-2H3-HexB and 1F10-RBL-2H3-HexB, was assayed. We measured the release of hexosaminidase due to interactions between multimeric antigen and sur-


Biomedicines **2022**, 10, 422 6 of 12

face-bound IgE, as described in the Materials and Methods (Figure 1b). A read of fluorescence was taken at different times of incubation, selecting 60 min as the optimal assay time (Table 1).

| Cell                  | Specific Exocytosis (RFU) | Background Release (RFU) |
|-----------------------|---------------------------|--------------------------|
| Mock                  | $6424 \pm 584$            | $2056 \pm 30$            |
| RBL-2H3-HexB          | $14,149 \pm 309$          | $2908 \pm 77$            |
| 1C4-RBL-2H3-HexB      | 39,185 ±1024              | $6302 \pm 132$           |
| 1F10-RBL-2H3-HexB     | $17,378 \pm 362$          | $4394 \pm 83$            |
| 1B7-RBL-2H3-HexB      | $28,554 \pm 790$          | $2611 \pm 61$            |
| Without cells (blank) | $1203 \pm 42$             | 1176 ± 47                |

RFU, relative fluorescence units.

The above results indicate that overexpression of HexB into RBL-2H3 cells produces a functional enzyme, as measured with 4MU-NGlc, that is stored inside the granules and is specifically released by exocytosis through cross-linking of IgE receptor. In the stably transfected population, RBL-2H3-HexB, the specific signal is increased 2.2 times with respect to mock transfectants, while background increased 1.4 times, that is, specific signal increased more than background release and this indicates that the uses of transfected cells as sensors is better than the use of untransfected cells. The specific signal to the background (S/B) ratio of untransfected RBL-2H3 cells in this experiment was 5.9, while the S/B ratio for RBL-2H3-HexB cells was 7.5. Clones selected by limiting dilution from the RBL-2H3-HexB increased specific release, but clones 1C4 and 1F10 also increased background release having a signal-to-background ratio of 7.4 and 5. However, the signal-to-background ratio of clone 1B7 was 19.1, which is an increment of 11.6 times of RBL-HexB ratio and 13.2 times of RBL untransfected cells ratio (Figure 2a). The clone 1B7-RBL-2H3-HexB, with the best signal-to-background ratio, was chosen as our working clone.

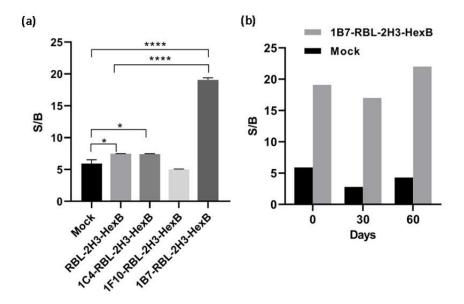

**Figure 2.** Activity of hexosaminidase released by transfected cells. (a) Signal-to-background ratio of mock transfectants, RBL-2H3 cells transfected with HexB and RBL-2H3-HexB clones. Means and SD of triplicates are plotted. (b) Stability of 1B7-RBL-2H3-HexB cells activity in the time. Representative data for one experiment are shown. \* p < 0.1, \*\*\*\* p < 0.0001.


Even more important than the signal reached is the fact that production and release of hexosaminidase by the 1B7-RBL-2H3-HexB cells obtained is extremely regular while

Biomedicines **2022**, 10, 422 7 of 12

that of RBL-2H3 untransfected cells has a very strong variability over time. Both, RBL-2H3 and 1B7-RBL-2H3-HexB, cells were cultured for 2 months, and exocytosis was measured as above every 1 month. The signal-to-background ratio of RBL-2H3 cells was 5.9 (month 0), 2.8 (month 1) and 4.3 (month 2) while the signal-to-background ratio of 1B7-RBL-2H3-HexB cells was 19.1 (month 0), 17 (month 1) and 22 (month 2). A strong natural variability was observed in RBL-2H3 cells with a reduction of 63% of S/B ratio between month 0 and month 1. But 1B7-RBL-HexB cells behaved much better and the S/B ratio, while still variable as corresponds to live cells, was more stable and the maximal variability was 22% between month 1 and month 2 (Figure2b).

# 3.2. Development of Calcium Mobilization Assay for CCR6h Receptor Coupled to $\beta$ -Hexosaminidase Readout


To study CCR6 activation by CCL20 with  $\beta$ -hexosaminidase enzymatic assay, 1B7-RBL-2H3-HexB cells were transfected with an expression plasmid containing viral gly-cosylation sequence (VGS) and c-myc tag fused in frame upstream of CCR6h (Figure 3a). The VGS sequence was used to enhance CCR6h surface expression [19]. Transfected cells were kept in antibiotic selection with neomycin for two weeks, after which cells were tested in flow cytometry. The cells 1B7-RBL-2H3-HexB-CCR6h had 69% c-myc and 99% NGFR expression (Figure 3b).

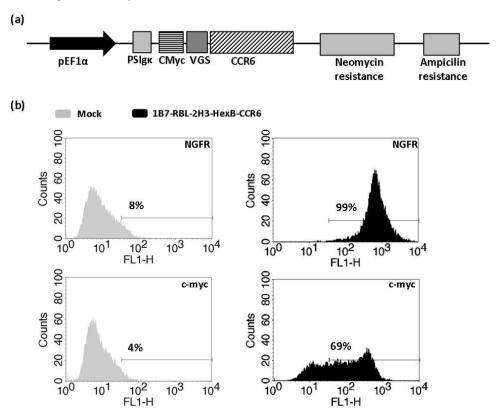

**Figure 3.** Obtaining of 1B7-RBL-2H3-HexB-CCR6h cells. (a) Vector construction. Plasmid vector with neomycin resistance cassette used to stably express CCR6h fused to c-myc on 1B7-RBL-2H3-HexB cells. (b) Flow cytometry. Expression level of stably transfected cells with HexB and CCR6h. One representative analysis from three independent experiments is shown. The lines represent the marker of the positive signal.

The cells were assayed in a functional calcium mobilization assay, based on HexB cells activity, incubating the 1B7-RBL-2H3-HexB-CCR6h cells with serial dilutions of CCL20 to measure specific signal of CCR6 receptor (Figure 4a). 1B7-RBL-2H3-HexB cells without CCR6h were used as negative control. A concentration-response curve was realized showing an EC50 of 10.72 nM (Figure 4b).


Biomedicines **2022**, 10, 422 8 of 12

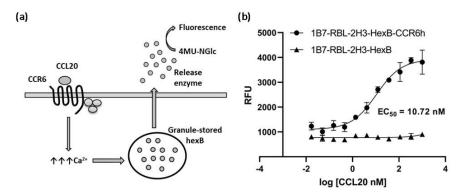

**Figure 4.** Calcium mobilization assay based on cell sensor. (a) Illustration of calcium mobilization assay for CCR6h. The arrows represent an increment of calcium. (b) Concentration-response curve of CCL20 in 1B7-RBL-2H3-HexB-CCR6h cells and 1B7-RBL-2H3-HexB cells. Means and SD of triplicates are plotted. RFU: relative fluorescence units.

# 3.3. Optimization and Performance of HTS Assay

Various experimental conditions were tested to optimize the assay for high-throughput screening (HTS). We found that the times of incubation tested (60, 70, 80 and 90 min) did not affect the EC50 value of the ligand CCL20 and the S/B ratio reached a plateau at 60 min of incubation (Figure 5a,b). Cells dispensing methods were examined by running the calcium mobilization assay with electronic multichannel pipette or manually plated cells. The EC50 value was the double when the cells were dispensed with electronic pipette and the signal-to-background ratio was better than the ratio obtained with manually plated cells (Figure 5c,d). Therefore, the final assay conditions for HTS were determined as follows: the time of ligand incubation was 60 min and the cells were dispensed with electronic pipette.

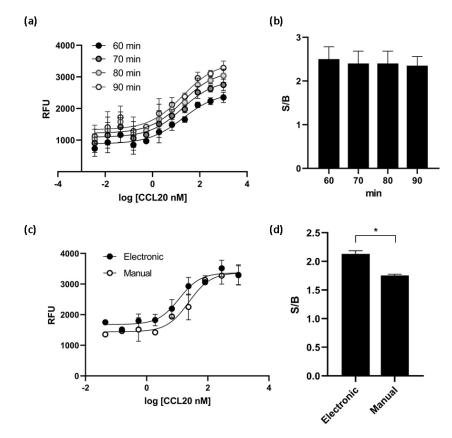


**Figure 5.** CCR6h calcium mobilization assay optimization. At different incubation times, the EC<sub>50</sub> values obtained (**a**) and the S/B ratio (**b**) were almost the same. The cell dispensing method affect

Biomedicines **2022**, 10, 422 9 of 12

to dose response curve (c) and the S/B ratio was significantly higher with electronic method (d). \* p < 0.1. Means and SD of triplicates are plotted. RFU: relative fluorescence units.

To see if this cell-based assay was suitable for HTS, Z' factor and signal window were calculated. Using the fully optimized assay, we generated Z'-factor and signal window of assay running 48 positive controls and 48 negative controls. The average of the positive controls was 2164 RFU with a standard deviation of 107.1 RFU and the average of the negative controls was 848 RF with a standard deviation of 54.9 RFU. As shown in Figure 6a, the Z' factor was 0.63 and the signal window 7.75. A Z'-factor  $\geq$  0.5 and a signal window  $\geq$  2 are appropriate for HTS assay; therefore, our assay is deemed to be highly robust and reproducible, and hence suitable for HTS applications.

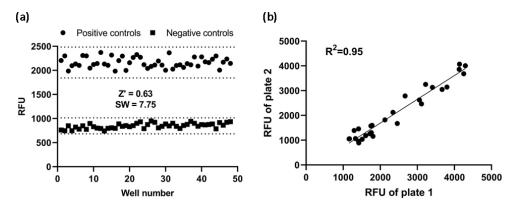

**Figure 6.** Assay performance. (a) Z' factor and signal window (SW) determination. At the fully optimized assay, 48 replicates of positive and negative controls were studied. Dotted lines indicate mean  $\pm$  3SD of 48 data points. (b) Reproducibility. The same concentration of CCL20 was added to corresponding wells of two plates. Reproducibility of data from duplicates plates was investigated with a linear regression analysis. RFU: relative fluorescence units.

Furthermore, to investigate reproducibility between plates, the corresponding wells from two plates were treated with the same concentration of CCL20 and the data were investigated with linear regression analysis. The correlation coefficient was 0.95, showing a high degree of reproducibility between duplicate plates (Figure 6b).

# 4. Discussion

GPCRs represent the largest family of druggable targets due to their numerous physiological and pathological roles, together with their potential for therapeutic intervention via using small molecules as regulators [1]. An important focus of the pharmaceutical industry is the development of robust, reliable, nonradioactive, homogenous, cost effective and easily adapted to a microtiter plate format (96-, 384-, or 1536-well) for robotic automation assays [22]. Although in the past decades signaling-dependent cell-based functional assays have been made, biased GPCR signaling complicates drug discovery efforts. The development of biased ligands as therapeutics heralds an era of increased drug efficacy with reduced drug side effects [23].

The chemokine receptor CCR6 is a GPCR that plays an important role in immunity and is involved in numerous diseases such as inflammatory bowel disease, psoriasis, multiple sclerosis and rheumatoid arthritis [24]. This GPCR is also involved in non-autoimmune diseases such as cancer or atherosclerosis. Thus, inhibition of CCR6 signaling would be an attractive strategy for the treatment of various diseases and a small molecule inhibitor of CCR6 could be useful for both in vitro and in vivo pathophysiological studies [25]. In this manner, the availability of a simple, rapid and robust assay to monitor CCR6h activity would expedite the search for drugs against CCR6.

There are some functional assays to screen compounds that bind to CCR6 receptor based on different biological measurement corresponding to different signaling pathways.

Biomedicines **2022**, 10, 422

Commercially available CCR6 assays include those that measure of cAMP accumulation (cAMP Hunter<sup>TM</sup> eXpress assay), flash increment of Ca<sup>2+</sup> by FLIPR technology or  $\beta$ -arrestin recruitment (PathHunter<sup>®</sup> eXpress CCR6 assay) in different cells: CHO-K1, cell lines like Chem-1 with endogenous high levels of promiscuous  $G_{\alpha 15}$  that promotes coupling  $G_{\alpha 15}$  receptors to an increase in intracellular calcium and CHO-K1/U2OS cell lines, respectively. But to the date, around sixteen CCR6-CCL20 inhibitors have been investigated and no therapeutic agents targeting CCR6 have progressed into clinical evaluation [12]; therefore, it seems to be necessary to develop reliable new technologies for screening of CCR6h.

 $\beta$ -hexosaminidase (HexB) assay is useful for detecting mast cell degranulation upon a short stimulation with drugs or chemical agents [26]. This protein is used as reporter for granule secretion but has been traditionally considered a low sensitivity reporter with a low signal-to-background ratio and a strong variability between experiments due to a large variation over time in the amount of enzyme stored in the granules [27].

In the present investigation,  $\beta$ -hexosaminidase was overexpressed inside the granules of hematopoietic cells, RBL-2H3, reaching an improved signal-to-background ratio assay and a stable production of HexB by transfected cells. While untransfected cells have strong variability, cells transfected with HexB gene were more stable. Thus, such mast cells overexpressing HexB become sensitive cell-based sensors with low variability to measure exocytosis.

To ensure that 1B7-RBL-2H3-HexB cells expressed CCR6 receptor on their surface, the receptor was cloned into the cells. In these cells, when CCL20 activates CCR6 there is a release of HexB, encoded by the gene that was previously overexpressed as describe above, due to the increase in intracellular calcium.

Calcium mediated degranulation of RBL-2H3 cells is induced by GPCR signaling pathways and is associated with  $Ca^{2+}$  mobilization [28]. The release of hexosaminidase contained in the granules was detected by a substrate that measures extracellular  $\beta$ -hexosaminidase activity due to the activation of CCR6 by CCL20 ligand calcium signaling. The untransfected cells 1B7-RBL-2H3-HexB, without expression of CCR6 in a natural manner similar to their parental cell RBL-2H3, had a very low level of hexosaminidase secretion under ligand stimulation, demonstrating that the release of hexosaminidase in transfected cells is specific of CCL20 activation. Therefore, production of  $Ca^{2+}$  in this reporter cell line, the subsequent releasing of HexB and the generation of fluorescent signal can be attributed to the activation of a particular signaling pathway, initiated by the activation of calcium mobilization signaling of CCR6 receptor by CCL20 ligand. Then, like IgE receptor, at least some GPCRs induce degranulation upon activation in RBL-2H3.

The VGS sequence has been reported [19] to enhance surface expression of GPCRs and, therefore, was used to promote CCR6 expression on the plasma membrane. In the sTable 1B7-RBL-2H3-HexB-CCR6 cells the response to CCL20 was potent and had a robust signal with an EC $_{50}$  of 10.72 nM.

The assay was then optimized, and it tested the time of ligand incubation and the method of cell dispensing. The results indicate that the assay performance is well within 60 min of CCL20 incubation and electronic cell dispensing.

With the optimized conditions, we calculated some parameters commonly used as tools to measure assay performance during validation of an assay for HTS. The Z' factor is a useful tool used to assess the robustness of assay for screening [20] and the signal window (SW) provides a degree of separation between signals [21]. Although both parameters measure that assay signal adjusted for assay variability, is described that Z' has better precision properties than SW to measure assay performance [29]. Nevertheless, we calculated both parameters for a better support to the method developed. In general, a Z' value greater than or equal than 0.5 suggests that separation band is large and that the assay is robust enough for HTS and a signal window above 2 is recommended. Therefore, the assay was deemed robust and reproducible with a Z' value of 0.63 and a SW value of 7.75. Moreover, the correlation of two plates of assay confirms the excellent reproducibility.

Biomedicines **2022**, 10, 422

An ideal GPCR screening technology should be simple, nonradioactive, with high signal-to-noise ratio, homogeneous, with minimal reagents additions and be amenable to a microtiter plate format to facilitate robotic automation [16]. We used the HexB stored inside granules of hematopoietic cells to detect CCL20 induced exocytosis on CCR6 overexpressing cells and obtain an assay with a high signal-to-background ratio and low inter-assay variability. The assay performed very well, had a low baseline response, was reproducible and was suitable for high throughput screens. The sensor of the present investigation is faster than sensors based on inducible promoters, no lysis is needed for release of reporters, no washing or stop steps are needed thus increasing throughput, is stable and high signal to background is obtained for a robust assay with low variability between inter-assay experiments.

In theory, this calcium mobilization assay could be used to identify other ligands as well, including antagonist, inverse agonist or allosteric modulators, although further studies are needed to confirm this hypothesis. Moreover, these applications could be expanded to other GPCRs, provided their activation can also be coupled to an increase in calcium release.

Thus, we established a cell-based HTS assay system for monitoring calcium mobilization due to CCR6h activation by CCL20 using  $\beta$ -hexosaminidase activity. This technology represents a valuable tool in screening, counter screening and SAR characterization of compounds that modulate CCR6 activity.

#### 5. Conclusions

In the present report, we generated a cell-based calcium mobilization assay for CCR6h receptor. The calcium mobilization assay established here, using measurement of  $\beta$ -hexosaminidase induced by the activation of 1B7-RBL-2H3-HexB-CCR6 cells calcium signaling by CCL20, can be reliably used to screen and study modulators of CCR6h.

This assay exhibited good values of Z' and signal window and- thus can be applied in high-throughput screening formats. It has been demonstrated that the assay described here is robust, time saving and cost effective, generating highly reproducible results. Moreover, the assay could be easily extended to study modulators of other GPCRs. In conclusion, our approach provides an alternative strategy for ligand screening of CCR6h.

**Author Contributions:** Conceptualization, E.P.-R. and F.I.G.-M.; methodology, S.G.-M.; software, T.G.-M.; validation, I.T., E.P.-R. and F.I.G.-M.; formal analysis, J.C.-V.; investigation, V.L.-G. and S.G.-M.; resources, G.M.-P. and S.G.-M.; data curation, T.G.-M. and S.G.-M.; writing—original draft preparation, S.G.-M.; writing—review and editing, E.P.-R. and F.I.G.-M.; visualization, I.T.; supervision, F.I.G.-M.; project administration, E.P.-R.; funding acquisition, F.I.G.-M. All authors have read and agreed to the published version of the manuscript.

**Funding:** This research was funded by Ministry of Economic Affairs and Digital Transformation (MINECO) of Spain, grant number PCB-010000-2010-23.

**Data Availability Statement:** The data that support the findings of this study are available from the corresponding author upon reasonable request.

**Conflicts of Interest:** The authors declare no conflict of interest. The funders had no role in the design of the study; in the collection, analyses, or interpretation of data; in the writing of the manuscript, or in the decision to publish the results.

#### References

- 1. Zhang, R.; Xie, X. Tools for GPCR drug discovery. Acta Pharmacol. Sin. 2012, 33, 372–384. https://doi.org/10.1038/aps.2011.173.
- 2. Nelson, A.L.; Dhimolea, E.; Reichert, J.M. Development trends for human monoclonal antibody therapeutics. *Nat. Rev. Drug-Discov.* **2010**, *9*, 767–774. https://doi.org/10.1038/nrd3229.
- 3. Wedegaertner, P.B.; Wilson, P.T.; Bourne, H.R. Lipid modifications of trimeric G proteins. J. Biol. Chem. 1995, 270, 503–506. https://doi.org/10.1074/jbc.270.2.503 10.
- 4. Neves, S.R.; Ram, P.T.; Iyengar, R. G protein pathways. Science 2002, 296, 1636–1639. https://doi.org/10.1126/science.1071550 11.
- 5. Strathmann, M.; Simon, M.I. G protein diversity: A distinct class of alpha subunits is present in vertebrates and invertebrates. *Proc. Natl. Acad. Sci. USA* **1990**, *87*, 9113–9117. https://doi.org/10.1073/pnas.87.23.9113 12.

Biomedicines **2022**, 10, 422

 Wilkie, T.M.; Scherle, P.A.; Strathmann, M.P.; Slepak, V.Z.; Simon, M.I. Characterization of G-protein alpha subunits in the Gq class: Expression in murine tissues and in stromal and hematopoietic cell lines *Proc. Natl. Acad. Sci. USA* 1991, 88, 10049–10053. https://doi.org/10.1073/pnas.88.22.10049.

- 7. Woszczek, G.; Fuerst, E. Ca<sup>2+</sup> mobilization assays in GPCR drug discovery. *Methods Mol. Biol.* **2015**, 1272, 79–89. https://doi.org/10.1007/978-1-4939-2336-6\_6.
- 8. Zhu, T.; Fang, L.; Xie, X. Development of a universal high-throughput calcium assay for G-protein-coupled receptors with promiscuous G-protein Gα15/16. *Acta Pharmacol. Sin.* **2008**, 29, 507–516. https://doi.org/10.1111/j.1745-7254.2008.00775.x.
- 9. Liu, K.; Southall, N.; Titus, S.A.; Inglese, J.; Eskay, R.L.; Shinn, P.; Austin, C.P.; Heilig, M.N.; Zheng, W. A multiplex calcium assay for identification of GPCR agonists and antagonists. *Assay Drug Dev. Technol.* **2010**, *8*, 362–374. https://doi.org/10.1089/ADT.2009.0245.
- Ma, Q.; Ye, L.; Liu, H.; Shi, Y.; Zhou, N. An overview of Ca<sup>2+</sup> mobilization assays in GPCR drug discovery. *Expert Opin. Drug Discov.* 2017, 12, 511–523. https://doi.org/10.1080/17460441.2017.1303473.
- 11. Ranasinghe, R.; Eri, R. Pleiotropic Immune Functions of Chemokine Receptor 6 in Health and Disease. *Medicines* **2018**, *5*, 69. https://doi.org/10.3390/medicines5030069.
- 12. Ranasinghe, R.; Eri, R. Modulation of the CCR6-CCl20 axis: A potential therapeutic target in inflammation and cancer. *Medicina (Lithuania)* **2018**, *54*, 88. https://doi.org/10.3390/medicina54050088.
- 13. Puri, N.; Roche, P.A. Mast cells possess distinct secretory granule subsets whose exocytosis is regulated by different SNARE isoforms. *Proc. Natl. Acad. Sci. USA* **2008**, *105*, 2580–2585. https://doi.org/10.1073/PNAS.0707854105.
- 14. Gómez-Melero, S.; García-Maceira, F.I.; García-Maceira, T.; Luna-Guerrero, V.; Montero-Peñalvo, G.; Túnez-Fiñana, I.; Paz-Rojas, E. Amino terminal recognition by a CCR6 chemokine receptor antibody blocks CCL20 signaling and IL-17 expression via β-arrestin. *BMC Biotechnol.* **2021**, 21, 41. https://doi.org/10.1186/S12896-021-00699-2.
- 15. Heydorn, A.; Ward, R.J.; Jorgensen, R.; Rosenkilde, M.M.; Frimurer, T.M.; Milligan, G.; Kostenis, E. Identification of a Novel Site within G Protein α Subunits Important for Specificity of Receptor-G Protein Interaction. *Mol. Pharmacol.* **2004**, *66*, 250–259. https://doi.org/10.1124/mol.66.2.250.
- 16. Kostenis, E.; Martini, L.; Ellis, J.; Waldhoer, M.; Heydorn, A.; Rosenkilde, M.M.; Norregaard, P.K.; Jorgensen, R.; Whistler, J.L.; Milligan, G. A Highly Conserved Glycine within Linker I and the Extreme C Terminus of G Protein α Subunits Interact Cooperatively in Switching G Protein-Coupled Receptor-to-Effector Specificity. *J. Pharmacol. Experim. Ther.* 2005, 313, 78–87. https://doi.org/10.1124/jpet.104.080424.
- 17. Offermanns, S.; Simon, M.I.  $G\alpha$ 15 and  $G\alpha$ 16 Couple a Wide Variety of Receptors to Phospholipase C. *J. Biol. Chem.* **1995**, 270, 15175–15180. https://doi.org/10.1074/jbc.270.25.15175.
- 18. Matsubara, M.; Masaki, S.; Ohmori, K.; Karasawa, A.; Hasegawa, K. Differential regulation of IL-4 expression and degranulation by anti-allergic olopatadine in rat basophilic leukemia (RBL-2H3) cells. *Biochem. Pharmacol.* **2004**, *67*, 1315–1326. https://doi.org/10.1016/J.BCP.2003.12.008.
- 19. Rojas, E.P.; Maceira, F.I.G.; Guerrero, V.I.L.; Peñalvo, M.G.M.; Maceira, T.G.; Martínez, J.A.M.; Morales, A.M.M. GPCR with Improved Cell Surface Expression. U.S. Patent 9,499,584, W02013113370, 8 August 2013.
- Zhang, J.H.; Chung, T.D.; Oldenburg, K.R.J. A Simple Statistical Parameter for Use in Evaluation and Validation of High Throughput Screening Assays. J. Biomol. Screen 1999, 4, 67–73. https://doi.org/10.1177/108705719900400206.
- 21. Lister, M.D.; Sittampalam, G.S.; Iversen, P.W.; Boadt, J.A.; Kahl, S.D.; Bright, S.; Janzen, W.P. Design of Signal Windows in High Throughput Screening Assays for Drug Discovery, J. Biomol. Screen. 1997, 2, 159–169. https://doi.org/10.1177/108705719700200306.
- 22. Thomsen, W.; Frazer, J.; Unett, D. Functional assays for screening GPCR targets. *Curr. Opin. Biotechnol.* **2005**, *16*, 655–665. https://doi.org/10.1016/J.COPBIO.2005.10.008.
- 23. Smith, J.S.; Lefkowitz, R.J.; Rajagopal, S. Biased signalling: From simple switches to allosteric microprocessors. *Nat. Rev. Drug Discov.* **2018**, *17*, 243–260. https://doi.org/10.1038/nrd.2017.229.
- 24. Meitei, H.T.; Jadhav, N.; Lal, G. CCR6-CCL20 axis as a therapeutic target for autoimmune diseases. *Autoimmun. Rev.* **2021**, 20, 102846. https://doi.org/10.1016/J.AUTREV.2021.102846.
- 25. Tawaraishi, T.; Sakauchi, N.; Hidaka, K.; Yoshikawa, K.; Okui, T.; Kuno, H.; Chisaki, I.; Aso, K. Identification of a novel series of potent and selective CCR6 inhibitors as biological probes. *Bioorg. Med. Chem. Lett.* **2018**, 28, 3067–3072. https://doi.org/10.1016/J.BMCL.2018.07.042.
- 26. Huang, L.; Pi, J.; Wu, J.; Zhou, H.; Cai, J.; Li, T.; Liu, L. A rapid and sensitive assay based on particle analysis for cell degranulation detection in basophils and mast cells. *Pharmacol. Res.* **2016**, *111*, 374–383. https://doi.org/10.1016/J.PHRS.2016.05.033.
- 27. Paz-Rojas, E.; de Gracia Montero-Peñalvo, M.; Luna-Guerrero, V.I.; Maceira, F.I.G.; Morales-Martínez, J.A.; Maceira, T.G.; Aragón-Gómez, A.B.; Quesada-Molina, A.; Marquez-Morales, A.M. Novel Ultrasensitive Cell based Sensors and Use Thereof. U.S. Patent 9,347,942, W02012013204, 2 February 2012.
- 28. Chen, Y.C.; Chang, Y.C.; Chang, H.A.; Lin, Y.Ś.; Tsao, C.W.; Shen, M.R.; Chiu, W.T. Differential Ca<sup>2+</sup> mobilization and mast cell degranulation by FcεRI- and GPCR-mediated signaling. *Cell Calcium* **2017**, *67*, 31–39. https://doi.org/10.1016/J.CECA.2017.08.002.
- 29. Iversen, P.W.; Eastwood, B.J.; Sittampalam, G.S.; Cox, K.L. A comparison of assay performance measures in screening assays: Signal window, Z' factor, and assay variability ratio. *J. Biomol. Screen.* **2006**, 11, 247–252. https://doi.org/10.1177/1087057105285610.