



Feasibility of Rehabilitation Training With a Newly Developed Wearable Robot for Patients With Limited Mobility

著者	Kubota Shigeki, Nakata Yoshio, Eguchi Kiyoshi, Kawamoto Hiroaki, Kamibayashi Kiyotaka, Sakane Masataka, Sankai Yoshiyuki, Ochiai Naoyuki
journal or publication title	Archives of physical medicine and rehabilitation
volume	94
number	6
page range	1080-1087
year	2013-06
権利	(C) 2013 by the American Congress of Rehabilitation Medicine
URL	http://hdl.handle.net/2241/119563

doi: 10.1016/j.apmr.2012.12.020

Running Title: A feasibility study of HAL Rehabilitation

Article Title: Feasibility of rehabilitation training with a newly developed wearable robot for patients with limited mobility

Grant Support: This study was supported by the "Center for Cybernetics Research (CCR) - World Leading Human-Assistive Technology Supporting a Long-Lived and Healthy Society" granted through "Funding Program for World-Leading Innovative R&D on Science and Technology (FIRST Program)" initiated by the Council for Science and Technology Policy (CSTP).

COI: A commercial party having a direct financial interest in the results of the research supporting this article has conferred or will confer a financial benefit on 1 or more of the authors. Yoshiyuki Sankai is CEO of Cyberdyne Inc, Ibaraki, Japan. Hiroaki Kawamoto is a stockholder of the company. Cyberdyne is the manufacturer of the robot suit HAL. This study was proposed by the authors. Cyberdyne was not directly involved in the study design; collection, analysis, or interpretation of data; writing the report; or the decision to submit the paper for publication.

Financial Disclosure: No commercial party having a direct financial interest in the results of the research supporting this article has or will confer a benefit on the authors or on any organization with which the authors are associated (Kubota, Nakata, Eguchi, Kamibayashi, Sakane, Ochiai).

Feasibility of rehabilitation training with a newly developed wearable robot for patients with limited mobility

ABSTRACT

Objectives: To investigate the feasibility of rehabilitation training with a new wearable robot

Design: Before-after clinical intervention

Setting: University Hospital and private rehabilitation facilities

Participants: A convenience sample of 38 patients with limited mobility. The underlying diseases were stroke (n = 12), spinal cord injuries (n = 8), musculoskeletal diseases (n = 4), and other diseases (n = 14).

Interventions: The patients received 90-minute training with a wearable robot twice per week for 8 weeks (16 sessions).

Main Outcome Measures: Functional ambulation was assessed with the 10-m walk test (10MWT) and the timed-up and go (TUG) test, and balance ability was assessed with the Berg balance scale (BBS). Both assessments were performed at the baseline and after the rehabilitation.

Results: Thirty-two patients completed 16 sessions of the training with the wearable robot. The results of the 10MWT included significant improvements in gait speed, number of steps, and cadence. Although improvements were observed, as measured with the TUG test and BBS, the results were not statistically significant. No serious adverse events were observed during the training.

Conclusions: Eight weeks of rehabilitative training with the wearable robot (16 sessions of 90 minutes) could be performed safely and effectively, even many years after the subjects received their diagnosis.

Key Words: Rehabilitation, Robotics, Feasibility Studies, Orthopedic Equipment, Mobility
Limitation

ABBREVIATIONS

HAL: Hybrid assistive limb

10MWT: 10-m walk test

TUG: Timed up-and-go

BBS: Berg balance scale

SCI: Spinal cord injury

MADS: Musculoskeletal Ambulation Disability Symptom Complex

CVC: Cybernic Voluntary Control

CAC: Cybernic Autonomous Control

2

3 Rehabilitation robotics emerged in the 1980s with the aim of using robotic technology to assist
4 people with movement dysfunction.¹ Robotic devices have recently been developed for use in
5 clinical settings. Tefertiller et al. reviewed 30 articles (14 randomized controlled trials, 16 non-
6 randomized controlled trials) that examined the effects of locomotor training with robotic
7 assistance in patients following stroke, spinal cord injury (SCI), multiple sclerosis, traumatic
8 brain injury, and Parkinson's disease. The review supports the conclusion that locomotor training
9 with robotic assistance is beneficial for improving walking function in individuals following
10 stroke and SCI.² The development of main gait training machines followed. These machines
11 either involve an exoskeleton robotic device (e.g., Lokomat, LOPES exoskeleton robot)³⁻⁴ or a
12 robotic device with foot-driven plates (e.g., Gait Trainer GT I, Haptic Walker).⁵⁻⁶ The
13 exoskeleton robotic device is equipped with programmable drives or passive elements that flex
14 the knees and hips during the swing phase, whereas with the other type of robotic device, the feet
15 are placed on foot plates, whose trajectories simulate the stance and swing phases. Other than
16 robotic gait training and conventional therapy, another treatment approach involves treadmill
17 training with partial body weight support (BWSTT).⁷ However, this approach requires
18 considerable involvement of a physical therapist, and, generally, 3 therapists are required to
19 induce movement of the paretic leg during the swing phase and to shift the patient's weight onto
20 the stance limb.

21

22 The potentially positive common benefits of robotic gait training are that it involves repeatedly
23 undergoing sufficient and accurate training for a prolonged period. Lokomat (Hocoma,
24 Volketswil, Switzerland) is the first robotic-driven gait orthosis with electromechanical drives to
25 assist the walking movements of gait-impaired patients on a treadmill by supporting the body

26 weight.⁸⁻⁹ Husemann et al. compared a Lokomat group that received 30 minutes of robotic
27 training with a control group that received 30 minutes of conventional physiotherapy¹⁰. After 4
28 weeks of therapy, although there was no significant difference in walking ability between the
29 groups, the walking ability in both groups as expressed by functional ambulation classification
30 was significantly improved. The researchers reported that the Lokomat group demonstrated an
31 advantage for robotic training over conventional physiotherapy in the improvement of gait
32 abnormality and body tissue composition.¹⁰ However, in a recent randomized controlled study
33 that compared robot-assisted locomotor training with therapist-assisted locomotor training in
34 chronic stroke patients, the results indicated that greater improvements in speed and single limb
35 stance time on the impaired leg were observed in subjects who received therapist-assisted
36 locomotor training.¹¹ Thus, the usefulness of robot-assisted rehabilitation is controversial.

37
38 The robot suit hybrid assistive limb (HAL)^{12-15,a} is a new wearable robot that has a hybrid
39 control system comprised of two subsystems: Cybernic Voluntary Control (CVC) and Cybernic
40 Autonomous Control (CAC) (fig 1). The HAL suit has power units and force-pressure sensors in
41 the shoes. The power units consist of angular sensors and actuators on bilateral hip and knee
42 joints. Muscle action potentials are detected through the electrodes on the anterior and posterior
43 surface of the wearer's thigh. These various biological signals are processed by a computer. The
44 HAL suit can support the wearer's motion by adjusting the level and timing of the assistive
45 torque provided to each joint according to the surface muscle action potential as well as the
46 pressure sensors. The HAL suit can enhance the wearer's motion through the wearer's muscle
47 action potential; thus, the HAL suit can appear as an actual motion. Therefore, if the wearer's
48 muscle action potential varies, the wearer's motion varies, too. The HAL training, using muscle

49 activity, has the potential to intensify the feedback by evoking by an appropriate motion more
50 strongly than standard robot training. Thus, after HAL training, patients with limited mobility
51 will improve their walking abilities (gait speed, number of steps, cadence, or ability to transfer).

52

53 Few studies have been conducted to clarify the feasibility of rehabilitation with HAL. Only 1
54 preliminary study has reported on the short-term effects of HAL on the walking pattern of stroke
55 patients.¹⁶ The purpose of the present study was to investigate the feasibility of 16-session (8-
56 week) HAL rehabilitation training for patients with limited mobility.

57

58 **METHODS**

59 **Study design**

60

61 A quasi-experimental study was utilized, with measurements before and after the clinical
62 intervention. The target population included patients with limitations in their walking (no matter
63 the diagnosis, the time since the diagnosis and diagnosis age). The protocol of this study was
64 approved by the Institutional Review Board of the University of Tsukuba Hospital and was
65 registered with the UMIN Clinical Trials Registry (UMIN000002969). The clinical intervention
66 was conducted at the University of Tsukuba Hospital and Cyberdyne, Inc. in Japan between
67 January 2010 and March 2012. The patients included in this study were volunteers recruited
68 through local newspaper advertisements or outpatients at the University of Tsukuba Hospital.
69 They were informed about the aim and design of this study, and they subsequently provided
70 written, informed consent. Informed consent was also obtained from the patient's guardian if the
71 patient was younger than 20 years old.

72

73

74 The inclusion criteria were (1) musculoskeletal ambulation disability symptom complex
75 (MADS) or the underlying disorders of MADS, which is a condition newly defined in 2006 by
76 Japanese medical societies;¹⁷ (2) requiring physical assistance or assistive devices in at least one
77 of the following daily activities: standing up, sitting down, and walking; (3) ability to understand
78 an explanation of the study and to express consent or refusal; (4) body size that can fit in the
79 robotic suit HAL (height range of 145 - 180 cm, and maximal body weight of 80 kg); and (5)
80 ability to undergo usual physical and occupational therapies. The exclusion criteria were the
81 following: (1) inadequately controlled cardiovascular disorders; (2) inadequately controlled
82 respiratory disorders; (3) intellectual impairments that limit the ability to understand instructions;
83 (4) moderate to severe articular disorders, including contracture in the lower extremities; (5)
84 moderate to severe involuntary movements, ataxia, or impairments of postural reflex in the trunk
85 or the lower extremities; and (6) severe spasticity in the lower extremities.

86

87 **Participants**

88

89 Thirty-eight patients (25 men, 13 women) were enrolled in this study (24 outpatients, 14
90 volunteers through advertisements). The mean age of the 38 patients was 53.2 ± 17.8 y, with a
91 range of 18–81 y. Table 1 summarizes their clinical characteristics. Their underlying diseases
92 were stroke (10 men, 2 women), SCI (6 men, 2 women), musculoskeletal diseases (2 men, 2
93 women), and other diseases (Parkinson's disease, gonadotropin-dependent myopathy, limb-
94 girdle muscular dystrophy, inclusion body myositis, traumatic brain injury, disuse syndrome,

95 secondary to malignant lymphoma, cerebral palsy, sequelae of poliomyelitis, and hypoxic
96 ischemic encephalopathy; 7 men, 7 women). Twenty patients were able to ambulate
97 independently without any help (n = 9) or with several assistive devices (t-cane, bilateral
98 crutches, or lateral crutch) (n = 11). Eleven patients were able to ambulate with several assistive
99 devices and under supervision. Three patients required human assistance to ambulate at least 10
100 m (cases 33, 34, and 38), and the remaining 4 patients were unable to ambulate even with
101 assistive devices and human assistance (cases 8, 15, 17, and 27). All of the stroke and SCI
102 patients were in chronic stages.

103

104 **Training Program**

105

106 HAL training was administered twice per week for 8 weeks (16 sessions). The 90-minute
107 training sessions consisted of single-leg motion, a standing and sitting exercise, and walking on
108 the ground with HAL. For safety reasons, a walking device (All-in-One Walking Trainer^b) with a
109 harness was used. Treadmill training with mild body-weight support (Unweighing System^c) was
110 also used for some patients. The HAL suit has a hybrid control system comprising the CVC and
111 CAC. The CVC mode of the HAL suit can support the patient's voluntary motion according to
112 the voluntary muscle activity and the assistive torque provided to each joint. The CAC mode
113 provides physical support autonomously, based on output from force-pressure sensors in the
114 shoes. This study mainly used the CVC mode, which allows the operator to adjust the degree of
115 physical support to the patient's comfort and gradually reduce support as training progresses.

116

117 **Outcome Measures**

118

119 The feasibility of HAL rehabilitation with HAL was assessed by the number of completers and
120 the amount of time or the number of therapists needed to implement training. Patients were asked
121 to report adverse events during the training period.

122

123 The primary outcomes were functional ambulation and balance ability. Functional ambulation
124 was assessed with a 10-m walk test (10MWT) and a timed up-and-go (TUG) test. In the 10MWT,
125 patients were instructed to walk without wearing HAL on a flat surface at their self-selected
126 comfortable pace. Patients began to walk before they reached the starting line of the 10 m
127 distance so that they could accelerate and attain a stable speed before the test. To calculate gait
128 speed (m/s) as a primary outcome, the 10-m walking time was measured using a handheld
129 stopwatch. In addition, the number of steps between the start and finish line was counted, and
130 patient cadence was calculated from the walking time and number of steps. Patients were
131 allowed to use their assistive device and/or lower limb orthosis as necessary. Each patient used
132 the same assistive device and/or orthosis during the pre- and post-intervention measurements.
133 Therapists closely attended the patients during the 10MWT but did not provide physical
134 assistance. For each measurement, the 10MWT was performed twice. The faster time of two
135 trials was selected for analysis. In the TUG test, the following actions were timed: standing up
136 from a standard-height chair, walking 3 m, returning to the chair, and sitting down without HAL.
137 Two trials (each turning clockwise and counterclockwise) were carried out for each measurement.
138 Balance ability was assessed with the Berg balance scale (BBS), consisting of 14 tasks, as
139 detailed by Berg et al.¹⁸ Each task was scored on a scale ranging from 0 to 4 points (0 indicates

140 inability to complete), and the total score was used as the index of balance ability. All primary
141 outcomes were assessed at baseline and after completion of the 16 training sessions.

142

143 **Statistical Analysis**

144

145 All parametric data are expressed as means with standard deviations. Paired *t* tests were used to
146 evaluate differences between the baseline measurements and outcomes after the 16 sessions.

147 Unpaired *t* tests were used to evaluate the differences in characteristics of those who completed
148 16 sessions and those who did not. An effect-size calculation (Cohen *d*) was used to assess the
149 effect of the training. Pearson correlation coefficients were used to assess the relationship among
150 outcome measures. Data were analyzed using IBM SPSS statistics 18 software^d, with the alphas
151 level set at 5%.

152

153 **RESULTS**

154

155 A typical 90-minute HAL training session proceeded as follows: assessment of blood pressure,
156 resting heart rate, and walking pattern (10 min); preparation of electrodes and putting on the
157 HAL suit (5 min); computer set-up (5 min); HAL training (60 min, including resting time during
158 computer operation); taking off the HAL suit and the electrodes (5 min); and reassessment of
159 walking pattern (5 min). The net walking time was approximately 20 min. Typically, 2 therapists
160 implemented the training: one supported the patient and the other operated the computer. All
161 therapists and related staff had participated in a 3-h training workshop conducted by the
162 manufacturer to learn how to operate the HAL system.

163

164 Of the 38 patients (25 men, 13 women), 32 (21 men, 11 women) completed all 16 training
165 sessions. The mean age of the 32 patients was 53.2 ± 17.3 y, with a range of 18–81 y. There was
166 no statistically significant difference in age between those who completed training and those who
167 did not (54.0 ± 19.8 y). It took 10.0 ± 3.1 weeks (range, 8–21 weeks) to complete 16 sessions. Of
168 the 6 patients who did not complete the 16 sessions, 2 (cases 15 and 21) dropped out for medical
169 reasons, and 4 (cases 1, 2, 29, and 35) dropped out for personal reasons (difficulty visiting the
170 hospital). One medical reason for dropout was low back pain that developed during the first
171 training session (case 21); the patient withdrew consent at the third session. The other medical
172 reason for dropout was a relapse (after the second session) of neuropathic pain due to SCI (case
173 15); the patient withdrew consent at the fifth session. There were no serious training-related
174 adverse events. One stroke patient (case 7) had knee pain (patellar tendinitis) at home after the
175 15th session but was able to complete the 16th session after 1 month of rest. Another patient with
176 inclusion body myositis (case 31) developed knee pain at home after an early session but was
177 able to complete 16 sessions.

178

179 **Outcome Measures**

180

181 Functional ambulation was not assessed for 5 patients at baseline because 3 were unable to
182 ambulate with any assistance (cases 8, 17, and 27), and the other 2 patients needed considerable
183 human assistance to ambulate (cases 34 and 38). The other 27 patients presented significant
184 improvements ($P < 0.05$) in gait speed, number of steps, and cadence after the 16-session HAL
185 training (10MWT, table 2). Improvements in gait speed, number of steps, and cadence are

186 defined as an increase, a decrease, and an increase in the respective parameters. The mean
187 improvements and effect sizes (Cohen *d*) in gait speed, number of steps, and cadence were 0.09
188 ± 0.11 m/s ($d = 0.82$), 3.0 ± 4.9 steps ($d = 0.61$), and 6.8 ± 7.1 steps/min ($d = 0.96$), respectively.
189 Improvements in gait speed, steps, and cadence were observed in 25, 18, and 25 patients,
190 respectively (figs 2–4). Worsened gait speed and cadence were observed in 2 patients (cases 28
191 and 30). In the regards to the number of steps, we observed no change in 8 patients (cases 3, 5,
192 16, 25, 28, 30, 33, and 37) and increased steps in 1 (case 20). Correlation coefficients for gait
193 speed with number of steps and with cadence were $r = 0.30$ (not significant) and $r = 0.73$ ($P <$
194 0.01), respectively. The effect sizes for gait speed in stroke patients ($n = 9$), SCI patients ($n = 6$),
195 musculoskeletal disease patients ($n = 3$), and patients with other diseases ($n = 9$) were 1.41, 0.78,
196 2.43, and 0.63, respectively. The results of the TUG test ($n = 26$; the patient in case 10 was
197 unable to perform the test) and the BBS ($n = 32$) indicated improvement after the 16 training
198 sessions, but these improvements were not statistically significant. The mean decrease (Cohen *d*)
199 in the TUG test was 6.4 ± 16.4 s ($d = 0.39$). Twenty-one of 26 patients were faster after training,
200 and 5 patients were slower (cases 5, 13, 30, 31, and 36) (fig 5). The mean increase (Cohen *d*) in
201 BBS was 1.9 ± 5.5 ($d = 0.35$). Nineteen of 32 patients had higher scores compared to baseline;
202 no change was observed in 6 (cases 12, 17, 23, 27, 36, and 37), and 7 had lower scores (cases 11,
203 16, 26, 30, 31, 32, and 34) (fig 6).

204

205 **DISCUSSION**

206

207 We investigated the feasibility of rehabilitation using a robot suit HAL. We demonstrated that
208 HAL rehabilitation could be implemented safely and effectively. Although a few patients

209 developed lumbar or knee pain during the training, no serious training-related adverse events
210 occurred. Significant improvements in gait speed, number of steps, and cadence were observed,
211 as assessed by the 10MWT. Improved TUG test and BBS results were also observed, but
212 because of the small sample size of this pilot study, these improvements were not statistically
213 significant. Overall, our results suggest that HAL rehabilitation has the potential to improve
214 ambulation in patients with limited mobility.

215
216 Two patients (cases 15 and 21) dropped out for medical reasons. One developed lumbar pain
217 (case 21), and 1 experienced a relapse of neuropathic pain due to SCI (case 15). Although it is
218 unclear whether there was a causal relationship between HAL training and the pain that
219 developed, the lumbar pain in case 21 had been persistent before the HAL training and even after
220 the training ended, and the neuropathic pain in case 15 followed a previous pattern of symptom
221 flares associated with seasonal change. Therefore, it is likely that HAL training did not directly
222 cause the pain that developed in these 2 cases. Two other patients complained of knee pain
223 during the training period, but this pain was not severe, and the patients were able to complete
224 the training. Although, once again, direct causality is unclear, safe implementation of HAL
225 rehabilitation requires adequate caution on the part of therapists and self-awareness on the part of
226 patients who have lumbar and knee pain. Regarding feasibility, approximately 10 min was
227 required for 2 to 3 therapists to put electrodes and the HAL suit on or off the patient. This
228 procedure is a slight inconvenience to address but not a major obstacle to HAL rehabilitation.

229
230 Significant improvements in functional ambulation were observed, and the effect sizes (Cohen d)
231 for gait speed, number of steps, and cadence were 0.82, 0.61, and 0.96, respectively. The

232 correlation coefficient for gait speed with cadence was higher than that of gait speed with steps (r
233 = 0.73 vs. $r = 0.30$). Therefore, the improvement in gait speed with HAL training was mainly
234 brought about by improvement in cadence. That is, HAL training improved stride frequency
235 more than stride length. This finding agrees with a previous robotic training study.¹⁹ The effect
236 sizes for the TUG test and BBS were smaller than that effect sizes for the 10MWT. This result
237 seems to occur because the TUG test and BBS involve complicated motions such as moving
238 from sitting to standing, walking and returning, reaching forward, and alternating feet on each
239 step. The effect sizes for gait speed in 9 stroke patients and in 6 SCI patients were large (1.41
240 and 0.78), respectively. Therefore, training effectiveness in stroke and SCI patients can be
241 expected. The effect size in 3 patients with musculoskeletal diseases was also large (2.43), but
242 the number of patients was small. Therefore, further studies are needed. In this study, we
243 recruited patients with a wide range of stroke and SCI severities. Future studies should examine
244 the influence of the severity of stroke and SCI on the effectiveness of HAL rehabilitation.

245
246 Many recent studies have reported the efficacy of robot-assisted rehabilitation. It is very difficult
247 to directly compare these studies and our study, due to differences in diseases, severity and
248 duration of the disorder, robotic features, methods of intervention, and outcome measures.²⁰ Wirz
249 et al. reported that after locomotor training with Lokomat, the 10MWT gait speed of 20 patients
250 with chronic incomplete SCI increased by 0.11 ± 0.10 m/s ($d = 1.10$).²¹ The number of SCI
251 patients in our study was limited to 6, but our results also indicate the efficacy of HAL
252 rehabilitation for these patients ($d = 0.78$). Hornby et al. reported that after robotic-assisted
253 locomotor training, the gait speed in chronic stroke patients increased by 0.07 ± 0.07 m/s ($d =$
254 1.0).¹¹ Our results also indicate the efficacy of HAL rehabilitation for 9 chronic stroke patients (d

255 = 1.41). We conjectured that the mechanism of this recovery of functional ambulation was due to
256 changes in plasticity in the spinal cord and supraspinal centers. Appropriate sensory inputs, such
257 as maximum weight loading, facilitating proper trunk posture, and hip extension, are essential for
258 maximizing functional recovery.²² Our experience with HAL indicates that the HAL-induced
259 motion might evoke the sensory input, which has a favorable feedback effect on the central
260 nervous system for a recovery of locomotor function. In addition, even if a patient's condition
261 were too severe for medical therapists to provide adequate rehabilitation training, HAL might
262 still make adequate training possible. HAL is a robotic device with potential rehabilitation
263 applications that are dependent on the physical support it can provide.

264

265 **Study Limitations**

266

267 This study was not a randomized controlled trial and could not compare the efficacy of HAL
268 training with conventional rehabilitation. Second, long-term efficacy was not assessed after HAL
269 training. Third, this study could not exclude observer bias and subject bias because the same staff
270 implemented assessment and training, and approximately half of the patients were recruited
271 through local newspaper advertisements. Finally, the statistical power was low because of the
272 small number of patients with each disease.

273

274 **CONCLUSIONS**

275

276 This quasi-experimental study revealed the feasibility of HAL training for rehabilitating patients
277 with limited mobility. This study has shown that it is possible to manage 8 weeks of

278 rehabilitation with HAL training (16 sessions of 90 minutes) safely and effectively, even with
279 persons who received their diagnosis many years ago. After HAL training, significant
280 improvements in gait speed, number of steps, and cadence were observed. Although
281 improvements were observed in the TUG test and BBS, they were not statistically significant.
282 There were no serious adverse events. Further studies are needed to compare the effectiveness of
283 HAL training and conventional rehabilitation.

284

285

286 **References**

287

- 288 1. van Vliet P, Wing AM. A new challenge: robotics in the rehabilitation of the neurologically
289 motor impaired. *Phys Ther* 1991;71:39-47.
- 290 2. Tefertiller C, Pharo B, Evans N, Winchester P. Efficacy of rehabilitation robotics for walking
291 training in neurological disorders: a review. *J Rehabil Res Dev* 2011;48:387-416.
- 292 3. Colombo G, Joerg M, Schreier R, Dietz V. Treadmill training of paraplegic patients using a
293 robotic orthosis. *J Rehabil Res Dev* 2000;37:693-700.
- 294 4. Veneman JF, Kruidhof R, Hekman EE, Ekkelenkamp R, Van Asseldonk EH,
295 van der Kooij H. Design and evaluation of the LOPES exoskeleton robot for interactive gait
296 rehabilitation. *IEEE Transactions Neural Systems Rehabilitation Engineering* 2007;15:379-
297 86.
- 298 5. Hesse S, Uhlenbrock D, Werner C, Bardeleben A. A mechanized gait trainer for restoring gait
299 in nonambulatory subjects. *Arch Phys Med Rehabil* 2000;81:1158-61.
- 300 6. Schmidt H, Werner C, Bernhardt R, Hesse S, Kruger J. Gait rehabilitation machines based on
301 programmable footplates. *J Neuroeng Rehabil* 2007;4:2.
- 302 7. Visintin M, Barbeau H, Korner-Bitensky N, Mayo NE. A new approach to retrain gait in
303 stroke patients through body weight support and treadmill stimulation. *Stroke* 1998;29:1122-
304 8.
- 305 8. Lünenburger L, Colombo G, Riener R, Dietz V. Biofeedback in gait training with the robotic
306 orthosis Lokomat. *IEEE Engineering Medicine Biology Society* 2004;7:4888-91.
- 307 9. Neckel N, Wisman W, Hidler J. Limb alignment and kinematics inside a lokomat robotic
308 orthosis. *IEEE Engineering Medicine Biology Society* 2006;1:2698-701.

- 309 10. Husemann B, Müller F, Krewer C, Heller S, Koenig E. Effects of locomotion training with
310 assistance of a robot-driven gait orthosis in hemiparetic patients after stroke: a randomized
311 controlled pilot study. *Stroke* 2007;38:349-54.
- 312 11. Hornby TG, Campbell DD, Kahn JH, Demott T, Moore JL, Roth HR. Enhanced gait-related
313 improvements after therapist- versus robotic-assisted locomotor training in subjects with
314 chronic stroke: a randomized controlled study. *Stroke* 2008;39:1786-92.
- 315 12. Kawamoto H, Sankai Y. Power assist method based on phase sequence and muscle force
316 condition for HAL. *Adv Robot* 2005;19:717-34.
- 317 13. Lee S, Sankai Y. Virtual impedance adjustment in unconstrained motion for an exoskeletal
318 robot assisting the lower limb. *Abv Robot* 2005;19:773-95.
- 319 14. Suzuki K, Gouji M, Kawamoto H, Hasegawa Y, Sankai Y. Intention-based walking support
320 for paraplegia patients with robot suit HAL. *Adv Robot* 2007;21:1441-69.
- 321 15. Tsukahara A, Kawanishi R, Hasegawa Y, Sankai Y. Sit-to-stand and stand-to-sit transfer
322 support for complete paraplegic patients with robot suit HAL. *Adv Robot* 2010;24:1615-38.
- 323 16. Maeshima S, Osawa A, Nishio D, Hirano Y, Takeda K, Kigawa H, Sankai Y. Efficacy of a
324 hybrid assistive limb in post-stroke hemiplegic patients: a preliminary report. *BMC Neurol*
325 2011;11:116.
- 326 17. Michikawa T, Nishiwaki Y, Takebayashi T, Toyama Y. One-leg standing test for elderly
327 populations. *J Orthop Sci* 2009;14:5:675-85.
- 328 18. Berg K, Wood-Dauphinee S, Williams JI. The balance scale: reliability assessment with
329 elderly residents and patients with an acute stroke. *Scand J Rehabil Med* 1995;27:27-36.
- 330 19. Nooijen CF, Ter Hoeve N, Field-Fote EC. Gait quality is improved by locomotor training in
331 individuals with SCI regardless of training approach. *J Neuroeng Rehabil* 2009;6:36.

- 332 20. Hachisuka K. Robot-aided training rehabilitation. *Brain Nerve* 2010;62:2:133-40 (in
333 Japanese).
- 334 21. Wirz M, Zemon DH, Rupp R, Scheel A, Colombo G, et al. Effectiveness of automated
335 locomotor training in patients with chronic incomplete spinal cord injury: a multicenter trial.
336 *Arch Phys Med Rehabil* 2005;86:672-80.
- 337 22. Barbeau H. Locomotor training in neurorehabilitation: emerging rehabilitation concepts.
338 *Neurorehabil Neural Repair* 2003;17:3-11.
- 339

340 **Suppliers**

341 a. Cyberdyne Inc, D25-1, Gakuen Minami, Tsukuba, Ibaraki, Japan 305-0818.

342 b. ROPOX A/S, 221 Ringstedgade, Naestved, Denmark 4700.

343 c. Biodex Medical Systems Inc, 20 Ramsay Rd, Shirley, NY 11967.

344 d. SPSS, Inc, 233 S Wacker Dr, 11th Fl, Chicago, IL 60606.

345

346 **FIGURE LEGENDS**

347 Figure 1. The robot suit HAL.

348 Figure 2. Change in 10MWT gait speed for 27 patients after HAL training.

349 Figure 3. Change in number of steps during 10MWT for 27 patients after HAL training.

350 Figure 4. Change in 10MWT cadence for 27 patients after HAL training.

351 Figure 5. Change in TUG test results for 26 patients after HAL training.

352 Figure 6. Change in BBS score for 32 patients after HAL training.



Figure 1.

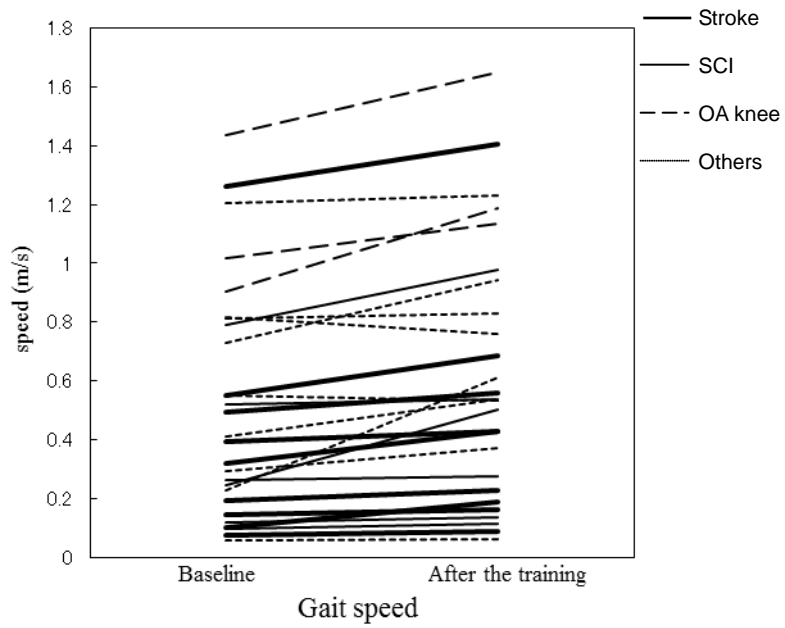


Figure 2.

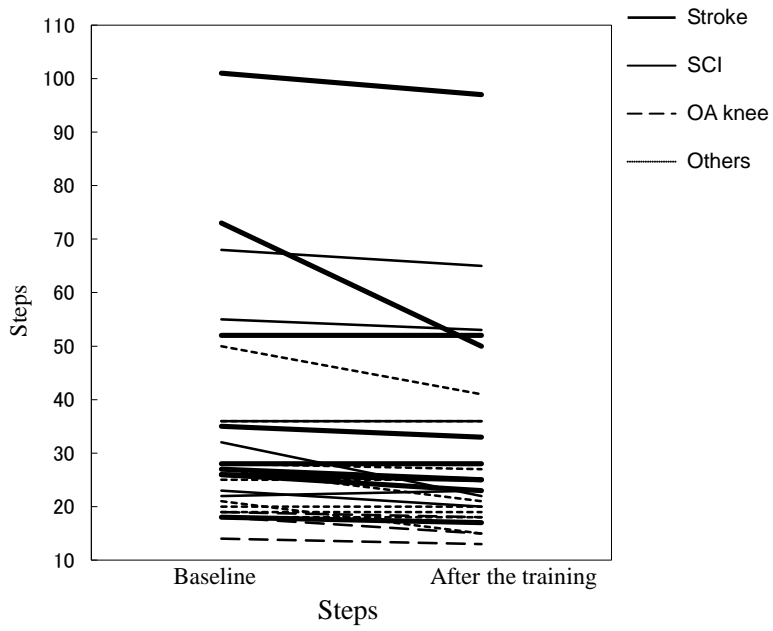


Figure 3.

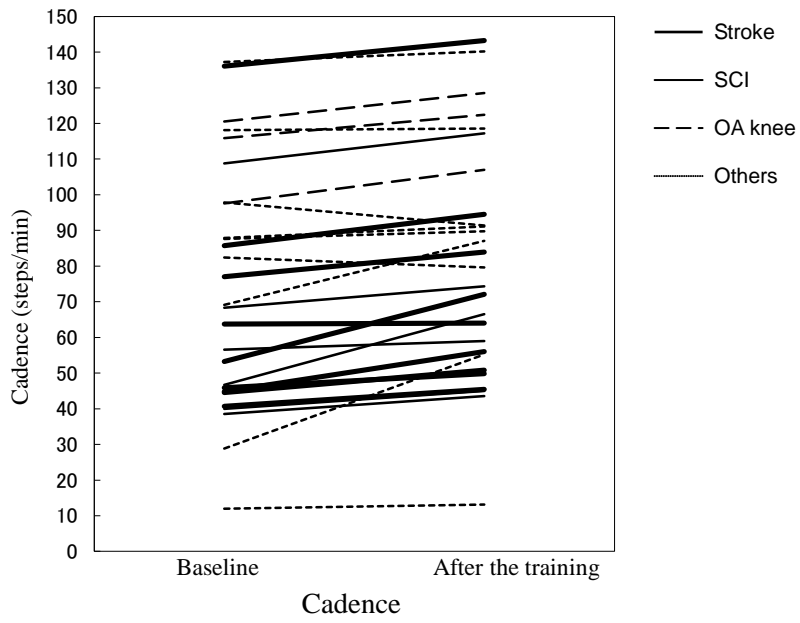


Figure 4.

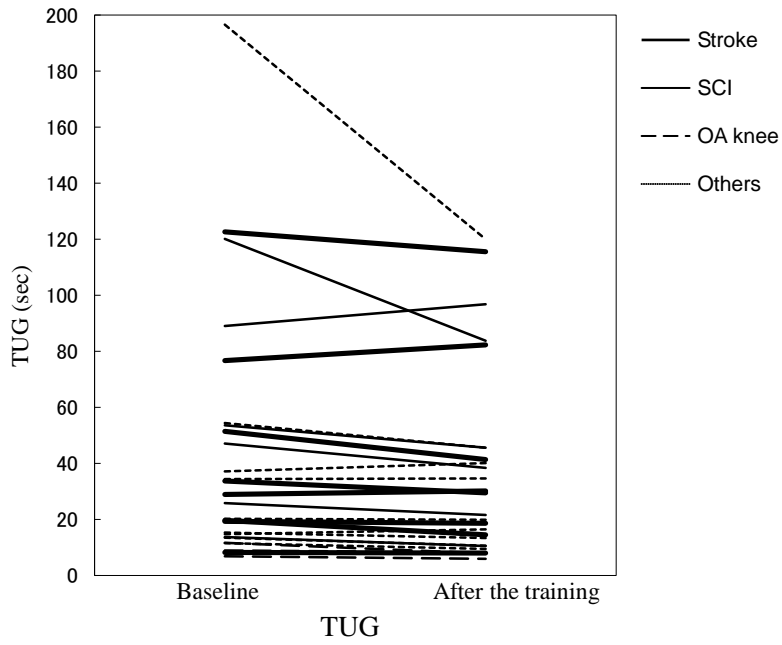


Figure 5.

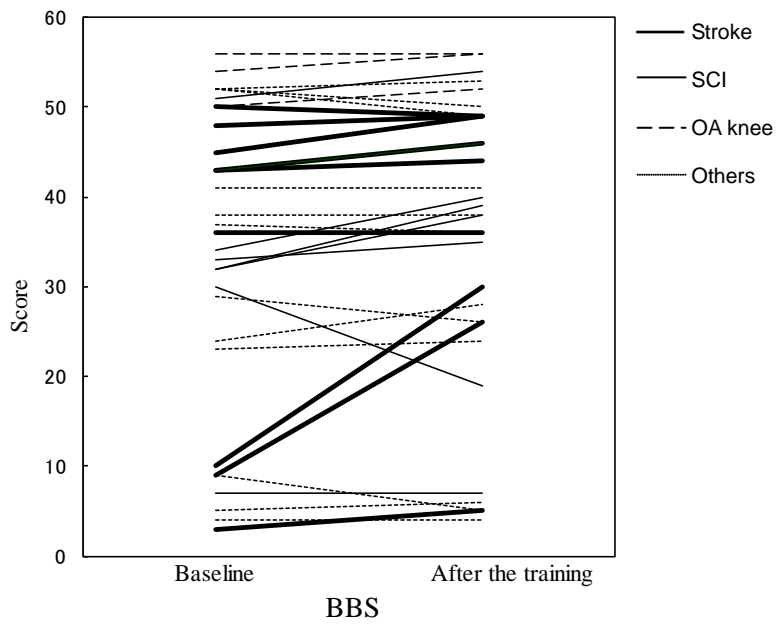


Figure 6.

Table 1: Clinical characteristics of patients

Case No.	Age (y)	Sex	Diagnosis	Paralysis type	Duration since disease	Ambulation	Assistive device	Orthosis	Training	Duration of training (wk)	Adverse events
1	69	M	stroke (cerebral infarcts)	paraplegia	15y	independently	t-cane	AFO	dropout (personal reason)	ND	Nothing
2	61	M	stroke (cerebral hemorrhage)	paraplegia	14y8mo	independently	t-cane	AFO	dropout (personal reason)	ND	Nothing
3	65	M	stroke (cerebral hemorrhage)	hemiplegia	2y2mo	supervision	quad-cane	AFO	complete	8	Nothing
4	37	F	stroke (cerebral hemorrhage)	quadriplegia	16y	independently	NA	AFO	complete	8	Nothing
5	72	M	stroke (cerebral infarcts)	hemiplegia	2y9mo	supervision	t-cane	AFO	complete	8	Nothing
6	54	M	stroke (cerebral hemorrhage)	hemiplegia	1y1mo	supervision	t-cane	NA	complete	8	Nothing
7	63	F	stroke (cerebral hemorrhage)	hemiplegia	1y6mo	independently	t-cane	AFO	complete	15	knee pain (patellar tendinitis)
8	52	M	stroke (cerebral hemorrhage)	ataxia	2y2mo	NA	NA	NA	complete	12	Nothing
9	74	M	stroke (cerebral infarcts)	hemiplegia	3y4mo	independently	t-cane	AFO	complete	9	Nothing
10	53	M	stroke (subarachnoid hemorrhage, cerebral infarcts)	hemiplegia	ND	supervision	pick up walker	KAFO	complete	9	Nothing
11	18	M	stroke (moyamoya disease)	hemiplegia	11y	independently	NA	AFO	complete	21	Nothing
12	64	M	stroke (cerebral hemorrhage)	hemiplegia	1y	supervision	t-cane	AFO	complete	8	Nothing
13	58	F	SCI (incomplete)	quadriplegia	3y3mo	supervision	lateral crutch	KAFO	complete	8	Nothing
14	69	M	SCI (incomplete)	quadriplegia	1y3mo	supervision	pick up walker	AFO	complete	8	Nothing
15	43	M	SCI (incomplete)	paraplegia	3y3mo	NA	NA	KAFO	dropout (medical reason)	ND	neuropathic pain after SCI
16	59	M	SCI (spina bifida)	paraplegia	6y4mo	supervision	t-cane	NA	complete	8	Nothing
17	31	M	SCI (complete)	paraplegia	3y	NA	NA	NA	complete	10	Nothing
18	64	F	SCI (incomplete)	quadriplegia	2y	independently	t-cane	AFO	complete	9	Nothing
19	54	M	SCI (central cervical cord injury)	quadriplegia	5y	supervision	t-cane	NA	complete	12	Nothing
20	47	M	SCI (spinal dural arteriovenous fistula)	paraplegia	1y1mo	independently	bilateral crutch	AFO	complete	8	Nothing
21	74	F	musculoskeletal disease (cervical spondylotic myelopathy)	quadriplegia	ND	independently	bilateral crutch	NA	dropout (medical reason)	ND	low back pain
22	81	F	musculoskeletal disease (OA Knee)	NA	ND	independently	NA	NA	complete	10	Nothing
23	44	M	musculoskeletal disease (OA Knee)	NA	ND	independently	NA	NA	complete	11	Nothing
24	74	M	musculoskeletal disease (OA Knee)	NA	ND	independently	NA	NA	complete	10	Nothing
25	62	M	Parkinson's disease	NA	8y	independently	NA	NA	complete	11	Nothing
26	72	F	Parkinson's disease	NA	7y8mo	independently	NA	NA	complete	9	Nothing
27	36	M	gonadotropin-dependent myopathy	paraplegia	19y	NA	NA	NA	complete	8	Nothing
28	52	F	limb-girdle muscular dystrophy	quadriplegia	24y	supervision	t-cane	NA	complete	9	Nothing
29	57	F	muscular dystrophy	NA	44y	independently	NA	NA	dropout (personal reason)	ND	Nothing
30	67	M	limb-girdle muscular dystrophy	NA	28y	independently	t-cane	NA	complete	8	Nothing
31	73	M	inclusion body myositis	NA	10y	independently	t-cane	NA	complete	10	knee pain
32	24	M	traumatic brain injury	quadriplegia	17y1mo	supervision	walker	NA	complete	8	Nothing
33	19	F	traumatic brain injury	quadriplegia	6y2mo	assistance	pick up walker	KAFO	complete	8	Nothing
34	29	F	traumatic brain injury	quadriplegia	10y7mo	assistance	pick up walker	KAFO	complete	9	Nothing
35	20	M	disuse syndrome, secondary to malignant lymphoma	NA	3y9mo	independently	t-cane	NA	dropout (personal reason)	ND	Nothing
36	31	F	cerebral palsy	quadriplegia	30y10mo	independently	lateral crutch	NA	complete	10	Nothing
37	55	M	sequelae of poliomyelitis	paraplegia	54y	independently	lateral crutch	NA	complete	19	Nothing
38	48	F	hypoxic ischemic encephalopathy	quadriplegia	2y	assistance	NA	NA	complete	12	Nothing

NOTE.

Abbreviations: AFO, ankle-foot orthosis; F, female; KAFO, knee-ankle-foot orthosis; M, male; NA, not applicable; ND, no data; OA, osteoarthritis

Table 2: Functional ambulation and balance ability at the baseline and after the 16-session HAL training

Outcome measurements	Baseline	After the training	Difference Mean (95% CI)	P-value	n
Speed (m/sec)	0.52 ± 0.40	0.61 ± 0.43	0.09 (0.05 to 0.14)	p < 0.001	27
10 MWT Steps	34.0 ± 20.4	31.0 ± 18.8	-3.0 (-4.9 to -1.0)	p < 0.001	27
Cadence (steps/min)	74.3 ± 34.1	81.1 ± 32.9	6.8 (4.0 to 9.6)	p < 0.001	27
TUG (sec)	43.7 ± 45.0	37.3 ± 34.1	-6.4 (-13.0 to 0.2)	0.057	26
BBS	33.6 ± 16.9	35.5 ± 16.3	1.9 (-0.1 to 3.9)	0.059	32

NOTE. Values expressed as mean ± SD.

Abbreviation: TUG, Timed-Up and Go. BBS, Berg Balance Scale.