

University of Rhode Island

DigitalCommons@URI

Open Access Master's Theses

1987

Gas Chromatographic/Infared Spectral Analysis Using a Computer Integrated Circular Variable Filter Spectrometer

Mark Ahmadjian

University of Rhode Island

Follow this and additional works at: <https://digitalcommons.uri.edu/theses>

Recommended Citation

Ahmadjian, Mark, "Gas Chromatographic/Infared Spectral Analysis Using a Computer Integrated Circular Variable Filter Spectrometer" (1987). *Open Access Master's Theses*. Paper 943.
<https://digitalcommons.uri.edu/theses/943>

This Thesis is brought to you for free and open access by DigitalCommons@URI. It has been accepted for inclusion in Open Access Master's Theses by an authorized administrator of DigitalCommons@URI. For more information, please contact digitalcommons@etal.uri.edu.

GAS CHROMATOGRAPHIC/INFRARED SPECTRAL ANALYSIS USING A COMPUTER
INTERFACED CIRCULAR VARIABLE FILTER SPECTROMETER

BY

MARK AHMADJIAN

A THESIS SUBMITTED IN PARTIAL FULFILLMENT OF THE
REQUIREMENTS FOR THE DEGREE OF

MASTER OF SCIENCE

IN

CHEMISTRY

UNIVERSITY OF RHODE ISLAND

1987

MASTER OF SCIENCE THESIS






OF

MARK AHMADJIAN

APPROVED:

THESIS COMMITTEE

MAJOR PROFESSOR

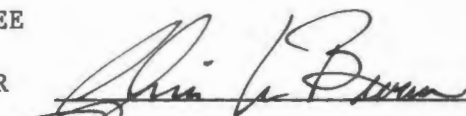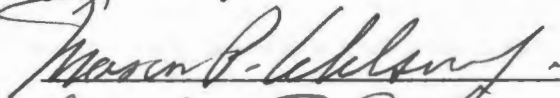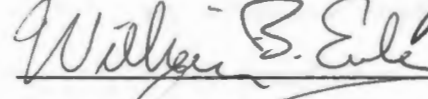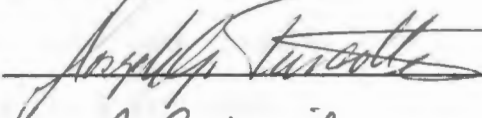DEAN OF THE GRADUATE SCHOOL

ABSTRACT

Gas chromatography/infrared spectroscopy (GC/IR) is a powerful analytical technique for "on-the-fly" analysis of multi-component mixtures. The gas chromatograph separates the mixture into single components which are then analyzed by the spectrometer as they elute off the gc column and into the spectrometer. Standard instrumentation usually involves a gas chromatograph, a fourier transform infrared spectrometer (FT-IR), and a computer system. The end products are a reconstructed gas chromatogram and a number of infrared spectra for each of the mixture components. These pieces of information often allow the identification of an unknown mixture. A GC/FT-IR/Computer System, however, is an elaborate and expensive laboratory instrument which requires extensive care, maintenance, and a well trained operator. It is not designed for field use and is seldom deployed outside of a laboratory. The objective of this research therefore was to develop a simple and inexpensive GC/IR/Computer System which would be rugged, easy to use, and readily adaptable for field use.

To do this, we interfaced a fast scanning circular variable filter (CVF) infrared spectrometer to a mini-computer and then coupled the spectrometer/computer to a gas chromatograph. The system was configured such that the same end products generated by a GC/FT-IR/Computer System, i.e. reconstructed gas chromatograms and infrared spectra, are also generated by this GC/CVF-IR/ Computer System. A five component mixture of hazardous class type compounds with closely related boiling points was analyzed to demonstrate the feasibility of this new system.

ACKNOWLEDGEMENTS

I would like to thank my major professor, Dr. Chris W. Brown, for his guidance, teaching, and most of all his friendship throughout the years.

My thanks to the U.S. Air Force for providing me with the opportunity to finish this degree and special thanks to Dr. Randall E. Murphy, Director, Infrared Technology Division, AF Geophysics Laboratory, for his support and encouragement. /

Mark Ahmadjian

PREFACE

This thesis was prepared in accordance with the standard thesis plan.

TABLE OF CONTENTS

	<u>Page</u>
ABSTRACT.....	ii
ACKNOWLEDGEMENT.....	iii
PREFACE.....	iv
LIST OF FIGURES.....	vi
INTRODUCTION.....	1
HISTORY.....	2
EXPERIMENTAL SECTION.....	3
RESULTS AND DISCUSSION.....	9
CONCLUSION.....	25
REFERENCES.....	27
BIBLIOGRAPHY.....	28
APPENDIX.....	30

LIST OF FIGURES

<u>Figure</u>		<u>Page</u>
1.	Block diagram of the GC/IR computer system.....	4
2.	Beckman IR-102 optical diagram.....	5
3.	Sequence of events for the GC/IR computer analysis..	7
4.	TWIRP (Total Wavelength IR Plot) A reconstructed gas chromatogram.....	12
5.	Spectral plots corresponding to the TWIRP peak in figure 4.....	12
6.	TWIRP for a five component mixture consisting of (1) methanol (2) acetone, (3) chloroform, (4) carbon tetrachloride, and (5) benzene in a carbon disulfide solution (CS ₂ vented to atmosphere).....	13
7.	Infrared spectra for the five component mixture corresponding to the TWIRP peak in figure 6. Scan #1 is the instrument background while scans 14, 22, 41, 47, and 55 correspond to the TWIRP peaks #1, #2, #3, #4, and #5 respectively.....	16
8.	Infrared spectrum for scan #1 showing instrument background. Absorptions are due to ambient CO ₂ and H ₂ O.....	17
9.	Infrared spectrum for scan #14 which corresponds to the first TWIRP peak of figure 6 and is due to methanol.....	18
10.	Infrared spectrum for scan #22 which corresponds to the second TWIRP peak of figure 6 and is due to acetone.....	19
11.	Infrared spectrum for scan #41 which corresponds to the third TWIRP peak of figure 6 and is due to chloroform.....	20
12.	Infrared spectrum for scan #44 which corresponds to the fourth TWIRP peak of figure 6 and is due to carbon tetrachloride.....	21

LIST OF FIGURES

<u>Figure</u>		<u>Page</u>
13.	Infrared spectrum for scan #55 which corresponds to the fifth TWIRP peak of figure 6 and is due to benzene.....	22
14.	Comparision of a gas chromatogram (flame ionization detector) to the TWIRP reconstructed gas chromatogram of a five component mixture of methanol, acetone, chloroform, carbon tetrachloride, and benzene.....	24
15.	Computer processing of Beckman IR-102 spectra demonstrating signal averaging and smoothing techniques..	26

INTRODUCTION

Since the 1950s infrared spectroscopy has been used to identify components separated by gas chromatography. Early efforts¹⁻¹⁰ involved grating or prism type spectrometers which, however, lacked the wavelength scanning speed to adequately measure chemical compounds eluting from gas chromatographs in real time, i.e. on-the-fly. It was usually necessary to trap the components and hold them in a trap or matrix until an infrared spectrum of the sample was obtained. This was time consuming and did not provide for generation of reconstructed gas chromatograms. Not until Cooley and Tukey (11) developed their algorithm for fast processing of interferograms and the ready availability of laboratory computers did actual GC/IR "on-the-fly" analysis become practical. For the first time a gas chromatograph could be coupled directly to a Fourier Transform - Infrared instrument (FT-IR) for analysis of multi-component mixtures with results plotted out by computer in the form of reconstructed gas chromatograms with corresponding spectra for the individual gas chromatographic peaks. The ability of the FT-IR system to rapidly acquire interferograms over a broad spectral region revolutionized the analysis of mixtures by gas chromatography. Today there are a number of GC/FT-IR instruments available (Nicolet, Perkin-Elmer, Digilab, IBM, and others); however, they are expensive, complicated (moving mirrors, air bearings, lasers, extensive data processing requirements), and are found only in a limited number of laboratories. In addition, they are not considered to be field type instruments. The objective of this research was to

develop a simple, inexpensive, and transportable GC/IR system using a circular variable filter (CVF) infrared spectrometer controlled by a small mini-computer. This combination could provide facilities with an inexpensive and easy to use survey type instrument for fast characterization of multi-component mixtures.

To do this we used a Beckman IR-102 fast scanning CVF infrared spectrometer to measure gas chromatographic effluents by passing the volatile fractions through a gold plated light pipe in the optical path of the spectrometer while it scanned the 2.5 - 14.5 μm region.

History

Analysis of gas chromatographic fraction with an IR-102 is not new. Beckman Instruments developed this instrument in 1965 for this purpose¹². At that time though there were few laboratory computers in routine use and spectra were recorded by a fast oscillographic recorder which necessitated manually and laboriously examining many feet of heat sensitive chart paper to pick out the infrared spectrum of interest. Needless to say this technique was not accepted by the analytical community and only a few IR-102's were ever sold. There is no evidence at Beckman Instruments or in the literature that this type of fast scanning CVF spectrometer was ever interfaced to a computer to perform the same tasks as a GC/FT-IR (however, the U.S. Air Force did use computerized CVF spectrometers for chemical analysis of the earth's upper atmosphere from sounding rockets and high altitude balloons).

This paper will show how a GC/IR-102 can be interfaced to a Data General Corporation NOVA 3/12 mini-computer and developed the software to perform GC/IR analysis on-the-fly. The system is demonstrated on a one component solution and then on a five component mixture of closely boiling compounds; acetone b.p. 56 C, chloroform b.p. 61 C, methanol b.p. 64 C, carbon tetrachloride b.p. 76 C, and benzene b.p. 80 C. The GC/IR not only provides a reconstructed gas chromatogram comparable to standard flame ionization detection (FID) gas chromatograms but, also generates infrared spectra for each of the chromatographic peaks.

EXPERIMENTAL SECTION

Instrumentation

Figure 1 is a block diagram of the GC/IR computer system. A HP-7100 gas chromatograph with a Supelco 10% SP-2250, 100/120 Supelcoport, 1/8"x10' column was directly coupled to the Beckman IR-102 gold plated 10 cm long light pipe. The interface line was heated to 150° C and was thermostatically controlled. 100% of the g.c. effluent was allowed to pass through the sample cell. The flow rate was varied from 10-15 cc/minute. The optical path and components of the IR-102 were used as designed except it was no longer necessary to use the trap and hold option valves for the sample cell. Figure 2 shows the optical layout. A Nernst glower was used as the IR source and chopped at 120 Hz. The beam is directed through the heated sample cell. A torroidal mirror collected the radiation and reflected it onto a flat mirror and through the CVF wheel,


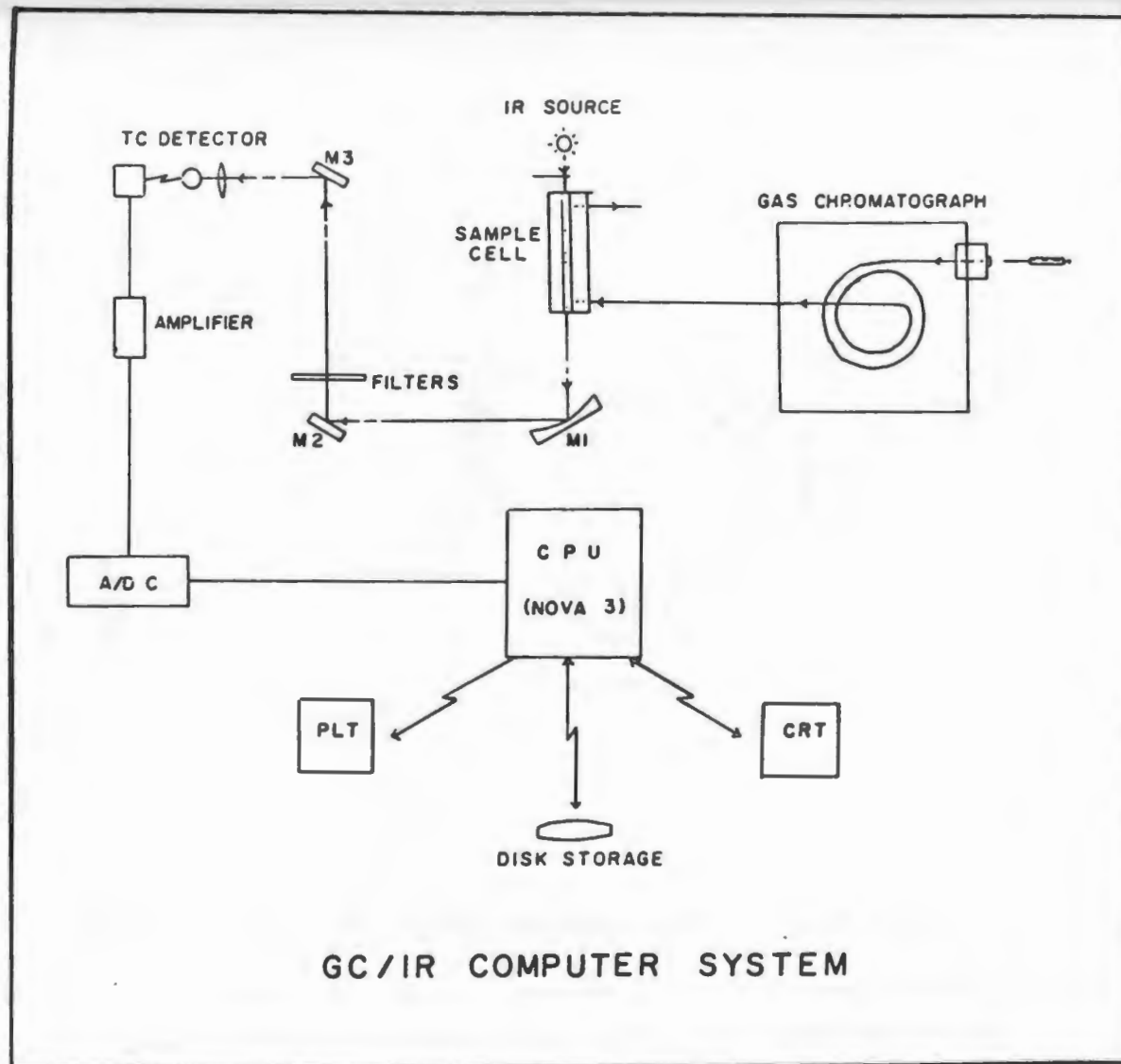


Figure 1. Block diagram of the GC/IR computer system

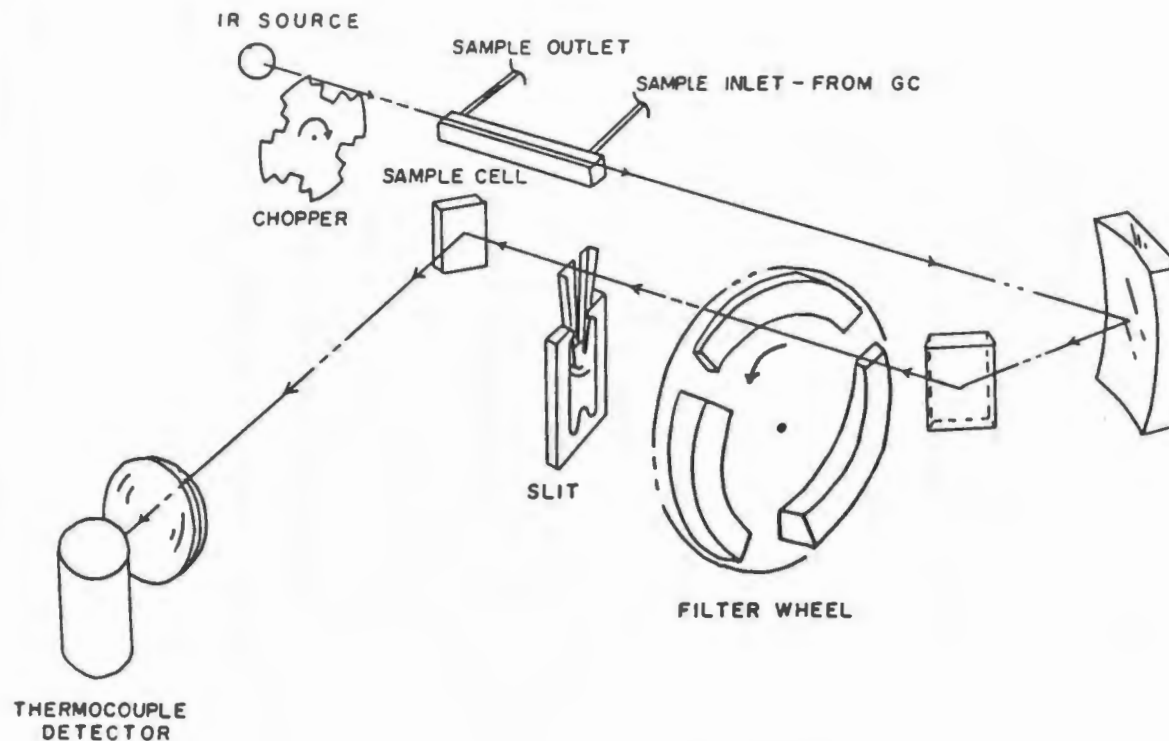

OPTICAL DIAG. BECKMAN IR-102 FAST SCANNING IR

Figure 2. Beckman IR-102 optical diagram (redrawn from Beckman user's manual)

a slit, and onto a thermocouple detector. The filter wheel consists of three wedge shaped segments covering the infrared wavelength range of 2.5 - 4.5 μm , 4.4 - 8.0 μm , and 7.9 - 14.5 μm . The analog output signal from the detector is passed through a signal conditioning channel where it is inverted and amplified to 0 to +10 VDC. The signal is then sampled by the DGC NOVA 3/12 16 bit analog digital converter. The sampling rate is determined by an internal oscillator, 1000 Hz clock, under software control. The digitized signal is then stored in contiguous data blocks on a 5 Megabit cartridge disk subsequent to any additional processing. The data can then be processed and displayed in a variety of ways. Usually a Total Wavelength IR Plot (TWIRP), i.e. the reconstructed gas chromatogram, is first generated to show how many components were in the mixture. The individual spectra for each of the reconstructed g.c. peaks can then be displayed as a function of scan number. The spectra can also be signal averaged, smoothed, and subtracted. Figure 3 shows a typical sequence of events for this analysis.

Samples

Mixtures were prepared by mixing equal volumes of reagent grade acetone, chloroform, methanol, carbon tetrachloride, and benzene in carbon disulfide.

Software

Eight computer programs were written for this system. By using

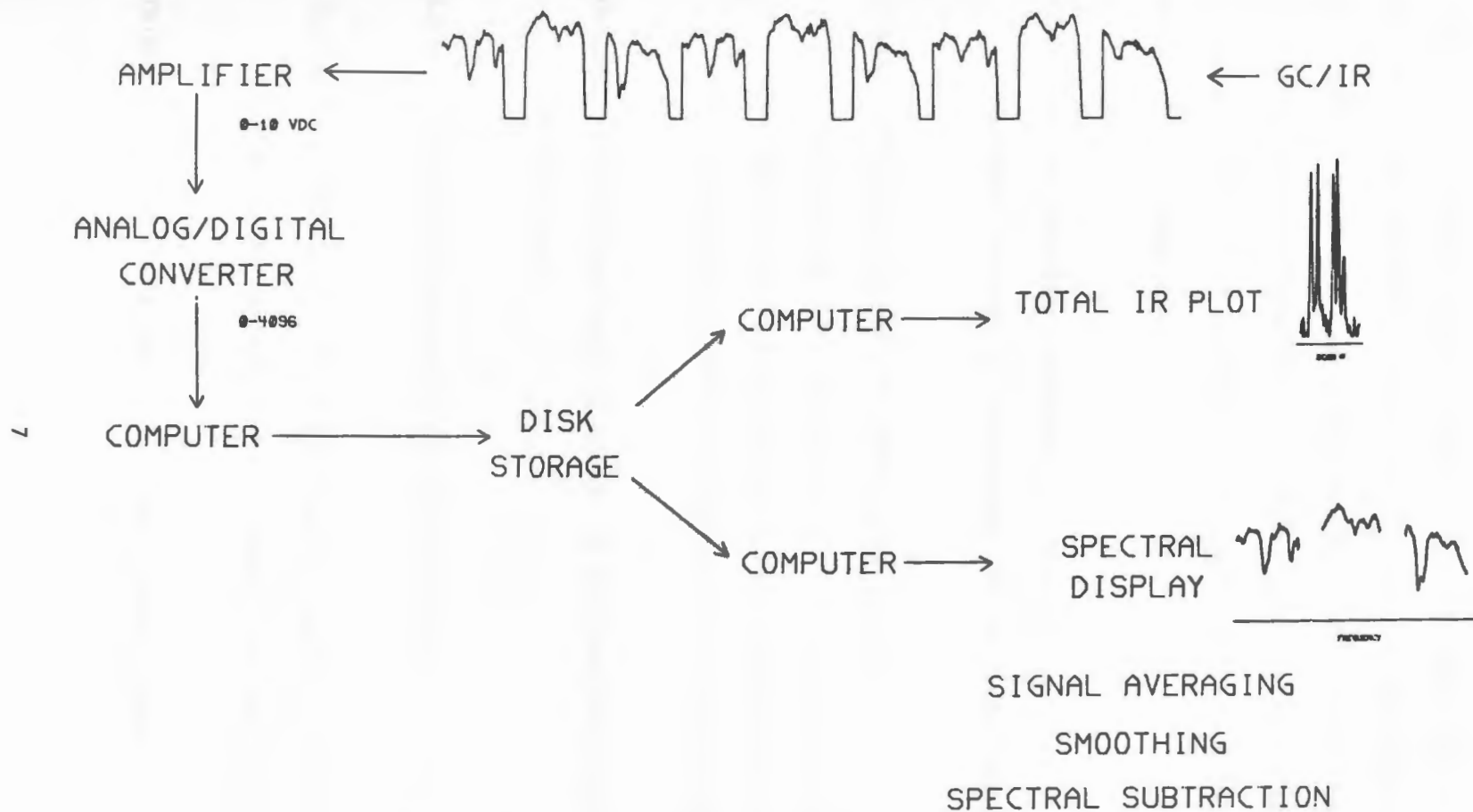


Figure 3. Sequence of events for the GC/IR computer analysis

individual programs we were able to develop, test, and debug the hardware/software one phase at a time. The individual programs are brought into memory by a main controller program. When not in use, they reside in one of the user directories on disk. The following is a list of the programs with a brief description of each. Copies of the programs are attached as an appendix.

GCIR: This is the main controller program, the executive, and allows the other programs to be brought into use via swap commands.

GCIRAD: This program and its subroutines samples the IR-102 analog detector signal through a multiplexer, performs an analog to digital conversion via the analog digital converter (ADC), and stores the digitized data into an optimally sized contiguous data file on disk.

GCIRSUB: A subroutine used by GCIRAD to initiate data acquisition and store data on disk.

ADGCIR: An assembly language subroutine to sample the ADC.

GCIRCAL: Calibrates all the individually acquired IR spectra using the 0 VDC signal between filter elements one and three.

TWIRP: Plots out the reconstructed gas chromatograms.

GCIRSPEC: Plots out individual spectra corresponding to peaks from the reconstructed gas chromatograms.

GCIRSUM: Signal averages any number of IR spectra.

RESULTS AND DISCUSSION

We want to show that this GC/IR system using a computer interfaced IR-102 infrared circular variable filter spectrometer is able to analyze a multi-component mixture in a manner similar to that of a sophisticated GC/FT-IR/computer system. As in any research endeavor, we started off with a simple feasibility demonstration which involved the measurement of a one component mixture. In this case we used acetone in carbon disulfide. Carbon disulfide was chosen as the solvent because our laboratory was using it on a routine basis for IR analysis of hydrocarbon solutions and because carbon disulfide has a large IR "window" and would therefore present minimal spectral interferences for low boiling compounds. To prevent gross CS₂ contamination of the gold plated IR sample cell, the solvent peak was vented via a bypass valve to an exhaust hood. Once the CS₂ had eluted from the g.c. column (time as determined from previous flame ionization detection measurements) the g.c. column effluent was rerouted to the Beckman IR-102 sample cell and the analog signal from the IR-102 detector amplifier was digitized via program GCIRAD. After sixty full wavelength scans, data acquisition was terminated and the IR data was transformed from contiguous data

data blocks to random data files, calibrated, converted to absorbance units, and plotted out in the form of a reconstructed gas chromatogram (figure 4). We call this a Total Wavelength IR Plot (TWIRP). The x-axis is plotted as a function scan number (same as time), i.e. each cycle of the circular variable filter wheel from 2.5 to 14.5 μm is one scan. In this case sixty scans are plotted. The y-axis represents the integrated signal of each single scan, i.e. the summation of each data point over the entire spectral region of one scan; $\sum A_i/\text{scan}$. For example, if each scan covered 2.5 to 3.0 μm and a data point was taken at every 0.1 μm , the result would be:

$$\text{scan \#1} = A_{2.5} + A_{2.6} + A_{2.7} + A_{2.8} + A_{2.9} + A_{3.0}$$

$$\text{scan \#2} = A_{2.5} + A_{2.6} + A_{2.7} + A_{2.8} + A_{2.9} + A_{3.0}$$

$$\text{scan \#3} = A_{2.5} + A_{2.6} + A_{2.7} + A_{2.8} + A_{2.9} + A_{3.0}$$

.
.
.
.
.

In this actual case the following would be represented:

$$\text{scan \#xxx} = A_{2.5} + A_{2.6} + A_{2.7} + \dots\dots\dots A_{14.5}$$

or in other words:

$$14.5$$

$$\text{scan \#xxx} = \sum A_i$$

$$i=2.5$$

The $\{A_i\}$'s are then normalized and plotted as a function of increasing scan number and we end up with a reconstructed gas chromatogram or a TWIRP. From this TWIRP, the number of compounds in the injected mixture can be easily seen. Of course the compounds need to have IR active groups to cause a signal on the IR detector. In this initial case, after sixty scans, there is one TWIRP peak starting at scan #25. Figure 5 shows the IR spectra associated with the scan numbers 22,24,25,26, and 28 plotted as a function of IR intensity vs wavenumber. Scan #22 is instrument background and since this is a single beam spectrometer the spectrum is due to ambient CO_2 and H_2O . At scan #24 however a distinctly different spectrum is seen and its intensity continues to increase at scan #25. At scan #28 the sample has completely eluted from the g.c. column thru the cell and all that is left is ambient background. These two figures would give two significant pieces of information for identification of an unknown sample. From having previously characterized the SP-2250 column, knowing the g.c. oven temperature is at 150°C , and seeing the sample elute after 24 scans (144seconds), this would identify it as a low boiling compound. Then by examining the IR spectrum and looking up the wavelengths of the absorbance bands a preliminary identification of the eluting compound would be made. Further chemical analysis would confirm that finding. Similar runs were performed for acetone, methanol, carbon tetrachloride, and benzene.

Having thus demonstrated the GC/IR system on one component solutions

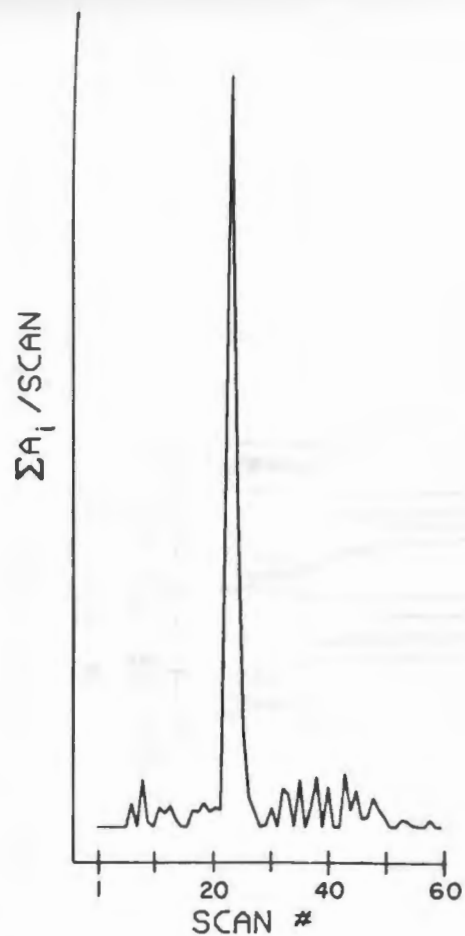


Figure 4. TWIRP (Total Wavelength IR Plot)
A reconstructed gas chromatogram

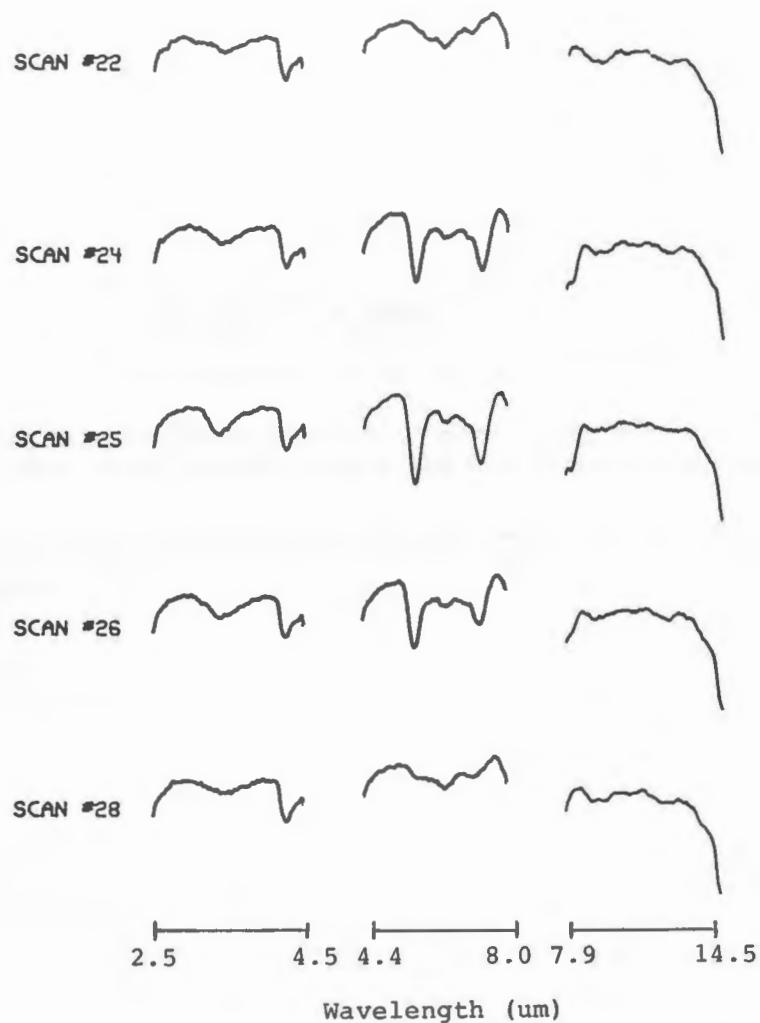


Figure 5. Spectral plots corresponding to the TWIRP peak
in figure 4

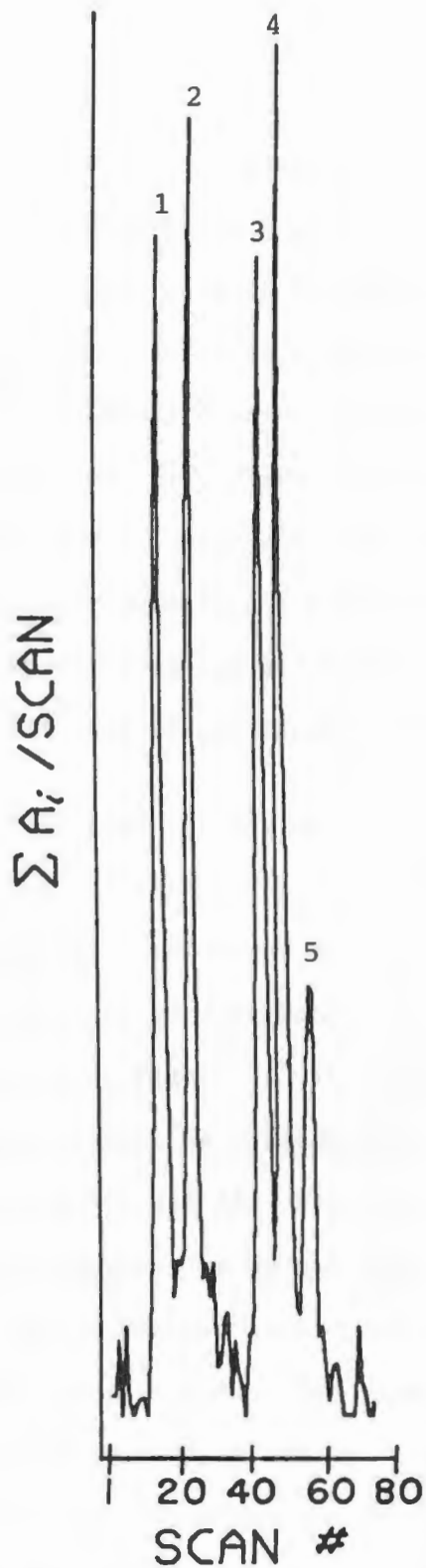


Figure 6. TWIRP for a five component mixture consisting of (1) methanol, (2) acetone, (3) chloroform, (4) carbon tetrachloride, and (5) benzene in a carbon disulfide solution (CS_2 vented to atmosphere).

the next step was the analysis of a multi-component mixture. Five compounds closely related in boiling points were chosen to provide a test as to whether we could separate and obtain spectra on-the-fly of rapidly eluting compounds. The five components, acetone b.p. 56°C, chloroform b.p. 61°C, methanol b.p. 64°C, carbon tetrachloride b.p. 76°C, and benzene b.p. 80°C were mixed in equal volumes into a carbon disulfide solution. This solution not only provided a stringent test as far as separation and identification but was also typical of compounds found in hazardous waste sites and which must be analyzed in an actual facility prior to the start of any cleanup operations.

The same procedure used for the one component solution was used for this five component analysis. The CS₂ was vented and the scan and acquisition initiated. Spectral scans were recorded, stored on disk, calibrated, converted to absorbance units, and the TWIRP displayed on the plotter as shown in figure 6. It is immediately evident that we have a five component mixture as evidenced by the reconstructed g.c. peaks at scan #14,22,41,47, and 55. Even though we had a known solution of equal amounts of compounds, we do not necessarily end up with equal intensity peaks. This is because each component has different IR active groups and absorbance coefficients. For example, CS₂ has very little IR absorbance in the 2.5 to 14.5 um region compared to acetone which has strong IR absorbance. Therefore the IR-102 detector would not "see" very much CS₂ relative to acetone and the resultant integrated IR signal would be very small. The same is true for all compounds and since the TWIRP is normalized to the strongest integrated signal, the

reconstructed g.c. peaks are all relative in intensity.

Again as done previously the spectra associated with each reconstructed g.c. peak are plotted out as a function of % transmission vs wavelength (figure 7). Scan #1 is the instrument background and the subsequent scans correspond to the five g.c. peaks. Figures 8, 9, 10, 11, 12, and 13 are the spectra displayed for more detailed analysis.

Scan #1: Instrument background. Absorptions correspond to atmospheric H_2O and CO_2 .

Scan #14: Represents TWIRP peak #1 which is due to methanol.

Scan #22: Represents TWIRP peak #2 which is due to acetone.

Scan #41: Represents TWIRP peak #3 which is due to chloroform.

Scan #47: Represents TWIRP peak #4 which is due to carbon tetrachloride.

Scan #55: Represents TWIRP peak #5 which is due to benzene.

For a real world application we would not know apriori what each IR spectrum is due to, but there is significant information contained in the GC/IR TWIRP and spectra to allow an educated selection of possible candidates. In scan #14 there is strong absorption at 3.0 um which is very

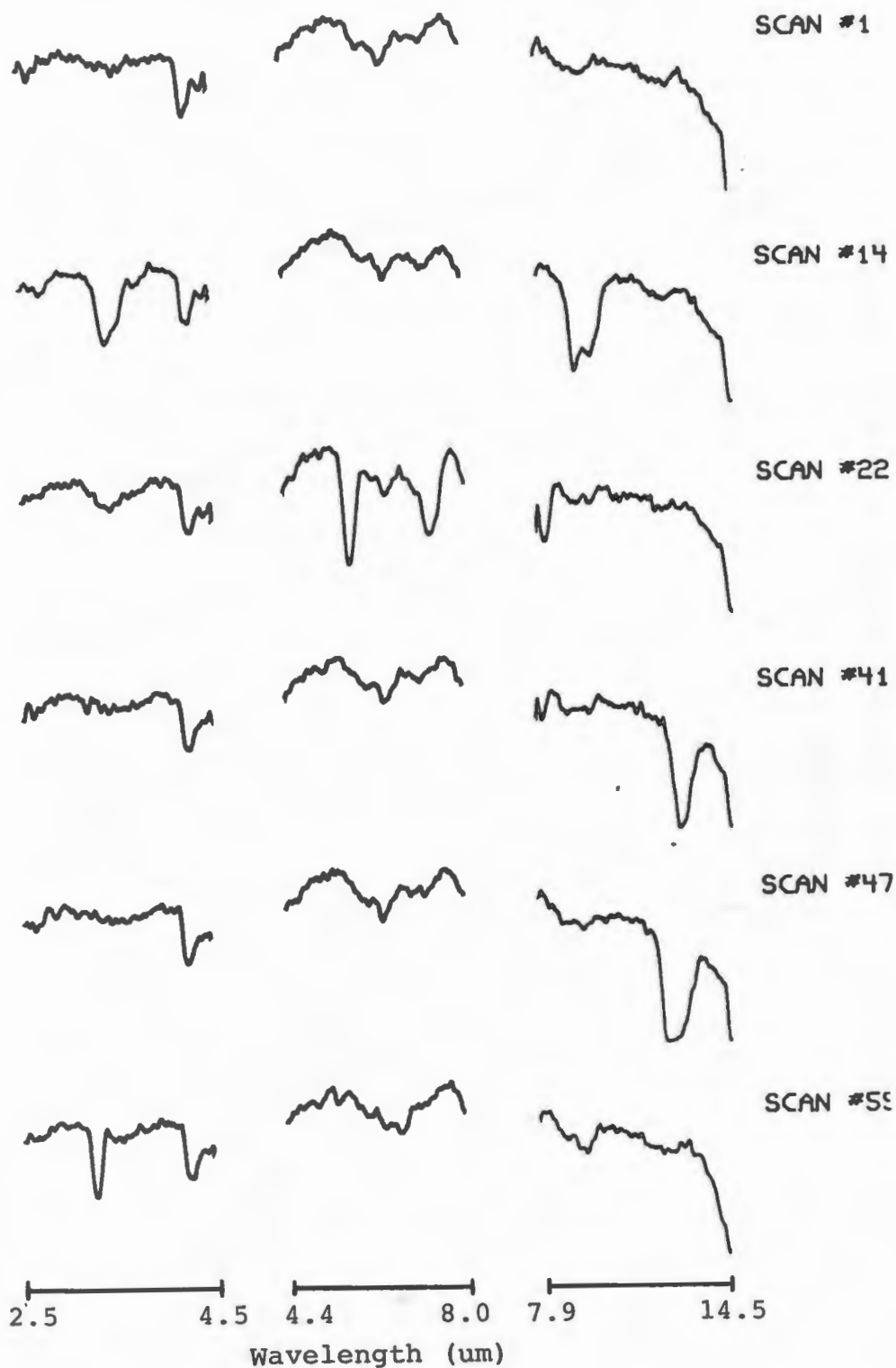


Figure 7. Infrared spectra for the five component mixture corresponding to the TWIRP in figure 6. Scan #1 is the instrument background while scans 14, 22, 41, 47, and 55 correspond to the TWIRP peaks #1, #2, #3, #4, and #5 respectively.

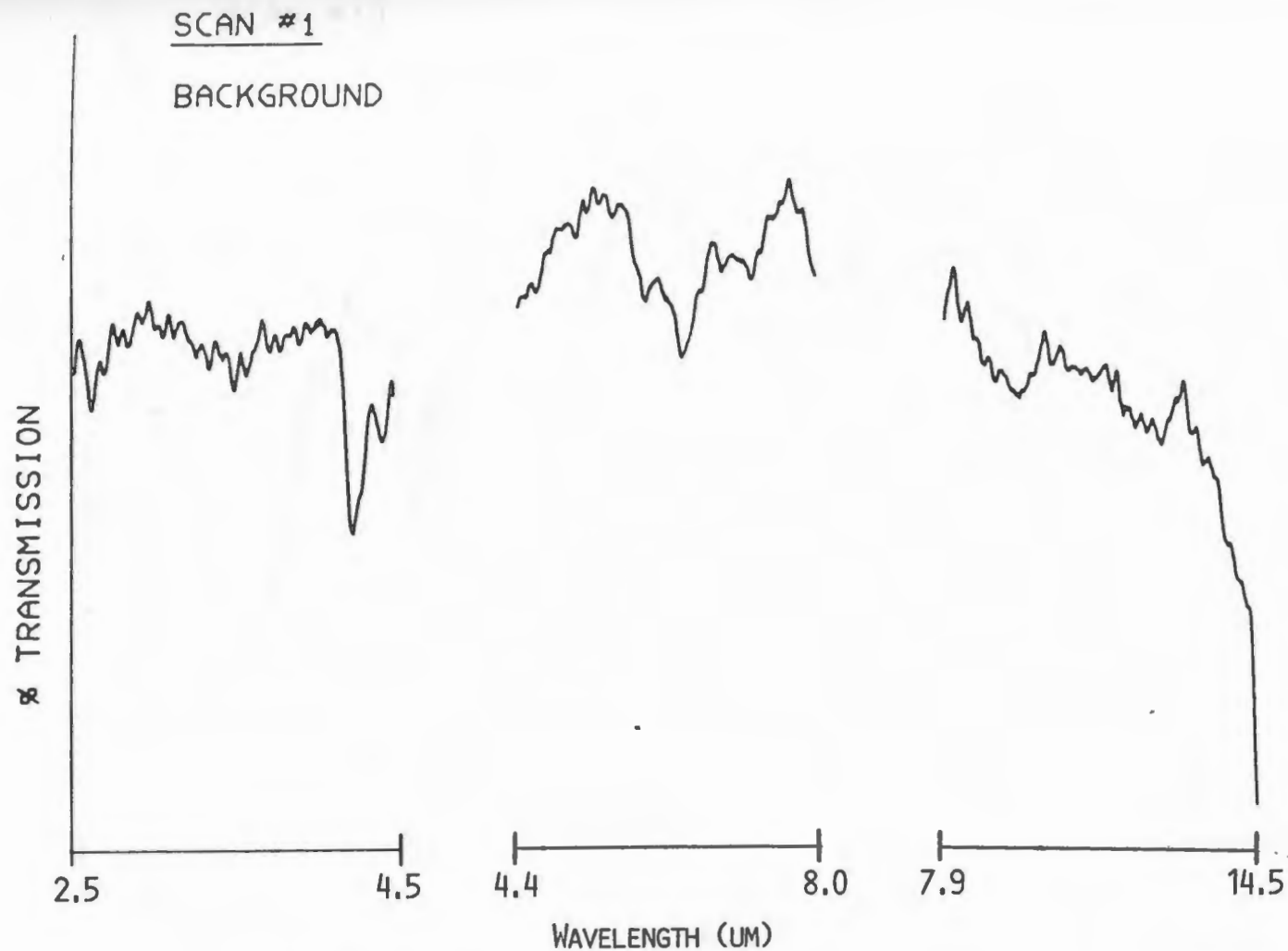


Figure 8. Infrared spectrum for scan #1 showing instrument background. Absorptions are due to ambient CO₂ and H₂O.

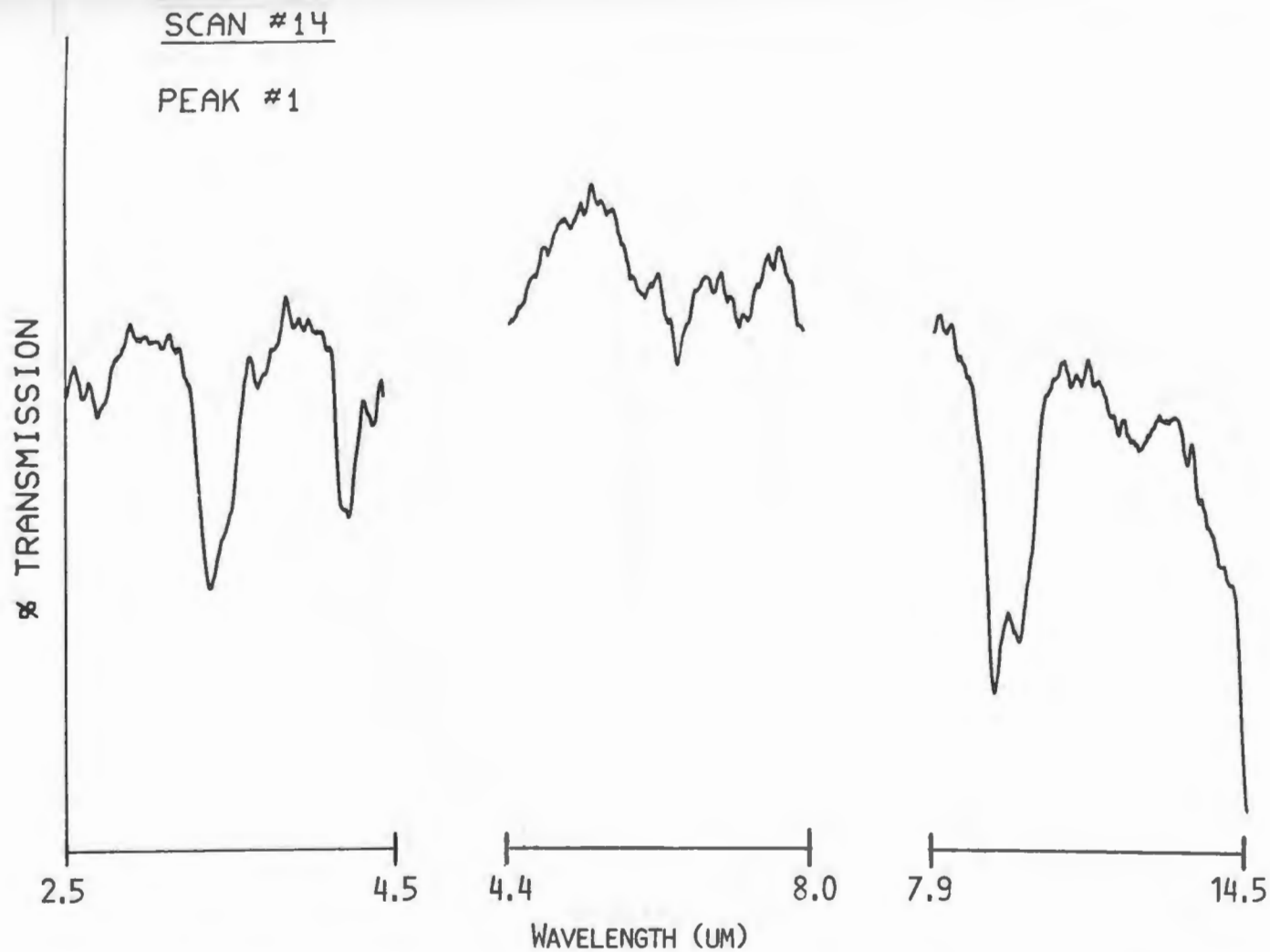


Figure 9. Infrared spectrum for scan #14 which corresponds to the first TWIRP peak of figure 6 and is due to methanol.

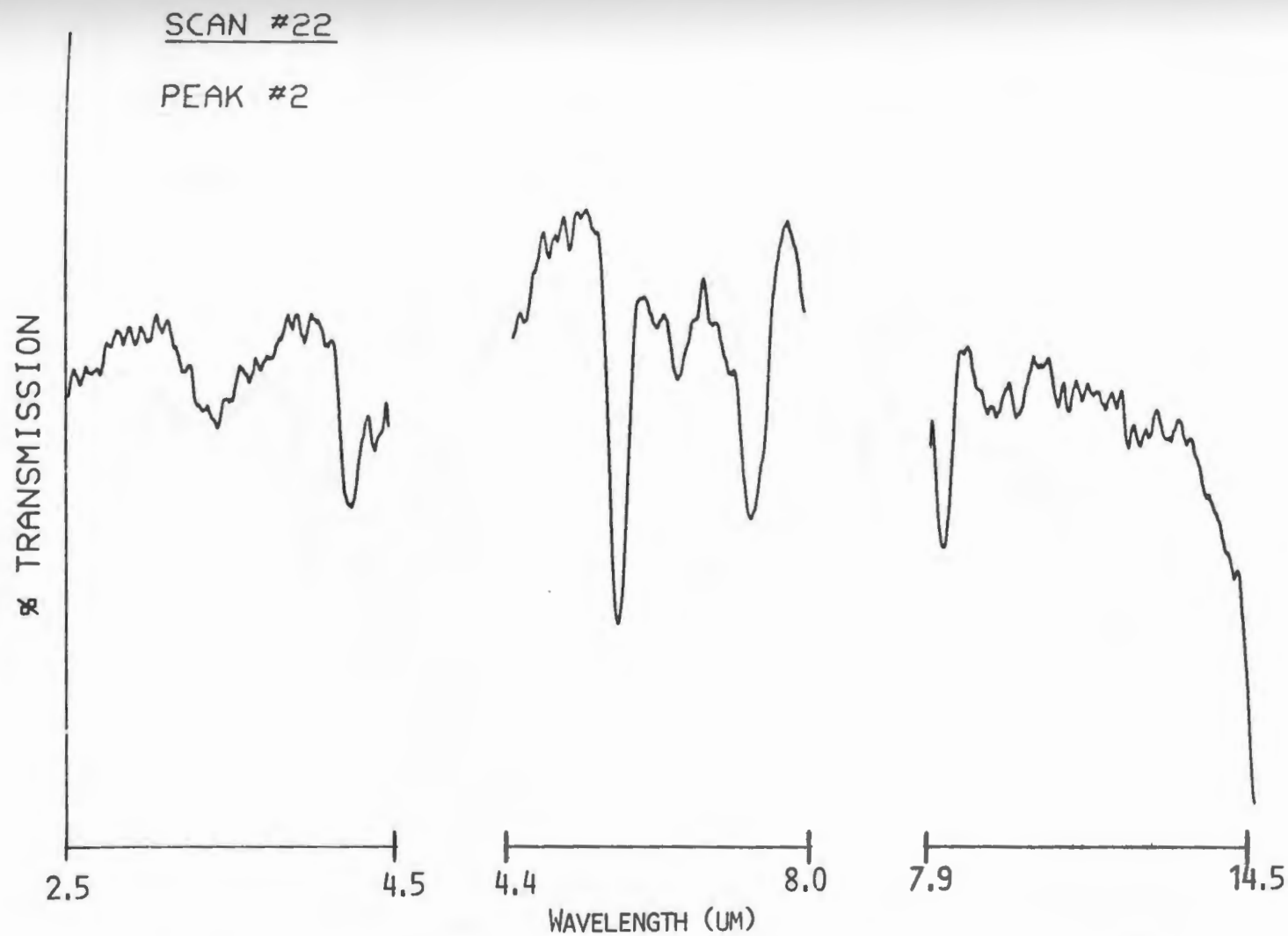


Figure 10. Infrared spectrum for scan #22 which corresponds to the second TWIRP peak of figure 6 and is due to acetone.

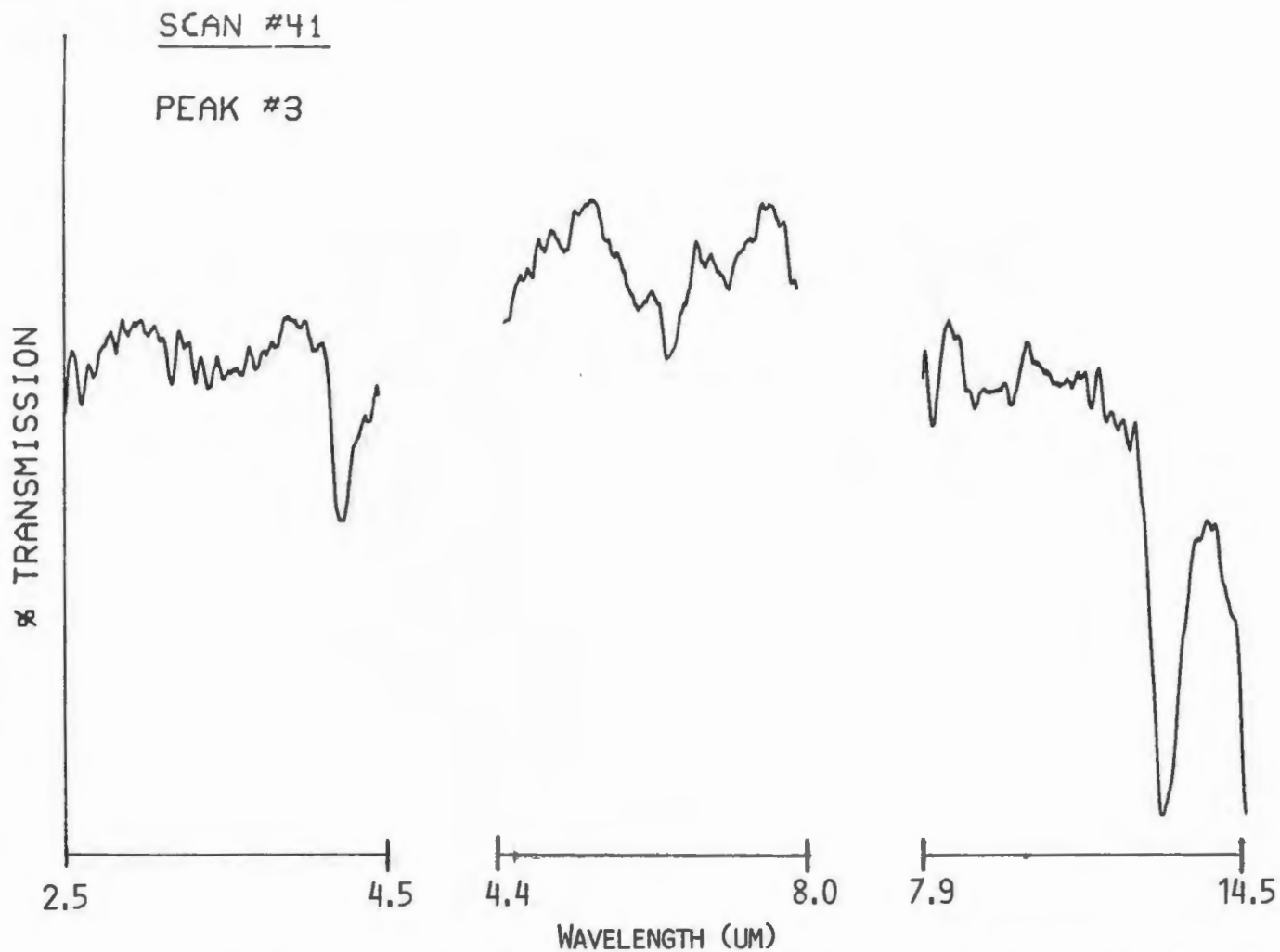


Figure 11. Infrared spectrum for scan #41 which corresponds to the third TWIRP peak of figure 6 and is due to chloroform.

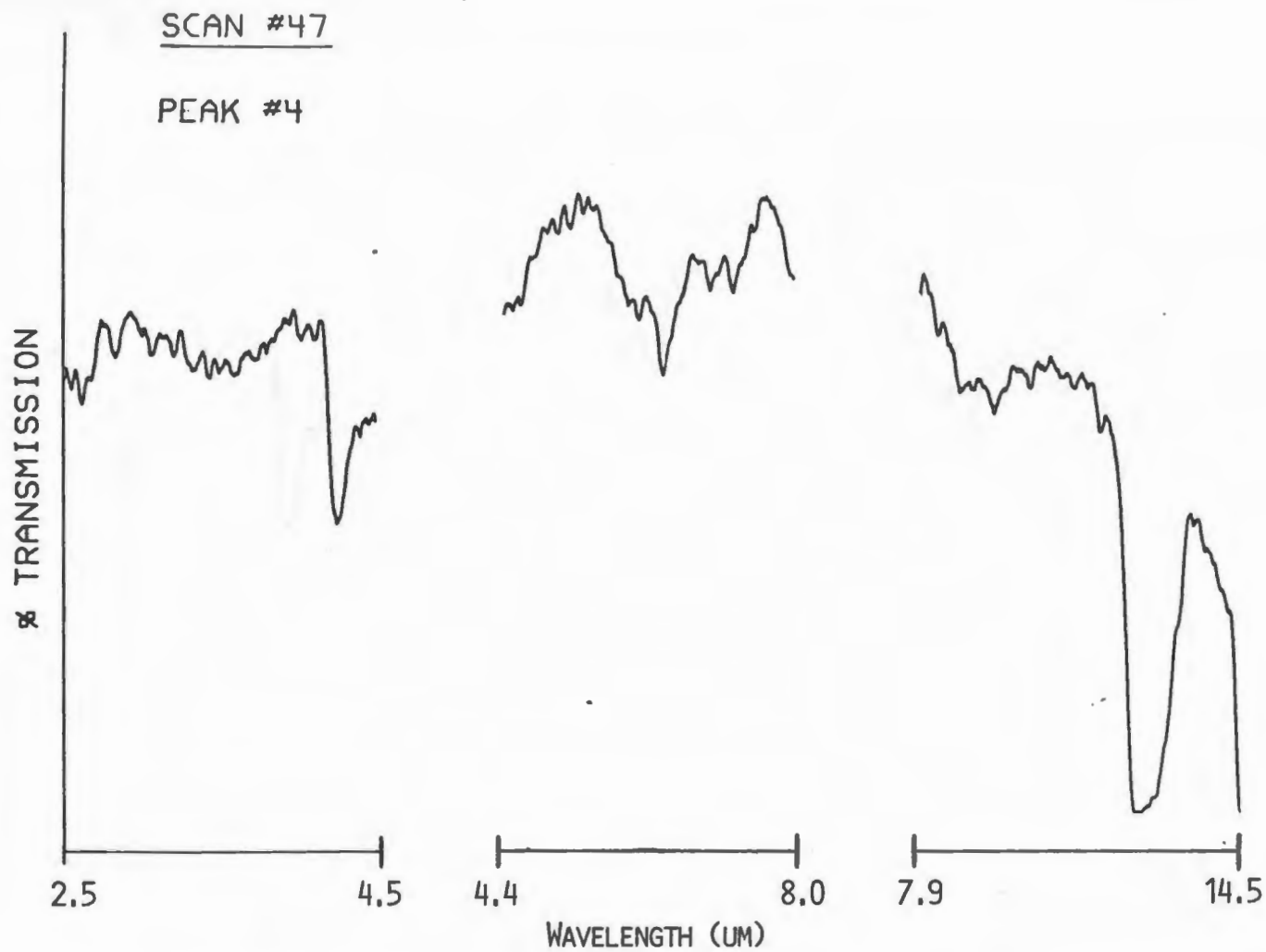


Figure 12. Infrared spectrum for scan #44 which corresponds to the fourth TWIRP peak of figure 6 and is due to carbon tetrachloride.

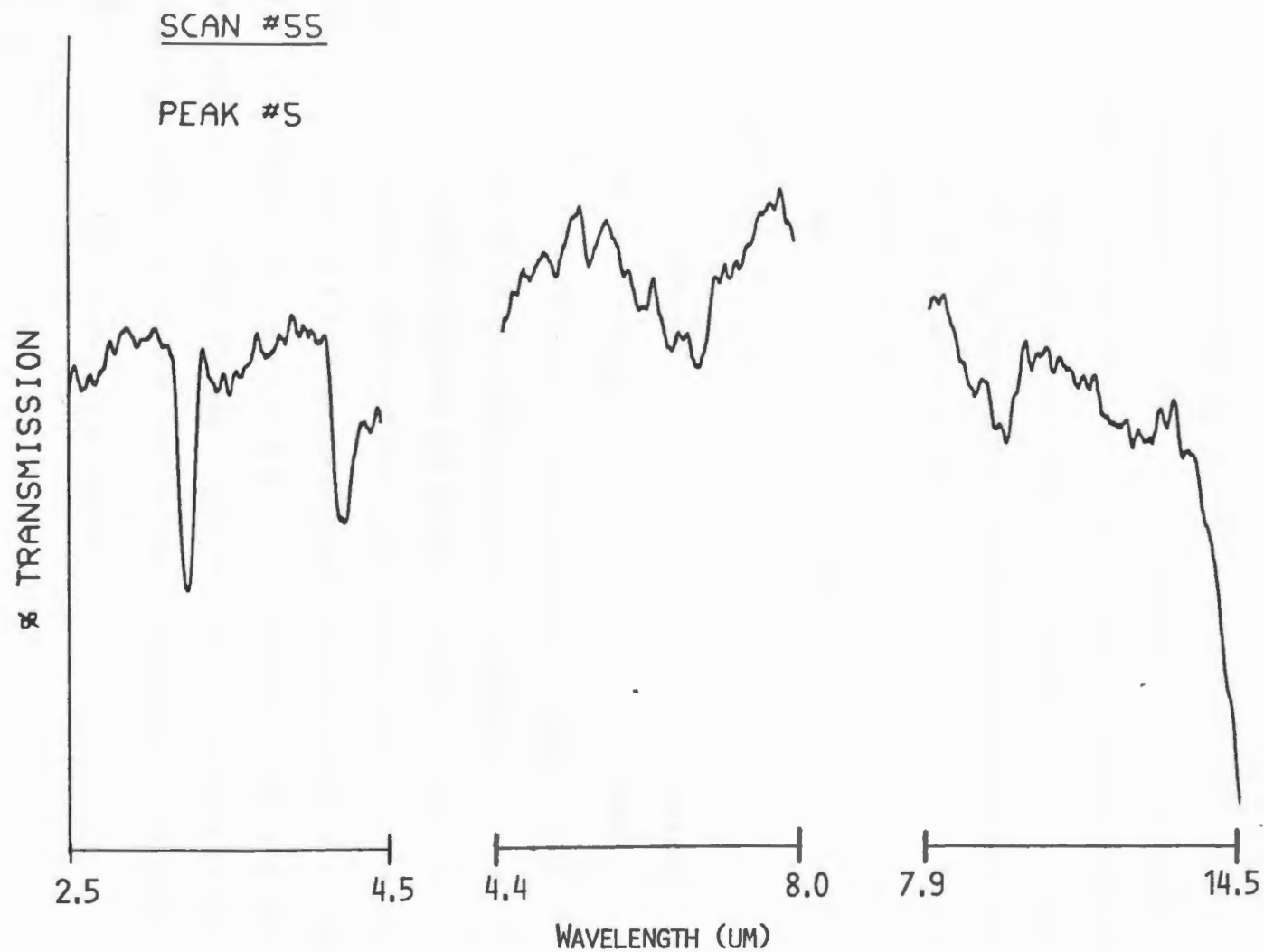


Figure 13. Infrared spectrum for scan #55 which corresponds to the fifth TWIRP peak of figure 6 and is due to benzene.

characteristic of an -OH stretch such as in an alcohol. This coupled with the fact that it is a low boiling point compound narrows the search for an unknown compound considerably. For scan #55 an absorption at 3.2 μm signifies an aromatic compound. Thus, based on very little interpretation the analyst would know the unknown mixture contains five low boiling compounds and that one is an alcohol and one is an aromatic. If only a gc was available, the analyst would know there were five compounds in the mixture and that they were low boiling. The additional IR data provides a great amount of valuable information and significantly narrows the choice of possible candidates for positive identification.

During this investigation we also conducted a comparison between the GC/IR TWIRP and a standard gc flame ionization detector (FID). The same five component mixture used previously was analyzed with the gc column reconnected to the FID as in normal gc operation. Figure 14 shows the comparison between the normal gc output and the TWIRP from the GC/IR. Both plots are similar in that they identify a five component mixture and just as the integrated IR detector signal for each particular compound is different so it is for the FID where the electron rich aromatic ring of benzene gives a much stronger signal than the other compounds. The FID analysis confirms the GC/IR results in regards to the number of components in the mixture.

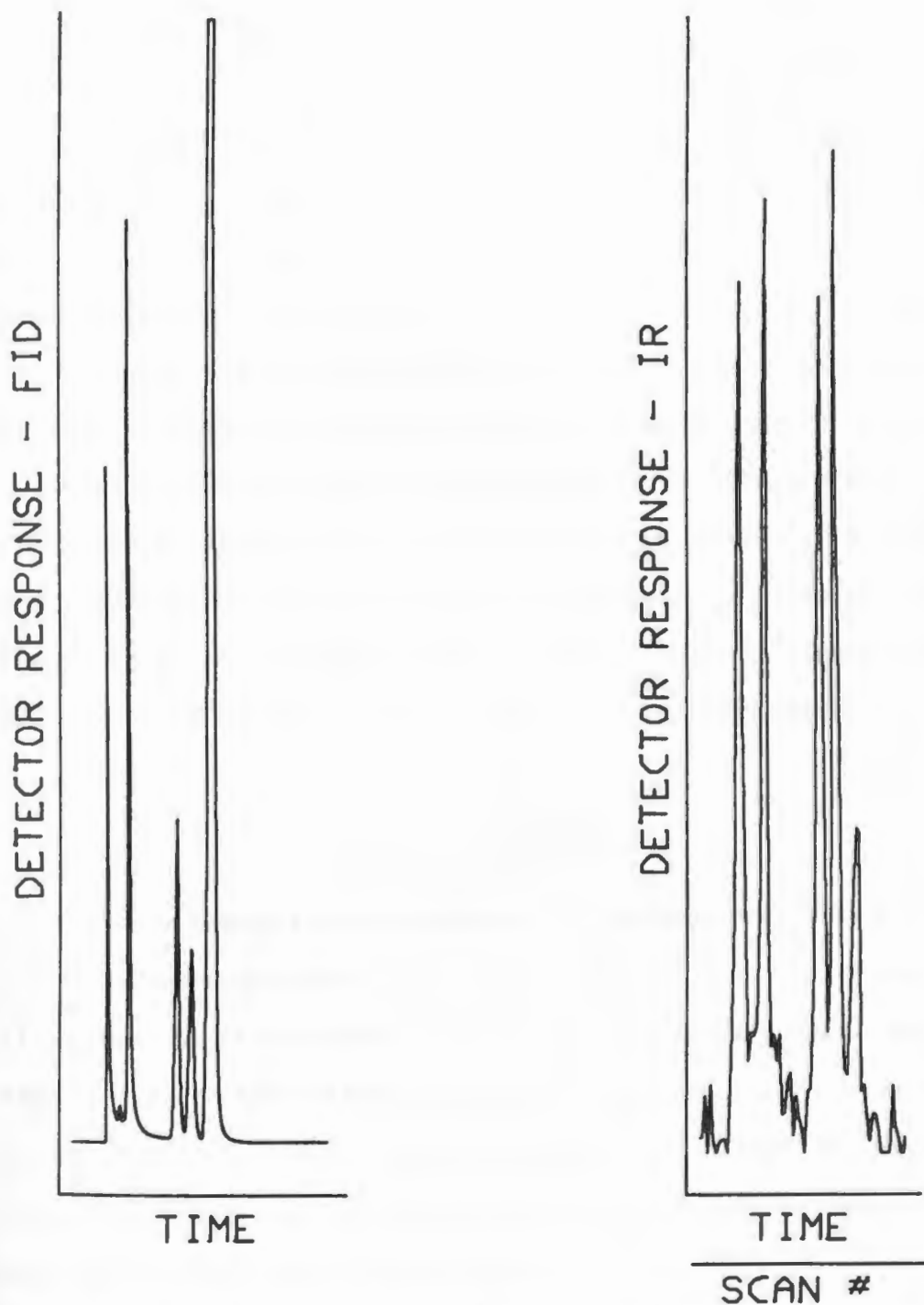


Figure 14. Comparison of a gas chromatogram (flame ionization detector) to the TWIRP reconstructed gas chromatogram of a five component mixture of acetone, chloroform, methanol, carbon tetrachloride, and benzene.

While checking the reproducibility of the CVF scan and data acquisition we decided to try signal averaging the spectra to decrease the noise in a plotted spectrum. Figure 15 shows the original spectrum and an averaged spectrum where four spectra were added together and then signal averaged. As can be seen there is a noticable reduction in noise (the original spectra were made noisy on purpose by running at a higher than normal amplifier gain). Also a three point smooth technique was used and again this resulted in noise reduction. As in any case when spectra are digitized and stored in memory a wide variety of signal processing options are possible including expansion, subtraction, and spectral searches to match an unknown to a known. In fact a laboratory which does routine analysis would maintain a spectral library of known compounds as well as a listing of gc elution times for a given gc column at known temperatures and flow rates.

Conclusion

We have demonstrated a simple GC/IR/Computer system to analyze multi-component mixtures "on-the-fly". The system is demonstrated on a five component mixture with closely related boiling points. The resulting reconstructed gas chromatograms and associated IR spectra provide sufficient information to identify each of the five components. This GC/IR consists of a rugged circular variable filter spectrometer and is ideally suited for field use.

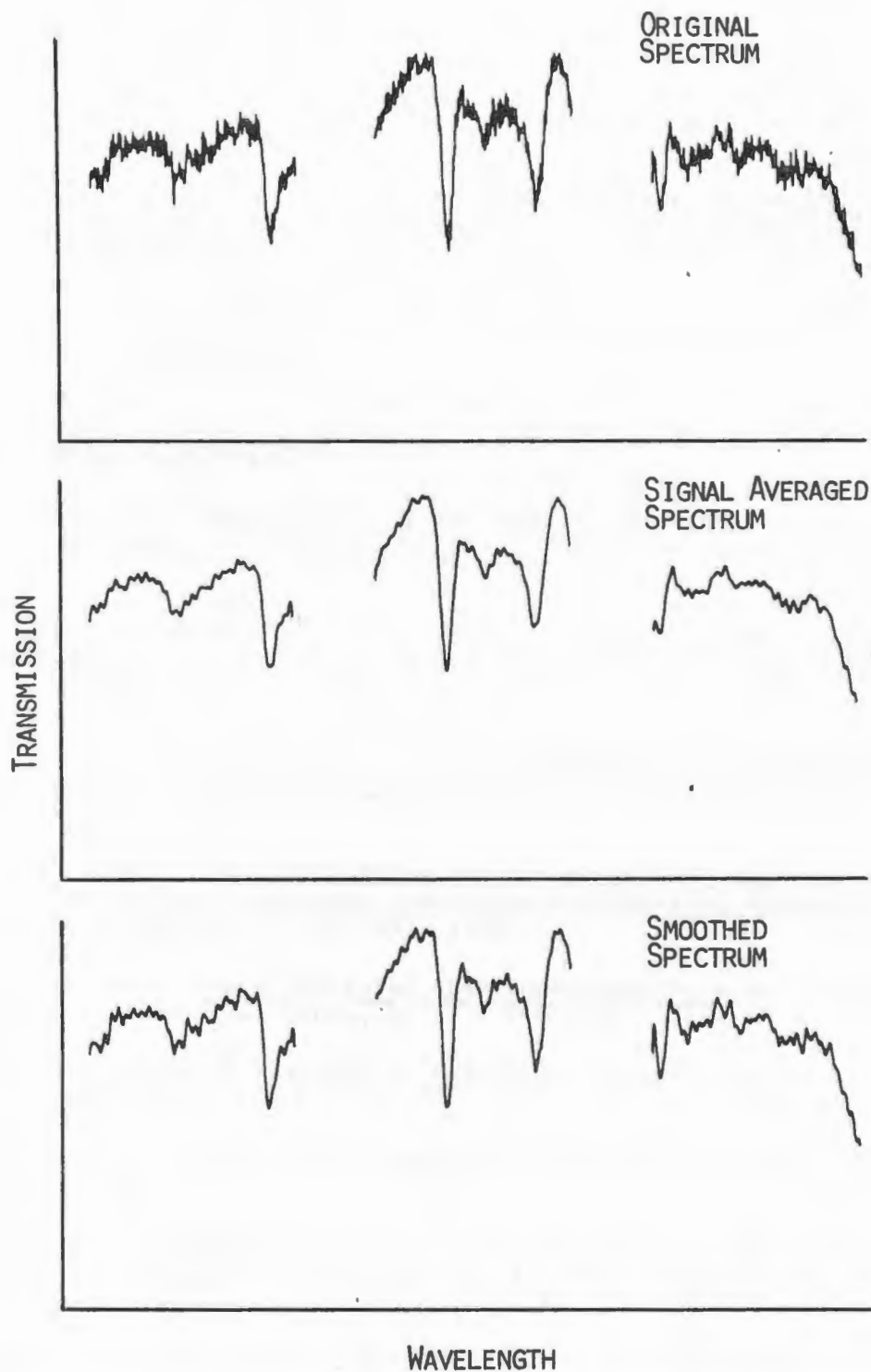


Figure 15. Computer processing of Beckman IR-102 spectra demonstrating signal averaging and smoothing techniques.

REFERENCES

1. J.G. Grasselli and M.K. Snavely, Analysis of Organic Reaction Products by Combined Infrared - Gas Chromatography Techniques, Applied Spectroscopy, (6), 90-4, 1952.
2. H.E. Bellis and E.J. Slowinski, Jr., Application of Vapor Chromatography to Infrared Spectroscopy of Liquids, J. Chem. Phys., (25), 794, 1956.
3. D.M.W. Anderson and J.L. Duncan, Identification and Quantitative Estimation of Gas-Chromatography Fractions by Vapor-Phase Infrared Spectroscopy, Chem. & Ind., 662, 1958.
4. D.M.W. Anderson, Applications of Infrared Spectroscopy: The Identification and Determination of Gas-Chromatographic Fractions, Analyst, (84), 50-5, 1959.
5. W.L. Senn, Jr. and H.V. Drushel, Organic Composition Analysis by Combination Gas-Chromatography with Infrared Spectroscopy, Anal. Chem. Acta, (25), 328-33, 1961.
6. P.A. Wilks, Jr. and R.A. Brown, Construction and Performance of an Infrared Chromatographic Fraction Analyzer, Analytical Chemistry, 36(10), 1896-9, 1964.
7. K.H. Kubeczka, Simple Arrangement for the Collection of Very Small, Gas Chromatographically Separated Substances for Infrared Spectroscopy, Naturwissenschaften, 52(14), 492, 1965.
8. G.P. Nicolais, Use of Infrared Spectrophotometry in Gas Chromatography: New Designs, Kem. Teollisuus, 23(12), 1095-102, 1966.
9. M.J.D. Low, Rapid IR Analysis of Gas-Chromatography Peaks, Chem. Commun., (12), 371-2, 1966.
10. N.W.R. Daniels, GC-IR (Gas-Chromatography-Infrared) Microanalysis, Column, 2(1), 2-5, 1967.
11. H. Bober and K. Buerner, Optimum Operation Conditions for the Combination of a Gas Chromatograph with the Fast Scanning Spectrophotometer IR-102, Fresenius Z. Anal. Chem., 238(1), 1-19, 1968.
12. R.J. Bell, Introductory Fourier Transform Spectroscopy, Academic Press, 1972.

BIBLIOGRAPHY

1. Ahmadjian, M. and Brown, C.W., IR/GC with an Inexpensive Non-Dispersive Spectrometer, The Pittsburgh Conference on Analytical Chemistry and Applied Spectroscopy, 1980.
2. Anderson, D.M.W., Applications of Infrared Spectroscopy: The Identification and Determination of Gas-Chromatographic Fractions, Analyst, (84), 50-5, 1959.
3. Anderson, D.M.W. and Duncan, J.L., Identification and Quantitative Estimation of Gas-Chromatography Fractions by Vapor-Phase Infrared Spectroscopy, Chem. & Ind., 662, 1958.
4. Azarraga, L.V., Gold Coating of Glass Tubes for Gas Chromatography Fourier Transform Infrared Spectroscopy "Light Pipe" Gas Cells, Appl. Spectry., 34(2), 224-5, 1980.
5. Azarraga, L.V. and McCall, Infrared Fourier Transform Spectrometry of Gas Chromatography Effluents, EPA-660/2-73-034, 1974.
6. Bellis, H.E., and Slowinski, E.J. Jr., Application of Vapor Chromatography to Infrared Spectroscopy of Liquids, J. Chem. Phys., (25), 794, 1956.
7. Bober, H. and Buerner, K., Optimum Operation Conditions for the Combination of a Gas Chromatograph with the Fast Scanning Spectrophotometer IR-102, Fresenius Z. Anal. Chem, 238(1), 1-19, 1968.
8. Brown, R.A., Kelliher, J.M., James, M., Heigl, J.J., Warren, C.W., and Charles, W., Rapid Scan Infrared Spectrometer for Operation with Support Coated Open Tubular or Packed Column Gas Chromatographs, Anal. Chem., 43(3), 353-8, 1971.
9. Daniels, N.W.R., GC-IR (Gas-Chromatography-Infrared) Microanalysis, Column, 2(1), 2-5, 1967. Fresenius Z. Anal. Chem, 238(1), 1-19, 1968.
10. Dijkstra, G., Coupling of a Gas Chromatograph with an IR Spectrometer, Lab. Instrum., 10(1), 19-20,23, 1974.
11. Grasselli, J.G. and Snavely, M.K., Analysis of Organic Reaction Products by Combined Infrared - Gas Chromatography Techniques, Applied Spectroscopy, (6), 90-4, 1952.
12. Kizer, K.L., GC/IR Analysis, Amer. Lab., 5(6), 40-2, 44-5, 1973.

13. Krakow, B., Continual Analysis of Gas Chromatographic Effluents by Repetitive Infrared Scanning, Anal. Chem., 41(6), 815-19, 1969.
14. Kubeczka, K.H., Simple Arrangement for the Collection of Very Small, Gas Chromatographically Separated Substances for Infrared Spectroscopy, Naturwissenschaften, 52(14), 492, 1965.
15. Low, M.J.D., Rapid IR Analysis of Gas-Chromatography Peaks, Chem. Commun., (12), 371-2, 1966.
16. Nicolais, G.P., Use of Infrared Spectrophotometry in Gas Chromatography: New Designs, Kem. Teollisuus, 23(12), 1095-102, 1966.
17. Pendry, A., Automated GLC/IR Analysis, Spectrovision, 19, 5-8, 1968.
18. Penzias, G.J., and Boyle, M.J., Combining Infrared and Gas Chromatography for On-Line Analysis, Amer. Lab., (Oct), 53-9, 1973.
19. Senn, W.L. Jr. and Drushel, H.V., Organic Composition Analysis by Combination Gas-Chromatography with Infrared Spectroscopy, Anal. Chem. Acta, (25), 328-33, 1961.
20. Shaps, R. and O. Dannhauser, GC-IR Coupling System, Chem-Anlagen Verfahren, (9), 94, 101, 103, 1976.
21. Skoog, D.A. and West, D.M., Principles of Instrumental Analysis, 1980.
22. Wall, D.L. and Mantz, A.W., High Sensitivity Infrared Spectroscopy of Gas Chromatographic Peaks, Appl. Spectry, 31(6), 552-60, 1977.
23. Wilkins, C.L., Giss, G.N., Brissey, G.M., and Steiner, S., Direct Linked Gas Chromatography-Fourier Transform Infrared-Mass Spectrometer System, Anal. Chem., 53, 113-17, 1981.
24. Wilks, P.A. Jr., Gas Chromatography and Infrared Spectroscopy: A Potent Analytical Team, Amer. Lab., (March), 14-20, 1969.
25. Wilks, P.A. Jr. and Brown, R.A., Construction and Performance of an Infrared Chromatographic Fraction Analyzer, Analytical Chemistry, 36(10), 1896-9, 1964.

APPENDIX

Program Listings and Descriptions

These computer programs were written for use on a Data General Corporation (DGC) NOVA 3/12 mini-computer system with an analog digital converter. The programs were written in Fortran and DGC assembly language.

Program: GCIR.FR

Purpose: This is the main controller program, i.e. executor, for the GC/IR/Computer System.

Method: GCIR.FR allows specific run time programs to be executed via swap commands.

Input: The user needs to specify at the prompt which real time programs are to be run. They are:

1. Data acquisition of the infrared data from the IR-102 signal conditioning electronics.
2. Wavelength calibration of the spectral data.
3. Signal averaging.
4. Plotting of the reconstructed gas chromatograms, TWIRPs, and spectra.

Output: User prompts for the above functions.

THIS IS THE MAIN CONTROLLER PROGRAM FOR THE GCIR COMPUTER SYSTEM. THE RUN TIME PROGRAMS ARE BROUGHT INTO EXECUTION VIA SWAP COMMANDS. THE FOLLOWING PROGRAMS ALONG WITH THEIR SUBROUTINES AND LIBRARIES ARE USED:

GCIRAD.FR : TO ACQUIRE DATA AND STORE THIS DATA IN
CONTIGOUS BLOCKS ON DISK FOR SUBSEQUENT PROCESS
AND DISPLAY.
GCIRSUB.FR
ADGCIR.SR
FORT.LB

GCIRCALIB.FR : TO CALIBRATE THE SPECTRAL DATA ON A FREQUENCY
SCALE
THE CALIBRATION DATA IS THE STORED
FOR SUBSEQUENT USE BY THE PLOTTER PROGRAM.
NPLOT.SR
FORT.LB

GCIRAVG.FR : TO SIGNAL AVERAGE SPECTRA
FORT.LB

GCIRPLOT.FR : TO PLOT OUT SPECTRA
PLOTTER4NS.LB
FORT.LB

```
20      CALL SWAP('GCIRCALIB.SV',IER)
        GO TO 100

30      CALL SWAP('GCIRAVG.SV',IER)
        GO TO 100

40      CALL SWAP('GCIRPLOT.SV',IER)
        GO TO 100

100     CONTINUE

        WRITE(10,110)
110     FORMAT(' RUN AGAIN?  00=NO, 01=YES  ')
        READ(11,120) JA
120     FORMAT(I2)
        IF(JA.EQ.1) GO TO 1
        STOP
        END
```

```

1      CONTINUE

      TYPE'DATA ACQUISITION      =01 '
      TYPE'CALIBRATION           =02 '
      TYPE'SIGNAL AVERAGING      =03 '
      TYPE'PLOTting              =04 '

      ACCEPT'? ',J

      GO TO(10,20,30,40),J

10     CALL SWAP('GCIRAD.SV',IER)
      GO TO 100

```

Program: GCIRAD.FR

Purpose: This program and its subroutines acquire data from the Beckman
IR-102 spectrometer.

Method: GCIRAD.FR and its subroutines samples the analog signal from the
IR-102 signal conditioning electronics, converts the analog data
to digital via the DGC ADC, and stores the digitized data on disk
for subsequent processing and display.

Input: The user needs to specify a file name into which to store the data
and an ADC sampling rate.

Output: Output are contiguous disk file containing the digitized spectral
data from the GC/IR.

THIS PRG AND ITS SUBROUTINES ACQUIRE DATA FROM THE GCIR AT EVERY CLOCK TICK(IF TIME DELAY=1) A DATA POINT IS TAKEN THE IR DETECTOR ELECTRONICS AND THIS DATA IS STORED IN BLOCKS.

SUBROUTINES: GCIRSUB.FR
ADGCIR.SR

```

        DIMENSION JSPEC(10), JDAT(8192)
        WRITE(10,2)
2       FORMAT(' SPECTRUM ID? ')

        READ(11,3) JSPEC(1)
3       FORMAT(S14)

        ACCEPT'# OF DISK TRANSFERS TO BE EFFECTED? ',N
        ACCEPT'# OF DATA PTS PER DISK TRANSFER? ',I3
        ACCEPT'# OF CLOCK PULSES BETWEEN DATA POINTS? ',I4

        IF(I3.GT.256) GO TO 5

C       I6=# OF DATA BLOCKS/SCAN

        I6=1
        GO TO 30

5       D3=I3/256,
        J3=D3
        R3=D3-I3
        IF(R3) 20,10,20
10      I6=J3
        GO TO 30
20      I6=J3+1
30      CONTINUE

C       IFILE=# OF DATA BLOCKS PER ENTIRE RUN
        IFILE=I6*(N+3)

C       I7=SIZE OF DIMENSION DATA ARRAY
        I7=256*I6

C       SET UP DATA FILE ON DISK
        CALL DFILW(JSPEC,IER)
        TYPE'DFILW IER =' ,IER

```

```
CALL CFILW(JSPEC,3,IFILE,IER)
IF(IER.NE.1) TYPE'ERROR ERROR ERROR ERROR'
TYPE 'CFILW IER =' ,IER

CALL OPEN(0,JSPEC,2,IER)
TYPE'OPEN  IER =' ,IER

CALL GCIRSUB(JDAT,I3,I4,I6,I7,N)

CONTINUE
STOP
END
```

Program: GCIRSUB.FR

Purpose: This is a fortran subroutine used by GCIRAD.FR to initiate data acquisition and write the data onto disk.

Method: Based on the inputs from GCIRAD.FR specifying the file name and the length of time for data acquisition this program will run an assembly language subroutine, ADGCIR.SR, to sample the ADC and write the data out to a disk.

Output: Contiguous data blocks on disk.

```

C      SUBROUTINE FOR GCIRAD.FR (ACQUIRING DATA FROM THE GCIR AND
C      STORING THIS DATA IN CONTIGUOUS DATA BLOCKS ON DISK).
C
C      *****
C
C      SUBROUTINE GCIRSUB (JDAT, I3,I4,I6,I7,N)
C      DIMENSION JDAT(I7), INFO(256)
C      COMMON IDAT(8192)
C      EQUIVALENCE (IDAT(1),IARRAY)
C      DO 90 I=1,I7
90      JDAT(I)=0
C      INFO(1)=I3
C      INFO(6)=I4
C      INFO(11)=I6
C      INFO(16)=I7
C      INFO(21)=N
C
C      WRITE(10,100)
100      FORMAT(' HIT RETURN TO START DATA ACQUISITION ')
C      READ(11,110) KR
110      FORMAT(I2)
C      J=2
C      LR=0
C      DO 140 IL=1,N
115      CONTINUE
C
C      CALL ADGCIR(IARRAY,I4,I3)
C
C      DO 130 I=1,I3
130      JDAT(I)=IDAT(I)
C
C      CALL WRBLK(0,J,JDAT,I6,IER)
C      IF (IER.NE.1) TYPE'ERROR ERROR ERROR ERROR'
C
C      J=J+I6
C
C      INFO(21)=IL
C      CALL WRBLK(0,1,INFO,1,IER)
140      LR=LR+1
C      CALL WRBLK(0,1,INFO,1,IER)
C      TYPE IER
C      TYPE 'SET UP FOR CALIBRATION RUN? 00-NO, 01-YES'
C      READ(11,110) JCL
C      IF (JCL.EQ.0) GO TO 150
C      PAUSE HIT RETURN TO START CALIBRATION RUN
C      GO TO 115
150      CONTINUE
C      RETURN
C      STOP
C      END

```

Program: ADGCIR.SR

Purpose: This assembly language subroutine runs the ADC and stores the digitized data into RAM.

Method: ADGCIR.SR samples the ADC channel #0 (multiplexed), stores the data into an accumulator, and then loads the accumulator into a RAM array for subsequent reading and writing by GCIRSUB.FR

Inputs: No user inputs are required.

Output: Digitized infrared spectral data in RAM storage.


```

        .TITL ADGCIR
        .NREL
        .TXIM 1
        .ENT ADGCIR
        .EXTU
        .EXTN .I

TTY=10
CLK=14
ADCONV=21
FS.=3
T.=-167

        FS.
ADGCIR: JSR @.FARL
        STA 0,AC0
        STA 1,AC1
        STA 2,AC2
        LDA 0,T.+0,3
        STA 0,ARRAY
        LDA 0,@T.+1,3    ;GET DELAY
        STA 0,STIME
        LDA 0,@T.+2,3    ;GET # OF DATA POINTS
        STA 0,COUNT
GO:      LDA 0,CNTP
        STA 0,CNTEMP
REPT:    LDA 1,STIME
        .SYSTEM
        .DELAY
        JSR ERR
        LDA 0,C5
        DOA 0,ADCONV
        INTDS
        NIOS ADCONV
        SKPDN ADCONV
        JMP .-1
        DICC 0,ADCONV
        INTEN
        LDA 2,ARRAY
        STA 0,0,2
        INC 2,2
        STA 2,ARRAY
        SUBZL 0,0
        DSZ COUNT
        JMP GO
RET:     LDA 0,AC0
        LDA 1,AC1
        LDA 2,AC2
        JSR @.FRET

```

Program: GCIRCAL.FR

Purpose: This program transforms the digitized infrared data into spectra.

Method: GCIRCAL.FR calibrates the infrared spectra on a wavelength scale using the 0 VDC signal between the circular variable filter elements one and three.

Inputs: The user needs to specify the file name the digitized data is stored in and a new file name to write the calibrated data out to.

Output: Output is a random data file containing calibrated spectra.

THIS PRG IS TO DETERMINE THE ARRAY BOUNDARIES OF INDIVIDUAL
SPECTRA ACQUIRED FROM THE GCIR SYSTEM VIA PRG GCIR...
AND TO STORE THE INDIVIDUAL SPECTRA IN A RANDOM DATA FILE

```
DIMENSION JSPEC(10),INFO(256),JDAT(8192),IPT(256)
DIMENSION YDAT(500),KI(5),KDAT(1051),TDAT(1050)
WRITE(10,2)
FORMAT(' SPECTRUM ID?  ')
READ(11,3) JSPEC(1)
FORMAT(S14)
```

```
CALL OPEN(0,JSPEC,2,IER)
TYPE'OPEN  IER =' ,IER
```

```
CALL RDBLK(0,1,INFO,1,IER)
TYPE'RDBLK IER =' ,IER
```

```
I3=INFO(1)
I6=INFO(11)
I7=INFO(16)
N=INFO(21)
TYPE'# OF DATA PTS/TRANSFER =' ,I3
TYPE'# OF DISK  TRANSFERS    =' ,N
```

```
TYPE'NEW FILE NAME?  '
READ(11,3) JSPEC(1)
CALL DFILW(JSPEC,IER)
TYPE 'DFILW IER =' ,IER
CALL CFILW(JSPEC,2,IER)
TYPE'CFILW IER =' ,IER
CALL OPEN(1,JSPEC,2,IER,2)
TYPE 'OPEN  IER =' ,IER
ACCEPT'SIGNAL AVG? 00=NO, 01=YES  ',KSIQ
IP=5
IC=0
DO 300 JS=1,N
I4=2+(I6*JS)-I6
```

```
CALL RDBLK(0,I4,JDAT,I6,IER)
TYPE'RDBLK IER =' ,IER
```

```
IL=1
```

```
DO 10 I=1,13
IF(JDAT(I).GT.300) GO TO 10
```

```
K=I+160
DO 20 J=I,K
IF(JDAT(J).LE.300) GO TO 20
I=J
GO TO 10
CONTINUE
```

CS: 6
CL: 2
CH: 3
AC0: 0
AC1: 0
AC2: 0
DTIME: 300.
STIME: 0
COUNT: 0
CNTF: 1.
CNTEMP: 0
ARRAY: 0
ERR: .SYSTEM
.BREAK
.END

Program: TWIRP.FR

Purpose: This program will plot the reconstucted gas chromatograms.

Method: TWIRP.FR inputs the infrared spectral data from the random files, sums the intensities/spectrum, and plots the sums as a function of scan number.

Inputs: The user needs to specify the file name, noise filtering if desired, and expansion factors for the plotter.

Output: Output is a Total Wavelength IR Plot, TWIRP, i.e. a gas chromatogram.

```

DO 30 J=I,I3
IF(JDAT(J).GE.1000) GO TO 25
30  CONTINUE
GO TO 11

25  IF((J+1050).GT.I3) GO TO 11

IPT(IL)=J
I=J
IL=IL+1

10  CONTINUE
11  CONTINUE

IL=IL-1
IF(KSIG.EQ.1) IL=IL-1
DO 100 I=1,IL
IPT1=IPT(I)
IPT2=IPT(I)+1050

L=1
DO 200 K=IPT1,IPT2
KDAT(L)=JDAT(K)
200  L=L+1
IF(KSIG.EQ.0) GO TO 202

IPT1=IPT(I+1)
IPT2=IPT(I+1)+1050

L=1

DO 201 K=IPT1,IPT2
KDAT(L)=(KDAT(L)+JDAT(K))/2
201  L=L+1
202  CONTINUE

IC=IC+1
YSUM=0.0

DO 210 K=1,1050
TDAT(K)=KDAT(K)
TDAT(K)=TDAT(K)*10./4096.+1
TDAT(K)=(ALOG10(TDAT(K)))*10.
TDAT(K)=10.-TDAT(K)
210  YSUM=YSUM+TDAT(K)

YDAT(IC)=YSUM

TYPE' IC =',IC

100  CONTINUE
300  CONTINUE

```

```

      JV=1
      DO 9 J=1,IC
      IF(YDAT(JV).GE.YDAT(J)) GO TO 9
      JV=J
9      CONTINUE

      XF=4096/YDAT(JV)

      DO 8 J=1,IC
8      KDAT(J)=YDAT(J)*XF
      KI(1)=IC

      CALL WRTR(1,0,KI,5,IER)
      TYPE 'IER =',IER

      CALL WRTR(1,5,KDAT,IC,IER)
      TYPE 'IER =',IER

      TYPE KDAT
      CALL RESET
      STOP
      END
R

```

Program: GCIRSPEC.FR

Purpose: This program plots out individual spectra corresponding
from the TWIRP.

Method: GCIRSPEC.FR obtains the calibrated spectral data from the
data file and plots it as a function of wavelength.

Inputs: The user needs to specify the file name, the spectrum number,
plotter expansion factors.

Output: Infrared spectra

THIS PRG IS TO DETERMINE THE ARRAY BOUNDARIES OF INDIVIDUAL
SPECTRA ACQUIRED FROM THE GCIR SYSTEM VIA PRG GCIR....
AND TO DETERMINE THE TOTAL WAVELENGTH IR PLOT IN A MANNER
ANALOGOUS TO THAT USED FOR THE GC/MS SYSTEM PRG TOTIUM.FR

```

        DIMENSION JSPEC(10),KI(5),KDAT(1050),JDAT(1050),YDAT
        DIMENSION ZDAT(1000),FLTY(300),FLTX(300)
        CALL INITIAL(6,200,10.5,30.0)
1      CONTINUE
        CALL PLOT(0.0,2.0,-3)
        TI=8.0/4096
        WRITE(10,2)
2      FORMAT(' SPECTRUM ID?  ')
        READ(11,3) JSPEC(1)
3      FORMAT(S14)

        CALL OPEN(0,JSPEC,2,IER,2)
        TYPE'OPEN  IER =',IER
        CALL RDRW(0,0,KI,5,IER)
        TYPE'RDRW  IER =',IER

        IC=KI(1)
        TYPE'IC =',IC
        CALL RDRW(0,5,KDAT,IC,IER)
        TYPE 'WRTR IER =',IER
        GO TO 401
        IP=KI(1)-1
        L=0
        IPT=5

        TYPE'# OF SPECTRA =',IP
        ACCEPT'SIGNAL AVG?  00=NO, 01=YES'  ,ISIG
        DO 300 I=1,IP
        CALL RDRW(0,IPT,JDAT,1050,IER)
        TYPE'RDRW  IER =',IER

        IF(ISIG.EQ.0) GO TO 4
        IPT1=IPT+1050
        CALL RDRW(0,IPT1,KDAT,1050,IER)
        TYPE 'RDRW  IER =',IER
4      YSUM=0.0

        DO 200 J=1,1050
        JDAT(J)=4096-JDAT(J)
        IF(ISIG.EQ.0) GO TO 200
        KDAT(J)=4096-KDAT(J)
        JDAT(J)=(JDAT(J)+KDAT(J))/2
200     YSUM=YSUM+JDAT(J)

```

```

      L=L+1
      ZDAT(L)=YSUM
      IFT=IFT+1050
300    CONTINUE

401    CONTINUE
      L=IC
      DO 402 I=1,L
402    ZDAT(I)=KDAT(I)
      JM=1
      DO 7 I=1,L
      IF(ZDAT(JM),LE,ZDAT(I)) GO TO 7
      JM=I
7      CONTINUE
      ZMIN=ZDAT(JM)
      TYPE'MIN VALUE =',ZMIN

      DO 8 I=1,L
8      YDAT(I)=ZDAT(I)-ZMIN
C      SCALE DATA
      JV=1
      DO 9 IK=1,L
      IF(YDAT(JV),GE,YDAT(IK)) GO TO 9
      JV=IK
9      CONTINUE

      XF=8.0/YDAT(JV)
      DO 400 IK=1,L
      YDAT(IK)=YDAT(IK)*XF
400    CONTINUE

      ACCEPT'LOW LEVEL DIGITAL FILTER IN PLOTTER UNITS?
      DO 500 IK=1,L
      YDAT(IK)=YDAT(IK)-XL
      IF(YDAT(IK),LT,0.0) YDAT(IK)=0.0
500    CONTINUE

      ACCEPT'Y-EXPANSION FACTOR FOR PLOTTER? ',XY
      ACCEPT'X-AXIS INCREMENT PER SCAN? ',XX
      DO 600 IK=1,L
      YDAT(IK)=YDAT(IK)*XY
      IF(YDAT(IK),GT,8.0) YDAT(IK)=8.0

      YDAT(IK)=8.-YDAT(IK)+1.0
      YDAT(IK)=(ALOG10(YDAT(IK)))*10.
      YDAT(IK)=10.-YDAT(IK)

      YSUM=YDAT(IK)
      XFLT=IK*XX
      FLTX(IK)=XFLT
600    FLTY(IK)=YSUM
      CALL LINE(FLTX,FLTY,L,0,1)

```

```

CALL PLOT(XPLT,-.25,3)
CALL PLOT(0.0,-.25,-2)

DO 700 1=1,N
X=XX
CALL PLOT(X,.1,3)
CALL PLOT(X,-.1,2)
700 CONTINUE
CALL PLOT(0.0,-1.75,-3)
GO TO 401
CALL RESET
GO TO 1

STOP
END

```

```

C      THIS PRG IS TO PLOT INDIVIDUAL SPECTRA
C
C      *****
DIMENSION JSPEC(10),INFO(256),JDAT(1280)
DIMENSION XAXIS(4096),YDAT(4096),KI(5)
CALL INITIAL(6,200,10.0,15.0)
WRITE(10,2)
2      FORMAT(' SPECTRUM ID?  ')
      READ(11,3) JSPEC(1)
3      FORMAT(S14)

      CALL OPEN(0,JSPEC,2,IER,2)
      TYPE'OPEN IER =' ,IER

      CALL PLOT(0.0,0.0,3)
1      ACCEPT'SPECTRUM # TO BE PLOTTED?' ,IP
      ACCEPT'Y-AXIS IN INCHES?  ' ,TY
      ACCEPT'X-AXIS/1050 PTS ?  ' ,XAX
      XY=XAX/1050
      TI=TY/4096

      IP=IP-1
      KF=IP*5+1
      CALL RDBLK(0,KF,JDAT,5,IER)
      TYPE'RDBLK IER =' ,IER
      DO 7 J=1,1050
      XAXIS(J)=J*XY
7      YDAT(J)=JDAT(J)*TI

      DO 8 I=5,1046
      YDAT(I-4)=(-YDAT(I-4)*21.+YDAT(I-3)*14.+YDAT(I-2)
1      *54.+YDAT(I)*59.+YDAT(I+1)*54.+YDAT(I+2)*39.+YDAT
2      YDAT(I+4)*21.)/231.
      IF(YDAT(I-4).LT.0.) YDAT(I-4)=0.
8      CONTINUE

      X=XAXIS(1)
      Y=YDAT(1)

      CALL LINE(XAXIS,YDAT,1050,0,1)

      X=1070*XY
      CALL PLOT(X,0.0,-1)

      GO TO 1
      CALL RESET
      STOP
      END

```

Program: GCIRSUM.FR

Purpose: This program will signal average any number of infrared spectra.

Method: GCIRSUM.FR sums corresponding spectral points for the specified spectra and finds the average.

Inputs: The user needs to specify the file name and the spectra numbers to be signal averaged.

Output: Signal averaged spectra.

THIS PRG IS TO DETERMINE THE ARRAY BOUNDARIES OF INDIVIDUAL
SPECTRA ACQUIRED FROM THE GCIR SYSTEM VIA PRG GCIR...
AND TO STORE THE INDIVIDUAL SPECTRA IN A RANDOM DATA FILE

```

DIMENSION JSPEC(10),INFO(256),JDAT(8192),IPT(256)
DIMENSION YDAT(500),KI(5),KDAT(1051)
WRITE(10,2)
2  FORMAT(' SPECTRUM ID?  ')
  READ(11,3) JSPEC(1)
3  FORMAT(S14)

CALL OPEN(0,JSPEC,2,IER)
TYPE'OPEN  IER =',IER

CALL RDBLK(0,1,INFO,1,IER)
TYPE'RDBLK IER =',IER

I3=INFO(1)
I6=INFO(11)
I7=INFO(16)
N=INFO(21)
TYPE'# OF DATA PTS/TRANSFER =',I3
TYPE'# OF DISK  TRANSFERS   =',N

TYPE'NEW FILE NAME?  '
READ(11,3) JSPEC(1)
CALL DFILW(JSPEC,IER)
TYPE 'DFILW IER =',IER
CALL CFILW(JSPEC,2,IER)
TYPE'CFILW IER =',IER
CALL OPEN(1,JSPEC,2,IER,2)
TYPE 'OPEN  IER =',IER
ACCEPT'SIGNAL AVG? 00=NO, 01=YES  ',KSIQ
IP=5
IC=0
DO 300 JS=1,N
I4=2+(I6*JS)-I6

CALL RDBLK(0,I4,JDAT,I6,IER)
TYPE'RDBLK IER =',IER

IL=1

DO 10 I=1,I3
IF(JDAT(I).GT.300) GO TO 10

```

```

      K=I+150
      DO 20 J=I,K
      IF(JDAT(J).LE.300) GO TO 20
      I=J
      GO TO 10
20    CONTINUE

      DO 30 J=I,I3
      IF(JDAT(J).GE.1000) GO TO 25
30    CONTINUE
      GO TO 11

25    IF((J+1050).GT.I3) GO TO 11

      IFT(IL)=J
      I=J
      IL=IL+1

10    CONTINUE
11    CONTINUE

      IL=IL-1
      IF(KSIG.EQ.1) IL=IL-1
      DO 100 I=1,IL
      IFT1=IFT(I)
      IFT2=IFT(I)+1050

      L=1
      DO 200 K=IFT1,IFT2
      KDAT(L)=JDAT(K)
200    L=L+1
      IF(KSIG.EQ.0) GO TO 202

      IFT1=IFT(I+1)
      IFT2=IFT(I+1)+1050

      L=1

      DO 201 K=IFT1,IFT2
      KDAT(L)=(KDAT(L)+JDAT(K))/2
201    L=L+1
202    CONTINUE

      IC=IC+1
      YSUM=0.0

      DO 210 K=1,1050
      KDAT(K)=4096-KDAT(K)
210    YSUM=YSUM+KDAT(K)

      YDAT(IC)=YSUM

      TYPE' IC =' ,IC

100    CONTINUE
300    CONTINUE

```

```

      JV=1
      DO 9 J=1,IC
      IF(YDAT(JV).GE.YDAT(J)) GO TO 9
      JV=J
9     CONTINUE

      XF=4096/YDAT(JV)

      DO 8 J=1,IC
8     KDAT(J)=YDAT(J)*XF
      KI(1)=IC

      CALL WRTR(1,0,KI,5,IER)
      TYPE 'IER =',IER

      CALL WRTR(1,5,KDAT,IC,IER)
      TYPE 'IER =',IER

      TYPE KDAT
      CALL RESET
      STOP
      END

```