

2020

A near-infrared probe for non-invasively monitoring cerebrospinal fluid flow by F-18-positron emitting tomography and fluorescence

H. Guo

H. Kommidi

C. C. Lekaye

J. Koutcher

M. S. Judenhofer

See next page for additional authors

Follow this and additional works at: <https://academicworks.medicine.hofstra.edu/publications>

 Part of the [Otolaryngology Commons](#)

Recommended Citation

Guo H, Kommidi H, Lekaye CC, Koutcher J, Judenhofer MS, Cherry SR, Wu AP, Akin O, Souweidane MM, Ting R, . A near-infrared probe for non-invasively monitoring cerebrospinal fluid flow by F-18-positron emitting tomography and fluorescence. . 2020 Jan 01; 10(1):Article 7764 [p.]. Available from: <https://academicworks.medicine.hofstra.edu/publications/7764>. Free full text article.

This Article is brought to you for free and open access by Donald and Barbara Zucker School of Medicine Academic Works. It has been accepted for inclusion in Journal Articles by an authorized administrator of Donald and Barbara Zucker School of Medicine Academic Works. For more information, please contact academicworks@hofstra.edu.

Authors

H. Guo, H. Kommidi, C. C. Lekaye, J. Koutcher, M. S. Judenhofer, S. R. Cherry, A. P. Wu, O. Akin, M. M. Souweidane, R. Ting, and +2 additional authors

ORIGINAL RESEARCH

Open Access



A near-infrared probe for non-invasively monitoring cerebrospinal fluid flow by ^{18}F -positron emitting tomography and fluorescence

Hua Guo^{1,2,3†}, Harikrishna Kommidi^{3†}, Carl C. Lekaye⁴, Jason Koutcher⁴, Martin S. Judenhofer⁵, Simon R. Cherry⁵, Amy P. Wu⁶, Oguz Akin⁷, Mark M. Souweidane⁸, Omer Aras⁷, Zhaohui Zhu^{1,2*} and Richard Ting^{3*} 

Abstract

Purpose: Knowing the precise flow of cerebrospinal fluid (CSF) is important in the management of multiple neurological diseases. Technology for non-invasively quantifying CSF flow would allow for precise localization of injury and assist in evaluating the viability of certain devices placed in the central nervous system (CNS).

Methods: We describe a near-infrared fluorescent dye for accurately monitoring CSF flow by positron emission tomography (PET) and fluorescence. IR-783, a commercially available near-infrared dye, was chemically modified and radiolabeled with fluorine-18 to give [^{18}F]-IR783-AMBF₃. [^{18}F]-IR783-AMBF₃ was intrathecally injected into the rat models with normal and aberrant CSF flow and evaluated by the fluorescence and PET/MRI or PET/CT imaging modes.

Results: IR783-AMBF₃ was clearly distributed in CSF-containing volumes by PET and fluorescence. We compared IR783-AMBF₃ (fluorescent at 778/793 nm, ex/em) to a shorter-wavelength, fluorescein equivalent (fluorescent at 495/511 nm, ex/em). IR783-AMBF₃ was superior for its ability to image through blood (hemorrhage) and for imaging CSF-flow, through-skin, in subdural-run lumboperitoneal shunts. IR783-AMBF₃ was safe under the tested dosage both in vitro and in vivo.

Conclusion: The superior imaging properties of IR783-AMBF₃ could lead to enhanced accuracy in the treatment of patients and would assist surgeons in non-invasively diagnosing diseases of the CNS.

Keywords: Positron emission tomography, Fluorescence, Non-invasively, CSF flow, CNS disease

Introduction

Advances in molecular imaging have enhanced our ability to non-invasively track molecules in patients. Yet, no current single imaging modality is ideal, as each individual imaging modality is limited by unique temporal, spatial, and depth (through-tissue) resolutions [1]. One popular strategy in the development of more universal

contrast agents is to combine different imaging modalities with synergistic properties, where one modality will complement the resolution shortcomings of another imaging modality [2]. One good combination is positron emission tomography (PET) and fluorescence imaging. Both ^{18}F -PET and fluorescence imaging allow the tracking of small molecules at sub-nanomolar concentration, making them among the most sensitive, and therefore also the potentially least toxic, of available imaging modalities [3–6]. Moreover, while PET emissions can be non-invasively quantitated through a patient's body, fluorescence imaging is ideal for visualizing submicron,

* Correspondence: 13611093752@163.com; rct2001@med.cornell.edu

[†]Hua Guo and Harikrishna Kommidi contributed equally to this work.

¹Department of Nuclear Medicine, Chinese Academy of Medical Sciences and Peking Union Medical College Hospital, Beijing 100730, China

³Department of Radiology, Molecular Imaging Innovations Institute (MI3), Weill Cornell Medical College, New York, NY 10065, USA

Full list of author information is available at the end of the article

histologic structure in superficial tissue. Therefore, by combining ^{18}F -PET and fluorescent imaging modalities within a single molecule, we can transcend the molecular resolution limits of any one imaging modality [2, 7–9].

Damage to the central nervous system (CNS) can significantly reduce a patient's quality of life, especially if not promptly addressed [10, 11]. Unfortunately, CNS injuries can be difficult to identify, especially in patients who are unconscious or who have spontaneous or deep-tissue cerebrospinal fluid (CSF) leakage, where the anatomical CSF leak source is not obvious. A popular strategy for indicating CNS damage is to introduce a contrast agent into the CNS [12–17], but tracking small molecules deep within the CNS represents a particular challenge for multiple reasons: (1) Toxicity to structures within the CNS result in particularly morbid prognoses; therefore, it is crucial that new contrast be useful at low concentrations to ensure non-toxicity [18]. (2) Dense tissue and the thick bone that shield the CNS can scatter both exciting and emitting photons, thereby reducing the quality of certain imaging modalities significantly. (3) Media such as CSF can move a molecule, thereby quickly reducing the usefulness of slower image acquisition strategies that involve faster scanning or signal averaging [13, 15, 19].

To meet this challenge, we previously developed a fluorescein-based PET/fluorescent probe (Fc-AMBF₃) for imaging within the CNS [20, 21]. To make this probe compatible with the most current FDA-approved intra-operative robotic systems, which are designed to visualize fluorescent near-infrared dyes (> 600 nm) [22], we further develop a near-infrared dye as a PET/fluorescent probe. In this study, we report the synthesis of a near-infrared dye, IR-783 derivative, IR783-AMBF₃, a.k.a. (((2-((E)-6'-((E)-2-(3,3-dimethyl-1-(4-sulfobutyl)-3H-indol-1-ium-2-yl)vinyl)-2'-(2-((E)-3,3-dimethyl-1-(4-sulfobutyl)indolin-2-ylidene)ethylidene)-2',3',4',5'-tetrahydro-[1,1'-biphenyl]-4-carboxamido)ethyl)dimethylammonium methyltrifluoroborate. IR783-AMBF₃ was radiolabeled via isotopic exchange radiolabeling with fluorine-18 to give a near-infrared PET/fluorescent probe. We explored the utility of IR783-AMBF₃ in multiple rat models bearing normal and aberrant CSF flow. We compared IR783-AMBF₃ with previously published Fc-AMBF₃. The ability to image CSF flow in deep tissue and through blood suggests that IR783-AMBF₃ can be used as a tool for diagnosing CNS diseases that involve CSF flow by both PET and fluorescence imaging.

Materials and methods

General synthetic methods

Commercial chemicals and solvents were obtained from suppliers and used as purchased. Endotoxin-free and mycoplasma-free, isotonic sterile 1 × PBS (pH 7.4) was

obtained from Corning (cellgro®, Product #21-040-CM). Mass spectra were measured with a Waters Acquity H class HPLC/SQD2 mass spectrometer (UPLC–MS) using Acquity UPLC 1.7 μm C18 100 Å, 2.1 × 50 mm column (186002350), with a 4.0 min 10–90% H₂O/acetonitrile (ACN) (0.05% TFA) gradient and a flow rate of 0.4 mL/min. Preparative HPLC was performed on an Agilent 1260 series HPLC equipped on a Phenomenex Luna C18(2)100 Å, 250 cm × 21.20 mm i.d. 10 μm reverse phase column (00G4253-P0 AX), with a 40 min 10–80% H₂O/ACN (0.05% TFA) gradient and a flow rate of 12 mL/min. Nuclear magnetic resonance (NMR) spectra were recorded in a deuterated solvent with a 500-MHz Bruker spectrometer. Fluoride-18 concentration was done in a Thermo Scientific 5-mL vial (React-Vial no. 13223). Purity of radiolabeling (≥ 95%) was verified on a Varian reverse phase HPLC, using a Waters Sunfire™ C18 3.5 μm, 4.6 mm × 50 mm column (186002551), an attached radiodetector, and a 10–90% H₂O/ACN (0.05% TFA) elution gradient with a flow rate of 2 mL/min.

Chemical synthesis


Reagents and conditions used to synthesize [^{18}F]-IR783-AMBF₃ are described in Scheme 1. IR-783 was coupled to 4-carboxy phenylboronic acid via a Suzuki reaction catalyzed with a trace amount of Pd(PPh₃)₄. The obtained acid, **1**, was then reacted with an equivalent of *N,N*-dimethyl ethylamine in the presence of 1-ethyl-3-(3-dimethylaminopropyl)carbodiimide hydrochloride (EDCI) and HOBT for 5 h at room temperature. This gave the IR-783-tertiary amine, **2**. Alkylation of **2** with bromo methylboronate, and potassium hydrogen fluoride workup gave the desired ammonium trifluoroborate (IR783-AMBF₃), **3**, which was isolated by preparative HPLC. [^{19}F]-**3** was stable at physiological pH. Less than 1% defluoridation of [^{19}F]-**3** was observed following 7 days of room temperature incubation with fetal bovine serum. The detailed synthesis is described in supporting information.

Cell culture

The murine gliosarcoma cell line (9 L/lacZ), the murine endothelioma cell line (bEnd.3), and the human microvascular endothelial cell line (HMEC-1) were purchased from ATCC and cultured in DMEM medium, supplemented with 10% (FBS, Seradigm, USA) and 100 U/mL penicillin/streptomycin (Gibco, USA), at 37 °C in a humidified incubator. Cell cytotoxicity assessment methodology is described in supporting information.

A note on IR783-AMBF₃ nomenclature due to the use of a PET-active multimodality imaging probe

In this manuscript, the notation “[^{18}F]-” on “[^{18}F]-IR783-AMBF₃” indicates that a radioactive mixture of fluoride-18 (<1%) and fluoride-19 (>99%) containing IR783-

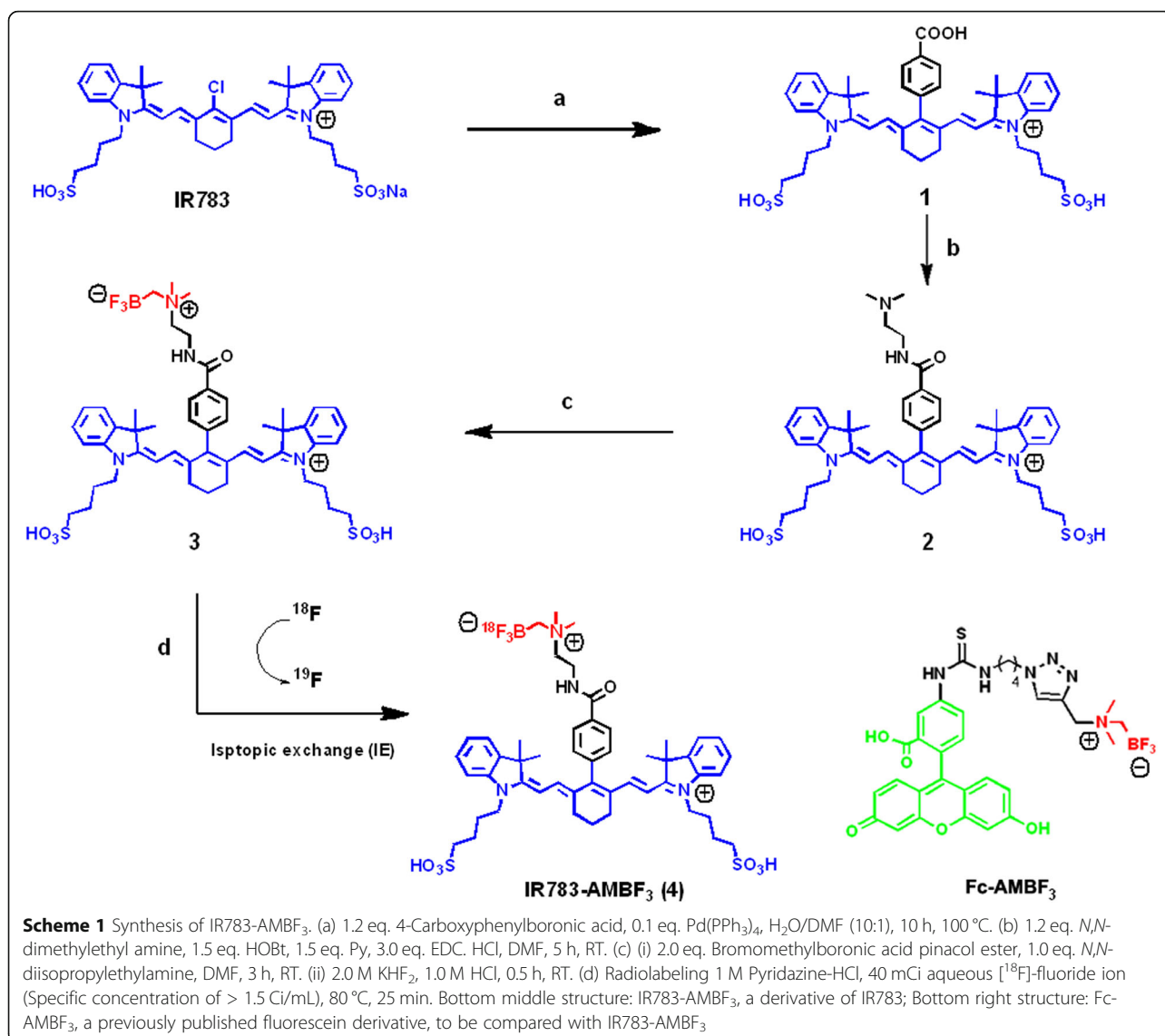

AMBF₃ is used. [¹⁸F]-IR783-AMBF₃ is visible on both fluorescent devices and PET/CT or PET/MRI devices. [¹⁸F]-IR783-AMBF₃ and IR783-AMBF₃ are electronically identical. The two molecules differ by only one neutron and the mixture is radioactive. The notation, IR783-AMBF₃, without the “[¹⁸F]-”, denotes a formulation containing only fluorine-19 atoms. This formulation is not radioactive, non-PET active, but is fluorescent. Non-radioactive IR783-AMBF₃ (as indicated by the lack of “[¹⁸F]”) is safer to handle and can be manipulated in facilities without radiation hazard protection standards.

In vitro fluorescent imaging

Red blood cells (RBCs) isolation was performed as described previously [23]. RBCs were collected from a Balb/c mouse and suspended in sterilized PBS-containing heparin (Product#27602, Fresenius Kabi

USA). Solutions of 1.5×10^9 RBCs were mixed with different concentrations of IR783-AMBF₃ (6.25–50 μM) or 25 μM of IR783-AMBF₃ and mixed with different RBC counts (1.5×10^6 to 1.5×10^9) in an Eppendorf tube. Fluorescent images were measured using an In-Vivo imager (Bruker, Billerica, MA, USA) at excitation/emission acquisition settings of 760/830 nm. Fc-AMBF₃ at the similar concentrations was used as a control. The excitation/emission acquisition settings for Fc-AMBF₃ was 450/535 nm.

CSF-rhinorrhea fluorescent imaging

Male Sprague-Dawley rats (6–7 weeks, 200–250 g) were purchased from Charles River Laboratories (Wilmington, MA, USA). Intrathecal catheterization was performed on rats by inserting a catheter (PE10 (0.26 mm I.D., 0.60 mm O.D), Instech Laboratories Inc., Plymouth Meeting,

PA) into the subarachnoid space of spinal cord between L5-L6, as described previously [20]. A 22-G needle was passed through the left nostril to pierce the cribriform plate to induce CSF leakage from the olfactory bulb into the paranasal sinus in spaces that also contained blood. Following CSF leak creation, 10 μM IR783-AMBF₃ (300 μL , 2.9 μg) was intrathecally injected into the rat, and CSF-rhinorrhea was collected on a gauze for fluorescent imaging. In control experiments, equimolar quantities of Fc-AMBF₃ were used as a control.

PET/MRI Scanning

A custom-built MRI compatible PET scanner with high-gain Silicon photomultiplier (SiPM) array detectors and clinical PET processing electronics (Cardinal, Siemens Molecular Imaging, Knoxville, TN) was used for the PET/MRI scan. This second generation PET/MRI scanner is considerably more thermal stable and less MRI sensitive than the first generation scanner [24]. The rat was anesthetized by 2% isoflurane and its respiration was monitored by a physiological monitoring system (SA Instruments, Stony Brook, NY). An intrathecal catheter line was placed before the animal was placed in the PET scanner. Two minutes after the start of PET data acquisition, 300 μL of 100 μCi [¹⁸F]-IR783-AMBF₃ (2.9 μg) was infused into the catheter with an auto-injector at 10 $\mu\text{L}/\text{min}$ for 30 min. ¹⁸F-PET data was continuously acquired for the same duration. After the PET scan, the animal was transported in the same cradle, in the same position, to the adjacent MRI scanner under anesthesia. Rat head and upper body MRI was performed on the 7-T MRI scanner using a rat brain coil for detection and a volume coil for RF excitation (Bruker Biospin Corp., Billerica, MA). After scout images were taken, 2D sagittal and coronal T2-weighted rapid acquisition with relaxation enhancement (RARE) images of the head and upper body were acquired. The FOV was 8 \times 5 cm with a data matrix of 384 \times 256, 20 1-mm slices, TR is 2.2 s and TE 50 ms, 8 averages, and a total acquisition of 10 min. PET data processing was achieved with home-written ordered subset expectation-maximization (OSEM) algorithm using the filtered back projection reconstruction method [25].

PET/CT Imaging

PET/CT scanning was performed on a Siemens Inveon PET/CT scanner (Siemens, Malvern, PA, USA). The rats ($n = 3$), bearing an intrathecal catheter (PE10), were anesthetized with 2.0–2.5% isoflurane. Then a 10 min CT/2 h PET scan was obtained, and 300 μL of 100 μCi [¹⁸F]-IR783-AMBF₃ was injected through catheter within 1 min during the PET phase. After scanning, the rats were sacrificed with CO₂ overdose followed by bilateral

thoracotomy. PET/CT was processed with Amide v1.0.4 and Inveon Research Workplace.

Statistical analysis

Values were shown as mean \pm SD for all experiments. Unpaired 2-tailed Student's *t* test was performed to assess statistical significance of the results and *P* values less than 0.05 were considered significant.

Results

Synthesis and Characterization of IR783-AMBF₃

We incorporated an alkylammoniummethyltrifluoroborate trap (AMBF₃) onto IR-783 [26, 27] to give an IR-783 derivative, IR783-AMBF₃. The synthesis of IR783-AMBF₃ is described in Scheme 1, and the characterizations of IR783-AMBF₃ by ¹H NMR, ¹³C NMR, ¹⁹F NMR, UPLC-MS, and HRMS are reported in supporting information. A mass of 11 mg was obtained starting with 100 mg of IR-783 in a 3-step chemical synthesis. The AMBF₃ and its linker did not alter IR-783's fluorescent properties, as the optical properties of IR783-AMBF₃ were similar to IR-783 (pH 7.4, 1 \times PBS, Supporting Table 1). At physiological pH (1 \times PBS), IR783-AMBF₃ had an extinction coefficient of 196,000 M⁻¹ cm⁻¹ ($\epsilon_{\text{max}} = 779$ nm) and an excitation maximum $\lambda = 778$ nm, with Stokes shift of 15 nm, while IR-783 had an excitation maximum $\lambda = 780$ nm, with Stokes shift of 22 nm (Supporting Table 1). IR783-AMBF₃ had a quantum yield of $\phi = 0.186$, similar to that of IR-783 ($\phi = 0.183$) (Supporting Table 1). The detailed fluorescence spectrum of IR783-AMBF₃ and IR-783 are shown in the supporting information.

Radiochemistry


The radiolabeling of IR783-AMBF₃ (3, Scheme 1) was performed in 30 min. In a typical synthesis, 9.2 mCi of [¹⁸F]-IR783-AMBF₃ (4) (molar activity of 185 mCi/ μmol , RCY = 23%, decay uncorrected) was obtained starting with 40.5 mCi of [¹⁸F]-fluoride ion (Scheme 1, d). Isotope exchange was undertaken in one step under aqueous, acidic pH conditions (pH = 2.0, pyridazine-HCl buffer, 10 μL) [28] and proceeded quickly (10–15 min) at high temperatures (90–100 °C). Unreacted [¹⁸F]-fluoride ion was removed by passing the [¹⁸F]-IR783-AMBF₃ (4) reaction mixture through a prewashed (2 mL ethanol followed by 5 mL deionized water) C18 cartridge (Waters no. 186005125). Unreacted [¹⁸F]-fluoride ion was flushed from the cartridge with 20 mL water. [¹⁸F]-IR783-AMBF₃, bound to the cartridge, was eluted with a 3.0 mM HCl in ethanol (300 L, 99%). All steps involving the cartridge were performed using a syringe pump set to deliver at a 40 mL/h flow rate. Resulting [¹⁸F]-IR783-AMBF₃ in acidic ethanol was immediately diluted 10- to 20-fold with 1 mM 1 \times PBS and filtered through a 0.22 μm filter. The reverse-phase HPLC of the resulting

pH 7.4 filtrate showed high purity of $[^{18}\text{F}]\text{-IR783-AMBF}_3$. The resulting filtrate was injected intrathecally.

PET and fluorescent visualization of IR783-AMBF₃ in the CSF-containing cisterns

To determine $[^{18}\text{F}]\text{-IR783-AMBF}_3$ utility in imaging CSF flow, we first evaluated $[^{18}\text{F}]\text{-IR783-AMBF}_3$ contrast distribution in the CNS following intrathecal administration (between the L5-L6 vertebrae) on PET and fluorescent imaging devices. Following L5-L6 intrathecal introduction, $[^{18}\text{F}]\text{-IR783-AMBF}_3$ presence was observed in the ambient cisterns, pineal recess, pituitary recess, and olfactory bulb both by PET/MRI (Fig. 1a) and PET/CT (Fig. 1b) imaging, which distributed the same way as a fluorescein-conjugated AMBF₃ (Fc-AMBF₃, Scheme 1)

[20, 21]. A dynamic video of intrathecal injection of $[^{18}\text{F}]\text{-IR783-AMBF}_3$ PET/MRI is included in Supporting Video 1. It should be noted that there was no contaminating free $[^{18}\text{F}]\text{-fluoride}$ ion after radiolabeling or before injection. Defluorination of IR783-AMBF₃ was not observed. $[^{18}\text{F}]\text{-fluoride}$ ion accumulation at the bone was not observed by PET/MRI or PET/CT. This lack of bone uptake demonstrates reliable synthesis, IR783-AMBF₃ radiolabeling, and molecular stability after intrathecal injection. Following PET/CT scanning, rats were sacrificed by CO₂ overdose and rat brains were collected and analyzed by ex vivo fluorescent imaging and high-magnification fluorescent histological analysis (Fig. 1c). Fluorescent signal was visible in the basal ambient cisterns, which was consistent with the PET images.

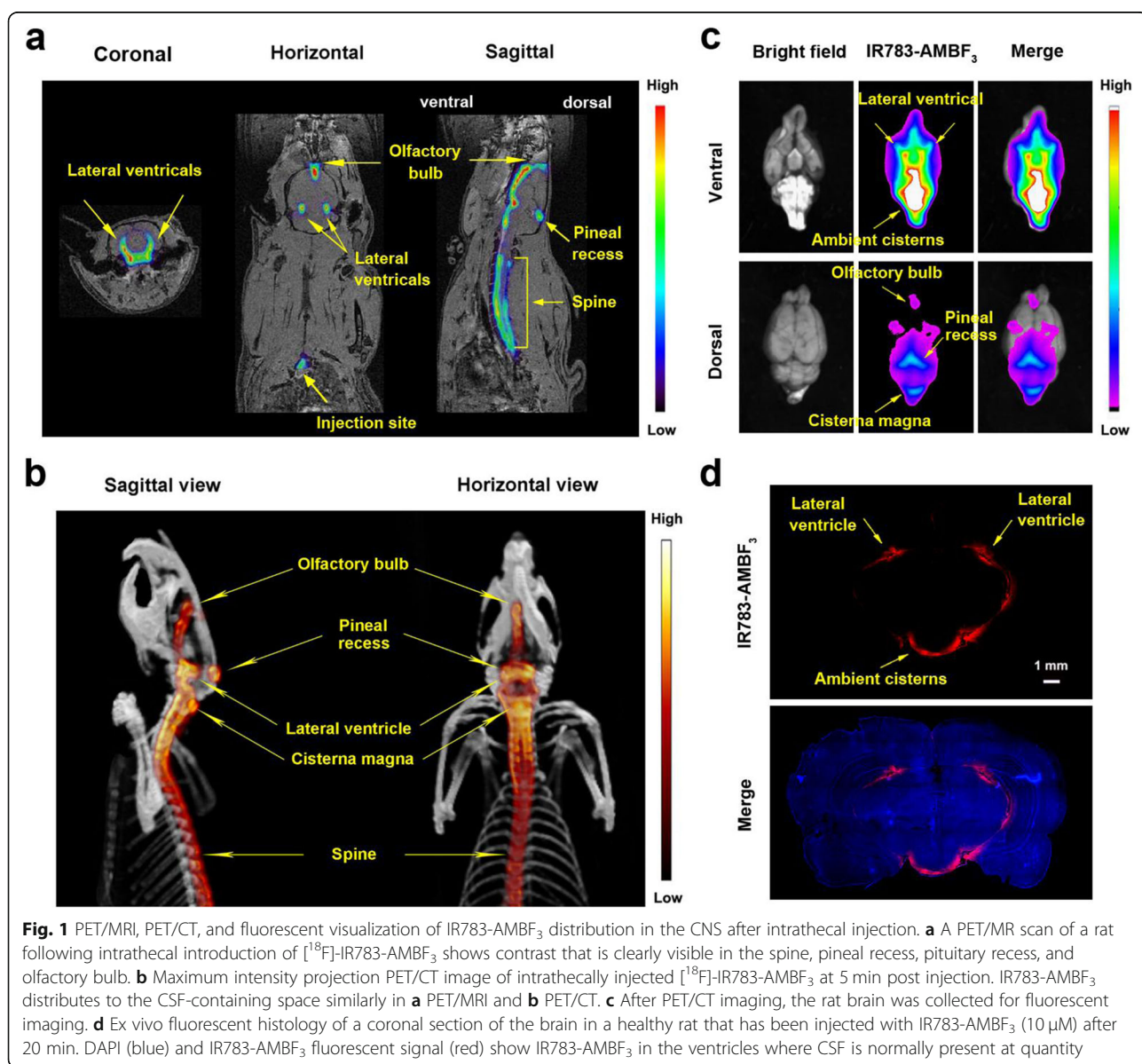

Notably, two lateral ventricles were observable through the unprocessed brain in the IR783-AMBF₃ fluorescent channel. In our previous publication using Fc-AMBF₃ [20, 21], these ventricles were difficult to visualize without sectioning the brain due to the strong absorption of exciting and emitted short wavelength photons by overlying brain tissue [2]. Fluorescent histology of a coronal section (- 19.11 mm, ant. ac) showed IR783-AMBF₃ in the ventricles where CSF is normally present (Fig. 1d). These data imply that newly designed IR783-AMBF₃ can serve as a multimodal contrast agent for delineating CSF-containing space by both PET and fluorescent imaging.

IR783-AMBF₃ was superior to Fc-AMBF₃ in diagnosing CSF leak when hemorrhage was present

In a prior study, we reported that Fc-AMBF₃, a fluorescein PET/Fluorescent agent, is useful for visualizing a paranasal-sinus CSF leak in rats by PET and fluorescence [20]. We concluded that Fc-AMBF₃ (excitation/emission = 495/511 nm) can be used to guide surgical CSF leak repair in endoscopic procedures in real time or to inform on the qualitative status of a CSF tear [29, 30].

However, CSF-bearing rhinorrhea collected from rats with a CSF leak was only fluorescently visible when the CSF was not contaminated with blood. A drawback to Fc-AMBF₃ is that the absorption of red blood cells (RBCs) can quench Fc-AMBF₃ fluorescence significantly (Fig. 2a, b, bottom) (Fig. S1b). Fortunately, IR783-AMBF₃ fluorescence was not quenched by RBCs (Fig. 2a, b, top), even when mixed with an increasing number of RBCs (Fig. S1a). To mimic the potential clinical ability of IR783-AMBF₃ in confirming CSF leak, in the presence of contaminating blood, an anterior skull base CSF leak in a rat with significant accompanying hemorrhage was created. A 22-G needle was inserted into left nostril of the deceased rats to induce the CSF leakage. As shown in Fig. 2c, CSF that was present in rhinorrhea collected from nostrils and month of rats bearing intrathecal IR783-AMBF₃ could not be visualized by the naked eye when it was mixed with blood. However, the CSF leak was able to be imaged and identified by IR783-AMBF₃ fluorescence even in the absence of massive bleeding (top). On the contrary, Fc-AMBF₃ fluorescence was not visible in a rat rhinorrhea containing significant quantities of blood (bottom). The above data demonstrate

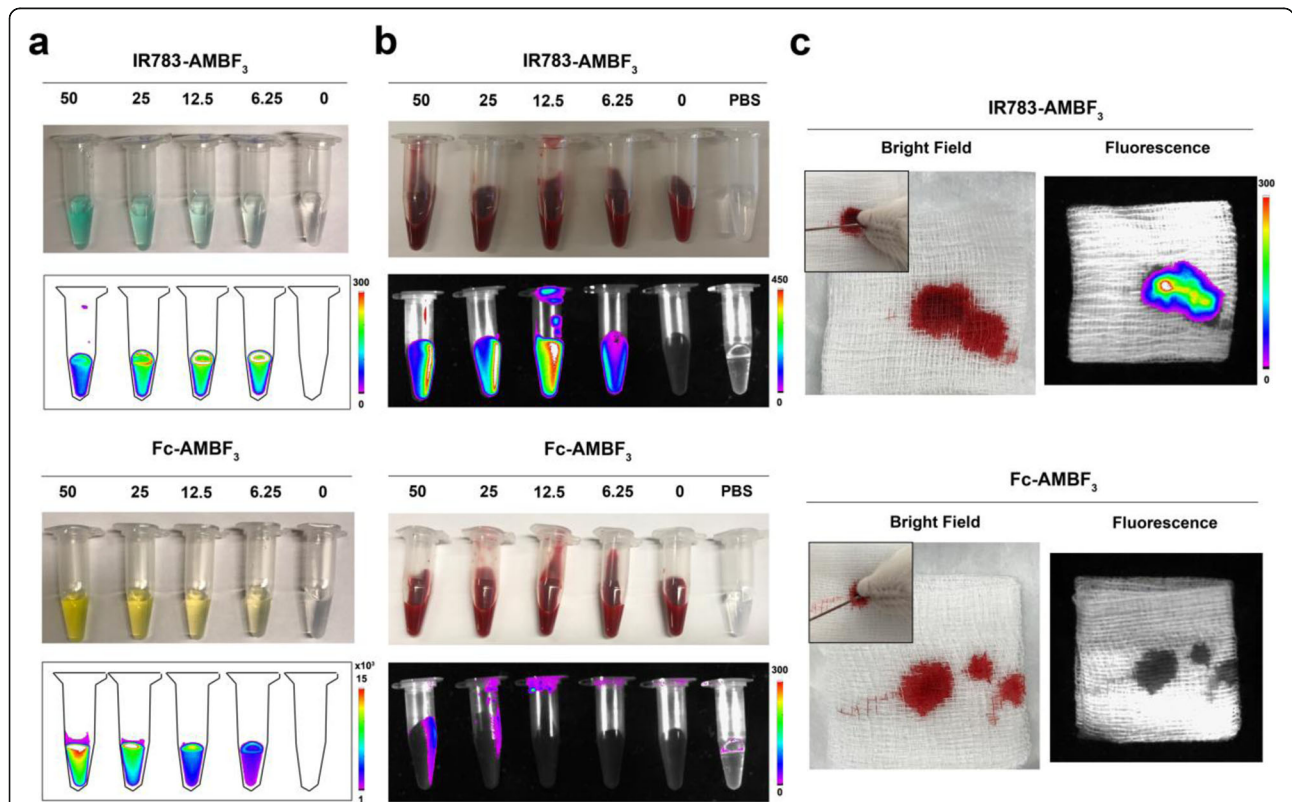


Fig. 2 RBC-mediated quenching of Fc-AMBF₃ but not IR783-AMBF₃ fluorescence. Fluorescence imaging of different concentrations of IR783-AMBF₃ or Fc-AMBF₃ in PBS, **a**, or when mixed with 1.5x10⁹ RBC/per tube, **b**. When mixed with blood, Fc-AMBF₃ fluorescence is not visible, while IR783-AMBF₃ fluorescence is clearly visible. Fluorescence imaging of blood contaminated CSF-containing rhinorrhea following the creation of a CSF leak using a 22G needle **c**. Excitation/emission filters for Fc-AMBF₃ and IR783-AMBF₃ were set at 450/535 nm and 760/830 nm, respectively. Note that in a CSF leak where RBC is not present, Fc-AMBF₃ is visible by fluorescence [20]

that IR783-AMBF₃ can be used in the identification of CSF leak in particularly traumatic situations: e.g. to assist surgeons in the rapid, accurate diagnosis and treatment of CSF lesions in the emergency setting.

IR783-AMBF₃ fluorescence was visible through skin in a rat bearing a lumboperitoneal shunt

To investigate IR783-AMBF₃ utility in observing CSF flow through a lumboperitoneal (LP) shunt, an LP shunt was placed in rats [21]. One end of the catheter (PE60) was inserted into the CSF space between L5 and L6 of a rat, and the outlet of this catheter was run subcutaneously along the abdomen, before being directed in the abdominal cavity (Fig. 3a). IR783-AMBF₃ was injected through the cisterna magna [21] and its flow was visualized in the fluorescent mode. IR783-AMBF₃ fluorescence was observable in vivo, through the skin (and some bone) and in the shunt, spine, and brain, demonstrating that IR783-AMBF₃ correctly distributed in the CNS via normal CSF flow (Fig. 3b). In control experiments performed with Fc-AMBF₃, no significant fluorescent signal in the shunt, spine, or brain could be observed through the skin (Fig. 3c). Most importantly, IR783-AMBF₃ clearance from CNS into the peritoneal cavity could be monitored through the skin by fluorescence

imaging (Fig. S2). The fluorescent intensity of IR783-AMBF₃ in the shunt decreased over time. The fluorescent intensity in the brain and peritoneal cavity first increased as IR783-AMBF₃ diffused into these structures with CSF flow, and then decreased as IR783-AMBF₃ cleared from the brain of rats (Fig. S3). This model is relevant in the clinical setting as ventriculoperitoneal shunts are run in patients, subdurally, from ventricles in the brain to the peritoneum. It was noted that the fluorescence signal in the kidney also increased after intrathecal administration of IR783-AMBF₃ and then decreased over time (Fig. S3). These fluorescence imaging data suggest that intrathecally injected [¹⁸F]-IR783-AMBF₃ clears from the CSF.

IR783-AMBF₃ was safe

The safety profile of IR783-AMBF₃ was evaluated in different cell lines and in healthy rats. IR783-AMBF₃ did not affect proliferation at concentrations as high as 50 μM following the incubation of IR783-AMBF₃ with cells for 24 h, suggesting that IR783-AMBF₃ does not induce cytotoxicity at concentrations under 50 μM (Fig. 4a). To evaluate the intrathecal toxicity of IR783-AMBF₃, rats (*n* = 4) were intrathecally injected with 45 nmol (15-fold higher than the imaging dose), and rat weight was monitored every

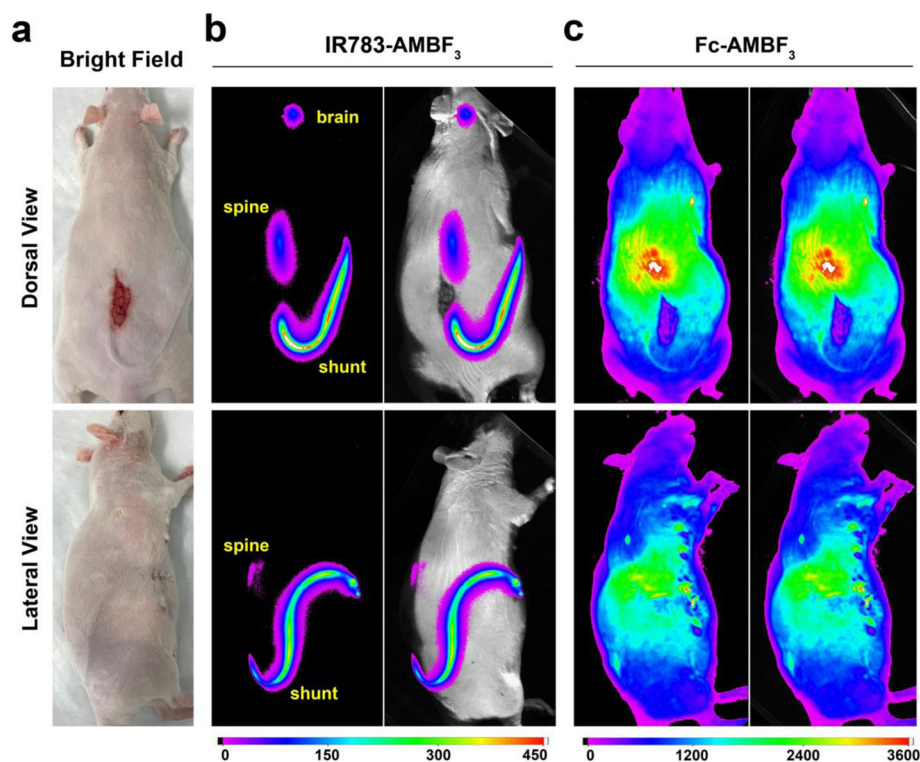


Fig. 3 Non-radioactive IR783-AMBF₃ can be visualized through deep tissue in the fluorescence mode. **a** Bright field images of rats (*n* = 2) with LP shunts that run superficially, under the skin of a shaved rat. **b** In vivo fluorescent images of the shunt in a rat after filling the shunt with IR783-AMBF₃ (10 μM). **c** (Control) Comparative in vivo fluorescent imaging of a shunt filled with Fc-AMBF₃ placed in a rat (1 mM, L5-L6). IR783-AMBF₃ fluorescence is visible in the shunt, spine, and brain, and through skin, bone, and a viable shunt

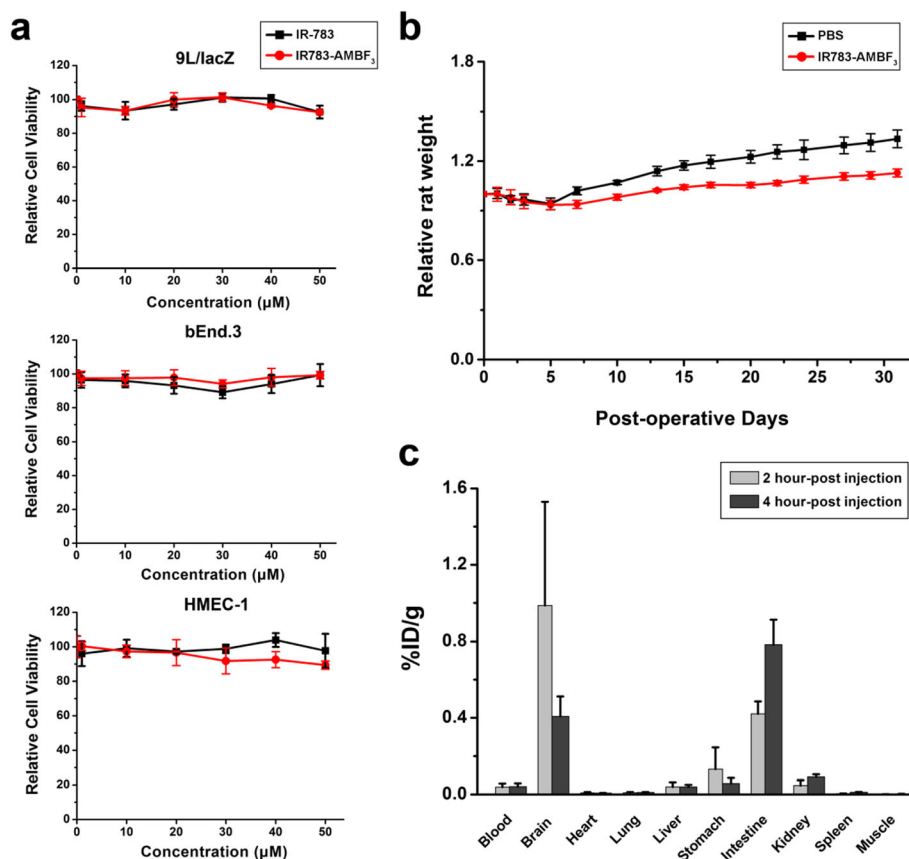


Fig. 4 Safety studies of IR783-AMBF₃. **a** Effect of IR783-AMBF₃ on murine gliosarcoma (9L/lacZ), murine endothelioma (bEnd.3) and human microvascular endothelial (HMEC-1) cell viability. Cells were treated with different concentrations of IR-783 (black line) or IR783-AMBF₃ (red line) for 24 h. **b** Relative mean weight of rats after intrathecal injection of IR783-AMBF₃ at 150 μM ($n = 4$). Error bars are \pm SD. **c** Biodistribution data collected in different tissues shows that intrathecally injected IR783-AMBF₃ clears from the CSF and blood quickly and is ultimately excreted and through the renal and hepatic systems ($n = 3$). Error bars are \pm SEM

other day for up to 1 month; normal recovery following intrathecal catheterization surgery (Fig. 4b) and post-surgical weight gain were observed in the rats. To confirm that IR783-AMBF₃ clears entirely from the CSF, rats were injected intrathecally with IR783-AMBF₃ and their organs were collected at 2 h and 4 h post-injection ($n = 3$) for γ -scintillated biodistribution data collection (Fig. 4c). Biodistribution data showed that IR783-AMBF₃ was cleared from the brain and spine, into the intestines and kidneys, suggesting that all intrathecally injected IR783-AMBF₃ was cleared from the CNS, and excreted through hepatic and renal routes.

Discussion

The direct radiochemical introduction of an ¹⁸F atom onto a cyanine dye using traditional carbon fluoride bond formation strategies requires non-trivial or multi-step approaches. These requirements demand automated synthesis, greater quantities of starting radioactivity, and lower chemical and radiochemical yields. Heptamethine cyanines are particularly unstable in basic conditions.

We describe a convenient synthesis of a PET-emitting heptamethine cyanine using aqueous, acidic ¹⁸F-isotopic exchange radiolabeling of IR783-AMBF₃ [31, 32]. IR-783 that is modified to bear a 4-carboxy phenylboronic acid provides a handle for further modification via Suzuki coupling (C–C bond formation). Isotopic exchange acidic radiolabeling gave the high-yield transfer of [¹⁸F]-fluoride ion onto [¹⁸F]-IR783-AMBF₃. The resulting radiolabeled, near-infrared wavelength fluorescent, dye-mixture contains [¹⁸F]-IR783-AMBF₃ for visualization by PET (both PET/CT and PET/MRI) and non-radioactive, non-PET visible IR783-AMBF₃ for fluorescence imaging (Fig. 1).

Rational design of IR783-AMBF₃

The use of AMBF₃ (dimethylammonium methyltrifluoroborate) ¹⁸F-isotopic exchange technology allows us to employ a stable, fluoride-19 containing, IR783-AMBF₃ radiochemical precursor to synthesize radioactive [¹⁸F]-IR783-AMBF₃ [32, 33]. Because of the AMBF₃ isotopic exchange technology employed, IR783-AMBF₃ is also

the stable radiochemical precursor to ^{18}F -IR783-AMBF₃, ^{18}F -IR783-AMBF₃ and IR783-AMBF₃ are electronically identical. i.e. No biological system can accurately distinguish between a molecule of ^{18}F -IR783-AMBF₃ and IR783-AMBF₃. This means that non-radioactive, fluoride-19 containing, IR783-AMBF₃ can substitute for ^{18}F -IR783-AMBF₃ in fluorescent-only studies or in toxicology studies (but not dosimetry studies). Additionally, IR783-AMBF₃ is hydrophilic and is low in molecular weight (MW < 1000). This hydrophilicity, conveyed by the zwitterionic AMBF₃, will prevent IR783-AMBF₃ from adhering to the meninges lining the CSF pathways, thus ensuring that IR783-AMBF₃ is cleared from the CSF and the body after a reasonable time, post-scan. Finally, IR783-AMBF₃ has a low molecular weight that resembles the molecular weights of both indocyanine green (ICG) and fluorescein, two current FDA approved agents for visualizing vascular flow in the fluorescent mode. In this study, when injected into the CSF, IR783-AMBF₃ allowed the CSF to be imaged in a rat model of CSF leak (Fig. 2) as well as allowed CSF flow to be imaged through a patent lumbo-peritoneal shunt (Fig. 3). In vitro or in vivo, IR783-AMBF₃ did not show any toxicity and was rapidly (< 2 h) cleared from the blood (Fig. 4).

Potential applications of IR783-AMBF₃

IR783-AMBF₃ is designed to roughly resemble the FDA approved agent indocyanine green (ICG), which is used in angiography and lymphangiography by ophthalmologists, interventional radiologists, and surgeons (sentinel node removal). IR783-AMBF₃ shares a molecular weight and a hydrophobicity that is similar to ICG. These properties were incorporated into the design of IR783-AMBF₃, in the hopes that IR783-AMBF₃ could serve as a substitute to ICG in future studies. For example, IR783-AMBF₃ potentially could allow for guided surgery on FDA-approved intraoperative fluorescent robotic devices [22], and to assist in the non-invasive through-skin evaluation of devices placed in the CNS [34]. In addition to being visible on the same fluorescent devices used to visualize ICG, IR783-AMBF₃ is also useful in PET scanning and will allow the molecularly-coherent-visualization of flow on PET/CT, PET/MRI and fluorescent devices. In CSF applications, drug and IR783-AMBF₃ co-delivery could be used to know exact CSF flow rates to allow for more accurate drug delivery in treating leptomeningeal disease [21], or to precisely localize injury sustained by the central nervous system (CSF leak).

Advantages of IR783-AMBF₃ vs. other modalities

CSF leaks can have a variety of appearances depending on their cause and rate leakage. Leak localization, characterization, and mitigation are very important to

aid both diagnosis and treatment. As a result, a number of imaging tests have been used to detect these leaks, each with its own strengths and weaknesses [35]. CT myelography and cisternography are the commonly used modalities in diagnosing CSF leaks [36]. MR cisternography has the advantage of being non-invasive and of not involving any ionizing radiation [35]. However, MR cisternography typically does not localize the exact leak site and suffers from artifacts. In addition, it can be difficult to distinguish CSF from sinusitis [37]. MR contrast myelography is superior to non-contrast MR; however, the introduction of high-dose gadolinium into the intrathecal space may bring up potential neurotoxicity in the CNS [18, 38, 39] and the long-term safety of intrathecal gadolinium is not conclusively known. Intrathecal gadolinium is not approved by the U.S. FDA currently. Another alternative approach for diagnosing CSF leak is ^{111}In -diethylene triamine pentaacetic acid (DTPA) intrathecal injection followed by single photon emission computed tomography (SPECT) imaging [40, 41]. The advantage of ^{111}In -SPECT is its prolonged monitoring capability—up to or greater than 72 h for the diagnosis of intermittent CSF leak [40]. However, this exam gives relatively low spatial resolution and sensitivity. On the contrary, IR783-AMBF₃ PET/MR is potentially a more sensitive and quantitative means of imaging CSF leak, since PET can produce higher resolution images than its ^{131}I , ^{111}In or $^{99\text{m}}\text{Tc}$ SPECT contrast counterparts [42]. Low-mass injections (< 1 mg) of IR783-AMBF₃ PET/CT are unlikely to have the same issues associated with the introduction of intrathecal gadolinium MR contrast in the CNS. IR783-AMBF₃ is superior to standalone PET due to its near-infrared fluorescence properties. IR783-AMBF₃ can immediately identify the precise anatomical location of CSF tears in patients with severe brain trauma in the emergency room using FDA-approved fluorescent imaging systems. There is no current universal approach for imaging CSF leaks most efficiently, and a single, superior imaging modality currently does not exist [1].

Advantages of near-infrared fluorescent IR783-AMBF₃ vs. fluorescent Fc-AMBF₃

We previously established the use of Fc-AMBF₃ to image CSF flow [21]. Like Fc-AMBF₃, IR783-AMBF₃ can be used to resolve CSF flow deep within the CNS in the PET mode. In this study, we observed three significant advantages of IR783-AMBF₃ over Fc-AMBF₃: (1) IR783-AMBF₃ can be observed in shunts, through the skin, and in subdural-run lumbo-peritoneal shunts, whereas Fc-AMBF₃ fluorescence cannot be observed through the skin of rats. Fc-AMBF₃ fluorescence is only visible in superficially run LP shunts (Fig. 3, Fig. S3). It is impractical to run a LP

shunt superficially in a patient, as LP shunts that are superficially run in patients introduce a significant likelihood of infection and shunt dislocation. (2) IR783-AMBF₃ fluorescence is visible when mixed with blood, while Fc-AMBF₃ fluorescence is not visible in the presence of blood (Fig. 2). In scheduled surgeries, the presence of blood is highly regulated and is therefore generally not an issue. However, IR783-AMBF₃ would be superior in an emergency room setting where a significant presence of blood can obscure trauma involving the CSF. (3) The greatest advantage of IR783-AMBF₃ is its current compatibility with currently FDA-approved intraoperative robotics which is already equipped to visualize IR783-AMBF₃.

Disadvantages of IR783-AMBF₃ vs. Fc-AMBF₃

IR783-AMBF₃ (excitation/emission = 778/793 nm) requires more complicated instrumentation for its fluorescent observation. The wavelength of photon emitted by IR783-AMBF₃ is not visible by the naked eye; therefore, a CCD/CMOS camera fitted with an optical surface is required to collect and report on the emitted, fluorescent photon. Additionally, a surgeon can only visualize IR783-AMBF₃ on a digital screen. Unlike IR783-AMBF₃, the fluorescent emission of Fc-AMBF₃ (excitation/emission = 495/511 nm) is visible with inexpensive black-light illumination, by the naked eye, and without an emission filter; therefore, surgeons can observe Fc-AMBF₃ directly in the operating field within a patient without costly equipment. Another drawback to IR783-AMBF₃ is that the safety of IR783-AMBF₃ in the CSF is unknown. Fc-AMBF₃ is attractive for clinical trials as it is molecularly similar to fluorescein. Fluorescein is currently used and is known to be safe at < 50 mg doses for the purposes of CSF leak repair. The safety of IR783-AMBF₃ is less understood; however, the noted advantages of IR783-AMBF₃ may make its utility in CSF imaging more attractive vs. Fc-AMBF₃.

Conclusion

A new synthesis and ¹⁸F-radiolabeling of a dual PET/near-infrared fluorescence probe, IR783-AMBF₃, is reported. Intrathecal injections of IR783-AMBF₃ are tolerated in rats under 45 nmol dosage, and it can rapidly clear from the body within 2 h. IR783-AMBF₃ can be used to non-invasively and accurately diagnose and monitor CSFL and CSF flow through lumboperitoneal shunts by both PET imaging and near-infrared fluorescence imaging. The ability to image CSF flow in deep tissue and in the presence of blood allows IR783-AMBF₃ to serve as a powerful dual-modality agent for diagnosing both CSF-related diseases and imaging vascular flow.

Supplementary information

Supplementary information accompanies this paper at <https://doi.org/10.1186/s13550-020-0609-3>.

Additional file 1. A dynamic [¹⁸F]-IR783-AMBF₃ PET/MRI of a rat showing IR783-AMBF₃-CSF flow into the brain (0–30 min) by a 300 μL, 10 μL/min [¹⁸F]-IR783-AMBF₃ lumbar (L5–L6) infusion.

Additional file 2. Synthesis and chemical characterization of IR783-AMBF₃, supporting materials, methods, and supporting figures.

Abbreviations

9L/lacZ: Murine gliosarcoma cell; ACN: Acetonitrile; AMBF₃: alkylammoniumethyltrifluoroborate; bEnd.3: murine endothelioma cell; CNS: Central nervous system; CSF: Cerebrospinal fluid; CSFL: Cerebrospinal fluid leak; CT: Computed tomography; DMF: Dimethyl formamide; DMSO: Dimethyl sulfoxide; DTPA: Diethylene triamine pentaacetic acid; EDCI: 1-ethyl-3-(3-dimethylaminopropyl)carbodiimide hydrochloride; eq: Equivalents; EtOH: Ethanol; HCl: Hydrochloric acid; HMEC-1: human microvascular endothelial cell; HOBt: Hydroxybenzotriazole; HPLC: High performance liquid chromatography; ICG: Indocyanine green; ID: Inside diameter; KHF₂: Potassium hydrogen fluoride; L5–L6: The space between fifth and sixth lumbar vertebrae; LP: Lumboperitoneal; MeOH: Methanol; MRI: Magnetic resonance imaging; MS: Mass spectrometer; NMR: Nuclear magnetic resonance; OD: Outside diameter; PBS: Phosphate-buffered saline; PE: Polyethylene; PET: Positron emission tomography; RBCs: Red blood cells; SPECT: Single photon emission computed tomography; TFA: Trifluoroacetic acid

Acknowledgements

Not applicable.

Authors' contributions

HG and HK contributed equally. RT and ZZ conceived and designed the experiments. HG conducted most of the in vitro and in vivo experiments. HK designed and synthesized the compound. CL, JK, MJ, and SR performed the PET/MRI experiments and analysis. AW, OA, MS, and OA designed the animal experiments and discussed part sections of the manuscript with corresponding authors. RT, ZZ, HG, and HK analyzed the experimental data and wrote the manuscript with input from all co-authors. All authors read and approved the final manuscript.

Funding

This work was funded by the National Institute of Biomedical Imaging and Bioengineering K99/R00 (EB013904), the NIH/NCI Cancer Center Support Grant P30 (CA008748), the US Department of Defense Congressionally Directed Research Projects Special Emphasis Grant on Pediatric Brain Tumors (CA160373), and an NIH Grant RC2 (CA148971).

Availability of data and materials

Please contact the authors for data or material requests.

Ethics approval and consent to participate

All procedures conducted in studies were approved by the Weill Cornell Medical Center Institutional Animal Care and Use Committee (no. 2014-0030) and were consistent with the recommendations of the American Veterinary Medical Association and the National Institutes of Health Guide for the Care and Use of Laboratory Animals, and complied with the ARRIVE guidelines.

Consent for publication

Not applicable.

Competing interests

The authors declare that they have no competing interests.

Author details

¹Department of Nuclear Medicine, Chinese Academy of Medical Sciences and Peking Union Medical College Hospital, Beijing 100730, China. ²Beijing Key Laboratory of Molecular Targeted Diagnosis and Therapy in Nuclear Medicine, Beijing 100730, China. ³Department of Radiology, Molecular Imaging Innovations Institute (MI3), Weill Cornell Medical College, New York,

NY 10065, USA. ⁴Department of Medical Physics, Memorial Sloan Kettering Cancer Center, New York, NY 10065, USA. ⁵Department of Biomedical Engineering, University of California at Davis, Davis, CA 95616, USA. ⁶Department of Otolaryngology – Head & Neck Surgery, Northwell Health, Hofstra Northwell School of Medicine, New York, NY 10075, USA. ⁷Department of Radiology, Memorial Sloan Kettering Cancer Center, New York, NY 10065, USA. ⁸Department of Neurological Surgery, Weill Cornell Medicine, New York, NY 10065, USA.

Received: 28 November 2019 Accepted: 14 February 2020

Published online: 16 April 2020

References

- Tsien RY. Imaging imaging's future. *Nat Rev Mol Cell Biol.* 2003;Suppl:SS16-21.
- An FF, Chan M, Kommidhi H, Ting R. Dual PET and near-infrared fluorescence imaging probes as tools for imaging in oncology. *AJR Am J Roentgenol.* 2016;207:266–73. <https://doi.org/10.2214/AJR.16.16181>.
- Disselhorst JA, Brom M, Laverman P, Slump CH, Boerman OC, Oyen WJ, et al. Image-quality assessment for several positron emitters using the NEMA NU 4-2008 standards in the Siemens Inveon small-animal PET scanner. *J Nucl Med.* 2010;51:610–7. <https://doi.org/10.2967/jnumed.109.068858>.
- Nguyen QT, Tsien RY. Fluorescence-guided surgery with live molecular navigation—a new cutting edge. *Nat Rev Cancer.* 2013;13:653–62. <https://doi.org/10.1038/nrc3566>.
- Paulus A, Desai P, Carney B, Carlucci G, Reiner T, Brand C, et al. Development of a clickable bimodal fluorescent/PET probe for in vivo imaging. *EJNMMI Res.* 2015;5:120. <https://doi.org/10.1186/s13550-015-0120-4>.
- Rosenthal EL, Warram JM, De Boer E, Chung TK, Korb ML, Brandwein-Gensler M, et al. Safety and tumor specificity of cetuximab-IRDye800 for surgical navigation in head and neck cancer. *Clin Cancer Res.* 2015;21:3658–66.
- Cheng L, Kamkaew A, Sun H, Jiang D, Valdovinos HF, Gong H, et al. Dual-modality positron emission tomography/optical image-guided photodynamic cancer therapy with chlorin e6-containing nanomicelles. *ACS Nano.* 2016;10:7721–30.
- Guo H, Hari Krishna K, Vedvyas Y, McCloskey JE, Zhang W, Chen N, et al. A fluorescent,[18F]-positron-emitting agent for imaging PMSA allows genetic reporting in adoptively-transferred, genetically-modified cells. *ACS Chem Biol.* 2019.
- Carlucci G, Carney B, Brand C, Kossatz S, Irwin CP, Carlin SD, et al. Dual-modality optical/PET imaging of PARP1 in glioblastoma. *Mol Imaging Biol.* 2015;17:848–55. <https://doi.org/10.1007/s11307-015-0858-0>.
- Khazim R, Dannawi Z, Spacey K, Khazim M, Lennon S, Reda A, et al. Incidence and treatment of delayed symptoms of CSF leak following lumbar spinal surgery. *Eur Spine J.* 2015;24:2069–76.
- Schievink WI, Maya MM, Chu RM, Moser FG. False localizing sign of cervico-thoracic CSF leak in spontaneous intracranial hypotension. *Neurology.* 2015; 84:2445–8.
- Jiang D, Dalong N, Li S, Ehlerding E, Sun T, Cao T, et al. Framework DNA nanocages alleviate ischemic stroke via intrathecal injection. *J Nucl Med.* 2019;60:335.
- Iliff JJ, Lee H, Yu M, Feng T, Logan J, Nedergaard M, et al. Brain-wide pathway for waste clearance captured by contrast-enhanced MRI. *J Clin Invest.* 2013;123:1299–309. <https://doi.org/10.1172/JCI67677>.
- Kwon S, Janssen CF, Velasquez FC, Sevic-Muraca EM. Fluorescence imaging of lymphatic outflow of cerebrospinal fluid in mice. *J Immunol Methods.* 2017;449:37–43.
- Strazielle N, Ghersi-Egea JF. Physiology of blood-brain interfaces in relation to brain disposition of small compounds and macromolecules. *Mol Pharm.* 2013;10:1473–91. <https://doi.org/10.1021/mp300518e>.
- Wolf DA, Hesterman JY, Sullivan JM, Orcutt KD, Silva MD, Lobo M, et al. Dynamic dual-isotope molecular imaging elucidates principles for optimizing intrathecal drug delivery. *JCI insight.* 2016;1.
- George N, Gean EG, Nandi A, Frolov B, Zaidi E, Lee H, et al. Advances in CNS imaging agents: focus on PET and SPECT tracers in experimental and clinical use. *CNS Drugs.* 2015;29:313–30.
- Sese D, Blaha T, Khurana R, Wachsmann A, Al-ali F, Murray T, et al. A leaky situation: gadolinium contrast induced neurotoxicity from intrathecal contrast. D45 CRITICAL CARE CASE REPORTS: I (DON'T) WANT TO BE SEDATED-NEUROCRITICAL CARE, SEDATION, AND DELIRIUM: American Thoracic Society; 2018. p. A6888-A.
- Starke RM, Dumont AS. Intraoperative imaging and assessment of cerebral blood flow in cerebrovascular surgery: hybrid operating rooms, intraoperative angiography and magnetic resonance imaging, Doppler ultrasound, cerebral blood flow probes, endoscopic assistance, indocyanine green videography, and laser speckle contrast imaging. *World Neurosurg.* 2014;82:e693–6. <https://doi.org/10.1016/j.wneu.2013.10.045>.
- Kommidhi H, Guo H, Chen N, Kim D, He B, Wu AP, et al. An [18F]-positron-emitting, fluorescent, cerebrospinal fluid probe for imaging damage to the brain and spine. *Theranostics.* 2017;7:2377–91. <https://doi.org/10.7150/thno.19408>.
- Guo H, Kommidhi H, Maachani UB, Voronina JC, Zhang W, Magge RS, et al. An [(18)F]-positron emitting fluorophore allows safe evaluation of small molecule distribution in the CSF, CSF fistulas, and CNS device placement. *Mol Pharm.* 2019;16:3636–46. <https://doi.org/10.1021/acs.molpharmaceut.9b00485>.
- Sajedi S, Sabet H, Choi HS. Intraoperative biophotonic imaging systems for image-guided interventions. *Nanophotonics.* 2018;8:99–116.
- Wang Y, An F-F, Chan M, Friedman B, Rodriguez EA, Tsien RY, et al. 18F-positron-emitting/fluorescent labeled erythrocytes allow imaging of internal hemorrhage in a murine intracranial hemorrhage model. *J Cereb Blood Flow Metab.* 2017;0271678X16682510. doi:10.1177/0271678X16682510.
- Catana C, Wu Y, Judenhofer MS, Qi J, Pichler BJ, Cherry SR. Simultaneous acquisition of multislice PET and MR images: initial results with a MR-compatible PET scanner. *J Nucl Med.* 2006;47:1968–76.
- Wu Y, Catana C, Farrell R, Dokhale PA, Shah KS, Qi J, et al. PET performance evaluation of an MR-compatible PET insert. *IEEE Trans Nucl Sci.* 2009;56:574–80. <https://doi.org/10.1109/TNS.2009.2015448>.
- Liu Z, Lin KS, Benard F, Pourghiasian M, Kiesewetter DO, Perrin DM, et al. One-step (18)F labeling of biomolecules using organotrifluoroborates. *Nat Protoc.* 2015;10:1423–32. <https://doi.org/10.1038/nprot.2015.090>.
- Liu Z, Pourghiasian M, Radtke MA, Lau J, Pan J, Dias GM, et al. An organotrifluoroborate for broadly applicable one-step F-labeling. *Angew Chem Int Ed Engl.* 2014. <https://doi.org/10.1002/anie.201406258>.
- Liu Z, Li Y, Lozada J, Wong MQ, Greene J, Lin K-S, et al. Kit-like 18F-labeling of RGD-19F-Aryltrifluoroborate in high yield and at extraordinarily high specific activity with preliminary in vivo tumor imaging. *Nucl Med Biol.* 2013;40:841–9. <https://doi.org/10.1016/j.nucmedbio.2013.05.002>.
- Seth R, Rajasekaran K, Benninger MS, Batra PS. The utility of intrathecal fluorescein in cerebrospinal fluid leak repair. *Otolaryngol Head Neck Surg.* 2010;143:626–32. <https://doi.org/10.1016/j.otohns.2010.07.011>.
- Banu MA, Kim J-H, Shin BJ, Woodworth GF, Anand VK, Schwartz TH. Low-dose intrathecal fluorescein and etiology-based graft choice in endoscopic endonasal closure of CSF leaks. *Clin Neurol Neurosurg.* 2014;116:28–34.
- Bernard-Gauthier V, Bailey JJ, Liu Z, Wangler B, Wangler C, Jurkschat K, et al. From unorthodox to established: the current status of (18)F-trifluoroborate- and (18)F-SiFA-based radiopharmaceuticals in PET nuclear imaging. *Bioconjug Chem.* 2016;27:267–79. <https://doi.org/10.1021/acs.bioconjugchem.5b00560>.
- Perrin DM. [18F]-organotrifluoroborates as radioprotective groups for PET Imaging: from design principles to preclinical applications. *Acc Chem Res.* 2016. <https://doi.org/10.1021/acs.accounts.5b00398>.
- Schirmacher E, Wängler B, Cypriak M, Bradtmöller G, Schäfer M, Eisenhut M, et al. Synthesis of p-(di-tert-butyl [18F] fluorosilyl) benzaldehyde ([18F] SiFA-A) with high specific activity by isotopic exchange: a convenient labeling synthon for the 18F-labeling of N-amino-oxy derivatized peptides. *Bioconjug Chem.* 2007;18:2085–9.
- Yadav YR, Parihar V, Sinha M. Lumbar peritoneal shunt. *Neurol India.* 2010; 58:179–84. <https://doi.org/10.4103/0028-3886.63778>.
- Kranz PG, Luetmer PH, Diehn FE, Amrhein TJ, Tanpitukpongse TP, Gray L. Myelographic techniques for the detection of spinal CSF leaks in spontaneous intracranial hypotension. *Am J Roentgenol.* 2016;206:8–19.
- Kranz P, Gray L, Taylor J. CT-guided epidural blood patching of directly observed or potential leak sites for the targeted treatment of spontaneous intracranial hypotension. *Am J Neuroradiol.* 2011;32:832–8.
- Melamed J, Stone B, Beaucher WN. Spontaneous cerebrospinal fluid rhinorrhea with associated sinusitis and allergic rhinitis. *Allergy Proc.* 1994; 15:197–200.
- Park K-W, Im S-B, Kim B-T, Hwang S-C, Park J-S, Shin W-H. Neurotoxic manifestations of an overdose intrathecal injection of gadopentetate dimeglumine. *J Korean Med Sci.* 2010;25:505–8.

39. Feng X, Xia Q, Yuan L, Yang X, Wang K. Impaired mitochondrial function and oxidative stress in rat cortical neurons: implications for gadolinium-induced neurotoxicity. *Neurotoxicology*. 2010;31:391–8.
40. Mokri B. Radioisotope cisternography in spontaneous CSF leaks: interpretations and misinterpretations. *Headache: The Journal of Head and Face Pain*. 2014;54:1358–68.
41. Howard BA, Gray L, Isaacs RE, Borges-Neto S. Definitive diagnosis of cerebrospinal fluid leak into the pleural space using 111In-DTPA cisternography. *Clin Nucl Med*. 2015;40:220–3.
42. Rahmim A, Zaidi H. PET versus SPECT: strengths, limitations and challenges. *Nucl Med Commun*. 2008;29:193–207. <https://doi.org/10.1097/MNM.0b013e3282f3a515>.

Publisher's Note

Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Submit your manuscript to a SpringerOpen[®] journal and benefit from:

- ▶ Convenient online submission
- ▶ Rigorous peer review
- ▶ Open access: articles freely available online
- ▶ High visibility within the field
- ▶ Retaining the copyright to your article

Submit your next manuscript at ▶ [springeropen.com](https://www.springeropen.com)
