

Publication Year	2019
Acceptance in OA@INAF	2022-07-15T09:40:26Z
Title	ESPRESSO: the new VLT workhorse for high-precision astronomy
Authors	DI MARCANTONIO, Paolo
Handle	http://hdl.handle.net/20.500.12386/32507






ESPRESSO: the new VLT workhorse for high-precision astronomy

P. Di Marcantonio

on behalf of

ESPRESSO Team

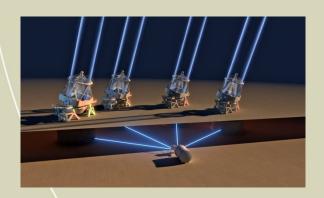





Outline of the talk





- the ESPRESSO (Echelle SPectrograph for Rocky Exoplanets and Stable Spectral Observations) instrument
- ESPRESSO at the VLT
- ESPRESSO Data Flow and its final deployment
- ESPRESSO operations
- ESPRESSO "science" software in a nutshell

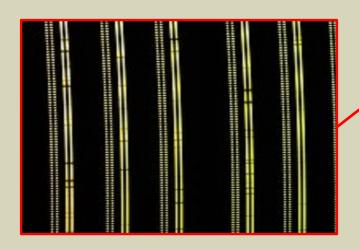

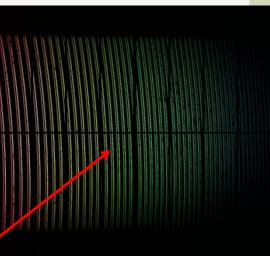



The ESPRESSO "facility"









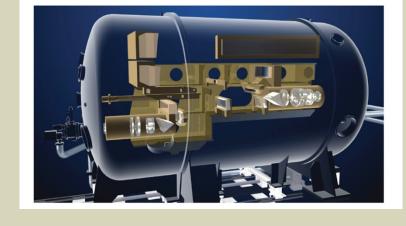

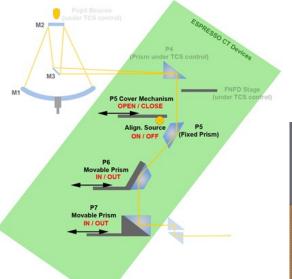

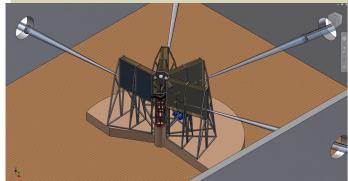

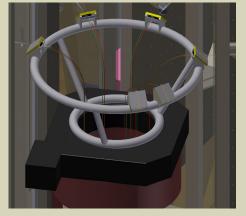

Galactic Archaelogy in the Gaia Era

Instrument characteristics





Parameter	singleUHR	singleHR	multiMR
Wavelengths		Blue arm: 380 – 520 n Red arm: 520 – 780 n	
Spectral coverage		Full	
Spectra format	Echelle, up to 4	spectra per order (2 fibe	ers, 2 spectra / fiber)
Resolving power	225'000	134'000	59'000
Aperture on sky	0.5 arcsec	1.0 arcsec	4x1.0 arcsec
Spectral sampling	2.5 pixels	4.5 pixels	10 pixels
Spatial sampling	9 pixels	18 pixels	44 pixels
Available binning	1x1	1x1 or 2x1	4x2 or 8x4
Sky/Simultaneous reference		Yes (mutually exclusive	e)
Instrumental RV precision	<10 cm/s	<10 cm/s	~1 m/s


ESPRESSO science








measurement of highprecision RVs of G, K and M dwarfs to search for rocky planets inside the HZ

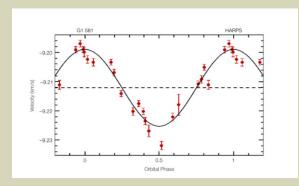

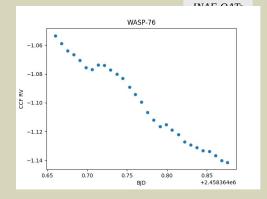

investigate possibly variability of fundamental constants



> study chemical composition of stars in local galaxies, IGM etc.

the ESPRESSO instrument (FDR movie)

Optics AF-OATs Fiber 100 um Equivalent fiber 100x300 um Pupil slicer Pupil 20mm Pupil 60x20 mm Equ. fiber (2X) 200x600 um FIBER EXIT 3X anamorphic collimator Fiber image (2X) 60x120 um Equ. fiber (2X) 200x400 um Pupil 300x200 mm (on echelle grating) Pupil 120x120 mm (on cross-disperser) DETECTOR 1.5X anamorphic cross disperser Figure 78. Paraxial image and pupil sizes along the spectrograph optics 905 nm 141 434 nm 151 405 nm

Figure 107. Blue spectral format. The box on the left represents the CCD area.


380 nm

Scrambler








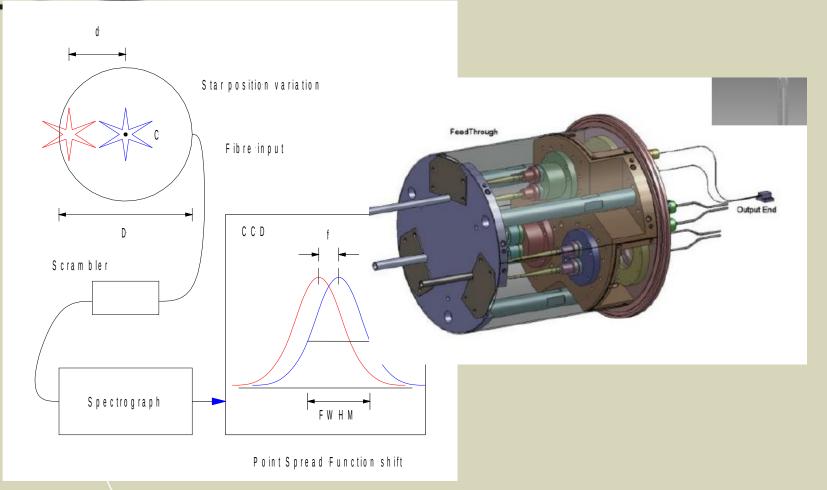

ESPRESSO CCDs





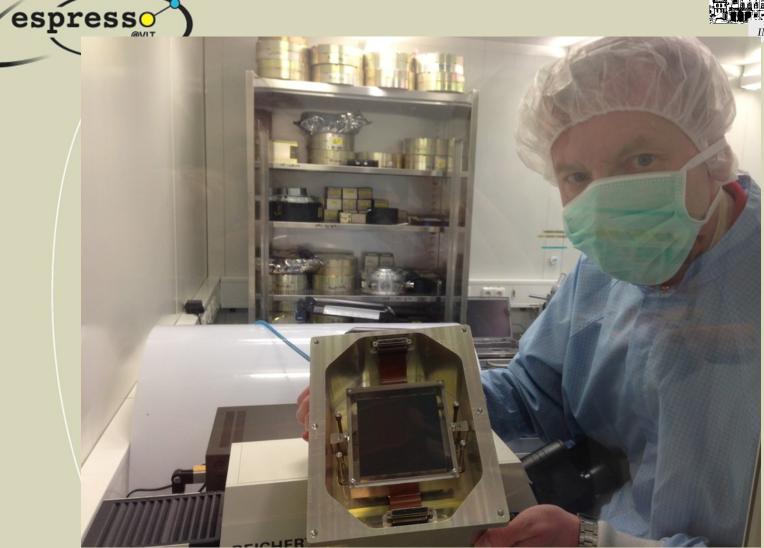

Laser Frequency Comb








Astro-comb: ~450 lines per order

Th-Ar: ~ 150 lines per order





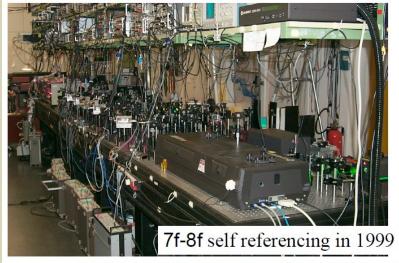

espr Astro-comb





The consortium





Institutes of four countries (1 PI or Co-PI in each)

- Switzerland: Univ. Geneva (Lead, F. Pepe), Univ. Bern
- Italy: INAF Trieste (S. Cristiani, P. Molaro, P. Di Marcantonio), INAF Brera
- Portugal: Univ. Porto (N. Santos), Univ. Lisbon
- Spain: IAC (R. Rebolo Lopez)

Associated partner (Representative in Executive Board)

ESO (H. Dekker)

ESPRESSO @ VLT





ESPRESSO milestones:

- ✓ PDR: end of November 2011
- ✓ FDR: May 2013
 - ✓ CT First light: 25 Sept 2016 on UT4
- ✓ PAE: Q2 Q3 2017, shipment on Aug 2017
 - ✓ ESPRESSO single-UT first light: 27 Nov 2017 on UT1
 - ✓ ESPRESSO 4-UT first light: 3 Feb 2018
- Delivery to community: Oct 2018 (single-UT mode only)
- Start of GTO: September 2018

AIT @ Paranal









AIT @ Paranal



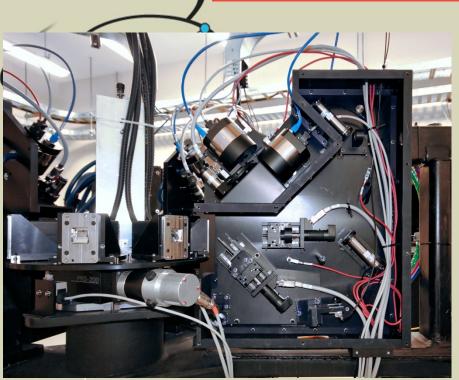



First light(s)

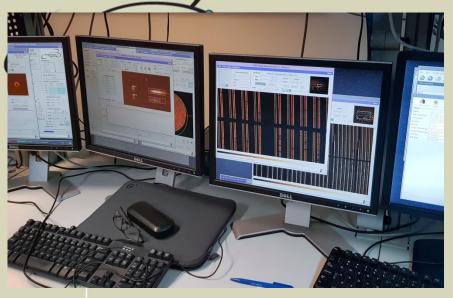






ESPRESSO Data Flow System







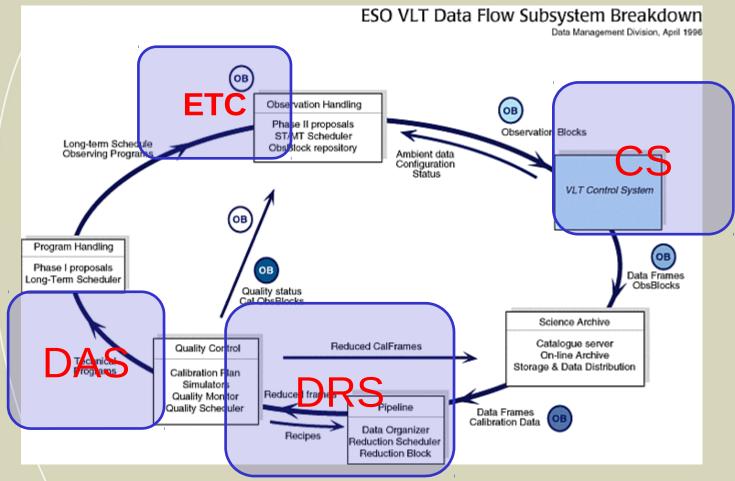


ESPRESSO Data Flow deployment





Data Flow model for Espresso (OHS-DHS)

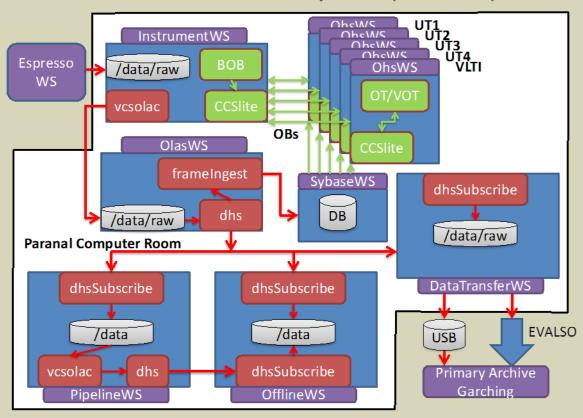

(courtesy C. Guirao)

Exposure Time Calculator






+ <u>E</u> \$+ 0 +							INAF-OATS
					Optical Echelle	elle Spectroscopy Mode Version P103.4	IIVII -O/IIS
Farget Input Fl	lux Distrib	ıtion				IQ: arcsec FWHM at the airmass and wavelength of observation (to be used for the OB constraint section).	et)
Template Spectrum	A0V (Pickles) ~					
O MARCS Stellar Mo	del Teff=4000 io	g(g)=-0.5 [Fe/H]= (0 M= 1	Redshift z = 0.00	Target Ma	Magnitude T	
O Upload Spectrum	Select				V V = 7	Instrument Setup	
O Blackbody	Temperature	: K			Magnitudes	des are given	
O Power Law	Index:	FΩ) ochinden			Telescope feed: 1 UT V	
O Emission Line	Lambda: Flux: FWHM:	nm 10 ⁻¹⁶ er nm	gs/s/cm²	(per arcsec² for	extended source	Resolution/detector mode: singleHR11 Single UT, high spectral/RV resolution, 1x1 binning and fast readout Exposure Time: 30.000 s	
Spatial Distribution:	Point Source						
Sky Conditions Override almanac Moon FLI:			pical fixe	d sky model par	ameters except l	Results pt Moon phi	
Almanac						Tables: ☑ Toggie All / No Tables Graphs: ☐ Toggie All / No Graphs	
						☑ Spectral Format ☐ Input Spectrum	
U 1 08/01/2019 12:50:54 Time step forward/back +							
Paranal Local C	1.53951389 Civil Time LCT: idereal Time LST:	= 09:56:54	Day!		-	□ Obj	
Paranal Local S	idereal Time LST :	= 15:26:28	Day.			□ Sky	
Target coordin	nates Name or ID:	Aldebaran / HIP20	109 / M33	search		□ Maximum Intensity	
α 04:00:00	δ 00:0		toggle unit			" S/N	
UT rise culm. Hour Ar	ngle azimuth altitu	de zenith distance	airmass	FLI moon/sur phase ang	n moon/target le separation		
Target 01:30 11:26: 07:31	28 199°.5 -64°					38-	Submit Reset
11:08 Moon 18:02	105°.0 9°.	target too low		0.06 "dark"	100°.7		
01:02 09:54 Sun 16:44				27°.9		10	
Sun 16:44 23:42	100°.0 37° day					https://www.eso.org/observing/etc/b	in/gen/form2INS N
An adva	anced almanac is ava	ilable in the ESO SI	kyCalc sky	model calculator		AME=ESPRESSO+INS.MODE=spe	octro
Seeing/Image Quality:						AIVIE-ESPRESSOTINS.IVIODE-SPE	<u> </u>
Seeing: 1.00	arcsec FWHM in	ı V-band at zeniti	h (use thi	is value in the pr	oposal)		
Probabilit	ty 87% of realisin	g the seeing ≤ 1 c	arcsec				

Operations





INAF-OATs

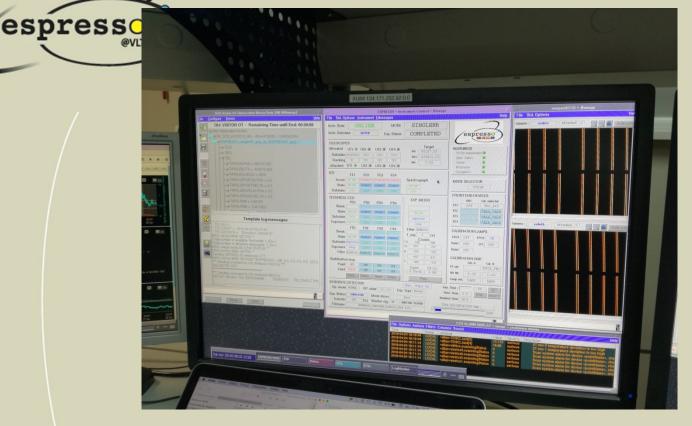

During the observing run, the Observation Block (OB) to be executed is loaded on the Observation Handling (OH) workstation and then sent to the BOB (Broker for Observation Block) tool on the IWS. BOB reads the contents of the OB and executes one by one the templates specified in there. Each template consists in general of a sequence of commands determining completely a scientific exposure. It allows presetting the telescope, to setup the instrument and detector so that the required scientific observation could be properly performed.

Acquisition phase






INAF-OATs

During the acquisition phase, after telescope presetting, a <u>correction</u> is necessary in order to center the object on the fibre entrance.

A carefully designed and tested algorithm is able to find automatically the most luminous object in the field (which usually coincides with the object under interest due to the small instrument field-ofview), computes the required correction and send it to the telescope.

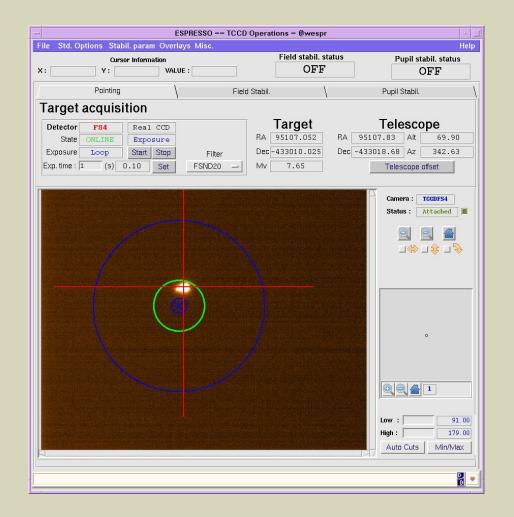

Stabilization phase








INAF-OATs

The purpose of the (field) stabilization is to maintain the object centered on the fiber during the whole duration of the scientific exposure. The TCCD(s), one per front end, acquires the image continuous loop up to 7 Hz; auxiliary software two in real-time, processes. analyze the obtained image and command the low-level piezo tip-tilt devices such to try keep the object well to centered on the fiber.

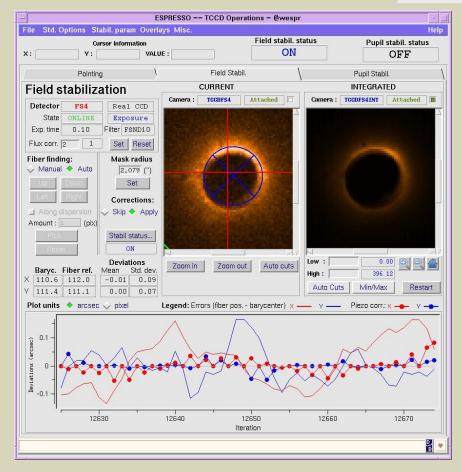


ESPRESSO exposure meter






The signal-to-noise ratio of the ongoing scientific exposure is continuously monitored by the <u>ESPRESSO</u> exposure meter (EM). In ESPRESSO, the EM is actually a <u>small spectrograph</u> in itself and a dedicated process continuously computes the counts in three spectra channels from images obtained by a dedicated TCCD.

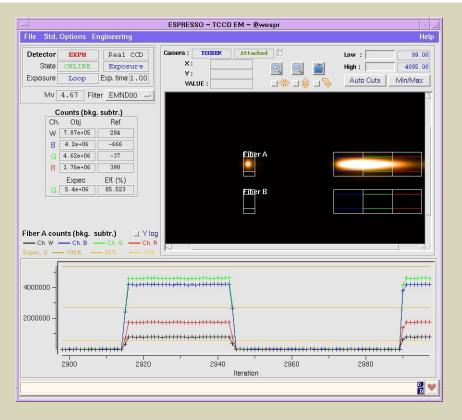

ESPRESSO degraded mode / 4 -UT





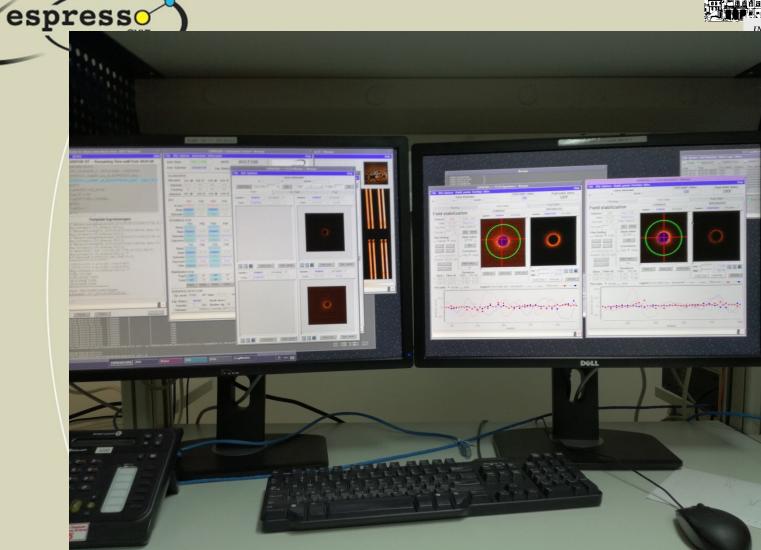

Science software





Thanks to its fixed spectral format and long-term stability, ESPRESSO has been conceived, starting already from the preliminary design phases, as "<u>truly science-grade products generating machine</u>".




To this purpose, in addition to the standard **Data Reduction Software (DRS)** package, a dedicated ESPRESSO **Data Analysis Software (DAS)** package has been developed with the aim to extract scientific information as soon as reduced data are available.

Pipeline, *esorex*, *esoreflex*





INAF-OATs

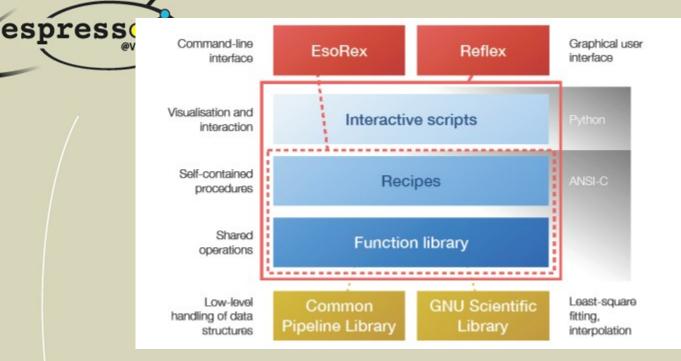

<u>@ Paranal</u>: pipeline as quick-look tools to do a real-time coarse assessment of the quality of calibrations and science observations; runs unsupervised with the help of data organizer (DO) and Reduction Block Scheduler (RBS) via OCA rules;

<u>ESPRESSO specific</u>: each science image can be reduced using the wavelength calibration taken few hours or even minutes before - <u>virtual product functionality</u>


@ Garching: more thorough data processing and evaluation - raw calibration frames into master calibrations; generate quality control parameters to monitor the instrument and detector performance

Data Reduction Cascade








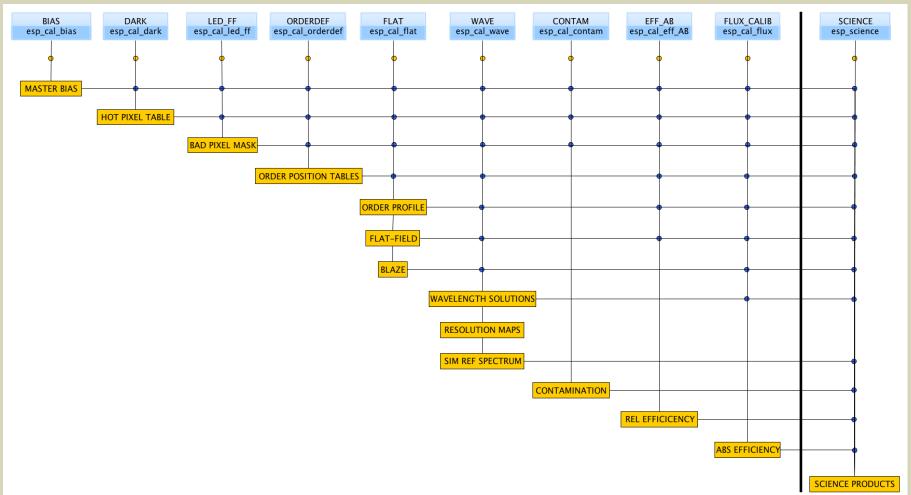


DRS example results





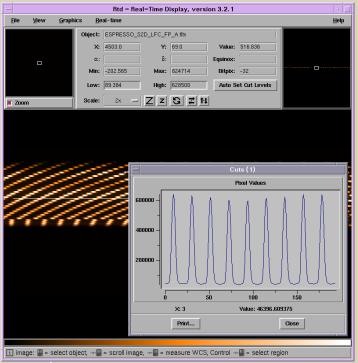

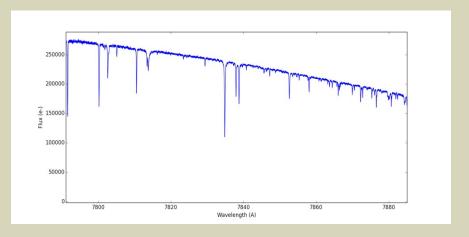

Wavelength coverage starts at 378.25 nm for the bluest order of fiber A and ends at 525.1 nm for the reddest order on the blue chip. On the red chip the shortest wavelength is 523.8 nm while the longest wavelength is 788.4 nm.

Single two-dimensional (S2D) spectrum extracted from the LFC exposure with one order shown by cuts. Regular spacing and the almost constant flux level could be easily noticed.

Additional info:

http://eso.org/sci/software/pipelines/espresso/espresso-pipe-recipes.html/

ESPRESSO Data Analysis





INAF-OATs

The ESPRESSO DAS comprises a total of 13 recipes. It is split into <u>four branches</u>: one for the <u>analysis of QSO spectra</u> and three for the <u>analysis of star spectra</u>.

Some of the recipes are common to more branches while others are specific to the treated spectra.

Each DAS branch is managed by a dedicated *Reflex workflow*.

All the recipes are available off-line at the telescope on a dedicated workstation and the users have the possibility to install the whole package also on <u>his/her computer</u>.

Operation	Branch	Name
Co-addition of spectra	star, quasar	espda_coadd_spec
Creation of a spectral mask	star, quasar	espda_mask_spec
Creation of a list of absorption lines	star, quasar	espda_create_linelist
Continuum fitting of a 1D spectrum	quasar	espda_fit_qsocont
Identification of the absorption systems	quasar	espda_iden_syst
Voigt-profile fitting of absorption lines	star, quasar	espda_fit_line
Measurement of the EW of spectral lines	star	espda_compu_eqwidth
Computation of the stellar parameters	star	espda_compu_starpar
Continuum fitting of single orders	star	espda_fit_starcont
Comparison with a synthetic spectrum	star	espda_synth_spec
Computation of the radial velocity using a synth. spec.	star	espda_rv_synth
Computation of the radial velocity	star	espda_compu_radvel
Computation of the stellar activity indexes	star	espda_compu_rhk

DAS workflow example






INAF-OATs



→ ♦ ♦ ♦

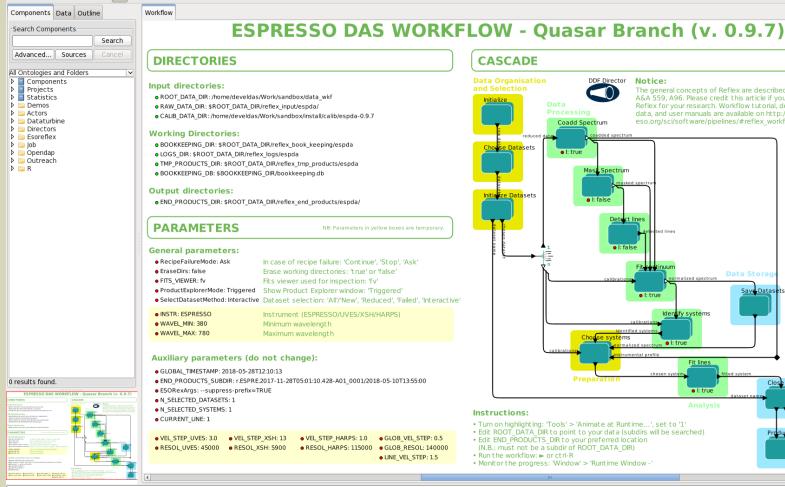

CASCADE DDF Director Notice: The general concepts of Reflex are described in A&A 559, A96. Please credit this article if you use Reflex for your research. Workflow tutorial, demo data, and user manuals are available on http://www. eso.org/sci/software/pipelines/#reflex workflows. Coadd Spectrum e Datasets • I: true e Datasets I: false Identify systems Fit lines Close Systems I: true Instructions: • Turn on highlighting: 'Tools' > 'Animate at Runtime...', set to '1' • Edit ROOT_DATA_DIR to point to your data (subdirs will be searched) • Edit END PRODUCTS DIR to your preferred location (N.B.: must not be a subdir of ROOT DATA DIR)

execution finished: 2517 ms. Memory: 681984K Free: 219687K (32%)

• Run the workflow: ► or ctrl-R


Monitor the progress: 'Window' > 'Runtime Window -

DAS quasar branch results





INAF-OATs

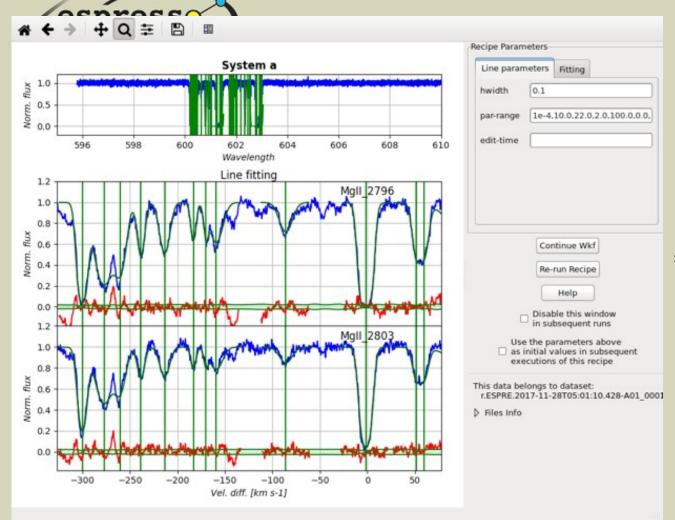


An absorption system (MgII) detected and fitted on the spectrum of quasar HE0515–4414 (observed during the comm 1b run, rebinned at 1 km/s). Blue: normalized spectrum; green: best-fit composite Voigt profile; red: residuals.


DAS star branch results






INAF-OATs

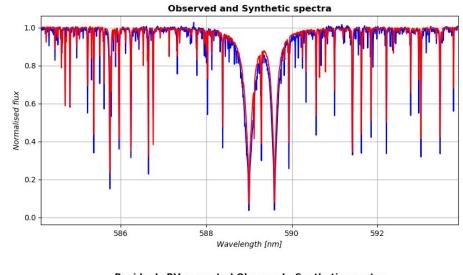

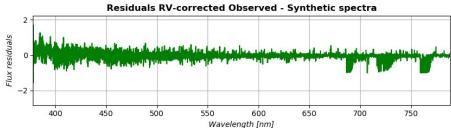

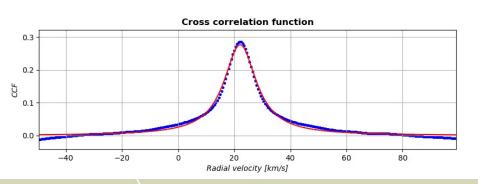

Top: plot of the output for the RV standard star HD 35854; **blue:** the observed normalized merged spectrum corrected for RV compared with the synthetic spectrum (**red**). **Center:** Observed - synthetic flux residuals. **Bottom:** cross-correlation function fitted with a Lorentzian function


First results I







INAF-OATs

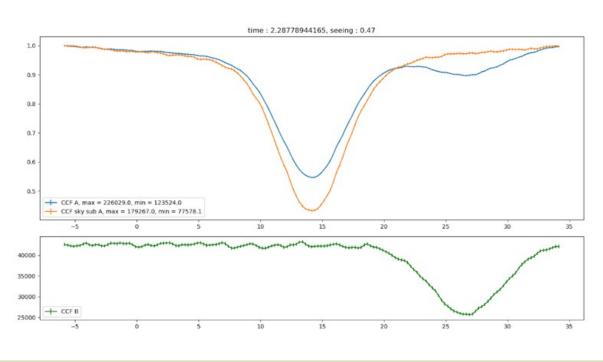


Cross correlation function (CCF) of a RV standard close to the Moon. The **blue line** is the CCF of the star spectrum and the RV dip corresponding to the Moon RV is clearly visible. The **orange line** is the CCF of the star spectrum after subtraction of the sky spectrum. **Bottom panel** show instead the CCF of the sky spectrum alone.

First results II



INAF-OATs

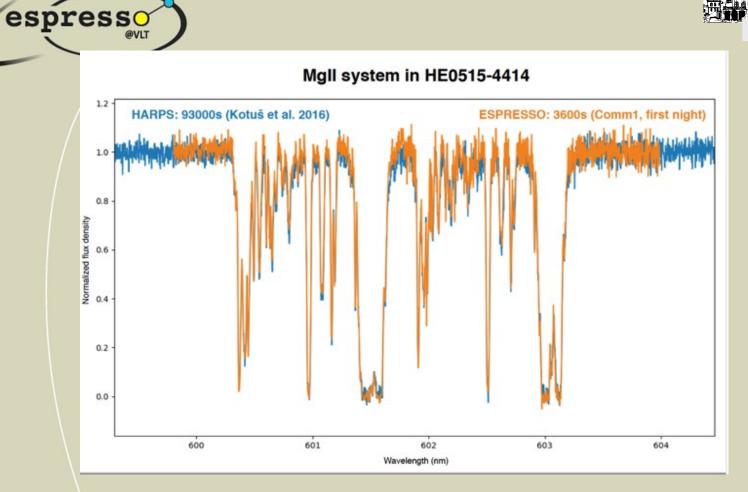


A comparison of a portion of HE0515-44 taken in 3600s the equivalent of 18 HARPS exposures for a total 93000 sec.

First results III





INAF-OATs

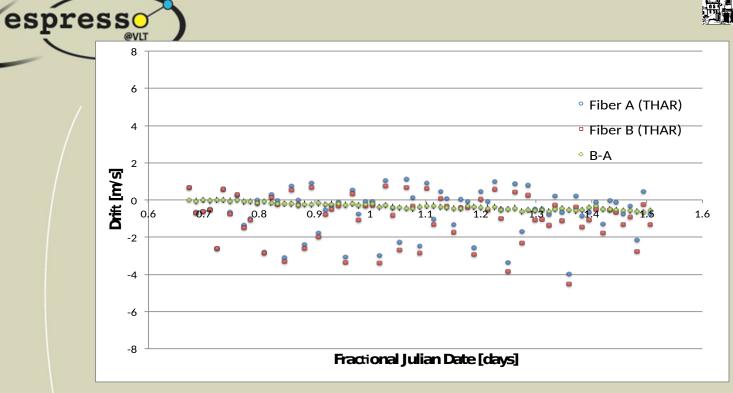

The image shows a drift sequence of more than 20 hours duration. The instrument drifted during this time by less than 1 m s⁻¹. On fiber A and B individually a scatter of about 1 m s⁻¹ rms is observed, which we consider to be due to the detector itself (self-heating due to fast read-out). The proof is that this scatter is perfectly common-mode. In fact, the differential drift B-A remained small for the whole duration. The 20-hours dispersion is about 18 cm s⁻¹, but if we consider the last 10 hours only, the dispersion is 9.4 cm s⁻¹, comparable with the photon noise and compliant with the requirements!

GTO - *Fundamental constants*






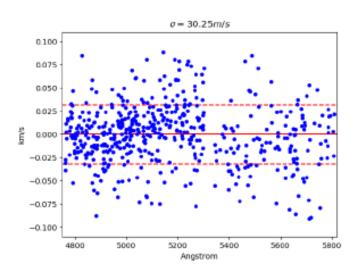


Figure 1 – Radial velocity differences between the solar lines as measured in the ESPRESSO 1UT Single 2x1 mode spectrum of Iris (6-Dec-2017) and the CERES-HARPS solar spectrum calibrated with the Laser Frequency Comb from ⁹. The gap at 530 nm is due to the HARPS gap.




Future: HIRES@ELT

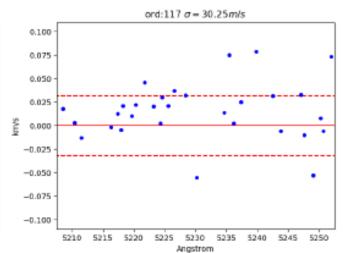

HIRES layout





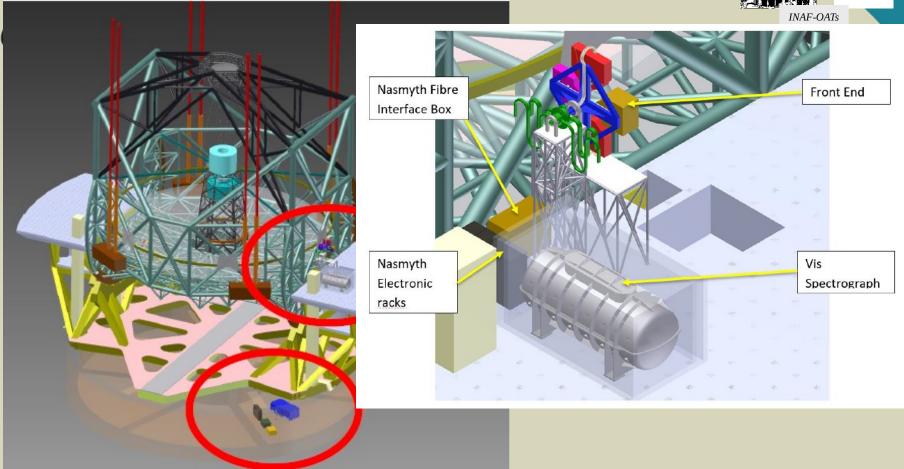







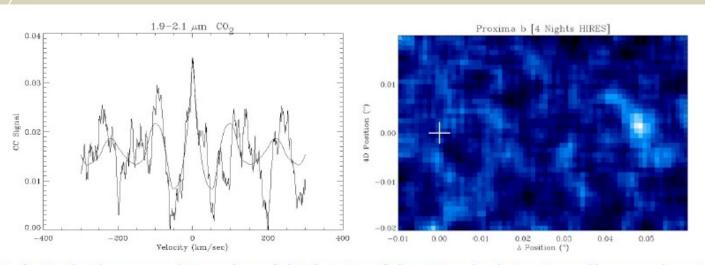


Figure 1 Left: Simulated cross-correlation values of the detection of the CO2 molecule in a Venus-like atmosphere in Trappist-1b or c, combining data of 4 transits with HIRES. Right: Reflected light cross-correlation signal of the direct surroundings of Proxima Cen.







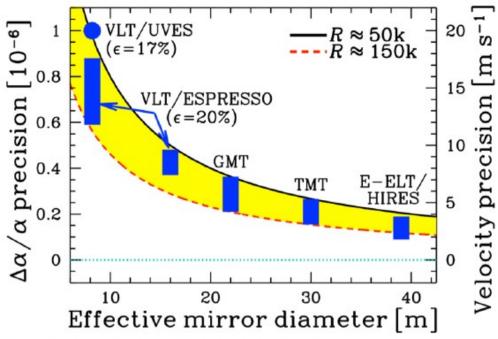


Figure 5 Expected statistical precision on variations in the fine-structure constant, $\Delta\alpha/\alpha$, achievable with future high-resolution spectrographs as a function of telescope diameter. The equivalent velocity precision is also shown (assuming a typical variety of metal-ion transitions). The length of the bars indicates the range of precision expected for different spectral resolutions available on those facilities. Two modes of operation for VLT/ESPRESSO are shown, its single-telescope mode (8-m effective diameter), with a range of resolving powers, and its anticipated four-telescope mode (16-m effective diameter), with R up to \approx 70,000.

Conclusions





In this talk I have presented ESPRESSO, its DFS, challenges and instrument peculiarities together with some first "results".

- <u>ESPRESSO</u>, the next generation ESO VLT high-resolution ultra-stable spectrograph, <u>is ready for operations</u>.
- After the very successful commissioning periods where instrument performance have been assessed and validated ESO gave green light for starting routine observations from October 2018 in single UT mode.
- 4 UT mode will be offered in 2019 after a final science verification phase.