

Investigating interfacial behavior of PDMS using an advanced JKR apparatus

Citation for published version (APA):

Vaenkatesan, V., Vellinga, W. P., Dekkers, E. C. A., Jeu, de, W. H., & Meijer, H. E. H. (2005). *Investigating interfacial behavior of PDMS using an advanced JKR apparatus*. Poster session presented at Mate Poster Award 2005 : 10th Annual Poster Contest.

Document status and date:

Published: 01/01/2005

Document Version:

Publisher's PDF, also known as Version of Record (includes final page, issue and volume numbers)

Please check the document version of this publication:

- A submitted manuscript is the version of the article upon submission and before peer-review. There can be important differences between the submitted version and the official published version of record. People interested in the research are advised to contact the author for the final version of the publication, or visit the DOI to the publisher's website.
- The final author version and the galley proof are versions of the publication after peer review.
- The final published version features the final layout of the paper including the volume, issue and page numbers.

[Link to publication](#)

General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal.

If the publication is distributed under the terms of Article 25fa of the Dutch Copyright Act, indicated by the "Taverne" license above, please follow below link for the End User Agreement:

www.tue.nl/taverne

Take down policy

If you believe that this document breaches copyright please contact us at:

openaccess@tue.nl

providing details and we will investigate your claim.

Investigating Interfacial Behavior of PDMS using an Advanced JKR Apparatus

V. Vaenkatesan,¹ Z. Li,¹ W. P. Vellinga,² E. C. A. Dekkers,³ W. de Jeu,¹ H. E. H. Meijer²

¹ Dept. of Chemistry and Chemical Engineering, ² Dept. of Mechanical Engineering, ³ GTD- Technical University of Eindhoven

Introduction

Johnson Kendall and Roberts (JKR) method is a technique to study the adhesive properties of soft polymeric systems that allows the determination of the work of adhesion (W) at an interface by means of contact mechanics. The evolution of W as a function of time for a contact is a natural extension of this subject and is expected to impact the adherence and frictional behavior of a contact as well. We developed an apparatus that is able to carry out JKR analysis, and also measure friction forces on a single contact. Additionally, the placement of the apparatus on a confocal microscope enables one to study the volume affected by the contact pressure along with the contact area.

Adhesion Studies

The JKR model allows a direct estimation of W & elastic modulus (E) of two elastic bodies in contact by the measurement of applied load (P), the resulting displacement (δ) and the radius of deformation (a) at the interface under equilibrium conditions.

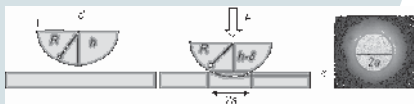


Fig 1. Schematic representation of the geometrical changes when an elastic semi-sphere and a sheet are brought in contact under a load P along with an image of the resulting contact area.

For the case of an elastic semi-sphere in contact with a flat sheet, the following equations describe the JKR model:

$$\delta = \frac{a^2}{3R} + \frac{2P}{3aK} \quad (1)$$

$$a^3 = \frac{R}{K} \left[P + 3\pi WR + \sqrt{6\pi WRP + (3\pi WR)^2} \right] \quad (2)$$

$$\frac{1}{K} = \frac{3}{4} \left[\frac{1-\nu_1^2}{E_1} + \frac{1-\nu_2^2}{E_2} \right]$$

ν - Poisson ratio
 R - radius of curvature of lens

Instrument

A double-leaf spring unit moved using the motors is used as the force probe. The deflection of the leaf-spring on contact of samples is proportional to the applied force. The system is mounted on a microscope to record the contact area.

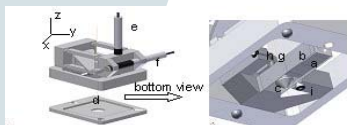


Fig 2. Schematic of the JKR instrument: a, b – vertical & horizontal leaf-spring unit; c – sample holder (lens); d – sample holder (sheet); e, f – motors; g, h – capacitive sensor for vertical and horizontal force detection (resolution ~ 15 μ N); i – capacitive sensor for indentation (resolution ~ 10 nm) measurement.

Materials

Influence of molecular mass (M_c) between crosslinks in poly(dimethylsiloxane) (PDMS) on its elastic modulus, work of adhesion and its friction dynamics has been investigated.

/ department of mechanical engineering
PO Box 513, 5600 MB Eindhoven, the Netherlands

Results

In a typical experiment, PDMS samples were mounted on the sample holder and a and δ were recorded as a function P . The behavior of $P - \delta$ and $P - a$ curves with changing M_c is shown in Fig. 3. Values obtained from fitting of the curves to Eq. (1) & (2) is shown in Table 1.

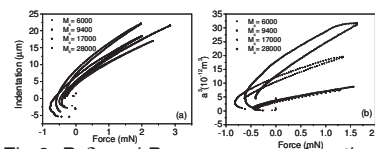


Fig. 3. $P - \delta$ and $P - a$ curves as a function of M_c .

M_c (g.mol ⁻¹)	E_δ (Eq. 1) (MPa)	E_L (Eq. 2) (MPa)	W_L (mJ/m ²)
6000	0.94±0.10	1.16±0.14	45±3
9400	0.91±0.05	1.03±0.05	45±3
17000	0.84±0.02	0.93±0.03	44±2
28000	0.60±0.11	0.71±0.14	42±3

Table 1. Shows the influence of M_c on E & W of PDMS.

Friction Studies

Influence of M_c and sliding speed on friction dynamics (lens sliding on sheet) of PDMS was investigated and is shown in Fig. 4.

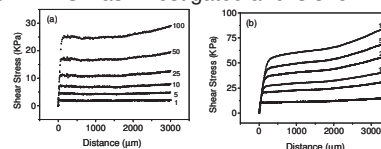


Fig. 4. Friction behavior of PDMS with M_c (a) 6000 (b) 28000 : lens was slid on sheet at various speeds (μ m/s) under a normal stress of 20 KPa.

Confocal Imaging

Fig. 5 shows a composite of confocal images recorded along the z-axis of PDMS lens in contact with a flat surface and surrounded by a fluorescent dye. The image gives an impression of the lens shape above the contact area and also the shape of the meniscus of the dye solution.

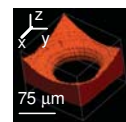


Fig. 5. Confocal image of PDMS lens in contact with a flat surface.

Conclusions

- The instrument is an excellent tool for performing adhesion and friction studies on a single contact of soft materials in a dynamic mode.² Imaging of the sample interface and volume during such processes is also possible.
- Novelty of the instrument allows accurate measurement of δ and determine E using Eq. (1) independent of Eq. (2).
- The E of PDMS decreases with increasing M_c though W remains largely unaffected. Adherence increases with increasing M_c (indicated by hysteresis behavior in Fig. 3).³
- Both adherence and friction are correlated in PDMS increasing with increasing M_c .³

References

[1] K. R. Shull, *Mat. Sci. & Eng. R* **36**, 1, 2002.
[2] V. Vaenkatesan et. al., (submitted to *Langmuir*, 2005).
[3] A. Galliano et. al., *J. of Colloids & Interfaces* **265**, 372, 2003.