JAAS

PAPER

Check for updates

Cite this: DOI: 10.1039/c9ja00288j

Received 20th August 2019 Accepted 27th November 2019

DOI: 10.1039/c9ja00288j

rsc.li/jaas

Introduction

Rhenium (Re, atomic mass 186.207) is one of the least abundant chemical elements present on Earth, with concentrations in most materials ranging from 10^{-12} g g⁻¹ (ppt) to 10^{-9} g g⁻¹ (ppb),¹⁻³ with 10^{-6} g g⁻¹ (ppm) only in found in specific phases.⁴ Several characteristics make Re a useful element for tracking geochemical processes. Notably, its redox-dependent solubility makes it a target for tracking oxidation reactions during chemical weathering^{1,5-7} and reconstructing redox cycling in lake and oceanic sediments.⁸

Rhenium oxidation states range from -1 to +7, with the most abundant being +7, +6 and +4.° Rhenium is a soluble element in oxic conditions and mostly present in surface waters in the form of the oxyanion Re^{VII}O₄⁻.^{1,2,10} In seawater, Re has a conservative behaviour with an average present-day concentration around 7.5 ppt.^{2,3} Under anoxic conditions, Re^{VII} is reduced to Re^{IV}, becomes insoluble and is removed from water either through complexation of Re^{IV} with organic matter and/or incorporated into sulphides.^{11–13} Among the redox sensitive metals, Re has the largest enrichment factor in authigenic phases of anoxic



Mathieu Dellinger, 🗅 ** Robert G. Hilton* and Geoffrey M. Nowell 🕑 b

Rhenium (Re) is a trace element whose redox chemistry makes it an ideal candidate to trace a range of geochemical processes. In particular, fractionation of its isotopes ¹⁸⁷Re (62.6% abundance) and ¹⁸⁵Re (37.4%) may be used to improve our understanding of redox reactions during weathering, both in the modern day and in geological archives. Published methods for measurement of Re isotopic composition are limited by the requirements of Re mass to reach a desirable precision, making the analysis of many geological materials unfeasible at present. Here we develop new methods which allow us to measure Re isotope ratios (reported as δ^{187} Re) with improved precision: $\pm 0.10\%$ (2σ) for a mass of Re of ~ 1 ng to $\pm 0.03\%$ (2σ) for a mass of Re of >10 ng. This is possible due to the combination of a modified column chemistry procedure and the use of 10^{13} Ω amplifiers for measurement *via* multicollector inductively coupled plasma mass spectrometry (MC-ICP-MS). For river water samples (with Re concentrations typically $\sim 10^{-12}$ g g⁻¹) we design a field-based pre-concentration of Re that can be used with large volumes of filtered water (5–20 L) shortly after sample collection to provide abundant Re for isotope analysis. As a result of these developments we provide new measurements of δ^{187} Re in standards reference materials (δ^{187} Re values range from $-0.06 \pm 0.07\%$ to $+0.19 \pm 0.05\%$) and a seawater standard (δ^{187} Re = $+0.10 \pm 0.04\%$), providing impetus for further exploration of the Re isotope system.

sediments relative to the detrital background.⁸ As such, high concentrations of Re are observed in Phanerozoic sedimentary rocks formed under reducing conditions, such as black shales (average 150 ppb), relative to the upper continental crust (~0.3 ppb) and most Precambrian sedimentary rocks (10–30 ppb).^{8,14} Earth's core represents the main Re reservoir on Earth (230 ppb) and the abundance of Re in Bulk Silicate Earth (BSE, 0.35 ppb) is about two orders of magnitude lower than in carbonaceous chondrite (about 40 ppb).^{15,16} During magmatic processes, Re behaves as a moderately incompatible element during melting and differentiation processes.^{15,17,18}

View Article Online

Rhenium is composed of two naturally abundant isotopes ¹⁸⁷Re (62.6%) and ¹⁸⁵Re (37.4%). The ¹⁸⁷Re isotope is radioactive and undergoes β -decay to ¹⁸⁷Os with a half-life of 4.35 \times 10¹⁰ years. Following the pioneering work of Miller et al., (2009),19 two studies have measured the Re isotopic composition expressed as δ^{187} Re, where δ^{187} Re = ((187 Re/ 185 Re)/ $\left(^{187}\text{Re}/^{185}\text{Re}\right)_{std}~-~1\right)~\times~1000,$ in shales and weathered sedimentary rocks9 and iron meteorites.20 These studies have documented a total of ~0.80% variability in δ^{187} Re values, including $\sim 0.30\%$ across a chemical weathering profile in soil developed on a Devonian black shale9 and 0.40% between iron meteorites.²⁰ In addition, an *ab initio* calculation⁹ shows that oxidized Re^{VII} species are usually enriched in ¹⁸⁷Re relative to more reduced Re^{IV}. The fractionation factor between Re^{IV} and $\operatorname{Re}^{\operatorname{VII}}(10^3 \times \log e(\alpha_{\operatorname{Re}^{\operatorname{IV}}-\operatorname{Re}^{\operatorname{VII}}}))$ could potentially vary between -0.80% and +0.50%, depending on the degree of thiolation of

^aDepartment of Geography, Durham University, DH1 3LE Durham, UK. E-mail: mathieu.dellinger@durham.ac.uk

^bDepartment of Earth Sciences, Durham University, DH1 3LE Durham, UK

[†] Electronic supplementary information (ESI) available. See DOI: 10.1039/c9ja00288j

the Re^{VII} species.⁹ Overall, the existing measurements and calculations suggest that Re isotope ratios have strong potential to help track redox processes.

The previous measurements of the Re isotopic composition of materials^{19,20} were made by multicollector inductively coupled plasma mass spectrometry (MC-ICP-MS - Thermo Fisher Scientific Neptune, fitted with $10^{11} \Omega$ amplifiers) with the combination of tungsten (W) doping and standard-sample bracketing to correct for instrumental mass bias. With that set-up, a precision of $\pm 0.10\%$ (2 σ) was achieved for a concentration of 10 ppb Re. Based on that pioneering work, many sample types were thus considered too low in Re concentration (e.g. river water, felsic rocks) to permit δ^{187} Re measurement. Here, we revisit the Re isotopic analysis by MC-ICP-MS with the aim to establish a method to measure Re isotopic composition in low abundance samples which are relevant to understanding crucial redox processes on Earth. These include grey shales with lower organic carbon contents than black shales (and Re < 1 ppb),^{6,7} felsic igneous rocks (Re < 2 ppb),^{17,21} weathered rocks and soils (Re \sim ppt-ppb range),^{6,7,22} river waters and seawater (Re \sim ppt).^{1,2} To do this we: (i) make use of recent mass spectrometry hardware advances (notably the development of $10^{13} \Omega$ amplifiers²³) and very low uptake rate nebulizer (CF 35) which allows us to run smaller volumes at higher concentrations; (ii) explore and refine sample preparation and purification by column chemistry for solids; (iii) propose a new field preconcentration method for river water samples; and (iv) measure and report a range of natural low Re abundance standard reference materials. We recommend repeated column chemistry procedures (loading and elution) to purify samples. The methods allow for the measurement of δ^{187} Re values to a precision of better than $\pm 0.05\%$ (2 σ) for Re mass > 3 ng, unlocking the vast majority of natural samples for the analysis of Re isotopic composition.

Experimental methods

Reference materials and samples

In this study, we used a wide range of SRM (Standard Reference Materials), including basalts (BCR-2, BCR-1, BHVO-2, BIR-1), sediments (MAG-1, SCO-1, NIST 1646a), granites (JG-2), peridotite (GP13), serpentinite (UB-N), andesite (AGV-1), diabase (TBD-1) and dolerite (DNC-1). We also used river samples from different locations (UK, Canada, France), a black shale and a carbonaceous chondrite sample (see ESI†).

Dissolution of solid samples

The low Re concentration of most rock and sediment samples mean that a mass of at least ~0.5 g is necessary for a precise δ^{187} Re. To ensure complete re-dissolution of fluorides, a maximum aliquot mass of 0.5 g per beaker was dissolved for each sample. A volume of 3 mL 29 M HF and 3 mL of 16 M HNO₃ was added to each sample in PFA Teflon beakers and heated at 120 °C for at least 24 h. Following initial digestion, samples were evaporated to dryness at 80 °C and then re-dissolved in aqua regia to destroy fluorides, heated at 120 °C for 24 h before

This article is licensed under a Creative Commons Attribution-NonCommercial 3.0 Unported Licence. Dpen Access Article. Published on 07 January 2020. Downloaded on 1/8/2020 9:42:56 AM

HCl. Care was taken to ensure complete dissolved in 1 M HCl. Care was taken to ensure complete dissolution and a volume of about 20 to 40 mL of 1 M HCl was required to achieve this. For samples containing refractory organic matter, a black residue was treated in 16 M HNO₃ and aqua regia for several days at 160 °C. Ultimately, only the most refractory organic matter (*e.g.* graphite) was left un-dissolved and removed from the solution by centrifugation and pipetting of the supernatant.

We tested for the potential loss of Re during heating and evaporation for 16 M HNO₃, 10 M HCl, H₂O, aqua regia and HClO₄, at temperatures of 80, 120, 150 and 190 °C. We found no loss of Re and recovery of 100% for all reagents except for HClO₄, for which there are significant evaporative losses of Re (>85%) at temperatures higher than 150 °C.

Chemical separation for solid samples

The chemical separation procedure is modified from those previously described.^{19,24} Polypropylene columns (inner diameter of 7.1 mm) are filled with 1 mL of AG1-X8 resin (200-400 mesh). A frit was added to the top of the resin-bed to avoid remobilization of the resin during addition of reagents. The resin is cleaned with 30 mL of 8 M HNO3, and conditioned with 5 mL of 1 M HCl. Samples are loaded on the column in 1 M HCl. At low acid concentration, Re is strongly bounded to this resin,²⁵ with the partition coefficient of Re with the AG1-X8 resin ($K_{d Re}$ = [Re]_{resin}/[Re]_{solution}) > 100 for HCl and HNO₃ concentrations < 1 mol L^{-1} . Before collecting the Re fraction, elution of the sample matrix is achieved in 3 steps: (i) addition of 10 mL of 1 M HCl, (ii) addition of 15 mL of 0.5 M HNO3 and (iii) 1.5 mL of 4 M HNO₃ (Table 1). The Re fraction is eluted with 12.5 mL of 4 M HNO₃ ($K_{\rm d, Re} \sim 8$ –10). The Re fraction is then evaporated at temperature of 100 to 120 °C to complete dryness and refluxed in 16 M HNO3 at 120 to 150 °C for at least 24 h. This full procedure is repeated two more times (i.e. three identical column chemistry steps) to purify the Re fraction from the residual matrix. Samples are refluxed in 16 M HNO3 after each column step to destroy resin-derived organic residues. A small aliquot (2% of the sample Re mass) is taken after 1 column pass to measure the Re concentration and after 3 column passes to check the total yield of the separation.

Chemical separation for dissolved samples

For dissolved Re samples (*e.g.* filtered river water and seawater), a volume of 1 to 20 L is usually necessary to recover enough Re for accurate δ^{187} Re (>0.4 ng of Re). Such volumes present significant handling and shipping issues. To overcome this, we developed a method for pre-concentrating Re in the field. River water samples are filtered at 0.2 µm (through PES Millipore filters) and collected in 20 L sterile double-lined beverage bags that are commercially available. The sample mass is measured and the bags are then connected (using a custom fitting) to a BioRad Econo-Pac® 20 mL column filled with 2 to 4 mL of AG1-X8 resin, and passed through the column in a field laboratory (typically over a 6–10 h period). The Re-loaded columns are then capped, packaged and returned to Durham. Upon

Table 1 Column chemistry procedure

Step	Acid	Volume (mL)	Elements eluted
Resin AG1-X8; 1 mL			
Precleaning	$HNO_3 8 N$ (twice)	15	
Precleaning	H ₂ O	5	
Equilibration	HCl 1 N	5	
Introduction	HCl 1 N	2-50	Bulk elements
Cleaning step 1	HCl 1 N	10	Residual bulk elements + Mo, Nb, Ti, Ta, Pb, Ag, W
Cleaning step 2	HNO ₃ 0.5 N	15	Zn, Cd and Residual Pb, Ta, Nb
Cleaning step 3	HNO ₃ 4 N	1.5	Residual Ta, Nb, Ag
Elution of Re	HNO ₃ 4 N	12.5	Re

^{*a*} Bulk elements are elements for which more than 85% of the elements is not fixed to the resin. It includes Li, Be, Na, Mg, Al, K, Ca, V, Cr, Mn, Fe, Co, Ni, Cu, Se, Sr, Y

return to the laboratory, 50 mL of 1 M HCl is passed through the resin to remove some of the matrix (and if required, to collect the SO₄ fraction for sulfur isotope analysis²⁶). Finally, the Re is collected with 30 to 40 mL of 8 M HNO₃. Using a higher HNO₃ molarity allows better recovery of the proportion of Re bound to the resin. Organic molecules from dissolved organic carbon also have strong affinity with this resin (the resin colour becomes darker for water samples with high dissolved organic carbon) and are partly eluted with Re. Therefore, to destroy these organics, the Re residue is re-dissolved in 16 M HNO₃ and refluxed for 5–7 days at 150 °C. These samples are then passed through the three steps of column chemistry separation as described for solid samples.

Mass spectrometry

Following the method developed by Miller et al., (2009),¹⁹ rhenium isotopic ratios were measured by MC-ICP MS (ThermoFisher Scientific NetunePlus without jet interface) at the Arthur Holmes Isotope Geology Laboratory, Dept. Earth Sciences, University of Durham. We use a combination of standard-sample bracketing and external normalization to tungsten (W). Following separation of Re, samples were redissolved in 0.5 M HNO3 and a small aliquot (2% of the sample mass) was taken to measure the Re concentration by MC-ICP MS. Once the mass of Re in the remaining sample aliquot is determined, a known amount of pure W solution (ROMIL PrimAg Mono-Component Reference Solution) was added to each sample, so that the ratio (in $g g^{-1}$) between W and Re in the sample is 20. Samples were run in 0.5 M HNO₃ and introduced using a either PFA-50 or CF 35 nebuliser together with a micro-cyclonic spray chamber. For the CF 35 nebuliser, the measured uptake rate was ${\sim}37\,\mu L~min^{-1},$ giving a sensitivity of 0.26 V (¹⁸⁷Re) for a 10 ppb Re solution. In this study, we have not tested whether the use of other introduction systems (Aridus, Apex) could improve the sensitivity of the measurement. The standard NEPTUNE sample and H skimmer cones were used. Solution uptake time is 50 to 60 seconds with a 70 second wash time. Complete measurement of a sample (2 sample replicate analysis and 3 standards) takes \sim 50 minutes, consuming ~ 0.65 mL. The cup configuration allows

simultaneous measurement of Re (¹⁸⁵Re and ¹⁸⁷Re), W (¹⁸²W, ¹⁸⁴W and ¹⁸⁶W) and monitoring of Os (¹⁸⁸Os, ¹⁸⁹Os and ¹⁹⁰Os) isotopes for possible isobaric interference (see ESI†).

Analysis were performed in low mass resolution mode (MR ~400) and consisted, like previous methods,19 of 25 cycles with an integration time of 16.77 s each. Whether this method is optimal for accuracy has not been fully assessed here and remains to be investigated. Hereafter, the term "analysis" is used to refer only to the 25 cycles measured on a given sample aliquot solution (one or several analysis) whereas the term "measurement" refers to the complete procedure (including digestion for solids, preparation and dilution). As Re concentrations are low in geological materials, our study used $10^{13} \Omega$ amplifiers for detectors measuring ¹⁸⁷Re and ¹⁸⁵Re isotopes and $10^{11} \Omega$ amplifiers for other detectors. The $10^{13} \Omega$ amplifiers allow improved precision over $10^{11} \Omega$ amplifiers on ion beams less than 0.52 V (ref. 23) (see section "Result and discussion"). Samples were typically measured at least twice ("replicate analysis" of the same purified Re solution) or more, depending upon the mass of Re available, with a bracketing standard in between (e.g. standard-sample-standard-sample-standard). All analytical sessions were carried out automatically using an SC-µ DX autosampler to ensure that the measurements and washouts followed a strict time sequence; necessary for interpolating and correcting for instrumental mass bias (IMB) using a standardsample bracketing analysis method (see below).

In plasma mass spectrometry IMB is the tendency of measured isotope ratios to be biased in favour of the heavier isotope.²⁷ To a first approximation IMB is only mass dependant in MC-ICP-MS, although in reality IMB can drift slightly over time or can change as a result of instabilities in the plasma or of non-spectral interferences/matrix load on the plasma (due to residual elements).²⁷ Rhenium has only two isotopes so IMB cannot be corrected for internally. There are two alternative approaches to correct for Re mass bias, both of which we employ in this study. It can be corrected for by doping the Re with an element of similar mass but with no isobaric masses (*e.g.* Ir or W)^{20,28} and normalizing to a known ratio of that element (a method referred to as "external normalization", EN). Similar to previous studies,^{19,20,29} we doped the Re with tungsten

and corrected for Re mass bias using an exponential law using the measured $^{186}W/^{184}W$. The related equations are:

$$\left(\frac{^{187}\text{Re}}{^{187}\text{Re}}\right)_{\text{corrected}} = \left(\frac{^{187}\text{Re}}{^{187}\text{Re}}\right)_{\text{measured}} \times \left(\frac{m_{187}}{m_{185}}\right)^{\beta_{\text{w}}}$$
(1)

with:

$$\beta_{\rm w} = \frac{\ln\left(\frac{\left(\frac{186}{184}{\rm W}\right)_{\rm true}}{\left(\frac{186}{184}{\rm W}\right)_{\rm measured}}\right)}{\ln\left(\frac{m_{186}}{m_{184}}\right)}$$
(2)

where $({}^{187}\text{Re}/{}^{185}\text{Re})_{\text{measured}}$ and $({}^{187}\text{Re}/{}^{185}\text{Re})_{\text{corrected}}$ are the measured and W-corrected Re isotope ratio, m_{187} and m_{185} the molar masses of each Re isotopes, $({}^{186}\text{W}/{}^{184}\text{W})_{\text{true}}$ and $({}^{186}\text{W}/{}^{184}\text{W})_{\text{measured}}$ are the real and measured W isotope ratios, while m_{186} and m_{184} are the molar masses of each W isotope. This method assumes that W and Re share the same degree of mass bias in the instrument (*i.e.* $\beta_{\text{W}} = \beta_{\text{Re}})^{29}$ and allows us to reduce both the uncertainty arising from instrumental mass bias and mitigate the influence of matrix effects.

The measured delta value expressed as below:

$$\delta^{187} \mathbf{R} \mathbf{e} = \left[\frac{(^{187} \mathbf{R} \mathbf{e} / ^{185} \mathbf{R} \mathbf{e})_{\text{sample}}}{\left((^{187} \mathbf{R} \mathbf{e} / ^{185} \mathbf{R} \mathbf{e})_{\text{std}-1} + (^{187} \mathbf{R} \mathbf{e} / ^{185} \mathbf{R} \mathbf{e})_{\text{std}-2} \right) \times 0.5} - 1 \right] \times 1000$$
(3)

Where $({}^{187}\text{Re}/{}^{185}\text{Re})_{\text{sample}}$ and $({}^{187}\text{Re}/{}^{185}\text{Re})_{\text{std}}$ can be either the W-corrected or measured ratio. Herein, we refer to the notation " $\delta^{187}\text{Re}$ " as the $\delta^{187}\text{Re}$ determined by combining W-corrected (External Normalisation) and Standard-Sample bracketing ("EN–SSB method"). We refer to " $\delta^{187}\text{Re}^{\text{SSB}}$ " as the Re isotopic composition determined only by SSB. Two consecutive replicate analysis of the same sample solution allow calculation of three Re isotope delta values. Following Miller *et al.* (2009),¹⁹ we use the Re standard reference material NIST SRM 989 as the normalising standard. Other pure Re standards ("HReO₄" and "NIST 3141")¹⁹ and in-house standard solution "DURH-Re-1" are also used as secondary standards.

Results and discussion

Yield of the column chemistry separation and Re concentration measurements

We evaluate the yield of our column chemistry procedure by measuring the Re concentration of standard reference materials (SRMs) and a pure Re solution after a single pass of the column chemistry. In general, we find a good agreement (within $\pm 10\%$) between our measured Re concentration values and those from the literature measured by isotope dilution^{30–32} (Table 2). Measured yields range from 92% to 124% with a median value of 96% (N = 44 digestions on 10 different SRM) with no systematic differences between rock types. This shows that HF– HNO₃ digestion and column chemistry can be used to determine Re concentration in solid material with only slightly less precision than by isotope dilution techniques. For the pure Re solution, we find yields ranging between 85% and 105%, with an average value of 98.9% (N = 10). This indicates that additional matrix in rock samples does not detrimentally lower the Re yield.

Efficiency of the separation

Our column chemistry procedure aims to remove inorganic and organic matrices that can create spectral and non-spectral effects and impact the accuracy of the Re isotope measurement. A large number of elements (Li, Be, Na, Mg, Al, K, Ca, V, Cr, Mn, Fe, Co, Ni, Cu, Se, Sr, Y, Cs, Nd, Hf, Th, U) are not retained by the resin (Fig. 1) and are eluted almost immediately following the initial loading step (>85%). Some elements are partly eluted during the loading step and during the initial 1 M HCl elution step (Ti, As, Zr, Mo, Sb, W). Zinc is fully retained by the resin until the 0.5 M HNO₃ step (Fig. 1) at which point it is quantitatively removed along with most of the remaining Cd (70%). About 60 to 80% of the Nb is eluted during the sample loading and 1 M HCl steps, with the remainder eluted during the 0.5 M and 4 M HNO3 steps. Thallium is the element with the highest proportion eluted during the Re elution step (about 30% of the original thallium). Importantly for Re isotope measurement protocol, where samples are spiked with a W solution that is isotopically homogenous, W originating in the sample is almost quantitatively removed prior to the Re elution peak (Fig. 1). Compared to previous methods,19,20 the use of the 1 M HCl step prior to the 0.5 M HNO₃ allows for a more efficient removal of most chemical elements. The Re blank of the whole method (including digestion) range from 0 to 14 pg with a mean value of 3 ± 4 pg ($n = 9, 1\sigma$).

Fractionation of Re isotopes during the elution

Previous studies^{19,20} have shown that Re isotopes are fractionated during column chemistry and that a minimum yield of 80% is required for preventing resolvable fractionation of Re isotopes on the column. We measure δ^{187} Re as a function of the yield of the procedure for a pure Re solution (Fig. 2) and find a similar behaviour as previous work.^{19,20}

Field pre-concentration of Re from water samples

The Re recovery of field pre-concentration columns were checked for all water samples. For a sample set of 55 water samples, we find recovery yields ranging from 43% to 118% with a median value of 95.7%, with 48 of 55 samples having recovery yield between 75 and 120% (considering an uncertainty of $\pm 15\%$). Possible fractionation of Re isotopes due to incomplete recovery of the Re pre-concentrated on the resin has been investigated by recovering various amount of pure Re solution on the 3 mL of BioRad resin used for pre-concentration (Fig. 2; filled squares). For Re recovery >30%, no Re isotopic fractionation compared to the 1 mL column is probably due to the different aspect ratio of the column used for the field and/or the different eluent concentration (here 8 M HNO₃).

Table 2	Results for Re concentration measurements
---------	---

	This study				Literature			
SRM name	Mass digested per sample (mg)	[Re] (ppb)	±2SD (in %)	N measurements	[Re] (ppb)	Uncertainty (in %, 2σ)	Ν	Reference
BCR-2	60 to 514	11.56	6.4	8	12.60	15.9	4	Jochum <i>et al.</i> , $(2016)^{32}$
BCR-1	480 to 653	0.81	9.6	9	0.84		30	Jochum <i>et al.</i> , $(2016)^{32}$
UB-N	466 to 698	0.20	14.3	9	0.21	4.8	14	Meisel and Moser (2004) ³
MAG-1	172 to 559	3.65	7.3	6	3.91		1	Meisel and Moser (2004) ³
BHVO-2	498 to 617	0.61	12.8	5	0.54	5.3	6	Jochum <i>et al.</i> , $(2016)^{32}$
BIR-1	208	0.70		1	0.65		2	Jochum <i>et al.</i> , $(2016)^{32}$
JG-2	379	0.015		1	0.016		1	Imai <i>et al.</i> , $(1995)^{33}$
GP-13	200	0.30		1	0.32	14.2	4	Meisel and Moser (2004) ³
AGV-1	197	0.38		1	0.38		3	Jochum <i>et al.</i> , $(2016)^{32}$
TDB-1	497	0.99		1	0.79	6.0	7	Meisel and Moser (2004) ³
DNC-1	197	0.85		1	0.91		1	Meisel and Moser (2004) ³
SCO-1	130	0.97		1	1.01	12.4	5	Meisel and Moser $(2004)^3$
SCO-1	591	0.04		1				
NIST 1646a	515	1.99		1				

Ti

Zn Zr

Nb Mo Ag Cd Ta W Re Pb

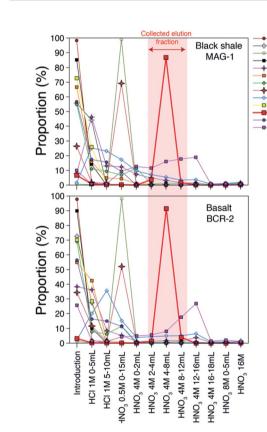


Fig. 1 Elution curves of chemical elements and Re for two standard materials. Element concentrations were measured by quadrupole ICP-MS. Bulk elements are shown in Table 1.

Evaluation of analytical errors

Here we discuss the different sources of error^{34,35} during Re isotopic measurement by MC ICP-MS. We refer to the "internal error" as the standard error (SE, in 2σ) of the 25 cycles of 16.77 second-integration periods (corresponding to one "analysis"). The "intermediate error" is the standard deviation (SD, in 2σ)

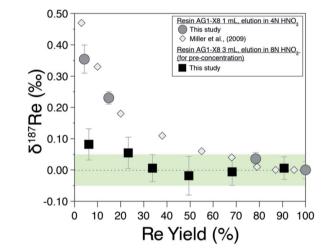


Fig. 2 Rhenium isotope composition of pure Re solution (NIST-SRM989) as a function of the proportion of recovery from (i) field-based columns (black squares) and (ii) laboratory columns (grey circles). White diamonds are from Miller *et al.*, (2009).¹⁹ The green area is the $\pm 0.05\%$ long-term external reproducibility of pure Re solution at 5 ppb Re concentration.

of repeated analysis of the same solution, (here 2 sample analysis and 3 bracketing standards, giving 3 measured δ^{187} Re). This within-run-precision will be influenced by mass bias instability over tens of minutes to an hour. Finally, the "external error" corresponds to the standard deviation (SD, in 2σ) of several measurements of a given solid or dissolved material (including digestion for solids, preparation and dilution) over an extended period of time. This between-run-precision will be influenced by day-to-day differences in instrumental plasma conditions and by sample preparation steps for samples.

There are two major sources of noise that influence the internal error: counting statistics and the Johnson-Nyquist noise.^{27,34,36} The internal error for a given signal intensity

corresponds to the sum of these two sources of error. The counting statistics can be predicted following the Poisson distribution.^{27,34,35} The relative standard error (RSE) in the isotope ratio due to counting statistics ($\sigma_{\text{counting statistics}}$) is:

$$\sigma_{\text{counting statistics}} = \sqrt{\frac{n_{187} + n_{185}}{n_{187} \times n_{185}}} \tag{4}$$

Where n_{187} and n_{185} represent the total number of ions collected for each Re isotope for a given period of time, in this case analysis. The latter are calculated as:

$$n_{18x} = \frac{V_{18x}}{(R \times \text{Gain})} \times t \times 6.241 \times 10^{18}$$
(5)

Where *R* is the resistor amplifier (here $10^{13} \Omega$), "Gain" the gain (0.01), " V_{18x} " the voltage for ^{18x}Re, "t" is the total integration time of the measurement (25 cycles of 16.77 seconds) and 6.241 \times 10¹⁸ is the conversion factor for ampere to number of ions. An intensity of 0.1 volt measured on a $10^{13} \Omega$ amplifier, with a gain of 0.01, for 1 second corresponds to 10^{-12} A, which gives a number of ions $n = 6.24 \times 10^6$ cps.V⁻¹. The Johnson–Nyquist noise corresponds to the fluctuation of the current in the resistor due to the thermal agitation of electrons:

$$V_{\rm jn} = \sqrt{\frac{4 \times k \times R \times T}{t}} \tag{6}$$

Where "k" is the Boltzman constant, R is the amplifier resistor value, T is the amplifier temperature in Kelvin, which for the Neptune is 319 K. Assuming that both detectors have the same amplifier, the Johnson noise error^{34,35} is:

$$\sigma_{\text{Johnson noise}} = \sqrt{\left(\frac{V_{\text{jn}} \times \text{Gain}}{V_{187}}\right)^2 + \left(\frac{V_{\text{jn}} \times \text{Gain}}{V_{185}}\right)^2} \qquad (7)$$

We can see from eqn (7) that an increase of the resistor amplifier from 10^{11} to 10^{13} Ω results into an increase of the Johnson-Nyquist noise by a factor of 10, whereas at the same time the gain increases by a factor 100, which implies that the signal/noise ratio is improved by a factor of 10 using $10^{13} \Omega$ amplifiers.²³ Finally, the total internal error, corresponding to the sum of the counting statistics noise and the Johnson noise is calculated^{34,35} as:

$$\sigma_{\rm internal\ error} = \sqrt{\sigma_{\rm Johnson\ noise}^2 + \sigma_{\rm Counting\ statistics}^2}$$
 (8)

These theoretical calculations show that for a signal intensity higher than 0.5 V, the only significant source of uncertainty is the counting statistics. In this case, the SE is the same for the $10^{11} \Omega$ and $10^{13} \Omega$ amplifiers (Fig. 3). Below 0.5 V, the role of Johnson-Nyquist noise increases for the $10^{11} \Omega$ amplifier but remains negligible for the $10^{13} \Omega$ amplifier (Fig. 3). The difference between the SE (2σ) of $10^{11} \Omega$ and $10^{13} \Omega$ amplifier is small $(<0.01^{\circ}_{00})$ for ¹⁸⁷Re intensity > 0.3 V, but significant (0.04°_{00}) for 187 Re intensity of 0.1 V, and very high (>0.10%) for intensity < 0.1 V. Our measurements of pure Re solutions agree well with theoretical predictions and confirm the lower SE at a given intensity for $10^{13} \Omega$ relative to $10^{11} \Omega$ amplifiers. No difference in the measured internal error was observed between pure Re solutions and samples processed through column chemistry. Hence, the use of $10^{13} \Omega$ amplifiers acts to significantly improve the internal error for small beam intensities, *i.e.* < 0.2 V. The internal error for 1 ng of Re, is less than 0.10% for $10^{13} \Omega$ amplifiers but higher than 0.20% for $10^{11} \Omega$ amplifiers.

We also evaluated the variability of the intermediate error as a function of the ¹⁸⁷Re beam intensity (Fig. 3B). We observe no significant difference between pure Re solutions and samples (which include column chemistry). The average values of intermediate errors increase with decreasing 187Re beam intensity. The average intermediate errors for a given ¹⁸⁷Re beam intensity are similar to internal errors, showing that the

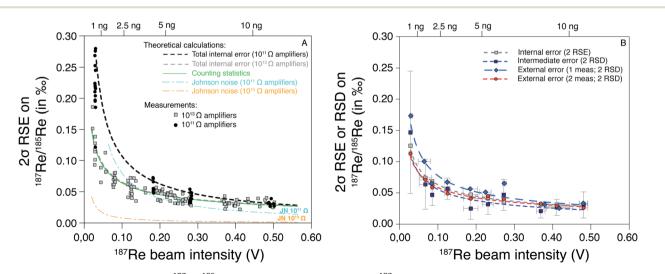


Fig. 3 Standard error and deviation on ¹⁸⁷Re/¹⁸⁵Re ratio (as 2*a*) as a function of the ¹⁸⁷Re beam intensity (in V). In (A) are the predicted trends (see main text) and data measured with 10^{11} and 10^{13} Ω amplifiers on both Re isotopes. (B) Average intermediate and external error calculated every 0.05 V for a HReO₄ solution measured over the course of three years (see Table 3). The relationship between the beam intensity (I^{157Re}) and the mass of Re required (m^{Re} , in ng) for two repeated analysis of each sample is $m^{\text{Re}} = 24.866 \times I^{187\text{Re}}$

8

		2 replicate analysis	lysis					1 analysis			
Standard solution	Range of ¹⁸⁷ Re intensity (V)	Average ¹⁸⁷ Re δ^{187} intensity (V) (%)	δ ¹⁸⁷ Re _{SRM989} (%)	Internal error (±2 SE)	$ \begin{array}{llllllllllllllllllllllllllllllllllll$	External error Number of (±2 SD) measureme	Number of measurements	Average ¹⁸⁷ Re intensity (V)	δ ¹⁸⁷ Re _{SRM989} (%)	External error Number of (±2 SD) measureme	Number of measurements
$HReO_4$	0.45 - 0.50	0.48	0.22	0.03	0.03	0.03	4	0.48	0.21	0.03	9
	0.40 - 0.45	0.41	0.23	0.03	0.03	0.03	5	0.42	0.23	0.03	7
	0.35 - 0.40	0.37	0.22	0.03	0.02	0.03	5	0.38	0.23	0.03	14
	0.25 - 0.30	0.27	0.23	0.04	0.05	0.04	8	0.27	0.23	0.07	17
	0.20 - 0.25	0.23	0.22	0.04	0.03	0.05	11	0.22	0.21	0.05	23
	0.15 - 0.20	0.19	0.20	0.05	0.03	0.04	4	0.18	0.20	0.06	7
	0.10 - 0.15	0.12	0.23	0.06	0.06	0.05	7	0.12	0.23	0.07	16
	0.08 - 0.10	0.09	0.21	0.07	0.05	0.06	5	0.08	0.21	0.07	12
	0.04 - 0.08	0.07	0.23	0.08	0.06	0.07	5	0.06	0.22	0.10	8
	0.02 - 0.04	0.03	0.21	0.13	0.15	0.11	8	0.03	0.18	0.17	14
SRM 3141	0.48 - 0.52	0.49	0.28	0.03	0.03	0.03	13	0.49	0.28	0.04	26
DURH-Re-1 1.08-1.10	1.08 - 1.10	1.09	0.44	0.02	0.02	0.02	2	1.09	0.44	0.03	10

internal errors can explain the variability between consecutive replicate analysis. To evaluate the long-term reproducibility (i.e. external error) associated to our Re isotope measurements, we repeatedly measured a pure Re solution (HReO₄) at various ¹⁸⁷Re beam intensity over the course of three years (Table 3). We calculated the external error as 2 SD of: (i) of individual analysis of HReO₄; and (ii) averages of consecutive replicate analysis of HReO₄ solutions. No significant difference is observed between the internal, intermediate and external error for replicate analysis (Fig. 3B). However, the external error is higher than the internal error on individual analysis. This shows that, on average, the internal errors can explain the variability between consecutive replicates, and between long-term repeated measurements. There is a tradeoff between running a sample once (1 analysis) at higher concentration (with issues to consider, e.g. potential for short-term blockage or instability) and running a sample multiple times (several analysis) at lower concentration. Here we choose to analyse each sample twice consecutively and use the relationship between the external error and the ¹⁸⁷Re beam intensity (red curve on the Fig. 3B) to determine the "overall uncertainty" for each sample.

Influence of W/Re ratio and intensity matching

Previous study¹⁹ doped samples to W/Re ratio of 2. However, in our case, since we measured Re isotopes at lower Re concentration with $10^{13} \Omega$ amplifiers on Re isotope cups and $10^{11} \Omega$ amplifiers on W isotope cups, we required a higher W/Re ratio. A W/Re ratio of 2 would result in an internal error on the ¹⁸⁶W/¹⁸⁴W higher than that for the ¹⁸⁷Re/¹⁸⁵Re ratio and this would be propagated onto the latter through the mass bias correction. Instead, we used a W/Re ratio of 20, so that even at very low Re concentration, the internal error on the ¹⁸⁶W/¹⁸⁴W is at least three times lower than the internal error on the ¹⁸⁷Re/¹⁸⁵Re. In addition, having a high W/Re ratio means that the contribution of residual W from the sample is lower. To ensure that this W/Re ratio does not create problems (in terms of abundance sensitivity or hydrides), we tested different W/Re ratios (1, 5, 10, 20 and 40) for two different Re concentrations (1 ppb and 5 ppb) and found no differences in precision and accuracy for W/Re ratios ranging from 5 to 40 (see ESI⁺).

Differences in the concentrations between the sample and bracketing standard can potentially influence isotopic measurements. To test this, we measured pure rhenium solutions prepared with distinct Re/W ratio relative to the Re/W of the bracketing standard. We observe that δ^{187} Re values are offset by more than 0.05% for a ratio mismatch of more than 40% (Fig. 4). This shows that is crucial to match within 10% the Re/W of the measured sample with that of the bracketing standard. The reasons for this effect may be due to different contributions of hydrides or W backgrounds.

Matrix effects and accuracy of the measurements

Residual inorganic and/or organic matrices following column separation can impact the accuracy of isotope measurement.^{27,29} We compared the measured δ^{187} Re after one, two, three and four repeated column separations (Fig. 5). This test illustrates

This article is licensed under a Creative Commons Attribution-NonCommercial 3.0 Unported Licence.

Dpen Access Article. Published on 07 January 2020. Downloaded on 1/8/2020 9:42:56 AM

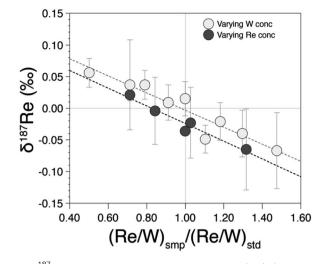


Fig. 4 δ^{187} Re as a function of the ratio between the (Re/W) ratio of the measured test standard relative to the (Re/W) ratio of the bracketing standard.

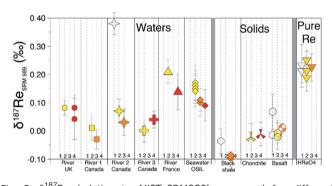


Fig. 5 δ^{187} Re (relative to NIST-SRM989) measured for different samples after repeated column separations (1 column pass = white, 2 = yellow, 3 = orange, 4 = red).

that the measured δ^{187} Re can change as a function of number of column passes but how it does so is somewhat dependent on sample type. For BCR-2 (basalt), the δ^{187} Re values are lower for two column passes relative to one, but similar for two and three column passes. For a chondrite (Allende), we also observe no difference in δ^{187} Re between the second and fourth column passes on the same aliquot. For a seawater reference material (OSIL Atlantic seawater standard), the δ^{187} Re for three and four column passes are identical, but lower than for two columns passes. For one river water (River 2 Canada), the δ^{187} Re is significantly higher for 1 column pass relative to the value measured after 2 or more passes. For others (River 1 and 2 Canada), there is a small decrease of the δ^{187} Re value after three passes relative to two, although this is not observed for every river sample (*e.g.* River UK and River 3 Canada).

Overall, these tests reveal two things: (i) in general, the δ^{187} Re values are unchanged after three column chemistry separations; (ii) the materials which continue to show a change in δ^{187} Re after more than two repeated separations are materials that originally contain organic matter (river water and

shales). To test whether the amount of remaining inorganic matrix can affect the accuracy of δ^{187} Re measurement, we performed a series of tests where a Re standard solution, doped with various elements (Al, Fe, Mg, Nb, Zr, Ta, Mo, Zn, U, Hf) at different concentrations, covering a range of value of X/Re (where X is the doped element), was measured. We observe that over the range of concentrations of dopant tested, there is no effect, within uncertainty, on the accuracy of the δ^{187} Re value (see ESI[†]). In addition, we used a granite SRM (USGS "G-2") with a very low Re concentration (11 ppt) that we doped with a known mass of ¹⁸⁷Re-rich standard solution DURH-1 (+0.44 \pm 0.02%) Table 3) so that the contribution of Re from the sample is negligible (<0.6%, see ESI[†]). This modified granite SRM was passed through column chemistry (3 times) and measured as a normal sample and returned the same value as pure DURH-1 (Fig. 6). These tests show that our column procedure is optimal for removing effects caused by an inorganic matrix.

Organics that influence isotopic measurement can originate from: (i) the resin or (ii) the sample. Matrix effects caused by resin degradation have been documented for other isotope systems.^{37,38} In this study, we systematically treat the column residue with concentrated HNO₃ at 130–150 °C for 24 h before isotope measurement. Measurements of pure Re solutions processed and non-processed through column chemistry return identical δ^{187} Re values (Table S4, ESI†) suggesting that any residual organic molecules derived from the resin do not affect the measurement. We note though that for some of these tests, abrupt changes in the ¹⁸⁶W/¹⁸⁴W mass bias are observed, leading to inaccurate δ^{187} Re^{SSB} and (δ^{187} Re^{SSB}– δ^{187} Re) values up to +0.35% but correct δ^{187} Re. This shows that matrix effects caused by residual resin-derived organics are adequately corrected for using W normalization.

Resin AG1-X8 is a cationic resin for which organic molecules (negatively charged) have a high affinity. Adsorption of organic molecules to the resin is evidenced by a change in the colour (orange to dark) of the resin. Some of the organics are eluted in

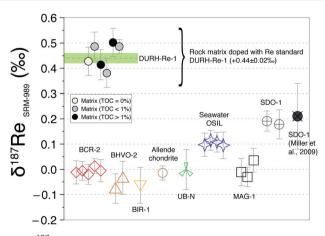


Fig. 6 δ^{187} Re measured for different SRMs and matrix tests doped with standard DURH-Re-1. Each point corresponds to one measurement (*i.e.* one digestion). The green line is the average δ^{187} Re of DURH-Re-1. The error bar for each measurement corresponds to the overall uncertainty (which is a function of the ¹⁸⁷Re beam intensity).

Paper

4 M HNO₃, at the same time as Re, observed as a yellowish colour of the solution and the size/colour of the residue. These organics could potentially create interferences on the Re isotopic measurement and influence the elution volume of Re (*e.g.* competition for binding sites on the resin). To assess the role of organic matter we used two samples (New Zealand soils with 1 and 6% organic carbon) that have a very low Re concentration (<35 ppt)⁷ and several shale-type Re-free matrix (by mixing volume cuts before and after Re elution peak, see ESI†) and doped them with Re solution DURH-Re-1. Following three column chemistry separations, we measured δ^{187} Re respectively which are in agreement within uncertainty with the value of the DURH-Re-1 standard (Fig. 6). Collectively, this shows that our column procedure successfully remove effects caused by sample-derived organic matrix.

Application to reference materials and assessment of data quality

The δ^{187} Re values of several reference materials (BCR-2, BHVO-2, UB-N, BIR-1, MAG-1, SDO-1, OSIL Atlantic Seawater) have been determined (Table 4). The variability in δ^{187} Re values defined by these standards is about 0.28% (Fig. 6). The intermediate and external error of the SRM are comparable to those of pure Re solutions, indicating that digestion and column chemistry does not introduce additional uncertainty on the δ^{187} Re measurement relative to pure Re solution. Since there are very few reference materials that have been measured for δ^{187} Re, it is of prime importance to assess the accuracy of the measurements. Our analysis of NIST 3141 and HReO₄ return δ^{187} Re values respectively of 0.28 \pm 0.03‰ and 0.22 \pm 0.03‰, in agreement with published values.¹⁹ The δ^{187} Re value we obtain for the SDO-1 standard, δ^{187} Re = 0.19 \pm 0.03%, also matches previously published value.19 These results are encouraging for accurate analysis, however in future a wider range of materials will need to be crosscalibrated between laboratories for a more thorough assessment.

Perspectives for the measurement of $\delta^{187} \mbox{Re}$ in geological samples

The first Re isotope study¹⁹ achieved precision of $\pm 0.10\%$ for 10 ng of Re. From this, the authors concluded that Re isotope

 Table 4
 Results of standard reference materials measured for Re isotopes

measurements would be limited only to Re-rich samples (such as black shales or sulfide ores). Here, our external reproducibility is $\pm 0.04\%$ for 5 ng of Re (Fig. 3) based on the use of 10^{13} Ω amplifiers in wet plasma conditions (micro-cyclonic chamber). There are potential alternative methods that could be used to deliver this improvement, or further refine it. A desolvator (Aridus, Apex) combined with a jet sample cone and X-skimmer cone could lead to a further significant increase of the sensitivity for Re. However, there are concerns regarding the Re beam stability in dry plasma conditions and the possibility of Re isotope fractionation on the desolvator membrane.³⁹ Future studies should assess whether these additional changes to the sample introduction and/ or interface could allow further reduction in the volume of sample required for making precise Re isotopic measurements.

Irrespective of any refinements to the sample introduction and interface, the modified column chemistry and analytical advances we have made allow a large range of geological materials to be measured at high precision (Fig. 7). For instance, for most organic rich shales, we estimate that only ~0.1 g of sample is required for a precision of $\pm 0.04\%$. For igneous rock samples

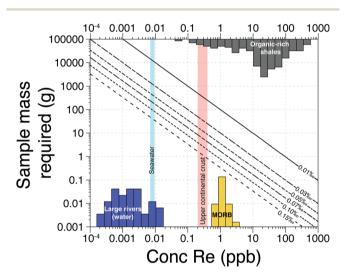


Fig. 7 Sample mass required (in g) as a function of the Re concentration of geological materials for different measurement precisions. Rhenium concentration data from published work.^{1,2,14,21}

	δ^{187} Re _{SRM989}	Overall uncertainty	External error	Number of	Intermedia error $(\pm 2$	error	Mass Re per analysis
Reference material	(‰)	$(\pm 2 \text{ SD})$	$(\pm 2 \text{ SD})$	measurements	SD)	$(\pm 2 \text{ SE})$	(ng)
OSIL (Atlantic seawater)	0.10	0.04	0.02	4	0.05	0.05	3.9
MAG-1 (marine mud)	0.00	0.05	0.07	3	0.05	0.05	4.2 - 5.5
BHVO-2 (basalt)	-0.06	0.06	0.07	2		0.06	2-2.5
BIR-1(basalt)	-0.06	0.07		1		0.06	0.9
BCR-2 (basalt)	-0.01	0.04	0.02	5	0.05	0.04	4.9-12.2
CV3 Allende (chondrite)	-0.01	0.03		1	0.01	0.04	21
SDO-1 (black shale)	0.19	0.05	0.03	2	0.01		3-12
UB-N (serpetinite)	0.00	0.08		1	0.04	0.08	1.3

(e.g. Mid-Ocean Ridge Basalts, MORB), we estimate that \sim 1–5 g would be required for a precision of $\pm 0.04\%_{oo}$. For river waters, a large volume of sample (~1 to 20 L) is still required because most rivers have dissolved Re concentrations in the ppt range (Fig. 7). Our field pre-concentration method addresses this. Instead of transporting many kilograms of filtered water, we are able to concentrate \sim 5–10 ng of Re onto a few mL of resin and subsequently recover Re for precise isotopic analysis. The required precision for any isotopic analysis will depend on the range of δ^{187} Re values expected in the samples undergoing analysis. So far, the range of measured δ^{187} Re in geological materials is $\sim 0.80\%$, a range much larger than the analytical precision. Along with the method developments described here which make the vast majority of natural samples available for analysis (Fig. 7), these calculations provide impetus for further exploration of the Re isotope system.

Conclusions

Here we provide new methods to measure the isotopic composition of Re in a wide range of solid and liquid samples. The main points are:

- An improved precision via MC-ICP-MS afforded by the use of $10^{13}~\Omega$ amplifiers on ^{185}Re and $^{187}\text{Re}.$

- The use of W doping, matching W/Re ratios of samples and standards, and standard-sample bracketing allows for the correction of instrumental mass bias and some matrix affects.

- The importance of column chemistry separation, and here we provide evidence that three repeated separations are necessary.

- The use of a field-based pre-concentration of Re for river water samples, which makes the analyses of river waters feasible.

Combining these approaches, we are able to report an improved precision on $\delta^{187}Re$ measurements for a lower required mass of Re $(\pm 0.10\%$ for a mass of Re of ~ 1 ng to $\pm 0.03\%$ for a mass of Re of >10 ng) compared to previous studies. 19,20 We have analysed a range of solid geological reference materials, and a seawater standard, and observe variability of 0.28%. The approaches here open up studies of Re isotopic fractionation, in particular opening up the potential of Re isotopes as a proxy for oxidative weathering, and for reconstruction of past redox cycling.

Conflicts of interest

There are no conflicts to declare.

Acknowledgements

This research was funded by the European Research Council (ERC) Starting Grant "ROC-CO₂" (project 678779) awarded to Robert Hilton and a COFUND Junior Research Fellowship supported by the European Union's Seventh Framework Programme to Mathieu Dellinger. We thank: Bernard Peucker-Ehrenbrink and Alex Dickson for sharing Re isotope standards; Martin West for laboratory assistance; Josh West, Kevin

Burton, Julie Prytulak, Ambre Luguet for discussions prior to submission.

Notes and references

- 1 C. A. Miller, B. Peucker-Ehrenbrink, B. D. Walker and F. Marcantonio, *Geochim. Cosmochim. Acta*, 2011, 75, 7146–7179.
- 2 D. Colodner, J. Sachs, G. Ravizza, K. Turekian, J. Edmond and E. Boyle, *Earth Planet. Sci. Lett.*, 1993, **117**, 205–221.
- 3 A. D. Anbar, R. A. Creaser, D. A. Papanastassiou and G. J. Wasserburg, *Geochim. Cosmochim. Acta*, 1992, **56**, 4099–4103.
- 4 J. Golden, M. McMillan, R. T. Downs, G. Hystad, I. Goldstein, H. J. Stein, A. Zimmerman, D. A. Sverjensky, J. T. Armstrong and R. M. Hazen, *Earth Planet. Sci. Lett.*, 2013, 366, 1–5.
- 5 T. K. Dalai, S. K. Singh, J. R. Trivedi and S. Krishnaswami, Geochim. Cosmochim. Acta, 2002, 66, 29–43.
- 6 R. G. Hilton, J. Gaillardet, D. Calmels and J.-L. Birck, *Earth Planet. Sci. Lett.*, 2014, **403**, 27–36.
- 7 K. Horan, R. G. Hilton, D. Selby, C. J. Ottley, D. R. Gröcke, M. Hicks and K. W. Burton, *Sci. Adv.*, 2017, 3, e1701107.
- 8 A. I. Sheen, B. Kendall, C. T. Reinhard, R. A. Creaser, T. W. Lyons, A. Bekker, S. W. Poulton and A. D. Anbar, *Geochim. Cosmochim. Acta*, 2018, 227, 75–95.
- 9 C. A. Miller, B. Peucker-Ehrenbrink and E. A. Schauble, *Earth Planet. Sci. Lett.*, 2015, **430**, 339–348.
- 10 J. Crusius, S. Calvert, T. Pedersen and D. Sage, *Earth Planet. Sci. Lett.*, 1996, 145, 65–78.
- 11 J. L. Morford, W. R. Martin and C. M. Carney, *Chem. Geol.*, 2012, **324–325**, 73–86.
- 12 D. Selby, R. A. Creaser and M. G. Fowler, *Geochim. Cosmochim. Acta*, 2007, **71**, 378–386.
- 13 F. Mahdaoui, L. Reisberg, R. Michels, Y. Hautevelle, Y. Poirier and J.-P. Girard, *Chem. Geol.*, 2013, **358**, 90–100.
- 14 A. Dubin and B. Peucker-Ehrenbrink, *Chem. Geol.*, 2015, **403**, 111–120.
- 15 R. J. Walker, Geochemical Perspectives, 2016, 5, 1–2.
- 16 K. Lodders, Astrophys. J., 2003, 591, 1220.
- 17 E. H. Hauri and S. R. Hart, Chem. Geol., 1997, 139, 185-205.
- 18 J. W. Morgan, J. Geophys. Res.: Solid Earth, 1986, 91, 12375– 12387.
- 19 C. A. Miller, B. Peucker-Ehrenbrink and L. Ball, *J. Anal. At. Spectrom.*, 2009, **24**, 1069–1078.
- 20 R. Liu, L. Hu and M. Humayun, *Meteorit. Planet. Sci.*, 2017, 52, 479–492.
- 21 A. Gannoun, K. W. Burton, I. J. Parkinson, O. Alard,
 P. Schiano and L. E. Thomas, *Earth Planet. Sci. Lett.*, 2007,
 259, 541–556.
- 22 L. A. Jaffe, B. Peucker-Ehrenbrink and S. T. Petsch, *Earth Planet. Sci. Lett.*, 2002, **198**, 339–353.
- 23 T. Breton, N. S. Lloyd, A. Trinquier, C. Bouman and J. B. Schwieters, *Procedia Earth Planet. Sci.*, 2015, 13, 240– 243.
- 24 Z. Chu, Y. Yan, Z. Chen, J. Guo, Y. Yang, C. Li and Y. Zhang, *Geostand. Geoanal. Res.*, 2015, **39**, 151–169.

- 25 J. P. Faris and R. F. Buchanan, *Anal. Chem.*, 1964, **36**, 1157–1158.
- 26 E. L. Gendre, E. Martin, B. Villemant, P. Cartigny and N. Assayag, *Rapid Commun. Mass Spectrom.*, 2017, **31**, 137– 144.
- 27 F. Albarède and B. Beard, *Rev. Mineral. Geochem.*, 2004, 55, 113–152.
- 28 J. Day, D. G. Pearson and G. M. Nowell, in *Plasma Source Mass Spectrometry: applications and emerging technologies*, 2003, pp. 374–390.
- 29 A. Poirier and R. Doucelance, *Geostand. Geoanal. Res.*, 2009, 33, 195–204.
- 30 T. Meisel and J. Moser, *Geostand. Geoanal. Res.*, 2004, 28, 233–250.
- 31 T. Meisel and M. F. Horan, *Rev. Mineral. Geochem.*, 2016, **81**, 89–106.

- 32 K. P. Jochum, U. Weis, B. Schwager, B. Stoll, S. A. Wilson,
 G. H. Haug, M. O. Andreae and J. Enzweiler, *Geostand. Geoanal. Res.*, 2016, 40, 333–350.
- 33 N. Imai, S. Terashima, S. Itoh and A. Ando, *Geostand. Newsl.*, 1995, **19**(2), 135–213.
- 34 S. G. John and J. F. Adkins, Mar. Chem., 2010, 119, 65-76.
- 35 G. Paris, A. L. Sessions, A. V. Subhas and J. F. Adkins, *Chem. Geol.*, 2013, 345, 50–61.
- 36 J. Liu and D. G. Pearson, Chem. Geol., 2014, 363, 301-311.
- 37 A. J. Pietruszka and A. D. Reznik, *Int. J. Mass Spectrom.*, 2008, 270, 23–30.
- 38 F. Kurzweil, C. Münker, J. Tusch and R. Schoenberg, *Chem. Geol.*, 2018, **476**, 407–417.
- 39 M. Pfeifer, N. S. Lloyd, S. T. M. Peters, F. Wombacher, B. M. Elfers, T. Schulz and C. Münker, J. Anal. At. Spectrom., 2017, 32, 130–143.