

Aalborg Universitet

Long-term storage of clinical samples in CyMol® medium for PNA- FISH® and culturing from the eSwab™ system

Larsen, Lone Heimann; Xu, Yijuan; Pedersen, Malene Schibler; Schønheyder, Henrik Carl; Thomsen, Trine Rolighed

Publication date: 2014

Document Version Early version, also known as pre-print

Link to publication from Aalborg University

Citation for published version (APA):

Larsen, L. H., Xu, Y., Pedersen, M. S., Schønheyder, H. C., & Thomsen, T. R. (2014). Long-term storage of clinical samples in CyMol® medium for PNA- FISH® and culturing from the eSwab™ system. Poster session presented at 24th European Congress of Clinical Microbiology and Infectious Diseases, ECCMID 2014, Barcelona, Spain.

General rights

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- ? Users may download and print one copy of any publication from the public portal for the purpose of private study or research. ? You may not further distribute the material or use it for any profit-making activity or commercial gain ? You may freely distribute the URL identifying the publication in the public portal ?

Take down policy

If you believe that this document breaches copyright please contact us at vbn@aub.aau.dk providing details, and we will remove access to the work immediately and investigate your claim.



Long-term storage of clinical samples in CyMol® medium for PNA- FISH® and culturing from the eSwabTM system

preserve samples at -80°C for study by PNA-FISH® for at least 12

months and we expect storage for even longer periods to be feasible.

We estimated the effective detection limit to be in the order of >10²

CFU/mL (updated from the abstract). Both the morphology and

The eSwab™ was a convenient system for documenting a broad

intensity of staining with nucleic acid and PNA probes were distinct.

range of bacterial pathogens associated with foreign body infections.

Lone Heimann Larsen^{1,2}, Yijuan Xu², Malene S. Pedersen², <u>Henrik C. Schønheyder^{1,3}</u>, Trine Rolighed Thomsen^{2,4} and the PRIS study group

¹Department of Clinical Microbiology, Aalborg Univervity Hospital, Aalborg, Denmark. ²Center for Microbial Communities, Department of Biotechnology, Chemistry and Environmantal Engineering, Aalborg University, Aalborg, Denmark. ³Department of Clinical Medicin, Aalborg University, Aalborg, Denmark. ⁴Danish Technological Inistitue, Aarhus, Denmark. Contact information: **lone.l@rn.dk**

Objectives: The diversity of bacteria reported in foreign body infections is steadily growing, and culture-independent methods have become a valuable supplement to established culture methods. Therefore, sampling and preservation of specimens for molecular analysis have become an important issue. We report here experience gained with different specimen types obtained in parallel both preand peroperatively in patients enrolled prospectively in a clinical study of prosthesis-related problems (www.joint-prosthesis-infection-pain.dk).

Methods: Parallel sampling for both culture-dependent and independent analyses were done over a period of two years. Specimens included tissue biopsies, bone biopsies, joint fluid, and eSwabsTM (Copan, Italy) taken from the prosthesis in situ, as well as sonication fluid (if prosthesis components were removed). Transfer to the laboratory was direct (ambient temperature) and storage temperature for culture specimens was +4°C (max 24 h), and -80°C for specimens for culture-independent methods (until batchwise analysis). Specimens for peptide nucleic acid-fluorescence in situ hybridization (PNA-FISH) analysis were stored in CyMol® medium (Copan, Italy) for up to one year. Direct visualization of microorganisms followed a previously published PNA-FISH® (AdvanDx, USA) protocol except for a substitution of the fixation step with filtration of 200 μ L sample through a 0.22 μ m white polycarbonate filter (prewashed with dH₂O, GE Water & Process Technologies, USA) in order to fix and concentrate the sample.

Results: With broad range and specific PNA-FISH® probes, we demonstrated bacteria with a bright signal and morphology comparable to the isolates obtained by culture of parallel specimens within 24 h. The detection limit for PNA-FISH® were >10^2 CFU/mL (estimated from the colony forming units cultured, updated from abstract). With the eSwab TM system we detected a broad range of Gram-positive taxa including Staphylococcus spp., Streptococcus spp., Enterococcus spp., and Corynebacterium spp. by culture and 16S rRNA gene amplicon sequencing (table 1).

Direct culturing 16S rDNA sequencing eSwabTM stored in modified Amies modified Amies medium with 20% glycerol at -80 °C Staphylococcus aureus Staphylococcus epidermidis Staphylococcus lugdunensis Staphylococcus caprae Enterococcus faecalis Streptococcus dysgalactiae

Streptococcus agalactiae

Corynebacterium jeikeium

Corynebacterium striatum

Escherichia coli

Propionibacterium acnes

Table 1: Bacterial species recovered by $eSwab^{TM}$ (Copan, Italia). $eSwab^{TM}$ used for culturing were transported in modified Amies medium and cultured within 24 h. eSwabsTM for sequencing were stored at -80 °C in modified Amies medium with 20% glycerol until analysed batchwise.


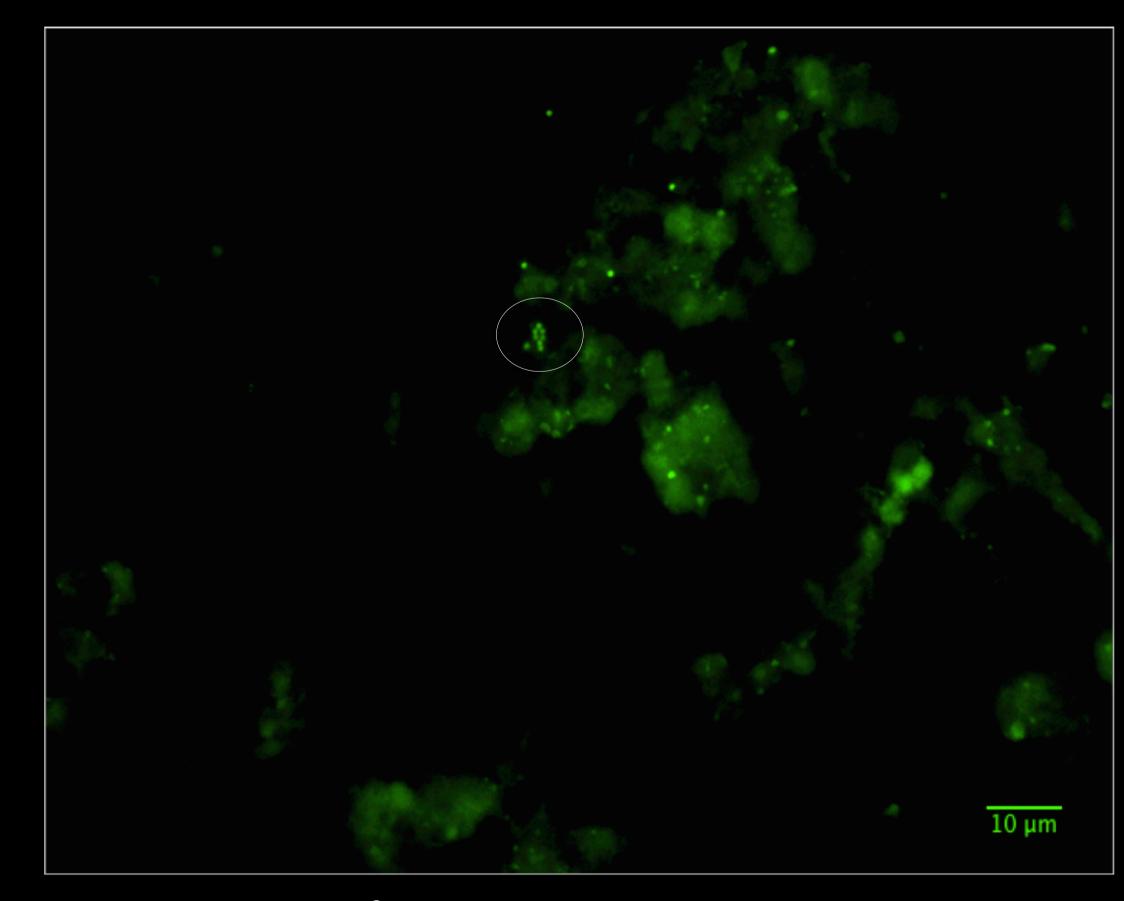


Figure 2: Filter PNA-FISH® on sonication fluid from knee prosthetic component. Staphylococcus aureus; green. The culture and sequencing were positive for S. aureus. The PNA-FISH confirmed the culturing and sequencing results illustrating green bacteria in clusters (eg. circled).

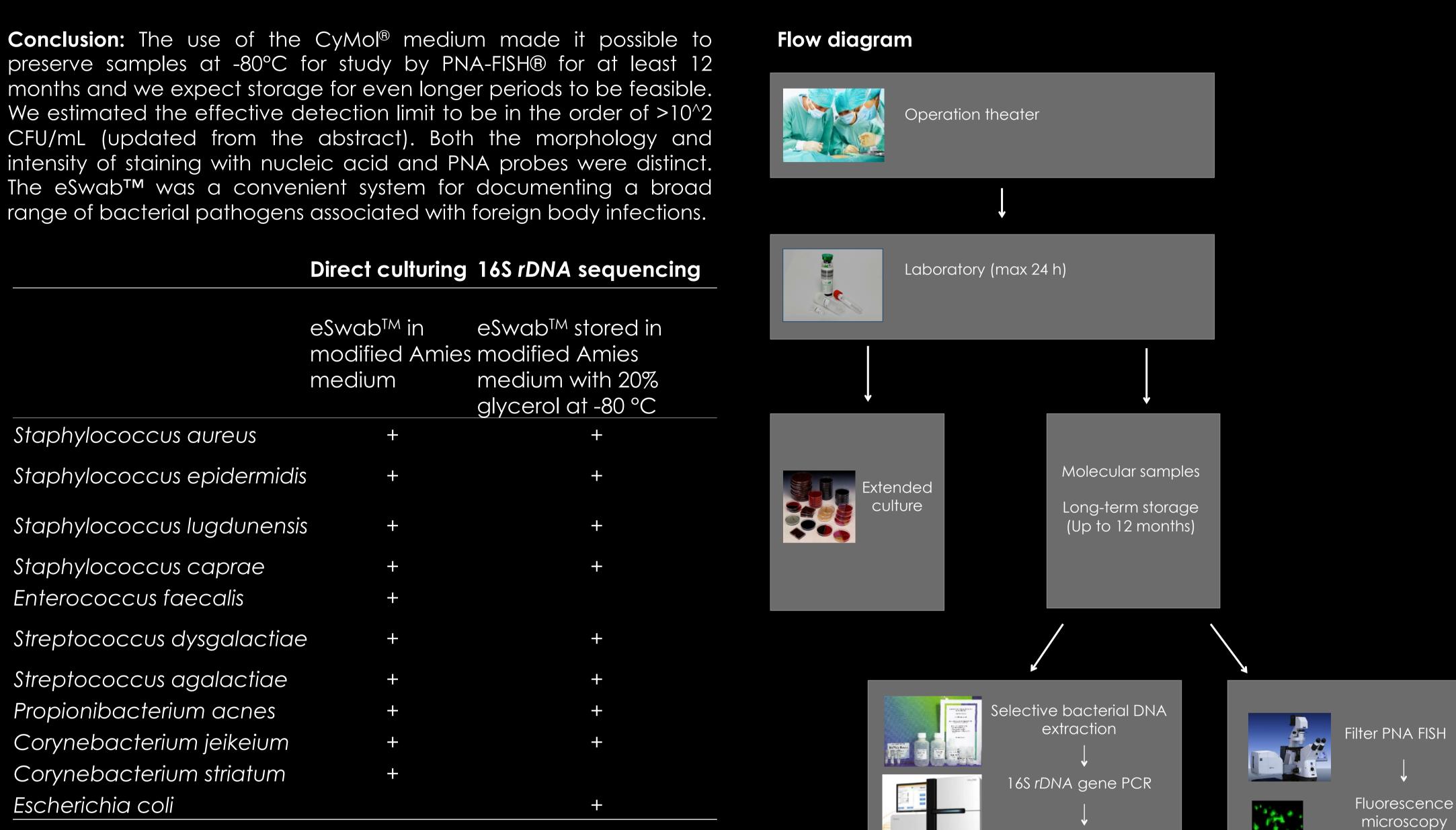

Amplicon sequencing

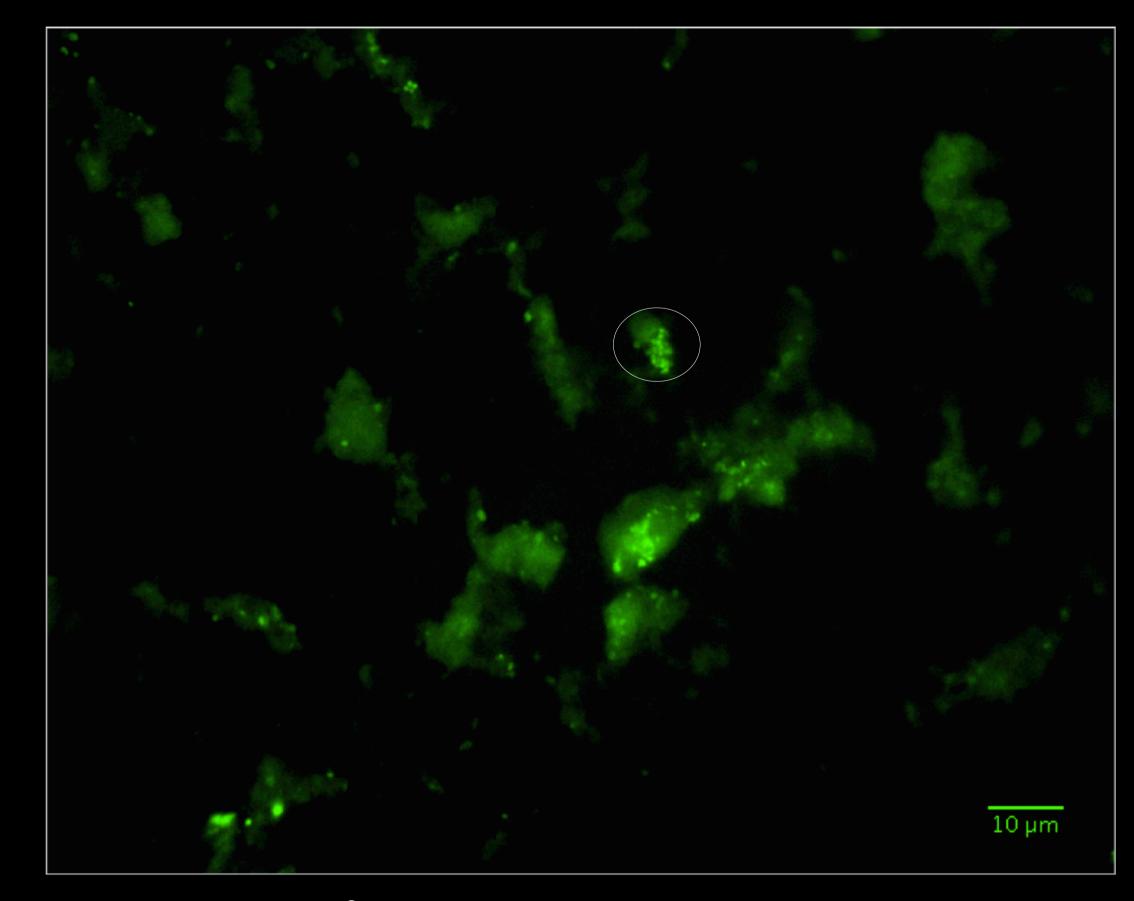


Figure 3: Filter PNA-FISH® on sonication fluid from knee prosthetic component (same case as figure 2). Staphylococcus aureus; green. The culture and sequencing were positive for S. aureus. The PNA-FISH confirmed the culturing and sequencing results illustrating green bacteria in clusters (eg. circled).

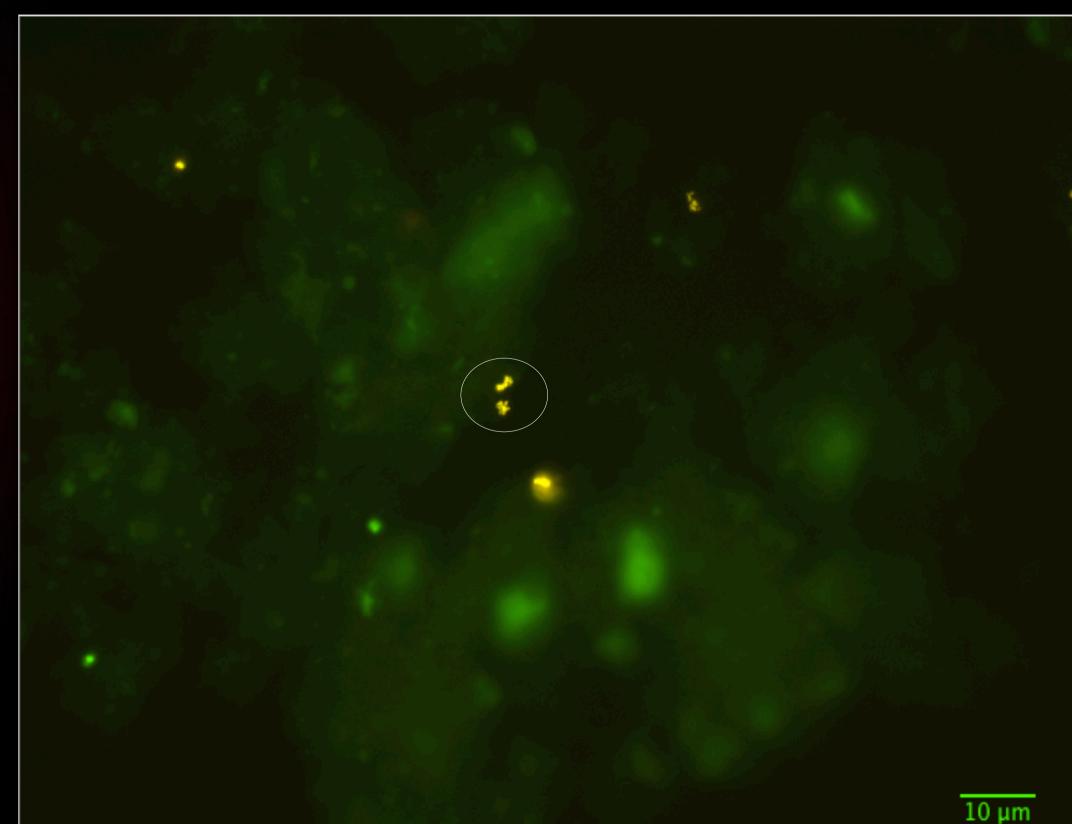


Figure 1: Filter PNA-FISH® on sonication fluid from hip prosthetic component. Universal bacterial probe (BacUni; green, Enterococcus faecalis; red. The culture and sequencing were positive for E. faecalis. The PNA-FISH confirmed the culturing and sequencing results illustrating yellow bacteria in clusters (eg. circled).



DANISH

INSTITUTE

TECHNOLOGICAL



AALBORG UNIVERSITY

DENMARK

