



**University of
Zurich**^{UZH}

**Zurich Open Repository and
Archive**

University of Zurich
University Library
Strickhofstrasse 39
CH-8057 Zurich
www.zora.uzh.ch

Year: 2024

Illumination matters part IV: blackout and whiteout in flexible ureteroscopy - first report on a phenomenon observed by PEARLS

Kwok, Jia-Lun ; De Coninck, Vincent ; Panthier, Frédéric ; Kamkoum, Hatem ; Pauchard, Felipe ; Shrestha, Anil ; Gauhar, Vineet ; Brachlow, Jan ; Schmid, Florian Alexander ; Poyet, Cédric ; Eberli, Daniel ; Traxer, Olivier ; Keller, Etienne Xavier

DOI: <https://doi.org/10.1007/s00345-024-04988-1>

Posted at the Zurich Open Repository and Archive, University of Zurich

ZORA URL: <https://doi.org/10.5167/uzh-262488>

Journal Article

Published Version



The following work is licensed under a Creative Commons: Attribution 4.0 International (CC BY 4.0) License.

Originally published at:

Kwok, Jia-Lun; De Coninck, Vincent; Panthier, Frédéric; Kamkoum, Hatem; Pauchard, Felipe; Shrestha, Anil; Gauhar, Vineet; Brachlow, Jan; Schmid, Florian Alexander; Poyet, Cédric; Eberli, Daniel; Traxer, Olivier; Keller, Etienne Xavier (2024). Illumination matters part IV: blackout and whiteout in flexible ureteroscopy - first report on a phenomenon observed by PEARLS. *World Journal of Urology*, 42(1):294.

DOI: <https://doi.org/10.1007/s00345-024-04988-1>



Illumination matters part IV: blackout and whiteout in flexible ureteroscopy - first report on a phenomenon observed by PEARLS

Jia-Lun Kwok^{1,2} · Vincent De Coninck^{3,4,5} · Frédéric Panthier^{3,6} · Hatem Kamkour^{3,7} · Felipe Pauchard^{3,8} · Anil Shrestha^{3,9} · Vineet Gauhar^{3,10} · Jan Brachlow^{3,11} · Florian Alexander Schmid¹ · Cédric Poyet¹ · Daniel Eberli¹ · Olivier Traxer^{3,6} · Etienne Xavier Keller^{1,3,4}

Received: 7 February 2024 / Accepted: 7 April 2024
© The Author(s) 2024

Abstract

Purpose To date, no study has evaluated effects of varying brightness settings on image quality from flexible ureteroscopes submerged in saline. The aim was to evaluate blackout and whiteout occurrences in an in-vitro kidney calyx model.

Material and methods We evaluated a series of contemporary flexible ureteroscopes including the Storz Flex-Xc and Flex-X2s, Olympus V3 and P7, Pusen 7.5F and 9.2F, as well as OTU WiScope using a 3D-printed enclosed pink in-vitro kidney calyx model submerged in saline. Endoscopic images were captured with ureteroscope tip placed at 5 mm, 10 mm and 20 mm distances. The complete range of brightness settings and video capture modes were evaluated for each scope. Distribution of brightness on a grayscale histogram of images was analyzed (scale range 0 to 255). Blackout and whiteout were defined as median histogram ranges from 0 to 35 and 220 to 255, respectively (monitor image too dark or too bright for the human eye, respectively).

Results Blackout occurred with the P7, Pusen 7.5F, 9.2F and WiScope at all distances, and V3 at 20 mm - with lowest brightness settings. Whiteout occurred with Flex-X2s, V3 and P7 at 5 mm and 10 mm, as well as with V3 and P7 at 20 mm - mostly with highest brightness settings. The Flex-Xc had neither blackout nor whiteout at all settings and distances.

Conclusion Blackout or whiteout of images is an undesirable property that was found for several scopes, possibly impacting diagnostic and therapeutic purposes during ureteroscopy. These observations form a guide to impact a urologist's choice of instruments and settings.

Keywords Flexible ureteroscopy · Blackout · Whiteout · Stone disease · Image

✉ Etienne Xavier Keller
etiennexavier.keller@usz.ch

¹ Department of Urology, University Hospital Zurich, Frauenklinikstrasse 10, CH-8091 Zurich, Switzerland

² Department of Urology, Tan Tock Seng Hospital, Singapore, Singapore

³ Progressive Endourological Association for Research and Leading Solutions (PEARLS), Paris, France

⁴ Endourology & Urolithiasis Working Group, Young Academic Urologists (YAU), Arnhem, The Netherlands

⁵ Department of Urology, AZ Klina, Brasschaat, Belgium

⁶ GRC N°20, Groupe de Recherche Clinique sur la Lithiase Urinaire, Hôpital Tenon, Sorbonne Université, 75020 Paris, France

⁷ Hamad Medical Corporation, Doha, Qatar

⁸ Urology Department, Hospital Naval Almirante Nef, 2520000, Viña del Mar, Chile

⁹ Department of Urology, National Academy of Medical Sciences, Bir Hospital and B&B Hospital, Gwarko Lalitpur, Nepal

¹⁰ Department of Urology, Ng Teng Fong General Hospital, Singapore, Singapore

¹¹ Zentrum Für Urologie Winterthur, Winterthur, Switzerland

Introduction

Ureteroscopy treatment has been increasing in various countries worldwide over the last 2 decades [1–3]. Quality of optics is an important part of evaluation of flexible ureteroscopes [4]. Imaging properties of flexible ureteroscopes have been extensively evaluated - mostly in air [4], and recently in saline [5]. Visibility has been prior evaluated [4], but no study to date has evaluated the effects of varying brightness settings on image quality from ureteroscopes submerged in saline - the medium used for most endoscopic interventions [6–23]. The authors have noticed that at low brightness settings for some scopes, there is blackout - where the endoscopic image is too dark to visualize the endoscopic field well. At higher brightness settings, the authors have noticed whiteout - where the endoscopic image is too bright to visualize the endoscopic field well. This phenomenon has not been objectively evaluated and reported before. The aim of the present study was to evaluate blackout and whiteout occurrences in an *in-vitro* kidney calyx model with various brightness settings and brightness modes.

Material and methods

A range of ureteroscopes accessible at University Hospital Zurich were evaluated, including the Flex-Xc and Flex-X2s (Karl Storz SE & Co. KG, Tuttlingen, Germany), the URF-P7 and URF-V3 (Olympus, Center Valley, PA, USA), and single-use scopes Uscope 7.5F PU3033A, Uscope 9.2F PU3022A (Zhuhai Pusen Medical Technology Co. Ltd. Guangdong, China) and WiScope (OTU Medical Inc, CA, USA). To simulate routine clinical conditions, the single-use scopes (Pusen 7.5F, Pusen 9.2F and OTU WiScope) were unused from sealed sterile packages. The reusable scopes (Storz and Olympus scopes) had all been processed and decontaminated post clinical use, with no record of prior usage numbers.

For the Storz Flex-X2s, the Power LED 175 light source (unit usage < 100 h) was utilized with a corresponding 230 cm and 3.5 mm fiberoptic cable, together with the IMAGE1 S HX-P HDTV 1-Chip pendular camera (Karl Storz SE & Co. KG, Tuttlingen, Germany). With the Olympus URF-P7 and URF-V3, the VISERA elite CLV-S190 light source (Xenon short-arc lamp used < 100 h) was paired with a WA03310A 300 cm and 4.3 mm fiberoptic light cable, and the CH S190 08 LB camera head (Olympus, Center Valley, PA, USA) were used. Fiberoptic cables used were entirely new.

A 3D-printed pink kidney calyx model was used to hold the tested ureteroscope at fixed distances from the concave

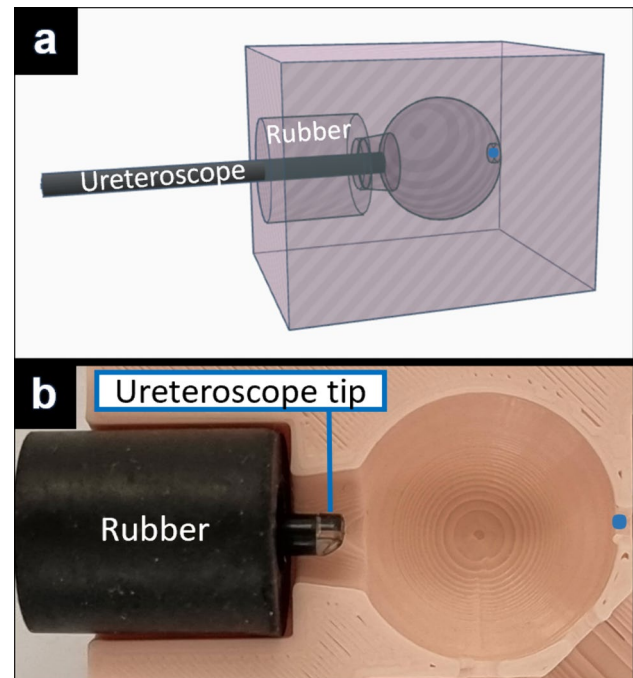


Fig. 1 Experimental setup **a** 3D modelling of the experimental setup with a ureteroscope inserted into a pink kidney calyx model with a spherical cavity. **b** Cross section of the actual 3D-printed pink model cut in half to display the position of the ureteroscope in relation to the spherical cavity. The black rubber cylinder was used to ensure solid tightening of the ureteroscopes within the model during measurements. The tip of the OTU WiScope is visible and looking out 4 mm from the black rubber for a 20 mm distance between tip of the scope and the center of the concave cavity surface (marked with a blue dot on this figure). The scope was inserted further for capturing images at 5, 10 mm distances (not shown in this figure)

surface of the kidney calyx model in a dark room (Fig. 1). The model consisted of a closed spherical cavity to replicate the human kidney calyx. Pink was chosen for the kidney calyx model to simulate the color of human urothelial mucosa. The size of the target field and distance from the light sensor was determined with reference to dimensions of models constructed on data from endocasts [24] used in studies investigating laser lithotripsy [25, 26], to reflect *in vivo* settings. The ureteroscope was maintained in a non-deflected position, with the center of the scope view aligned to the center of the opposite concave pink surface. Endoscopic images were captured with the tip of the ureteroscope at 5 mm, 10 mm and 20 mm distances from the inner surface of the sphere - mimicking situations frequently found in clinical routine. All experiments were performed in saline (NaCl 0.9%) to replicate the usual operative medium in ureteroscopy.

The complete range of brightness settings and brightness modes were evaluated for each scope - where available. These included: First, video image brightness setting (Storz Flex-Xc and Flex-X2s). Second, light stack or processing

unit brightness setting (all scopes, except the Storz Flex-Xc). Lastly, brightness mode - Manual vs Auto (Storz Flex-X2s, Olympus V3 and P7). Depending on the scope, some of these adjustment options were not present. This is further described:

For the Storz Flex-Xc, brightness adjustment was available only via the video image brightness menu buttons on the scope handle, and not on the light stack. This was tested at settings of 1/3/5 out of 5. This brightness adjustment function was found not to affect the actual ureteroscope tip light source brightness in previous studies, but only of the projected video image [5, 27, 28]. Similarly, the Storz Flex-X2s had buttons on the camera head allowing for separate brightness adjustment of the video image (tested at settings 1/3/5 out of 5). Assessment of the Flex-Xc and Flex-X2s was done in “Standard” mode with no additional post processing provided by the IMAGE1 STORZ Professional Image Enhancement System (e.g. CLARA/CHROMA).

For all scopes except the Storz Flex-Xc, brightness settings could be changed on the light stack or processing unit, and was tested at the lowest, middle and highest of the available range in each brightness mode and video image brightness setting.

Additionally, the Storz Flex-X2s, Olympus V3 and P7 allowed both automatic and manual modes on the light stack. For the Storz Flex-X2s in manual mode, testing was done with the light stack brightness at 0/50/100% settings. With the automatic mode, the Storz Flex-X2s video image menu interestingly allowed only for setting a specific automatic mode maximum (max) brightness value (50–100%), which when adjusted, was also concurrently reflected on the light stack brightness bar setting. After the ureteroscope was inserted into the model, the light brightness bar on the light stack then adjusts to an “automatic” brightness setting within that maximum value. Images were captured for max 50%/max 75%/max 100% automatic settings on the video image menu. In contrast, for the Olympus V3 and P7 brightness automatic modes, the light stack interestingly allowed setting a fixed brightness value (0–100%), with no subsequent auto adjustment reflected by the light stack brightness bar. The two Olympus scopes were tested at 0/50/100% light stack brightness settings for both auto and manual modes.

For the Pusen 7.5F and 9.2F scopes, testing was done at processing unit light brightness settings of 0/4/7 out of 7, as provided by the manufacturer. For the WiScope, this was tested at 0/50/100% brightness setting, as provided by the manufacturer.

For each different combination of light brightness setting, video capture mode, and distance, a set of five repeated images were captured on the ureteroscope processing unit. Each individual combination is referred to as “situation” from here on. The ureteroscope was withdrawn out of the model and reinserted before each repeated image capture.

Statistical analysis

All captured images were transferred and analyzed for their histogram grayscale pixel count distribution of brightness using the software *ImageJ* (Version 1.53t RRID: SCR_003070) [29]. This resulted in 256 pixel count measurements for each of the 256 individual histogram grayscale values (histogram scale 0–255). 0 on the scale represents the darkest on the grayscale, and 255 the brightest.

For each situation, the mean pixel count for each 256 individual histogram grayscale values (0, 1, 2, 255) was calculated from the five values of the five repeated captured images, giving a mean histogram grayscale distribution for each individual situation.

The overall histogram median value for each situation was defined as the median grayscale value based on the above calculated pixel counts (Fig. 2). All analyses were performed with GraphPad Prism 10.0.1 (GraphPad Software, La Jolla CA, USA).

Blackout was defined as median histogram range from 0 to 35, and whiteout from range of 220 to 255 (image on monitor too dark or too bright for the human eye, respectively). Cut-offs values were determined by consensus between two authors, J-LK and EXK after review of all captured images. An example of images on a grayscale can be found in Fig. 2, and an example of endoscopic image median histogram values at various brightness settings and modes can be found in Fig. 3.

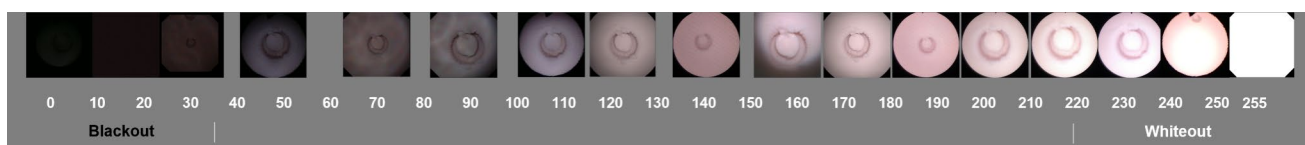


Fig. 2 Distribution of endoscopic images over grayscale histograms. Images shown are representative images within the set of five endoscopic images captured for various scopes, at various brightness and distances from the target, arranged in order of the overall histogram

median grayscale value for that set of five images. This is meant to illustrate the appearance of blackout and whiteout on the histogram scale (range 0–255). Blackout was defined as 0–35 and whiteout 220–255 on the histogram scale

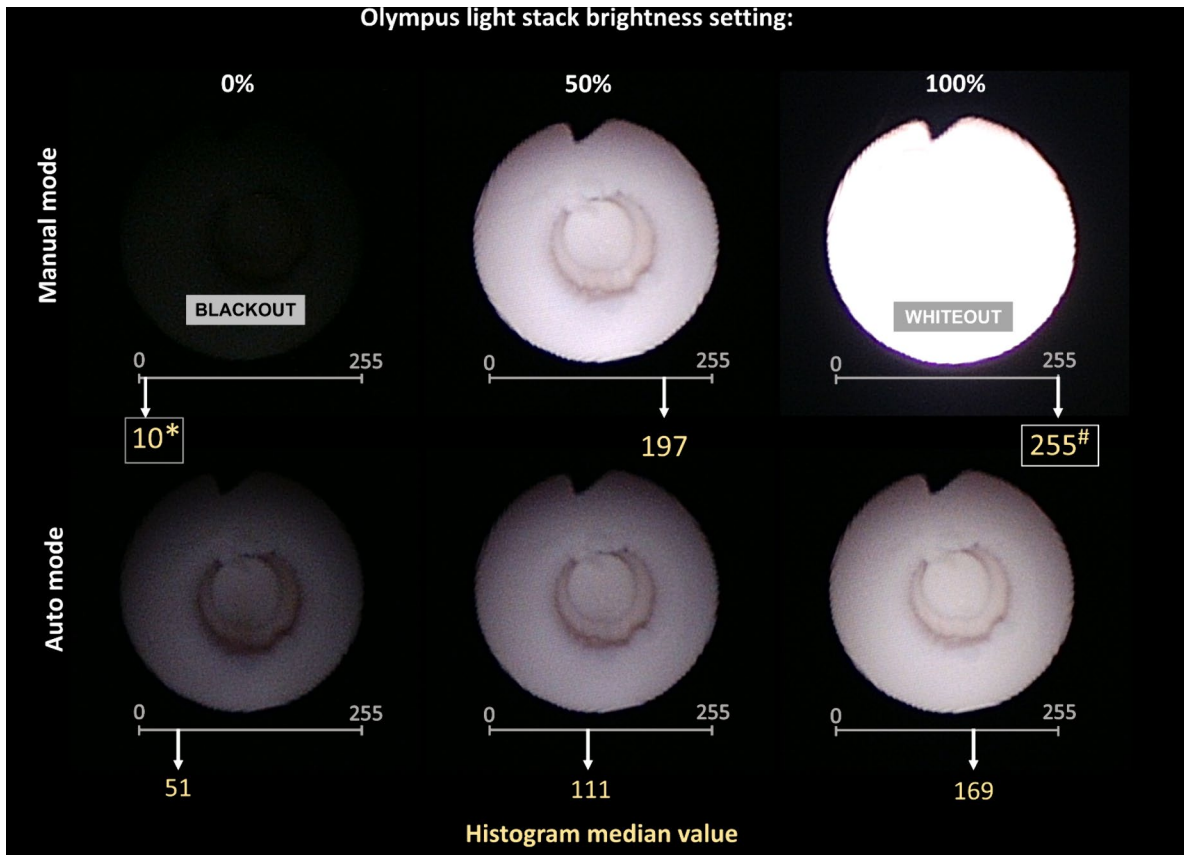


Fig. 3 Blackout and whiteout in flexible ureteroscopy. Example images are from the Olympus P7 at 5 mm distance. These are representative images within the set of five endoscopic images captured for each situation. The median histogram value is based on a set of five

endoscopic images. The histogram median values are for each brightness setting and mode on the histogram scale of 0 to 255. *Histogram value meeting blackout criteria. # Histogram value meeting whiteout criteria

Results

Blackout

Blackout occurred with the P7, Pusen 7.5F and 9.2F, WiScope and at all evaluated distance settings (5 mm, 10 mm, 20 mm) - all with lowest brightness settings. For the V3, this was only found at the furthest 20 mm distance (Table 1).

For scopes with adjustable brightness modes, blackout occurred only in the manual mode and not in automatic mode (V3 and P7). The Flex-X2s did not have blackout even in manual mode.

Interestingly, the Flex-Xc and Flex-X2s had no blackout even at the furthest 20 mm distance and lowest brightness mode setting available.

Whiteout

Whiteout occurred with the Flex-X2s, V3 and P7 at 5 mm and 10 mm distance, as well as at 20 mm distance with

the V3 and P7 - mostly with 50–100% brightness settings (Table 1).

For scopes with adjustable brightness mode, whiteout happened for the Flex-X2s in both manual and automatic mode, mostly when the video image capture setting was five out of five. For the V3 and P7, whiteout happened exclusively in the manual mode at all tested distances (5, 10, 20 mm).

No blackout or whiteout

The Flex-Xc was the only scope that had neither blackout nor whiteout over all brightness settings and brightness modes.

Discussion

Our findings reveal that within the range of brightness settings available, some scopes can blackout or whiteout, while others do not. This phenomenon mainly occurred when light

Table 1 Blackout and whiteout characteristics of flexible ureteroscopes at various distances from target

Scope		Blackout/whiteout /None	Distance (5,10, 20mm)	Brightness mode	Light Stack / Processing unit Brightness Setting	Video image brightness settings	Median Histogram Value (IQR)		
Reusable	Storz Flex-Xc	None at all tested settings	5, 10, 20mm	-	-	1, 3, 5 of 5	-		
	Storz Flex-X2s	Whiteout	5mm	Manual	0%	5 of 5	232 (212 – 246)		
					50%	5 of 5	235 (216 – 248)		
					100%	5 of 5	238 (219 – 250)		
				Auto	Max 50%	5 of 5	233 (213 – 246)		
					Max 75%	5 of 5	230 (212 – 242)		
					Max 100%	5 of 5	232 (213 – 244)		
			10mm	Manual	100%	5 of 5	226 (213 – 235)		
					50%	5 of 5	221 (209 – 232)		
				Auto	Max 50%	5 of 5	220 (207 – 232)		
					Max 75%	5 of 5	221 (208 – 234)		
	Olympus V3	Blackout	20mm	Manual	0%	-	33 (29 – 36)		
			Whiteout	5mm	Manual	50%	-	255 (255 – 255)	
		100%				-	255 (255 – 255)		
		10mm		Manual	50%	-	255 (252 – 255)		
					100%	-	255 (255 – 255)		
		20mm		Manual	100%	-	255 (255 – 255)		
		Olympus P7		Blackout	5mm	Manual	0%	-	10 (8 – 14)
					10mm	Manual	0%	-	6 (6 – 7)
					20mm	Manual	0%	-	6 (6 – 6)
Whiteout			5mm		Manual	100%	-	255 (255 – 255)	
	10mm		Manual	100%	-	255 (255 – 255)			
	20mm		Manual	100%	-	246 (228 – 255)			
	Single-use		Pusen 7.5F	Blackout	5mm	Not adjustable	0 of 7	-	17 (16 – 17)
10mm					Not adjustable	0 of 7	-	17 (16 – 17)	
20mm		Not adjustable			0 of 7	-	14 (14 – 15)		
Pusen 9.2F		Blackout	5mm	Not adjustable	0 of 7	-	34 (33 – 35)		
			10mm	Not adjustable	0 of 7	-	34 (33 – 35)		
			20mm	Not adjustable	0 of 7	-	34 (33 – 35)		
OTU WiScope		Blackout	5mm	Not adjustable	0%	-	20 (15 – 23)		
			10mm	Not adjustable	0%	-	18 (14 – 22)		
			20mm	Not adjustable	0%	-	18 (14 – 22)		

Settings explored: (1) lowest, mid and highest brightness settings of the available light stack or processing unit for all ureteroscopes, (2) Auto/manual brightness mode for the Storz Flex-X2s, Olympus V3, Olympus P7, (3) Video image brightness setting for Storz Flex-Xc and Flex-X2s

source brightness or image brightness settings were at their extreme low or high ends. These undesirable settings may occur during clinical routine, whether consciously or inadvertently set by operating room staff or members.

Blackout of scopes is particularly undesirable, as the surgeon would like to maintain a decent level of vision even

at the lowest brightness setting of the scope. The only rare situation where blackout is desired is in a setting of an endoscope meeting procedure, typically referred to as a “cut to the light” procedure [30]. In this special situation, one would need a scope that is capable of switching off its light source to artificially create a blackout situation. From there, the

light source of the second endoscope can be searched for. All single-use scopes in our study had blackout at all distances with the lowest brightness settings. It would be of interest to test more single-use scopes to see if this characteristic is found in other single-use scopes.

Generally, whiteout happens during higher light stack brightness settings for the Flex-X2s, V3 and P7. This happens as well during highest settings on the Flex-X2s video image brightness settings. Based on these findings, it may be prudent for most scopes to start at 50% brightness light stack brightness setting and 3 out of 5 video image brightness setting (when available), before adjusting for the particular in vivo situation. It is also important to be aware that some scopes like the Flex-X2s, depending on the camera head used, may have both light brightness settings on the light stack and on the camera head video image menu.

The automatic mode of the modern ureteroscopic light stack helps to adjust the brightness of the image. Based on our findings of whiteout happening mostly in the manual mode for scopes with adjustable modes, the authors would recommend keeping light mode on automatic for most clinical applications. It is also interesting that we were able to adjust brightness settings in the auto mode for the Flex-X2s, V3 and P7, influencing brightness in an “automatic” mode.

The present study is not devoid of limitations. First, the authors arbitrarily used median histogram value cut-offs to determine blackout and whiteout criteria. The strength of this methodology is that it is based on histogram values which provide an objective way to qualitatively determine blackout and whiteout in images from scopes. Second, this is an in-vitro attempt to reflect in vivo use of ureteroscopes. Still images were used instead of videos, with the methodology dictating an objective assessment via still image histogram values. If videos were used, this would be technically challenging to compute for each frame of the video. To ensure reliability of the results, a set of five images was captured for each individual combination of brightness setting, brightness mode and distance. Environmental factors may play a role in affecting the clinical translation of our study findings, and thus interpretation of the data must, therefore, be taken with care. The effect of blood clots, stone dust and urine on blackout and whiteout characteristics will need to be further evaluated. Third, other situations not evaluated in this study include that of situations in the ureter with the enclosed working cavity being a hollow tube rather than spherical in nature. This warrants further evaluation.

Conclusions

Blackout or whiteout of images is an undesirable property that was found for several scopes, possibly impacting diagnostic and therapeutic purposes during ureteroscopy. These

observations form a guide for urologists which may impact their choice of instruments and corresponding settings.

Author contributions Kwok: Protocol/project development, data collection or management, data analysis, manuscript writing/editing. De Coninck: Protocol/project development, data analysis, manuscript writing/editing. Panthier: Protocol/project development, data analysis, manuscript writing/editing. Kamkoum: Protocol/project development, data analysis, manuscript writing/editing. Pauchard: Protocol/project development, data analysis, manuscript writing/editing. Shrestha: Protocol/project development, data analysis, manuscript writing/editing. Gauhar: Protocol/project development, data analysis, manuscript writing/editing. Brachlow: Protocol/project development, data analysis, manuscript writing/editing. Schmid: Protocol/project development, data analysis, manuscript writing/editing. Poyet: Data analysis, manuscript writing/editing. Eberli: Data analysis, manuscript writing/editing. Traxer: Research concept, data analysis, manuscript writing/editing. Keller: Research concept, protocol/project development, data collection or management, data analysis, manuscript writing/editing.

Funding Open access funding provided by University of Zurich.

Data availability On request to corresponding author for raw data on the experimental setup.

Declarations

Conflict of interests Vincent De Coninck is a consultant for Axess Vision Technology, BD Bard and Coloplast, and has no specific conflicts relevant to this study. Frédéric Panthier is a consultant for Dornier Medtech and has no specific conflicts relevant to this study. Olivier Traxer is a consultant for Coloplast, Karl Storz, Rocamed, Quanta Systems and IPG Medical and has no specific conflicts relevant to this study. Etienne Xavier Keller is a speaker and/or consultant for Coloplast, Olympus and Boston Scientific, and has no specific conflicts of interest relevant to this work. All other authors have no relevant financial or non-financial interests to disclose.

Ethics approval Not applicable for this in-vitro study.

Open Access This article is licensed under a Creative Commons Attribution 4.0 International License, which permits use, sharing, adaptation, distribution and reproduction in any medium or format, as long as you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons licence, and indicate if changes were made. The images or other third party material in this article are included in the article's Creative Commons licence, unless indicated otherwise in a credit line to the material. If material is not included in the article's Creative Commons licence and your intended use is not permitted by statutory regulation or exceeds the permitted use, you will need to obtain permission directly from the copyright holder. To view a copy of this licence, visit <http://creativecommons.org/licenses/by/4.0/>.

References

- Geraghty RM, Jones P, Somani BK (2017) Worldwide trends of urinary stone disease treatment over the last two decades: a systematic review. *J Endourol* 31(6):547–556. <https://doi.org/10.1089/end.2016.0895>
- Monga M, Murphy M, Paranjpe R, Cutone B, Eisner B (2023) Prevalence of stone disease and procedure trends in the United

- States. *Urology* 176:63–68. <https://doi.org/10.1016/j.urology.2023.03.040>
3. Heers H, Stay D, Wiesmann T, Hofmann R (2022) Urolithiasis in Germany: trends from the national DRG database. *Urol Int* 106(6):589–595. <https://doi.org/10.1159/000520372>
 4. Hendriks N, Henderickx MMEL, Schout BMA, Baard J, van Etten-Jamaludin FS, Beerlage HP et al (2021) How to evaluate a flexible ureterorenoscope? Systematic mapping of existing evaluation methods. *BJU Int* 128(4):408–423. <https://doi.org/10.1111/bju.15544>
 5. Kwok J-L, De Coninck V, Corrales M, Sierra A, Panthier F, Ventimiglia E, et al. Illumination matters Part I: Comparative analysis of light sources and illumination in flexible ureteroscopy – fundamental findings from a PEARLS analysis. Undergoing peer review. Not yet published.
 6. Al Qahtani SH, Abdelhamied MH, AlMuhrij AH, Al Rawashada MY, Al Askar AM, Abdelhamid AM et al (2020) Prospective comparative study between the effect of CIDEX® OPA and STERRAD NX on the durability of digital flexible ureteroscope. *World J Urol* 38(3):775–781. <https://doi.org/10.1007/s00345-019-02800-z>
 7. Al-Qahtani SM, Geavlette BP, De Medina SGD, Traxer OP (2011) The new Olympus digital flexible ureteroscope (URF-V): Initial experience. *Urol Ann* 3(3):133–137. <https://doi.org/10.4103/0974-7796.84963>
 8. Cho SY, Lee JY, Shin DG, Seo IY, Yoo S, Park HK (2018) Evaluation of performance parameters of the disposable flexible ureterorenoscope (LITHOVUE) in patients with renal stones: a prospective, observational, single-arm multicenter study. *Sci Rep*. <https://doi.org/10.1038/s41598-018-28247-7>
 9. Doizi S, Kamphuis G, Giusti G, Andreassen KH, Knoll T, Osther PJ et al (2017) First clinical evaluation of a new single-use flexible ureteroscope (LithoVue™): a European prospective multicentric feasibility study. *World J Urol* 35(5):809–818. <https://doi.org/10.1007/s00345-016-1936-x>
 10. Dragos LB, Somani BK, Sener ET, Buttice S, Proietti S, Ploumidis A et al (2017) Which flexible ureteroscopes (digital vs fiber-optic) can easily reach the difficult lower pole calices and have better end-tip deflection In vitro study on K-box A PETRA evaluation. *J Endourol* 31(7):630–7. <https://doi.org/10.1089/end.2017.0109>
 11. Emiliani E, Mercadé A, Millan F, Sánchez-Martín F, Konstantinidis CA, Angerri O (2018) First clinical evaluation of the new single-use flexible and semirigid pusen ureteroscopes. *Central Eur J Urol* 71(2):208–213. <https://doi.org/10.5173/cej.2018.1620>
 12. Emiliani E, Talso M, Baghdadi M, Barreiro A, Orosa A, Servià P et al (2017) Evaluation of the spies™ modalities image quality. *Int Braz J Urol* 43(3):476–480. <https://doi.org/10.1590/S1677-5538.IBJU.2016.0324>
 13. Johnston TJ, Baard J, de la Rosette J, Doizi S, Giusti G, Knoll T et al (2018) A clinical evaluation of the new digital single-use flexible ureteroscope (UscopePU3022): an international prospective multicentered study. *Central Eur J Urol* 71(4):453–461. <https://doi.org/10.5173/cej.2018.1787>
 14. Kam J, Yuminaga Y, Beattie K, Ling KY, Arianayagam M, Canagasingham B et al (2019) Single use versus reusable digital flexible ureteroscopes: a prospective comparative study. *Int J Urol* 26(10):999–1005. <https://doi.org/10.1111/iju.14091>
 15. Kim DH, Shin JH, Choi SY, You D, Kim CS, Park HK (2018) In vitro, in vivo, and clinical tests of a novel flexible ureteroscope for the diagnosis and treatment of kidney and ureteral diseases. *Investig Clin Urol* 59(5):328–334. <https://doi.org/10.4111/icu.2018.59.5.328>
 16. Legemate JD, Kamphuis GM, Freund JE, Baard J, Zanetti SP, Catellani M et al (2019) Durability of flexible ureteroscopes: a prospective evaluation of longevity, the factors that affect it, and damage mechanisms. *Eur Urol Focus* 5(6):1105–1111. <https://doi.org/10.1016/j.euf.2018.03.001>
 17. Lusch A, Abdelshehid C, Hidas G, Osann KE, Okhunov Z, McDougall E et al (2013) In vitro and in vivo comparison of optics and performance of a distal sensor ureteroscope versus a standard fiberoptic ureteroscope. *J Endourol* 27(7):896–902. <https://doi.org/10.1089/end.2013.0003>
 18. Multescu R, Geavlete B, Geavlete P (2013) A new era: performance and limitations of the latest models of flexible ureteroscopes. *Urology* 82(6):1236–1239. <https://doi.org/10.1016/j.urology.2013.07.022>
 19. Proietti S, Dragos L, Molina W, Doizi S, Giusti G, Traxer O (2016) Comparison of new single-use digital flexible ureteroscope versus nondisposable fiber optic and digital ureteroscope in a cadaveric model. *J Endourol* 30(6):655–659. <https://doi.org/10.1089/end.2016.0051>
 20. Schlager D, Hein S, Obaid MA, Wilhelm K, Miernik A, Schoenthaler M (2017) Performance of single-use flexorvue vs reusable boavision ureteroscope for visualization of calices and stone extraction in an artificial kidney model. *J Endourol* 31(11):1139–1144. <https://doi.org/10.1089/end.2017.0454>
 21. Schlager D, Obaid MA, Hein S, Wilhelm K, Schönthaler M, Gratzke C et al (2020) Current disposable ureteroscopes: performance and limitations in a standardized kidney model. *J Endourol* 34(10):1015–1020. <https://doi.org/10.1089/end.2020.0185>
 22. Talso M, Proietti S, Emiliani E, Gallioli A, Dragos L, Orosa A et al (2018) Comparison of flexible ureterorenoscope quality of vision: an in vitro study. *J Endourol* 32(6):523–528. <https://doi.org/10.1089/end.2017.0838>
 23. Winship B, Wollin D, Carlos E, Li J, Preminger GM, Lipkin ME (2019) Avoiding a lemon: performance consistency of single-use ureteroscopes. *J Endourol* 33(2):127–131. <https://doi.org/10.1089/end.2018.0805>
 24. Marroig B, Favorito LA, Fortes MA, Sampaio FJB (2015) Lower pole anatomy and mid-renal-zone classification applied to flexible ureteroscopy: experimental study using human three-dimensional endocasts. *Surg Radiol Anat* 37(10):1243–1249. <https://doi.org/10.1007/s00276-015-1503-y>
 25. Aldoukhi AH, Roberts WW, Hall TL, Teichman JMH, Ghani KR (2018) Understanding the popcorn effect during holmium laser lithotripsy for dusting. *urology* 122:52–57. <https://doi.org/10.1016/j.urology.2018.08.031>
 26. Aldoukhi AH, Hall TL, Ghani KR, Roberts WW (2021) Strike rate: analysis of laser fiber to stone distance during different modes of laser lithotripsy. *J Endourol* 35(3):355–360. <https://doi.org/10.1089/end.2020.0298>
 27. Kwok J-L, Panthier F, De Coninck V, Ventimiglia E, Barghouthy Y, Danilovic A, et al (2024) Illumination matters Part II: Advanced comparative analysis of flexible ureteroscopes in a kidney model by PEARLS. *World J Urol*. <https://doi.org/10.1007/s00345-024-04987-2>
 28. Kwok J-L, Ventimiglia E, De Coninck V, Sierra A, Panthier F, Corrales M, et al (2024) Illumination matters Part III: Impact of light obstruction on illuminance from flexible ureteroscopes - a comparative PEARLS analysis. *World J Urol* 42(1):188. <https://doi.org/10.1007/s00345-024-04910-9>
 29. Rasband WS. U. S. National Institutes of Health Bethesda. Maryland, USA1997–2018. p. ImageJ.
 30. Shiozaki K, Izaki H, Fukuta K, Fukawa T, Takahashi M, Kanayama H (2022) Laser endoureterotomy with cut-to-the-light technique for complete ureteral obstruction—a case report. *Urol Case Rep*. <https://doi.org/10.1016/j.eucr.2022.102036>