



University of Pennsylvania
ScholarlyCommons

Tool Data

Browse by Type

5-7-2016

Reactive Ion Etch (RIE) of Silicon Dioxide (SiO₂) with Trifluoromethane and Oxygen (CHF₃/O₂)

Meredith Metzler

Singh Center for Nanotechnology, metzlerm@seas.upenn.edu

Follow this and additional works at: http://repository.upenn.edu/scn_tooldata

 Part of the [Nanoscience and Nanotechnology Commons](#)

Metzler, Meredith, "Reactive Ion Etch (RIE) of Silicon Dioxide (SiO₂) with Trifluoromethane and Oxygen (CHF₃/O₂)", *Tool Data*. Paper 38.

http://repository.upenn.edu/scn_tooldata/38

This paper is posted at ScholarlyCommons. http://repository.upenn.edu/scn_tooldata/38

For more information, please contact repository@pobox.upenn.edu.


Reactive Ion Etch (RIE) of Silicon Dioxide (SiO₂) with Trifluoromethane and Oxygen (CHF₃/O₂)

Summary/Description

This report discusses the CHF₃/O₂ etch process of SiO₂ using the Oxford 80 Plus RIE.

Disciplines

Nanoscience and Nanotechnology

	Reactive Ion Etch (RIE)	Document No:
	Silicon Dioxide (SiO ₂) with CHF ₃ / O ₂	Revision:
	Oxford 80 Plus	Author: Meredith Metzler

1. Introduction

The purpose of this document is to examine the etch properties of the Oxford 80 Plus RIE system.

2. Baseline Recipe

Units:

Gas flow rate: standard cubic centimeters per minute (sccm)

Pressure: millitorr (mT)

Temperature: degrees Celsius (C)

High frequency (RF) power: Watts (W)

Step 1: Pump to 5e-04 Torr, “Pump to Pressure” checked

Step 2: Etch Step

Trifluoromethane (CHF₃) flow rate: 100 sccm

Oxygen (O₂) flow rate: 4 sccm

Pressure: 50 mT

RF Power: 150 W

Capacitor starting points: Capacitor #1: 60 %, Capacitor #2: 80 %

Time set point is hh:mm:ss (hours:minutes:seconds)*

Temperature: 15 C

Step 3: Pump to 5e-04 Torr, “Pump to Pressure” checked

*notes for Step 2: The time set point for the etch step should be kept below 10 minutes due to thermal issues and to avoid resist burning. If a longer time is needed for a thicker film then the system should be vented prior to running the process again.

3. Etch Characteristics

Film thickness is measured using a Filmetrics F50 optical interferometer which is equipped with a motorized stage allowing for the collection of full wafer maps. See the following link for more information about this instrument:

<http://www.filmetrics.com/thicknessmeasurement/f50>

The film being etched is PECVD SiO₂ deposited on 100 mm, <100> orientation, wafers that are 525 ± 25 micron thick.

Figure 1 below shows a screen capture image of a “Difference Map” from the Filmetrics software with 115 data points and a 5 mm edge exclusion. The standard SiO₂ material file supplied in the software is used for these measurements. This is data from a 3 minute etch displayed as a “difference map” in the software that is already averaged to display the etch rate in nm/min.

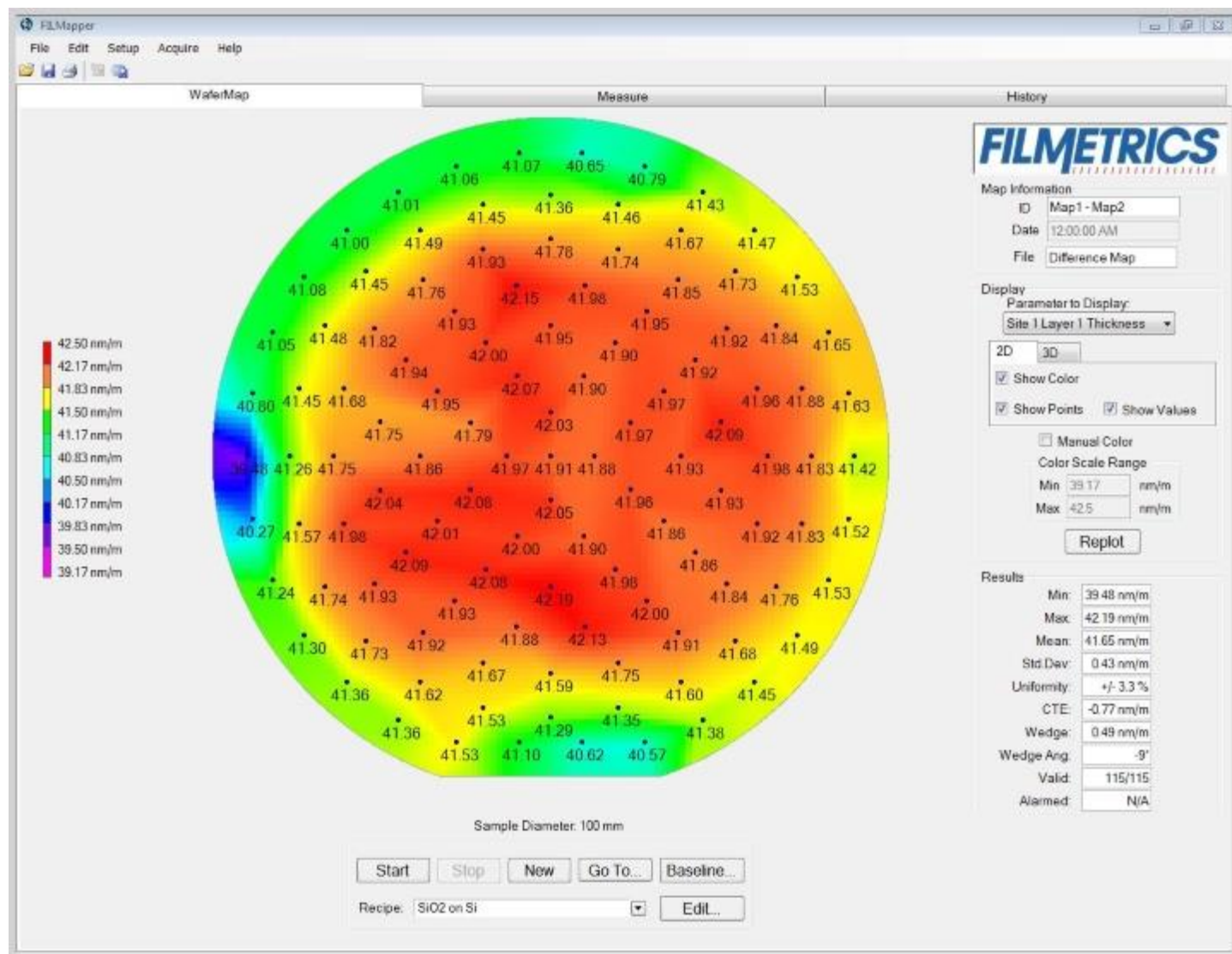


Figure 1. Wafer map showing the results for a 3 minute SiO₂ etch using CHF₃ and O₂ showing 42 nm/min etch rate with a standard deviation of 0.43 nm and a uniformity across the wafer of ± 3.3%.