

ELECTRONIC SUPPORTING INFORMATION

Theranostic Body Fluid Cleansing: Rationally designed Magnetic Particles enable Capturing and Detection of Bacterial Pathogens

*M. Lattuada, Q. Ren, F. Zuber, M. Galli, N. Bohmer, M.T. Matter, A. Wichser, S. Bertazzo, G.B. Pier and I.K. Herrmann**

**Department Materials Meet Life,*

Swiss Federal Laboratories for Materials Science and Technology (Empa),

Lerchenfeldstrasse 5,

9014 St. Gallen, Switzerland.

email: inge.herrmann@empa.ch; phone: +41 58 765 71 53

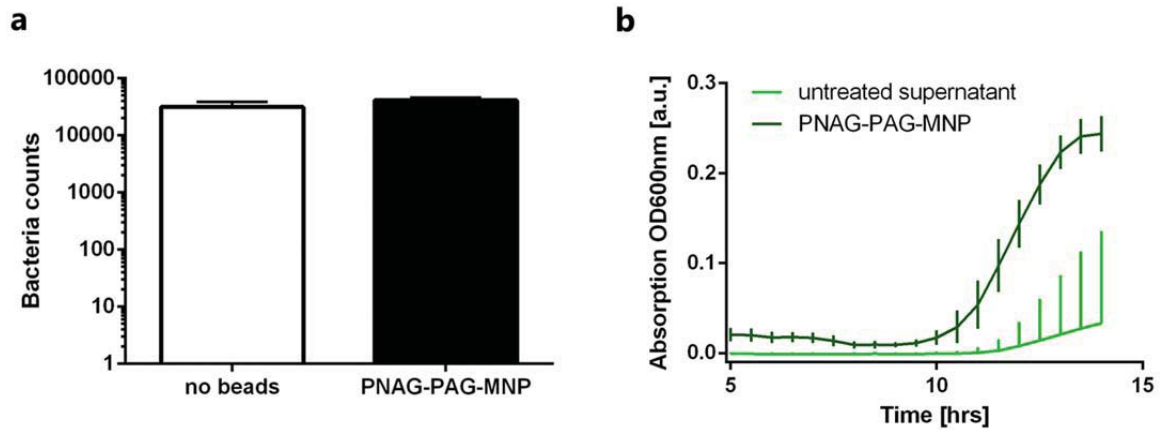


Figure S1: Bacteria samples incubated with PNAG-PAG-MNP show normal growth rate compared to controls (no beads)(**a**). Measurements of re-growth of *S. aureus* in the untreated supernatant (no beads) and on collected beads show enrichment of bacteria on the beads (**b**).

Bacteria viability quantification with proliferation test (Re-growth). 100 μ l of the resuspended beads was mixed with 900 μ l fresh 30 % TSB (9 g/l) containing 0.25% glucose. 200 μ l of the mixture was added to a microplate well (3 repeats) and the plate was closed with a sticking transparent foil (Diversified Biotech - Breathe-Easy, T093.1, Carl Roth GmbH+Co. KG) and incubated in a Synergy HT Multi-Detection Microplate Reader (BioTek®) at 37 °C and medium shaking (18Hz - preset value of the instrument). Optical density at 600 nm was measured every 30 min for 14 hours. For each sample, two independent experiments were performed.