

LITERATURE REVIEW

Zirconia: Substitute for Metal Ceramics

¹C Gyan Kumar, ²DP Shruthi, ³K Sounder Raj, ⁴D Kalpana, ⁵G Harish**ABSTRACT**

In the search for the ultimate esthetic restorative material, many all-ceramic systems have been proposed. Dental research is nowadays directed toward metal-free prosthetic restorations in order to improve an esthetical outcome of crown and bridge restorations. Zirconia is a polycrystalline ceramic without a glassy phase and exists in several forms. Zirconia cores for fixed partial dentures (FPD) on anterior and posterior teeth and on implants are now available. This article is a review of zirconia crowns as replacement for conventional metal-ceramic crowns.

Keywords: Zirconia, Metal ceramic, Cosmetic dentistry, Esthetics, Mechanical properties.

How to cite this article: Kumar CG, Shruthi DP, Raj KS, Kalpana D, Harish G. Zirconia: Substitute for Metal Ceramics. *J Orofac Res* 2014;4(4):209-212.

Source of support: Nil

Conflict of interest: None declared

INTRODUCTION

Natural look of soft tissue in contact with fixed partial dentures (FPD) is influenced by two factors: mucosal thickness and typology of restorative material. Metal free restorations allow preserving soft tissue color more similar to the natural one than porcelain fused to metal restorations.

Meta-analyses have demonstrated good long term clinical results for the conventional fixed dental prosthesis (FDP) with a metal framework.¹⁻³ As dentistry has evolved, the demand for metal free materials with increased translucency (Figs 1 and 2) that mimic the natural dentition has arisen.⁴

This has led to the development of several different ceramics that are esthetically pleasing and biocompatible.⁵⁻⁷ What is still being established is which materials are acceptable not only for their favorable optical properties, but also for their favorable mechanical properties, adequate clinical function and longevity.

Zirconia (ZrO₂) is a ceramic material with adequate mechanical properties and is a crystalline dioxide of zirconium. Its mechanical properties are very similar to those of metals and its color is similar to tooth color.⁸ In 1975, Garvie proposed a model to rationalize the good mechanical properties of Zirconia, by virtue of which it has been called 'ceramic steel'.⁹

Zirconia crystals can be organized in three different patterns: monoclinic (M), cubic (C) and tetragonal (T). By mixing zirconia (ZrO₂) with other metallic oxides, such as MgO, CaO, or Y₂O₃, great molecular stability can be obtained.¹⁰ Yttrium-stabilized zirconia, also known as tetragonal zirconia polycrystal (TZP), is presently the most studied combination.¹¹ Zirconia stabilized with Y₂O₃ has the best properties for these applications.

When a stress occurs on a zirconia (ZrO₂) surface, a crystalline modification opposes the propagation of cracks. Its resistance to traction can be as high as 900 to 1200 MPa and its compression resistance is about 2000 MPa. Cyclical stresses are also tolerated well by this material. Applying an intermittent force of 28 kN to zirconia substrates, Cales found that some 50 billion cycles were necessary to break the samples, but with a force in excess of 90 kN structural failure of the samples occurred after just 15 cycles.¹² Surface treatments can modify the physical properties of zirconia. Exposure to wetness for an extended period of time can have a detrimental effect on its properties.¹³ This phenomenon is known as zirconia ageing. Moreover, also surface grinding can reduce toughness.¹⁴ Kosmac confirmed this observation and reported a lower mean strength and reliability of zirconium oxide after grinding.¹⁵

Basic Properties

The first proposal of the use of zirconium oxide for medical purposes was made in 1969 and concerned orthopedic application. Orthopedic research focused on the mechanical behavior of zirconia, on its wear, and on its integration with bone and muscle. Since 1990, *in vitro* studies

¹Postgraduate Student, ²Private Practitioner, ^{3,4}Professor and Head, ⁵Reader

^{1,3}Department of Prosthodontics and Implantology, Government Dental College and Research Institute, Victoria Hospital Campus, Bengaluru, Karnataka, India

²No. 790, 20th Cross, 18th Main, Banashankari Stage II Bengaluru, Karnataka, India

^{4,5}Department of Prosthodontics, Dayanand Sagar College of Dental Sciences, Bengaluru, Karnataka, India

Corresponding Author: C Gyan Kumar, Postgraduate Student, Department of Prosthodontics and Implantology Government Dental College and Research Institute, Victoria Hospital Campus, Fort, Bengaluru, Karnataka, India, e-mail: gyan_kumarc@yahoo.co.in

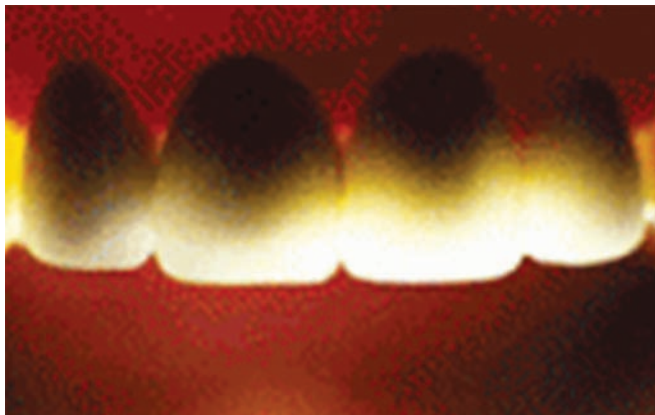


Fig. 1: Metal ceramic restorations

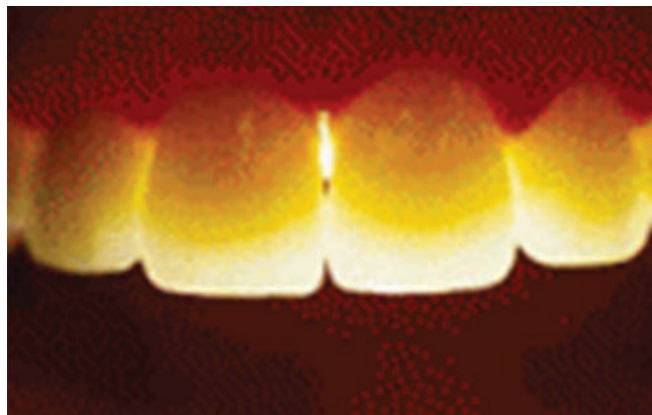


Fig. 2: Zirconium based restorations

have also been performed in order to obtain information about cellular behavior toward zirconia.¹⁶ *In vitro* evaluation confirmed that zirconia (ZrO_2) is not cytotoxic.¹⁷⁻¹⁹ Zirconium oxide creates less flogistic reaction in tissue than other restorative materials, such as titanium.²⁰ This result was also confirmed by a study about peri-implant soft tissue around zirconia healing caps in comparison with that around titanium ones.²¹ Inflammatory infiltrate, microvessel density, and vascular endothelial growth factor expression was found to be higher around the titanium caps than around the zirconia (ZrO_2) ones. Also, the level of bacterial products, measured with nitric oxide synthase, was higher on titanium than on zirconium oxide. Zirconia can up- or down-regulate expressions of some genes, so that zirconia can be regarded as a self-regulatory material that can modify turnover of the extracellular matrix.²²

Clinical Applications

Raigrodski analyzed different all-ceramic systems and concluded that reinforced ceramics can only be used to replace anterior teeth with single crown restorations or maximum with three units FPD. On the other hand, zirconia (ZrO_2) restorations have a wider application field.

Other ceramic technologies only allow the construction of structures that are resistant to chewing stresses on anterior teeth. On the contrary, zirconia-ceramic FPD can also be used on molars.²³

Tinschert compared lifetime of different metal-free core for FPD and reported that zirconia-ceramic with alumina oxide had the highest initial and most favorable long-term strength.²⁴ Connecting surface area of the FPD must be at least 6.25 mm^2 .²⁵ For this reason, ceramic FPD should only be used when the distance between the interproximal papilla and the marginal ridge is close to 4 mm. In a comparison between 3, 4 and 5 unit zirconia FPD and minimal connecting surface resulted, respectively, 2.7 mm^2 , 4.0 mm^2 and 4.9 mm^2 .²⁶ Height of abutment is fundamental to obtain zirconia (ZrO_2) frameworks with correct shape and dimension in order to ensure mechanical resistance of restoration. This aspect must be carefully considered when realizing a metal-free FPD.

Zirconia restorations have found their indications for FPD supported by teeth or implants. Single tooth restorations and FPD with a single pontic element are possible on both anterior and posterior elements because of the mechanical reliability of this material.²⁷⁻²⁹ It is possible to use juxtagingival marginal preparations and various



Fig. 3: Prepared discolored non-vital teeth



Fig. 4: Cemented zirconia-ceramic crowns

finishing lines to obtain a good esthetic (Figs 3 and 4).³⁰ Fixed partial denture extension is nowadays a limitation in using zirconia-ceramic restorations. Although some manufacturer allows obtaining also full arch restorations, five units FPD are reported to be as maximal possible.³¹

Zirconia opacity is very useful in adverse clinical situations, for example, for masking of dischromic abutment teeth.

Radiopacity can aid evaluation during radiographic controls. Zirconia frameworks are realized by using computer-aided design/manufacturing (CAD/CAM) technology. As the development of zirconia crown using CAD/CAM technology, the usage of full zirconia crown is gradually increased.

When compared to the layering technique which is fabricated with zirconia coping and veneering porcelain, full zirconia crown shows higher strength and easier laboratory procedure. To prevent mechanical fracture of the full zirconia crown, the thickness of zirconia crown and proper sintering process should be considered.³² Nearly all (87.5%) conventionally veneered crowns failed already during chewing simulation, whereas crowns with CAD/CAM manufactured veneers were nonsensitive to artificial ageing. Crowns veneered with lithium disilicate ceramic displayed ultimate loads to failure of about 1600 N.³³

Although many types of zirconia-containing ceramic systems are currently available only three are used to date in dentistry. These are yttrium cation-doped tetragonal zirconia polycrystals (3Y-TZP), magnesium cation-doped partially stabilized zirconia (Mg-PSZ) and zirconia-toughened alumina (ZTA).³⁴

CONCLUSION

Although clinical long-term evaluations are a critical requirement to conclude that zirconia has good reliability for dental use, biological, mechanical, and clinical studies published to date seem to indicate that zirconia (ZrO₂) restorations are both well tolerated and sufficiently resistant. Ceramic bonding, luting procedures, ageing and wear of zirconia abutment should be evaluated in order to guide the adequate use of zirconia as prosthetic restorative material. Patient selection, coupled with adequate clinical and technical protocols, are imperative in order to obtain good performance of these restorations. The CAD/CAM production of veneers for restorations with zirconia framework is a promising way to reduce failures originating from material fatigue.

REFERENCES

- Creugers NH, Kayser AF, van't Hof MA. A meta-analysis of durability data on conventional fixed bridges. *Community Dent Oral Epidemiol* 1994 Dec;22(6):448-452.
- Scurria MS, Bader JD, Shugars DA. Meta analysis of fixed partial denture survival: prostheses and abutments. *J Prosthet Dent* 1998 Apr;79(4):459-464.
- Tan K, Pjetursson BE, Lang NP, Chan ES. A systematic review of the survival and complication rates of fixed partial dentures (FPDs) after an observation period of at least 5 years. *Clin Oral Implants Res* 2004 Dec;15(6):654-666.
- Raptis NV, Michalakis KX, Hirayama H. Optical behavior of current ceramic systems. *Int J Perio Restor Dent* 2006 Feb;26(1):31-41.
- McLean JW, Hughes TH. The reinforcement of dental porcelain with ceramic oxides. *Br Dent J* 1965 Sep 21;119(6):251-267.
- Raigrodski AJ, Chiche GJ. All-ceramic fixed partial dentures, part I: in vitro studies. *J Esthet Restor Dent* 2002;14(3):188-191.
- Raigrodski AJ, Chiche GJ, Swift EJ Jr. All-ceramic fixed partial dentures, part III: clinical studies. *J Esthet Restor Dent* 2002;14(5):313-319.
- Manicone P, Iommetti P, Raffaelli L. An overview of zirconia ceramics: basic properties and clinical applications. *J Dent* 2007;35(11):819-826.
- Garvie RC, Hannink RH, Pascoe RT. Ceramic steel? *Nature* 1975 Dec 25;258:703-704.
- Piconi C, Maccauro G. Zirconia as a ceramic biomaterial. *Biomaterials* 1999 Jan;20(1):1-25.
- Gupta TK, Bechtold JH, Kuznickie RC, Cadoff LH, Rossing BR. Stabilization of tetragonal phase in polycrystalline zirconia. *J Materials Sci* 1977;12(12):2421-2426.
- Cales B, Stefani Y. Mechanical properties and surface analysis of retrieved zirconia femoral hip joint heads after an implantation time of 2 to 3 years. *J Materials Sci Materials in Med* 1994;5:376-380.
- Swab JJ. Low temperature degradation of Y-TZP materials. *J Materials Sci* 1991 Dec 15;26(24):6706-6714.
- Luthardt RG, Holzhueter M, Sandkuhl O, Herold V, Schnapp JD, Kuhlisch E, et al. Reliability and properties of ground Y-TZP-zirconia ceramics. *J Dent Res* 2002 Jul;81(7):487-491.
- Kosmac T, Oblak C, Jevnikar P, Funduk N, Marion L. The effect of surface grinding and sandblasting on flexural strength and reliability of Y-TZP zirconia ceramic. *Dent Materials* 1999 Nov;15(9):426-433.
- Bukat A, Fassina P, Greco F, Piantelli S, Piconi C, Zaghini C. GSP-made Ca-PSZ for Biomedical Applications. 3rd Symposium on ceramics in medicine 1990.
- Dion I, Bordenave L, Levebre F. Physicochemistry and cytotoxicity of ceramics. *J Materials Sci Materials in Med* 1994 Jan;5(1):18-24.
- Torricelli P, Verne E, Brovarone CV, Appendino P, Rustichelli F, Krajewski A, et al. Biological glass coating on ceramic materials: in vitro evaluation using primary osteoblast cultures from healthy and osteopenic rat bone. *Biomaterials* 2001 Sep;22(18):2535-2543.
- Lohmann CH, Dean DD, Koster G, Casasola D, Buchhorn GH, Fink U, et al. Ceramic and PMMA particles differentially affect osteoblast phenotype. *Biomaterials* 2002 Apr;23(8):1855-1863.
- Warashina H, Sakano S, Kitamura S, Yamauchi KI, Yamaguchi J, Ishiguro N, et al. Biological reaction to alumina, zirconia, titanium and polyethylene particles implanted onto murine calvaria. *Biomaterials* 2003 Sep;24(21):3655-3661.

21. Degidi M, Artese L, Scarano A, Perrotti V, Gehrke P, Piattelli A. Inflammatory infiltrate, microvessel density, nitric oxide synthase expression, vascular endothelial growth factor expression, and proliferative activity in peri-implant soft tissues around titanium and zirconium oxide healing caps. *J Periodontol* 2006 Jan;77(1):73-80.
22. Carinci F, Pezzetti F, Volinia S, Francioso F, Arcelli D, Farina E, et al. Zirconium oxide: analysis of MG63 osteoblast-like cell response by means of a microarray technology. *Biomaterials* 2004;25:215-228.
23. Raigrodski AJ. Contemporary materials and technologies for all-ceramic fixed partial dentures: a review of the literature. *J Prosthet Dentist* 2004 Dec;92(6):557-562.
24. Tinschert J, Natt G, Mohrbotter N, Spiekermann H, Schulze KA. Lifetime of alumina- and zirconia ceramics used for crown and bridge restorations. *J Biomedical Materials Res Part B Applied Biomaterials* 2007 Feb;80(2):317-321.
25. Oh WS, Anusavice KJ. Effect of connector design on the fracture resistance of all-ceramic fixed partial dentures. *J Prosthet Dentist* 2002 May;87(5):536-542.
26. Studart AR, Filser F, Kocher P, Gauckler LJ. Fatigue of zirconia under cyclic loading in water and its implications for the design of dental bridges. *Dent Materials* 2007;23:106-114.
27. Gamborena I, Blatz MB. A clinical guidelines to predictable esthetics with Zirconium oxide ceramic restorations. *Quintessence of Dent Technol* 2006;29:11-23.
28. Piwowarczyk A, Ottl P, Lauer HC, Kuretzky T. A clinical report and overview of scientific studies and clinical procedures conducted on the 3M ESPE Lava all-ceramic system. *J Prosthodont* 2005 Mar;14(1):39-45.
29. Potiket N, Chiche G, Finger IM. *In vitro* fracture strength of teeth restored with different all-ceramic crown systems. *J Prosthetic Dentist* 2004 Nov;92(5):491-495.
30. Coli P, Karlsson S. Fit of a new pressure-sintered zirconium dioxide coping. *Int J Prosthodont* 2004;17(1):59-64.
31. Larsson C, Vult von Steyern P, Sunzel B, Nilner K. All ceramic two- to five-unit implant-supported reconstructions. A randomized, prospective clinical trial. *Swedish Dent J* 2006;30(2):45-53.
32. Jang GW, Kim HS, et al. Fracture strength and mechanism of dental ceramic crown with zirconia thickness. *J Proeng* 2011;10:1556-1560.
33. Schmitter M, Mueller D, Ruer S. Chipping behavior of all ceramic crowns with zirconia framework and CAD/CAM manufactured veneer. *J Dent* 2011;40(2):154-162.
34. Denry I, Kelly J. State of the art of zirconia for dental applications. *J Dent* 2008;24(3):299-305.