Tree Genetics & Genomes (2009) 5:211–223 DOI 10.1007/s11295-008-0176-7

ORIGINAL PAPER

Development and test of 21 multiplex PCRs composed of SSRs spanning most of the apple genome

A. Patocchi · F. Fernández-Fernández · K. Evans ·

D. Gobbin · F. Rezzonico · A. Boudichevskaia ·

F. Dunemann · M. Stankiewicz-Kosyl ·

F. Mathis-Jeanneteau · C. E. Durel · L. Gianfranceschi ·

F. Costa · C. Toller · V. Cova · D. Mott · M. Komjanc ·

E. Barbaro · L. Kodde · E. Rikkerink · C. Gessler ·

W. E. van de Weg

Received: 29 February 2008 / Revised: 29 May 2008 / Accepted: 11 June 2008 / Published online: 29 August 2008 © Springer-Verlag 2008

Abstract A series of 21 multiplex (MP) polymerase chain reactions containing simple sequence repeat (SSR) markers spanning most of the apple genome has been developed. Eighty-eight SSR markers, well distributed over all 17 linkage groups (LGs), have been selected. Eighty-four of them were included in 21 different MPs while four could not be included in any MPs. The 21 MPs were then used to

Electronic supplementary material The online version of this article (doi:10.1007/s11295-008-0176-7) contains supplementary material, which is available to authorized users.

Communicated by E. Dirlewanger

A. Patocchi (⊠) · D. Gobbin · F. Rezzonico · C. Gessler
Plant Pathology, Institute of Integrative Biology (IBZ),
ETH Zurich,
8092 Zurich, Switzerland
e-mail: andrea.patocchi@acw.admin.ch

F. Fernández-Fernández · K. Evans East Malling Research, New Road, East Malling, Kent ME19 6BJ England, UK

A. Boudichevskaia · F. Dunemann Bundesforschungsinstitut für Kulturpflanzen, Julius Kühn-Institut (JKI), Pillnitzer Platz 3a, 01326 Dresden, Germany

M. Stankiewicz-Kosyl Laboratory of Basic Research in Horticulture, Faculty of Horticulture and Landscape Architecture, Warsaw Agricultural University (WAU), ul. Nowoursynowska 166, 02-787 Warsaw, Poland genotype approximately 2,000 DNA samples from the European High-quality Disease-Resistant Apples for a Sustainable agriculture project. Two SSRs (CH01d03 and NZAL08) were discarded at an early stage as they did not produce stable amplifications in the MPs, while the scoring of the multilocus (ML) SSR Hi07d11 and CN44794 was too complex for large-scale genotyping. The testing of the

F. Mathis-Jeanneteau · C. E. Durel UMR1259 Genetics and Horticulture (GenHort), Institut National de la Recherche Agronomique (INRA), BP 60057, 49071 Beaucouzé, France

L. Gianfranceschi Department of Biomolecular Sciences and Biotechnology, University of Milan, Via Celoria 26, 20133 Milan, Italy

F. Costa Department of Fruit Tree and Woody Plant Sciences, University of Bologna, 40127 Bologna, Italy

C. Toller · V. Cova · D. Mott · M. Komjanc · E. Barbaro Istituto Agrario di San Michele all'Adige, 38010 San Michele all'Adige, Trento, Italy

L. Kodde · W. E. van de WegDepartment of Biodiversity and Breeding,Plant Research International,P.O. Box 16, 6700 AA Wageningen, The Netherlands

remaining 80 SSRs over a large number of different genotypes allowed: (1) a better estimation of their level of polymorphism; as well as of (2) the size range of the alleles amplified; (3) the identification of additional unmapped loci of some ML SSRs; (4) the development of methods to assign alleles to the different loci of ML SSRs and (5) conditions at which an SSR previously described as ML would amplify alleles of a single locus to be determined. These data resulted in the selection of 75 SSRs out of the 80 that are well suited and recommended for large genotyping projects.

Keywords SSR · Multiplex PCR · Genotyping · Malus

Introduction

Recently, the interest in developing a set of well-scattered highly polymorphic molecular markers spanning the whole apple genome has increased. Such a set of markers are necessary for the generation of genetic maps to be used for the accurate identification of genes and quantitative trait loci (QTLs) by association studies as well as for "fast breeding" strategies which identify progeny plants with the lowest proportion of undesired genome inherited from a wild apple (Volz et al. 2007). A subset of these selected markers can also be used for the genome scanning approach (GSA; Patocchi and Gessler 2003) or used to fingerprint cultivars and selections.

The development of such a set of markers was necessitated by the European Project High-quality Disease-Resistant Apples for a Sustainable agriculture (HiDRAS), established in 2002 (Gianfranceschi and Soglio 2004), which aimed to identify genetic loci controlling apple fruit quality with the objective of increasing the fruit quality of disease-resistant apples and therefore their acceptability and dissemination, leading to a reduction in the use of fungicides. HiDRAS is based on an innovative approach of association of genotype with phenotype called pedigreebased analysis (PBA). PBA includes the analysis of the

E. Rikkerink

The Horticulture and Food Research Institute of New Zealand Ltd, Mt Albert Research Center, Private Bag, 92169 Auckland, New Zealand

Present address: A. Patocchi · F. Rezzonico Agroscope Changins-Wädenswil Research station, Plant protection, Phytopathology, Schloss, B. O. 185, 8820 Wädenswil, Switzerland segregation of specific chromosomal regions in genetically related cultivars, breeding selections, and small progenies through highly polymorphic codominant markers. This "identity by descent" approach allows the detection of QTLs (Bink et al. 2002; van de Weg et al. 2004).

To increase the efficiency and reduce the cost of this type of genetic study, this set of markers must be polymerase chain reaction (PCR)-based and suitable for organization into multiplexes (MPs). Simple sequence repeat (SSR) markers were our markers of choice as they meet these requirements and because they are also generally highly polymorphic, ideal for map alignments, and highly transferable between laboratories.

More than 300 SSRs have been developed in apple (Guilford et al. 1997; Gianfranceschi et al. 1998; Hokanson et al. 1998; Liebhard et al. 2002; Hemmat et al. 2003; Vinatzer et al. 2004, Silfverberg-Dilworth et al. 2006, Khan et al. 2007, Igarashi et al. 2008; Van Dyk 2008). In addition, the use of pear SSRs can be explored as it has been demonstrated that they can generally be transferred to the apple and vice versa (Yamamoto et al. 2002; Pierantoni et al. 2004; Silfverberg-Dilworth et al. 2006) due to the high level of synteny between the two species. Nearly all apple SSRs have been mapped; the majority of them being mapped in the "reference map of the apple" (Liebhard et al. 2003; Silfverberg-Dilworth et al. 2006). The position within the apple genome of the few that were not mapped in the reference map can be inferred from map alignments using SSR markers in common as references.

Silfverberg-Dilworth et al. (2006) estimated that approximately 100 single-locus markers are required to span the apple genome with an intramarker distance of approximately 15-20 cM and a maximal 10-cM distance from the linkage group ends. The known distribution of the SSRs in the apple genome (reference map) as well as small-scale studies of their level of polymorphism (Liebhard et al. 2002; Silfverberg-Dilworth et al. 2006) allowed the selection of a set of 84 SSRs that covers about 85% of the apple genome (Silfverberg-Dilworth et al. 2006). This set of SSRs was composed of: 75 single-locus or presumed single-locus SSRs, seven multilocus SSRs (CH01d03, CH03g12, CH04c06, CH04g09, CH04h02 Hi07d11, and HB03-SSR) where the position of all the loci is known, and two multilocus SSRs (CN493139 and Hi23g12) where the position of only one locus was known. For one out of the seven multilocus SSRs (CH03g12), both loci (LG1 and LG3) are in interesting positions (i.e., located at extremities of a linkage group). For the other six SSRs, only one of the loci was of main interest as their second locus maps close to a single-locus SSR, which was chosen in preference to the multilocus SSRs to represent their genomic region. For 16 regions of the genome, no SSRs were identified by Silfverberg-Dilworth et al. (2006).

In this paper, we present an extension of this core set of SSRs spanning most of the apple genome, the development of multiplex PCRs for high throughput application of these SSRs, and finally an evaluation of their performance carried out by screening them over about approximately 2,000 DNA samples of the HiDRAS project.

Materials and methods

Plant material

The HiDRAS plant material genotyped with the selected SSRs is composed of 27 genetically related progenies of variable sizes (on average, 51 individuals with a range of 26 to 98), for a total of 1,373 progeny plants. In addition, another 359 cultivars and breeding lines were analyzed, many of which were represented by replicated DNA samples from different countries (Table S1). These genotypes were part of the pedigrees of the populations analyzed or were so-called common progenitors (genotypes present at different partner sites) or new founders (e.g., heritage cultivars). Including internal controls and repeated samples, 1,997 samples were genotyped.

DNA extraction, quantification, and standardization

DNA was extracted using the Qiagen DNeasy (Qiagen, Hilden, Germany) kit by the owner of the plant material and sent to the Swiss Federal Institute of Technology Zürich. DNA concentration was determined with the Quant-iT PicoGreen dsDNA Assay Kit (Invitrogen, Basel, Switzerland), on a SPECTRAFluor Plus microplate reader (Tecan, Männedorf, Switzerland). DNA concentration of each sample was adjusted to 5 ng/µl. DNA samples were reorganized into 96well plates using a Genesis RSP Liquid-Handling robot (Tecan, Switzerland). On each plate, samples of the genotypes "Prima" and "Fiesta" were added in reference positions H6 and H12. Aliquots were then prepared using a Tecan Genesis RSP Liquid-Handling robot (Tecan, Switzerland) and shipped to the six partners responsible for the genotyping, defined in this paper as "genotypers".

SSRs spanning the whole apple genome

The set of genome-spanning SSRs examined is composed of 88 SSRs developed by different groups: Liebhard et al (2002, 54 SSRs) Silfverberg- Dilworth et al. (2006, 24 SSRs), Yamamoto et al. (2002; NH009 on LG13; NH029 on LG9), Hokanson et al. (1998; GD103 on LG5; GD147 on LG13) Guilford et al. (1997; NZ02b01 on LG15), Vinatzer et al. (2004; CH-Vf1 on LG1), Broggini et al. (submitted manuscript; HB11-SSR on LG1 and HB03-SSR on LG6),

Rikkerink (unpublished; NZAL08 on LG3, NZEST67774 on LG4). The primer sequences of the latter two SSRs are: NZAL08-F GGC ACA AGC ACA AGG AAA CA, NZAL08-R GTT TGA GCC AGT CCA TTT TTC CCT AT for SSR NZAL08 and NZEST67774-F CCC GAC TGA CTG AAC CTT TT, NZEST67774-R GTT TCC GTG GAA GTG GAG TGA AG for SSR NZEST67774 (Tm of 60°C).

Assembly of multiplex PCRs for ABI sequencers

The assembly of MP PCRs of the 61 SSRs for genotyping on ABI sequencers was performed at Plant Research International (PRI) on an ABI PRISM[®] 377 DNA Sequencer (Applied Biosystems). As far as possible, SSRs that mapped to the same LG were assigned to the same MP. Different combinations of primer concentrations were tested until all SSRs of a MP produced similar quantities of amplicons. SSR markers giving unsatisfactory results were moved to other MP PCRs. SSR markers with primers labeled with the same fluorochromes were assigned to the same MP only if the ranges of their allele length were well separated. Further optimization for local conditions was performed at the site of each genotyper, who all used a capillary system. Protocols of the MPs at each genotyping site are presented in Table 1 and Table S1 (full information).

Assembly of multiplex PCRs for LI-COR sequencer

Twenty-seven SSRs were assigned to Julius Kuehn-Institute, the only partner working with an automatic dual-laser DNA sequencer LI-COR 4200 with two different wavelengths (IRD 700–IRD 800). As with the SSRs used on ABI sequencers, the SSRs analyzed on LI-COR were first pooled by linkage group. Protocols for each MP were obtained using the Qiagen Multiplex PCR kit. Due to different intensities of the alleles of different loci as well as to differences in the labeling intensity of the IRD-labeled primers, the primer concentrations had to be specifically adapted for each MP PCR. Protocols for the MP reactions are presented in Table 1.

Nomenclature of the multiplex PCR

The names of the MP PCRs were generated as follows. The first two letters of each MP are "Hi" indicating that the MPs have been developed within the frame of the HiDRAS project. The next number indicates the linkage group represented by the highest number of SSRs in the MP. If two linkage groups are represented by an identical number of SSRs, the number of both linkage groups is indicated separated by a dash. If the resulting name was identical to that of another MP, the letters "a" and "b" were added at the end of the name.

$ \begin{array}{cccccccccccccccccccccccccccccccccccc$		SSKs	ΓG	Fragment analysis	For primer conc (µM)	Rev primer conc (µM)	Labeling	DNA (ng)	Cycling conditions	Vol (µl)
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$	i1/6	Hi02c07	1	ABI 3100	0.08	0.08	PET	10	94°C 5', 10 cycles: 94°C 30",	10
$ \begin{array}{c ccccccccccccccccccccccccccccccccccc$		CH05a05	9		0.08	0.08	VIC		58–53°C 45" (dropping 0.5°C per	
		CH-Vfl	1		0.08	0.08	NED		cycle) and 72°C 1', 25 cycles: 94°C 30",	
		CH03d12	9		0.08	0.08	6FAM		53°C 45" and 72°C 1'; 72°C 15"	
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$	i2 ^a	CN493139SSR	2/2/5	ABI 3100	0.033	0.033	NED	10	94°C 2.5'; 33 cycles: 94°C 30",	10
$ \begin{array}{llllllllllllllllllllllllllllllllllll$		CH03d01	2		0.042	0.042	HEX		55°C 30" and 72°C 1'; 72°C 5'	
$ \begin{array}{llllllllllllllllllllllllllllllllllll$		CH02f06	2		0.045	0.045	6-FAM			
$ \begin{array}{lclcrcr} CH0367 & 3 & LLCOR 4200 (dual laser) & 0065 & 0.175 & RD 300 & 0 & 94°C 5.3 $ $ $ $ $ $ $ $ $ $ $ $ $ $ $ $ $ $ $$		CH05e03	2		0.068	0.068	HEX			
$ \begin{array}{llllllllllllllllllllllllllllllllllll$	i3	CH03g07	З	LI-COR 4200 (dual laser)	0.065	0.175	IRD 800	10	94°C 5', 38 cycles: 94°C 30", 60°C 30"	10
$ \begin{array}{llllllllllllllllllllllllllllllllllll$		Hi03d06	б		0.05	0.15	IRD 700		and 72°C 1', 72°C 10'	
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$		AU223657SSR	З		0.03	0.1	IRD 800			
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$	i4a ^b	CH02c02b	4	ABI 310	0.19	0.19	NED	10	94°C 2.5'; 5 cycles: 94°C 30", 65–61°C	15
$ \begin{array}{lcccccccccccccccccccccccccccccccccccc$		CH04e03	5		0.13	0.13	NED		for 1' (dropping 1°C per cycle) and 72°C	
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$		CH04e02	4		0.17	0.17	6-FAM		1', 30 cycles: 94°C 30", 60°C 1', 72°C 1'; 72°C 10'	
H23602 4 0.12 0.28 FET 0.12 0.28 FET NZ02001 15 0.12 0.28 6FAM 0.12 0.28 6FAM NZDS101 15 0.12 0.28 6FAM 0.12 0.28 6FAM NZEST6774 4 0.12 0.12 0.28 6FAM 0.12 0.28 6FAM NZEST6774 4 1 1 1 N Doyloging for per cycle) and 72°C NZEST6774 4 1 1 1 NED 0.12 0.28 6FAM NZEST6774 4 1 1 1 NED 0.12 0.28 6FAM NZEST677 5 ABI 3100 0.02 0.125 IRD 700 10 Identical to Hi3 H0402 5/7 ABI 3100 0.03 0.125 IRD 700 10 Identical to Hi3 H0402 5/7 ABI 3100 0.03 0.125 IRD 700 10 Identical to Hi3	i4h	CH02g01	13	ABI 3100	0.12	0.28	NED	10	Identical to Hi1/6	10
NZ02b01 15 0.12 0.28 6FAM CH05d02 4 0.12 0.28 VIC H122f12 5 ABI 3100 1.25 1.25 HEX 10 NZST6774 4 0.12 0.28 VIC 30° (dropping 1°C per cycle) and 72°C NZST6774 4 1 1 NED 94°C 5(.10 cycles: 94°C 30°, 55°C 30° and 72°C NZST6774 5 ABI 3100 0.02 0.1 RD 700 10 94°C 5(.10 cycles: 94°C 30°, 55°C 30° and 72°C Hi04a02 5?? ABI 3100 0.02 0.1 RD 700 10 17.25°C 55°C 30° and 72°C CH03a06 5 ABI 3100 0.02 0.12 RD 700 10 140°C 55°C 30° 55°C 30° 55°C 30° Hi03a10 7 2 125 RD 700 10 140°C 55°C 50° 55°C 30° 55°C 30° Hi03a10 7 10 0.02 0.125 RD 700 10 140°C 55°C 50° 50° 55°C 50° 50°C 50° 55°C 50° 55°C 50° Hi03a10 7 10 10 10 <td></td> <td>Hi23g02</td> <td>4</td> <td></td> <td>0.12</td> <td>0.28</td> <td>PET</td> <td></td> <td></td> <td></td>		Hi23g02	4		0.12	0.28	PET			
$ \begin{array}{lcccccccccccccccccccccccccccccccccccc$		NZ02b01	15		0.12	0.28	6FAM			
H:2f12 5 ABI 3100 1.25 1.25 HEX 10 94°C 5', 10 cycles: 94°C 30', 65–55°C NZEST67774 4 NED 30' (dropping 1°C per cycle) and 72°C H:04a08 5 LL-COR 4200 (dual laser) 0.02 0.1 RD 700 1', 25 cycles: 94°C 30', 55°C 30" and 72°C H:04a08 5 LL-COR 4200 (dual laser) 0.02 0.1 RD 700 10 1/25°C 50" and 72°C H:04a02 5/? ABI 3100 0.015 0.075 RD 700 10 1/25°C 50" and 72°C CH03a09 5 ABI 3100 0.02 0.1 RD 700 10 1/25°C 50" and 72°C CH03a00 5 ABI 3100 0.02 0.1 RD 700 10 1/25°C 50" and 72°C CH02b12 5 ABI 3100 0.020 0.040 VIC 10 94°C 25'. 33 cycles: 94°C 30", 56°C 30" MS0cg03 10 77°C 1', 72°C 5' RD 700 10 94°C 25'. 33 cycles: 94°C 30", 56°C 30" MS0cg03 10 0.030 0.09 0.009 NED		CH05d02	4		0.12	0.28	VIC			
NZEST6774 4 1 1 NED 30" (dropping 1°C per cycle) and 72°C Hi04a08 5 LL-COR 4200 (dual laser) 0.02 0.1 RD 700 1; 25 cycles: 94°C 30" and 72°C Hi04a08 5 LL-COR 4200 (dual laser) 0.02 0.1 RD 700 1 72°C 1; 72°C 5' Hi04a02 5/? ABI 3100 0.03 0.125 RD 700 10 Identical to H13 CH03a09 5 ABI 3100 0.03 0.125 RD 700 10 Identical to H13 CH02b12 5 ABI 3100 0.040 0.040 VIC 10 Identical to H13 K103a10 7 7 U 0.090 NED 0.040 0.040 NED 0.050 0.050 NED 0.075 SG 7.33 cycles: 94°C 30", 56°C 30" SG 7.30" Md 72°C 1', 72°C 5' NG 73°C 1', 72°C 5' NG 73°C 1', 72°C 5' NG 73°C 1', 72°C 5'	i4/5	Hi22f12	5	ABI 3100	1.25	1.25	HEX	10	94°C 5', 10 cycles: 94°C 30", 65–55°C	20
Hi04a08 5 Ll-COR 4200 (dual laser) 0.02 0.1 RD 700 17, 25 cycles: 94°C 30°, 55°C 30° and 72°C 1; 72°C 5° CH03a09 5 0.015 0.015 0.075 RD 700 10 Identical to Hi3 CH03a09 5 0.015 0.015 0.075 RD 700 10 Identical to Hi3 CH05e06 5 ABI 3100 0.03 0.125 RD 800 72°C 1; 72°C 5° CH05b05 5 ABI 3100 0.030 0.090 PET and 72°C 1; 72°C 5° MS06g03 10 7 0.090 0.090 0.090 PET and 72°C 1; 72°C 5° MS06g03 10 7 0.090 0.090 0.090 PET and 72°C 1; 72°C 5° MS06g03 10 0.018 6.FAM and 72°C 1; 72°C 5° and 72°C 1; 72°C 5° Art 72°C 1; 72°C 5° CH02b03b 10 0.005 0.018 6.FAM and 72°C 1; 72°C 5° Art 72°C 1; 72°C 5° CH02b03b 10 0.018 0.018 6.FAM and 72°C 1; 72°C 5° Art 72°C 1; 72°C 5° CH04660 10 0.018 <		NZEST67774	4		1	1	NED		30" (dropping 1°C per cycle) and 72°C	
$ \begin{array}{lcccccccccccccccccccccccccccccccccccc$									1', 25 cycles: 94°C 30", 55°C 30" and 72°C 1': 72°C 5'	
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$	i5	Hi04a08	5	LI-COR 4200 (dual laser)	0.02	0.1	IRD 700	10	Identical to Hi3	10
$ \begin{array}{llllllllllllllllllllllllllllllllllll$		CH03a09	5	·	0.015	0.075	IRD 700			
$ \begin{array}{lcccccccccccccccccccccccccccccccccccc$		Hi04d02	5/?		0.03	0.125	IRD 800			
$ \begin{array}{cccccccccccccccccccccccccccccccccccc$	i5/10	CH05e06	5	ABI 3100	0.040	0.040	VIC	10	94°C 2.5', 33 cycles: 94°C 30", 56°C 30"	10
Hi03a10 7 0.090 0.090 NED MS06g03 10 0.018 6-FAM MS06g03 10 0.018 6-FAM CH02b03b 10 0.018 0.018 6-FAM CH04e05 7/? LI-COR 4200 (dual laser) 0.02 0.125 IRD 800 10 CH04e05 7/? LI-COR 4200 (dual laser) 0.02 0.125 IRD 700 CH05c06 16 0.035 0.125 IRD 700 0.035 0.125 IRD 700 CH05c06 16 0.015 0.125 IRD 700 0.0125 IRD 700 CH01c06 8 LI-COR 4200 (dual laser) 0.015 0.125 IRD 700 CH01f09 8 0.015 0.125 IRD 700 0.016 Identical to Hi3 CH01h10 8 0.015 0.125 IRD 700 0.016 Identical to Hi3		CH02b12	5		0.090	0.090	PET		and 72°C 1', 72°C 5'	
MS06g03 10 0.018 0.018 6-FAM CH02b03b 10 0.065 0.065 6-FAM CH02b03b 10 0.065 0.065 6-FAM CH04e05 7/? L1-COR 4200 (dual laser) 0.02 0.125 IRD 800 10 CH04e05 7/? L1-COR 4200 (dual laser) 0.02 0.125 IRD 700 CH05c06 16 0.155 0.125 IRD 700 CH05c06 8 L1-COR 4200 (dual laser) 0.015 IRD 700 CH01c06 8 L1-COR 4200 (dual laser) 0.015 IRD 700 CH01c06 8 L1-COR 4200 (dual laser) 0.015 IRD 700 CH01c06 8 L1-COR 4200 (dual laser) 0.015 IRD 700 CH01c06 8 0.125 IRD 700 Identical to Hi3 CH01c06 8 0.125 IRD 700 Identical to Hi3		Hi03a10	7		0.090	0.090	NED			
CH02b03b 10 0.065 0.065 6-FAM CH04e05 7/? Ll-COR 4200 (dual laser) 0.02 0.125 IRD 800 10 Identical to Hi3 CH04e05 7/? Ll-COR 4200 (dual laser) 0.02 0.125 IRD 700 CH05e06 16 0.15 0.125 IRD 700 CH05a04 16 0.035 0.125 IRD 700 CH01e06 8 Ll-COR 4200 (dual laser) 0.015 0.1 CH01f09 8 0.1 0.015 0.1 IRD 700 CH01h10 8 0.1 0.225 IRD 700 10 Identical to Hi3		MS06g03	10		0.018	0.018	6-FAM			
CH04e05 7/? LL-COR 4200 (dual laser) 0.02 0.125 IRD 800 10 Identical to Hi3 CN444794SSR 7/7 0.15 0.325 IRD 700 10 Identical to Hi3 CN444794SSR 7/7 0.15 0.325 IRD 700 16 16 16 16 17 18 10 16 16 16 16 17 17 17 17 17 17 17 17 17 17 17 17 17 17 18 10 16 16 16 16 16 16 17 17 17 17 17 16 16 <td< td=""><td></td><td>CH02b03b</td><td>10</td><td></td><td>0.065</td><td>0.065</td><td>6-FAM</td><td></td><td></td><td></td></td<>		CH02b03b	10		0.065	0.065	6-FAM			
CN44794SR 7/7 0.15 0.325 IRD 700 CH05c06 16 0.035 0.125 IRD 700 CH05a04 16 0.4 0.125 IRD 700 CH01c06 8 LI-COR 4200 (dual laser) 0.015 0.1 IRD 800 10 Identical to Hi3 CH01h09 8 0.1 0.225 IRD 700 0.125 IRD 700 CH01h10 8 0.1 0.125 IRD 700 10 Identical to Hi3	i7/16	CH04e05	i/L	LI-COR 4200 (dual laser)	0.02	0.125	IRD 800	10	Identical to Hi3	10
CH05c06 16 0.035 0.125 IRD 700 CH05a04 16 0.4 0.125 IRD 700 CH01c06 8 L1-COR 4200 (dual laser) 0.015 0.1 IRD 800 10 Identical to Hi3 CH01h09 8 0.1 0.225 IRD 700 0.125 IRD 800 10 CH01h10 8 0.03 0.125 IRD 800 10 Identical to Hi3		CN444794SSR	L/L		0.15	0.325	IRD 700			
CH05a04 16 0.4 0.125 IRD 700 CH01c06 8 LI-COR 4200 (dual laser) 0.015 0.1 IRD 800 10 Identical to Hi3 CH01f09 8 0.1 0.225 IRD 700 0.125 IRD 700 CH01h10 8 0.1 0.125 IRD 800 10 Identical to Hi3		CH05c06	16		0.035	0.125	IRD 700			
CH01c06 8 LI-COR 4200 (dual laser) 0.015 0.1 IRD 800 10 Identical to Hi3 CH01f09 8 0.1 0.225 IRD 700 CH01h10 8 0.03 0.125 IRD 800		CH05a04	16		0.4	0.125	IRD 700			
8 0.1 0.225 IRD 8 0.03 0.125 IRD	i8	CH01c06	8	LI-COR 4200 (dual laser)	0.015	0.1	IRD 800	10	Identical to Hi3	10
0.03 0.125 IRD		CH01f09	8		0.1	0.225	IRD 700			
		CH01h10	8		0.03	0.125	IRD 800			

Table 1 PCR protocols of the multiplexes used to screen the HiDRAS germplasm

10	20	10	10	10	15	10	10	20	10 10
Identical to Hi3	Identical to Hi5/10	95°C 15', 30 cycles: 94°C 30", 60°C 1'30" and 72°C 1', 60°C 15', 72°C 15'	94°C 5'; 37 cycles: 94°C 30", 60°C 30" and 72°C 45"; 72°C 8'	Identical to Hi3	Identical to Hi1/6 Identical to Hi2/10	Identical to Hi1/6	Identical to Hi3	94°C 2.5', 33 cycles: 94°C 30″, 65°C 30″, 72°C 1'; 72°C 5'	Identical to Hi3 Identical to Hi11
10	10	10	10	10	10	10	10	Ś	10 10
IRD 700 IRD 800 IRD 700 IRD 800	NED 6-FAM HEX HEY	TAMRA HEX FLUO FLUO HFX	NED HEX HEX NED 6FAM	IRD 800 IRD 700 IRD 700	VIC NED 6FAM PET	of Alm NED VIC 6FAM PET	IRD 700 IRD 800 IRD 800	HEX 6-FAM NED NED	IRD 700 IRD 700 IRD 800 FLUO HEX FLUO
0.1 0.075 0.1 0.075	1 0.875 1.75	0.2 0.2 0.0 0.2 0.0 0.0 0.0 0.0 0.0 0.0	0.093 0.240 0.200 0.187 0.130	0.125 0.125 0.125 0.125	0.17 0.17 0.17 0.17 0.17	0.09 0.28 0.28 0.28	0.1 0.075 0.075	0.15 0.15 0.12 0.12	0.1 0.1 0.15 0.2 0.2 0.2
0.02 0.015 0.03 0.015	1 0.875 1.75	0.2 2 2 2 0	0.02 0.240 0.200 0.187 0.187	0.03 0.02 0.02	0.17 0.17 0.17 0.17 0.17	0.00 0.12 0.12 0.12 0.12	0.04 0.03 0.015	0.15 0.15 0.12 0.12	0.015 0.02 0.04 0.2 0.2 0.2
LI-COR 4200 (dual laser)	ABI 3100	ABI 3100	ABI 3100	LI-COR 4200 (dual laser)	ABI 3100	ABI 3100	LI-COR 4200 (dual laser)	ABI 3100	LI-COR 4200 (dual laser) ABI 3100
6 6 6 6	10 10 10/5	2/01 11 11 11	10/17 12 4/12 12 12	1 2 2 2	113 113 113	0 5 1 1	$\begin{array}{c}1\\1\\4\\1\end{array}$	15 15 15	16 16 16 17 17 17
CH01f03b Hi05e07 CH01h02 CN444542SSR	CH03d11 CH02b07 CH02c11 CH04c00	CH04h02 CH04h02 H107d11 CH02d08 CH04a12 H116d07	CH05404 CH05404 CH01403 CH04402 CH04g04 CH01501	CH01g12 CH03c02 Hi07f01	Hi04g05 CH03e03 CH05h05 GD147 HP03 SSP ^d	CH05704 GD103 NH009b CH05g08	CH01g05 CH04c07 MDAJ761SSR	CH02c09 CH01d08 CH02d11 Hi03g06	CH02a03 Hi04e04 CH04f10 AT000174SSR Hi03c05 Hi02f12
Hi9	Hi10 ^c	Hill	Hi12a	Hi12b	Hi13a	Hi13b	Hi14	Hi15	Hi16 Hi17

J	
-	
e	
Tabl	
F	

Table 1 (continued)	ontinued)								
MP name	SSRs	ΓG	Fragment analysis	For primer conc (µM)		Labeling	DNA (ng)	Rev primer Labeling DNA (ng) Cycling conditions conc (µM)	Vol (µl)
	CH05d08	9/17		0.2	0.2	HEX			
	Hi07h02	17		0.2	0.2	FLUO			
	CH04g07	11		0.2	0.2	TAMRA			
^a CH02g01 v ^b CH05d02 ; ^c PCRs perfe ^d PCRs perfe	¹ CH02g01 worked in this multiplex at PRL CH05d02 and Hi23g02 worked in this mu PCRs performed separately for CH04g09 PCRs performed separately. Two microlite	ltiplex at PI ed in this r or CH04g09 Two microli	^a CH02g01 worked in this multiplex at PRI, CH02g01 moved to MP Hi4b for germplasm screen ^b CH05d02 and Hi23g02 worked in this multiplex at PRI; both moved to MP Hi4b for germplasm screen ^e PCRs performed separately for CH04g09 ^d PCRs performed separately. Two microliters HB03AT loaded with 1 μ l MPHi13	b for germplasr MP Hi4b for g MPHi13	n screen ermplasm screei	- г			

SSR, MP name name of the multiplex PCR, SSRs SSR markers present in the multiplex, LG linkage group on which the

reaction, Cycling conditions PCR profile used for the specific multiplex, Vol final volume of the PCR

SSRs have been mapped, Fragment analysis indicates the instrument used to separate the amplicons, For and Rev primer conc (µM) concentration of the different primers in the multiplex PCRs,

amount of DNA used in the PCR

the primer, DNA (ng)

label

5

Labeling fluorochrome used

eaction.

? Unknown map position of the second locus of the presumed multilocus \$

🖉 Springer

Division of multiplexes

Multiplexes were divided over the six genotyper sites, so that each multiplex was applied at a single site. In this way, we avoided inconsistencies in allele sizing that easily occur when data for the same SSR come from different laboratories. ABI capillary platforms were used at five of the six sites, while a LI-COR was used at the sixth site.

Fragment analysis

For fragment analysis performed on the ABI PRISM® 3100 or 310 DNA capillary sequencer (Applied Biosystems), 1 µl of PCR product was mixed with 10 µl deionized formamide and 0.2-µl 500-LIZ or 500-ROX ladder (Applied Biosystems). Chromatographs were generated using Genescan 3.7 software and SSR fragment lengths were scored with Genotyper 3.6 (Applied Biosystems).

For fragment analysis performed on the LI-COR sequencer, 2 µl of the multiplex PCR reactions was diluted with 15-µl formamide loading buffer and denatured for 3 min at 95°C. A 0.5-µl diluted reaction was loaded on the gel. As a sizing standard, a 50-350 DNA ladder composed of 14 IRD-labeled fragments (LI-COR) was loaded separately in every ninth lane (five lanes per gel total). The detection of alleles of SSRs analyzed on LI-COR 4200 was performed using the fragment analysis software SAGA (LI-COR). Analysis was initially performed automatically; however, for a large part of the data, a manual scoring was necessary.

Verification of the consistencies of allele sizing

Following verification of nonsignificant shift of the allele sizes of the genotypes of reference "Prima" and "Fiesta" between plates, the allele sizes were rounded to the most appropriate integer according to the length of the repeated motif (2 or 3nt repeat). Rounded values were verified for consistency of the data within and among pedigrees as described in Patocchi et al. (2007). In brief, the procedure consisted of firstly checking if a progeny of a cross inherited an allowed combination of alleles from its two parents using an excel macro called "Gob validator" (D. Gobbin, unpublished). Then, consistency among all the connected HiDRAS pedigrees was verified using the statistical software FlexQTL[™] (www.flexqtl.nl). Then, the parts of the pedigrees generating inconsistencies were visualized using the software Pedimap (Voorrips 2007). Finally, the electropherograms of the plants involved in the inconsistencies were checked, and, if justified, the data scores were corrected. If scores were changed, another check with the software FlexQTL[™] was performed. This cycle of the validation procedure was performed until no improvement of data was possible.

Results

Modification of the set of SSRs spanning the apple genome

A core set of 84 SSRs spanning about 85% of the apple genome was already proposed by Silfverberg-Dilworth et al. (2006). One SSR in this set (CH03b06, LG15) showed difficulties in multiplexing. As NZ02b01 (SSR of the core set of Silfverberg-Dilworth et al. 2006) mapped at only 3 cM from CH03b06, this SSR was not replaced. Five SSRs were

added as they mapped on unrepresented regions: CH02b12 (LG5, Liebhard et al. 2002), GD103 (LG5) and GD147 (LG13; Hokanson et al. 1998), and NZAL08 (LG3) and NZEST67774 (LG4; Rikkerink, unpublished).

The current core set is now composed of 88 SSRs representing 89 regions (Fig. 1). For 13 regions, no SSR has yet been identified. As 101 single-locus SSRs are necessary to span the whole apple genome (to meet the criteria explained in the "Introduction"), the present coverage of the genome is approximately 89%.

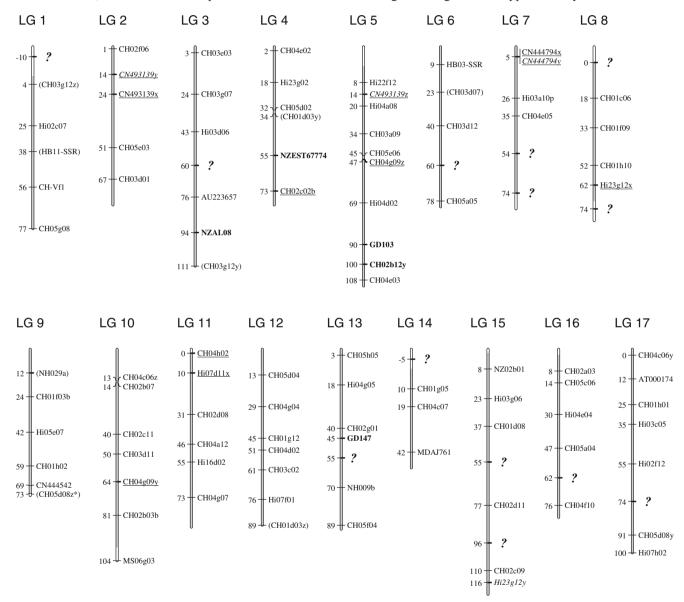


Fig. 1 Coverage of the apple genome with the selected set of SSRs. Five SSRs (in *bold*) have been added to the set selected by Silfverberg-Dilworth et al. (2006); one SSR (CH03b06, LG15) has been removed. Regions for which no SSR is currently available are indicated by *question marks*. SSRs of the set which could not be used to screen the germplasm set (see text) are indicated in *brackets*. Loci of the five SSRs for which scoring was judged too complex for large-scale applications (CN493139, CH04g09, CH04h02, Hi07d11, and

CN444794) or loci judged not sufficiently polymorphic (CH02c02b and Hi23g12) are *underlined*. These SSRs need to be replaced as soon as alternatives are available. Additional loci of the SSRs CN493139 (LG2 and 5) and Hi23g12 (LG15) identified in this study are indicated with *italics*. *Asterisks*, alleles of locus z (LG9) of multilocus SSR CH05d08 (LG9/17) were not amplified under the conditions used in this study (see Table 2)

Multiplex PCRs

Twenty-one MP PCRs containing in total 82 out of the 88 selected SSR markers have been developed. Thirteen MP PCRs, containing in total 56 SSR markers, were developed and tested on ABI sequencers (ABI-MPs), while eight MP PCRs, containing 27 SSRs, were developed and tested with a LI-COR sequencer (Li-MPs; Table 1 and Table S2). Four SSRs, HB11-SSR, CH03g12, CH03d07, and NH029 did not perform well in any MP and have not been further tested. SSR NZAL08 performed well under the conditions at PRI but failed under the conditions at the site of the genotyper (therefore, also considered not included in a MP). In addition, SSR CH01d03 (MP Hi12a) did not give stable amplifications and has not been further considered.

The highest number of SSR markers (six) placed in MPs was achieved by MPs Hi12a and Hi17. In two MPs, Hi10 and Hi13a, one SSR could not be placed directly in the MP (CH04g09 in Hi10 and HB06-SSR in Hi13); however, following an independent PCR, the products could be added to the MP for fragment analysis.

MPs Hi13b and Hi4b (four SSRs each) allowed the highest efficiency as their amplicons could frequently be combined for fragment analysis. Their SSRs generally amplified alleles in well-separated allele size ranges. For the screening of the segregating populations with these two MPs, it was first ascertained whether the alleles of the parents overlapped and, if this was not the case, the PCR products of both MPs were mixed together before performing fragment analysis. Therefore, within a single run, alleles of eight loci could be separated.

Three out of the eight LI-COR-MPs contained four SSRs, while the other five MPs each contained three SSRs. Using the Qiagen Multiplex PCR kit, it was possible to establish a robust protocol for each multiplex. None of the LI-COR SSR markers had to be amplified on its own.

Multilocus SSRs and assignment of alleles to the different loci

Although single-locus SSRs were preferred to multilocus SSRs (ML) in the selection of the set, for some genomic regions, the only available SSRs were multilocus (CH04c06, CN493139, CH04g09, HB03-SSR, GD103, CH04h02 Hi07d11, and Hi23g12) or were presumed to be multilocus (CH02b12, CH05d08, Hi04d02, CH04e05, CN444794). Therefore, it was necessary to assign the alleles to the different loci in order for these SSRs to be used for the "identity by descent" approach.

SSRs CH02b12, CH04e05, Hi04d02 (all presumed ML), and CH05d08 (previously proven to be ML) only amplified one locus in the conditions used in this study and therefore no splitting of the data was necessary for these SSRs. For the other SSRs, CH04c06, HB03-SSR, GD103, and Hi23g12, a simple procedure was developed to assign the alleles to the different loci (Table 2) and this was applicable also for large-scale genotyping efforts. For SSRs CN493139, CH04g09, CH04h02, CN444794, and Hi07d11, adequate

Table 2 Description of the procedure applied to identify the alleles of the different loci amplified by multilocus and presumed multilocus SSRs

SSR	MP	LGs	Procedure
CH02b12	Hi5/10	5/?	Under the conditions used only LG5 alleles amplified
CH04e05	Hi7/16	7/?	Under the conditions used only LG7 alleles amplified
CH05d08	Hi17	9/17	Under the conditions used only LG17 alleles amplified
Hi04d02	Hi5	5/?	Under the conditions used only LG5 alleles amplified
GD103	Hi13b	5/10	Simple, alleles >90 bp belong to locus z (LG5), alleles <90 bp to locus y (LG10)
HB03-SSR	Hi13a	6/16	Simple, alleles >350 bp belong to locus y (LG6), alleles <350 bp to locus x (LG16)
Hi23g12	Hi8	8/15	Simple, the alleles 221 and 224 bp belong to locus x (LG 8), alleles >233 bp belong to locus y (LG15)
CH04c06	Hi12	10/17	Simple, intensity difference, only two alleles found for locus z on LG10 (157 and 171 bp); however, allele ranges overlap; therefore, this SSR is not fully adequate for high throughput genotyping
CH04g09	Hi10	10/5	Complicated, overlapping size ranges; only the z locus can be scored and distinguished from the y locus due to great differences in amplification efficiency
CH04h02	Hi11	11/11	Allele assignment is relatively easy, but adequate scoring is difficult due to a combination of overlapping size ranges and large differences in amplification efficiency; use possible only in some segregating populations and using flanking markers (mapping)
CN444794	Hi7/16	7/7	Complicated, use possible only in some segregating populations
CN493139	Hi2	2/2/5	Complicated, use possible only in some segregating populations
Hi07d11	Hi11	11/?	Complicated, too many unspecific amplicons for use

? Presumed multilocus SSR (position of the second locus unknown), MP name of the multiplex, LGs linkage groups carrying the loci amplified by the multilocus SSR, Procedure method used to assign the alleles to the different loci

scoring and assignment of alleles to the corresponding loci was too complex for large-scale genotyping efforts due to a combination of overlapping allele ranges and differences in amplification efficiency between the loci. For the first three SSRs, the splitting of the data was only possible in certain crosses and circumstances, and this was considered too time-consuming for CN444794 and Hi07d11 and therefore not undertaken in this study. Data could therefore be produced for 80 SSRs (85 loci).

The assignment of alleles to the different loci allowed the mapping of additional loci from three ML SSRs (CN493139, Hi23g12, and CN444794) that were not identified by Silfverberg-Dilworth et al. (2006; Fig. 1). SSR CN493139 amplified three loci; the first additional locus (y) is on LG2 at 10 cM proximal from the known locus x of CN493139, while the second additional locus (z) was identified on LG5 between SSR markers Hi22f12 and Hi04a08. The second locus of SSR Hi23g12 has been mapped at the distal end of LG15 at approximately 6 cM from marker CH02c09. The second locus of CN444794 is a copy of the first locus at exactly the same position (top of LG7; Fig. 1).

Degree of polymorphism and size ranges of alleles of the tested SSRs

Screening the 80 SSRs over a much wider germplasm collection than the eight to nine cultivars tested during the development of the SSRs allowed a better estimation of their degree of polymorphism (Table 3). The average number of alleles per locus (85 loci) was 13.3. The most polymorphic single-locus SSRs were CH03d12 (24 alleles), CH04f10 (22 alleles), and Hi05e07 (20 alleles). The single-locus SSR with the lowest level of polymorphism was CH02c02b (four alleles). The ML SSR Hi23g12 was selected for its position on LG8. Regrettably, this locus proved to have a low level of polymorphism (three alleles), which was not evident in Silfververg-Dilworth et al. (2006).

The large-scale testing of the selected SSRs also allowed a better estimation of the allele size range that can be expected from each specific SSR (or each locus in case of ML SSR; Table 3). Direct comparisons between our results and those of Liebhard et al. (2002) and Silfverberg-Dilworth et al. (2006) are feasible, as all their reference cultivars except for "Nova Easygro" were tested by us. Such comparisons showed various differences in relative allele sizes, e.g., a difference of 10 nt between two alleles observed by us was previously reported as a 12-nt difference. As such differences can easily occur due to different fragment analysis and sizing methods and as they have little substantial impact on the estimation of the degree of polymorphism and size ranges of alleles, we did not elaborate on this.

Discussion

Set of SSRs

Out of the initial 88 SSRs that had been chosen for the core set of genome-covering SSR, 82 were extensively analyzed. This analysis allowed the identification of the 75 SSR markers that are well suited for large genotyping projects. These SSRs are well distributed within the apple genome, give robust amplifications in the MPs, and are sufficiently polymorphic and, for the few ML SSRs, the identification of the alleles belonging to the different loci is simple. We suggest that seven of the 82 tested SSRs need to be replaced as soon as better options become available (CH02c02b, Hi23g12, CN493139, CH04g09, CH04h02, Hi07d11, and CN444794). SSRs Hi23g12 (locus on LG8) and CH02c02b were discovered to have an extremely low level of polymorphism. Hi23g12 was included in the set because of its strategic position at the distal part of LG8 where no other SSRs were available, even though the presumed ML SSR Hi23g12 was already suspected to have low polymorphism (five alleles in total for two loci identified in nine genotypes, Silfverberg-Dilworth et al. 2006). The fact that the single-locus SSR CH02c02b amplifies only four different alleles in this study is somewhat surprising because the same SSR was reported to amplify five different alleles from as few as eight cultivars (Liebhard et al. 2002). This difference is probably caused by a combination of suboptimal PCR conditions and underloading of the capillaries in the current study, resulting in inability to detect the larger alleles.

SSRs CN493139, CN444794, CH04g09, CH04h02, and Hi07d11 are all ML SSRs with alleles of the different loci having similar sizes, making the assignment of alleles to the loci labor-intensive and, moreover, making the scoring of individual alleles prone to errors. All these SSRs were chosen because of their strategic position and lack of better SSRs in terms of quality of the amplifications and expected degree of polymorphism. The use of these five ML SSRs in large genotyping projects is not advised.

For full coverage of the apple genome, around 101 singlelocus SSRs are required. Currently, no SSRs are available for 13 loci (Fig. 1). Including the five ML SSRs found unsuitable within this project, new SSRs for 18 loci have to be found (Fig. 1). A new series of apple SSRs have been recently developed and mapped (Igarashi et al. 2008; Van Dyk 2008) some of which are located in the current gaps (e.g., LG15). New pear SSRs have also been developed (Fernández-Fernández et al. 2006; Inoue et al. 2007) and others are in development (Yamamoto T. personal communication). Once these SSRs have been mapped, making use of the synteny between apple and pear, it should be ascertained whether any of these SSRs could fill the gaps.

						,			,																		
Locus name	ΓG	No of	No of	Range (bp) ^b		Allele size ^{b,c}																					I
		alleles ^a	alleles ^b																								
AT000174-SSR	17	9	14	181-211	llun	181	185	187	189			195 19 7	7 199	9 201	203	205	211										1
AU223657-SSR	3	9	8	219–233	llun	219	221	223	225	229 2	231 23	233															
CH01c06	8	9	15	150-190	150	152	154	156	158			164 168				186	188	190									
CH01d08	15	9	16	244–299	244	246	248	250	254		259 26	261 265	5 273	3 277	279	285	287	295	299								
CH01f03b	6	7	17	137-193	137	139	141	143	145		,	159 161				177	179	181	183	193							
CH01f09	8	7	16	113-162	113	115	121	123	125	127 1	129 13	131 133	3 135	5 137		143	153	160	162								
CH01g05	14	9	19	138-186	lluu	138	140	142	144	146 1	150 15	152 154	4 156	5 158		164	166	168	172	174	184	186					
CH01g12	12	9	18	106 - 186	106	108	011	112	132	134 1	136 13	138 140	0 146	5 150		154	156	178	182	184	186						
CH01h01	17	9	13	110-145	110	118	120	122	124	126 1	128 13	130 132	2 134			145											
CH01h02	6	9	13	226–254	lluu	226	230	232	236		242 24	244 246	6 248			254											
CH01h10	8	5	15	92–221	92	98	100	102	104				4 116			134	194	221									
CH02a03	16	7	15	100 - 170	100	126	136	140	144	146 1	148 IS	150 152		4 156	164	166	168	170									
CH02b03b	10	8	13	75-107	null	75	79	81	89		93 9	95 97	7 99	101 €	103	107											
CH02b07	10	7	15	98-140	98	100	102	104	106	108 1	110 11	112 116	6 118	3 120	122	124	126	140									
CH02b12y	5	13	11	118-146	null	118	128	130	132	134 1	138 14	140 142	2 144	4 146													
CH02c02b	4	5	4	112-122	lluu	112	911	122																			
CH02c09	15	9	7	238-262	238	244	246	248	250		260																
CH02c11	10	7	14	202-233	202	208	210	212		216 2		222 22				231	233										
CH02d08	11	7	15	217-269	lluu	217	221	223	225	227 2	229 23	231 233	3 235	5 253		257	261	269									
CH02d11	15	7	14	118-156	118	122	124	126				134 136			152	154	156										
CH02f06	2	7	13	137-164	null	137	139	143	145			154 156				164											
CH02g01	13	5	15	184–236	llun	184	194	196	200		204 21	212 214	4 216	5 218	222	224	230	236									
CH03a09	5	9	12	121-145	121	125	127	129	131	133 1	135 13	137 139	9 141	1 143	145												
CH03c02	12	5	12	106-162	106	114	120	122	124	126 1	128 13	132 134	4 138	8 142	162												
CH03d01	2	7	11	94-118	lluu	94	98	104	106	108 1	110 11	112 114	4 116	6 118													
CH03d11	10	9	12	108-176	llun	108	110	114	911			122 124			176												
CH03d12	9	7	24	96-171	llun	96	102	104			110 11	112 114	4 116	5 118	120	126	128	130	138	140	142	146 1	148 1	150 1:	152 15	154 10	171
CH03e03	б	9	11	188-218	llun	188	191	198				206 208															
CH03g07	3	5	14	119-179	llun	119	121	123	125							177	179										
CH04a12	11	8	13	165-197	llun	165	171	173	175							197											
CH04c06y	17	8	12	167-194	167	175	177	179	181	184 1	186 18	187 188	8 190	0 192	194												
CH04c06z	10		7	158-171	158	121																					
CH04c07	14	8	16	96–134	lluu	96	98	100	106		_	112 114	4 116	5 118	120	122	130	132	134								
CH04d02	12	Э	7	123-151	lluu	123	125	131	135		151																
CH04e02	4	9	6	141 - 168	lluu	141	143	145	152			160 16	8														
CH04e03	5	11	19	182-236	lluu	182	187	189	191	193 1	195 19	197 19	9 201	1 203		207	209	210	212	214	224	236					
CH04e05	7	8	17	175-227	175	189	197	199	201	203 2	205 20				217	219	221	223	225	227							
CH04f10	16	6	22	166-316	lluu	166	172	174	186	188 1	194 22	224 226	6 230	0 234		238	240	242	244	246	248	250 2	252 2	264 3	316		
CH04g04	12	5	14	155-211	155	161	163	173	175		179 18	181 183	3 187	7 191		195	211										
CH04g07	Π	6	15	153-214	153	155	164	166	168	170 1					182	185	197	214									
CH04g09y	10	11	9	138-154	llun	138	143	145	152	154																	
CH04g09z	5		16	135-178	llun	135	137	144	146	148 l	150 15	152 155	5 161	1 167	169	172	174	176	178								

Table 3 SSR allele number and size identified within HiDRAS plant material, separated by locus

4				
254				
249				
246				
233				
221				232
214				230
210		223		169 228 254
209		211	172	2266
207	215	191	150 162 422	151 215 240 240 2222 2222 2220
1 206	213 126	207 187 162 162	148 153 154 154 154 154	149 213 250 255 255 255 256 279 270 270 270 271 271 271 271 271 272 273 274 274 274 275 276 277 278 279 279 270 271 271 272 273 274 274 274 274 274 274 274 274 274 274 </td
3 204	5 197 7 124	3 205 3 205 5 158 7 158	4 146 133 9 152 5 418	147 2147 2248 2251 2251 2251 2251 2251 2252 2252 2253 2254 2254 2255 2254 2254 2255 2255 2255 2255 2255 2255 2255 2255 2255 2354 2355 244 255
) 203		 223 203 203 182 181 181 	8 144 141 131 131 150 150	145 246 2246 2246 2273 2274 2274
7 199			4 188 0 142 1 164 1 164 1 129 8 412 8 412 8 150	1 143 9 2013 9 244 0 224 1 225 2 216 2 216 3 240 6 188 6 188 0 162 0 162 1 162
5 197			8 174 4 140 9 121 9 121 6 148 6 148	 217 219 133 139 141 133 139 141 133 149 240 240 240 241 214 216 215 215 255 254 260 214 215 215 216 214 216 215 214 216 214 216 214 216 214 216 214 216 214 216 216
1 195		• • • • • • • • • • • • • • • • • • • •	4 168 5 134 7 19 7 19 8 382 8 382 9 9 1 295	 13.9 13.9 13.9 13.9 13.9 13.9 13.9 13.9 13.9 13.4 13.4
161 6				213 214 131 133 214 214 224 226 223 226 223 226 223 237 223 247 252 247 252 247 252 247 252 247 252 247 252 247 252 247 252 254 252 254 252 ert genol
185 189	298 167 169 216 220 110 112	• • • • • • • •	,	227 2.15 127 131 131 132 132 238 2.12 236 2.28 236 2.28 255 2.29 164 167 164 167 154 167 152 154 178 158 152 2.22 154 2.26 232 2.24 170 178 152 154 152 2.22 248 2.50 different
				0
			142 1: 124 1: 124 1: 149 1: 135 1: 135 1: 135 1: 137 1: 13	[5] · · · · · · · · · · · · · · · · · · ·
			140 1. 122 1. 122 1. 122 1. 144 1. 144 1. 144 1. 1107 1. 133 1. 1	119/11 119/11 119/11 119/12 119/12 119/12 119/12 119/12 119/12 119/12 119/12 119/12 12 12 12 12 12 12 12 12 12 12 12 12 1
175 1			· · · · · · · · · · · · · · · · · · ·	 194 197 119 119 119 118 119 118 118 118 118 118 119 119 119 119 114 114
			-	null null null null 186 186 186 186 186 186 180 133 133 231 231 231 231 232 222 222 222
				10 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0
161–298	151–215 202–257 84–130	173–223 173–223 122–146 166–197 120–162 158–170 148–205 149–185	128-188 112-150 1136-164 134-149 134-149 123-172 103-153 103-153 353-422 353-422 353-422 106-150 129-149	115-169 115-169 115-250 176-250 176-250 180-254 180-232 180-232 180-233 180-213 197-213 220-263 231-245 231-245 231-224 231-224 155-192 231-224 231-224 231-224 231-224 231-2260 231-2260 231-2260 232-2260 222-20
_				ypes t
31	16 13 16	13 18 17 17 15 8 8 8 12 12 12	12 12 12 12 12 12 12 12 12 12 12 12 12 1	10 11 12 13 13 14 15 15 12 12 12 13 13 13 13 11 10 10 11 11 12 13 13 14 15 15 15 15 15 15 15 15 15 15
14	× • × •	4 5 6 8 0 1 0 4 4 4 4 4 4 4 4 4 4 4 4 4 4 4 4 4	n.a 88 8 6 6 6 6 6 6 6 6 6 6 6 6 6 6 6 6 6	° ° ° ° ° ° ° ° ° ° ° ° ° ° ° ° ° ° °
				8 diff
11	16 6 16	4 17 13 13 13 13	1 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0	10 11 15 13 13 14 14 15 15 16 16 17 11 17 11 17 18 18 18 18 18 18 18 18 18 18
CH04h02	CH05a04 CH05a05 CH05c06	CH05404 CH05408 CH05408 CH05603 CH05606 CH05f04 CH05g08 CH05b05 CH05b05	CH-Vf1 CN444542-SSR CN493139-SSRx CN493139-SSRz CN493139-SSRz GD103 GD147 HB03-SSRy Hi0267 Hi0267 Hi02612 Hi02612 Hi02612 Hi02612	nuoscos10100-10101Hi03d0615516115-515null115Hi03d0615516115-515null175Hi04d025316175-215null176Hi04d025316175-240null176Hi04d0213816175-240null176Hi04d021381613214-246218Hi04d021381613214-246218Hi04d02111015214-246218283Hi07f0112513180-218180183Hi07f011258126-167null226Hi07f021158126-167null226Hi122f1255231-245231233Hi23g12y1555231-245231233Hi23g12y18813152-124101226MDA1761-SSR182231-245231233Hi23g12y1883221-224101155MDA1761-SSR18913155-192null155MDA1761-SSR1823135-192null155NS06g0310913155-192null155N202b0115na10212-236212214N205b0115

 $\underline{\textcircled{O}}$ Springer

Multiplex PCRs

Very high levels of efficiency in genotyping have been obtained by developing PCR multiplexes. As 1,997 DNA samples were genotyped with 82 SSRs organized in 21 MPs, 121,817 PCRs and fragment analyses were saved from a total of 163,754. This allowed a considerable (approximately 75%) saving on labor, time, and consumables. The efficiency of the current MP can still be further increased. Where the current ABI-MP were based on a three-dye system, now up to four different dyes can be used. The increased knowledge about size ranges also allows the number of SSRs within MPs to be increased, as SSRs with well-distinct allele ranges can be included using the same dye. This approach has already systematically successfully applied with the LI-COR-MPs where only two different dyes can be used, as well as with some ABI-MPs.

ABI-platforms proved to be much more useful for high through SSR-genotyping than the LI-COR system, both in terms of number of SSR per MP as in consistency in allele sizing and scoring. The information on compatibility of SSRs from the LI-COR system could be used to develop ABI-MPs for these SSRs, for laboratories with the ABI system only.

As far as possible, SSRs mapping on the same LG have been placed in the same or in a maximum of two MPs. This approach was generally successful, as only the SSRs mapping on LG4, LG5, and LG13 had to be placed in more than two MPs. SSRs organized in this way are very useful in view of data validation during the genotyping process as well as for some applications, e.g., to add an SSR backbone to a specific LG of a genetic map based on poorly transferable (AFLP) markers where a QTL has been identified or for genome scanning approaches aiming at the identification of markers associated with a specific trait.

Multilocus SSRs in the set

For more than half of the ML or presumed ML SSRs, either no second locus was amplified (four cases) or it was possible to find a simple way to distinguish the alleles belonging to the different loci (four cases, Table 2). This demonstrates that ML and presumed ML SSRs need not be excluded *a priori* from large genotyping projects; however, as large genotyping projects by definition require high throughput, the splitting of the data must be as robust as possible. If this is not the case, those SSRs should be substituted as soon as alternative single-locus markers are available.

Degree of polymorphism of the tested SSRs

Coart et al. (2003) found a rough correlation between the number of alleles identified in a small group of accessions

and the total number of alleles present in a wide germplasm collection. Liebhard et al. (2002) and Silfverberg-Dilworth et al. (2006) estimated the level of polymorphism of their SSRs on the basis of eight and nine guite unrelated cultivars, respectively. Enlarging the test to more than 350 genotypes showed that all the 34 single-locus SSRs amplifying at least seven different alleles in a small set of accessions amplified at least ten different alleles in the larger screen (Table 3). The average number of alleles amplified by these 34 SSRs is 15.5. The average number of alleles of the 36 single-locus SSRs amplifying less than seven alleles in the small set is 12.3, indicating that, among these SSRs, highly polymorphic SSRs (SSRs with at least ten alleles) are also present. Indeed, the probability of finding them would decrease only slightly, namely from 100% to about 72% (26 out of 36 SSRs). It may thus be wise to include such putatively low polymorphic SSRs in directed genotyping approaches if no alternative markers are available.

Conclusions

A large set of highly variable, easy-to-score SSRs has been identified that are well distributed over the whole apple genome. The SSRs have been organized in multiplex PCRs and used to generate, in a fast and cost-efficient way, solid genotyping data of around 2,000 apple accessions. These data will now be combined with different parameters of fruit quality (e.g., firmness, acidity) in a PBA to identify QTL for these traits.

Acknowledgments This publication was carried out with the financial support from the Commission of the European Communities (contract no. QLK5-CT-2002-01492), Directorate-General Research—Quality of Life and Management of Living Resources Program. It does not necessarily reflect the Commission's views and in no way anticipates its future policy in this area. Its content is the sole responsibility of the publishers.

References

- Bink MCAM, Uimari P, Sillanpää MJ, Janss LLG, Jansen RC (2002) Multiple QTL mapping in related plant populations via a pedigree-analysis approach. Theor Appl Genet 104:751–762
- Coart E, Vekemans X, Smulders MJM, Wagner I, Van Huylenbroeck J, Van Bockstaele E, Roldán-Ruiz I (2003) Genetic variation in the endangered Wild apple (*Malus sylvestris* (L.) Mill.) in Belgium as revealed by AFLP and microsatellite markers. Consequences for conservation. Mol Ecol 12:845–857
- Fernández-Fernández F, Harvey NG, James CM (2006) Isolation and characterization of polymorphic microsatellite markers from European pear (*Pyrus communis* L.). Mol Ecol Notes 6:1039–1041
- Gianfranceschi L, Seglias N, Tarchini R, Komjanc M, Gessler C (1998) Simple sequence repeats for the genetic analysis of apple. Theor Appl Genet 96:1069–1076

- Gianfranceschi L, Soglio V (2004) The European project HiDRAS: innovative multidisciplinary approaches to breeding high quality disease resistant apples. Acta Hort 663:327–330
- Guilford P, Prakash S, Zhu JM, Rikkerink E, Gardiner S, Bassett H, Forster R (1997) Microsatellites in *Malus x domestica* (apple): abundance, polymorphism and cultivar identification. Theor Appl Genet 94:249–254
- Hemmat M, Weeden NF, Brown SK (2003) Mapping and evaluation of *Malus* x *domestica* microsatellites in apple and pear. J Am Soc Hortic Sci 128:515–520
- Hokanson SC, Szewc-McFadden AK, Lamboy WF, McFerson JR (1998) Microsatellite (SSR) markers reveal genetic identities, genetic diversity and relationships in a *Malus x domestica* Borkh. core subset collection. Theor Appl Genet 97:671–683
- Igarashi M, Abe Y, Hatsuyama Y, Ueda T, Fukasawa-Akada T, Kon T, Kudo T, Sato T, Suzuki M (2008) Linkage maps of the apple (*Malus x domestica* Borkh.) cvs 'Ralls Janet' and 'Delicious' include newly development EST markers. Mol Breed (in press)
- Inoue E, Matsuki Y, Anzai H, Evans K (2007) Isolation and characterization of microsatellite markers in Japanese pear (*Pyrus pyrifolia* Nakai).. Mol Ecol Notes 7:445–447
- Khan MA, Durel CE, Duffy B, Drouet D, Kellerhals M, Gessler C, Patocchi A (2007) Development of Molecular markers linked to the 'Fiesta' linkage group 7 major QTL for fire blight resistance and their application for marker-assisted selection. Genome 50:568–577
- Liebhard R, Gianfranceschi L, Koller B, Ryder CD, Tarchini R, Van de Weg E, Gessler C (2002) Development and characterisation of 140 new microsatellites in apple (*Malus x domestica* Borkh.). Mol Breed 10:217–241
- Liebhard R, Koller B, Gianfranceschi L, Gessler C (2003) Creating a saturated reference map for the apple (*Malus x domestica* Borkh.) genome. Theor Appl Genet 106:1497–1508
- Patocchi A, Gessler C (2003) Genome scanning approach (GSA), a fast method for finding molecular markers associated to any trait. In: Proceedings of the plant and animal genomes XI conference (Abstract) San Diego, California, p 187
- Patocchi A, Fernandez-Fernandez F, Evans K, Silfverberg-Dilworth E, Matasci CL, Gobbin D, Rezzonico F, Boudichevskaia A,

Dunemann F, Stankiewicz-Kosyl M, Matisse F, Soglio V, Gianfranceschi L, Durel CE, Toller C, Cova V, Mott D, Komjanc M, Barbaro E, Costa F, Voorrips RE, Yamamoto T, Gessler C, Van de Weg WE (2007) Development of a set of apple SSRs markers spanning the apple genome, genotyping of HiDRAS plant material and validation of genotypic data. Acta Hort (in press)

- Pierantoni L, Cho KH, Shin IS, Chiodin RI, Tartarini S, Dondini L, Kang SJ, Sansavini S (2004) Characterisation and transferability of apple SSRs to two European pear F₁ populations. Theor Appl Genet 109:1519–1524
- Silfverberg-Dilworth E, Matasci CL, Van de Weg WE, Van Kaauwen MPW, Walser M, Kodde LP, Soglio V, Gianfranceschi L, Durel CE, Costa F, Yamamoto T, Koller B, Gessler C, Patocchi A (2006) Microsatellite markers spanning the apple (*Malus x domestica* Borkh.) genome. Tree Genet Genomes 2:202–224
- Van de Weg WE, Voorrips RE, Finkers R, Kodde LP, Jansen J, Bink MCAM (2004) Pedigree genotyping: a new pedigree-based approach of QTL identification and allele mining. Acta Hort 663:45–50
- Van Dyk MM (2008) Identification of quantitative trait loci controlling the requirement for chilling in vegetative budbreak in apple (*Malus x domestica* Borkh.). Ph.D. Thesis University of the Western Cape, South Africa
- Vinatzer BA, Patocchi A, Tartarini S, Gianfranceschi L, Sansavini S, Gessler C (2004) Isolation of two microsatellite markers from BAC clones of the Vf scab resistance region and molecular characterization of scab-resistant accessions in Malus germplasm. Plant Breed 123:321–326
- Voorrips RE (2007) Pedimap: software for visualization of genetic and phenotypic data in pedigrees. Plant Research International, Wageningen
- Volz R, Rikkerink E, Sustin P, Lawrence T, De Silva N, Bus V (2007) "Fast-breeding" in apple: a strategy to accelerate introgression of new traits into elite germplasm. Acta Hort (in press)
- Yamamoto T, Kimura T, Shoda M, Imai T, Saito T, Sawamura Y, Kotobuki K, Hayashi T, Matsuta N (2002) Genetic linkage maps constructed by using an interspecific cross between Japanese and European pears. Theor Appl Genet 106:9–18