



POLITECNICO DI TORINO
Repository ISTITUZIONALE

Data from cyclic tensile tests on sutured organs to evaluate creep behaviour, distraction, and residual thread strength

Original

Data from cyclic tensile tests on sutured organs to evaluate creep behaviour, distraction, and residual thread strength / Pascoletti, G.; Pressanto, M. C.; Putame, G.; Terzini, M.; Franceschini, G.; Zanetti, E. M.. - In: DATA IN BRIEF. - ISSN 2352-3409. - 30(2020), p. 105644.

Availability:

This version is available at: 11583/2827727 since: 2020-05-21T09:21:05Z

Publisher:

Elsevier Inc.

Published

DOI:10.1016/j.dib.2020.105644

Terms of use:

openAccess

This article is made available under terms and conditions as specified in the corresponding bibliographic description in the repository

Publisher copyright

(Article begins on next page)



Data Article

Data from cyclic tensile tests on sutured organs to evaluate creep behaviour, distraction, and residual thread strength



Giulia Pascoletti^a, Maria Chiara Pressanto^b, Giovanni Putame^{c,d},
Mara Terzini^{c,d}, Giordano Franceschini^a, Elisabetta M. Zanetti^{a,*}

^a Department of Engineering, University of Perugia, Italy

^b Dip. di Medicina Veterinaria – Sezione Chirurgia e Radiodiagnostica, University of Perugia, Italy

^c Polito^{BIO} MedLab, Politecnico di Torino, Italy

^d DIMEAS, Politecnico di Torino, Italy

ARTICLE INFO

Article history:

Received 21 April 2020

Accepted 24 April 2020

Available online 30 April 2020

Keywords:

Suture testing

Distraction

Failure

Creep

Neuropathy

Tissue relaxation

ABSTRACT

A number of applications in the surgical practice are based on tensile sutures aimed to keep soft tissues in place and compensate the exit of neuropathies, prolapses or general tissue relaxation. Long-term behaviour of these constructs need to be carefully examined in order to define tensile forces to be applied and to compare different suture anchors.

Data here reported refer to equine laryngoplasties, where a suitable loading system has been designed in order to be able to test sutures in-situ, applying known forces (“On-site testing of sutured organs: an experimental set up to cyclically tighten sutures” (Pascoletti et al., 2020 [1])). The loading protocol was made of two steps: in the first step, 3000 loading cycles have been performed; in the following step, a tensile test up to rupture was performed.

Cyclic load/displacement curves allow evaluating suture distraction, as a consequence of suture migration and/or soft tissues creep. Tensile curves allow evaluating the residual thread strength and its ultimate displacement.

* Corresponding author.

E-mail address: elisabetta.zanetti@unipg.it (E.M. Zanetti).

These data can provide a detailed insight of long-term suture behaviour and can be a reference to compare different threads and/or suture anchors.

© 2020 The Authors. Published by Elsevier Inc.

This is an open access article under the CC BY license.

(<http://creativecommons.org/licenses/by/4.0/>)

Specifications table

Subject	Biomechanics
Specific subject area	Sutures
Type of data	Raw Data in a Table
How data were acquired	Data were acquired performing mechanical tests. The loading machine is Instron Electropulse E3000. The samples underwent 3000 loading cycles, displacement-controlled (rate equal to 1 mm/s), between 30 N and 50 N force limits. Finally, the samples underwent a tensile test up to rupture (displacement rate equal to 1 mm/s). The suture was stretched in-situ with known loads, thanks to the system described in a previous work Pascoletti et al. [1]. Samples were equine laryngoplasties performed by an expert surgeon; they were chosen as representative of prosthetic sutures used to lift soft tissues.
Data format	Raw data concerning 8 samples.
Parameters for data collection	Data have been sampled at 100 Hz
Description of data collection	Data come from the output of the loading machine: force data were measured by a Dynacell biaxial dynamic load cell (axial load range ± 5 kN; accuracy equal to 0.10% in the 0–60 N load range here used, according to the most recent calibration curve); displacement was measured by an optical encoder.
Data source location	PolitoBIOMedLab Politecnico di Torino Torino (ITALY) 45°03'58.1"N 7°39'30.1"
Data accessibility	Repository name: Mendeley Data identification number: doi:10.17632/zpf2jrtgxs.1 Direct URL to data: https://data.mendeley.com/datasets/zpf2jrtgxs/1
Related research article	Giulia Pascoletti, Maria Chiara Pressanto, Giovanni Putame, Mara Terzini, Alberto L. Audenino, Elisabetta M. Zanetti On-site testing of sutured organs: an experimental set up to cyclically tighten sutures Journal of Mechanical Behaviour of Biomedical Materials In Press

Value of the data

- These data provide an in-sight into the long-term behaviour of tensile sutures.
- Information here shared can be useful to all surgeons working on soft tissues with localised loads, provided by tensioned sutures: aesthetic surgeons working with thread lifts [2] otorhinolaryngologists performing sling arytenoid adduction [3] or arytenoid abduction lateropexy in neonatal care [4]. Surgeons performing blepharoptosis [5,6], surgeons correcting various kinds of prolapses with sutures [7], and so on.
- These data can be compared to data obtained with different threads or different thread anchorages (suture buttons or suture anchors [8]).

- There is paucity of 'In vitro' data concerning suture distraction: most data in literature have been obtained with set up which cannot allow estimating the exact load being carried by the suture, unless the suture alone or the suture constrained on one single end has been tested. This experimental set up allows testing the suture 'on-site' so that the compliance and creep behaviour of both its end points is taken into account.

1. Data description

This article includes raw data (organised in 16 files and uploaded on Mendeley repository: <https://data.mendeley.com/datasets/zpf2jrtgxs/1>) and two figures. With reference to raw data, file names have been assigned in this way:

- TestxCyclic.txt: sample x (with x ranging from 1 to 8), undergoing cyclic loading,
- TestxRamp: sample x (with x ranging from 1 to 8), undergoing tensile test.

TestxCyclic.txt files are made of 7 columns:

1. Total Time (s)
2. Cycle Elapsed Time (s): time elapsed since the beginning of a new cycle
3. Step: '1' when the load is rising; '2' when the load is lowering
4. Loop1(1): number of performed cycles
5. Position (mm): absolute position
6. Force (N)
7. Displacement (Linear:Digital Position) (mm): relative displacement, taking zero-load position as a reference

TestxRamp.txt files are made of 5 columns:

1. Total Time (s)
2. Cycle Elapsed Time (s): time elapsed since the beginning of a new cycle
3. Position (mm): absolute position
4. Force (N)
5. Displacement (Linear:Digital Position) (mm): relative displacement, taking zero-load position as a reference.

Fig. 1 illustrates the loading set up: the organ is laid on its support; thread ends cross, they pass through lower pulleys and are knotted above an upper pulley, moved by the loading machine.

Fig. 2 reports force/displacement patterns for cyclic load (Fig. 2(a)) or for the final load-to-failure test (Fig. 2(b)).

2. Experimental design, materials, and methods

The objective of this research has been evaluating the distraction (or creep behaviour) and residual strength of a tensile suture undergoing 3000 cyclic tests.

3. Materials

Equine laryngoplasties have been used as a benchmark. This surgery is performed in case of recurrent laryngeal neuropathy and it is aimed to provide enough airflow for the horse. The prosthesis is actually a suture which is placed between the caudal aspect of the cricoid cartilage and the muscular process of the arytenoid cartilage to simulate the action of the cricoarytenoideus dorsalis muscle.

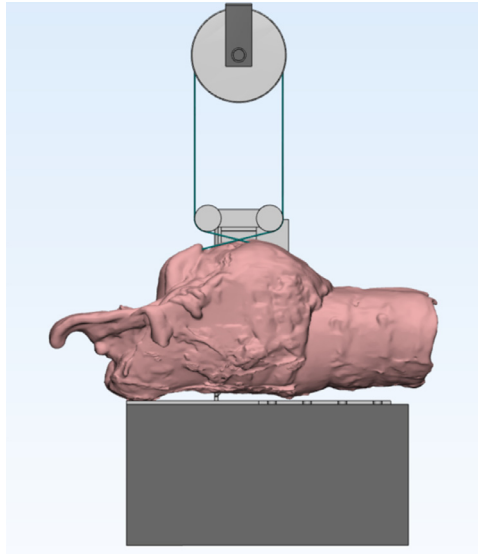


Fig. 1. Experimental set-up.

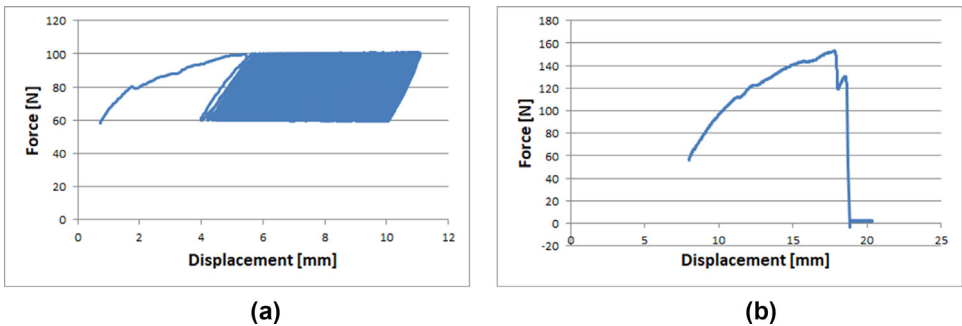


Fig. 2. Sample curves obtained during cyclic tests (a) or during the final tensile test (b).

Larynges samples were collected from a local abattoir. They belonged to 5 male and 3 female horses whose age ranged between 3 and 13 years. Laryngoplasties were performed by an experienced veterinary surgeon, who used Ethibond UPS 6 thread.

A suitable experimental set-up has been designed in order to test the suture on-site. It is depicted in Fig. 1.

The organ is securely fastened to the lower base of the loading machine. Lower pulleys allow stressing the suture along its axis, independently from the thread elongation. The load applied by the machine acts on the upper pulley and it is equal to twice the thread tension. Testing the suture on-site allows taking into account the local compliance of both suture ends at the respective anchorage point on the organ.

4. Methods

The loading protocol has been so established:

- Step 1: 3000 loading cycles, displacement controlled at 1 mm/s, between two load limits (30 N and 50 N).
- Step 2: tensile test, up to rupture, displacement controlled at 1 mm/s.

The value of load limits and displacement rate have been taken from physiological conditions, according to in vivo measurements [9].

Sample curves have been reported in Fig. 2.

Acknowledgments

This research was partially funded on the Project “Metodi di progettazione e verifica di dispositivi protesici e di sintesi ossea innovative”, fund “Ricerca di Base, 2018” by University of Perugia.

Conflict of interest

The authors declare that they have no known competing financial interests or personal relationships which have, or could be perceived to have, influenced the work reported in this article.

References

- [1] G. Pascoletti, C. Pressanto, G. Putame, M. Terzini, G. Franceschini, E.M. Zanetti, On-site testing of sutured organs: an experimental set up to cyclically tighten sutures, *J. Mech. Behav. Biomed. Mater.* (2020) In Press, doi:[10.1016/j.jmbbm.2020.103803](https://doi.org/10.1016/j.jmbbm.2020.103803).
- [2] S. Halepas, X.J. Chen, E.M. Ferneini, Thread-lift sutures: anatomy, technique, and review of current literature, *J. Oral Maxillofac. Surg.* (2020), doi:[10.1016/j.joms.2019.11.011](https://doi.org/10.1016/j.joms.2019.11.011).
- [3] M. Hess, D. Schroeder, K. Püschel, Sling arytenoid adduction, *Eur. Arch. Oto-Rhino-Laryngol.* 268 (2011) 1023–1028, doi:[10.1007/s00405-010-1429-2](https://doi.org/10.1007/s00405-010-1429-2).
- [4] V. Matievics, A. Bach, B. Sztano, Z. Bere, Z. Tobias, P.F. Castellanos, A.H. Mueller, L. Rovo, Functional outcomes of endoscopic arytenoid abduction lateropexy for unilateral vocal cord paralysis with dyspnea, *Eur. Arch. Oto-Rhino-Laryngol.* 274 (2017) 3703–3710, doi:[10.1007/s00405-017-4696-3](https://doi.org/10.1007/s00405-017-4696-3).
- [5] M. Bouazza, M. Elbelhadji, A. McHachi, L. Benhmidoune, A. Amraoui, Traitement du ptosis congénital par suspension au muscle frontal par fil de polypropylène: résultats d'une étude de 21 cas, *J. Fr. Ophtalmol.* 37 (2014) 520–525, doi:[10.1016/j.jfo.2014.01.014](https://doi.org/10.1016/j.jfo.2014.01.014).
- [6] R. Gazzola, E. Piozzi, L. Vaienti, F. Wilhelm Baruffaldi Preis, Therapeutic algorithm for congenital ptosis repair with levator resection and frontalis suspension: results and literature review, *Semin. Ophthalmol.* 33 (2018) 454–460, doi:[10.1080/08820538.2017.1297840](https://doi.org/10.1080/08820538.2017.1297840).
- [7] T.E. David, C.M. David, M. Lafreniere-Roula, C. Manlhiot, Long-term outcomes of chordal replacement with expanded polytetrafluoroethylene sutures to repair mitral leaflet prolapse, *J. Thorac. Cardiovasc. Surg.* (2019), doi:[10.1016/j.jtcvs.2019.08.006](https://doi.org/10.1016/j.jtcvs.2019.08.006).
- [8] L. Pirtea, O. Balint, C. Secoşan, D. Grigoraş, R. Iliina, Laparoscopic pectopexy with burch colposuspension for pelvic prolapse associated with stress urinary incontinence, *J. Minim. Invasive Gynecol.* (2019), doi:[10.1016/j.jmig.2019.10.022](https://doi.org/10.1016/j.jmig.2019.10.022).
- [9] T.H. Witte, J. Cheetham, L.V. Soderholm, L.M. Mitchell, N.G. Ducharme, Equine laryngoplasty sutures undergo increased loading during coughing and swallowing, *Vet. Surg.* 39 (2010) 949–956, doi:[10.1111/j.1532-950X.2010.00742.x](https://doi.org/10.1111/j.1532-950X.2010.00742.x).