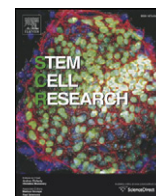


Contents lists available at [ScienceDirect](http://ScienceDirect.com)

Stem Cell Research

journal homepage: www.elsevier.com/locate/scr

Lab resource

Derivation of the clinical grade human embryonic stem cell line RCe015-A (RC-11)

P.A. De Sousa^{a,b,c,*}, B.J. Tye^a, K. Bruce^a, P. Dand^a, G. Russell^a, D.M. Collins^a, A. Greenshields^a, K. McDonald^a, H. Bradburn^a, A. Laurie^a, M. Canham^c, T. Kunath^c, J.M. Downie^a, M. Bateman^a, A. Courtney^a^a Roslin Cells Limited, Nine Edinburgh Bio-Quarter, 9 Little France Road, Edinburgh EH16 4UX, UK^b Centre for Clinical Brain Sciences, University of Edinburgh, UK^c MRC Centre for Regenerative Medicine, University of Edinburgh, UK

ARTICLE INFO

Article history:

Received 22 April 2016

Accepted 28 April 2016

Available online 7 May 2016

ABSTRACT

The human embryonic stem cell line RCe015-A (RC-11) was derived under quality assured compliance with UK regulation, European Union Directives and International guidance for tissue procurement, processing and storage according to Good Manufacturing Practice (GMP) standards. The cell line was derived from a fragmented cleavage stage embryo voluntarily donated as unsuitable or surplus to fertility requirements following informed consent. RCe015-A (RC-11) shows normal pluripotency marker expression and differentiation to the three germ layers in vitro and in vivo. It has a normal 46XX female karyotype and microsatellite PCR identity, HLA and blood group typing data are available.

© 2016 The Authors. Published by Elsevier B.V. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Resource table

Name of stem cell construct	RCe015-A
Alternative name	RC-11, RC11
Institution	Roslin Cells Ltd.
Person who created resource	B.J. Tye, K. Bruce, P. Dand, G. Russell, D.M. Collins, A. Greenshields, K. McDonald, H. Bradburn, A. Laurie, M.A. Canham
Contact person and email	Paul.desousa@roslincells.com ; Paul.desousa@ed.ac.uk Janet.downie@roslincells.com Aidan.courtney@roslincells.com Malcolm.bateman@roslinfoundation.com Tilo.kunath@ed.ac.uk
Date archived/stock date	11 March 2011 (seed bank)
Type of resource	Biological reagent: cell line
Sub-type	hESC, clinical grade
Origin	Cleavage stage embryo (Mitosis)
Key transcription factors	Oct4 (confirmed by flow cytometry),
Authentication	See Quality Control Certificate of Analysis (Fig. 1)
Link to related literature (direct URL links and full references)	N/A
Information in public databases	http://hpscreg.eu/cell-line/RCe015-A

Ethics

Informed consent obtained. Scotland A Research Ethics committee approval obtained (07/MRE00/56). Conducted under the UK Human Fertilisation and Embryology Authority licence no R0136 to centre 0202 and UK Human Tissue Authority (HTA) licensing number 22,631.

Resource details

RCe015-A (RC-11) was received as a fragmented cleavage stage embryo that was surplus to requirement or unsuitable for clinical use. Human embryonic stem cell (hESC) isolation, expansion and qualification was performed in a facilities whose specification, operation and monitoring complied with GMP standards enabling; i) a fully traceable procurement procedure with informed ethical consent which includes provision for commercial use, ii) detailed medical history and blood borne virus (BBV) screening of donors, and iii) compilation of a cell line history providing details on hESC manufacturing process and quality control testing regime.

Human ESC culture and processing was performed in a grade A tissue culture cabinet in a grade B clean room environment monitored for particulate and microbiological contamination during cell processing, in accordance with Rules and Guidance for Pharmaceutical Manufacturers and Distributors – The Orange Guide, compiled by the UK Medicines Healthcare Products Regulatory Authority (Go to: <https://www.gov.uk/guidance/good-manufacturing-practice-and-good-distribution-practice>). Accordingly, the facility was operating under a mature Quality

* Corresponding author at: Reader, Centres for Clinical Brain Sciences & Regenerative Medicine University of Edinburgh Chancellor's Building 49 Little France Crescent, Edinburgh, EH16 4SB, Scotland, UK; Chief Scientist, Executive Director, Founder Roslin Cells Ltd Scottish Centre for Regenerative Medicine, 5 Little France Drive, Edinburgh, EH16 4UU, Scotland, UK.

Management System, compliant with ISO9001:2008 standards. Further hESC derivation was performed under licensure from the UK HFEA (R0136 to centre 0202) and HTA (Licensing Number 22631).

The embryo was grown to blastocyst stage and the cell line was derived by whole embryo outgrowth on mitotically inactivated human dermal fibroblast (HDF) feeder cells. HDFs were derived and manufactured according to GMP and had been approved for clinical use by the Food and Drug Administration, USA. During derivation on HDFs, hESCs were grown in a xeno-free cell therapy grade media (XF KODMEM) supplemented with xeno-free human recombinant bFGF. HESCs were subsequently expanded in a GMP grade serum-free

medium (StemPro hESC Serum Free Medium) and xeno-free matrix (CellStart). The former contained bovine serum albumin (BSA) from a Transmissible Spongiform Encephalopathy (TSE)-free country of origin. The cell line was cryopreserved in a GMP compliant cryopreservation solution (CryoStor CS10).

By flow cytometry, RCe015-A (RC-11) expressed the pluripotency makers Oct-4, Tra-1-60 and SSEA-4 (90.0%, 87.2% and 94.9%, respectively), whereas low expression of the differentiation marker SSEA-1 (0.0%) was observed (Fig. 1, Fig. 2). Differentiation to the three germ layers, endoderm, ectoderm and mesoderm, was demonstrated using embryoid body formation in vitro, and expression of the germ layer markers

ROSLIN CELLS

Quality Control Test Certificate

Sample Point 2 Test Results

Certificate Number:	QCC-11-011	Version:	3
Grade:	CLINICAL		
Sample ID:	RC-11 P10A		

Assay	Test Method	Roslin Cells Assay Code	Date of Assay	Result
Mycoplasma Detection	RT-qPCR (SOP/QCP/22)	MYCO-10-003 MYCO-10-004	01Sep10 07Sep10	Not Detected**
Endotoxin Detection	Kinetic Chromogenic LAL (SOP/QCP/12)	ENDO-10-002 ENDO-10-003	31Aug10 01Sep10	4.14 EU/ml*** 0.921 EU/ml***
Viral Screening*	PCR (CMV,HTLV1,HIV1,HCV, HBV,EBV) (SOP/QCP/60)	N/A	09Mar11	Not Detected****
Karyotype*	G-banding (SOP/QCP/59)	N/A	22Mar11	46,XX****
Pluripotency / Differentiation	Flow Cytometry (SOP/QCP/25)	FLOW-10-002	01Sep10	% Positive
				SSEA-4 – 94.9
				Oct 3/4 – 90.0
				Tra-1-60 – 87.2
				SSEA-1 – 0
Microsatellite Genotyping*	PCR (SOP/QCP/6)	N/A	14Oct10	ID Obtained

*Subcontracted to a Third Party
**Refer to IR/0214
***Both results were generated for this line, refer to IR/0232
****Analysis performed on RC-11 P17A

Certificate Prepared by (QC):

Certificate Reviewed by (QC):

Date: 25oct11

Date: 21OCT11

Confidential

Page 1 of 2

Fig. 1. Quality Control Certificate of Analysis for RC-11 (RCe015-A) P10A seed lot.

ROSLIN CELLS

Quality Control Test Certificate
Sample Point 2 Test Results

Certificate Number:	QCC-11-011	Version:	3
Grade:	CLINICAL		
Sample ID:	RC-11 P10A		

Assay	Test Method	Roslin Cells Assay Code	Date of Assay	Result
HLA Typing*	PCR-SSO (SOP/QCP/62)	N/A	10Nov10	HLA Typed Class I and Class II
Blood Group Genotyping*	PCR (SOP/QCP/63)	N/A	14Oct10	ABO Genotype: O ⁺ O ⁺ (probable phenotype O ⁺)
Differentiation	Embryoid Body Formation (Endoderm, Ectoderm, Mesoderm) (SOP/QCP/7 & SOP/QCP/58)	EB-CUL-11-001 EB-STAIN-11-001	04Mar11 31Mar11	Endoderm – Detected**
				Ectoderm – Detected**
				Mesoderm – Detected**
Viability	Flow Cytometry (SOP/QCP/40)	VIAB-10-001	01Sep10	% Non-Viable: 10.3 ± 5.3

*Subcontracted to a Third Party
**Analysis performed on RC-11 P17A

Certificate Prepared by (QC): ASd Date: 20Oct11

Certificate Reviewed by (QC): Laura Kennedy Date: 21OCT11

Certificate Approved by (QA): Julie Hean Date: 23OCT2011

Confidential Page 2 of 2

Fig. 1 (continued).

α -fetoprotein, β -tubulin and muscle actin was observed (Fig. 3, top panel). In vivo teratoma formation yielded typical hESC derived teratomas. Histological examination of fixed and stained sections clearly showed generation of cell types from ectoderm and mesoderm lineages. Endoderm differentiation was also present, but the quality of the structures meant these could not be fully characterised (Fig. 3, bottom panels).

A microsatellite PCR profile has been obtained for the cell line, and HLA Class I and II typing is available (Table 1). Blood group genotyping gave the blood group O₁O₁, expected to give rise to blood group O⁺ (Table 1). The cell line is free from mycoplasma contamination as determined by RT-qPCR.

Verification and authentication

The cell line was analysed for genome stability by G-banding and showed a normal 46XX female genotype (Fig. 4). SNP genotyping was carried out using the Illumina HumanCytoSNP-12 v2.1 BeadChip and revealed a 3.6 Mb Copy-Neutral Loss of Heterozygosity (CN-LOH) on chromosome 2 (Chromosome 2p16.2-16.1) as described in Canham et al., 2015. While CN-LOH can be constitutional acquired CN-LOH can occur during oncogenic transformation however such changes are almost exclusively telomeric or greater than 25 Mb when interstitial (Kearney et al., 2011; Stephens et al., 2006). This region was cross-referenced to the Genomic Imprinting database (<http://www.>

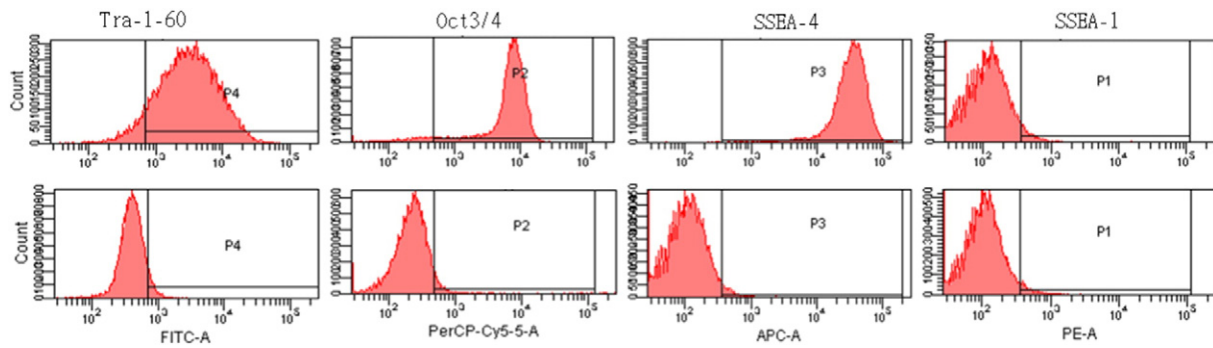


Fig. 2. RCE015-A (RC-11) was subjected to flow cytometry analysis for markers of pluripotency with specific antibody (top row) or isotype control (bottom row) as indicated above the histograms. Percentage staining is indicated in the Certificate of Analysis (Fig. 1).

geneimprint.com) (Falls et al., 1999) and was found not to reside in an imprinted region, although *CCDC85A* in the CN-LOH is predicted, but not validated, to be imprinted on this database.

Materials and methods

Ethics

Derivation of hESC from surplus to requirement and failed to fertilise/develop embryos was approved by The Scotland A Research

Ethics Committee and local ethics board at participating fertility clinics and conducted under licence no R0136 from the UK HFEA with informed donor consent. The processing and storage of hESC cells for human application was conducted under Licence Number 22631 from the UK Human Tissue Authority.

Cell culture

Fresh embryos were cultured in Sydney cleavage medium (Cook Medical, Hertfordshire, UK) until day 3 and Sydney blastocyst medium

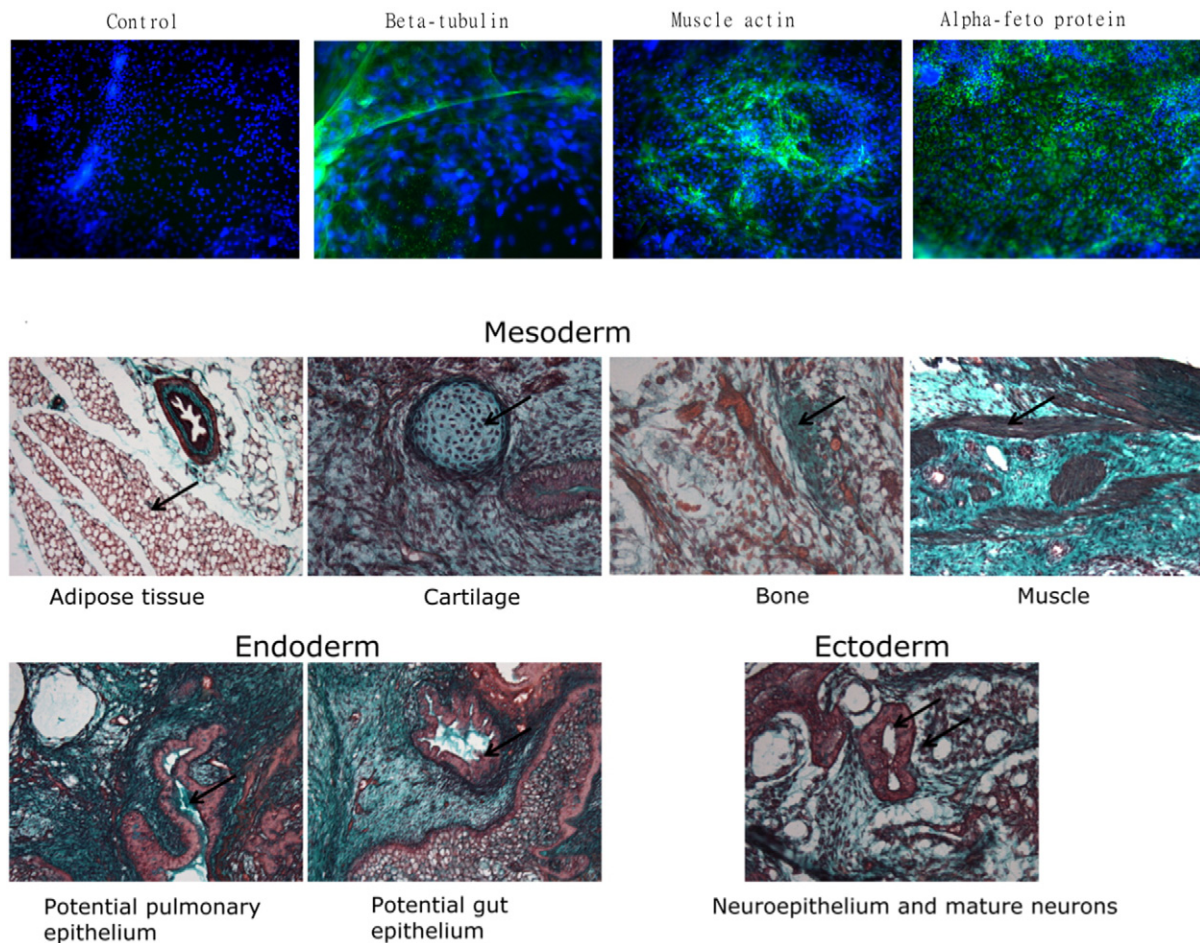


Fig. 3. RCE015-A (RC-11) can differentiate into the three germ layers in vitro and in vivo. Embryoid body mediated differentiation (top panel) resulted in expression of β -tubulin III (ectoderm), muscle actin (mesoderm), and α -fetoprotein (endoderm). Specific staining is shown in green, cell nuclei are counterstained with DAPI (blue). Histological sections of teratomas formed under the kidney capsule of mice yielded clear evidence for mesoderm and ectoderm differentiation and indicated that endoderm structures were also present.

Table 1

Microsatellite PCR, blood group and HLA tissue typing results for RCE015-A (RC-11).

<i>Microsatellite PCR results</i>							
D3S1358 1	D3S1358 2	vWA 1	vWA 2	D16S539 1	D16S539 2	D2S1338 1	D2S1338 2
15	16	16	17	11	13	18	19
Amelogenin 1	Amelogenin 2	D8S1179 1	D8S1179 2	D21S11 1	D21S11 2	D18S51 1	D18S51 2
X	X	13	14	29	31.2	18	20
D19S433 1	D19S433 2	THO1 1	THO1 2	FGA 1	FGA 2	CSF1PO 1	CSF1PO 2
14	14	6	9	23	23	10	12
D5S818 1	D5S818 2	D7S820 1	D7S820 2	D13S317 1	D13S317 2	TPOX 1	TPOX 2
11	12	10	13	10	12	8	11
<i>Blood group genotyping</i>							
RhD	RhC	Rhc	RhE	Rhe	Fy a	Fy b	Fy GATA
pos	pos	pos	pos	pos	neg	pos	neg
Jka	Jkb	K	k	M	N	S	s
pos	pos	neg	pos	pos	neg	neg	pos
Do a	Do b	ABO					
pos	neg	O1O1					
<i>HLA tissue typing</i>							
HLA class I type	HLA-A*01, A*24; B*07, B*08; Cw*07						
HLA class II type	HLA-DRB1*03, DRB1*15; DRB3*01; DRB5*01; DQB1*02, DQB1*06; DPB1*01, DPB1*04						
Comment	DRB1*03 is expressed serologically as DR17, Unable to exclude DPB1*0602, 8201.						

(Cook Medical) after day 3 of development. Embryos were cultured at 36.5–37.5 °C, 5.0 ± 0.5% CO₂, 5.0 ± 0.5% O₂ in drops under paraffin oil (Cook Medical) and transferred to fresh medium at least every 2–3 days.

By day 8 of development, embryos were placed in derivation conditions consisting of mitotically inactivated GMP grade neonatal human dermal fibroblasts (HDFs) (Forticell Biosciences, NJ, USA) on tissue



Fig. 4. RCE015-A (RC-11) was analysed by Giesma staining of 30 metaphase spreads and showed a normal 46XX female karyotype.

culture plastic in XF KODMEM medium (Knockout-DMEM, 15% KOSR-XF, 2 mM L-glutamine, 1% MEM Non essential amino acids, 2% XF Growth Factor Cocktail, 0.1 mM β -mercaptoethanol (ThermoFisher Scientific, Paisley, UK) supplemented with 80 ng/ml human bFGF (ThermoFisher Scientific). When available, cell therapy system quality reagents were used. Assisted hatching was performed by removing the zona pellucidae mechanically using Swemed cutting tools (Vitrolife, Göteborg, Sweden).

HDF cells were cultured in DMEM (Lonza, Slough, UK), 10% Pharma grade FCS (GE Healthcare (PAA), Buckinghamshire, UK) and 2 mM L-glutamine (ThermoFisher Scientific). HDF were mitotically inactivated using gamma irradiation at 50 Gy using a Gammacell Elite 1000 machine. For use as a feeder layer, irradiated HDFs were plated at 50,000 cells/cm² in XF KODMEM medium supplemented with 80 ng/ml human bFGF (ThermoFisher Scientific). Cells were cultured at 36.5–37.5 °C, 5.0 \pm 0.5% CO₂, 5.0 \pm 0.5% O₂ and 50% medium exchanged 6 days a week.

The established cell line was expanded and banked using CellStart matrix and Stempro hESC Serum Free Medium (cell therapy system quality reagents, ThermoFisher Scientific). This contained BSA from a TSE-free country of origin. Passaging was performed mechanically using an EZ passage tool (ThermoFisher Scientific). hESC lines were expanded to 25–30 wells of a 6-well plate and cryopreserved in 0.5–1 ml Cryostor CS10 (Biolife Solution, Washington, USA) using an EF600-107 controlled rate freezer (Grant Instruments, Cambridge, UK) before being stored in a –150 °C freezer (Panasonic Biomedical, Loughborough, UK).

Mycoplasma

In process mycoplasma detection was performed using Applied Biosystems PrepSEQ™ Mycoplasma Nucleic Acid Extraction Kit and MicroSEQ™ Mycoplasma Real-Time PCR Detection Kit (ThermoFisher Scientific (Applied Biosystems)) according to the manufacturer's instruction. European pharmacopoeia (EP) mycoplasma testing was carried out by Moredun Scientific Ltd. (Edinburgh, UK), under a quality and technical agreement.

Endotoxin

Endotoxin levels were determined using the Kinetic-QCL assay (Lonza) and an incubating plate reader (BioTek ELx808) according to the manufacturer's instructions. Briefly, an unknown sample was compared with a standard curve of known levels of control endotoxin. An assay was deemed valid if the coefficient of correlation, $r \geq 0.980$ and the CV (%) for the standard curve was $\leq 10\%$, and the reaction time of the negative control was greater than the reaction time of the lowest standard on the standard curve.

Flow cytometry

Pluripotency was determined using the Human and Mouse Pluripotent Stem Cell Analysis kit (BD, Oxford, UK). Oct 3/4 and SSEA-4 were included as pluripotency markers, and SSEA-1 as a differentiation marker. FITC conjugated Tra-1-60 (BD) was used as an additional pluripotency marker. Fixed and permeabilised cells were analysed using a FACS Aria flow cytometer (BD) or a Guava easyCyte flow cytometer (Millipore, Watford, UK). Percentage expression of each marker was compared to isotype control or unstained cells.

Viability

Viability was determined using the Guava ViaCount assay. Briefly, the Guava Viacount reagent (Millipore) containing a nuclear and a viability dye, was mixed with a single cell suspension, incubated for 5 min and analysed using the Guava easyCyte flow cytometer (Millipore).

Total cell count, viable cell count and percentage viable cells was obtained.

In vitro differentiation

hESCs were pre-treated for 1 h with 10 μ M ROCK inhibitor in Stempro hESC SFM (ThermoFisher Scientific) and embryoid bodies (EBs) generated in ultra low attachment plates (Corning) for 7 days before being transferred into EB medium (20% FBS (GE Healthcare (PAA)), 80% KO-DMEM 1 mM L-glutamine, 3.5 μ M β -mercaptoethanol, 1% non-essential amino acids (all ThermoFisher Scientific)), on glass slide tissue culture chambers (Nunc, ThermoFisher Scientific) coated with 0.2% gelatin (Sigma Aldrich, Dorset, UK) at 0.1 ml/cm² for 14 days.

Immunocytochemistry

Cells were fixed in 4% paraformaldehyde (ThermoFisher Scientific (Alfa Aesar)), permeabilised using 100% ethanol (ThermoFisher Scientific) and stained with AFP (1:500; Sigma), β -tubulin III (1:1000; Sigma) and muscle-specific actin (1:50; DAKO, Glostrup, Denmark) and secondary antibody anti mouse IgG-AlexaFluor 488 (1:200; Sigma). Images were acquired using a Zeiss S100 Axiovert fluorescence microscope or Nikon eC1 confocal microscope.

In vivo differentiation

The developmental potential of hESC lines to form teratomas consisting of tissues representative of all three germ layers was evaluated following transplantation under kidney capsule in NOD scid gamma mice. After three months, the animals were culled and assessed for teratoma formation. Teratomas were fixed in 4% paraformaldehyde, embedded in paraffin wax and serial sections of 7 μ m thickness were cut according to standard procedures. For histological assessment, the tissue sections were dewaxed, rehydrated and stained with Masson staining. Tissue sections were analysed using bright field and microscopy and digital images were recorded.

SNP genotyping and analysis

DNA samples were assayed using the Illumina HumanCytoSNP-12 v2.1 BeadChip. Genotyping data was initially assessed using GenomeStudio genotyping module (v1.94, Illumina). Karyostudio (v1.4, Illumina) was employed to perform automatic normalisation and to identify genomic aberrations utilising default settings of the built-in cnvPartition algorithm (3.07, Illumina) to generate B-allele frequency and smoothened Log R ratio plots for detected regions. These parameters are designed to detect CNVs greater than 75 kb and CN-LOH regions larger than 1 MB with a confidence value greater than 35. All identified regions were first cross-matched to the Database of Genomic Variants (DGV; <http://dgv.tcag.ca>) to identify naturally-occurring structural variations in the human. CNVs that were not identified on the DGV were then checked against a list of ES cell-associated culture adaptation genomic variants published by the International Stem Cell Initiative (Amps et al., 2011). See also Canham et al., 2015 for further details.

Genomic analysis and outsourced assays

All outsourced assays were carried out under a Quality and Technical Agreement. DNA was extracted using the QIAamp DNA Mini kit (Qiagen, Manchester, UK) according to manufacturer's recommendations and provided in recommended quantities to the service providers.

Microsatellite PCR, or Short Tandem Repeat analysis, was used to determine cell line identity and was carried out by Public Health England. A profile was obtained for the following core alleles: vWA,

D16S539, Amelogenin, TH01, CSF1PO, D5S818, D75820, D135317 and TPOX.

Human Leukocyte Antigen (HLA) tissue typing was carried out by the Scottish National Blood Transfusion Service.

Blood group genotyping was carried out by the Molecular Diagnostics laboratory at NHSBT.

Karyotype analysis was carried out by The Doctors Laboratory (London, UK) or the Western General Cytogenetics Laboratory (Edinburgh, UK). Live cells at 60–70% confluency were shipped in warm containers, fixed and analysed by standard G-banding analysis. For clinical grade lines, 30 spreads were analysed.

Viral screening for cytomegalovirus (CMV), Human T-cell lymphotropic virus (HTLV)-1, Human immunodeficiency virus (HIV)-1, Hepatitis C virus (HCV), Hepatitis B virus (HBV) and Epstein–Barr virus (EBV) was carried out by The Doctors Laboratory.

Acknowledgements

Research culminating in the derivation of this line was funded by a grant (PM07321) from Scottish Enterprise Economic Development

Agency to PDS, MB, and AC. Research culminating in SNP genotype analysis of the cell line was funded by a Medical Research Council RMRC project grant (MR/K017276/1) and The Cure Parkinson's Trust to TK and MAC.

We thank the transgenic service core facility at the Scottish Centre for Regenerative Medicine (Edinburgh) for the teratoma work.

References

- Amps, K., et al., 2011. Screening ethnically diverse hESC identifies a chromosome 20 minimal amplicon conferring a growth advantage. *Nat. Biotechnol.* 29, 1131–1144.
- Canham, M.A., et al., Nov 26 2015. The molecular karyotype of 25 clinical-grade human embryonic stem cell lines. *Sci. Rep.* 5, 17258.
- Falls, J.G., et al., 1999. Genomic imprinting: implications for human disease. *Am. J. Pathol.* 154, 635–647.
- Kearney, H.M., et al., 2011. Diagnostic implications of excessive homozygosity detected by SNP-based microarrays: consanguinity, uniparental disomy, and recessive single-gene mutations. *Clin. Lab. Med.* 31, 595–613– ix.
- Stephens, K., et al., 2006. Interstitial uniparental isodisomy at clustered breakpoint intervals is a frequent mechanism of NF1 inactivation in myeloid malignancies. *Blood* 108, 1684–1689.